dosis® BT 100

CE0197

VICTRIX AB
SVEAVÄGEN 82
S-113 59 STOCKHOLM
SVERIGE
TEL: +46-8-673 77 77
FAX: +46-8-673 77 75
Internet: www.dosis.se

SE BRUKSANVISNING
SK KÄYTTÖOHJE
DF BRUGSANVISNING
NO BRUKSANVISNING
GB INSTRUCTIONS
DE GEBRAUCHSANWEISUNG
NL GEBRUIKSAANWIJZING
FR MANUEL D'INSTRUCTIONS
ES MANUAL DE INSTRUCCIONES
IT MANUALE DI ISTRUZIONI

SVENSKA

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ENGLISH

Please read the instruction before using the device.

Normal body temperature in babies lies between 36.4 °C and 38.0 °C. Ideally you should take your baby's temperature when he/she is fully fit and write down the reading, so that in the event of illness you may determine straightaway whether your baby has a temperature or not. After the baby has had hot or cold drinks, wait 30 minutes before taking the measurement.

Instructions

Ensure an adult supervises the use of the thermometers. Note that oral measurements are usually 0,5 °C lower than the actual body temperature.

1. Take off the soother protective cap.
2. Press the ON/OFF button and you will hear a beep.
3. The last temperature reading appears for about 2 seconds on the LCD display. Then the "L °C" appears in the display.
4. The thermometer is ready for use when "°C" flashes in the display.
5. Place the soother in the baby's mouth and leave it there until the thermometer beeps four times. The display shows now the measured temperature and the "°C" lights continuously.
6. Remove the soother from the baby's mouth.
7. Press the ON/OFF button to switch off the thermometer. The device will automatically switch off 10 min. after the last reading.

The thermometer does not measure!

Press the ON/OFF button and you will hear a beep. Dip the soother in warm water until approx. 38°C appears in the display. Take the soother out of the water, wait until the thermometer beeps four times and then switch off the thermometer. To make a new measurement, press the ON/OFF button etc...

Changing the battery

Once the display fades or the thermometer cannot be switched on, the battery need to be changed. Take off the lid very gently. Use an insulated probe to remove the battery (a toothpick for example). Insert the new battery (LR41) with the positive side (+) facing upwards. Put the lid back in place.

Care

Switch off the thermometer before cleaning and never use any chemical solvents. The soother is made from non-poisonous silicone and plastic.

- Before each use check the nipple of the soothers for tears or cracks and pull it to ensure that it is firmly attached to the shield. Do not use it if a defect is found or if it is loose.
- After the reading, clean the soother in the water and disinfect with isopropyl alcohol.
- Do not clean the soother in boiling water, the LCD display is not waterproof.
- Do not stand on the soother, as it may break.
- Never expose the soother to direct sunlight or extreme heat, as it might damage the display.
- Store the soother out of reach of children; it is not a toy.
- Authorized specialists may only carry out repairs.

Warranty

The unit has been carefully checked for defects. If nevertheless you do have cause for complaint, please send us the device with your proof of purchase. We offer 1 year warranty from the date of purchase. We are not liable for damage arising from incorrect handling, improper use or wear and tear. We reserve the right to make technical modifications.

Notes on Disposal







Old units, marked with the symbol as illustrated, may not be disposed of in the household rubbish. Used rechargeable and non-rechargeable batteries, which are marked with one of the symbols illustrated, may not be disposed of in the household rubbish. You must take them to a collection point for old units, old batteries or special waste (enquire at your local authority) or the dealer from whom you bought them. These agencies will ensure environmentally friendly disposal.

Technical data

Temperature range: 32 °C–42 °C

Type of battery: LR 41, G 3-A, SR 41 or UCC 392

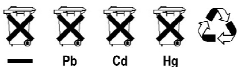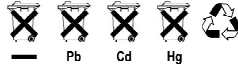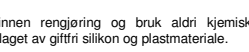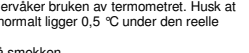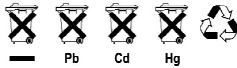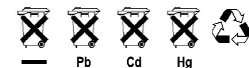Battery life: In excess of 200 hours when used continuously

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