FaceFX™

IMETEC

www.imetec.com

Tenacta Group S.p.A. Via Piemonte, 5/11 24052 AZZANO S. PAOLO (BG) ITALY

Tel. +39.035.688.111 Fax +39.035.320.149

Istruzioni per l'uso
Instrucciones de uso
Operating instructions
Navodila za uporabo
Instruções para o uso

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M1000574
0412 (MMVY)

SKIN TREATMENT APPLIANCE INSTRUCTIONS FOR USE

Dear customer, IMETEC thanks you for purchasing this product. We are certain you will appreciate the quality and reliability of this appliance, designed and manufactured with customer satisfaction in mind. These instructions for use have been drawn up in compliance with the European Standard EN 62079.

IMPORTANT! Instructions and warnings for a safe use.

Before using this appliance, carefully read the instructions and in particular the safety warnings, which must be complied with. Keep this manual and illustrative guide with the appliance for future consultation. Should you pass the appliance on to another user, make sure to also include this documentation.

Should certain parts of this user manual be difficult to understand or should doubts arise, contact the company before using the product, at the address indicated on the last page.

INDEX

Safety instructions	Page 36
Symbols	Page 40
Contraindications for use	Page 40
Possible adverse effects	Page 41
Description of appliance and accessories	Page 42
Treatment areas	Page 43
Treatment programme	Page 43
Use	Page 44
Maintenance	Page 46
Troubleshooting	Page 47
Disposal	Page 49
Assistance and warranty	Page 51
Illustrative guide	I
Technical data	II

SAFETY INSTRUCTIONS

- **This appliance cannot be used by persons under the age of 16 and can only be used by persons with reduced physical, sensorial or mental abilities or persons lacking experience if they have received**

prior instructions concerning its safe use and if they are informed about the dangers associated with this product. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

- If the appliance is used in the bathroom, remove the plug from the socket after use, as bringing it close to water is dangerous even if the appliance is switched off.**
- For enhanced safety, it is recommended to install a differential circuit breaker with a nominal operating differential current that does not exceed 30 mA (10 mA is ideal) in the electric circuit of the bathroom. Ask your installer for advice.**

IMPORTANT! Do not use this appliance near bath tubs, wash basins, showers or other water containers! Never place the appliance in water!

- After having removed the appliance from its packaging, verify its integrity based on the drawing and check for any damage caused during transport. If in doubt, do not use the appliance and contact the authorised service centre.
- Do not let children play with the packaging! Keep the plastic wrapping out of the reach of children; risk of suffocation!
- Before connecting the appliance, please make sure that the voltage indicated on the identification plate corresponds to that of the mains voltage. The identification plate is found on the appliance.
- This appliance must be used solely for the purpose it has been designed for, that is as a face and neck skin treatment appliance for domestic use. Any other use is improper and potentially dangerous. This appliance is not intended for commercial or professional use. Do not use it intensively.

- **⊗ DO NOT use the appliance with wet hands, damp feet or barefoot.**
- **⊗ DO NOT pull the power cord or the appliance itself to remove the plug from the socket.**
- **⊗ DO NOT expose the appliance to humidity or weather conditions (rain, sun).**
- **Always disconnect the appliance from the mains before cleaning or performing maintenance and when the appliance is not in use.**
- **If the appliance is faulty or malfunctions, switch it off and do not tamper with it. Any repairs must only be carried out by an authorised technical service centre.**
- **If the power cord is damaged, it must be replaced by an authorised technical service centre, so as to prevent all risks.**
- **This appliance is not suitable to be used in the presence of flammable anaesthetic mixtures with air or Oxygen or Nitrous Oxide.**
- **⊗ DO NOT use the appliance if the adaptor (8) is cracked or broken.**
- **Use the appliance only with the PART N°: UE120503GWAD01-P power supply provided.**
- **⊗ DO NOT obstruct the ventilation filters (5).**
- **⊗ DO NOT shake, drop or bang the appliance.**
- **The appliance is temperature-sensitive during use. Do not expose it to temperatures below 15 °C or above 35 °C during operation.**

- **The appliance cannot withstand extreme changes in temperature. When moved from a cold room (below 10 °C) to a warm one, wait for approximately 3 hours before using it.**
- **Do not cover the appliance and the adaptor while charging.**

SYMBOLS

Warning.

Double insulation.

Prohibited.

Only for indoor use.

Power supply connector polarity.

CONTRAINDICATIONS FOR PRODUCT USE

Some conditions can temporarily prevent the the use of the appliance. **DO NOT USE** under the following conditions:

- Dark brown or black marks, such as freckles, birth marks, moles or warts in the treatment area.
- Eczema, psoriasis, lesions, open wounds or active infections, such as herpes, in the treatment area. Wait for the affected area to heal before using the appliance.
- Abnormal skin conditions due to diabetes or other systemic or metabolic disorders.

- During or after Alpha Hydroxy Acid (AHA), Beta Hydroxy Acid (BHA), Retin-A®, topical retinoid or azelaic acid treatments.
- In the six months following Accutane® (isotretinoin) treatment.
- In the following three months following steroid treatment.
- In the event of past history of herpes in the treatment area, unless a doctor was consulted and prophylactic measures were taken before and after use.
- In the event of doubts on the possibility of using the appliance, consult your doctor or dermatologist.

POSSIBLE SIDE EFFECTS

Side effects and complications arising from the use of the appliance are rare if the instructions are complied with. However, every cosmetic procedure, including those intended for domestic use, involves a certain level of risk. This is why it is important to understand and accept the risks and complications that can occur while using energy systems intended for domestic use. Treatments with this appliance are generally well tolerated. Most users feel moderate heat transmitted by the appliance during use, without particular discomfort.

- **Minor skin irritation**
The skin can get slightly irritated after the use of the appliance, but it is a typical and temporary phenomenon. If the irritation does not disappear within three days, consult your doctor or dermatologist.
- **Skin redness**
The skin can become red after the treatment. The redness should disappear within a few hours.

DESCRIPTION OF THE APPLIANCE AND ACCESSORIES

Thus is a manual appliance with light-emitting diodes (LED) to be used for facial treatments intended for reducing periorbital wrinkles, reduce pore size and temporarily improve skin appearance.

1. Handpiece
2. LED treatment area
3. On /off switch 
4. Indicator light panel
 - A. Green:
 - Fixed - Ready for treatment
 - Blinking - Low charge
 - B. Yellow: System error – charge the unit and try again
 - C. Orange: Battery charging
 - D. Green: Battery charged
5. Ventilation filters
6. Charging base
7. Charging connection
8. AC/DC adaptor with power cord

For the appliance features, refer to the external packaging.
To be used only with the accessories supplied.

TREATMENT AREAS

The typical treatment areas for this appliance are: cheeks, forehead, neck and eye contours.

TREATMENT PROGRAMME

Results can vary from person to person, depending on the frequency of use. Some can see results sooner than others.

A single session normally lasts for 15-20 minutes, during which, the handpiece must be moved gently and slowly on the area to be treated. During this procedure, you should feel a slight heat penetrating deep into the skin.

Initial phase of the treatment for the first month (4 weeks), perform the treatment 3 times a week for 15-20 minutes per session.

Maintenance phase: use the appliance once a week for 15-20 minutes. The correct use of the appliance, in compliance with the instructions, will result in thin line improvement, pore size reduction and facial appearance improvement. The best results can be noticed 3-7 weeks after completing the recommended treatment series.

USE

PREPARATION

- Remove the charging base (6), handpiece (1), and adaptor (8) from the packaging.
- Fully charge the battery before using the appliance for the first time. This operation can take up to 6 hours.

CHARGING THE BATTERY

- Insert the handpiece (1) into the charging base (6).
- Insert the adaptor plug (8) into the charging connection (7).
- Insert the adaptor (8) into a socket:
- The orange - battery charging - indicator light (C) lights up to indicate that the appliance is charging.
- The appliance is fully charged and ready to be used when the orange - battery charging - indicator light (C) goes off and the green - battery charged - one (D) lights up.

The appliance must be charged when the orange - battery charging - indicator light (C) lights up again after approximately 25 minutes of use.

DO NOT cover the appliance and the adaptor while charging!

Note

Do not use the appliance while it is charging.

FIRST TREATMENT

The skin must be clean, dry and free from powder, gel or cream.

- Remove the handpiece (1) from the charging base (6).
- press the On/Off button (3) ① : the green indicator light (A) - appliance ready for treatment - starts blinking and then becomes steady.
- Place the treatment surface (2) on your skin: the LED will light up automatically.
- Gently press and move the handpiece (1) in small circles. The temperature of the treatment surface (2) will slowly increase.
- Perform the treatment for 5-7 minutes on every area.
- Turn the appliance off at the end of the session by pressing the on/off switch (3).

**DURING THE TREATMENT**

During treatment it is normal to experience and feel:

- Heat – you can feel a moderate heat. This is normal and should disappear approximately one hour after the session.
- Slight redness – after the treatment you can notice slight skin redness (mild erythema) over the treated area. This is normal and should fade within an hour or two.

- Fan noise – it is normal to hear the noise of a fan that is activated and deactivated during the treatment. This is a system that helps adjust the temperature during use.

AFTER THE TREATMENT

After every session, we recommend:

- Applying a small amount of moisturising cream on the treated areas, gently massaging it in circular movements.
- Cleaning the appliance and, in particular, the treatment surface (2) (see “Maintenance” chapter).
- Insert the handpiece (1) in the charging base (6) and charge the appliance (see chapter “Charging the battery”).

MAINTENANCE

CLEANING AND STORING

We recommend cleaning the appliance after every session:

- Switch the appliance off by pressing the ON/OFF (3) switch ①.
- Disconnect the adapter plug (8) from the socket.
- Leave the treatment surface (2) to cool down.
- Clean the appliance and, in particular, the treatment surface (2), using a clean, soft cloth, either dry or slightly damp.
- Store the clean appliance and all its components in a clean, dry and dust-free place at a temperature ranging between -10°C and 60°C.

ATTENTION:

- Never immerse the appliance or any of its parts in water!
- Do not clean the appliance or its components under the tap or in the dishwasher.

TROUBLESHOOTING

Problem	Possible cause	Solution
Not working	Flat rechargeable batteries	Charge the appliance
	The appliance is broken	Contact the Imetec customer service
It turns off abruptly	Flat rechargeable batteries	Charge the appliance
The temperature of the appliance increases while charging	This is normal	No action required
Orange indicator light (C) - battery charging - does not light up when the adaptor is connected to the charging connection (7) via the power cord (8)	the adaptor (8) is not plugged in the socket	Insert the adapter plug (8) into the socket
		Connect another appliance to the wall socket to see whether it works
		If the appliance still does not charge, contact the IMETEC customer service
	The adaptor (8) is not inserted properly into the socket and/or the power cord is not inserted properly into the charging connection (7)	Ensure that the adaptor (8) is inserted properly into the socket and/or that its power cord (8) is inserted properly into the charging connection (7).
	The handpiece is broken	Contact the Imetec customer service

The handpiece temperature increases during use	This is normal	No action required
The green indicator light (A) - ready for treatment - does not light up when the on/off button (3) is pressed	Flat batteries	Charge the appliance
	The batteries are charged. The green indicator light (D) lights up when the handpiece (1) is placed in the charging base (6)	Press the on/off switch (3) twice to turn the appliance off and then back on. If the green indicator light (A) - ready for treatment - does not light up again, contact the IMETEC customer service
The yellow indicator light (B) - system error - blinks	Appliance malfunction	Press the on/off switch (3) twice to turn the appliance off and then back on. If the yellow indicator light (B) - system error - lights up again, contact the IMETEC customer service
The green indicator light (A) - ready for treatment - blinks	Flat rechargeable batteries	Charge the appliance
The orange (C) - battery charging - and green (D) - battery charged - indicator lights light up alternatively	Appliance malfunction	Disconnect the adaptor (8) from the socket and press the on/off switch (3) twice to turn the appliance off and then back on. Connect the adaptor (8) to the wall socket. If the problem persists, contact the IMETEC Customer Service Centre
	Appliance malfunction	Press the on/off switch (3) twice to turn the appliance off and then back on. If the problem persists, contact the IMETEC Customer Service

DISPOSAL

The packaging is made up of recyclable material. Dispose of it in compliance with the environmental protection standards.

The appliance must be disposed of properly, in accordance with the provisions of European standard 2002/96/EC.

The recyclable materials contained in the appliance must be recovered so as to prevent environmental pollution.

For further information, contact the local disposal body or the appliance dealer.

BATTERIES DISPOSAL

This appliance contains Lithium Ion rechargeable batteries.


WARNING: Do not dispose of the batteries in domestic waste: these must be disposed of separately. The battery is to be disposed of safely.

- **The battery must be removed from the appliance before it is scrapped.**
- **Dispose of the batteries at a specific collection centre.**


REMOVING THE BATTERIES

- Remove the batteries only once they are flat and cannot be recharged
- Disassemble the batteries:


1. Use a screwdriver to loosen the screw on the lower panel of the handpiece.

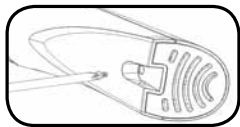
2. Lift and remove the handpiece panel.

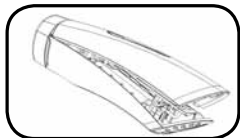
3. Push the three snaps and remove the electronic board.

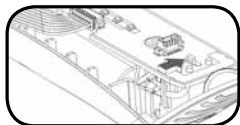4. Cut the wires that connect the rechargeable batteries to the handpiece.

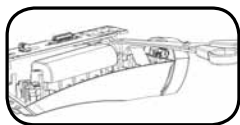5. Cut the wires that connect the electronic board to the lower handpiece panel.

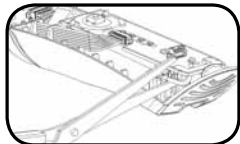Note: Do not plug the appliance back into the socket once the batteries have been removed.

ASSISTANCE AND WARRANTY

To request repairs or purchase spare parts, please contact IMETEC's authorised customer service centre on the Toll-Free Number below or via the website. The appliance is covered by the manufacturer's warranty.

For details, read the attached warranty sheet. Failure to comply with the instructions for use and maintenance provided in this manual will render the manufacturer's warranty null and void.

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