beurer

MG 156

**READ THIS MANUAL
COMPLETELY AND CAREFULLY
BEFORE USING THIS PRODUCT.**

Keep this manual in a safe location
for future reference.

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INDICATION FOR USE

This electrical massage device is intended to massage the back and thighs. Massages can be both relaxing and stimulating, and are often used in the treatment of stiff or tense muscles, pains and fatigue. This device massages with the aid of 3 vibration motors and offers a wide range of setting options. The massage seat cover has to be used when placed on a suitable seating surface. Please do not place a cover over it.

The device is only for **private use** and is not intended for medical or commercial purposes.

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1. Important safety notes

Signs and symbols

The following signs appear in the Safety Section (pages 2 to 3) and in this manual on pages 4 to 7.

READ AND UNDERSTAND THIS ENTIRE MANUAL, INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS, COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	DANGER DANGER indicates a hazardous situation which, if not avoided, will result in death or serious injury.
	WARNING WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.
	CAUTION CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.
	NOTICE NOTICE addresses practices not related to personal injury, such as product and/or property damage.

DANGER:

When using the AC-Adapter take precaution of all warnings and safety information in order to avoid hazardous situations.

- Do not use any AC-Adapter other than the one supplied with the unit.
- Always unplug this product immediately after use and make sure it is not accessible to children.
- Do not place or store product where it can fall or be pulled into a bathtub, sink, water, or other liquid. Risk of electric shock.
- Do not use this product while bathing, showering, washing dishes, or close to water sources of any kind. Risk of electric shock.
- Do not place device in or drop into water or other liquid. Do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry. Risk of electric shock.
- Do not reach for a product that has fallen into water or other liquid. If possible, safely unplug immediately.

diately. Do not touch a wet power cable or plug. Do not touch the power cable or plug with wet extremities.

- Never touch the unit, including the AC-Adapter, with wet hands while it is plugged in. This can cause an electric shock.
- Before use, ensure that there is no visible damage to the unit or accessories. Never operate the unit if it
 - has any damaged parts, e.g. damaged cord, plug or housing,
 - has been exposed to extreme humidity or to any liquids inside the housing,
 - has been dropped,
 - is not working properly.
- If the unit does not function correctly or if you notice any abnormalities, e.g. an unusual noise or smell, stop using the unit immediately. Switch off and unplug the unit immediately. When in doubt, do not use the unit and contact Beurer customer service. (Please see warranty for service contact).
- Do not open or disassemble the massage seat cover. Risk of electric shock. All disassembly and maintenance of the device must be performed by a qualified service technician. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact.
- Before you clean the unit, always switch it off, unplug it and let it cool off.
- This product should never be left unattended when plugged in.
- Do not pull the AC-Adapter from the power outlet by the cord. Always grasp the AC-Adapter housing and unplug it. Failure to do so could result in risk of electric shock, short circuiting, or fire.
- Do not carry this appliance by supply cord or use cord as handle.
- Do not pinch or bend the power cord or pull it over sharp objects. This can cause an electric shock.
- The unit and the AC-Adapter must be protected and stored away from sources of heat. (e.g. stoves, heating radiators)
- The AC-Adapter has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not modify the plug in any way. Do not attempt to manipulate this safety feature.
- Do not operate this device outdoors. Risk of electric shock.
- Do not use extension cords, multi-sockets and/or voltage converters when using the AC-Adapter.

WARNING:

- Using the unit is not a substitute for consultation with or treatment from a physician. Whenever you have any pain, illness or discomfort, always contact your physician first. If you have an emergency please call 911 immediately.
- If you start feeling any discomfort, pain, drowsiness or other change in your condition, stop using the unit immediately. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- Never use this product while in bed.
- This device is not intended to be used by persons with restricted physical (e.g. activity limitations, disabilities), sensory (e.g. insensitivity to pain) or mental disabilities (e.g. mental disorders) or persons lacking the required experience or knowledge for safe operation of the device, unless closely supervised by a person responsible for their safety. **NEVER LEAVE THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.**
- This product contains small parts that may present a choking hazard to children. Keep all small parts out of reach of children. The power cord also presents a strangulation hazard. **NEVER LEAVE CHILDREN OR THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.**
- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Use the unit only on human beings and for its intended purpose and in the manner described in these instructions for use. Any improper use can be dangerous!
- Do not operate the unit in the presence of flammable gases, oxygen and aerosol spray products. Non-compliance may result in e.g. fire, burns, explosions, etc.
- Do not use this unit while operating a vehicle or a machine.
- Before using the AC-Adapter, check the data tag sewn onto the side of the massager to ensure that the voltage and wattage indicated on the unit correspond to the voltage and wattage available.
- The massage seat cover must only be used with suitable Beurer accessories. The use of other accessories may damage the device and can impair its usage.
- Never drop or insert any object into any opening.
- This unit is not intended for commercial or clinical use, but only for private household use!
- Lay the power cord in such a way that no one can trip over it. There is a danger of injury due to e.g. stumbling, tripping, falling down.
- When cleaning, ensure that all residues are removed. Never use any substances that may be

poisonous upon contact with the skin or mucous membranes.

NOTICE:

- Protect the unit or the parts from strong impacts such as dropping the unit on the floor.
- Never open or attempt to repair the unit yourself. Failure to observe this regulation shall void the warranty.
- Ensure that no water penetrates inside the unit! Never immerse the unit in water or any other liquid. Never hold the complete device under running water to clean it. Otherwise liquid can enter it and cause damage.
- Do not put the device or accessories in a dishwasher, washing machine or tumble dryer.
- Store the device and accessories in a location protected against environmental influences, i.e. in a safe and dry location. Do not expose the device and accessories to direct sunlight, extreme hot or cold temperatures or humidity (e.g. in the bathroom).
- Protect the unit against dust, dirt and moisture. This may damage the device.
- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.
- Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.
- Portable and mobile RF communication devices can disrupt electrical medical equipment. Technical electromagnetic compatibility data is available on request from Beurer customer service. (Please see warranty for service contact)
- **NOTE:** The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

Contents

1. Important safety notes	2
⚠ DANGER	2
⚠ WARNING	3
NOTICE	3
2. Package Contents	4
3. Getting acquainted	4
4. Important information on the device	4
5. Description of device and accessories	5
6. Initial use	5
7. Operation	6
8. Cleaning	6
9. Storage and Caring for the device	7
10. Disposal	7
11. Troubleshooting	7
12. Technical specifications	8
13. Warranty and Service	8

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

The manufacturer shall not be held liable for damage or injuries caused by improper or incorrect use.

2. Package Contents

For details of the delivery scope, see description of the device and accessories on page 5.

- Massage seat cover
- Power supply unit
- These instructions for use

3. Getting acquainted

Dear valued customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air. Please read these instructions completely and carefully before using this product and keep them for later reference. Be sure to make them accessible to other users and observe the information they contain.

For product information, questions or customer service call toll free 1-800-536-0366 or email to info@beurer.com.

With kind regards,
Your Beurer team

Area of application

This device is solely intended for massaging the backs and legs of human beings. It cannot replace medical treatment. Do not use the massage device if any of the following warnings apply to you.


4. Important information on the device

⚠ WARNING:

- The massage seat cover must only be used with suitable Beurer accessories. The use of other accessories may damage the device and can impair its usage.
- Use of the massage device is prohibited in the following cases:
 - not with changes caused by disease or injury in the back region (e.g. disc problems, open wounds),
 - any health issues,
 - not on children,
 - not during pregnancy,
 - never while you sleep,
 - not on animals,
 - not during activities during which an unexpected reaction may be dangerous,
 - not while ingesting substances that reduce reaction times (e.g. painkiller medications, alcohol),
 - never longer than 20 minutes (risk of overheating) and allow it to cool off for at least 15 minutes before using it again.
- This device is not intended for use by persons (including children) with restricted physical, sensory or mental disabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety or they are given instructions on how to use the device.
- Children should be supervised to ensure that they do not use the device as a toy.
- The device has a relatively hot surface. Persons who are sensitive to heat must be careful when using the device.
- You should consult your doctor before using the massage device if:
 - you suffer from a serious illness or are recovering from an operation on the upper body, buttocks or thighs.
 - you are fitted with a cardiac pacemaker, implant or other aid,
 - Thromboses,
 - Diabetes,
 - in case of pains of unclarified origin.

The device is only intended for the purpose described in these operating instructions. The manufacturer cannot be held liable for damage occurring as a result of inappropriate or careless use.

5. Description of device and accessories

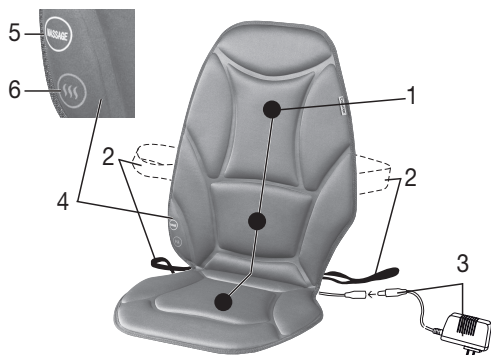- 1 3 massage motors
- 2 Fastening straps
- 3 Power supply unit with plug connector/AC-Adapter
- 4 Control panel
- 5 MESSAGE button
- 6 HEAT button

6. Initial use

WARNING:

- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Use the unit only on human beings and for its intended purpose and in the manner described in these instructions for use. Any improper use can be dangerous!

6.1 Before using the device for the first time

WARNING:

- Do not subject the device to any impact and do not drop it.
- Do not use the device if damage is evident to the device itself or any of its accessories.
- Connect the power supply unit to the seat cover using the connector.
- Route the cable safely so that there is no risk of anyone tripping over it.

- Remove all packaging materials. Ensure that the packaging materials are disposed of properly due to risk of suffocation.
- Position the seat cover on a suitable seating surface (chair, sofa, or similar). Ensure that the seat cover is supported over its entire area and has a sufficiently high backrest.
- In order to secure the seat cover, use the flexible tie straps if necessary.
- Connect the power supply unit to the seat cover using the connector.
- Plug in the power supply to the wall outlet.

6.2 Operating the device with the AC-Adapter

DANGER:

- Never touch the unit, including the AC-Adapter, with wet hands while it is plugged in. This can cause an electric shock.
- Do not pull the AC-Adapter from the power outlet by the cord. Always grasp the AC-Adapter housing and unplug it. Failure to do so could result in risk of electric shock, short circuiting, or fire.
- Do not pinch or bend the power cord or pull it over sharp objects. This can cause an electric shock.
- The AC-Adapter has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not modify the plug in any way. Do not attempt to manipulate this safety feature.
- Do not operate this device outdoors. Risk of electric shock.
- Do not use extension cords, multi-sockets and/or voltage converters when using the AC-Adapter.

WARNING:

- Before using the AC-Adapter, check the data tag sewn onto the side of the massager to ensure that the voltage and wattage indicated on the unit correspond to the voltage and wattage available.
- Lay the power cord in such a way that no one can trip over it. There is a danger of injury due to e.g. stumbling, tripping, falling down.

Plug the AC-Adapter (item 3 from picture above) only into a 120 Volt AC wall outlet only.

- Ensure that there is a power outlet near to where the device will be used.
- Insert the power supply into a suitable power outlet.
- Push the plug completely into the power outlet.
- To disconnect the massage seat cover from the power source after use, first switch off the device before removing the AC-Adapter (item 3 from picture above) from the outlet.

7. Operation

DANGER:

When using the AC-Adapter take precaution of all warnings and safety information in order to avoid hazardous situations.

- Do not use any AC-Adapter other than the one supplied with the unit.
- Always unplug this product immediately after use and make sure it is not accessible to children.
- Before use, ensure that there is no visible damage to the unit or accessories. Never operate the unit if it
 - has any damaged parts, e.g. damaged cord, plug or housing,
 - has been exposed to extreme humidity or to any liquids inside the housing,
 - has been dropped,
 - is not working properly.
- This product should never be left unattended when plugged in.
- Do not use this product while bathing, showering, washing dishes, or close to water sources of any kind. Risk of electric shock.
- Do not place device in or drop into water or other liquid. Do not allow any water to be sprayed onto the unit. The unit must be operated only when it is completely dry. Risk of electric shock.
- Do not reach for a product that has fallen into water or other liquid. If possible, safely unplug immediately. Do not touch a wet power cable or plug. Do not touch the power cable or plug with wet extremities.
- If the unit does not function correctly or if you notice any abnormalities, e.g. an unusual noise or smell, stop using the unit immediately. Switch off and unplug the unit immediately. When in doubt, do not use the unit and contact Beurer customer service. (Please see warranty for service contact).

WARNING:

- Using the unit is not a substitute for consultation with or treatment from a physician. Whenever you have any pain, illness or discomfort, always contact your physician first. If you have an emergency please call 911 immediately.
- If you start feeling any discomfort, pain, drowsiness or other change in your condition, stop using the unit immediately. For medical advice please contact a licensed physician. If you have an emergency please call 911 immediately.
- This device is not intended to be used by persons with restricted physical (e.g. activity limitations, disabilities), sensory (e.g. insensitivity to pain) or mental disabilities (e.g. mental disorders) or persons lacking the required experience or knowledge for safe operation of the device, unless closely su-

pervised by a person responsible for their safety. NEVER LEAVE THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.

- Never use this product while in bed.
- This product contains small parts that may present a choking hazard to children. Keep all small parts out of reach of children. The power cord also presents a strangulation hazard. NEVER LEAVE CHILDREN OR THOSE WHO REQUIRE CLOSE SUPERVISION UNATTENDED WITH THIS DEVICE.
- Use the unit only on human beings and for its intended purpose and in the manner described in these instructions for use. Any improper use can be dangerous!
- Do not use this unit while operating a vehicle or a machine.
- Never drop or insert any object into any opening.
- This unit is not intended for commercial or clinical use, but only for private household use!
- Lay the power cord in such a way that no one can trip over it. There is a danger of injury due to e.g. stumbling, tripping, falling down.
- Do not operate the unit in the presence of flammable gases, oxygen and aerosol spray products. Non-compliance may result in e.g. fire, burns, explosions, etc.

How to use

Before use, ensure the device and accessories are at room temperature.

- Sit comfortably on the seat cover. Take care that you sit centered and your body is in good even contact with the seat cover.
- Switch on the massage device with the MASSAGE button. The massage intensity is preset to level 1. Press the MASSAGE button again to increase the massage intensity to level 2. Press the MASSAGE button again to switch off the unit.
- The automatic switch-off function is activated by starting the massage and is preset to 20 minutes.
- Press the HEAT button to switch the heat function in the lumbar region on and off.
- Do not use the massage device for longer than 20 minutes. The massage can be finished at any time by pressing the MASSAGE button.
- After the massage is finished, disconnect the device from the power outlet.

8. Cleaning

DANGER:

- Before you clean the unit, always switch it off, unplug it and let it cool off.
- Never touch the unit, including the AC-Adapter, with wet hands while it is plugged in. This can cause an electric shock.

WARNING:

- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Do not use the device again until it is completely dry.
- When cleaning, ensure that all residues are removed. Never use any substances that may be poisonous upon contact with the skin or mucous membranes.

NOTICE:

- Ensure that no water penetrates inside the unit! Never immerse the unit in water or any other liquid. Never hold the complete device under running water to clean it. Otherwise liquid can enter it and cause damage.
- Do not put the device or accessories in a dishwasher, washing machine or tumble dryer.
- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.

How to clean the unit

The massage seat cover and accessories are intended for multiple use.

Before you clean the device, always switch it off, unplug it and let it cool off.

Then you can clean it with a cloth which is no more than lightly moistened. You can also use a mild detergent.

9. Storage and Caring for the device

DANGER:

- Always unplug this product immediately after use and make sure it is not accessible to children.
- Do not open or disassemble the massage seat cover. Risk of electric shock. All disassembly and maintenance of the device must be performed by a qualified service technician. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact.
- The unit and the AC-Adapter must be protected and stored away from sources of heat. (e.g. stoves, heating radiators).
- Do not place or store product where it can fall or be pulled into a bathtub, sink, water, or other liquid. Risk of electric shock.

WARNING:

- The massage seat cover must only be used with suitable Beurer accessories. The use of other accessories may damage the device and can impair it's usage.

NOTICE:

- Do not use any aggressive solvents, cleaning agents, detergents or any other strong chemicals to clean the device.
- Protect the unit or the parts from strong impacts such as dropping the unit on the floor.
- Never open or attempt to repair the unit yourself. Failure to observe this regulation shall void the warranty.
- Store the device and accessories in a location protected against environmental influences, i.e. in a safe and dry location. Do not expose the device and accessories to direct sunlight, extreme hot or cold temperatures or humidity (e.g. in the bathroom).
- Protect the unit against dust, dirt and moisture. This may damage the device.

Keep the massage device in a cool and dry place where it is out of the reach of children. We recommend storing the device in its original packaging. Store and transport the massage seat cover away from direct sunlight. If need in repair please contact customer service.

Avoid any contact with sharp or pointed objects which could damage the device or cover.

10. Disposal

NOTICE:

- Dispose of the device, components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.

If you have any questions, please contact the local authorities responsible for waste disposal.

11. Troubleshooting

NOTICE:

- Repair or replacement is our only responsibility and your only remedy under this written warranty.

Problem	Cause	Remedy
The device is very loud during massage.	The seat cover is not fully in contact with the base.	Ensure that the seat cover is supported over its entire area.
Massage stopped.	The preset maximum period of operation has elapsed.	Let the device cool down for at least 15 min. before using again. Switch the device on again.

Problem	Cause	Remedy
Message stopped.	The connection to the electrical outlet has been interrupted.	Check the plug connection and the contact of the AC-Adapter with the electrical outlet.
No message.	Appliance not connected to the power supply.	Make all connections and switch on the device.
	Device not switched on.	Press the Message button.
	Massage device defect. Does not work.	Switch off and unplug the unit immediately and contact Beurer customer service. (Please see warranty for service contact)

12. Technical specifications

The following symbols are used on the packaging and on the nameplate of the device and accessories.

Dimensions	38.19" x 17.72" x 1.18" (970 x 450 x 30 mm)
Device weight	21.16 oz (600 g)
Mains connection	100–240V ~; 50/60 Hz; 1.3 A

Subject to technical modifications.

Symbols Key

	WARNING! Read safety information to avoid serious injury.
	Device protection class 2
	Indoor use only!

Note:

If the device is used outside of the specifications, proper function is no longer guaranteed! We reserve the right to make technical changes to improve and further develop the product.

Note:

This device and its accessories comply with US standards and are subject to special safety measures in terms of electromagnetic tolerance. Note that portable and mobile RF communication equipment can affect this device. For detailed information, contact customer service at the address provided.

ELECTROMAGNETIC COMPATIBILITY

- The device complies with current specifications with regard to electromagnetic compatibility and is suitable for use in all premises, including those designed

for private residential purposes. The radio frequency emissions of the device are extremely low and in all probability do not cause any interference with other devices in the proximity.

- It is recommended that you do not place the device on top of or close to other electronic devices. Should you notice any interference with other electrical devices, move the device or connect it to a different socket.
- Radio equipment and TV interference may affect the operation of this device.

13. Warranty and Service

Limited Lifetime Warranty For Original Purchaser

Your Beurer Massage seat cover, Model MG 156, excluding all enclosed accessories as these require routine maintenance, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Product, Model MG 156, without charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide quotation for repair, respectively, replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; failures of the Massage seat cover, tubing, or accessories; parts or cleaning solutions; the attachment of unauthorized accessories; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss

of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER.

Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
Hallandale Beach, FL 33009
www.beurer.com

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