- . Manuale di istruzioni
PRESSOTERAPIA SEQUENZIALE
- . Instruction manual
SEQUENTIAL COMPRESSION THERAPY
- . Manual de instrucciones
PRESOTERAPIA SECUENCIAL
- . Manual de instruções
PRESSOTERAPIA SEQUENCIAL

Moretti S.p.A.
Via Bruxelles 3 - Melegnano
52022 Cavriglia (Arezzo)
Telefono +39 055 96 21 11

Fax +39 055 96 21 200
www.morettispa.com
info@morettispa.com

LA PRESSOTERAPIA

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Via _____ Località _____

COMPRESSION THERAPY

The sequential compression therapy system has been designed for lymphatic drainage. The sections in the sleeves are made of partially overlap so as to eliminate any gap between one another and reduce the possibility of stagnant fluid. The sections in the sleeves inflate and deflate sequentially. Compression starts in the distal chamber and moves on to the proximal one, so favouring fluid drainage. The pressure can be adjusted by means of a knob situated on the front of the pump. Two systems are available:

ADVANCE 1000 with a 3-outlet pump for full leg or full arm sleeves;

ADVANCE 2000 with a 2-outlet pump for half leg or half arm sleeves.

Both systems are lightweight and easy to operate.

INSTALLATION

BEFORE USE

1. Before use make sure that the feeding cable is not damaged or crushed.
2. Supply power to the pump and turn it on without connecting the tubes. Check that air comes out of the outlets of the tubes (caution: the air comes out alternately from the outlets which can be in 2 or 3 according to the model).
3. Check that the connecting tubes are not damaged or crushed.
4. If necessary, replace the damaged parts only with original ones.

INSTRUCTIONS

TURN THE APPLIANCE OFF

1. Place the pump on a flat surface or hang it where you can by means of the hooks.
2. Connect the tubes to the pump and to the arm/leg sleeve.
3. Check that the tube connecting points on the sleeve are not damaged.
4. Turn the pressure-adjusting knob to the minimum.
5. Put the sleeve on and fasten the zipper completely.
6. Switch the pump on. The green light indicates that the pump is working. Do not put on or take off the sleeve while the pump is working.
7. Adjust pressure turning the knob on the front of the appliance clockwise to the value indicated by your physician according to the following chart.

KNOB	1	2	3	4	5	6	7	8
mmHg	30	38	47	55	64	72	81	90

The chart indicates the millimetres of mercury as units of measure and pressure and the corresponding numbers on the pressure-adjusting knob.

TROUBLESHOOTING

IN CASE OF DYSFUNCTIONING CHECK THE FOLLOWING POINTS.

1. Check that the pump is connected to mains and that the green switch is on and lit.
2. Check that the tubes are properly connected to the sleeve and to the pump.
3. The sleeve does not inflate immediately. You have to wait for a few minutes before noticing the cyclical working of the pump.

CONTRAINDICATIONS

COMPRESSION THERAPY IS CONTRAINDICATED FOR PEOPLE SUFFERING FROM THE FOLLOWING PATHOLOGIES: Cardiopathy, suspected or ascertained thrombosis, arteriosclerosis or ischemia, dermatitis, gangrene, local infections.

CLEANING INSTRUCTIONS

Before cleaning the appliance, make sure it is not connected to the electric network. Do not clean the appliance while you are using it. To clean the pump, only use water and deterging soap. Avoid using solvents. You can wash the sleeves with liquid detergents in water not exceeding 50°C. Do not scrub strongly and never immerge the tube connectors into water. Do not dry-clean or iron. Sterilization can be carried out only with gas at a temperature never exceeding 50°C.


TECHNICAL FEATURES

	ADVANCE 1000	ADVANCE 2000
Power:	220V/50Hz	220V/50Hz
Pressure:	30 - 90mmHg	30 - 90mmHg
Air output:	5 liters	5 liters
Cycle time:	3/6 minutes	3/6 minutes
Feeding cable.:	2,5m	2,5m
Outlets:	3	2
Dimensions:	24x13x10cm	24x13x10cm

	LTM520	LTM510
Type:	Full leg	Full arm
Material:	Nylon+ PU	Nylon + PU
Dimensions:	95x35cm	77,5x30cm

	LTM515	LTM525
Type:	Half leg	Half arm
Material:	Nylon + PU	Nylon + PU
Dimensions:	62,5x35cm	52,5x27,5cm

INSTRUCTION MANUAL FOR SEQUENTIAL COMPRESSION THERAPY, TWO LEGS WITH A SEQUENTIAL PUMP

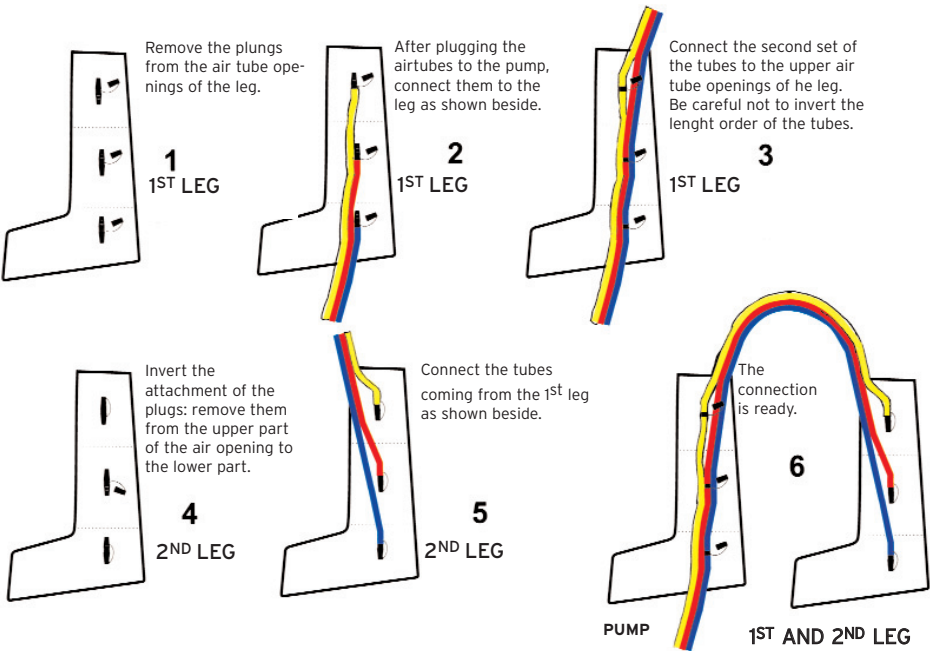WARRANTY

This device has **2-YEAR** warranty from the purchase date which is certified by the receipt or invoice. The distributor reserves the right either to repair or replace necessary parts to correct defects in the materials or workmanship without any charge. The warranty does not apply to damage resulting from failure to follow the operating instructions, accidents, abuse, alteration or disassembly by unauthorized individuals. In the case you need to get your device repaired and you fulfil the conditions of warranty, send the device to the nearest retailer together with the purchase receipt.

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