IC0908

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MASSAGE ACCESSORY
ACCESSOIRE DE MASSAGE

Thank you for purchasing this product. Please read this instruction manual carefully before operating it.

Please pay attention to the important safety warnings and keep this manual for future use.

Note: EBI reserves the right to revise the design and description of this product without notice.

CONTENTS

1. Safety warnings.....	2
2. Environment.....	2
3. Restrictions.....	2
4. Maintenance.....	3
5. Troubleshooting.....	3
6. Features and benefits.....	4
7. Components.....	5
8. Applications.....	6
9. Operation.....	7
10. Product specifications.....	7
11. Warranty.....	8
FRANÇAIS.....	9

1. SAFETY WARNINGS

- This product is not recommended for children.
- Disconnect the massage unit before cleaning.
- Please operate this product according to the instructions in this manual.
- This product is designed for indoor use only.
- Recommended use time is 30 minutes or less.
- Please avoid use if the massage unit is damaged and call your iComfort technician or dealer.
- Do not use wet hands to pull the plug.
- Do not spill liquids onto this product to avoid electric shock or cause damage to this product.
- In the event of a power failure disconnect the massager.

2. ENVIRONMENT

The iComfort shiatsu neck massager should be used only in well ventilated, normally dry and dust free environments.

3. RESTRICTIONS

Stop immediately the use of this massage unit if you feel unwell and consult your doctor.

4. MAINTENANCE

- This product should only be maintained by iComfort certified technicians otherwise, EBI reserves the right to void the massager warranty.
- Make sure the wall outlet is in good working order before use.
- If this massage unit will not be used for an extended period of time, store in a dry and dustless environment.
- Avoid product exposure to high temperatures and direct sunlight.
- Clean this product with a dry cloth. Do not use thinner, benzene or alcohol.
- We recommend to use the iComfort shiatsu neck massager intermittently and not for long time periods.

5. TROUBLESHOOTING

- It is normal for the motor to make noise while in use.
- If the function keys doesn't work normally, please check that the unit AC power cord are properly connected to the wall AC outlet and the ON/OFF switch of the unit is turned on.
- When the automatic program mode working time is over, the foot massager will power-off automatically.
- After several use, the shiatsu neck massager might power-off because of overheat. It can be used again after half an hour of cooling time.


6. FEATURES AND BENEFITS

Living the cubicle life? Spend your days staring at a monitor? Enjoy a true Shiatsu neck massage in the comfort of your own home. The adjustable kneading Shiatsu iComfort® ic0908 Neck Massager simulates the hands of a professional massage therapist. The mechanical kneading thumbs provide a rhythmic action that helps to provide a soothing relaxation feeling. The new iComfort® ic0908 Neck Massager even includes Infrared thermotherapy for an even more incredible massage experience. This massage therapist won't ever get tired!

- 2-Way Direction Selector.
- Deep Kneading Action.
- Simulates the hands of a professional massage therapist.
- Shiatsu massage action
- Can be used on the neck, arm or foot. You can lie back, relax and enjoy the benefits of a soothing massage.
- 2-way Direction Selector allows you to adjust kneading action inwards or outwards.
- Angled design for greater comfort lying down.
- Large handles make it easy to hold when you want to target specific areas.

7. COMPONENTS

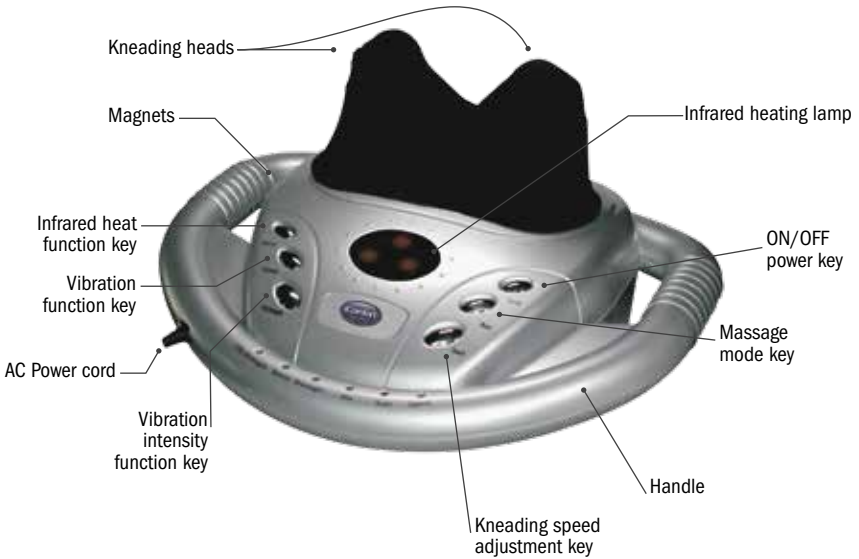8. APPLICATIONS

In addition to the massage of the neck, the iComfort infrared shiatsu neck massager may also be used to massage the following body parts:

NECK

SHOULDER

BACK

LUMBAR

ARM

SOLE OF FOOT

9. OPERATION

- Press the ON/OFF power key to start the massaging function; press it again to shut the massager off.
- **Massage mode key:** Press the Massage mode key to choose between three kneading massage actions and an automatic kneading massage mode. In automatic mode the speed and intensity cannot be adjusted. However, the infrared and vibration massage functions can be activated or deactivated.
- **Kneading speed adjustment key:** Pressing the kneading speed adjustment key multiple times will let you choose between fast, medium or slow massage speeds indicated by the red, orange and green LED light.
Note: When you power on the unit, it will automatically massage at medium speed.
- **Infrared heat function key:** Press the infrared heat function key once to start, press again to activate the infrared heating function intermittently. Press again to deactivate completely.
- **Vibration function key:** Press the vibration function key to activate the massage by vibration. Press the key again and the vibration will alternate left and right. Press again and the vibration feature will work intermittently. Press the key a fourth time and the feature will deactivate completely.
- **Vibration intensity function key:** Once the vibration feature is activated, this key will let you adjust vibration intensity between strong, medium and weak. indicated by the red, orange and green LED light.
Note: When you power on the massage by vibration feature, the unit will automatically massage by vibration at medium intensity.

10. PRODUCT SPECIFICATIONS

Rated Voltage: AC110V 60Hz

Rated power: 35 Watt

11. WARRANTY

LIMITED ONE (1) YEAR WARRANTY

This product is guaranteed against all manufactured faults and defects for a period of one (1) year from the date of purchase.

For repairs or product defects, please call our customer service department at : 1-866-967-7333. Please have your proof of purchase on hand. If we are not able to replace your defective part, we will send you a service technician for repairs.

Important notes: Any damages caused by an accident or an abusive usage is NOT covered under warranty. We also decline all responsibility concerning indirect or consecutive damages. In certain Canadian provinces, it is not permitted to limit indirect or consecutive damages. Therefore it is possible these limitations are not applicable to you. This warranty gives you specific legal rights. You may also have other legal rights which vary from one province to another.

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Pictures are for reference only, colors and models may vary without notice. Made in China. iComfort™ is a trademark of EBI Inc. Photographies pour référence seulement. Les couleurs et modèles sont sujet à varier sans préavis. Fabriqué en Chine iComfortMC est une marque de commerce de EBI inc