www.LaserComb.net

Patented in the USA. Patent Pending in 104 Countries.
HairMax LaserComb is a Registered Trademark of Lexington International, LLC.

AM10020-b
0105-5000

USER MANUAL

*Proudly Patented and Manufactured in the USA by
Lexington International, LLC*

The logo features the word "Lexington" in a simple, spaced-out, sans-serif font. Below it, "HairMAX" is written in a large, bold, italicized serif font, with the "X" being particularly large and stylized. Underneath "HairMAX" is the word "LaserComb" in a smaller, italicized serif font. A large, elegant, black swoosh arches over the text, starting from the left and ending on the right, with a small "TM" trademark symbol at its peak.

L e x i n g t o n
HairMAX
LaserComb™

HairMax LaserComb™

A Typical User Experience

As we all know, skin cells 'like' light. Your hair is no different. Both men and women can benefit from the laser light provided by the HairMax LaserComb.

Individual results vary, but within the first 12 to 16 weeks, you should notice early general improvements to the appearance of your hair. These improvements should include shinier and thicker looking hair.

Then, over the next few months, you may notice your hair gradually becoming fuller, stronger and denser looking. Over time, all users - men and women - usually notice progressively increasing benefits as they continue to use the HairMax LaserComb. Conversely, once the usage stops, our experience is that hair gradually regresses back to its original state.

Like brushing your teeth to maintain healthy teeth & gums, regular LaserComb usage is recommended. As little as three times a week for 10 to 15 minutes a session is sufficient to achieve results over time.

***Please visit
www.LaserComb.net
for more information***

Copyright 2005 HairMax LaserComb, LLC. All rights reserved.

HAIRMAX LASERCOMB



INSTRUCTIONS FOR USE

The LaserComb is a very effective product if used correctly. It is manufactured using precision components and should last a lifetime under normal usage. Please take time to read this entire instruction manual before operating the LaserComb.

Recommended usage of the LaserComb is 10 to 15 minutes per treatment, three times a week.

Studies conducted during 18 years of development of the LaserComb indicate that there is a maximum rate of return when the LaserComb is used between 10 and 15 minutes three times a week. While there aren't any negative side effects from overuse, the best results are consistently obtained by using the LaserComb according to directions.

It is recommended that the LaserComb be used on a clean head after washing. Your hair can be wet, towel dried or dry. The LaserComb should be used before applying any other topical treatments, gels, sprays or concealing products.

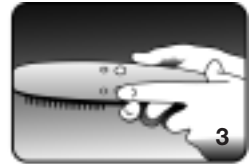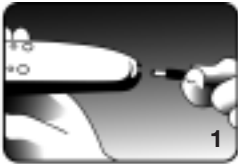1. Insert the power attachment into the handle end of the LaserComb.
2. Plug the power adaptor into any convenient wall electrical outlet.
3. Turn on the LaserComb by first switching on the 'Power' button followed by the 'Laser' button. Indicator lights for 'Power' and 'Laser' will show that the device is on. You can clearly see the laser light when the device is operating.
4. Place the LaserComb flat on your scalp using the two rows of teeth.

5. Move the LaserComb slowly using a 'glide' method by placing the LaserComb on a spot and leaving it there for 4 seconds then moving it 1/2 inch to the next spot. (The LaserComb features an audible beep tone every 4 seconds which can aid you in timing the movements). Move the LaserComb against the direction of the hair growth - front to back on the first pass (illustration 5) and bottom to top on the second pass (illustration 6). **IMPORTANT:** Do not lift the device off of the head when moving every 4 seconds. The device features a pair of comb teeth that part the hair so that the light can more effectively penetrate the desired areas.
6. Make two complete passes with the LaserComb over the entire head— first from front to back and then a second pass from bottom to top. With long hair, use an ordinary comb or your hand to lift the hair out of the way of the LaserComb.
7. After treatment turn off the power button, remove the wall plug and return the LaserComb and power supply to the case for storage.

Please Note:

- The total time required for treatment with the LaserComb should take between 10 to 15 minutes.
- It is important to move the LaserComb very slowly. Remember, you are bathing the hair with laser light. If you move it too fast, the treatment will not be as effective.

Audible Beep Tone for Proper Usage


The HairMax LaserComb is a very effective treatment for thicker, fuller, healthier looking hair when used correctly. In order to mark the time period between movements of the LaserComb, an audible tone has been added every 4 seconds. Now users can have an audible indication of when to move the LaserComb to the next area of treatment.

Instructions for Turning Beep Tone On and Off

To turn on the LaserComb and the beeping tone, push the 'Power' button (1) followed by the 'Laser' button (2). An audible beeping tone will be heard every 4 seconds with this normal activation of the LaserComb.

The tone is added to help you determine the time period between movements. Hold the LaserComb on one spot for 4 seconds, then move it approximately 1/2 inch and hold it still on the next spot. Move the LaserComb across the scalp in this fashion until you have covered the entire head twice.

If you prefer to use the LaserComb without the beep tone, press the 'Laser' button (2) a second time. You will hear a beep response and the audible tone will stop. To re-activate the tone, simply press the Laser button (2) again.

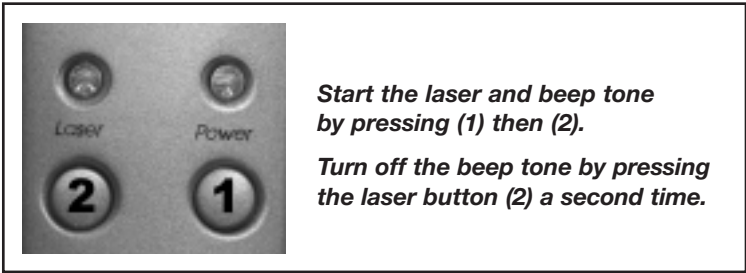
WARNINGS AND PRECAUTIONS

LASER SAFETY *Warnings and Cautions:*

1. Caution - use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.
2. Do not stare directly at the laser light or its reflection off a mirror-like surface (as this can act like a direct beam on the eye).
3. Do not shine the laser light in or around your eyes.
4. Keep out of reach of children. Do not allow minors to use lasers.
5. Never shine a laser at anyone.
6. The use of magnifying optical instruments with this product can cause eye danger if the laser is shined directly into the eye.
7. Only use directions, cleaning methods and operating procedures specified in this manual.

The following laser safety label can be found on this product:

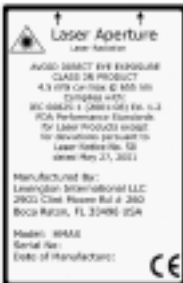Note: This appliance has no user-serviceable parts.

Location: On the underside of the handle.

ELECTRICAL SAFETY

Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

DANGER: Keep away from water. As with most electrical appliances, contact with water may result in an electrical shock.

To reduce the risk of hazard by electric shock:

- Always unplug the appliance from the power outlet immediately after using.
- Do not use while bathing.
- Do not place or store this appliance where it can fall or be pulled into a tub or sink.
- Do not expose this appliance to running water, drop into water or any other liquid.
- If the appliance falls into water, unplug it immediately from the power outlet. Do NOT reach into or touch the water.
- If using this appliance near water, the use of an electrical outlet protected by a Ground Fault Circuit Interrupter is recommended.

WARNING: To reduce the risk of burns, electrocution, fire or injury:

This appliance should never be left unattended when plugged into a power outlet.

Care should be taken when operating this appliance near children. Children are naturally curious and this appliance is not a toy and should not be played with by minors.

Use this appliance only for its intended use as described in this manual.

Never operate this appliance if it has a damaged cord or plug. If it is not working properly, or if it has been dropped, damaged or immersed in water, return the appliance to Lexington International, LLC for examination, service or replacement.

Keep the cord away from heated surfaces.

NOTICE: If any malfunction occurs, do not attempt to repair it yourself. This appliance has no user-serviceable parts.

POWER SUPPLY SPECIFICATIONS:

Output from transformer: 9V, minimum 500mA. Plug: 2.1mm x 5.5mm.

The power transformer supplied with the LaserComb is designed for use both in the United States and Internationally. This power supply will automatically adapt to all power supplies from 110 volt/60 cycle (US) to 240 volt/50 cycle. If outlets in your country do not conform to the plugs or adaptor supplied, an inexpensive adaptor can be obtained locally.

MAINTENANCE INSTRUCTIONS

Always unplug this appliance from the power outlet while cleaning it.

Warning! Do NOT use acetone or any other solvents on any part of this appliance. Use of acetone or other solvents to clean the window can immediately block most of the laser light from coming through.

To clean the laser window (from where the laser light emanates), gently clean with a soft, slightly moistened cloth.

If the cord becomes twisted, unplug and straighten it prior to use.


Do not put strain on cord if you need to untangle or straighten it.

Cleaning comb teeth

To keep the teeth of the combs clean, gently remove any hairs or dust with a slightly moistened towel.

Multi-user hygiene

For use by more than one person, a second set of combs have been included in the carrying case. Gently pull the combs down to install the second set of combs.

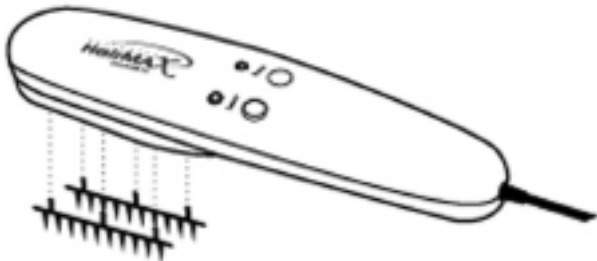For a more thorough cleaning, gently pull out hair combs and replace them using a little gentle manipulation.

SAVE THESE INSTRUCTIONS

WARRANTY

The HairMax LaserComb is warranted by Lexington International LLC ("Manufacturer") against defects in materials or workmanship for 12 months from date of purchase. If this product is determined to be defective, Lexington will repair or replace it, at its option, at no charge during the warranty period.

After the warranty period, the user must pay for all parts and labor charges. If the product is repaired or replaced, either within or outside the warranty period, Lexington reserves the right to provide equivalent parts in different colors to your original product.

This warranty does not cover damages caused by acts of God, misuse, negligence, accident, modification of, or to any part of the product. This warranty does not cover damage due to improper operation or maintenance, connection to improper voltage supply, or attempted repair by anyone other than Lexington. This warranty is void if product is tampered with.

To obtain warranty service, you must deliver the product freight prepaid, in either its original packaging or packaging affording an equal degree of protection, together with your name and address, to:

**Lexington International LLC
2650 N. Military Trail, Ste. 360
Boca Raton, FL 33431, USA**

Limitation of Liability; Indemnity: Except as expressly set forth above, Lexington International LLC shall not be liable for any special, indirect, or consequential damages that are directly or indirectly related to the use of, or the inability to use, the HairMax LaserComb, even if Lexington International LLC has been advised of the possibility of such damages whether in an action of contract, negligence or other tortious action. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

Governing Law: The usage of the HairMax LaserComb and the interpretation of its terms shall be governed by and construed in accordance with the laws of the State of Florida applicable to agreements negotiated, executed and performed wholly and solely within the State of Florida. The parties agree that the State of Florida's conflicts of laws rules shall not apply to this Agreement or any transaction performed by or under this Agreement.

Jurisdiction and Venue: The parties irrevocably submit and consent to the sole and exclusive jurisdiction and venue of the Florida state courts in and for Palm Beach County, Florida, U.S.A. and the Federal Courts in and for the Southern District of Florida, U.S.A. The parties agree not to raise the defense of forum non conveniens.

Authorized Representative in Europe (Regulatory Affairs Only): Emergo Europe, Molenstraat 15, 2513 BH The Hague, The Netherlands. Tel: (31) (0) 70 345-8570 Fax: (31) (0) 70-345-7299

Congratulations on becoming a HairMax LaserComb user!

For free customer support information providing useful insights into hair care, visit www.LaserComb.net

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