
Oral-B®

AdvancePower

D 9513 K Kids

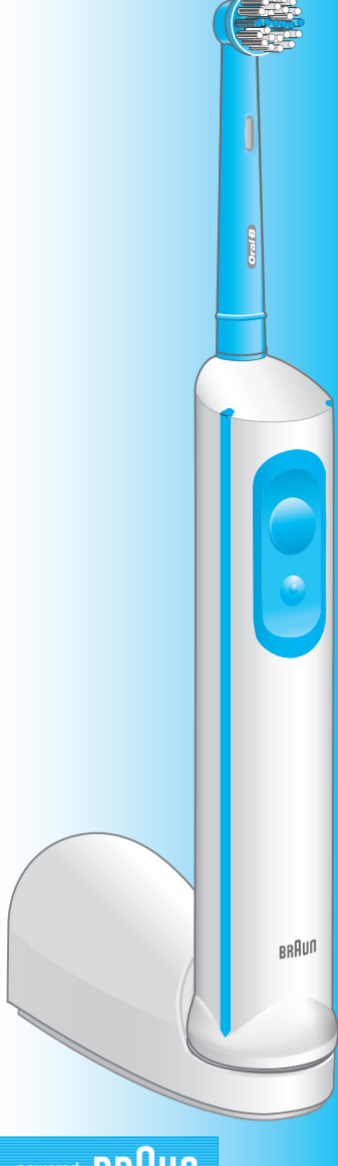powered by **BRAUN**

Type 4728

D	A	00 800 27 28 64 63
CH		08 44 - 88 40 10
UK		0 800 416 434
IRL		1 800 509 448
F		0810 365 855
B		0800 14 592
L		4 00 50 51
ES		901 11 61 84
PT		808 20 00 33
I		800 124 600
NL		0800-4 45 53 88
PL		0 801 127 286 0 801 1 BRAUN
CZ		221 804 335
SK		02/5710 1135
H		06-1-451-1256

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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English

Dear parents,
the Oral-B AdvancePower Kids has been developed with leading dental professionals to effectively remove plaque from children's teeth. We hope your child will enjoy his/her new Oral-B appliance.

This appliance is suitable for children 3 years or older. Please read the use instructions carefully and completely before showing your child how to use the appliance. We recommend that parents supervise their children when they start using the appliance in order to help them become familiar with it.

Important

- Periodically check the cord for damage. A damaged or non-functioning unit shall no longer be used. Take the charging unit (E) to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are given supervision by a person responsible for their safety.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.

Description

- A Brush head
- B On switch
- C Off switch
- D Handle
- E Charging unit

Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

Connecting and charging

The appliance is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. Handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approx. 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Using the toothbrush

Use the following steps to show and explain to your child how to use the appliance properly:

- Put a small amount of toothpaste onto the brush head. (We recommend fluoride toothpaste.)
- First guide the brush head to the teeth, only then turn on the handle with switch (B).
- Start with the lower teeth: Slowly guide the brush head from tooth to tooth;
 - first the outsides,
 - then the insides,
 - finally the chewing surfaces.Do not scrub or press too hard. The Oral-B AdvancePower Kids automatically provides the proper cleaning action.
- After one minute of operation, an incorporated timer signals by a short stutter. This tells your child to move on to the upper teeth and clean them – again tooth by tooth – until the next timer signal is given.
- Turn off the handle with switch (C) while the brush head is still in the mouth.

We recommend that you replace the brush head every three months.

The brush head is provided with a row of blue INDICATOR[®] bristles (1a).

With proper brushing twice a day for two minutes, the colour will disappear half-way down these bristles in a 3-month period (1b), provided toothpaste is used. If bristles begin to splay before the colour starts to recede, this indicates that your child exerts too much pressure on teeth and gums (1c).

Replacement brush heads EB10

Available at your retailer or at Oral-B Braun Service Centres.

Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on (2). Then switch off the handle and take off the brush head. Rinse both parts separately under running water (3) and wipe them dry before putting the handle back onto the charging unit.

Subject to change without notice.

This product conforms to the European

Directives EMC 2004/108/EC and Low

Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

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