Oral-B®

Vitality™

D 12513

D 12013

powered by **BRAUN**








Type 3709

- GB** If you should require further assistance, please call our Consumer Service line.
0 800 73 11 792 (free call)
- IRL** **Helpline**
1 800 509 448 (free call)
- F** Appelez le Service Consommateurs Oral-B :
0 810 365 855 (numéro azur)
- B** **02-711 92 11**
- LU** Vous avez des questions sur ce produit? Appelez:
4 00 50 51
- CZ** **Zelená linka**
00420-800-11-33-22
- SK** **Zelená linka**
00420-800-11-33-22
- H** **Információs vonal**
+36 (1) 801-3800
- TR** **Danışma Hattı**
0212-473-75-85
- RUS** **Телефонная помощь**
8-800-200-1444 (звонки бесплатны)
- UA** **Телефон допомоги**
38044 417 24 15

Internet:
www.oralb.com
www.braun.com

3-709-110/00/I-06/M
GB/F/PL/CZ/SK/H/TR/RUS/UA/Arab
Printed in Germany

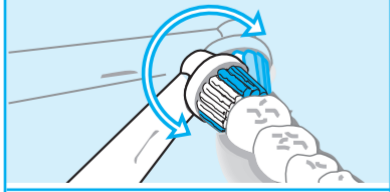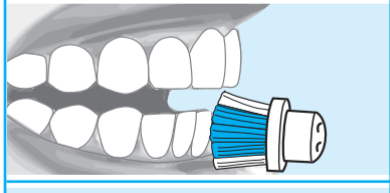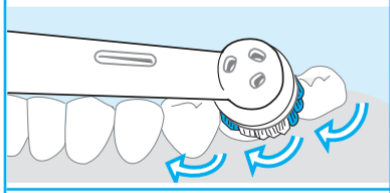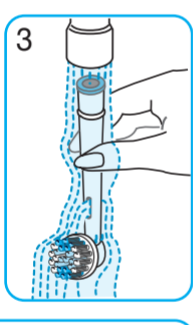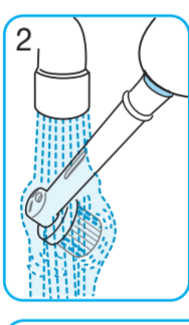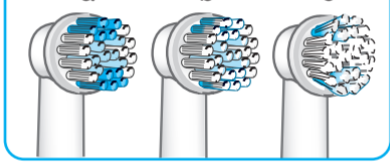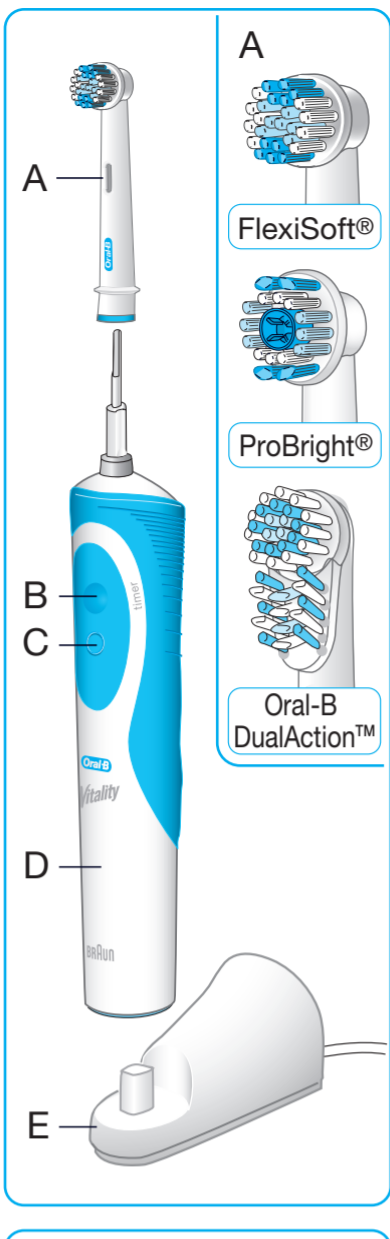English

Please read the use instructions carefully before first use.

Important

- Periodically check the cord for damage. A damaged or non-functioning unit shall no longer be used. Take the charging unit (E) to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- If the product is dropped, the brushhead should be replaced before the next use even if no damage is visible.

Description

- A Brushhead (depending on model) – interchangeable
- B On switch
- C Off switch
- D Handle
- E Charging unit

Charging

- Plug the charging unit (E) into an electrical outlet. Place the handle on the charging unit. Handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

Using the toothbrush

- Wet brushhead and apply toothpaste.
- Place brushhead on teeth, then switch on.
- Slowly guide the brushhead from tooth to tooth, first the outsides, then the insides, finally the chewing surfaces. Do not press too hard or scrub.
- After two minutes (brushing time recommended by dentists) the timer (D 12513 model only) will signal by a short stutter.
- During the first days of use, gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.
- Replace brushhead every 3 months, when bristles are worn (1c) or when colour on Indicator bristles (1b) has faded half-way down the bristles.

Cleaning

After use, rinse the brushhead under running water (2). Remove the brushhead and rinse brushhead and handle separately (3), then dry.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

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English

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brushheads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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Bahrain

Yaquby Stores
 Government road,
 Bab ALBahrein
 P.O. Box: 158, Manama
 ☎ 02-28887

Barbados

Dacosta Mannings Inc.
 Household & Furniture
 Pier Head
 Division Dacosta
 Mannings Building
 P.O. Box: 176
 Bridgetown
 ☎ 0431-8700

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Tehran Bouran Company
 145 Ghaem Magham
 Farahani
 Tehran
 ☎ 18312766

Ireland (Republic of)

Gavin's Electronics,
 83-84, Lower Camden
 Street,
 Dublin 2,
 ☎ 1800 509 448

Israel

S. Schestowitz Ltd.
 8, Shacham Str.
 49517 Tel-Aviv
 ☎ 1 800 335 959

Italia

Servizio Consumatori
 Braun
 Gillette Group Italy S.p.A.,
 Via G. B. Pirelli, 18,
 20124 Milano,
 ☎ 02 / 6678623

Jordan

Interbrands
 Wadi Al-Sir, Adwiat
 alhekma
 St. Mekadabi bld.
 Amman
 ☎ +692 582 75 67

Korea

Gillette Korea Ltd.
 144-27 Samsung-dong,
 Kangnam-ku,
 Seoul, Korea,
 ☎ 080-920-6000

Kuwait

Union Trading Company,
 Braun Service Center,
 P. O. Box 28 Safat,
 Safat Code 13001,
 Kuwait,
 ☎ 04 83 32 74

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Philippines

Gillette Philippines Inc.,
Corporate Corner
Commerce Avenues
20/F Tower 1,
IL Corporate Centre
1770 Muntinlupa city
☎ 02-771071 02-06/-16

Poland

Gillette Group Poland
Sp.Z o.o.,
ul. Domaniewska 41,
02-672 Warszawa,
☎ 22 548 88 88

Portugal

Grupo Gillette Portugal,
Lda.,
Braun Service,
Rua Tomás da Fonseca,
Torre G - 2º A,
1600-209 Lisboa,
☎ 808 2 000 33

Réunion

Dindar Confort,
Rte du Gymnase,
Boite Postale 278,
97940 St. Clotilde,
☎ 026 92 32 03

Romania

Gillette Romania srl.
Calea Floreasca
nr. 133-137 et 1, sect 1,
714011 Bucuresti
☎ 01-2319656

Russia

RTC Sovinservice,
Rusakovskaya 7,
107140 Moscow,
☎ (095) 264 41 61

Saudi Arabia

AL Naghi company
AL Madinah road
opposite to Fetihi center,
Al Forsan
P.O. Box: 269
21411 Jeddah
☎ 02 - 651 8670

Schweiz/Suisse/Svizzera

Telion AG,
Rütistrasse 26,
8952 Schlieren,
☎ 0844-88 40 10

Serbia/Montenegro

BG Elektronik,
Bulevar Kralja
Aleksandra 34,
11000 Beograd,
☎ 011 3240 030

Singapore (Republic of)

Beste (S) Pte. Ltd.,
No. 6 Tagore Drive,
03-04 Tagore Industrial
Building,
Singapore 787623,
☎ (65) 6552 2422

Slovakia

Techno Servis Bratislava
Trhová ul.
841 02 Bratislava IV
☎ 02 / 64 46 36 43

Slovenia

Lomas d.o.o.
Trzaska cesta 132
1000 Ljubljana
☎ 01-4233620

South Africa

Fixnet CC.
159 Queen Street,
Kensington South
P.O. Box: 751770
Johannesburg
☎ 6156765/6

St. Maarten

Rupchand Sons n.v.
(ram's),
P. O. Box 79, Philipsburg,
Netherlands Antilles,
☎ 52 29 31

St. Thomas

Boolchand's Ltd.,
31 Main Street,
P. O. Box 5667,
00803 St. Thomas,
US Virgin Islands,
☎ 340 776 0302

Suomi

Gillette Group Finland Oy,
P.O. Box 9,
Niittykatu 8, PL 9,
02200 Espoo,
☎ 09-45 28 71

Sverige

Gillette Group Sverige AB,
Dept. Sweden,
Stockholm Gillette
Räsundavägen 12,
Box 702, 16927 Solna,
☎ 020-21 33 21

Syria

Ahmed Hadaya Company
Shallan Houboubi
street, 1st.
P.O. Box: 33403
Damascus
☎ 113346976

Taiwan

Gillette Taiwan Inc.
Taiwan Branch office
7Fl, No. 99, Fu-Hsin North
Road, 105,
(105) Taipei
☎ 02-27197366917

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Yemen (Republic of)

Saba Stores for Trading,
26th September Street,
P. O. Box 5278,
Taiz,
☎ 4-25 23 88

For latest updates, please
refer to www.Braun.com.

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Gillette de Venezuela S.A.,