Oral-B®

OxyJet® 3D Center

OC 15 525 X

powered by **BRAUN**

Type 4715

English	4
Français	12
Español	20

Thank you for purchasing an
Oral-B OxyJet 3D Center.
If you have any questions, please call:
US residents 1-800-984-9668
Canadian residents 1 800 387 6657

Use instructions	6
Mode d'emploi	14
Instrucciones de uso	22

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USA


Warranty / Service	11
24 Month Performance Guarantee	27

Canada

Warranty / Service	19
Garantie / Service	19
24 Month Performance Guarantee	
Garantie de Rendement de 24 mois	28

Mexico

Garantía	26
Garantía de reembolso	29

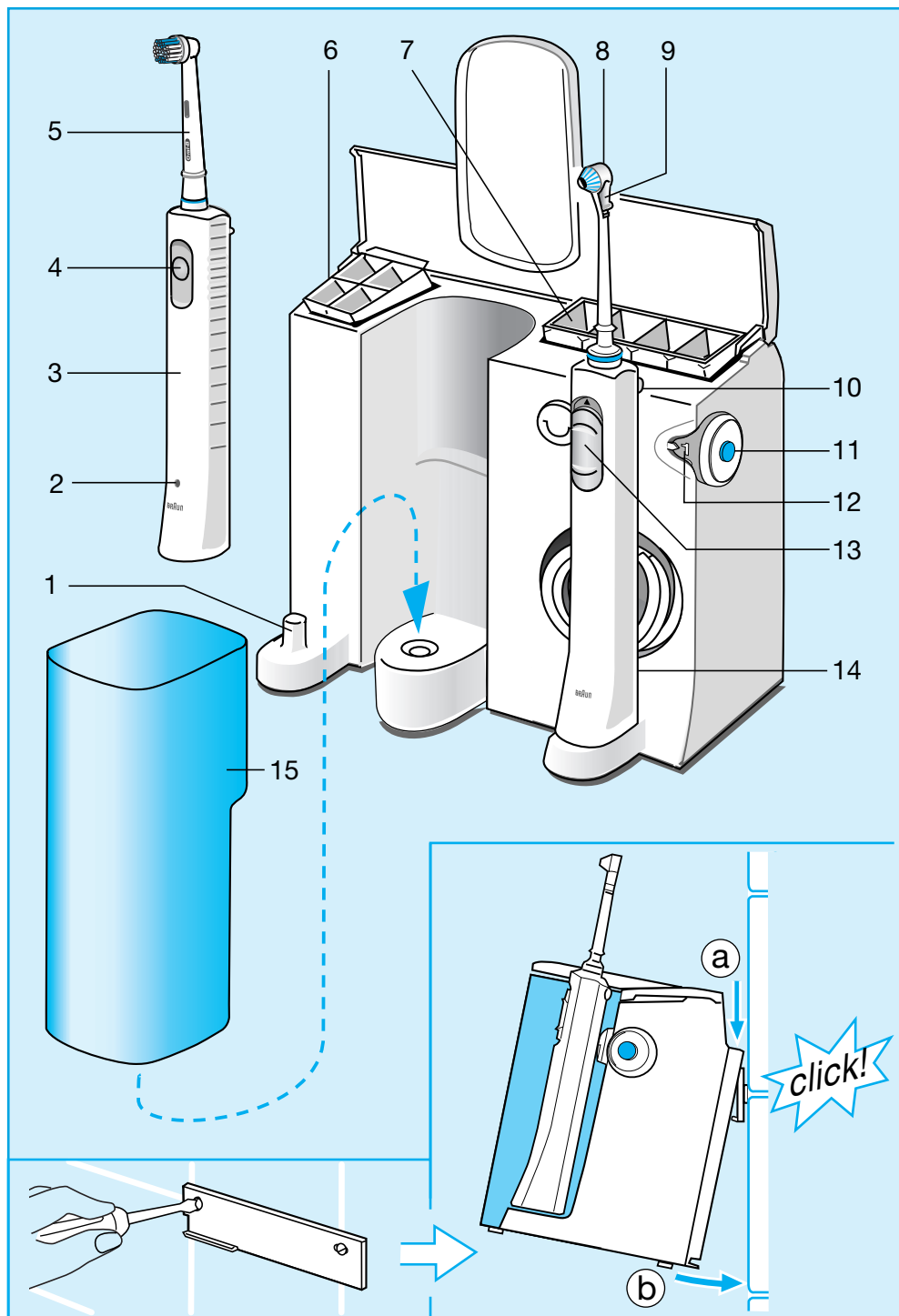IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFOR USING.

DANGER

To reduce the risk of electrocution:

1. Always unplug this appliance immediately after using or recharging batteries.
2. Do not use while bathing.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for an appliance that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. An appliance should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Remove the water container before filling. Fill with lukewarm water or appropriate oral rinse products only. Do not overfill.
7. Never drop or insert any object into any opening or hose.
8. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

**SAVE THESE
INSTRUCTIONS**

The Oral-B OxyJet® 3D Center has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth and to attack plaque bacteria.

The **Oral-B 3D pulsating toothbrush** adds a gentle pulsating movement to its oscillation to create a deep-cleaning 3D brushing action. The toothbrush should be used twice a day for at least two minutes, the minimum brushing time recommended by dental professionals.

The **Oral-B OxyJet®** irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria.

At the same time, it removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The plug should always be accessible.

Periodically check the cord for damage. Should this be the case, take the appliance to a Braun Service Center. A damaged or non-functioning unit should no longer be used.

Description

Toothbrush

- ① Charging base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ Switch
- ⑤ Brushhead
- ⑥ Brushhead compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet switch
- ⑩ Jet release button
- ⑪ On/off switch
- ⑫ Water pressure control dial
- ⑬ Sliding switch (water stop and go)
- ⑭ Irrigator handle
- ⑮ Water container (600 ml/20 fl. oz.)

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using the 3D pulsating toothbrush

Connecting and charging

Plug the appliance into an electrical outlet. Place the handle ③ on the charging base ①, the charge indicator will light up.

A full charge will take 16 hours and will provide 45 minutes operation time. For everyday use, the handle can be stored on the charging base to maintain it at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, the following procedure should be done every 6 months: Do not place the handle back onto the charging base in order to discharge it by regular use. Once discharged, recharge the handle on the charging base.

Operating the switch

Press the switch ④ to operate the toothbrush.

Your toothbrush features a reduced speed setting for sensitive areas of the mouth:

- press the switch to operate with normal speed,
- press a second time to operate with reduced speed,
- press a third time to switch the handle off.








Removing plaque

The toothbrush can be used with any leading brand of toothpaste.

To avoid splashing, guide the brushhead to your teeth before switching the handle on.

Guide the brushhead from tooth to tooth, following the curve of the gum and the shape of each tooth. Hold the brushhead in place for a few seconds before moving on to the next tooth. Brush the gums as well as the teeth. Do not press too hard nor scrub – simply let the brush do all the work.

During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

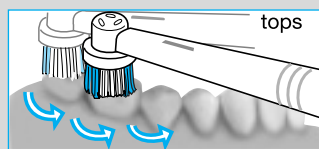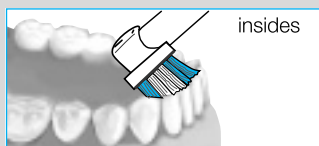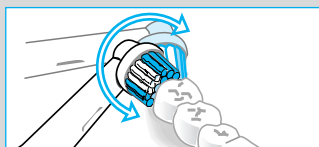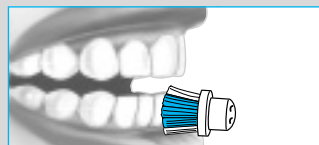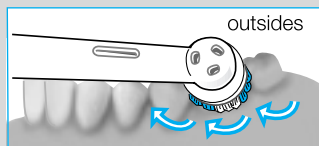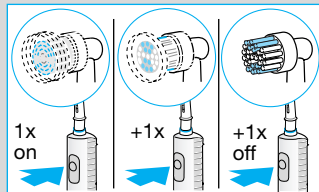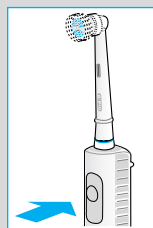
Pressure Sensor

For safety and optimum brushing results, the toothbrush has a pressure control system: If too much pressure is applied, the oscillating movement of the brushhead will continue, but the pulsation will stop. There is a noticeable change in feeling and sound when the pressure sensor is activated.


Memory timer

An incorporated timer memorises the elapsed brushing time, even when the handle is briefly switched off (max. 30 seconds). After two minutes – which is the minimum recommended brushing time – the timer signals by a short stutter.

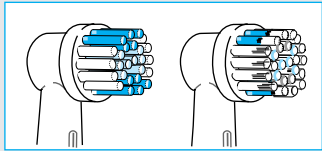INDICATOR® brushhead

Dentists recommend that you replace any toothbrush every 3 months, since worn bristles do not remove plaque as effectively as new ones. Therefore, the brushhead is provided with INDICATOR bristles. With proper brushing twice a day for two minutes using toothpaste, the color will disappear half-way down the bristles in a 3-month period.

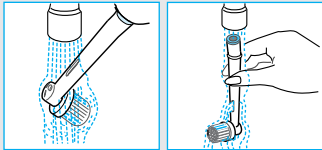After use

Rinse the brushhead thoroughly under running water for several seconds with the handle switched on. Switch off the handle and take off the brushhead. Clean both parts separately under running water, then wipe them dry.

The environment

For environment reasons, your toothbrush uses a «environmentally sound» nickel-hydride rechargeable battery. The battery does not contain any poisonous heavy metals that can damage the environment.

In the interest of recycling, please do not dispose of this appliance in the household waste at the end of its useful life. Instead dispose of it at an authorized Braun Service Center or at an appropriate collection site.

Using the OxyJet® Irrigator

Before using for the first time, fill the water container and allow the water to run through completely by holding the irrigator handle over a washbasin.

1. Attach one of the irrigator jets to the handle until it locks in place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch. (When doing so, the sliding switch on the handle should not be in its upper position.)
3. Select the water pressure with the pressure control dial (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over the washbasin and guide the jet to your teeth and gums. Slide the handle switch up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

You can choose between two different jet functions by sliding the jet switch down or up (before doing so, turn off the water stream by pushing the sliding switch down):










Rotating stream for general cleaning and gum massage (jet switch in lower position)


The rotation of the water stream and its micro-bubbles is designed for easy penetration below the gumline.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position)

Max. operating time: 15 minutes

Cooling off time: 2 hours

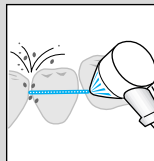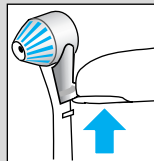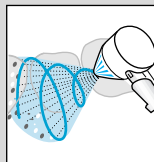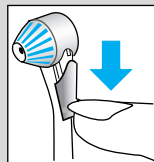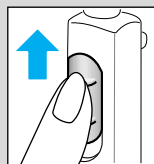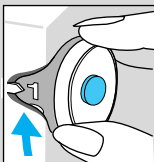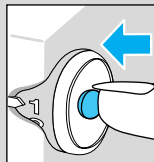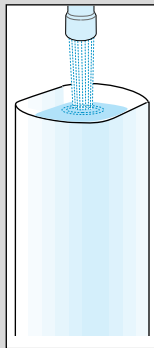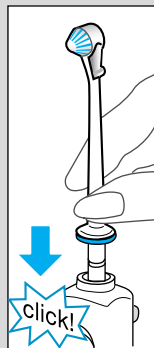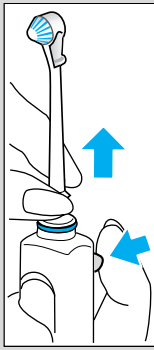After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch to turn off the irrigator. The irrigator will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container, the brushhead compartment and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

How to obtain service in the USA

Should this product require service, please call
1-800-BRAUN 11 (272-8611) to be referred to the Authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.
Pack the product well.
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Contact an authorized Braun Service Center.
Please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-892-8534.

To speak to a Braun Consumer Service Representative:

Please call 1-800-984-9668.

FOR USA ONLY

Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

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For Canada only

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service centre, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-387-6657 to be referred to the authorized Braun Service Center closest to you.

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To obtain original Braun replacement parts & accessories:

Contact an authorized Braun Service Center. Please call 1-800-387-6657 to be referred to the Authorized Braun Service Center closest to you.

To speak to a Braun Service Representative:

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24 Month Performance Guarantee

For USA only

Use the Oral-B OxyJet 3D Center daily until you are fully familiar with it. If you are not convinced that the Oral-B OxyJet 3D Center significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B OxyJet 3D Center within 24 months, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

24 Month Performance Guarantee

Braun Division,
The Gillette Company,
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B OxyJet 3D Center must be shipped in its original packaging and be postmarked no later than 24 months after the date of purchase.

Please allow 6–8 weeks for handling.

Name of Regular Toothbrush

Your Name _____

Street

City

State

Zip Code

24 Month Performance Guarantee

For Canada only

Use the Oral-B OxyJet 3D Center daily until you are fully familiar with it. If you are not convinced that the Oral-B OxyJet 3D Center significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B OxyJet 3D Center within 24 months, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

24 Month Performance Guarantee

Braun Products
Gillette Canada Company
4 Robert Speck Parkway
Mississauga, ON L4Z 4C5

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B OxyJet 3D Center must be shipped in its original packaging and be postmarked no later than 24 months after the date of purchase.

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Gillette Distribuidora, S.A. de C.V.