

Oral-B

**Oral-B** PROFESSIONALCARE™ 9000  
**TRIUMPH**™  
SMART TECHNOLOGY...BRILLIANT RESULTS™

Oral-B

BRAUN

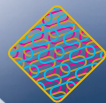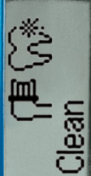**Smart**  
TECHNOLOGY™

powered  
by **BRAUN**

Oral-B

English	4
Español	15
Français	25







**Internet:**

[www.oralb.com](http://www.oralb.com)  
[www.braun.com](http://www.braun.com)

D 25.546.4  
D 25.526.4  
D 25.516.4

**Type 3731**  
**Modèle 3731**  
**Modelo 3731**
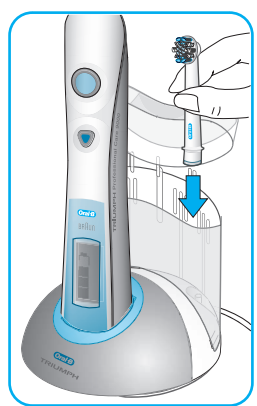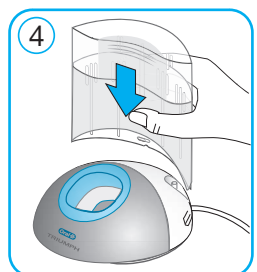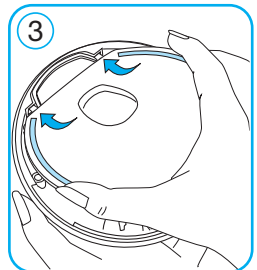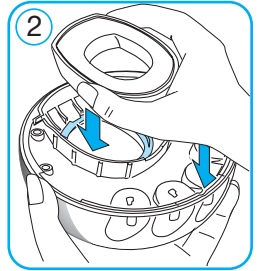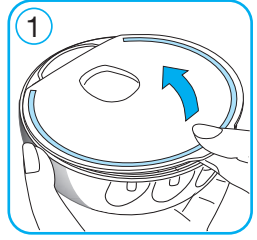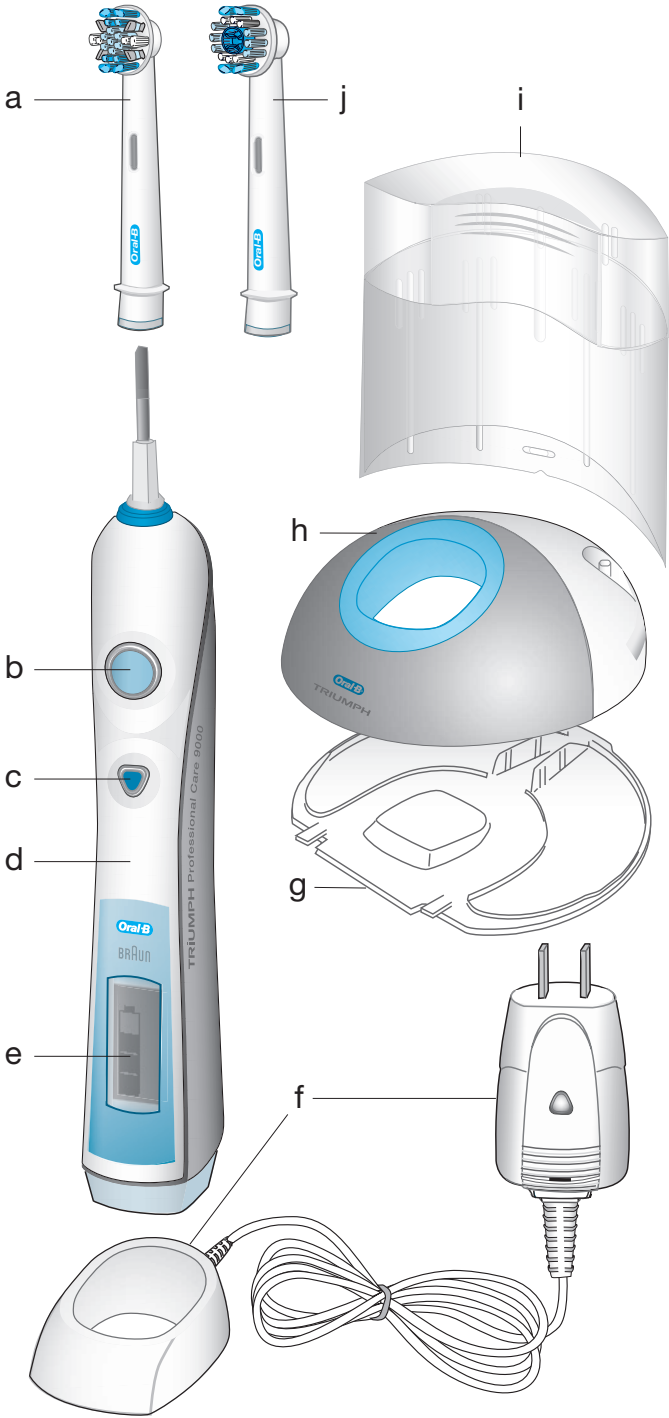
## Congratulations!

You have chosen the world's most technologically advanced rechargeable toothbrush. **Oral-B® Triumph™** will provide you with the superior\* oral care results you have been striving to achieve.

**Oral-B Triumph** has an innovative new brushing technology that is clinically proven to deliver cleaner and whiter teeth. Its advanced Oral-B® FlossAction™ brushhead features MicroPulse™ bristles that gently penetrate deep between teeth for a floss-like clean. It provides such a superior clean, you won't believe how great you feel about your teeth.

**Oral-B Triumph** is a significant advancement in brushing technology featuring the following:

- **Smart Technology**, an interactive feature that gives you feedback to maximize your brushing performance.
- **Innovative Brushheads** with new brushing technology provide complete cleaning for your teeth and outstanding gum care.
- **Different Brushing Modes** let you customize your brushing routine to address your specific needs.
- **An innovative SmartPlug™ charger** makes **Oral-B Triumph** convenient for traveling.

**Oral-B Triumph** is brought to you by Oral-B, the worldwide leader in toothbrushes and the brand more dentists use themselves.

If you have any questions about using the new **Oral-B Triumph** toothbrush, please visit our web site at [www.oralb.com](http://www.oralb.com) or call our toll free telephone number: 1-800-566-7252

Thank you for choosing  
**Oral-B Triumph.**

Sincerely,

Bruce Cleverly  
President, Oral-B

**Smart Technology ... Brilliant Results.™**

\* vs. brushing with Oral-B® FlexiSoft® brushhead

---

# IMPORTANT SAFEGUARDS

**Oral-B Triumph** has been carefully designed to offer you and your family a new and unique brushing experience that is both safe and effective.

However, when using electrical products, basic safety precautions should always be followed. **Please read carefully before use.**

## DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.

## WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has been dropped or damaged, if it is not working properly, or if the SmartPlug is damaged or has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
6. This appliance is provided with a SmartPlug with integrated Safety Extra Low Voltage power supply. Do not exchange or manipulate any part of it.

---

## SAVE THESE INSTRUCTIONS

## Description of Oral-B Triumph



- a Oral-B FlossAction brushhead
- b On/off button
- c Brushing mode button
- d Handle
- e SmartDisplay
- f SmartPlug charger (includes charging ring and SmartPlug)
- g Base plate
- h Base station
- i Protective cover
- j Oral-B® PowerPolisher™ brushhead (not with all models)

## Connecting and charging

**Oral-B Triumph** has a waterproof handle, is electrically safe and is designed for use in the bathroom.



To assemble the base station (h) as a charging unit, detach the base plate (g) first (①, page 3). Place the charging ring into the base station and attach the base plate as shown (② and ③, page 3). Then insert the protective cover (i) on the back of your base station (see diagram ④, page 3).

- Plug the SmartPlug into an electrical outlet and place the toothbrush handle (d) on the charging unit. Alternatively, **Oral-B Triumph** can be charged by placing the handle directly into the charging ring while the SmartPlug is plugged in.

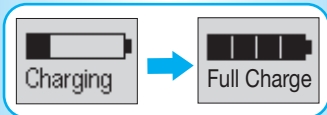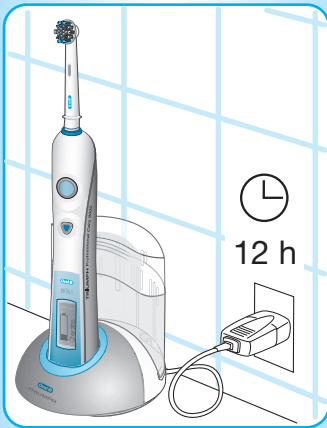- A battery symbol should appear in the SmartDisplay. The battery symbol will fill up while the handle is being charged and «Charging» will be displayed. When the handle is fully charged, the SmartDisplay shows «Full Charge». A full charge takes at least 12 hours and enables approximately two weeks of regular brushing.

Please note: If the battery symbol is not displayed immediately, continue charging. It will appear in 10 - 15 minutes.

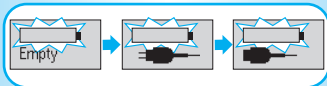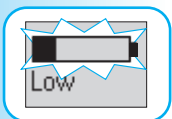
- When the rechargeable battery is low, the bars of the battery symbol will disappear and «Low» will be displayed. Your toothbrush will slow down and then stop. It is time to recharge **Oral-B Triumph**. When the battery is fully discharged, «Empty» is displayed. It will then need a 15-minute charge before you can use it for a 2-minute brushing.

- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

**Oral-B Triumph** is convenient for traveling because you only need the compact SmartPlug charger and an appropriate travel adapter for recharging in foreign countries. The SmartPlug automatically adapts to any voltage worldwide.

For electric specifications, see printing on the SmartPlug.

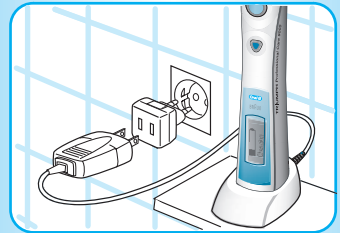
## Personalizing Oral-B Triumph

Before getting started, you can program your toothbrush so that it meets your specific needs. **Oral-B Triumph** has innovative Smart Technology that lets you choose the language and the timer that works best for you.


### Select your language

**Oral-B Triumph** comes with the English language activated. If you wish to switch to another language, proceed as follows: Press and hold the brushing mode button (c) until «English» is displayed. By briefly pressing the on/off button (b) you can scroll through 13 different languages. To confirm your desired language, briefly press the brushing mode button again. This will also lead you directly to the timer choice menu.
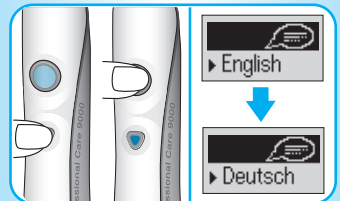
### Select your timer

The «2-Minute Timer» indicates by a long stuttering sound when the dentist recommended 2-minute brushing time has been reached. The «Professional Timer» helps you brush all four quadrants of your mouth equally. It signals with a short stuttering sound at 30-second intervals that it is time to move to the next quadrant of your mouth, followed by a long stuttering sound when you reach the recommended two minute mark.

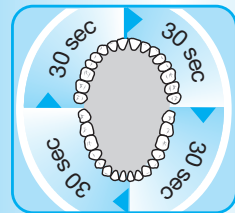
**Oral-B Triumph** comes with the «Professional Timer» activated. If you wish to switch to the «2-Minute Timer» after you have confirmed the language, simply press the on/off button to change the setting. Then press the brushing mode button to confirm your choice.

## Using Oral-B Triumph


Attach an Oral-B brushhead to the handle. You can use any kind of toothpaste. To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
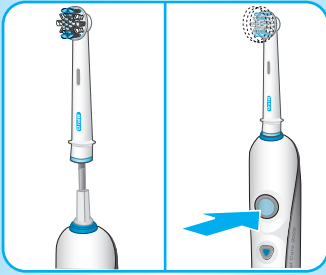
### Brushing technique

- Guide the brushhead slowly from tooth to tooth.
- Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

To help you get optimum brushing results, **Oral-B Triumph** has a pressure control feature to ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brushhead will continue but the pulsation will stop. In addition to feeling the change, you can also hear a different sound.
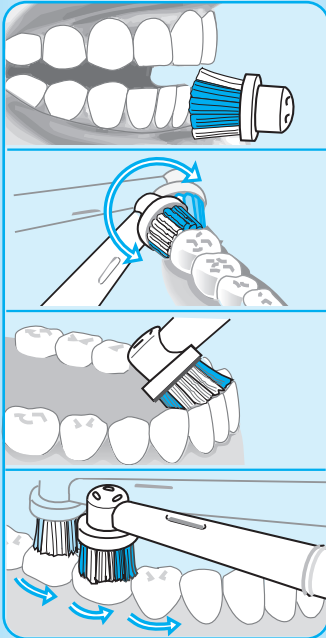
## Brushing modes

**Oral-B Triumph** offers different brushing modes for varying oral care needs:



- «Clean» – Superior whole mouth cleaning.
- «Soft» – Gentle, yet thorough cleaning for sensitive areas.
- «Massage» – Gentle stimulation of gums.
- «Polish» (deluxe model only) – All over polishing and whitening for occasional or everyday use.

How to switch between brushing modes:

- **Oral-B Triumph** automatically starts in the «Clean» mode when pressing the on/off button.
- To shift from the «Clean» to the «Soft» mode, press the brushing mode button once.
- If you wish to switch to the «Massage» mode, press the brushing mode button again and let the brushhead pulse against your gum line.
- The «Polish» mode (deluxe model only) is best used with an Oral-B PowerPolisher brushhead. Once the Oral-B PowerPolisher is attached, press the on/off button and your toothbrush automatically starts in the «Polish» mode (please see «Brushheads»).
- To return from any other mode to the «Clean» mode, press and hold the brushing mode button.
- To turn your brush off, press the on/off button.


During brushing, the SmartDisplay shows the brushing mode you have selected, followed by your actual brushing time. When you finish brushing, the total brushing time will be displayed for 30 seconds before the display switches back to the battery charge level. To directly get to the battery charge level display, you can clear the brushing time display by pushing the brushing mode button.

In the first few days of using **Oral-B Triumph**, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Soft» mode.

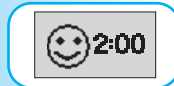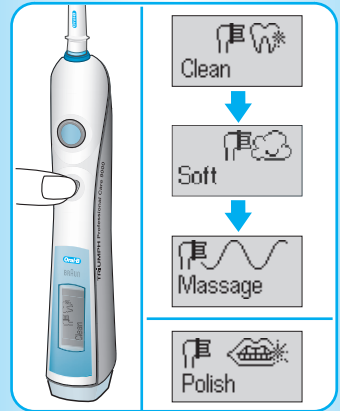
## Brushheads

The Oral-B FlossAction and Oral-B PowerPolisher brushheads contain a «smart chip» to communicate with the electronics inside the **Oral-B Triumph** handle. The «smart chip» works for you, automatically switching to the appropriate brushing mode for the attached brushhead. Nevertheless, you can select other brushing modes to customize your routine.
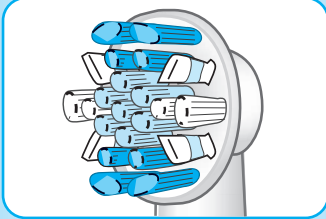
### Oral-B FlossAction brushhead

The Oral-B FlossAction brushhead is recommended for everyday deep cleaning. It features MicroPulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque. The brushhead's outer bristles surround each tooth, providing superior overall plaque removal.
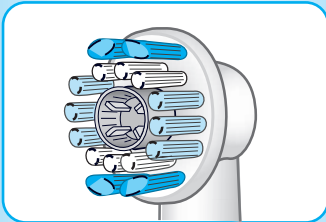
### Oral-B PowerPolisher brushhead (not with all models)

The Oral-B PowerPolisher brushhead is for occasional or everyday polishing. It features a dentist-inspired polishing cup to remove stains from coffee, tea and tobacco for effective brightening.




We do not recommend using the Oral-B PowerPolisher or Oral-B FlossAction brushhead with braces. The Oral-B Ortho brushhead is specifically designed to clean around braces and wires and can be used with the **Oral-B Triumph**.

Replace  
Brushhead

**Oral-B Triumph** has Smart Technology that helps you remember when to replace your brushhead to get optimum results. Dentists recommend that you change your brushhead every three months. As soon as your Oral-B FlossAction or Oral-B PowerPolisher brushhead has surpassed the recommended usage, the SmartDisplay shows «Replace Brushhead». Replace with a new brushhead to get optimum brushing results. If you notice unusual wear before receiving the «Replace Brushhead» message, replace your brushhead immediately.

## Cleaning recommendations

- After brushing, rinse your brushhead thoroughly under running water.
- Remove the brushhead from the handle and clean both parts separately under running water. Then wipe them dry before you put the handle on the charging unit.
- Before cleaning your charging unit, make sure you disassemble it first. The charging ring and the SmartPlug should never be placed in the dishwasher. However, all other parts (base station, base plate and protective cover) are dishwasher-safe.
- From time to time, you can also clean the charging ring and the SmartPlug with a damp cloth. The SmartPlug charger should never be placed in water.
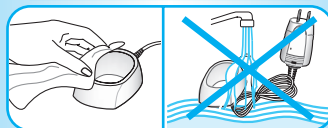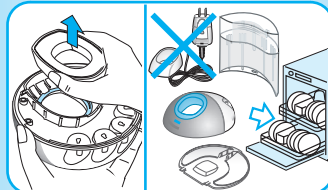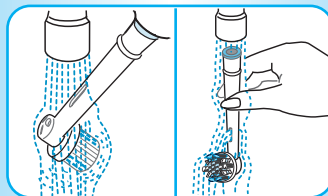
## Environmental notice


This product is provided with a nickel-hydride battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

### Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.
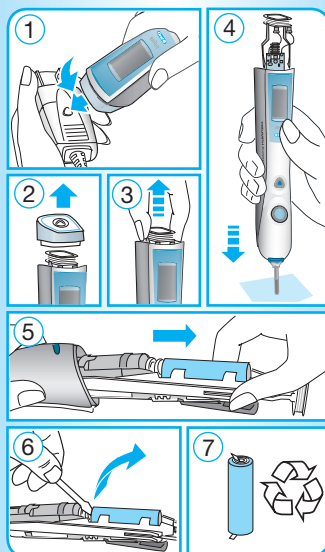
# FOR USA ONLY

## Limited 2-Year Warranty

Oral-B Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

## How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B Braun replacement parts and accessories:*

Contact an authorized Oral-B Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

*To order original Oral-B Braun replacement parts and accessories (including replacement brushheads) by phone:*

Please call 1-800-566-7252.

*To speak to an Oral-B Braun Consumer Service Representative:*

Please call 1-800-566-7252.

# FOR CANADA ONLY

## Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Oral-B Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Oral-B Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

## How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B Braun replacement parts & accessories:*

Contact an authorized Oral-B Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

*To speak to an Oral-B Braun Service Representative:*

Please call 1-800-566-7252.

# The 30-Day Challenge

Take the challenge! Experience the **Oral-B Triumph** difference!

**Oral-B Triumph** is the world's most technologically advanced rechargeable toothbrush for superior cleaning, natural whitening, gum care and freshening of your whole mouth.

No wonder that Oral-B is the brand more dentists use themselves worldwide. In fact, we are so convinced you'll fall in love with the new **Oral-B Triumph** toothbrush that we give you this risk-free offer.

Try **Oral-B Triumph** for 30 days. If you are not 100% satisfied, receive a full refund. No questions asked. Guaranteed.

For a full refund, please return your handle, refills, and charging unit with your original cash register receipt. For further details, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

## U.S. only:

Attention: 30-Day Challenge  
Oral-B Braun Division  
The Gillette Company  
Sobin Park and A Street, MS 1F  
South Boston, MA 02127

## Canada only:

Attention: 30-Day Challenge  
Oral-B Braun Products  
Gillette Canada Company  
4 Robert Speck Parkway  
Mississauga ON L4Z 4C5

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