Copyright 2013 © InnoEssentials All rights reserved.  
Print date: July 2014

InnoEssentials reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by InnoEssentials is believed to be accurate and reliable at the time of publication. However, InnoEssentials assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of InnoEssentials.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of InnoEssentials.  
Data is subject to change without notification.

InnoEssentials has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from InnoEssentials. Specifications are subject to change without notice.

InnoEssentials International Invention Works B.V. Donk 1B 2991 LE Barendrecht, The Netherlands

## Table of Contents

|                                                                            |           |
|----------------------------------------------------------------------------|-----------|
| <b>1. Before You Start</b> .....                                           | <b>2</b>  |
| 1.1. Silk'n Glide™ device description .....                                | 2         |
| 1.2. Package contents .....                                                | 2         |
| <b>2. Silk'n Glide™ Intended Use</b> .....                                 | <b>3</b>  |
| <b>3. Safety with Silk'n Glide™</b> .....                                  | <b>3</b>  |
| <b>4. Contraindications</b> .....                                          | <b>4</b>  |
| <b>5. Do not use Silk'n Glide™ if:</b> .....                               | <b>6</b>  |
| <b>6. Avoid using Silk'n Glide™ if any of the following applies:</b> ..... | <b>7</b>  |
| <b>7. Possible Side Effects using Silk'n Glide™</b> .....                  | <b>8</b>  |
| <b>8. Long Term Hair Removal the Silk'n Glide™ Way</b> .....               | <b>9</b>  |
| <b>9. Get Started!</b> .....                                               | <b>10</b> |
| 9.1. What to Expect when treating with Silk'n Glide™? .....                | 13        |
| 9.2. After treating with Silk'n Glide™ .....                               | 13        |
| <b>10. Maintenance of Silk'n Glide™</b> .....                              | <b>13</b> |
| 10.1. Cleaning Silk'n Glide™ .....                                         | 13        |
| <b>11. Troubleshooting</b> .....                                           | <b>14</b> |
| 11.1. "My Silk'n Glide™ does not start!" .....                             | 14        |
| 11.2. "A light pulse is not emitted when I press the PULSE SWITCH" .....   | 14        |
| <b>12. Customer Service</b> .....                                          | <b>15</b> |
| <b>13. Specifications</b> .....                                            | <b>15</b> |

Before using Silk'n Glide™ for the first time, please read this User Manual in its entirety. Pay particular attention to sections on device use procedures, device operation, and after-use procedures. We recommend you re-familiarize yourself with this User Manual before each use of Silk'n Glide™. Silk'n Glide™ is a powerful electrical device. As such, it should be used with special attention to safety. Please read all warnings and safety precautions before use, and strictly follow them when using Silk'n Glide™


## 1. Before You Start

**Silk'n Glide™ is a light-based device for permanent hair reduction designed for home-use.**

### 1.1 Silk'n Glide™ device description

The Silk'n Glide™ HAIR REMOVAL DEVICE has a TREATMENT SURFACE, a CONTROL PAD and a PULSE TRIGGER.

- ① Control pad
- ② Pulse Trigger
- ③ Skin Color Sensor
- ④ Treatment Surface
- ⑤ Central button
- ⑥ Energy level indicator lights
- ⑦ Ready/Skin tone warning indicator light
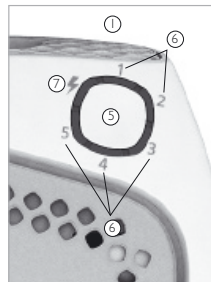
### 1.2. Package contents

Upon opening the Silk'n Glide™ package, you will find the following parts:

- Silk'n Glide™ Hair Removal Device
- Adaptor
- This User's Manual
- Silk'n warranty card

## 2. Silk'n Glide™ Intended Use

Silk'n Glide™ is intended for removal of unwanted hair. Silk'n Glide™ is also intended to effect long term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime.

Silk'n Glide™ may be used to remove unwanted body hair. Ideal body areas for Silk'n Glide™ use include the underarms, bikini line, arms, legs, face, back, shoulders and chest.

## 3. Safety with Silk'n Glide™

**With Silk'n Glide™ Safety Comes First.**

**HPL™ technology in Silk'n Glide™ - Superior safety with lower energy level:**

Home Pulsed Light™ technology is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy used in Silk'n Glide™ reduces its potential to cause harm or complications, and contributes to your overall safety.

**Silk'n Glide™ protects your skin:**

Light based hair removal is not suitable for naturally darker skin tones or tanned skin. Silk'n Glide™ comes with a built-in SKIN COLOR SENSOR that is designed to measure the complexion of the applied surface and enable application only on suitable skin complexions. This unique safety feature will not let you treat where your skin is too dark or too tanned.

Furthermore, the 3cm² TREATMENT SURFACE through which pulses of light are delivered is recessed inside the device. This enables Glide™ to protect your skin by avoiding direct contact between the TREATMENT SURFACE and the skin.

**Silk'n Glide™ protects your eyes:**

Silk'n Glide™ has a built-in safety feature for eye protection. It has been designed so that a light pulse can not be emitted when the TREATMENT SURFACE is facing open air. The safety feature enables treatment only when the TREATMENT SURFACE is in contact with the tissue.

## 4. Contraindications

### Important Safety Information – Read Before Use!

Silk'n Glide™ is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice you may also visit [www.silk'n.com](http://www.silk'n.com) or your local Silk'n™ domain.

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.

| Warning                                                                                                                                                                                                                                          | Possible results                                                                                                                                                                                                                                                                                                                                                                                                                                                                            | Safety Tip                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <p><b>Do not use Silk'n Glide™ on the naturally darkest skin complexion.</b> <small>Not Safe!</small></p>  <p>See hair and skin color chart on the package.</p> | <p>Treating dark skin with Silk'n Glide™ can result in adverse effects such as redness and skin discomfort</p>                                                                                                                                                                                                                                                                                                                                                                              | <p>A unique SKIN COLOR SENSOR is embedded in Silk'n Glide™ to measure the treated skin complexion at the beginning of each session and occasionally during the session. The SKIN COLOR SENSOR ensures that pulses will only be emitted on suitable skin tones.</p>                                                                                                                                                                                                                                                                                                                                                                                                                                   |
| <p><b>Do not use Silk'n Glide™ energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Do not expose treated areas to the sun.</b></p>                                                            | <p>Treating with Silk'n Glide™ energy levels 2-5 before or after any sun exposure can result in adverse effects such as redness and skin discomfort. Tanned skin particularly following sun exposure, contains large quantities of the pigment Melanin. This applies to all skin types and complexions, including those which don't seem to tan quickly. The presence of large quantities of Melanin exposes the skin to higher risk when using any method of light based hair removal.</p> | <p>If you are using any energy level higher than the lowest, avoid exposure to the sun for 4 weeks before and 2 weeks after your Silk'n Glide™ treatment! Sun exposure includes unprotected exposure to direct sunlight of over 15 minutes constantly, or unprotected exposure to diffused sunlight of over 1 hour constantly. To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen SPF 30 or higher; for 2 weeks after each hair removal session. Silk'n Glide's lowest energy level (1) is designed especially for safe use even after sun exposure. After sun exposure you may continue using Silk'n Glide™, but only at the lowest energy level.</p> |

| Warning                                                                                                                                                                             | Possible results                                                                                                                                                                               | Safety Tip                                                                                                                                                                                                                                                                      |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Never use Silk'n Glide™ around or near the eyes.                                                                                                                                    | Although, Silk'n Glide™ can be used for the treatment of facial hair, extra caution should be used on the face to avoid the eyes.                                                              | Silk'n Glide™ has a built-in safety feature for eye protection. It has been designed so that a light pulse can not be emitted when the TREATMENT SURFACE is facing open air. The Pulse Trigger enables treatment only when the TREATMENT SURFACE is in contact with the tissue. |
| Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, birth marks, moles or warts), nipples, genitals or lips.                           | Treating any area with the Silk'n Glide™ that is dark in color or has more pigment can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypopigmentation). | Cover area with material that will not absorb the light such as a white cloth or white medical tape.                                                                                                                                                                            |
| Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example) | Light based treatments can cause adverse reactions on already sensitive areas.                                                                                                                 | Wait for the effected area to heal before using Silk'n Glide™.                                                                                                                                                                                                                  |
| A treatment area that has a history of herpes outbreaks.                                                                                                                            | Light based treatments can cause adverse reactions on already sensitive areas.                                                                                                                 | Consult your physician to received preventative treatment before using Silk'n Glide™.                                                                                                                                                                                           |
| Do not treat with Silk'n Glide™ if you are pregnant or nursing.                                                                                                                     | Hormonal changes during pregnancy could increase skin sensitivity and the risk of adverse effects.                                                                                             |                                                                                                                                                                                                                                                                                 |
| Do not treat on areas where you may want hair to grow.                                                                                                                              | Results are permanent.                                                                                                                                                                         |                                                                                                                                                                                                                                                                                 |
| Keep Silk'n Glide™ away from water! Do not place or store Silk'n Glide™ where it can fall, be pushed or placed into a tub, sink or any other vessel containing water.               | This may cause severe electrocution.                                                                                                                                                           | Do not use Silk'n Glide™: while bathing or if it becomes damp or wet. Unplug immediately if it has fallen into water.                                                                                                                                                           |

| Warning                                                     | Possible results                                                                                                                                                                                                                             | Safety Tip                                                                                        |
|-------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------|
| Do not attempt to open or repair your Silk'n Glide™ device. | Opening Silk'n Glide™ may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury. It may also damage your device and will void your warranty. | Please contact Silk'n™ Customer Service if you have a broken or damaged device in need of repair. |

## 5. Do not use Silk'n Glide™ if:

|                                                                                                             |                                                                      |                                                            |
|-------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------|------------------------------------------------------------|
| The device is operated by any plugs, attachments, or accessories other than those provided with your Glide™ | the attachments or accessories are not recommended by InnoEssentials | Stop use immediately and contact Silk'n™ Customer Service. |
| it has a damaged cord or plug.                                                                              | Keep the power cord away from heated surfaces.                       | Stop use immediately and contact Silk'n™ Customer Service. |
| you see or smell smoke when it is in use.                                                                   |                                                                      | Stop use immediately and contact Silk'n™ Customer Service. |
| it is not working properly or if it appears damaged.                                                        |                                                                      | Stop use immediately and contact Silk'n™ Customer Service. |
| its SKIN COLOR SENSOR is cracked, or broken.                                                                |                                                                      | Stop use immediately and contact Silk'n™ Customer Service. |
| the outer shell is cracked or is coming apart.                                                              |                                                                      | Stop use immediately and contact Silk'n™ Customer Service. |

## 6. Avoid using Silk'n Glide™ if any of the following applies:

|                                                                                                                                |                                                               |
|--------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------|
| History of keloidal scar formation.                                                                                            | Consult with your physician prior to using the Silk'n Glide™. |
| If you suffer from epilepsy.                                                                                                   |                                                               |
| If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc. |                                                               |
| If you have a history of skin cancer or areas of potential skin malignancies.                                                  |                                                               |
| If you have received radiation therapy or chemotherapy treatments within the past 3 months.                                    |                                                               |
| If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.                  |                                                               |
| If your skin was treated with a physician-prescribed drug within the past 6 months, please consult with your physician.        |                                                               |

## 7. Possible side effects using Silk'n Glide™

When used according to the instructions, side effects and complications associated with use of Silk'n Glide™ are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

| Adverse Reaction                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    | Degree of Risk |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----------------|
| <b>Minor Skin Discomfort –</b><br>Although home pulsed light hair removal is generally very well-tolerated, most users do feel some mild discomfort during use, usually described as being a mild stinging sensation on the treated skin areas. The stinging sensation usually lasts during the time of the application itself or for a few minutes thereafter. Anything beyond this minor discomfort is abnormal and means that either you should not continue to use Silk'n Glide™ because you are unable to tolerate the hair removal application, or that the energy level setting is too high. | Minor          |
| <b>Skin Redness –</b><br>Your skin may become red right after using Silk'n Glide™ or within 24 hours of using Silk'n Glide™. Redness generally clears up within 24 hours. See your doctor if redness does not go away within 2 to 3 days.                                                                                                                                                                                                                                                                                                                                                           | Minor          |
| <b>Increased Sensitivity to Skin –</b><br>The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.                                                                                                                                                                                                                                                                                                                                                                                                                                                       | Minor          |
| <b>Pigment Changes –</b><br>Silk'n Glide™ targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is risk of temporary hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin tones. Usually discoloration or changes to skin pigment are temporary and permanent hyperpigmentation or hypopigmentation rarely occur.                                                        | Rare           |
| <b>Excessive Redness and Swelling –</b><br>In rare cases treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but one should avoid exposure to sun.                                                                                                                                                                                                                                                        | Rare           |

## 8. Long Term Hair Removal the Silk'n Glide™ Way

A typical full hair growth cycle may take 18-24 months. During this time multiple Glide™ sessions may be required in order to achieve permanent hair removal.

The efficiency of hair removal varies from person to person according to body area, hair color; and how Silk'n Glide™ is used.

Typical Silk'n Glide™ hair removal plan during a full hair growth cycle:

Treatments 1-4 – plan two weeks apart

Treatments 5-7 – plan four weeks apart

Treatments 8 + – treat as needed, until desired results are achieved.

### Expected Treatment Outcome for Hair Removal

| Fitzpatrick skin type | Hair Color           | Light Fluence [J/cm²] | Anatomic area of the body <sup>2</sup> | Average Number of Treatments <sup>1</sup> | Effectiveness (% of hair removal/regrowth) <sup>1</sup> |
|-----------------------|----------------------|-----------------------|----------------------------------------|-------------------------------------------|---------------------------------------------------------|
| I-IV                  | Light brown to brown | 3-4.5                 | Lower or upper legs                    | Session 10-12, approx. four weeks apart.  | 60% Hair reduction                                      |
|                       |                      |                       | Arms                                   |                                           | 50% Hair reduction                                      |
|                       |                      |                       | Bikini Line                            |                                           | 60% Hair reduction                                      |
|                       |                      |                       | Underarms                              |                                           | 55% Hair reduction                                      |
| I-IV                  | Brown to black       | 3-4.5                 | Lower or upper legs                    | Session 8-10, approx. four weeks apart.   | 70% Hair reduction                                      |
|                       |                      |                       | Arms                                   |                                           | 60% Hair reduction                                      |
|                       |                      |                       | Bikini Line                            |                                           | 70% Hair reduction                                      |
|                       |                      |                       | Underarms                              |                                           | 65% Hair reduction                                      |

1. Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

2. One cannot expect permanent hair removal in a single or even double treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

## 9. Get Started!

1. Remove Silk'n Glide™ HAIR REMOVAL DEVICE and other components from box.
2. Plug the ADAPTOR cord into the Silk'n Glide™ HAIR REMOVAL DEVICE socket.
3. Plug the ADAPTOR into an electrical outlet. All indicator lights in the CONTROL PAD will turn on and off repeatedly in a circular motion.

Your Silk'n Glide™ is now ready to start.

### **Treat with Silk'n Glide™ safely and properly!**

**Choose areas to treat that have not had recent exposure to the sun.**

4. The skin of the treatment areas should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. Do not wax, pluck or tweeze the hair out. Shaving is an important step to get the results you want! Best results will be reached when you shave 3 days prior to the treatment or when the hair is 1-2 mm - 3/32 inch.
5. Press the CENTRAL BUTTON. The device will turn on, the Energy Level 1 INDICATOR LIGHT will turn on and a fan sound will start.
6. Approximately 1 second after pressing the CENTRAL BUTTON, the READY/SKIN TONE WARNING INDICATOR LIGHT will turn on. The device is then ready for you to trigger the first pulse at the lowest energy level.
7. If this is your first treatment - Choose the energy level. Glide™ offers 5 levels of energy from 1 (the lowest) to 5 (the highest). To choose the right energy level for your first treatment, test each treatment area separately:

- Place the TREATMENT SURFACE on the skin and press the trigger to deliver one pulse. If you did not experience abnormal discomfort press the CENTRAL BUTTON again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.
- Continue testing at increased energy levels until you have reached the highest level that you feel comfortable with. If within one hour you experience no adverse effects, you can start your full treatment at that energy level.
- Repeat this test for each body part you intend to treat.

8. Set the energy level. To set an energy level higher or lower, press the CENTRAL BUTTON repeatedly to increase the energy level from 1 to 5 and from 5 directly back to 1, until the desired energy level is set. The ENERGY LEVEL INDICATOR LIGHTS will indicate the energy level setting. (For details see "Energy Level" chapter inside this user manual).

### Energy Level

Energy level determines the intensity of the Silk'n Glide™ light pulse delivered to your skin. As energy level increases, so does hair removal results as well as the risk of possible side effects and complications. Glide™ offers 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 ENERGY LEVEL INDICATOR LIGHTS.

The energy level set for the treatments determines the number of pulses the device can produce. Whenever Silk'n Glide™ is turned on its energy level will automatically be reset to the lowest energy level. Only one ENERGY LEVEL INDICATOR LIGHT will be on. To set the energy level, press the CENTRAL BUTTON repeatedly to increase the energy level from 1 to 5 and from 5 directly back to 1.

9. Apply the TREATMENT SURFACE to the skin, making sure the skin is spread evenly and smoothly. As soon as the TREATMENT SURFACE is in full contact with the skin, the READY/SKIN TONE WARNING INDICATOR LIGHT will start blinking slowly.
10. Press the PULSE TRIGGER. The device will first determine the color of your skin. If the color of the skin is light enough for safe application, the device will flash a pulse of light onto your skin, and the READY/SKIN TONE WARNING INDICATOR LIGHT will turn off. You will see a bright flash of light and simultaneously hear a subtle pop sound, which is a normal noise for the device. You will feel a mild sensation of warmth and tingling.  
Silk'n Glide™ will immediately recharge for the next pulse.  
After a few seconds (depending on the set energy level) the READY/SKIN TONE WARNING INDICATOR LIGHT will turn on again, and start blinking slowly when attached again to the skin.
11. Remove the TREATMENT SURFACE from the treated area of skin.
12. Choose a method of application for your session. Glide™ can be applied in two ways: "Pulsing" and "Gliding" (for more details see application methods frame below).
  - For the PULSING method, follow steps 9-11 above on another spot repeatedly until the entire area is fully covered.
  - For the GLIDING method, hold the PULSE TRIGGER down continuously, and slide the device smoothly along the treated area.

If your skin blisters or burns, STOP USE IMMEDIATELY!

### Methods of application

#### PULSING

PULSING is the best technique for working with higher energy levels as it enables better control of skin coverage. When pulsing, Glide™ pulses should be administered in rows, starting at one end of each row

and progressing sequentially towards the other end helping you avoid treating the same area more than once or overlapping skin areas. When pulsing, the Silk'n Glide™ TREATMENT SURFACE is designed to create temporary pressure marks on the treated area. In every repositioning of the TREATMENT SURFACE use these pressure marks to guide you for proper positioning of the next pulse, avoiding both gaps and overlaps between pulses.

Try to avoid overlapping pulses!

Do not treat the same area of skin more than once per hair removal session!  
Treating the same area of skin more than once per session increases the likelihood of adverse effects.

## GLIDING

GLIDING is the best way to work within the low energy levels and achieve a fast and more efficient session. When Gliding, hold the PULSE TRIGGER down continuously, and using a slow sliding motion while keeping contact with the skin, smoothly move the device along the treated area. In the lower energy levels the device takes only 1-1.5 seconds to recharge between pulses and provides a nonstop sequence of flashes.

If your skin blisters or burns, STOP USE IMMEDIATELY!

## The Skin Color Sensor

Light-based hair removal on darker skin tones can result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation).

To prevent such misuse, a unique SKIN COLOR SENSOR in Silk'n Glide™ measures the treated skin complexion at the beginning of each session and occasionally during the session. If the SKIN COLOR SENSOR detects a skin tone that is too dark for Silk'n Glide™ application, the device will automatically stop emitting pulses.

If you see no light pulse and the READY/SKIN TONE WARNING INDICATOR LIGHT is rapidly blinking, this is an indication that your skin tone, measured by the SKIN COLOR SENSOR, is too dark for safe application. Try using the device on a different body part or contact the Silk'n™ support.

### 9.1. What to Expect when treating with Silk'n Glide™

For many people, using Silk'n Glide™ may be their first experience with a light-based device designed for home use. Silk'n Glide™ is simple to use, and hair removal sessions go by quickly. During a Silk'n Glide™ session it is normal to experience and feel:

- **A Fan Noise** – The cooling fan in Silk'n Glide™ makes noise similar to a hairdryer. This is normal.
- **A Pop Sound with Each Pulse** – When a pulse of HPL™ light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.
- **A Sensation of Warmth and Tingling** – During each pulse of light it is normal to feel a mild sensation of warmth and tingling from the light energy.
- **Some Mild Red or Pink Color** – During and just after your Silk'n Glide™ session it is not uncommon to see some very mild, pink-like color of the skin. This is usually most noticeable around the hairs themselves. However if you see full redness of the skin, blistering or burns stop use of Silk'n Glide™ immediately.

### 9.2. After treating with Silk'n Glide™

1. When Glide™ session has been completed turn Glide™ off by giving the CENTRAL BUTTON an extended press. Be sure to remember the last energy level setting you used, as it will not be restored when turning Glide™ on again.
2. Unplug the ADAPTOR cord from the electrical outlet.
3. After each hair removal session it is recommended that you clean your Silk'n Glide™ device, especially the TREATMENT SURFACE (See: "Cleaning Silk'n Glide™").
4. After cleaning, it is recommended to store your Silk'n Glide™ device in its original box, and keep it away from water.

#### Skin care following hair removal session

Do not expose areas of skin treated in energy levels 2-5 to the sun. Be sure to carefully protect the treated skin with sunscreen, throughout the hair removal period and for at least 2 weeks following the last Silk'n Glide™ session.

## 10. Maintenance of Silk'n Glide™

### 10.1. Cleaning Silk'n Glide™

After each hair removal session, it is recommended to clean your Silk'n Glide™ device, and especially the TREATMENT SURFACE.

1. Unplug Silk'n Glide™ before cleaning.
  2. Use a dry, clean cloth and a specially formulated cleaning detergent for electronic equipment to gently wipe Silk'n Glide™ surface, and especially the TREATMENT SURFACE.
- Never immerse Silk'n Glide™ or any of its parts in water!

## 11. Troubleshooting

### 11.1. “My Silk’n Glide™ does not start.”

- Make sure the ADAPTOR is properly connected to the Glide™ device.
- Make sure the ADAPTOR is plugged into an electrical outlet on the wall.

### 11.2. “A light pulse is not emitted when I press the PULSE TRIGGER”

Make sure that the READY/SKIN TONE WARNING INDICATOR LIGHT is blinking slowly to ensure you have good contact with the skin and that the TREATMENT SURFACE is evenly and firmly pressed to the skin. For your safety, the PULSE TRIGGER will activate a pulse only if the TREATMENT SURFACE is firmly pressed against the skin.

- Check the SKIN TONE WARNING INDICATOR LIGHT. If it blinks rapidly, this is an indication that your skin tone, measured by the SKIN COLOR SENSOR, is too dark for safe application. Try using the device on a different body part or contact the Silk’n™ support.
- If the SKIN TONE WARNING INDICATOR LIGHT is constantly on, together with any number of indicator lights, you have a system error. If this problem persists, contact your local Glide™ Customer Service Center.
- If these problems persist, contact your local Silk’n™ Customer Service Center.
- Check if all lights are blinking together in a circle. If so, unplug the device from its power source. The device has reached its End of Life stage. This device is to be orderly disposed of. Deliver your old equipment to the nearest electronic waste collection point for safe and efficient disposal.

Do not attempt to open or repair your Silk’n Glide™ device. Only authorized Silk’n™ repair centers are permitted to perform repairs.

Opening Silk’n Glide™ may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury.

Trying to open Silk’n Glide™ may also damage the device and will void your warranty.

Please contact Silk’n™ Customer Service if you have a broken or damaged device in need or repair.

## 12. Customer Service

For more information about Silk'n™ products please enter your local Silk'n™ website, [www.silkn.eu](http://www.silkn.eu) or [www.silkn.com](http://www.silkn.com).

If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service:

For USA and Canada customers: 1-877-DO-SILKN / 1-877-367-4556, [contact@silkn.com](mailto:contact@silkn.com)

For customers in other countries: [info@silkn.com](mailto:info@silkn.com)

## 13. Specifications

|                                 |                                                                         |                |
|---------------------------------|-------------------------------------------------------------------------|----------------|
| Spot Size                       | 0.9cm x 3cm [2.7cm²]                                                    |                |
| Speed                           | energy level 1 - 1 pulse every 1 second                                 |                |
|                                 | energy level 2 - 1 pulse every 1.5 seconds                              |                |
|                                 | energy level 3 - 1 pulse after 2.3 seconds                              |                |
|                                 | energy level 4 - 1 pulse after 3 seconds                                |                |
|                                 | energy level 5 - 1 pulse after 3.4 seconds                              |                |
| Technology                      | Home Pulsed Light™                                                      |                |
| Max Energy Level                | Max 5J/cm²                                                              |                |
| Wavelength                      | 475-1200nm                                                              |                |
| Charging Time/Power Source      | Continuous operation                                                    |                |
| Electrical Requirements         | 100-240VAC, 2A                                                          |                |
| Time Needed to Treat Lower Legs | 10 minutes on average                                                   |                |
| Operation and Safety            | Skin Color Sensor seamlessly ensures use only on appropriate skin types |                |
| Package Size                    | 180x180x100cm                                                           |                |
| System weight                   | 200gr                                                                   |                |
| Temperature                     | Operating                                                               | 10°C to 35°C   |
|                                 | Storage                                                                 | -40 to +70°C   |
| Relative Humidity               | Operating                                                               | 30 to 75%rH    |
|                                 | Storage                                                                 | 10 to 90%rH    |
| Atmospheric pressure            | Operating                                                               | 700 to 1060hPa |
|                                 | Storage                                                                 | 500 to 1060hPa |
| Adaptor Rated Model             | KSA24A1200150HE in the UK KSA24A1200150HK   100-240V; 50-60Hz; 0.5A     |                |

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Tous droits réservés. Copyright 2013 © InnoEssentials
**Copyright © 2011 InnoEssentials Με επιφύλαξη παντός δικαιώματος.
Copyright 2013 © InnoEssentials Alle rettigheter reservert.
Copyright 2013 © InnoEssentials Kaikki oikeudet pidätetään.
Copyright 2013 © InnoEssentials Alla rättigheter förbehållna.
Copyright 2013 © InnoEssentials Alle rettigheder forbeholdes.