TARJETA DE EMERGENCIA:

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01/13

14464, 14474

Please read the following instructions and warnings carefully. Keep this instruction manual for future reference.

Adult assembly required.

For Technical Support Call
1-800-268-6237 or E-Mail at
customerservice@summerinfant.com.

For use with children from birth & up.

INSTRUCTION MANUAL

Complete Nursery Care Kit

WARNING:

DO NOT TIE MEDI-DISPENSER AROUND CHILD'S NECK AS IT PRESENTS A STRANGULATION DANGER.

CAUTION:

THE COMPONENTS OF THIS PRODUCT, EXCEPT FOR THE TEETHER, ARE NOT TOYS. DO NOT ALLOW CHILDREN TO PLAY WITH THEM.

BRUSH & COMB:

Soft brush bristles are gentle on baby’s scalp.Wide-toothed comb makes untangling hair easier. Handle design with soft grip texture is easy for toddlers to hold and use.

- To use:**
- Wash brush & comb before using.
 - Untangle hair from the roots out. Keep a hand on the roots while brushing or combing to avoid pulling on child’s scalp.
 - Encourage your toddler to try brushing his or her own hair with your assistance (approximately 3 years old).
- To clean:** Wash with mild soap and warm water. Rinse thoroughly; air dry.

NAIL CLIPPERS:

This standard size nail clipper’s handle rotates 360 degrees so you can safely trim nails from any angle. Soft over-molded handle is comfortable and helps you maintain a steady grip.

Note: Do not clip too close to baby’s skin. We recommend clipping baby’s nails while he or she is sleeping.

To clean: Sterilize with isopropyl rubbing alcohol and rinse in warm water. Wipe dry.

MEDICINE SYRINGE:

Easily dispense and measure medicine.

To use: Clean before first use. Push the inner half of the syringe all the way in and place the tip of the syringe into the liquid medicine. Pull the inner half out until the medicine reaches the correct marked dosage. Make sure your child is elevated to avoid choking and then place one finger inside their mouth and gently pull the cheek outward. Place the syringe toward the back of the cheek (not toward the back of the mouth to avoid gagging) and push down on the inner plunger to dispense the liquid.

To clean: Remove the inner plunger and wash both halves in warm, soapy water. Air dry. Top rack dishwasher safe. Do not place in boiling water.

5 EMERY BOARDS:

- Gently glide across baby’s fingernails to smooth them.

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NASAL ASPIRATOR:

To use:

- Hold baby upright.
 - While squeezing the aspirator bulb, gently place the tip of the aspirator into baby's nostril and tilt it to the side. Never point the tip directly in the center of the nostril. Be careful not to over-insert the aspirator.
 - Slowly release the bulb allowing the suction to draw the mucus out of baby's nostril.
 - Repeat with other nostril.
 - Repeat as necessary until baby's experiences relief from congestion.
 - To clean, put the aspirator under soapy, warm water. Squeeze and release the bulb several times to clean the inside of the bulb. Rinse with warm water. Squeeze out all excess water from the bulb.
-

PACIFIER MEDICINE DISPENSER:

The Medicine Dispenser provides parents with a convenient way to dispense medicine to their child. The unique design reduces the likelihood of spillage and helps ensure proper dosage. The soft, contoured applicator is comfortable for baby.

To use: Boil in water for 2 minutes before first use. Twist the clear cap counterclockwise and pivot open. Pour in the desired amount of liquid using the dosage marks. Holding the cap upright, pivot the pacifier shield down on the cap, and then twist clockwise to close. Place the Summer medicine dispenser into baby's mouth. Works best when the child is in a reclined position.

To clean: Immediately after use, open the cap and wash with hot soapy water. Rinse thoroughly. Dishwasher safe.

TEETHER:

To clean: Rinse under clean, warm water before first use and after each use.

Note: Do not boil.

FOREHEAD THERMOMETER & CASE:

This temperature indicator is reusable and comes with a convenient carrying case. The parent simply holds it in place on the child's forehead and gets a reading in less than 15 seconds.

EMERGENCY CARD:

- Keep important phone numbers and medical information handy for quick reference.
- Keep in a safe place. Familiarize family members, baby-sitters and care givers with the location of the card and the information on it.

STORAGE CASE:

- Keep items organized with the convenient storage case.

To clean: Surface wash with mild soap and warm water. Rinse thoroughly; air dry.

DIGITAL THERMOMETER:

For complete instructions for use and care, please see separate instruction manual enclosed.

FINGERBRUSH & CASE:

This unique item fits compactly into its own sanitary case. Whether a parent is massaging sore teething gums on the go, or maintaining a good oral brushing hygiene after eating out, this wonderful little product fills the need.

- Clean infant finger toothbrush before first use and after each use.
 - Place on parent's fingertip.
 - Gently clean and massage gums for emerging teeth.

GUM MASSAGER:

- Clean gum stimulator before first use and after each use.
- As your child's first tooth emerges: massage gums to ease tooth emergence and relieve discomfort.

We at Summer Infant stand behind all of our products. If you are not completely satisfied or have any questions, please contact our Customer Service Department at 1-800-268-6237 or e-mail us at customerservice@summerinfant.com.

CE 00434 Forehead thermometer strip manufactured by:
Hallcrest, Inc.
Glenview, IL USA

CE 0044 Digital thermometer manufactured by:
ONBO Electronic (Shenzhen) Co., Ltd.
China

Colors and styles may vary.

**Forehead thermometer strip and alcohol swabs
MADE IN USA.**

**All other products
MADE IN CHINA.**

Please retain information for future reference.

Summer Infant, Inc.
1275 Park East Drive
Woonsocket, RI 02895 USA
1-800-268-6237

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STANDARD DIGITAL THERMOMETER

14374, 14444, 14454, 14464, 14474

Read the instructions carefully before using this device.

IMPORTANT NOTES:

Performance of this device may be degraded if it is operated outside of stated temperature humidity range, stored outside of stated temperature/humidity range, dropped or subjected to mechanical shock, and/or patient temperature is below ambient temperature.

- High,prolonged fever requires medical attention.Be sure to contact your physician.
- For parental use.Do not leave child unattended with thermometer.
- Avoid exercise and drinking hot or cold beverages before taking a reading.
- Clean the probe tip with isopropyl alcohol or mild soap and warm water
Clean before first use and after each use.
- Thermometer is water resistant.
- Store thermometer in protective case when not in use.
Store at room temperature and out of direct sunlight.


HOW TO TAKE TEMPERATURE:

Keep in mind that taking a temperature can be trying for both parent and child.
When taking a reading be sure to remain calm and confident.The sound of your voice can be very soothing for an infant.Explaining what you are about to do and why can be helpful for a toddler or older child.

- Remove thermometer from protective case.Be sure to clean thermometer before first use and after each use.
- Push the ON-OFF button located on the large end of the thermometer.The display will read:“188.8 $\frac{\text{°F}}{\text{°C}}$ ”.This is a function check indicating the unit is working properly
If a reading is in memory it will be displayed briefly.When the °F or °C begins to flash,the thermometer will begin measuring temperature when exposed to body heat.
- The °F (or °C) symbol will continuously flash until completion of the temperature taking process.Please note that the “ ∇ ”indicates a low battery.
An “ERR”message in the window indicates your unit is not functioning properly. Do not expect an accurate temperature reading if this low battery indicator is present.


ORAL USE: (for children 4 years and older)

- Do not drink hot or cold fluids,exercise or perform other activities that will raise or lower temperature readings when compared to your child’s normal,average temperature.The mouth should remain closed up to 5 minutes before attempting a reading.
- Place the probe tip well under the tongue as indicated by the “★” (see Figure A) Instruct children not to bite thermometer.
- The peak temperature should be reached in approximately 60 seconds.Opening the mouth or improper placement of the probe tip could result in a longer time for a reading.
- Follow numbered instructions below.

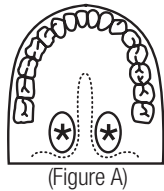UNDERARM (AXILLARY) USE:

- Wipe underarm with a dry towel.
- Place probe tip in armpit and bring arm down against the body. (Figure B)
The probe should be pointing upward in the same direction as your body if possible.
- The time required for an accurate underarm temperature may be 90 seconds or longer depending on placement of the probe and closure of the arm against the body.
- Important Note:An underarm reading will be approximately 1 degree lower than an oral reading.To calculate the oral equivalent add 1°F or 0.5°C to the underarm temperature reading.

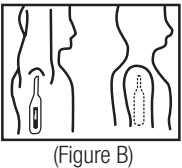RECTAL USE: (Infants and small children)

- Commonly used for babies,young children or when it is difficult to take an oral or underarm temperature.
- Lubricate probe tip with a water-soluble jelly or petroleum jelly.
- For easier insertion, position baby on his tummy, keeping his legs at a 45°angle to his back.It may be helpful to position the baby on your lap. Spread baby’s buttocks with your thumb and forefinger. Gently insert probe a MAXIMUM of 1/2 in or 1.3 cm into baby’s rectum and hold thermometer lightly in place. Do not force probe tip if you encounter resistance.
- The peak temperature should be reached in approximately 40 seconds.
- Follow numbered instructions below.
- Important Note: A rectal reading will be approximately 1 degree higher than an oral reading.
To calculate the oral equivalent subtract 1°F or 0.5°C from the rectal temperature reading.
- Once you have used the thermometer to take a rectal reading,do not take any oral readings,for hygienic reasons.

BATTERY INFORMATION:

When the “ ∇ ” appears in the lower right of the display, the battery is exhausted. Dispose of the thermometer as the battery is non-replaceable
Battery Contains Mercury. Do Not Put in Trash. Recycle or Manage as Hazardous Waste. Properly dispose of product according to your local regulations concerning disposal of this product.
The battery is well sealed inside the thermometer to reduce the risk of being swallowed by children. Keep away from small children and heat.

Storage: Room temperature 10 °C to 40 °C /50 °F to 104 °F. Humidity: 15% to 95% (noncondensing)

Operating temperature: -25 °C to 60°C/ -13°F to 140°F. Humidity: 15% to 95% (noncondensing)

WARRANTY

This Thermometer is guaranteed for a period of 24 months from date of purchase. This warranty does not cover batteries, or damage resulting from improper handling or use. Specific legal rights stipulated in your state may guide this warranty. If you have a legitimate claim under this warranty, please return the product to the place of purchase along with your sales receipt. If you have any questions or comments please contact our customer relations department.

Digital thermometer
manufactured by:
ONBO Electronic
(Shenzhen) Co., Ltd.
China

Type BF applied part

Batteries must be disposed of in
accordance with the locally applicable
regulations, not with domestic waste.

Distributed by:
Summer Infant, Inc.
1275 Park East Drive
Woonsocket, RI 02895 USA
1-800-268-6237
© 2005 Summer Infant, Inc.

Summer Infant Europe, LTD.
1st Floor North Wing, Focus 31
Cleveland Road
Hemel Hempstead, Herts
HP2 7BW UK
0144 250 5000

visit us at
www.summerinfant.com
Please retain information
for future reference.
Colors and styles may vary.
MADE IN CHINA.

RADIO AND TELEVISION INTERFERENCE

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to connect the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

You may also find helpful the following booklet, prepared by the FCC: “How to Identify and Resolve Radio-TV Interference Problems.” This booklet is available from the U.S. Government Printing Office, Washington D.C. 20402.

Changes and Modifications not expressly approved by Summer Infant, Inc. can void your authority to operate this equipment under Federal Communications Commissions rules.

This device complies with part 15 of FCC rules. Operation is subject to the following two conditions; (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

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Termómetro digital fabricado por: ONBO Electronic (Shenzhen) Co., Ltd. China