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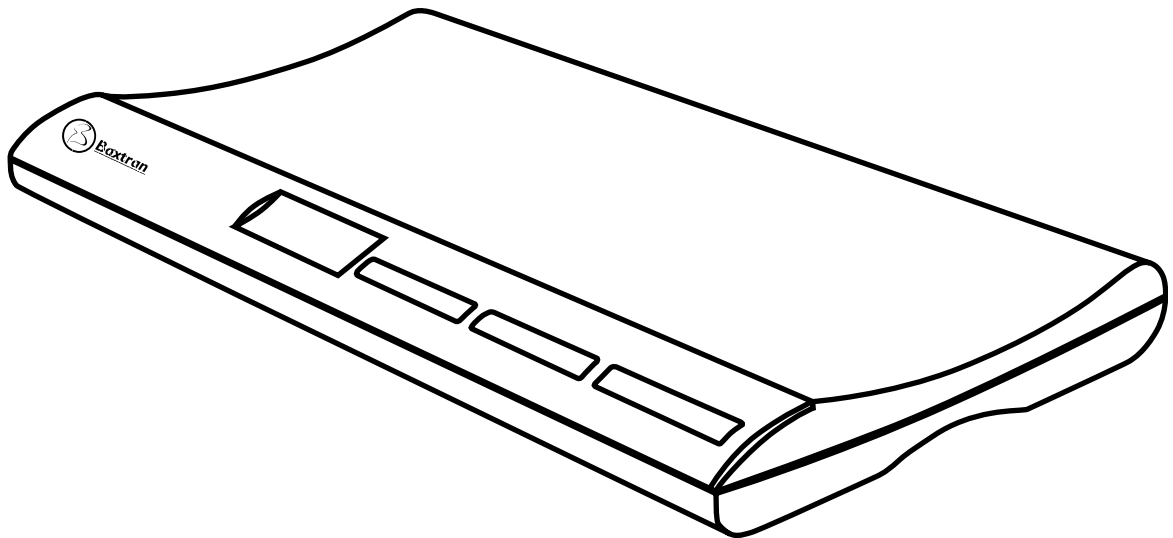
The electronic BAC baby scale has been designed to help you monitor the growth of your baby, accurately and reliably quantifying how much milk has been assimilated at each feed.

The TARE function allows you to weight the baby using a soft cloth or towel.

The scale come with a useful WEIGHT-BLOCK function that allows a steady weight to be displayed despite the child wriggling about.

Taking up to 20 kg, and with 10 g divisions, the scales will also help you monitor your child's progress in terms of weight gain during all stages of growth.

YOU ARE STRONGLY RECOMMENDED TO READ THE INSTRUCTIONS AND WARNINGS HEREIN BEFORE USING THE EQUIPMENT AND TO KEEP THE BOOKLET IN A SAFE PLACE.

TIPS FOR USE

Place the scale on a flat and hard surface. A soft and uneven floor (e.g. carpets, linoleum) could cause an alteration to the weight reading.

INSTALLING BATTERIES

Use four 1.5V alkaline batteries only.

Remove the batteries cover on the underside of the scale.

Insert the batteries and refit the cover.

Do not throw the batteries away together with your normal domestic rubbish. Take them to a used battery bin.

FUNCTION KEYS

ON/OFF

The ON/OFF key turns the scale on and off.

SELECT UNIT (SEL. UNITA')

Press this key to select the unit of measure you require: kg, lb or oz. These scale has a capacity of 20 kg, division 10 g. After turning the scale on, display shows 0.00 kg / 0.00 lb /0.00 oz.

HOLD (BLOCK WEIGHT FUNCTION)

Due to the baby moving, the weight shown on the display changes all the time. With the scale switched on, push the button "HOLD"; on bottom of the display left side you see a black round (●) icon flickering. Lay the baby on the scale and wait a few seconds till weight indication on the display flickers twice.

The display keeps weight indication for some seconds, then the scale switched off automatically.

TARE

If the baby is laid on a towel or cover, use the TARE function.

Turn the scale on with the ON/OFF/TARE key. The display shows "8888". Wait until the display shows 0.00 kg / 0.00 lb /0.00 oz, then place the towel or cover on the plate.

Press the ON/OFF/TARE key and wait until the display shows 0.00 kg / 0.00 lb /0.00 oz; the scale is now ready for use.

USING THE SCALE

Press the ON/OFF/TARE key. The display shows "8888". Wait until the display shows 0.00 kg / 0.00 lb /0.00 oz; the scale is now ready for use. Lay the baby on the plate and wait until the display shows its weight. Press the ON/OFF/TARE key to turn the scale off, otherwise the scale turn off automatically after 20 seconds without being used.

MAINTENANCE

- Take care of your scale, as it a precision instrument. Therefore, avoid exposing to extreme temperatures, damp, knocks, dust and direct sunlight. Do not drop or knock the scale.
- Keep the scale away from sources of heat.
- Do not jump on the scale.
- Clean the platform with a damp cloth.
- Clean the bottom of the scale with just a soft dry cloth.
- Never use chemical or abrasive products.
- Take care that no liquids penetrate into the scale, as it is not waterproof. Do not leave or use the scale in very damp surroundings, humidity above 85% or contact with water and other liquids could cause corrosion and effect the precision of the scale.

SYMBOLS

- "LO" indicates that the batteries are flat. Replace them.
- "ERR" shows the weight is too much.
- The scale is ready for use when the display shows "0.00"; if the instructions have not been followed correctly, the display shows "ERR" or "----".


TECHNICAL DATA

- Range of measurement: **0.5 - 22 kg**
- Division: **10 g**
- Battery: **4 alkaline batteries type 1.5V AAA LR03 series.**
- Working conditions: **0°C + 50°C RH 85%**
- Storage conditions: **0°C + 50°C RH 85%**
- Dimensions: **520 x 320 x 56 mm.**
- Accuracy and repeatability with tolerance level within: **5%.**

¡Attention! Read carefully the instructions for use.

DISPOSAL PROCEDURE (DIR. 2002/96/CE - RAEE)

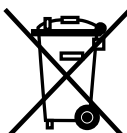The symbol on the bottom of the device indicates the separated collection of electric and electronic equipment. At the end of life of the device, do not dispose it as mixed solid municipal waste, but dispose it referring to a specific collection centre located in your area or returning it to the distributor, when buying a new device of the same type to be used with the same functions.

This procedure of separated collection of electric and electronic devices is carried out forecasting a European environmental policy aiming at safeguarding, protecting and improving environment quality, as well as avoiding potential effects on human health due to the presence of hazardous substances in such equipment or to an improper use of the same or of parts of the same.

Caution! The wrong disposal of electric and electronic equipment may involve sanctions.

WARRANTY

The warranty covers only any faults in manufacturing and does not cover accidental damage, incorrect use, negligence or misuse of the product. Use only the accessories provided; the use of accessories other than the approved ones can make your warranty void.

Never open the appliance; in the event of any opening or tampering the warranty will automatically become invalid.

The warranty does not cover either the parts subject to wear and tear or the batteries (if provided). When the warranty expires any technical service shall be charged. Information regarding technical service either under the terms of this warranty or to be paid for can be obtained by contacting giropes@giropes.com. No compensation of any kind is owed for repairs or substitutions of products under the terms of this guarantee. In the event of faults, go back to the retailer; DO NOT send the appliance directly to Giropes. Use in public places is not allowed.

The manufacturer shall not be liable for any damages to people, property or pets, either directly or indirectly, caused by the non compliance with the instructions laid down in the instructions manual, in particular the warnings regarding the installation, use and maintenance of the appliance.

It is Giropes policy to improve its products constantly; Giropes therefore reserves the right to change its products, either partially or totally, without previous notification and according to the manufacturing requirements, without either Giropes or its retailers being in any way liable for it.

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trade mark propriety of: