Oral-B®

ProfessionalCare™

8900 DLX OxyJet® Center

powered
by





BRAUN

English	4
Español	15
Français	24

Type 3719
Modelo 3719
Modèle 3719

3-719-135/00/II-05
USA/CDN

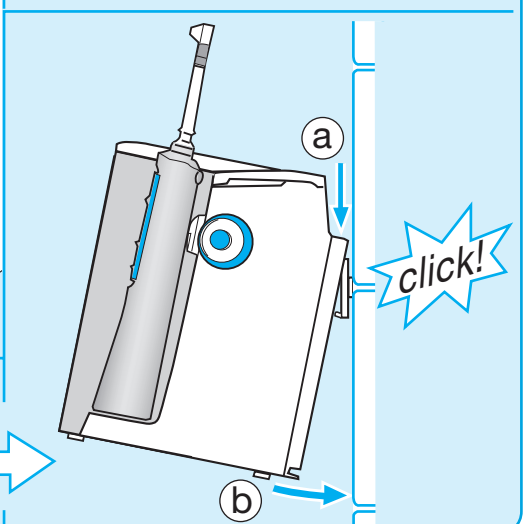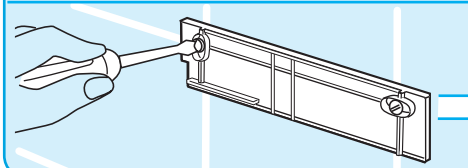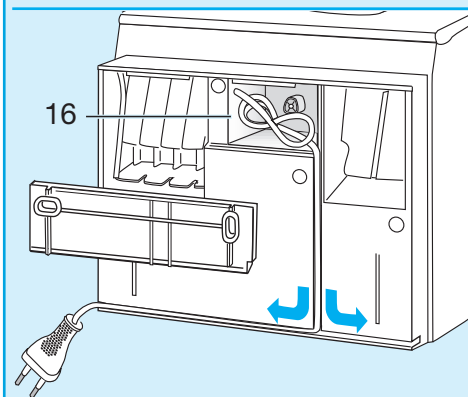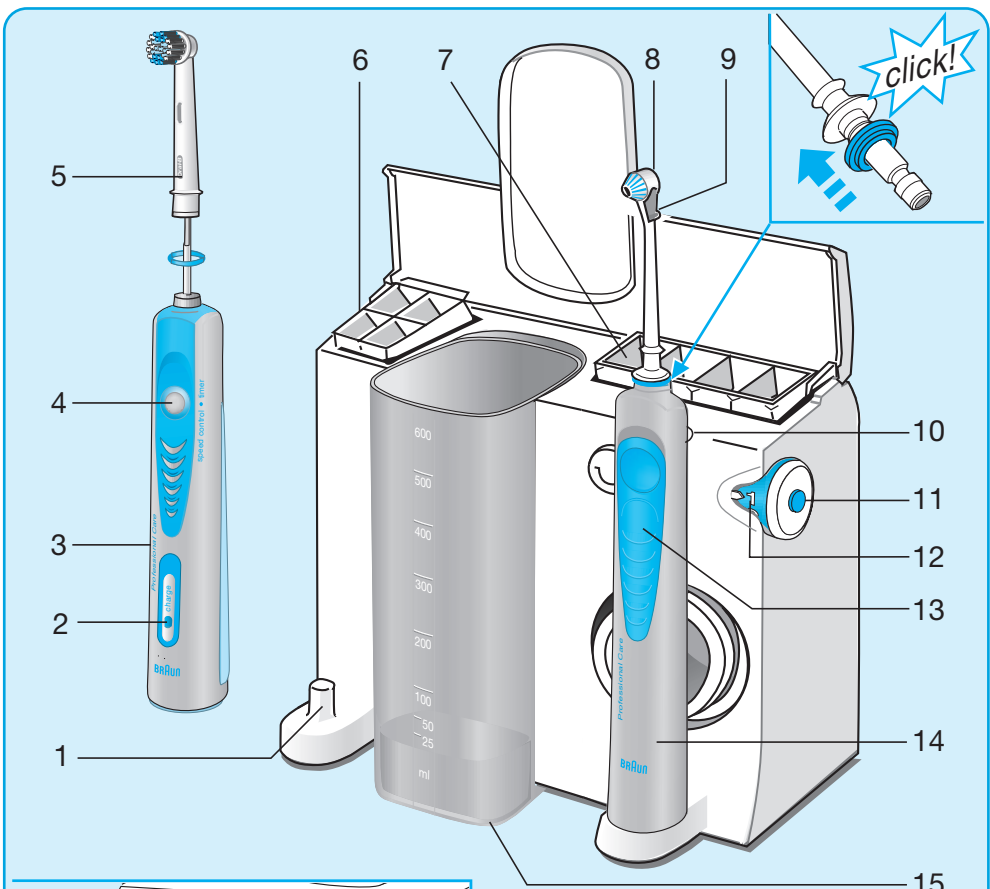IMPORTANT SAFEGUARDS

Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center has been carefully designed to offer you and your family a brushing experience that is both safe and effective.

However, when using electrical products, basic safety precautions should always be followed. **Please read carefully before use.**

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by or near children.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop or insert any object into any opening.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

SAVE THESE INSTRUCTIONS

The **Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center** is a dentist-inspired oral-care system for outstanding cleaning, gum care, whiter teeth and breath freshening.

The **Oral-B® ProfessionalCare™ 8900 DLX** toothbrush delivers outstanding performance:

- Patented 3D oscillating-pulsating technology – for outstanding cleaning and gum care
- Innovative accessories and brushheads to clean and polish teeth as well as freshen your tongue and mouth
- Customized Speed Control which offers a broad variety of speed levels to meet your individual brushing needs
- Premium handle design and ergonomics for greater comfort and control

The **Oral-B® ProfessionalCare™ OxyJet®** irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water. The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning appliance should no longer be used. Take it to a Braun Service Center.

Description

Toothbrush

- ① Charger base
- ② Charge indicator
- ③ Toothbrush handle (rechargeable)
- ④ On/off switch
- ⑤ Brushhead (see «Brushheads and accessories»)
- ⑥ Brushhead compartment

Irrigator

- ⑦ Jet compartment
- ⑧ Irrigator jet
- ⑨ Jet switch
- ⑩ Jet release button
- ⑪ On/off switch
- ⑫ Water pressure control dial
- ⑬ Handle switch (water stop and go)
- ⑭ Irrigator handle
- ⑮ Water container (20 fl. oz 600 ml)
- ⑯ Cord storage

Specifications

Power supply: see bottom of housing

Power intake: 24 W

Using your Oral-B® ProfessionalCare™ 8900 DLX toothbrush

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and can be used in the bathroom without hesitation.

- Plug the appliance into an electrical outlet and place the toothbrush handle ③ on the charger base ①.
- The blue charge indicator ② flashes while the handle is being charged. If you then use your brush and place it back on the charger base, the charging light will again flash to let you know it is charging.
- When the handle is fully charged, the charge indicator shines continuously. A full charge takes at least 16 hours and enables approximately two weeks of regular brushing.
- When the rechargeable battery is running low, your toothbrush will slow down, then stop.
- For everyday use, the toothbrush handle can be stored on the charger base to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, refrain from placing the handle on the charger base and fully discharge it through regular use at least every 6 months.

Customized Speed Control

Oral-B® ProfessionalCare™ 8900 DLX offers a variable range of speeds to best meet your individual needs. High speed is ideal for tooth and gumline cleaning. For cleaning sensitive teeth and gum tissue as well as precision use of the accessory (III), we recommend using a lower speed.

- Press the on/off button once to operate at maximum speed.
- Continue to hold the on/off button to shift from the highest speed to the Customized Speed Control mode.
- As long as you keep the on/off button pressed, the handle will continuously run through all speed levels.
- Release the on/off button to brush at your preferred speed level.
- Push the on/off button once more to turn the toothbrush off.
- If you then turn it on again you will automatically come back to the maximum speed level.

Brushing technique

- Attach one of the brushheads to the handle. They can be used with any kind of toothpaste.
- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Guide the brushhead slowly from tooth to tooth.
- Hold the brushhead in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Pressure Sensor

For optimum brushing results **Oral-B® ProfessionalCare™ 8900 DLX** has a pressure control feature to ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brushhead will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.



Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a better overall clean.

2-Minute Timer


A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

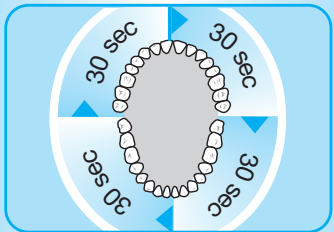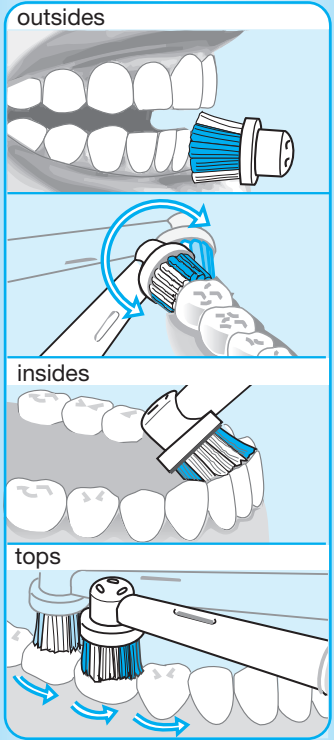
Brushheads and accessories

The **Oral-B® ProfessionalCare™ 8900 DLX** offers you a variety of brushheads and accessories, allowing you to select those which best meet your individual needs. They are all designed for daily use:

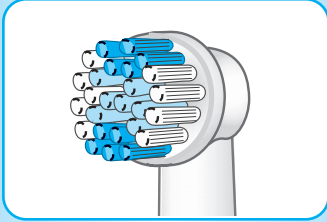I FlexiSoft® brushhead


The FlexiSoft® brushhead features our outstanding oscillating-pulsating action to effectively clean your teeth, along the gumline and in hard-to-reach areas such as behind teeth. The green FlexiSoft® bristles bend, enabling both gentle brushing and deeper reach of the blue Interdental Tips between teeth.

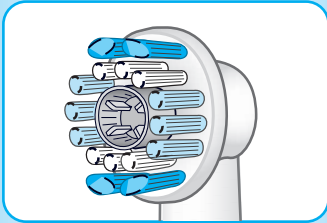
II Oral-B PowerPolisher™ brushhead

The Oral-B PowerPolisher™ brushhead is specifically designed to polish and naturally whiten your teeth by removing coffee, tea and tobacco stains. Although you can use any kind of toothpaste, a whitening paste is likely to result in a better whitening performance. We do not recommend using the Oral-B PowerPolisher™ brushhead with braces.

Apply the same brushing technique as described in the section «Brushing technique» for both brushheads.

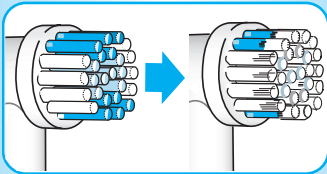Both brushheads feature blue INDICATOR® bristles to help you monitor brushhead replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway within 3 months, signaling it is time to change the brushhead.

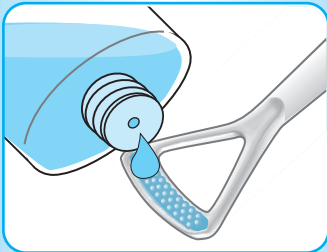III Tongue Freshener


The Tongue Freshener freshens breath by removing odor-causing bacteria from the tongue. It can be used every day – before or after brushing. You can also use mouthwash with the Tongue Freshener.

Attach it to the handle, turn the appliance on and drag it across the tongue from back to front. We recommend a low speed setting for enhanced control during use. Rinse after use. Replace the Tongue Freshener when materials become worn.

Oral-B® ProfessionalCare™ 8900 DLX replacement brushheads and accessories are available in stores that sell Oral-B products, at the Essentials.com web site and at Braun Service Centers.

After use

- Rinse your brushhead or accessory thoroughly under running water for several seconds with the handle switched on.
- Switch off the handle and remove the brushhead or accessory from the handle. Clean all parts separately under running water, then wipe them dry before you reassemble the brush and put the handle on the charger base.

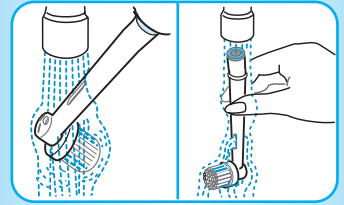
Environmental notice

This product is provided with a nickel-hydride battery, which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre or at approved recycling or disposing locations in accordance with local regulations.


Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.


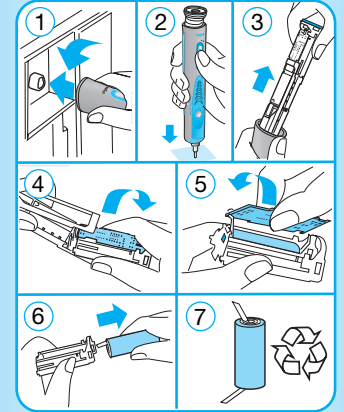
Using your Oral-B® ProfessionalCare™ OxyJet® irrigator

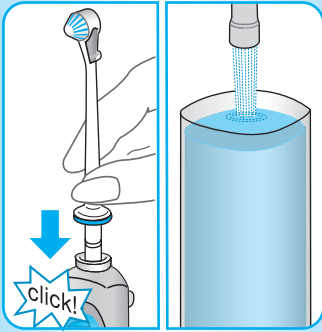Before using for the first time, fill the water container ⑮ and allow the water to run through completely by holding the irrigator handle ⑭ over a washbasin.


1. Attach one of the irrigator jets ⑯ to the handle until it locks into place. Fill the water container with water and add mouthwash, if desired.
2. Turn on the irrigator by pressing the on/off switch ⑩. (When doing so, the handle switch ⑪ should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial ⑫ (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch ⑬ up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.

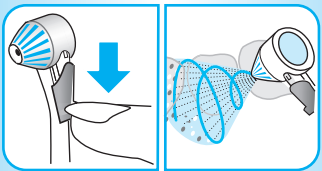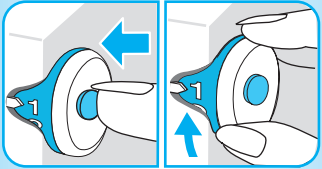You can choose between two different jet functions by sliding the jet switch ⑩ down or up (before doing so, turn off the water stream by sliding the handle switch down):


Rotating stream: for general cleaning and gum massage (jet switch in lower position). The rotation of the water stream and its micro-bubbles is designed for easy penetration.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position).

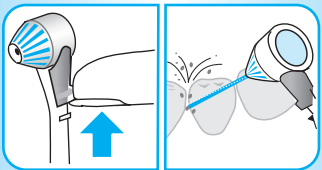Max. operating time: 15 minutes
Cooling off time: 2 hours

To ensure the jet functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

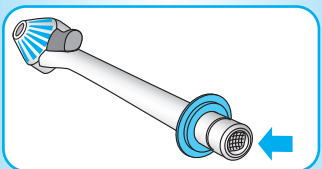After use

Always empty the water container entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.


Press the on/off switch to turn off the irrigator. It will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button [®]. The jets can be stored in the jet compartment. Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

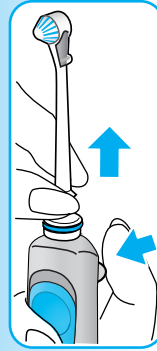FOR USA ONLY

Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Contact an authorized Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative:

Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts & accessories:

Contact an authorized Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

To speak to an Oral-B Service Representative:
Please call 1-800-566-7252.

The 30-Day Challenge

Take the challenge! Experience the **Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center** difference!

Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center is a complete oral care system for superior cleaning, natural whitening, gum care and freshening of your whole mouth.

No wonder that Oral-B is the brand more dentists use themselves worldwide. In fact, we are so convinced you'll fall in love with your new **Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center** that we give you this risk-free offer.

Try **Oral-B® ProfessionalCare™ 8900 DLX OxyJet® Center** for 30 days. If you are not 100% satisfied, receive a full refund. No questions asked. Guaranteed.

For a full refund, please return your handle, refills, and charging unit with your original cash register receipt to the address provided below. For further details, please call 1-800-566-7252 to be referred to the authorized Braun Service Center closest to you.

U.S. only:

Attention: 30 Day Challenge
Braun Division
The Gillette Company
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Canada only:

Attention: 30 Day Challenge
Braun Products
Gillette Canada Company
4 Robert Speck Parkway
Mississauga ON L4Z 4C5

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