Also try these other great Oral-B  
refills with this Vitality toothbrush!

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Dentists recommend changing your  
brushhead every 3 months.

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


D 12513
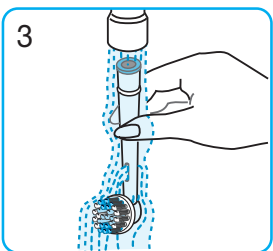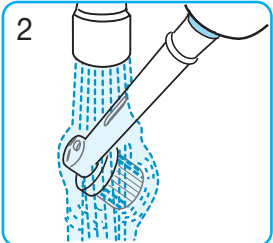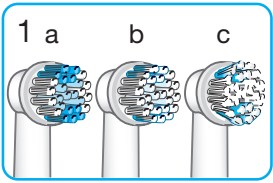
# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed.

Read all instructions before using.

## DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.

## WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. If any of these occur, please contact an Oral-B Braun Service Center and return the product for examination and repair.

# SAVE THESE INSTRUCTIONS

---

## Description

- A Brushhead (depending on model) – interchangeable\*
- B On/off switch
- C Handle
- D Charging unit

Please note: We do not recommend using the Oral-B Pro White brushhead with braces.

## Charging

- Plug the charging unit (D) into an electrical outlet. Place the handle on the charging unit. Handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

## Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

## Using the toothbrush

- Wet brushhead and apply toothpaste.
- Place brushhead on teeth, then switch on.
- Slowly guide the brushhead from tooth to tooth, first the outsides, then the insides, finally the chewing surfaces. Do not press too hard or scrub.
- After two minutes (brushing time recommended by dentists) the timer will signal by a short stutter.
- During the first days of use, gums may bleed slightly. Should bleeding persist after 2 weeks, contact your dentist.

\* All Oral-B replacement brushheads can be used on the Oral-B Vitality handle, except for Oral-B Sonic Complete and Oral-B CrossAction Power brushheads.

- Replace brushhead every 3 months, when bristles are worn (1c) or when color on Indicator bristles (1b) has faded half-way down the bristles.


## Cleaning

After use, rinse the brushhead under running water (2). Remove the brushhead and rinse brushhead and handle separately (3), then dry.

## Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center.

However, if you prefer to dispose of the battery yourself, open the handle as shown below:

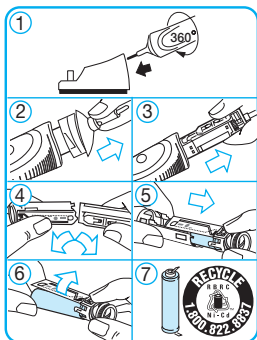Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

## FOR USA ONLY

### Limited 2-Year Warranty

Oral-B Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you

specific legal rights, and you may have other legal rights which vary from state to state.

## **How to obtain service in the USA**

Should this product require service, please call **1-800-566-7252** to be referred to the authorized Oral-B Braun Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B Braun replacement parts and accessories:*

Contact an authorized Oral-B Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

*To order original Oral-B Braun replacement parts and accessories (including replacement brushheads) by phone:*

Please call 1-800-566-7252.

*To speak to an Oral-B Braun Consumer Service Representative:*

Please call 1-800-566-7252.

## **FOR CANADA ONLY**

### **Limited 2-Year Warranty**

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair

or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

### **How to obtain service in Canada**

Should this product require service, please call **1-800-566-7252** to be referred to the authorized Oral-B Braun Service Center closest to you.

*To obtain service:*

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

*To obtain original Oral-B Braun replacement parts & accessories:*

Contact an authorized Oral-B Braun Service Center. Please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

*To speak to an Oral-B Braun Service Representative:*  
Please call 1-800-566-7252.

## 12 - Month Performance Guarantee

Improved oral care does not happen overnight. At Oral-B, we recognize that fact. That's why we offer you 12 months to adjust to the new Oral-B® Vitality™. With daily use, you will see how the Oral-B Vitality significantly reduces plaque. If you don't, send it back for a full refund.

That's our promise to you – if you are not happy with the results, then we aren't either. To help us provide you with your refund, please send us your Oral-B Vitality within 12 months of purchase along with your name and address and the original cash register receipt indicating the purchase price and date of purchase.

### **USA:**

Oral-B Performance  
Guarantee  
Oral-B Braun Division  
The Procter & Gamble  
Company  
Sobin Park and A Street, MS 1F  
South Boston, MA 02127

### **CANADA:**

Oral-B Performance  
Guarantee  
Braun Canada  
a division of Procter &  
Gamble Inc.  
Toronto, ON  
M5W 1C5

Sorry: No refunds will be made if your Oral-B Vitality was damaged by accident, neglect, unreasonable use or is lost in the mail.

The Oral-B Vitality must be sent in its original packaging and be postmarked no later than 12 months after the date of purchase.

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