BRAUN

Oral-B

D 9525

D 9525 T

D 9511

D 9513

D 9521

ultra
Plak Control
power toothbrush

Type 4728

Type 4725

Braun Infoline

English Thank you for purchasing a Braun Oral-B plaque remover. If you have any questions, please call:

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












Printed in Germany

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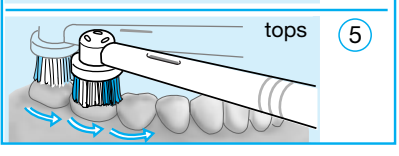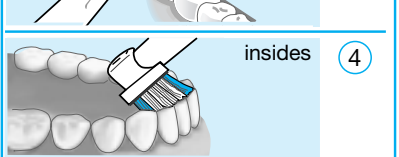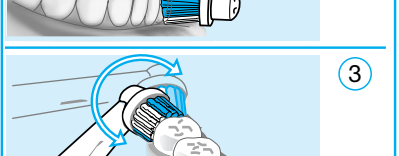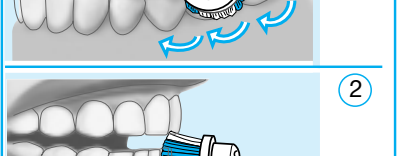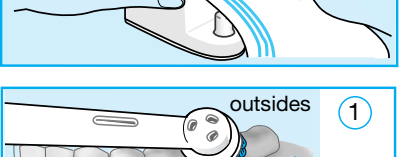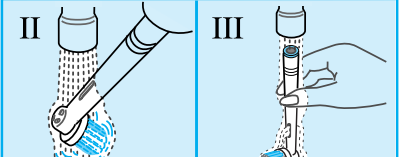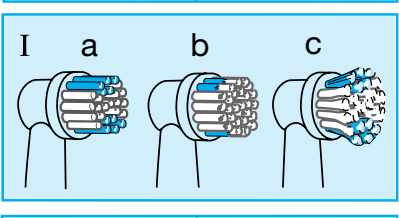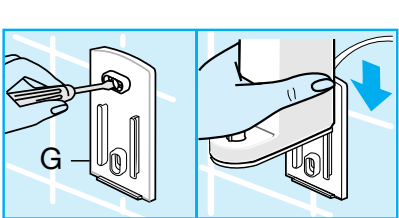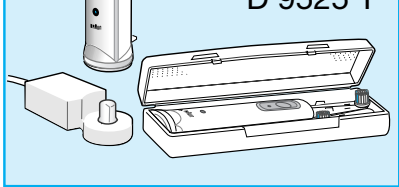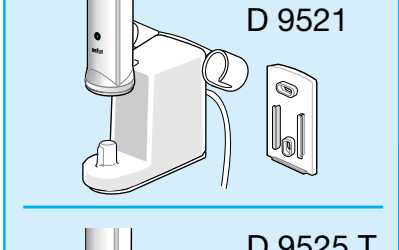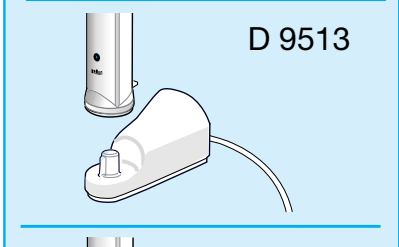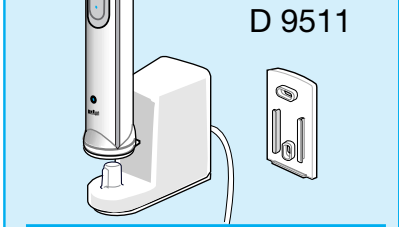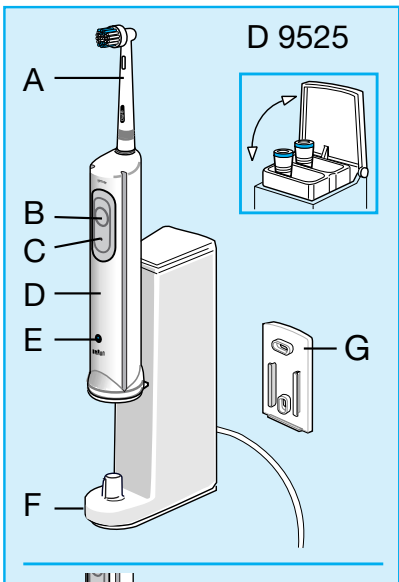IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER

To reduce the risk of electrocution:

1. Do not place or store charger where it can fall or be pulled into a tub or sink.
2. Do not place charger in or drop into water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never drop or insert any object into any opening.

SAVE THESE INSTRUCTIONS

The Braun Oral-B Plak Control Ultra has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth. It has been designed for daily use for the whole family.

Please read the use instructions carefully before first use.

Important

Periodically check the cord for damage. Should this be the case, take the charging unit (F) to a Braun Service Center.

A damaged or non-functioning unit should no longer be used.

Description

- A Brushhead
- B On switch
- C Off switch
- D Handle
- E Charge indicator
- F Charging unit
- G Wall holder

Specifications

Power supply: see bottom of charging unit

Handle voltage: 2.4 V

Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit.
- A full charge takes 16 hours and will provide an operation time of approx. 30 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Using the toothbrush

The toothbrush can be used with any leading brand of toothpaste.

- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Slowly guide the brushhead from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

Timer

An incorporated timer memorizes the elapsed brushing time. After two minutes the timer signals by a short stutter.

During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days.

Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months. Therefore, the brushhead is provided with INDICATOR® bristles (I.a). With proper brushing twice a day for two minutes, the color will disappear half-way down the bristles in a 3-month period (I.b), provided you use toothpaste.

If the bristles splay before the color recedes (I.c), you are exerting too much pressure on teeth and gums.

Replacement brushheads

This product can be used with all Braun Oral-B brushheads.

Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

From time to time, also clean the charging unit with a damp cloth (IV).

Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Center.

However, if you prefer to dispose of the battery yourself, please refer to the section «Battery removal».

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FOR USA ONLY

Limited Warranty

Braun warrants this Braun appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an Authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an Authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes, normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

To obtain service in the USA

Should your Braun product require service, please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A Take the product in to the authorized Braun Service Center of your choice or,
- B Ship the product to the authorized Braun Service Center of your choice.
 - Pack the product well.
 - Ship the product prepaid. Insurance is recommended.
 - Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Replacement parts and accessories are available at authorized Braun Service Centers.

Please call 1-800-BRAUN 11 (272-8611) to be referred to the Authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-892-8534.

To speak to a Braun Consumer Service Representative:

Please call 1-800-984-9668.

FOR USA ONLY

Nickel-Cadmium Battery Pack

SAVE THIS INFORMATION

PRODUCT CONTAINS A NICKEL-CADMIUM RECHARGEABLE BATTERY PACK. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

May not be disposed of in municipal waste. Recycle or dispose of in an environmentally sound manner.

RBRC

Ni-Cd

Ni-Cd

Call Braun 800-272-8611 for recycling information.

See next page for easy battery pack removal.

DISPOSE OF YOUR BATTERY PACK PROPERLY.


DO NOT PLACE USED BATTERIES IN YOUR REGULAR TRASH.

THIS NICKEL-CADMIUM BATTERY PACK MUST BE COLLECTED, RECYCLED OR DISPOSED OF IN AN ENVIRONMENTALLY SOUND MANNER.

The incineration, landfilling or mixing of nickel-cadmium batteries with the municipal solid waste stream is **PROHIBITED BY LAW** in most areas.

Return this battery pack to a federal or state approved nickel-cadmium battery recycler. This may be where you purchased the battery.

In MINNESOTA call 1-800-225-PRBA if further disposal information is required; or Contact your local waste management officials for other information regarding the environmentally sound collection, recycling and disposal of this battery pack.

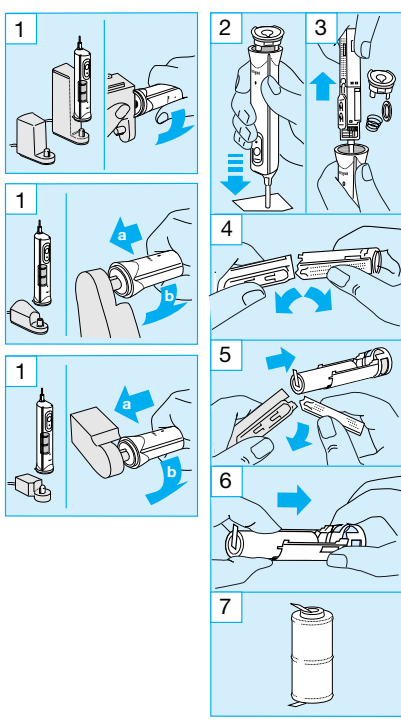English

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution: Opening the handle will destroy the appliance and invalidate the guarantee.

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PLAQUE GUARANTEE

For USA only

Use the Braun Oral-B Plak Control daily for at least three weeks until you are fully familiar with this mechanical plaque removal instrument. After that time, if you are not convinced that the Braun Oral-B Plak Control significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Braun Oral-B Plak Control within 30 days, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

Plaque Guarantee

Braun Division,
The Gillette Company,
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Sorry: No refunds will be made for any plaque remover damaged by accident, neglect or unreasonable use.

The Braun Oral-B Plak Control must be shipped in its original packaging and be postmarked no later than 30 days after the date of purchase.

Please allow 6–8 weeks for handling.

Name of Regular Toothbrush

Your Name

Street

City

State

Zip Code

PLAQUE GUARANTEE

For Canada only

Use the Braun Oral-B Plak Control daily for at least three weeks until you are fully familiar with this mechanical plaque removal instrument. After that time, if you are not convinced that the Braun Oral-B Plak Control significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Braun Oral-B Plak Control within 30 days, postage-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

Plaque Guarantee

Braun Canada Div.,
Gillette Canada Company,
110 Matheson Blvd. W. Suite 200,
Mississauga, Ontario L5R 3T4

Sorry: No refunds will be made for any plaque remover damaged by accident, neglect or unreasonable use.

The Braun Oral-B Plak Control must be shipped in its original packaging and be postmarked no later than 30 days after the date of purchase.

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For CANADA only

Limited 2-Year Warranty

In the event a Braun appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

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For CANADA only

How to obtain service in Canada

Should your Braun product require service, please call 1-800-387-6657 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.
Pack the product well.
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts & accessories:

Contact a Braun authorized Service Center. Please call 1-800-387-6657 to be referred to the authorized Braun Service Center closest to you.

To speak to a Braun Service Representative:

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