Oral-B®

AdvancePower

D 9535
D 9525
D 9521
D 9513
D 9511
D 9013

powered by **BRAUN**

Type 4728

English

Thank you for purchasing an Oral-B AdvancePower, which is designed and engineered by Braun. If you have any questions, please call:

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






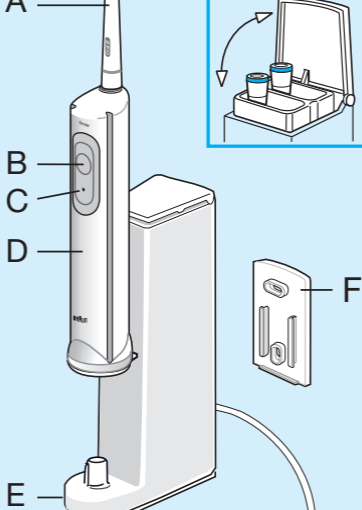
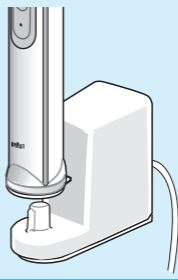
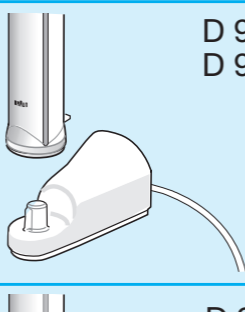
D 9521

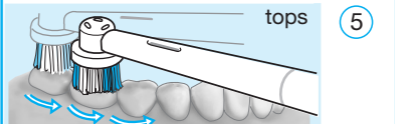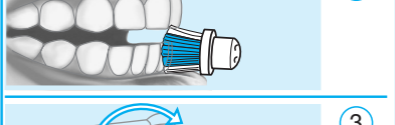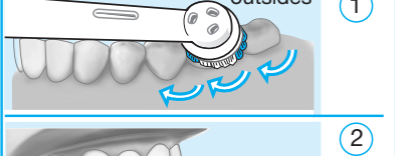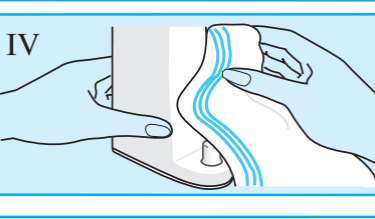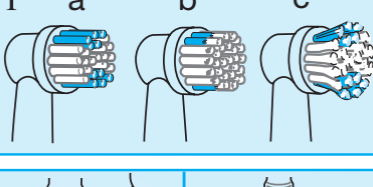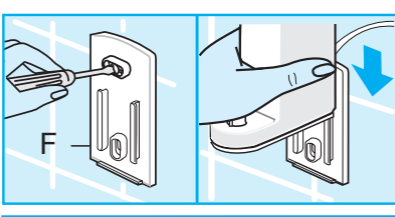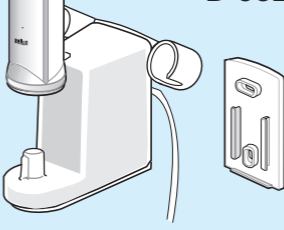IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER

To reduce the risk of electrocution:

1. Do not place or store charger where it can fall or be pulled into a tub or sink.
2. Do not place charger in or drop into water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never drop or insert any object into any opening.

SAVE THESE INSTRUCTIONS

The Oral-B AdvancePower has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth. It has been designed for daily use for the whole family.

Please read the use instructions carefully before first use.

Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Center.

A damaged or non-functioning unit should no longer be used.

Description

- A Brushhead
- B On switch
- C Off switch
- D Handle
- E Charging unit
- F Wall holder

Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit. Handle must be switched off.
- A full charge takes 16 hours and will provide an operation time of approximately 20 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

Using the toothbrush

The toothbrush can be used with any toothpaste.

- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Slowly guide the brushhead from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

During the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days.

Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months. Therefore, the brushhead is provided with INDICATOR® bristles (I.a).

With proper brushing twice a day for two minutes, the color will disappear half-way down the bristles in a 3-month period (I.b), provided you use toothpaste.

If the bristles splay before the color recedes (I.c), you are exerting too much pressure on teeth and gums.

Replacement brushheads

Available at retail stores or Braun Service Centers.

Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

From time to time, also clean the charging unit with a damp cloth (IV).

Environmental notice

This product contains nickel-cadmium rechargeable batteries. In the interest of protecting the environment, do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Center.

However, if you prefer to dispose of the battery yourself, please refer to the section «Battery removal».

Subject to change without notice.

FOR USA ONLY

Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

To obtain service in the USA

Should this product require service, please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A Take the product in to the authorized Braun Service Center of your choice or,
- B Ship the product to the authorized Braun Service Center of your choice.
 - Pack the product well.
 - Ship the product prepaid. Insurance is recommended.
 - Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:
Contact an authorized Braun Service Center. Please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:
Please call 1-800-892-8534.

To speak to a Braun Consumer Service Representative:
Please call 1-800-984-9668.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event a this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1 800 387-6657 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.
Pack the product well.
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts & accessories:

Contact a Braun authorized Service Center. Please call 1 800 387-6657 to be referred to the authorized Braun Service Center closest to you.

To speak to a Braun Service Representative:

Please call 1-800-387-6657.

12 Month Performance Guarantee

For USA only

Use the Oral-B AdvancePower daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B AdvancePower significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B AdvancePower within 12 months, postpaid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

12 Month Performance Guarantee

Braun Division,
The Gillette Company,
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B AdvancePower must be shipped in its original packaging and be postmarked no later than 12 months after the date of purchase.

Please allow 6–8 weeks for handling.

Name of Regular Toothbrush

Your Name

Street

City

StateZip Code

12 Month Performance Guarantee

For Canada only

Use the Oral-B AdvancePower daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B AdvancePower significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B AdvancePower within 12 months, postage-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

12 Month Performance Guarantee

Braun Products
Gillette Canada Company
4 Robert Speck Parkway
Mississauga, ON
L4Z 4C5

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B AdvancePower must be shipped in its original packaging and be postmarked no later than 12 months after the date of purchase.

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PACK. MUST BE RECYCLED

OR DISPOSED OF PROPERLY.

May not be disposed of in municipal waste.

Recycle or dispose of in an environmentally

sound manner.

RBRC

Ni-Cd

Ni-Cd

Call Braun 800-272-8611 for recycling

information.

See next page for easy battery pack

removal.

DISPOSE OF YOUR BATTERY PACK

PROPERLY.

DO NOT PLACE USED BATTERIES IN

YOUR REGULAR TRASH.

THIS NICKEL-CADMIUM BATTERY PACK

MUST BE COLLECTED, RECYCLED OR

DISPOSED OF IN AN ENVIRONMENTALLY

SOUND MANNER.

The incineration, landfilling or mixing

of nickel-cadmium batteries with the

municipal solid waste stream is

PROHIBITED BY LAW in most areas.

Return this battery pack to a federal or state

approved nickel-cadmium battery recycler.

This may be where you purchased the

battery.

In MINNESOTA call 1-800-225-PRBA if


further disposal information is required;

or Contact your local waste management

officials for other information regarding

the environmentally sound collection,

recycling and disposal of this battery pack.

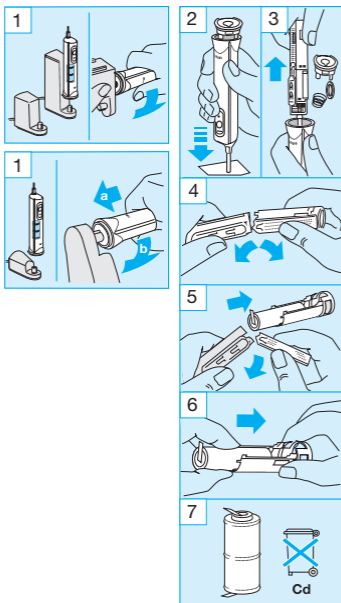English

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

Caution: Opening the handle will destroy the appliance and invalidate the guarantee.

Note: You will not be able to place new batteries in the handle.

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