15 35

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Lase Peroxide Flex⁽¹⁾ is a bleaching gel of Hydrogen Peroxide base, being able to assume the concentration of 15% or 35% of the active principle. The product can be used for the bleaching of vital and not-vital teeth and obligatorily with the gum protection by the application of photo-polymerizer gum barrier, Lase Protect⁽²⁾, or by the use of total isolation with rubber sheet.

It is a product of high bleaching power, that must be used in set with the Whitening Lase II⁽³⁾ or Whitening Lase Light Plus⁽⁴⁾ equipment or an equipment that emits blue light + infra red, with a density of minimum energy of 200mW/cm² and maximum of 350mW/cm².

ITEM THAT COMPOSE THE PRODUCT

Item	Patients		
	1	3	6
Phase 1 Bottle - Hydrogen Peroxide (H ₂ O ₂)	4,5 g	10 g	2 x 10 g
Phase 2 Bottle - Thicking	1,5 g	5 g	2 x 5 g
Phase 3 Bottle - Thinner	1,5 g	5 g	2 x 5 g
Bicarbonated Water Bottle - Neutralizing	5 mL	5 mL	5 mL

Table 01: Item that Compose the Product

The dentist will be able to form a gel with 15% or 35% of hydrogen peroxide. Follows below the table with the item amount that must be inserted to get the concentrations mentioned above.

Concentration	Phase 1 (Drops)	Phase 2 (Drops)	Phase 3 (Drops)
15 %	1	1	1
35 %	3	1	Not applicable

Table 02: Concentration Attainment of 15 % or 35 %

PRODUCTS COMPOSITION

- Mixture of Phases 1, 2 and 3 (15% hydrogen peroxide): hydrogen peroxide, thickening, nano-particled catalyser, neutralizing, coloring, chelator agent and purified water;
- Mixture of Phases 1 and 2 (35 % hydrogen peroxide): hydrogen peroxide, thickening, nano-particled catalyser, neutralizing, coloring, chelator agent and purified water;
- Bicarbonated Water: sodium bicarbonate, preservative and purified water.

USE AND ACTION PURPOSE

Lase Peroxide Flex⁽¹⁾ is indicated for processes of dental bleaching of vital and not vital teeth. This possesses in its formularization hydrogen peroxide, powerful whitening agent that acts oxidating organic pigments (and some not-organic) that are such in the enamel, as in the dentine.

As the bleaching process understands an initial stage of peroxide penetration in the dental structure, a time dependent peroxide concentration gradient is formed. The

pigmentations that concentrate in the dental enamel are removed with bigger rapidity than the existing ones in the dentine.

STORAGE AND TRANSPORT

Keep in place protected of direct light, heat and far from the children reach. Store the product between 10 °C and 25 °C.

CARES WITH THE PATIENT

Cares Before the Bleaching

- To make a good anamnesis, raising the aspects: odontological and medical history (if he is in treatment or use of medicines, etc.), age, habits, vices, teeth color, etc;
- The patient must be with the dental surface perfectly clean;
- The gum must be healthy, without inflammation and tartar;
- Prophylaxis and tartar removing, at least 3 days before the bleaching;
- In case that the patient uses bactericidal buccal rinses, to suspend at least 2 days before the bleaching;
- In case that the patient has spontaneous sensitivity in teeth, treat it before the bleaching;
- In case that the teeth to be bleached present gum retraction, isolate the area with the Lase Protect⁽²⁾;
- Do not leave the hydrogen peroxide gel reach the retraction, preventing sharp pain;
- Certify that the teeth to be bleached do not contain carieses or infiltrations in the restorations.

Cares During the Bleaching

- In case that the gel has contact with some soft tissue during the bleaching, the place will be white and causing

pain. In this in case, pass in the place a cotton absorbed of bicarbonated water (supplied in the kit);

- If some tooth that will be being bleached has breaking, or incisal consuming, protect with the gum barrier photopolimerizer Lase Protect⁽²⁾ before the application of peroxide preventing sensitivity;
- During the bleaching, if only one tooth present sensitivity, verify if is not having gel leakage for the gum causing the burning sensation. If it's occurring, the gum must be cleaned and be protected with the gum barrier again;
- If the patient presents pain in all the teeth, the bleaching will have to be interrupted. Request, then, that the patient makes mouth washing with fluorine, followed by the application of a desensitizing gel or if necessary to carry through a lasertherapy session tooth per tooth with the infra-red laser. The equipment Whitening Lase II⁽³⁾, Photon Lase III⁽⁵⁾ and Flash Lase III⁽⁶⁾ can be used for this function. The return to another bleaching section will have to occur only one week later;
- In the end of the bleaching if the patient presents sensitivity, apply a desensitizing gel. A laser of low power also could be used in these cases.

Cares After the Bleaching (Right after the bleaching)

- Do not smoke in the day of the bleaching;
- Do not ingest acid foods as: fruits or soft drinks;
- During 48 hours, do not ingest anything that contains colorings. Ex.: coffee, chocolate, cola soft drinks, red gravy, beetroot, red wine, etc;
- Do not use colorful dental cream or that contains sodium bicarbonate;
- Make mouths washing with fluorine.

CONTRA INDICATIONS

- Pregnant or nursing;
- Patient carriers or with predisposition to cancer;
- Patient with periodontal illnesses;
- Patient with incomplete root formation (Age < 17);
- Patient with allergic reaction to the product components;
- Patient with sensible teeth.

LIMITATIONS

Restorations, prosthesis and exposed cervix, are not open to bleaching.

PRODUCT EFFECTIVENESS AND SECURITY DESCRIPTION

Clinical Effectiveness and Security

The dental bleaching is an efficient, traditional, secure technique, of simple execution and that provides satisfaction to the patient. The final result is a natural and pleasant appearance. The same must be considered the first option of treatment in the cases of color alterations, for being a technique of low cost, with few collateral effects and does not imply in consuming of the dental structure. However, in more severe cases of color alterations, it can be associated to restoring treatments, as a previous step to these.

Toxicity and Inflammability

In normal conditions of use the product is perfectly safe and it is not inflammable. It must not be exposed to the extreme heat or flames.

Emissions

It has production of ions of fluorine, ions of potassium and oxygen with purpose foreseen in the treatment. It does not have emission of thermal or luminous energy, ionizing ra-

diation or emission of volatile gases.

COLLATERAL EFFECT

Cases of hypersensitivity after the treatment are told in literature. In these cases must be made use of a desensitizing gel or a lasertherapy must be carried through tooth per tooth with the infra-red laser. The equipment Whitening Lase II⁽³⁾, Photon Lase III⁽⁵⁾ and Flash Lase III⁽⁶⁾ can be used for this function.

GENERAL INSTRUCTIONS OF SECURITY

The manipulation of the product must be made exclusively by dentists. All the necessary precautions must be taken, in order to prevent the accidental contact of the product with skin, mucous or clothes.

The high oxidant power of the Lase Peroxide Flex⁽¹⁾ makes that any contact of the product with live tissues provokes a temporary irritation, with the sprouting of white spots, that normally disappears in a maximum period of 2 hours.

In the case of accidental contact with the skin or intra-verbal mucous, the region must be cleaned with water in abundance. In these cases, the action of peroxide can be neutralized covering the affected area with some drops of the bicarbonated water solution (supplied in the Bleaching Kit).

In case of accidental contact with the eyes, wash them immediately with water in abundance and look for a doctor.

DISCARDING INSTRUCTIONS

The discarding of the Lase Peroxide Flex⁽¹⁾, after its validity, is responsibility of the user and must attend the local legislation.

EQUIPMENT TO BE USED

The procedures described in this manual had been step by step elaborated for the use of the Dental Bleacher Lase Peroxide Flex⁽¹⁾ with systems of light generation that conjugate a matrix of LED's, with wave length of 460 nm to 480 nm, and with an emitting diode infra-red, with wave length between 808 nm - 830 nm. Other equipment emitting light, with wave length between 460 nm and 480 nm, also can be used, having always in mind that for being with a photo-activation process, a minimum level of Energy is necessary so that the process is unchained. The minimum power density recommended is of 200 mW/cm² and maximum of 350 mW/cm².

The procedure of photo-polimerization of the Lase Protect⁽²⁾ can be effected with common photo-polimerizers with light emission between 460 nm - 480 nm or using the systems Whitening Lase II⁽³⁾ or Whitening Lase Light Plus⁽⁴⁾.

The reduction of sensitivity procedure, by lasertherapy, must be carried through with the following equipment, which possess the infra-red laser: Whitening Lase II⁽³⁾, Photon Lase III⁽⁵⁾ and Flash Lase III⁽⁶⁾. The minimum fluency must be of 40 J/cm² and the maximum of 60 J/cm² or a minimum energy of 1,16 J and maximum of 1,74 J.

The professional will have to be intent, about the use of other light sources, in relation to the adequacy of the suggested products, considering the proper characteristics of these power plants to be used.

VALIDITY

The product has validity of 2 years from the manufacture date, which is printed in the primary packing (label) and secondary packing (external box).

HAVE ALWAYS IN MIND THE PREVENTION
OF ACCIDENTS AND THE SECURITY WHEN
MANIPULATING THE PRODUCT.

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After previous Prophylaxis (3 days before), register the initial coloration of the dental elements.

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In the option for the relative isolation, locate initially a labial retractor.

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
Dry the teeth and the region of the gum tissue for application of the Lase Protect⁽²⁾. After, remove the syringe cover and incase the application tip firmly.

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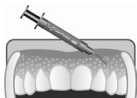Apply the Lase Protect⁽²⁾ in small layers covering about 0,5 cm of the surrounding gum tissue, stamping all the region to be worked, also root expositions.

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Photo-polimerize the Lase Protect⁽²⁾ per at least 30 seconds before the bleaching agent application.

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
Notes:

- 1)** *Agitate the bottle of Phase 2 before the use, to homogenise the product;*
- 2)** *17 drops from Phase 1 + 17 drops from Phase 2 + 17 drops from Phase 3, in case of 15%, are the sufficient for both the arches;*
- 3)** *39 drops from Phase 1 with 13 drops from Phase 2, in case of 35%, are the sufficient for both the arches.*

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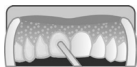1mm a 2mm de espessura sobre todo o esmalte vestibular.

Apply the gel on the initial surface of teeth. A gel layer must be applied with 1mm to 2mm of thickness on the initial enamel.

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Radiate with the Whitening Lase II⁽³⁾ or Whitening Lase Light Plus⁽⁴⁾ for:


15% - 2 minutes the jaw arc and after that, per 2 minutes the mandible arc (repeat 4 times this procedure) giving 16 minutes of gel exposition in each arc.

35% - 1 minute the jaw arc and after that, for 1 minute the mandible arc and to leave it resting 1 minute (repeat 3 times this procedure) giving 8 minutes of gel exposition in each arc.

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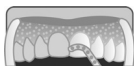Each gel application must have the duration of 16 minutes in case of 15% and 8 minutes in case of 35%. After this time remove the gel with the aspirating tip and clean the surface with a gauze.

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Repeat the procedures from 6 to 9 for until more two times depending on the gotten result.

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
After the ending of the third bleaching agent application, inhale with teething ring and wash with water in abundance.

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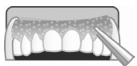Remove the gum barrier detaching it for front.

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
13 Polir os elementos dentais utilizando um disco impregnado umedecido.

Polish the dental elements using a humidified impregnated disc.

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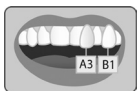Make a final color taking. Apply a desensitizing gel and if necessary to carry through lasertherapy.

NOTE: In case of more than a session be necessary, will have to be given a minimum interval of 7 days between the sessions.

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
Important:

- *The emitting light equipment (blue LED (visible) and infra red (invisible)) used must present an energy of, at least, 200 mW/cm² and maximum of 350 mW/cm² of irradiance. LED or LED/Laser System with low luminous energy are not indicated for this dental bleaching technique;*
- *The use of bicarbonate spurt immediately before the bleaching procedure can inhibit the action of the bleaching agent, since the composed related is exactly the neutralizing agent of the hydrogen peroxide solution that composes the kit;*
- *Must not use any desensitizing product right before the application of the Lase Peroxide Flex⁽¹⁾.*

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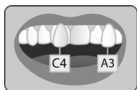After previous prophylaxis (3 days before), register the initial coloration of the tooth to be bleached and neighboring teeth.

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Promote the coronary opening and removal of the restoring material and residues. Remove gutta-percha from the pulp chamber up to 3 mm below the amelo-cementoblasts junction, entering in the root conduit.

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Confection the peripheral sealing with glass or zinc phosphate cement ionomer with minimum thickness of 2 mm.

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Apply the gum barrier for proximal vestibular and palatal/lingual with thickness of 0,5 cm protecting all gum tissue area. Photopolimerize the gum barrier for 30 seconds.

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
Effect the hygienic cleaning of the pulp chamber before preparing the gel. In this in case, the indicated concentration is of 35%. See gel preparation in Table 02.

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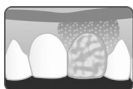Apply a layer of the bleaching gel Lase Peroxide Flex⁽¹⁾ with 2 mm of thickness for vestibular and internally the pulp chamber. The gel must be applied 3 times in the same session and in each application, to remain from 10 to 15 minutes in position.

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Radiate with the Whitening Lase II⁽³⁾ or Whitening Lase Light Plus⁽⁴⁾ per 3 minutes in the initial face, followed of application in the lingual/palatine for the same time. Repeat the irradiations for more two times.

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Inhale the gel locating for the initial face and inside the pulp chamber. Clean the surface with a gauze.

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Repeat the procedures from 6 to 8 for until more two times depending on the presented result.

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
After the third application, inhale the gel and after that, wash with water in abundance.

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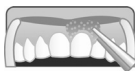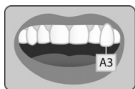
After the ending of the bleaching session, carry through the provisory restoration, applying delay dressing, if necessary.

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If the desired color after the first session is not reached, the procedure must be repeated after 7 days.

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The results presented here are only illustrative.

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