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Shiatsu Massage Device for the back, buttocks and thighs.

- Relieves muscle tension and relaxes your body
- Massage device with 4 heated massage heads moving up and down to stimulate the whole back + vibrations for buttocks
- Imitates the hands of a professional masseur
- Remote control for a comfortable use

The pack contains: (photo 1 page 1)

- 1 Back heated massager chair with vibration for buttocks (A)
- 1 controller (B)
- 1 AC adapter (C)
- 1 plug to connect with cigarette lighter (D)

When using this unit for the first time, read carefully the instructions and the following guide-lines in order to become familiar with the controller.

How to install the device

1. Plug the AC adapter to a power socket; connect the 12V cable to the controller.
2. Fasten the Back massager to your chair or seat with scratch straps
3. Sit down comfortably

Functions

This device performs kneading heating massage of the back and vibration for thighs to relieve the muscle tension.

Double built-in vibrating motors on the thighs level relaxes the body. Unique up-down device. Based on the properties of physical nature and characteristics of people. This product was developed to be able to imitate human hands to knead, press, roll and vibrate. Multi massage heads to massage in big area. Combination of kneading, heating and vibrations.

4 kneading heating heads have walking function: can roll and slide down to shiatsu and knead your body.

- Three massage variations for the full back, lower back or upper back.
- Heating option for the back
- Three vibration intensity levels for buttocks weak, mid, strong.

BE CAREFUL

Suitable to be used in an office, at home, in a car. (Engine shut down).

In case of using the Shiatsu massager in a mechanized vehicle, please, follow those cautions:

- Always use the device when the vehicle is park and stopped, turn off the engine.
- After use, pack the device before turn the engine on.

Safety Precautions

This unit is designed to be used by healthy adults. However, the following precautions should be taken when using the unit

- Do not use the device in high temperature or humid environment
- Avoid the breakage of the power cord
- Do not put the cable or plug on any heated surface
- Do not dismantle the device yourself. In case of any problem, apply to the Customer Service of Sport-Elec Institut
- Do not use the device under a pillow, a blanket, in case of any fire accident.
- Do not hit the unit with a sharp article.
- Do not use the device out of doors.
- Keep this unit out of reach of children.
- Do not connect the unit with any other power supply different from that one mentioned at the AC adapter.
- Use this appliance only with supplied adapter.
- Before undertaking any exercise program, including this unit, consult your family physician to ensure that you are capable of conducting a regular routine without putting undue stress on your body.
- If you have never exercised before, if you are pregnant, recently had surgery or suffer from any significant illness, including epilepsy, hemophilia, diabetes, cancer, hypertension, heart disorders, neuromuscular disorders and more, you should not conduct any exercise routine, including this unit, unless under the direction of your physician.
- Do not use the device while drowsy or sleeping.
- Do not use this device with any other electric equipment at the same time (e.g., another electrical muscle stimulation device).
- Do not use this unit and any other machine at the same time.
- The device is still under tension when the adapter is plugged into a power socket. Please, switch to "OFF" and unplug the device to stop use it.

This unit is neither a medical nor therapeutic device.

Warning

This product is designed to be used only while seated upright.

- Never use this product in bed
- Never use this unit while sleeping
- Never cover yourself, or the product, with a blanket or with any other material while using this product
- Do not use this product for more than 15 minutes at a time
- This product should never be used by individuals suffering from:
 - sensory deficiencies in the lower back
 - any physical ailment that would limit the user's capacity to operate control.

Never leave this product in use, or plugged in while unattended.

Any violation of the above safeguards may cause the unit to overheat, resulting in serious damage or fire

Start the device : (Picture 2 p.1)

Choose your program of massage: -

- 1- Start the device by pressing on the button "On / Off".
- 2- "Heating": Hit option for the back massage
- 3- «Top back»: to mass the top of the back and shoulders, from bottom to the top and from top to bottom.
- 4- «Full back»: to mass the whole back, from bottom to the top and from top to bottom, from renal area to shoulders.
- 5- «Bottom back»: to mass the bottom of the back and the renal area, from bottom to the top and from top to bottom.
- 6- «LED» Red leds to see massage heads movements along the body
- 7- In this program, you can also choose the intensity of the massage with 3 possibilities: "Low" = weak; "Mid" = middle; and "Hig" = strong.
- 8- «Stagger»: oscillating massage for buttocks.
- 9- «Continue»: vibrating massage continues for buttocks.
- 10- «Tap», massage by tapping for buttocks..

Every program lasts 15 minutes after which the device switches off automatically.

At the end of session :

At the end of your massage session, put the remote control in the case situated on the left of the shiatsu seat, then unplug the AC adapter from the power socket and unplug the 12V cable from the controller.

Finally store your Massage Device in the box and keep it in a dry place.


Specification:

Product name : Back heated massager Shiatsu

Ref.: JKW-832F

Adaptator: 100V - 240V ~ 50Hz
12V --- 4A

Dimensions : 123 x 57 x 5 cm
Net Weight : 5.8 Kg

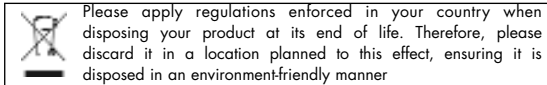WARRANTY

This appliance is guaranteed for 24 months from the date of purchase. During the period of the guarantee, the distributor will take responsibility for defects occur, the appliance will be repaired or exchanged as the distributor deems necessary. Damage arising from improper use or normal wear and defects, which do not materially affect the performance of the appliance, is not covered under this guarantee. The guarantee is not valid if unauthorized repairs have been attempted or if replacement parts used are not approved by the supplier. The guarantee is valid only if it is stamped and signed by the retailer at the time of the purchase. For all claims made during the period of the guarantee, send the appliance and the warranty card to your retailer. All other claims are excluded from the guarantee unless specifically prohibited by law.

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