

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following: **READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE.**

DANGER— To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2) Do not reach for an appliance that has fallen into water. Unplug immediately.
- 3) Do not use while bathing or in a shower.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place or drop into water or other liquid.

WARNING— To reduce the risk of burns, electrocution, fire, or injury to persons:

1. This product is for household use only.
2. Close supervision is necessary when any appliance is used by or near children or by invalids.
3. Unplug the appliance from outlet when not in use and before cleaning. Keep cord away from heated surfaces.
4. Do not operate this appliance with a damaged cord or plug or if the appliance malfunctions or has been dropped, dropped in water or damaged in any manner. Return appliance to nearest service facility for examination, electrical or mechanical adjustment or repair. No repairs should be attempted by the consumer.
5. Use appliance only for intended use as described in this instruction book. The use of attachments not recommended or sold by OBUSFORME by CONAIR may cause hazards.
6. Do not use outdoors or operate where aerosol (spray) products are being used, or where oxygen is being administered.
7. Turn switch to "OFF" position and unplug cord before cleaning the unit.
8. Never use while sleeping.
9. Do not overfill unit past "MAX" line inside unit. Use this appliance on level surface.
10. Do not heat paraffin wax in its plastic packaging and container.
11. Do not heat paraffin wax on stove or microwave. The paraffin unit is regulated at an even temperature to prevent scalding.
12. Do not operate unit without paraffin wax.
13. Do not pour melted paraffin wax from the tank.
14. Do not add water, lotions or any other additives to the wax.



15. Do not use if you have any skin irritations, rashes, cuts, burns, or warts.
16. Use only **OBUSFORME BY CONAIR Paraffin Spa** wax in this unit.
17. Do not soak or immerse unit in water. See **MAINTENANCE** for cleaning instructions.
18. Do not let the cord hang over the edge of the table or counter where it can be pulled on by children or tripped over.
19. Do not use extension cord with this appliance.
20. This appliance has no serviceable parts.

SAVE THESE INSTRUCTION OPERATING INSTRUCTIONS


1. Insert two “AA” Batteries into battery box on bottom of unit.
2. Place unit on a stable, level surface at a height that is convenient for treating the desired area.
3. Remove lid, wax container and grate from inside the unit. Discard all packaging on wax.
4. Replace only the wax into the unit. Place grate on top of the wax. As the wax melts, the grate will fall to the bottom of the unit. The grate is used to prevent your hand or foot from touching the bottom of the unit.
5. Place lid on unit. Lid should fit securely on unit.
6. Plug unit into electrical outlet (120 VAC 60 Hz). This appliance has a polarized plug (One blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature. (The indicator light will illuminate.)
7. **SETTING THE CLOCK**
 - A) Press and hold “SET CLOCK” button.
 - B) Press “Hour” and “Min” buttons to correct times; note AM/PM settings.
 - C) Release buttons; time is now displayed.
 - D) Repeat above steps to change the time.
- 8). **USING THE HIGH/LOW TEMPERATURE ADJUSTMENT**

Turn dial to “High” for a quicker heat-up and higher temperature. Turn dial to “Low” for a lower operating temperature. Use “low” for your first use, and adjust temperature control for subsequent uses.
- 9). 6 lbs. of wax will take approximately 5 hours to melt on the “HIGH” setting, and longer on lower settings. (See Treatment Instruction before dipping in wax). Do not place additional wax

**WARNING:
IF YOU THINK THE POWER IS OFF
WHEN THE SWITCH IS OFF,
YOU'RE WRONG.**

KEEP AWAY FROM WATER

Everybody knows that electricity and water are a dangerous combination. But did you know that an electrical appliance is still electrically alive even if the switch is off? If the plug is in, the power is on. So when you are not using your appliances, keep them unplugged.

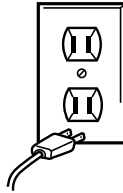

ALWAYS UNPLUG SMALL APPLIANCES

Public Service AD from Association of Home Appliance Manufacturers and Underwriters Laboratories, Inc.

into the unit until wax level drops at least 1”.

10) MANUAL OPERATION


Press “ON” button located on the left side of the control panel. Note that a light behind “ON” button is illuminated to indicate that you are operating in the manual mode. Wax will melt until shut off or until unit automatically shuts off after approximately 12 hours.

11) SETTING & USING THE AUTOMATIC TIMER

- A) Press & Hold “SET AUTO” Button.
- B) Press “Hour” & “Min” button to SET THE TIME YOU WANT TO USE THE WAXER.
- C) The SET TIME must be at least 5 hours from the current time. If you have programmed in too early a SET TIME, the display will show an error message (“ERR”).
- D) Press AUTO button, it will be illuminated to indicate that you are in the AUTO mode.
- E) Wax will be melted (or almost melted) at the time you have set. The base of the waxer will illuminate to indicate that the wax is ready to be checked for dipping use (See TREATMENT INSTRUCTIONS before dipping in wax).
- B) Unit will automatically shut off 12 hours after initial start.

TREATMENT INSTRUCTIONS

1. Place the plastic liners and universal mitts next to the unit for your convenience.
2. Wash hands or feet thoroughly. Check to make sure that you have no cuts or visible patches of sensitive skin.

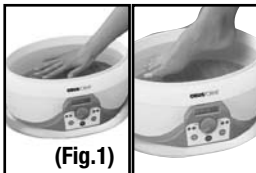

3. Remove any jewelry and roll up long sleeved shirts.
4. If you are sensitive to temperature, test the wax before beginning. Simply dab a bit onto the inside of your wrist. It should be very warm to the touch, but not uncomfortable. **DO NOT DIP YOUR HAND IF WAX IS TOO HOT.**
5. With fingers (or toes) slightly apart, immerse hand slowly in paraffin unit. **(Fig.1)**
6. Immediately pull your hand out, being sure to maintain the open-finger position. Wait 3 to 5 seconds until the paraffin loses that “wet” look and then dip again. Try to immerse your hand

LIMITED ONE YEAR WARRANTY

OBUSFORME by CONAIR will (at our option) repair or replace your appliance free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service centre, together with your sales slip and \$7.50 for postage and handling. In the absence of a purchase receipt, the warranty period will be 12 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12 MONTH DURATION OF THIS WARRANTY. IN NO EVENT SHALL OBUSFORME by CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

This warranty gives you specific legal rights, and you may also have other rights which vary from province to province. For service, please bring or send defective unit with sales slip to the service centre listed below.

Distributed by:

Conair Consumer Products Inc.

Obus Forme® and The Body Care Experts™ are registered trademarks owned by Integral Orthopedics, Inc. used under license.

Distributed by:

Conair Consumer Products Inc.

SERVICE CENTRE:

Conair Consumer Products Inc.

100 Conair Parkway,

Woodbridge, ON L4H 0L2

OBPB21C

For more information on any Conair Products call:

1-800-472-7606

Visit our Web site: www.conaircanada.ca

Consumer Call Centre E-mail: Consumer_Canada@Conair.com



IB-7906

07CN10377

to the same depth each time, so the lip of the wax “glove” thickens uniformly.

7. Repeat this process until 5 layers have been formed. Additional layers can be added to retain heat longer.
8. Slip your hand into the plastic liner and then insert it into the terry cloth mitt.
9. Sit down and relax for 10-15 minutes, while the warm paraffin wax continues to release heat to soothe and soften the hand.
10. Remove the terry cloth mitt and plastic liner once you no longer feel any heat. Wiggle your fingers slightly to loosen the wax. Then peel off the paraffin starting from the wrist and pulling it downward like a glove, into the plastic liner.
11. Discard the plastic liner and wax. See replacement section on how to order replacement wax blocks and liners.
12. Return lid to unit, switch the unit to “Off” and unplug. Allow unit and wax to cool before stowing unit away.
13. For the feet, do the same as the above, but simply place the unit on the floor.

MAINTENANCE

After 25 paraffin treatments or when visible evidence of debris occurs, the unit should be cleaned and the paraffin wax should be replaced.

1. Unplug unit and remove lid, grate and wax.
2. If wax is still hot, wait until it solidifies. Push on one end of the wax block to lift up.
DO NOT USE SHARP OBJECT. Remove wax and discard.
3. Clean inside of the unit with a light household cleaner and wipe dry.
4. Once unit is thoroughly clean and dry, place lid on unit and stow away. Do not add new paraffin wax until you are ready to use.

REPLACEMENT

OBUSFORME BY CONAIR Paraffin Spa replacement kit (model #PWR6) is available through the Conair service centre (1-800-472-7606) and retail stores carrying the OBPB21C unit. The kit includes: 2 lbs of wax and 20 disposable liners.

OBUSFORME®

THE BODYCARE EXPERTS™



Programmable Paraffin Spa

with Automatic timer so the wax is melted
when you want to use it

Model OBPB21C Instruction Booklet

For your safety and continued enjoyment of this products, always read the instruction book carefully before using.