

Digital thermometer Model Flex Temp Smart

- Digitales Fieberthermometer
Thermomètre digital
Termómetro digitale
Termómetro digital
Digitalt termometer
Digitaalinen kuumemittari
Digitalt termometer
Digitalt feberthermometeret
Digitale koortsthermometer
Ψηφιακό θερμοόμετρο
Dijital klinik termometresi
Цифровой термометр

الترمومتر الرقمي

- DE
FR
IT
ES
PT
DK
FI
NO
SE
NL
EL
TR
RU
AR

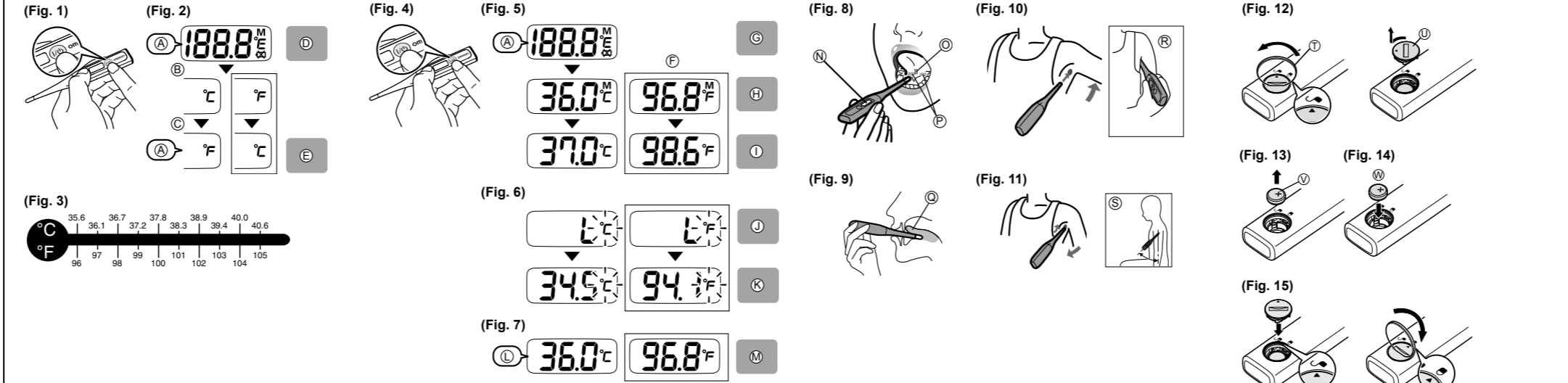


Table with 24 columns representing different languages: English, Deutsch, Français, Italiano, Español, Português, Dansk, Suomi, Norsk, Svenska, Nederlands, Ελληνικά, Türkçe, Русский, العربية. Each column contains instructions for using the thermometer in that language.

Digital Thermometer English

A. Sensor Tip
B. Display Window
C. I/O Button
D. Storage Case
Introduction
Thank you for purchasing the OMRON Digital Thermometer Flex Temp Smart.

Important Safety Notes
To assure the correct use of the product basic safety measures should always be followed including the precautions listed below.

Warning
A high or prolonged fever requires medical attention, especially for small children. Please contact your doctor.

Caution
Do not bite down on the thermometer. Doing so may lead to breakage and/or injury.

General safety precautions
Do not use the thermometer to measure anything other than human body temperature.

Switching between Fahrenheit/Centigrade (Celsius)
Press the I/O button for 7 seconds until a beep is emitted.

6. Display Indications and Troubleshooting

Table with 3 columns: Indicator, Cause, Remedy. Includes indicators for low/high temperature, air temperature display, and storage case.

7. Storage
Do not store the thermometer in the following types of places. Doing so may damage the thermometer.


8. Cleaning and Disinfecting
OMRON recommends to wipe the device with damp cloth. For disinfecting, 70% Ethanol or Isopropyl alcohol can be used.

9. Specifications
Product Type and Number: Digital Thermometer MC-343F-E/MC-343F-EB
Product Name: Flex Temp Smart

Applied part:
Type B
Notes:
The specifications may be changed without prior notice.

Read the Instructions in this manual carefully.
This device fulfills the provisions of the EC directive 93/42/EEC (Medical Device Directive) and the European Standard EN12470:2000.

Battery Replacement Indicator
Replace the battery when the battery replacement indicator appears when the thermometer is switched on.

Manufacturer 	OMRON HEALTHCARE Co., Ltd. 53, Kunotsubo, Terado-cho, Muko, Kyoto, 617-0002 JAPAN		
EU-representative <table border="1" data-bbox="39 200 177 260"> <tr> <td data-bbox="39 200 104 260">EC</td> <td data-bbox="104 200 177 260">REP</td> </tr> </table>	EC	REP	OMRON HEALTHCARE EUROPE B.V. Scorpius 33, 2132 LR Hoofddorp, THE NETHERLANDS www.omron-healthcare.com
EC	REP		
Production facility Produktionsstätte Site de production Stabilimento di produzione Oficina de producción Local de produção Produktionsanlæg Tuotantolaitos Produksjonsanlegg Produktionsenhet Productiefaciliteit Εργοστάσιο παραγωγής Üretim Tesisi Производственное подразделение مصنع الإنتاج	OMRON (DALIAN) CO., LTD. Dalian, CHINA		
Subsidiary Niederlassung Succursale Consociata Empresa filial Filial Datterselskab Tytäryhtiö Filial Dotterbolag Dochteronderneming Θυγατρική Yan Kuruluş Дочерняя компания الشركات التابعة	OMRON HEALTHCARE UK LTD. Opal Drive, Fox Milne, Milton Keynes, MK15 0DG, U.K. OMRON MEDIZINTECHNIK HANDELSGESELLSCHAFT mbH John-Deere-Str. 81a, 68163 Mannheim, GERMANY www.omron-medizintechnik.de OMRON SANTÉ FRANCE SAS 14, rue de Lisbonne, 93561 Rosny-sous-Bois Cedex, FRANCE Uniquement pour le marché français: OMRON Service Après Vente N° Vert 0 800 91 43 14 consommateurs@omron-sante.fr www.omron-sante.fr/contact		

Made in China	Produseret i Kina
Hergestellt in China	Tillverkad i Kina
Fabriqué en Chine	Geproduceerd in China
Prodotto in Cina	Κατασκευάζεται στην Κίνα
Fabricado en China	Çin'de Üretilmiştir
Fabricado na China	Сделано в Китае
Fremstillet i Kina	صُنِعَ في الصين
Valmistettu Kiinassa	