Oral-B®

Professional Care

D 16.525

D 16.513

powered
by

BRAUN

English	4	UK	0 800 416 434
Dansk	8	IRL	1 800 509 448
Norsk	12	DK	70 15 00 13
Svenska	16	N	22-72-88-10
Suomi	20	S	020 - 21 33 21
		FIN	020 377 877
		HK	852-25249377 (Audio Supplies Company Ltd.)

Internet:

www.oralb.com

www.braun.com

www.service.braun.com

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





Printed in Germany

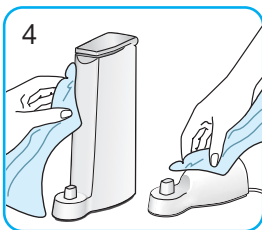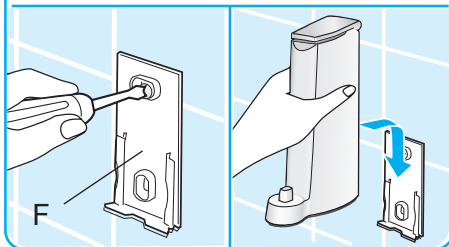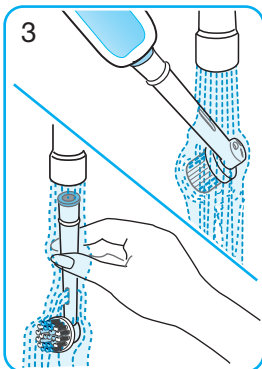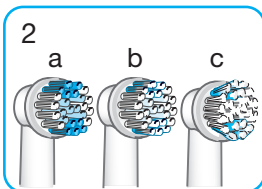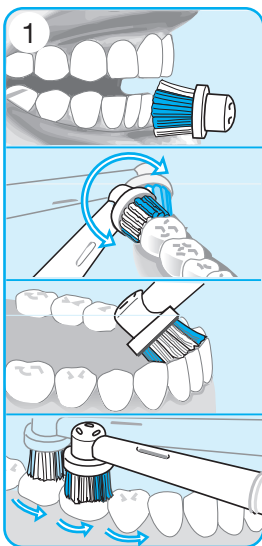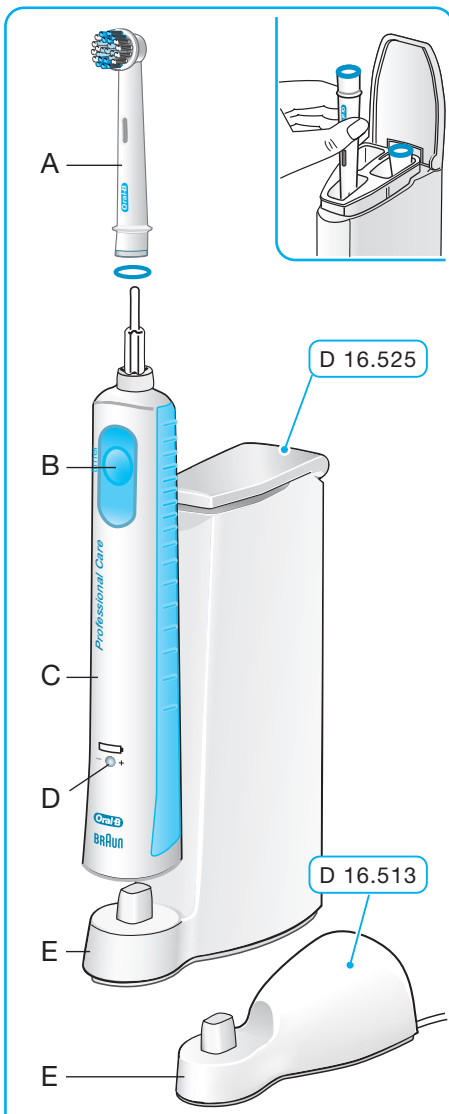English

Please read the use instructions carefully before first use.

Important

- Periodically check the cord for damage. A damaged or non-functioning unit shall no longer be used. Take the charging unit (E) to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age three.
- This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are given supervision by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.

Description

- A Brush head
- B On/off button
- C Handle
- D Charge indicator
- E Charging unit
- F Wall holder (D 16.525 only)

Specifications

Power supply: see bottom of charging unit
Handle voltage: 1.2 V

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and is designed for use in the bathroom.

- Plug the charging unit (E) into an electrical outlet and place the toothbrush handle (C) on the charging unit.
- The green charge indicator (D) light flashes while the handle is being charged. If you then use your brush and place it back on the charging unit, the green charging light will again flash to let you know it is charging.
- When the handle is fully charged, the green charge indicator light shines continuously. A full charge takes at least 12 hours and enables approximately one week of regular brushing (twice a day, 2 minutes).
- The red charge indicator (D) light starts blinking when the rechargeable battery is low. Now it is time to recharge your toothbrush.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle by regular use at least every 6 months.

Using the toothbrush

The toothbrush can be used with any kind of toothpaste.

- To avoid splashing, guide the brush head to your teeth before switching the appliance on.
- Slowly guide the brush head from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces (1).
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

Professional Timer

A short stuttering sound at 30-second intervals, which you can hear and feel, reminds you to brush the four quadrants of your mouth equally for a thorough consistent clean. A long stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

The timer memorises the elapsed brushing time, even when the handle is briefly switched off during brushing. Only when pausing longer than 30 sec. the timer resets to the initial starting point.

Pressure Sensor

For optimal brushing results your toothbrush has a pressure control feature to help you brush gently. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling it, you can also hear a different sound.

In the first days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months. The INDICATOR bristles help you monitor and improve your brushing performance. With proper brushing, twice a day for two minutes, the blue colour on the INDICATOR bristles will fade halfway approximately within 3 months, signaling it is time to change the brush head (2a-2c).

Replacement brush heads are available at your retailer or Oral-B Braun Service Centres.

Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on. Then switch off the handle and remove the brush head. Rinse both parts separately under running water (3) and wipe them dry. From time to time, also clean the charging unit with a damp cloth (4).

Subject to change without notice.

This product conforms to the European Directives EMC 2004/108/EC and Low Voltage 2006/95/EC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

The 30-Day Money Back Guarantee

Take the challenge! Experience the **Oral-B ProfessionalCare** difference! Satisfaction guaranteed or your money back.

Oral-B ProfessionalCare provides outstanding cleaning for clean teeth and healthy gums. In fact, we are so convinced you'll love your **Oral-B ProfessionalCare** that we give you this special offer.

Try **Oral-B ProfessionalCare** for 30 days starting from the day of purchase. If you are not 100% satisfied, receive a full refund. For a full refund, please send us your handle, refills, and charging unit along with your name, address, contact phone number and the original till receipt indicating purchase price and date.

United Kingdom / Ireland:

«Oral-B 30-Day Money Back Guarantee»
Dept PX415,
PO Box 10,
Diss, Norfolk,
IP98 1SP

This money back guarantee in no way affects your statutory rights.

Your **Oral-B ProfessionalCare** must be sent in its original packaging and be postmarked no later than 30 days after the date of purchase. Please allow 6–8 weeks for handling. Sorry: No refunds will be made for any **Oral-B ProfessionalCare** product damaged by accident, neglect, unreasonable use or lost in the post.

Guarantee

We grant 2 years limited guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

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