TANITA

Monitoring Your Health

Digital Lithium


Operating Instructions – HD-372

OPERATING INSTRUCTIONS

Digital Lithium Scale – HD-372

- ① Measuring Platform
- ② LCD Display
- ③ Battery Compartments on Bottom of Scale
- ④ Weight Mode Switch on Bottom of Scale

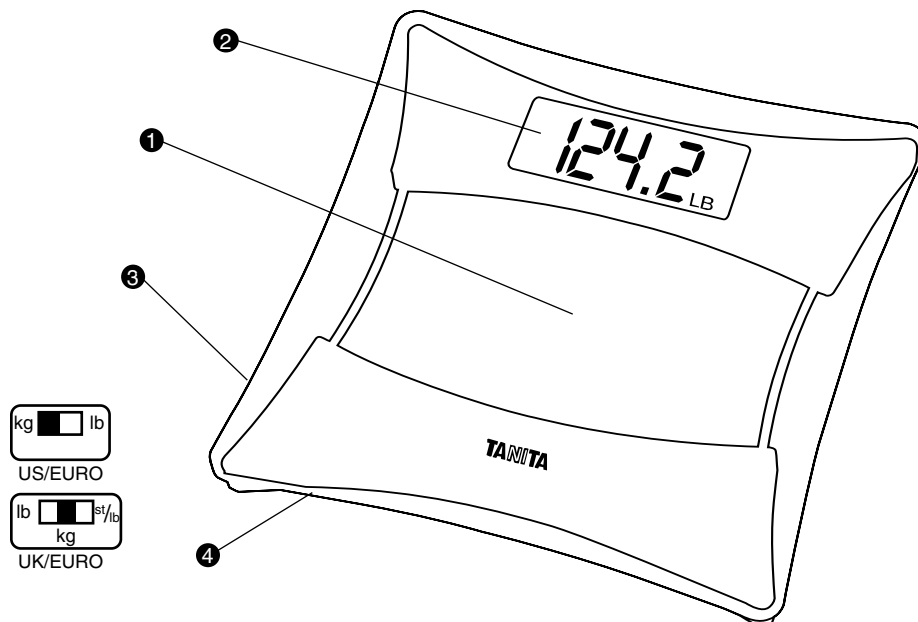Specifications

Maximum Capacity: 400 lb (180 kg)

Graduation: 0.2 lb (0.1 kg)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)
1 x 9 volt Alkaline Replaceable Battery (included)

UK Specifications

Maximum Capacity: 28 st 8 lb (180 kg) (400 lb)

Graduation: 1/4 lb (0.1 kg) (0.2 lb)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)
1 x 9 volt Alkaline Replaceable Battery (included)


Features

- Digital accuracy and reliability.
- Elegant, contemporary styling.
- Durable construction, non-skid surface.
- “L” message appears when the lithium battery needs replacing.
- Back light display using 9 volt battery

Before Use

Pull-off the protective tab from the battery.

Select your preferred weight mode using the Weight Mode Switch.

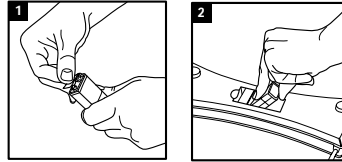Placing Scale

The scale should be used on a hard even surface for accurate weighing. Not recommended for use on carpeted surfaces. Placing scale on a very soft, deep pile rug may affect weighing accuracy.

Instructions for Use

Gently tap the center of the scale with your foot. When “00.0” shows on the display, step on the scale (Fig. 1).

Stand still while the LCD flashes “00.0” for 3 seconds. The LCD will display your weight (Fig. 2). If you step off from the scale, your weight will be displayed for 5 seconds then the scale will automatically shut off. If you remain on the scale, the weight will be displayed for 10 seconds then automatically shut off.


Note: If you step on the scale without tapping the platform, the display will show “Err” for 2 seconds then automatically shut off.

Once the scale displays “00.0”, if you do not step on the scale in 10 seconds, the scale will automatically shut off.

“Err” will appear on the display when the scale’s weight capacity is exceeded.

FIG. 1

FIG. 2

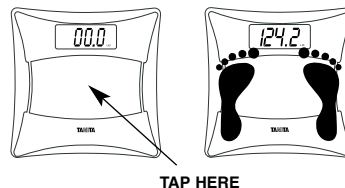
TAP HERE

Care and Maintenance

- Place the scale in an area free from heating equipment, high humidity, or extreme temperature changes.
- Never submerge the scale in water or spill chemical liquids on it. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impact or vibration to the unit.

Replacing the Batteries

Your scale is operated by one lithium battery for weighing. Lithium batteries can be typically used for many years before replacement is necessary. When the low battery sign, L, appears on the LCD, please replace with one new lithium battery. The blue backlight display is operated by one 9-volt alkaline battery. If the backlight does not work, please replace with new battery.

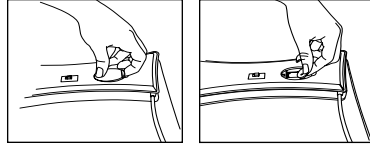Open the battery cover on the back of the measuring platform. Insert the new battery as indicated.

Note: Be sure that the polarity of the batteries is set properly. If the batteries are incorrectly positioned, battery fluid may leak and damage floors. If you do not intend to use this unit for a long period of time, it is advisable to remove the batteries before storage.

After Sales Service

Please contact the Tanita office nearest you regarding operation or trouble, request for repair, and/or placement of an order for parts for products.

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Tanita France S.A.
Tanita Health Equipment H.K. Ltd.
Tanita Europe GmbH
Tanita UK Ltd.