



mellerware®

P E R S O N A L

Retro Mechanical Bathroom Scale



complies with
IEC60355
SAFETY SPECIFICATIONS

m

CE
2 YEAR
GUARANTEE

Model:
20500

PLEASE READ ALL INSTRUCTIONS BEFORE USE

2 year guarantee






Dear customer,

Many thanks for purchasing a MELLERWARE product. Thanks to its technology, design and operation and the fact that it goes through strict quality controls - we feel sure you will be entirely satisfied.

- Place the scale on a hard and flat surface to ensure the greatest accuracy and reliability. Weighing on a carpet might affect accuracy of the scale.
- To ensure the most accurate readings, always weigh yourself on the same area of the scale platform and on the same floor surface. Uneven floors may affect readings.
- Weigh yourself without clothes and footwear, before meals and always at the same time of a day for a more accurate weigh. Early morning is a good time.
- Be sure the needle points to zero (0) before using it.

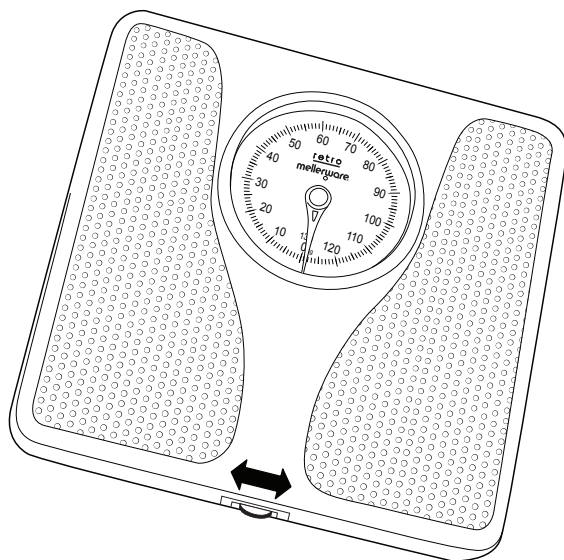

**Adjust the calibrating
dial left or right to set
the scale to zero (0).**

- Stand squarely and gently on the scale, without leaning forward or backward and distribute your weight evenly between your feet.
- Clean your scale with a damp cloth and do not allow water to get inside the scale.
- Do not use any chemical cleanser.
- Do not immerse your scale into any liquid, as this can damage the inner parts.
- Do not store the scale in an upright position.
- Treat your scale with care. It's a precision instrument, do not drop it or jump on it.
- This scale is only used in household for body weight measurement, not for commercial or other use.





creative housewares

27 Nourse Avenue Epping 2 Cape Town
P.O. Box 975 Eppindust 7475 South Africa
Tel: +27 21 534 8114 Fax: +27 21 534 8086
Service Centre Share Call: 086 111 5006
web: www.creativehousewares.co.za

