


**THE
SHARPER
IMAGE®**

Slim Scale™

EM150 Instructions

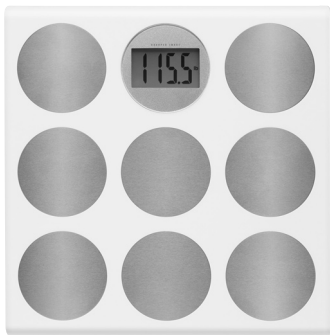

Read and save these instructions.

Dear Sharper Image Customer,

Thank you for your purchase of a Sharper Image product. We appreciate your business very much.

Our goal is that you will be completely satisfied with this product's quality and performance. We offer generous return privileges and excellent warranties on all our products.

Since 1977, The Sharper Image has shared the fun of discovering innovative products that make life better and more enjoyable. We've been in business a long time because we are an exciting place to shop and — just as importantly — we strive to deliver outstanding customer service that exceeds your expectations.

Enjoy your purchase!

Call 1-800-344-4444 for ordering and customer service, or visit
www.sharperimage.com, *our award-winning Website.*

Important Notes

- This scale is intended for use as a personal scale in a private household.
- This scale contains sensitive electronic components and must not be placed close to heat sources or in environments that are exposed to extreme humidity.
- Measurement results from various scales may differ slightly.
- The scale represents a sensitive electronic instrument. High-frequency devices may sometimes cause disturbance when used in the immediate vicinity of other products (for example, mobile phones, CB radio, remote controls and microwave ovens).
- If there are errors such as irregular or faulty readouts on the screen, remove the scale from the source of disturbance or deactivate the source of disturbance when you are using the scale.

Setup

1. Remove the insulation tab from the battery compartment.
2. Open the battery compartment and set the desired measurement unit by sliding the switch located in the compartment.

Setup

3. Place the scale on a solid, flat surface. To help prevent slipping, ensure that both the floor and the scale surface are dry.

Using the Scale

IMPORTANT: Never jump on the scale.

1. Gently tap your foot on the center of the scale. CAL appears on the LCD. CAL indicates that the scale is calibrated.
2. Stand with both feet on the scale once 0.0 appears on the LCD. Stand still, without shifting your weight.
3. During weight measurement, the displayed value increases until your weight has been determined.
4. Your weight is shown in 0.1kg/0.2lb. increments.

Warning Indicators

o-Ld: EEEE appears on the LCD if the weight exceeds the maximum weighing capacity of 180 kg/396 lbs./28 ST. Step off the scale immediately.

Warning Indicators

Lo: If Lo appears on the LCD, replace the batteries.
(See Replacing the Batteries.)

Err: The scale is unable to measure weight if Err appears on the display. Step off the scale, wait for a few seconds, then try again. Ensure you stand still.

Replacing the Batteries

1. Remove the battery-compartment cover.
2. Remove the batteries and replace with two new ones (CR2032 button cell). Observe proper +/– polarities.
3. Replace the battery-compartment cover, ensuring the cover snaps into place.

Storage

- Keep the scale in a clean, dry environment.
- Avoid using the batteries unnecessarily. For example, do not store any items on the scale.
- When storing the scale for a long period of time, remove the batteries.

Specifications

Maximum Load-Bearing Capacity	180 kg/396 lbs./28 ST
Minimum Measurable Weight	2 kg/4.4 lbs./4 ST
Selectable Measurement Units	kg/lbs./ST
Increments	0.1 kg/0.2 lbs.
Operating Temperature	32°F-104°F (0°C~ +40°C)
Power	2 CR2032 button cell (included)



This product has been tested to the relevant standards and complies with the requirements of the applicable EU directives.

Warranty

What Is Covered

This warranty covers any defects in materials or workmanship, with the exceptions stated below.

How Long Coverage Lasts

This warranty runs for 90 days from the date of original purchase.

What Is Not Covered

This warranty does not cover damage caused by misuse or use other than as intended and described in the product instruction manual, or loss or damage to batteries or removable parts.

What The Sharper Image Will Do

The Sharper Image warrants this product against defects in materials and workmanship. When repair is not practical, The Sharper Image will have the option to replace the defective product with one that is the same or has similar functionality.

How to Obtain Service

Return the product and receipt, along with a brief explanation of the problem, to:

T.S.I. - Returns Department
2901-A West 60th Street
Little Rock, AR 72209

If you have questions please feel free to contact our customer service representatives at 1-800-344-5555 or by email at care@web.sharperimage.com.

This warranty gives you specific legal rights, and you may have other legal rights that vary from state to state.

© 2006, 2007 Sharper Image Corporation

Sharper Image Corporation
350 The Embarcadero, Sixth Floor
San Francisco, CA 94105-1218

Shop online: **www.sharperimage.com**

Actual product may differ slightly from the product described herein.

V2 02/07