Professional Care™

7000 DLX series
série

English	4
Français	10
Español	15

Thank you for purchasing an Oral-B Professional Care™ which is designed and engineered by Braun.

If you have any questions, please call:

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
Type 4729

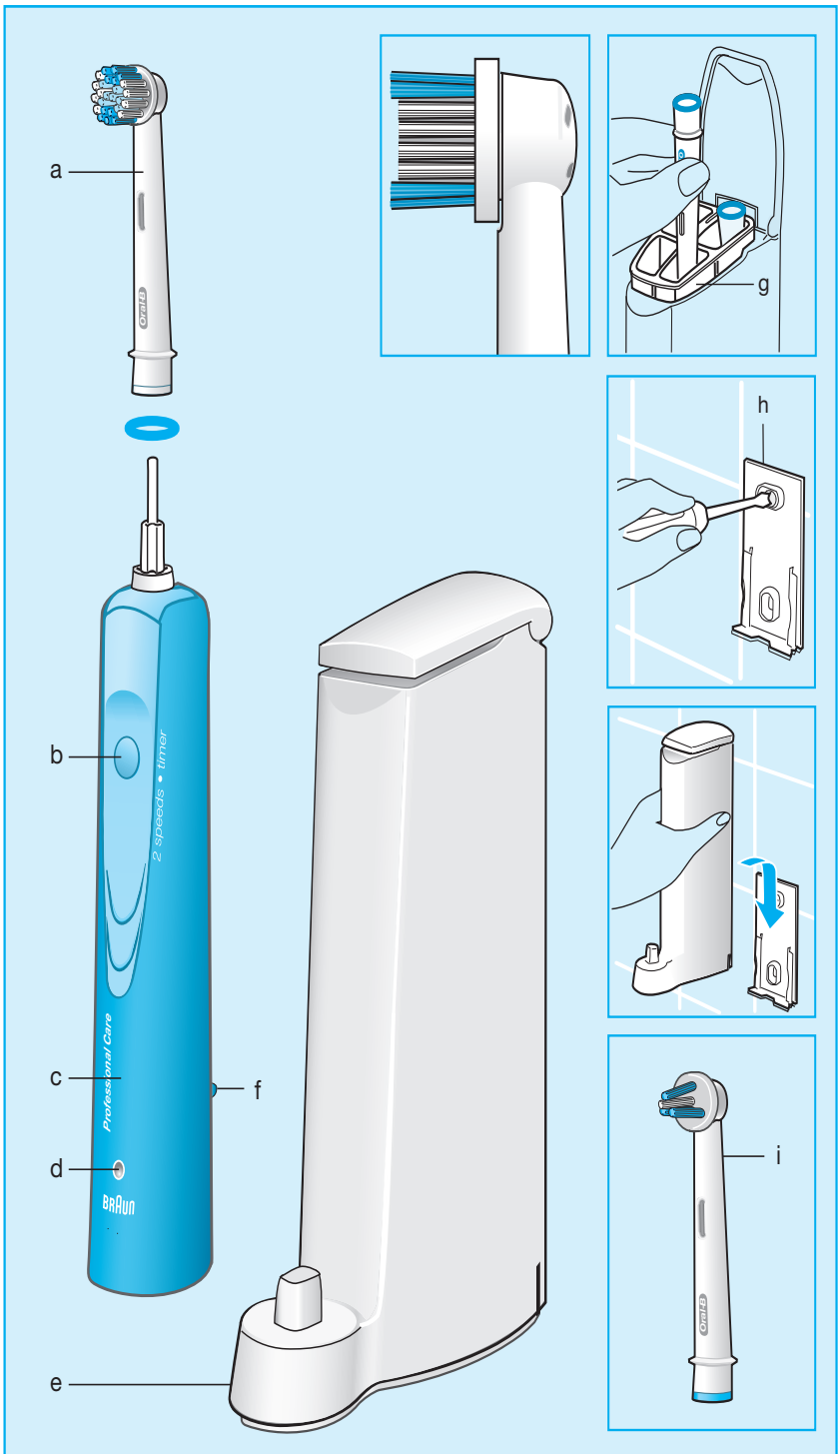IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

Read all instructions before using.

DANGER

To reduce the risk of electrocution:

1. Do not place or store charger where it can fall or be pulled into a tub or sink.
2. Do not place charger in or drop into water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if the charger has been dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never drop or insert any object into any opening.

SAVE THESE INSTRUCTIONS

The Oral-B Professional Care™ has been developed together with leading dental professionals to effectively remove plaque from all areas of your teeth. It adds a gentle pulsating action to its oscillation to create a deep-cleaning 3D brushing action. It has been designed for daily use for the whole family.

Please read the use instructions carefully before first use.

Important

Periodically check the cord for any damage. Should this be the case, take the charging unit (e) to a Braun Service Center. A damaged or non-functioning unit should no longer be used.

Description

- a Brushhead
- b Switch
- c Handle
- d Charge indicator
- e Charging unit
- f Roll-stop
- g Brushhead compartment
- h Wall holder
- i Power Tip®

Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V


Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit
- The charge indicator (d) flashes while the handle is being charged.
- When the handle is fully charged, the charge indicator lights permanently. A full charge takes at least 16 hours and provides an operation time of approximately 45 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

Battery maintenance



To maintain the maximum capacity of the rechargeable battery, you must discharge it at least once every 6 months. To do this, unplug the charging unit and allow the battery to run down with normal use.

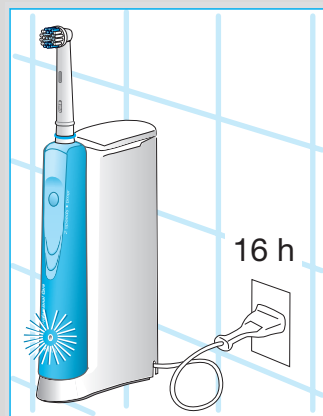Using the toothbrush

This product can be used with any toothpaste.

To avoid splashing, guide the brushhead to your teeth before switching the appliance on.

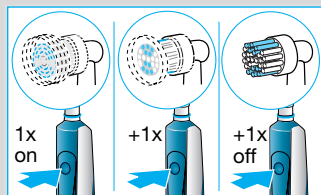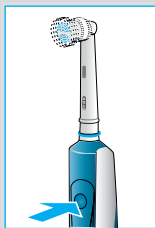Operating the switch






The toothbrush features two speed options, the lower speed being optimal for sensitive areas in the mouth.

- Press the switch to operate with normal speed,
- Press a second time to operate with reduced speed,
- Press a third time to switch the unit off.

Removing plaque

Guide the brushhead slowly from tooth to tooth. Hold the brushhead in place for a few seconds before moving on to the next tooth.

Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces. Do not press too hard or scrub, simply let the brush do all the work.

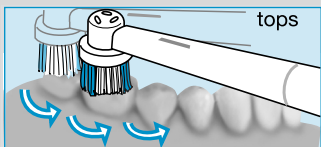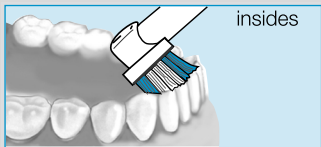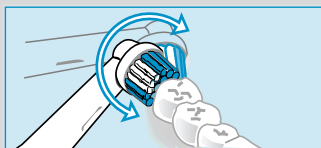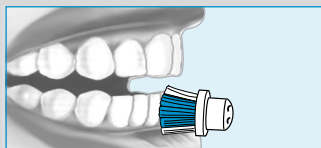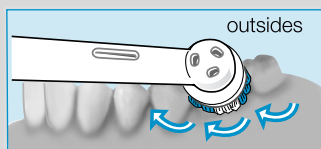Professional timer

With a signal you can hear and feel, the timer invites you to brush quadrant by quadrant. After every 30 seconds, it signals by a very short stutter that it is time to move to the next quadrant. After two minutes – which is the recommended brushing time – the timer signals by a longer stutter.

Pressure Sensor

For safety and optimum brushing results, the toothbrush has a pressure control system: If too much pressure is applied, the oscillating movement of the brushhead will continue, but the pulsation will stop. In addition to feeling it, you can also hear the difference.

INDICATOR® bristles

The INDICATOR® bristles help you monitor and improve your brushing performance. With proper brushing, twice a day for two minutes, the blue color on the INDICATOR® bristles will disappear half-way down within 3 months, signaling it is time to change the brushhead.

Power Tip®

The Power Tip® brushhead is designed for cleaning between teeth as well as around bridges, crowns and implants. To use, position brushhead at 45 degree angle to gumline. Slowly guide the brushhead along the gumline, pausing between teeth for optimal results.

Replacement brushheads

Available at retail stores or Braun Service Centers.

During the first few days of using the toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or hygienist.

Cleaning

- Rinse the brushhead thoroughly under running water for several seconds with the handle switched on.
- Switch off the appliance and remove the brushhead from the handle. Clean both parts separately under running water, then wipe them dry.
- From time to time, also clean the charging unit and the brushhead compartment with a damp cloth.
- For cleaning purposes, the brushhead compartment is removable.

The environment





For environmental reasons, your toothbrush uses an «environmentally sound» nickel-hydrate rechargeable battery. The battery does not contain any poisonous heavy metals that can damage the environment.

In the interest of recycling, please do not dispose of this appliance in the household waste at the end of its useful life. Instead dispose of it at an authorized Braun Service Center or at an appropriate collection site.

However, if you prefer to dispose of the battery yourself, please refer to the section «Battery removal» on page 20.

Subject to change without notice.

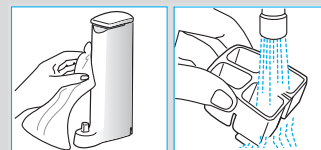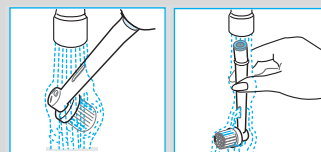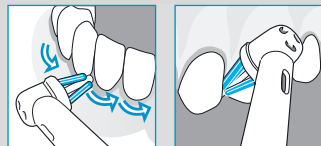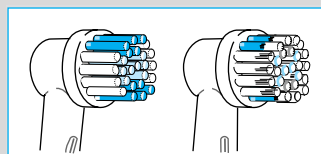FOR USA ONLY

Limited 2-Year Warranty

Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-BRAUN 11 (272-8611) to be referred to the Authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.
Pack the product well.
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts and accessories:

Contact an authorized Braun Service Center.
Please call 1-800-BRAUN 11 (272-8611) to be referred to the authorized Braun Service Center closest to you.

To order original Braun replacement parts and accessories (including replacement brushheads) by phone:

Please call 1-800-892-8534.

To speak to a Braun Consumer Service Representative:
Please call 1-800-984-9668.

For Canada only

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service centre, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages.

The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-387-6657 to be referred to the authorized Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Braun Service Center of your choice or,
- B. Ship the product to the authorized Braun Service Center of your choice.
Pack the product well.
Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Braun replacement parts & accessories:

Contact an authorized Braun Service Center.
Please call 1-800-387-6657 to be referred to the Authorized Braun Service Center closest to you.

To speak to a Braun Service Representative:
Please call 1-800-387-6657.

24 Month Performance Guarantee

For USA only

Use the Oral-B Professional Care™ daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B Professional Care™ significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B Professional Care™ within 24 months, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

24 Month Performance Guarantee

Braun Division,
The Gillette Company,
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B Professional Care™ must be shipped in its original packaging and be postmarked no later than 24 months after the date of purchase.

Please allow 6–8 weeks for handling.

Name of Regular Toothbrush

Your Name

Street

City

State

Zip Code

24 Month Performance Guarantee

For Canada only

Use the Oral-B Professional Care™ daily until you are fully familiar with this mechanical plaque removal instrument. If you are not convinced that the Oral-B Professional Care™ significantly reduces plaque, send it back for a full refund.

If not satisfied please send the Oral-B Professional Care™ within 24 months, post-paid and insured along with your name and address filled in below and the original sales slip indicating the purchase price and the date of purchase to:

24 Month Performance Guarantee

Braun Products
Gillette Canada Company
4 Robert Speck Parkway
Mississauga, ON
L4Z 4C5

Sorry: No refunds will be made for any appliance damaged by accident, neglect or unreasonable use.

The Oral-B Professional Care™ must be shipped in its original packaging and be postmarked no later than 24 months after the date of purchase.

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