



# KneeRover™ Evolution Scooter

## User Manual

### **Congratulations on your purchase of the KneeRover™ Evolution Scooter!**

***The versatile KneeRover™ Evolution seated scooter provides a new level of freedom to users recovering from knee or leg injuries. The large gel padded seat and convenient front leg stand allows the patient to rest in comfort while elevating the non-weight bearing leg. Increase your mobility without the use of painful crutches.***

### **Operating Instructions:**

This KneeRover™ scooter is designed to maximize mobility during your recovery from foot/ leg surgery or injury.

It is important to set the seat and handlebar at a comfortable and safe operating position. The best position for the handlebar is approximately at hip level, with the seat at a comfortable height for propelling the KneeRover™ scooter and resting your injured leg.

Ensure that the hand brake is fully engaged and locked before mounting the KneeRover™ scooter. Place your hands on the handgrips and while seated, place the heel of the injured leg on the leg rest. Now release the parking brake and propel yourself using your good leg. Allow the KneeRover™ scooter to glide, and as it slows down, push again.

### **Safety First:**

The KneeRover™ Evolution scooter is designed for easy steering and stopping and should be used to carefully move about at slow, controlled speeds. When turning the scooter, you should have the brake applied and be moving at a slow, safe, controlled speed. It is easy to tip the scooter if you are making a turn with too much speed. Use caution and balance and operate at your own risk.

With routine use, nuts and bolts on the steering column and brake mechanism can and likely will become loose. Periodically check nuts and bolts after use and re-tighten as appropriate.

# Evolution Set Up Instructions:


Please follow the step by step instructions on the following pages to properly assemble your KneeRover™ Evolution scooter. The diagram below identifies the key components of the scooter to aid with assembly.

**Weight Capacity = 300 pounds**



## Step 1: Install Steering Column

Insert the steering column tiller into the tube at the front of the frame (as shown in the left photo). Lower the steering column towards the rear by lifting the clamp lever up to release the clamp from the steering column and pull the lever forward to disengage (as shown in the right photo). Push the lever to the left and downward simultaneously to disengage the security pin and gently guide the steering column downward. Please ensure that the handle bars are aligned with the wheels before proceeding.

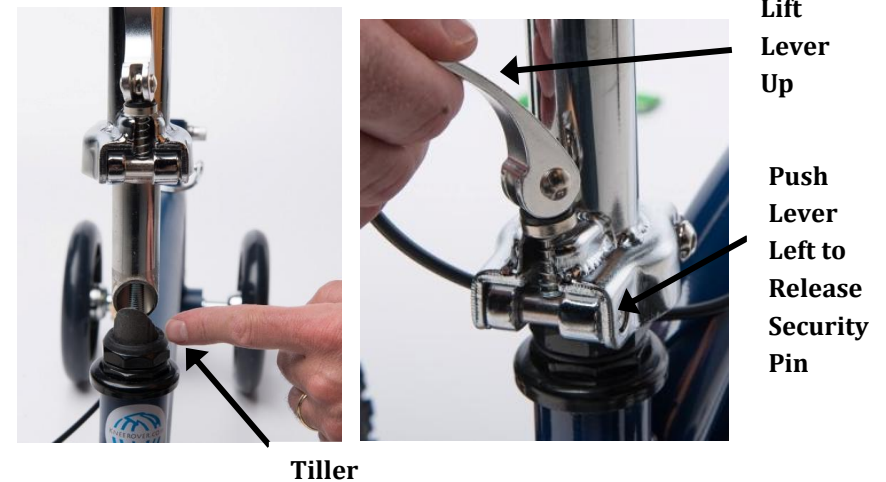

Next, take the provided allen wrench and tighten the exposed allen screw - turn clockwise until tight (as shown in the upper left photo on the next page). Raise the steering column, and again push the lever to the left and downward simultaneously to allow the security pin to lock into place. Raise the clamp lever upward and tighten by turning it clockwise. Finally, push clamp lever down firmly to secure the steering column (as shown in the upper right photo on the next page). Use the black knob to set the handlebars at a comfortable height for operation - typically hip level.

## Set Up Instructions Continued

### Tighten Allen Screw



### Step 3: Brake Operation

To engage the hand brake, squeeze the brake lever on handle bar.

The KneeRover™ Evolution scooter comes equipped with a brake lock / parking brake. Simply squeeze brake handle firmly and press the button to engage the lock. To release the lock, squeeze the brake handle again.



**Brake Lock**

### Step 2: Install Seat & Attach Basket

The KneeRover™ Evolution scooter is designed to be used for right or left leg injuries. To setup the seat and/or adjust the height, remove the locking pin and place the seat at the desired height then reinsert the pin. Pull the clamp lever downward so the end of the lever is pointing towards the ground to lock the seat in the desired position (rotate the lever counter clockwise to loosen and clockwise to tighten the pressure on the lever). Removable basket attaches to the steering column.



**Clamp  
Lever**

**Locking  
Pin**

### Ongoing Maintenance

With routine use, nuts and bolts on the steering column and brake assembly can become loose. To ensure safe operation of this KneeRover™ scooter it is important to periodically inspect the knee walker for loose nuts and bolts and re-tighten as appropriate.

## Who we are:

We are committed to providing innovative and affordable mobility solutions to get you moving as soon as possible after injury or surgery and for the duration of your recovery – free to go about your busy, active life – and have a little fun doing it. We are all about leaving the painful crutches behind and giving you the Freedom to Move.

The KneeRover™ Team is Driven to Change Your Outcome!

Your Satisfaction is Our Top Priority.

## Need Parts or Accessories?

Add an accessory like a cup holder or headlight to enhance your experience on your new KneeRover™ scooter. Keep your scooter operating in top shape with original manufacturer parts. Visit our exclusive online partner: [kneewalkershop.com](http://kneewalkershop.com), to explore our full range of parts and accessories.

## Need Help – Contact Us:

Phone: 1-877-737-9877

Web: [www.kneerover.com](http://www.kneerover.com)



## Limited Lifetime Warranty:

This warranty applies only to the original purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. Your KneeRover™ brand product is warranted to be free of defects in materials and workmanship for a period of one year from the date of purchase. The side frames and cross braces are warranted for the lifetime of the KneeRover™ scooter.

This warranty does not cover device failure due to owner misuse, negligence, or normal wear and tear. The warranty does not extend to non-durable components, which are subject to normal wear and need periodic replacement. These components include, but are not limited to: wheels, grips, brakes and upholstery.

The warranty shall not apply for damages or defects caused by unauthorized parts, attachments or repairs, negligence, accident, or any circumstances beyond the control of KneeRover™, such as improper operation, maintenance or storage.

If you have any questions regarding your KneeRover™ product or this warranty, please contact our Customer Care Team. We are available, for your convenience, toll free at 1-877-737-9877, Monday – Friday from 9 am to 5 pm Eastern.