

# RU+7

## User manual



---

## Warning


Dear user:

Thanks for using our company's products, in order to fully make use of the products, we sincerely suggests follow:

- ◆◆ Please read the instruction in detail and keep it well for further reading and consulting.
- ◆◆ Please follow the instruction to fit and operate the machine correctly.
- ◆◆ Please do not remove and change any accessories of the machine.
- ◆◆ It forbids others to open and remove the mainframe except the authorized people in the company.
- ◆◆ It forbids others to change and maintain the machine except the professional people.
- ◆◆ Please contact us if there is any trouble during operation.

**This product does not belong to medical, please don't use for disease treatment!**

## Parts Introduction

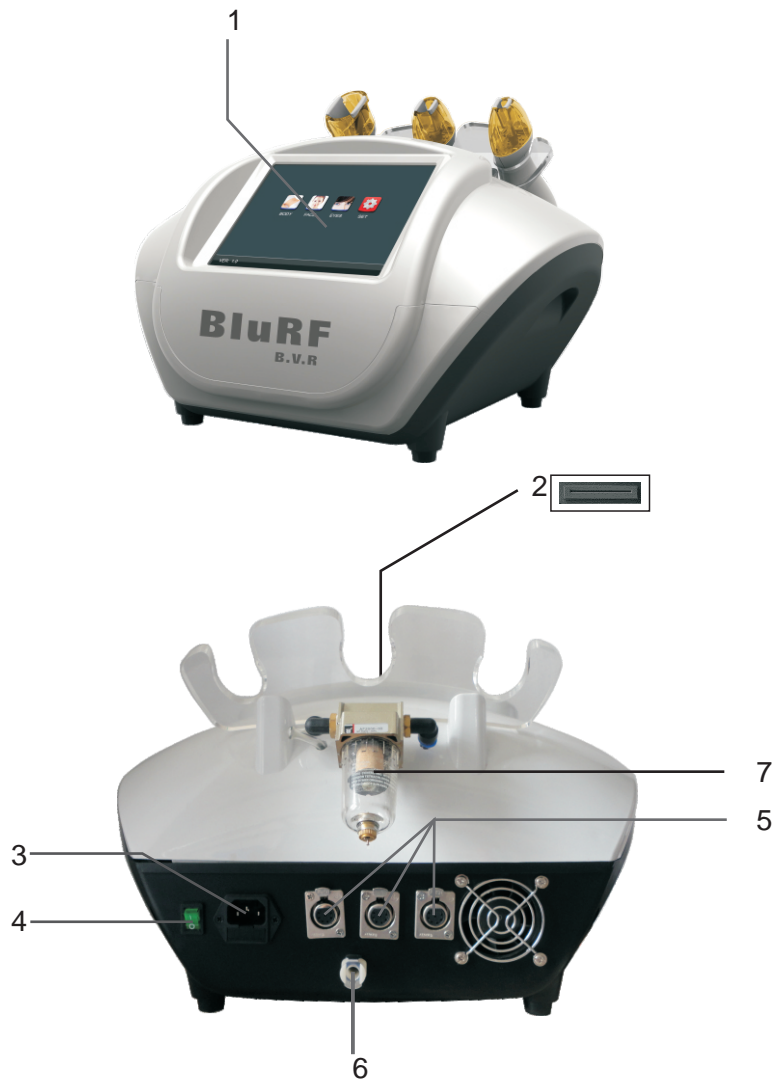

- 1.LCD screen
- 2.IC card socket
- 3.Power supply line jack
- 4.General power switch
- 5.Bipolar RF Vacuum work handle jack  
(eyes, face, body work handle can be automatically identified.)
- 6.Vacuum work handle public side jack

## Accessories



Bipolar RF vacuum work handle for eyes




Bipolar RF vacuum work handle for face



Bipolar RF vacuum work handle for body

## Basal Operation

- Put the accessories well and connect them to the host.
  - Insert Bipolar RF vacuum work handle for eyes, face and body to [5];
  - While using Bipolar RF vacuum work handle for eyes, face and body, please insert the transparent tube to [6];
- Insert power supply line to [3], turn on the General power switch [4], Insert the IC card to [2] correctly, the system will enter the interface as (Fig 1) shows.
- The system automatically recognizes the work handle.

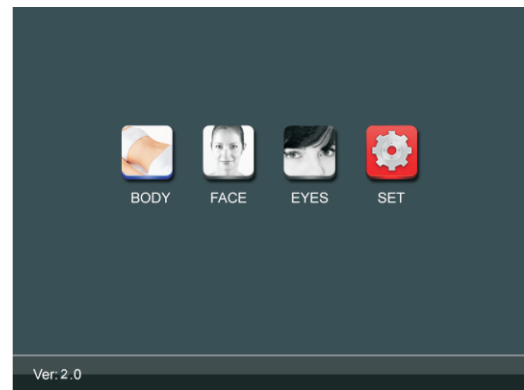


(Fig 1)

## Body treatment


- Click "BODY" as (Fig 1) shows, enter interface as (Fig 2).
- Click "🕒" to enter interface as (Fig 3) shows. Click "←" and "→" to set work time, 20 minutes for body treatment.
- Click "Vacuum" to enter interface as (Fig 4) shows.
  - Click "Mode" to set vacuum mode, Click "←" and "→" to regulate: P1~P6.
  - Click "Intensity" to set vacuum intensity, Click "←" and "→" to regulate: 10kPa~80kPa. When it is 0kPa means that close the vacuum function.
- Click "Bipolar RF" to enter interface as (Fig 5) shows. Click "←" and "→" to regulate RF output intensity: 1(weak)~10(strength). When it is 0 means that close the RF function.
- Click "PDT" to open or close the Photodynamic Therapy.
- Smear the essential oils in the treatment area, set work time, vacuum mode, vacuum intensity and bipolar RF intensity, Click "⏸", the machine start working. If there is too much fat, you can increase the output intensity.
- The machine will pause automatically if the setting time is over. Click "🕒" to return to (Fig 1).
- Clean the leftovers with hot towel, use warm towel to clean the work head.



(Fig 2)

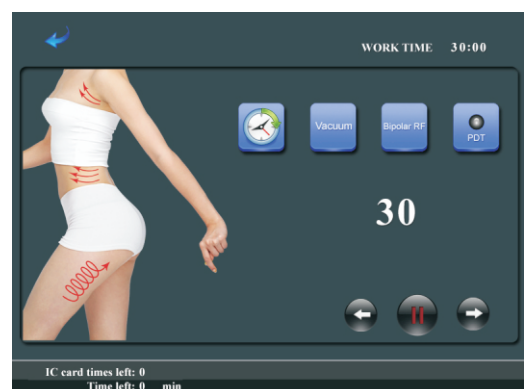

(Fig 3)

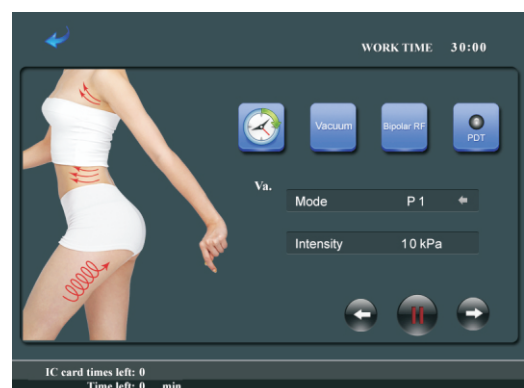

(Fig 4)

## Face treatment

1. Click "FACE" as (Fig 1) shows, enter interface as (Fig 6).
2. Set work time, 20 minutes for body treatment.
3. Click "Vacuum" to set vacuum mode and intensity.
  - (1) Click "Mode" to set vacuum mode: P1~P6.
  - (2) Click "Intensity" to set vacuum intensity: 10kPa~60kPa. When it is 0kPa means that close the vacuum function.
4. Click "Bipolar RF" to regulate RF output intensity:
  - 1(weak)~10(strength). When it is 0 means that close the RF function.
5. Click "PDT" to open or close the Photodynamic Therapy.
6. Smear the essential oils on skin around face, Click "⏸", the machine start working.
7. Beautician hold the handle and operating based on the direction of the arrow in the picture. from inside to outside, operating on both side, right and left, one side is over then to start the other side, and must reduce the intensity according to the less fat on forehead. if you want to pause during using, please click "⏪".  
NOTE :Do not let essential oils into eyes .
8. The machine will pause automatically if the setting time is over.
9. Clean the leftovers with hot towel, use warm towel to clean the work head.




(Fig 5)




(Fig 6)

## Eyes treatment

1. Click "FACE" as (Fig 1) shows, enter interface as (Fig 7).
2. Set work time, 15 minutes for body treatment.
3. Click "Vacuum" to set vacuum mode and intensity.
  - (1) Click "Mode" to set vacuum mode: P1~P6.
  - (2) Click "Intensity" to set vacuum intensity: 10kPa~50kPa. When it is 0kPa means that close the vacuum function.
4. Click "Bipolar RF" to regulate RF output intensity:
  - 1(weak)~10(strength). When it is 0 means that close the RF function.
5. Click "PDT" to open or close the Photodynamic Therapy.
6. Smear the essential oils on skin around eyes, Click "⏸", the machine start working.
7. Beautician hold the handle and operating based on the direction of the arrow in the picture. and must reduce the intensity according to the less fat around eyes.  
NOTE :Do not let essential oils into eyes .
8. The machine will pause automatically if the setting time is over.
9. Clean the leftovers with hot towel, use warm towel to clean the work head.

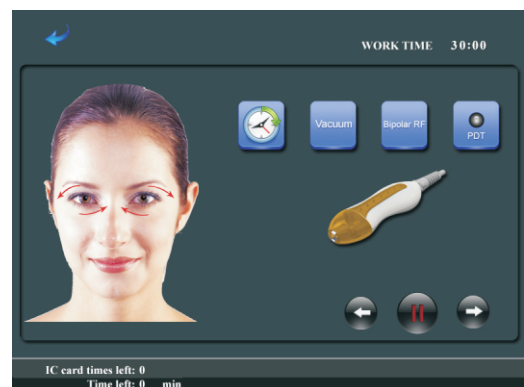

(Fig 7)

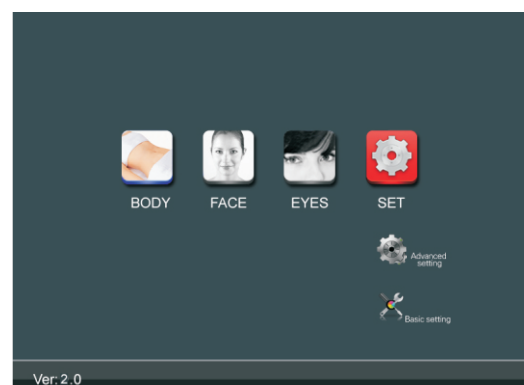

(Fig 8)

---

---

## Caution for using RF

- 1.Beautician should wear insulated gloves during operation.
- 2.Do not wear any aglet, both for the beautician and the users.
- 3.The electrode probe must be fully in contact with skin.
- 4.Do not make the RF electrode stay on any part without moving. Keep it moving on skin slowly.
- 5.Skin around eyes and forehead is very thin, the RF output intensity should be relatively weakened.

## Others

Click “set” as the(Fig 1) shows. Enter interface as (Fig 8):  
Click“ basic setting ”, enter the selecting language interface.

**Notes “advanced setting” is for manufacturer only.**

## Forbidden group

- 1.People with heart disease or high blood pressure, or configured cardiac pacemaker.
- 2.Patients with acute inflammation, asthma, deep vein thrombosis, thyroncus, cancer.
- 3.People with hemorrhagic disease, trauma or who is bleeding.
- 4.Pregnant women
- 5.Medical Plastic parts, or parts with artificial in filling inside.
- 6.People with metal inside the body(include metallic tooth)
- 7.Patient with skin inflammation or with edema.
- 8.People with an abnormal immune system.
- 9.Numb or insensitive to heat.

---

---

## Technical parameters

Input: AC110V/60Hz  AC220V/50Hz

Power:  $\leq 150\text{VA}$

Output RF: 2.6MHz

## Accessories list

1.Host	1set
2.Power supply line	1pc
3.User manual	1copy
4.Bipolar RF vacuum work handle for eyes	1pc
5.Bipolar RF vacuum work handle for face	1pc
6.Bipolar RF vacuum work handle for body	1pc
7.Accessory frame	1pc
8.Hexagonal screws	2pcs
9.Hexagon wrench	1pc
10.Sealing ring	15pcs