


XZ030

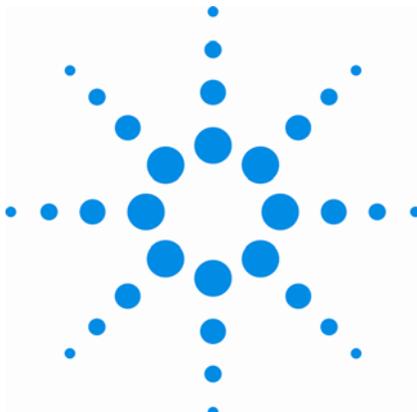

**AGILENT ACQIRIS
USER MANUAL
FAMILY OF
ANALYZERS**

Models covered:

AP235/AP240 with SSR

AP235/AP240 with Peak^{TDC}

AP101/AP201



June 2007

The information in this document is subject to change without notice and may not be construed in any way as a commitment by Agilent Technologies Inc. While Agilent Technologies makes every effort to ensure the accuracy and contents of the document it assumes no responsibility for any errors that may appear.

All software described in the document is furnished under license. The software may only be used and copied in accordance with the terms of license. Instrumentation firmware is thoroughly tested and thought to be functional but it is supplied "as is" with no warranty for specified performance. No responsibility is assumed for the use or the reliability of software, firmware or any equipment that is not supplied by Agilent Technologies Inc. or its affiliated companies.

Any versions of this manual which are supplied with a purchased product will be replaced at your request with the latest revision in electronic format. At Agilent Technologies we appreciate and encourage customer input. If you have a suggestion related to the content of this manual or the presentation of information, please contact your local Acqiris representative or Acqiris Technical Support (support@acqiris.com) or come visit our web site at <http://www.acqiris.com>.

Trademarks: product and company names listed are trademarks or trade names of their respective companies

Acqiris Headquarters:

**Agilent Technologies SA
12, chemin des Aulx
CH-1228 Plan-les-Ouates
Geneva
Switzerland**

Tel: +41 22 884 32 90

Fax: +41 22 884 32 99

Acqiris USA:

**Agilent Technologies Inc.
P.O. Box 2203
Monroe, NY 10949
USA**

Tel: 845 782 6544

Fax: 845 782 4745

Acqiris Asia-Pacific:

**Agilent Technologies
Australia Pty Ltd
347 Burwood Highway,
Forest Hill VIC 3131
Australia**

Tel: +61 3 9210 2890

Fax: +61 3 9210 5929

© Copyright Agilent Technologies Inc. June 2007

CONTENTS

1. OUT OF THE BOX	6
1.1. Message to the User	6
1.2. Using this Manual	6
1.3. Conventions Used in This Manual	6
1.4. Disclaimer and Safety	6
1.5. Warning Regarding Medical Use	7
1.6. Packaging and Handling	7
1.7. Warranty.....	7
1.8. Warranty and Repair Return Procedure, Assistance and Support	8
1.9. System Requirements.....	8
1.10. Transport & Shipping.....	9
1.11. Maintenance	9
1.12. Cleaning	9
1.13. Disposal and Recycling.....	9
2. INSTALLATION	10
2.1. Installing the Software under Windows	10
2.1.1. Warnings	10
2.1.2. Multiple Versions.....	10
2.1.3. Installation.....	10
2.2. Installing the Software for Linux	20
2.2.1. Kernel Mode Driver Compilation	21
2.2.2. Special cases.....	22
2.2.3. Environment variables for the Firmware.....	22
2.3. Installing the Hardware	22
2.4. After Restarting.....	23
2.4.1. Windows 2000	23
2.4.2. Windows XP	23
2.5. LabVIEW RT	25
2.6. Installing the IVI-COM/C Driver.....	25
2.7. Distribution for Windows 2000/XP and Linux	26
3. PRODUCT DESCRIPTION	27
3.1. Overview	27
3.2. Channel Input.....	29
3.2.1. Coupling & Impedance	29
3.2.2. Input Protection.....	29
3.2.3. Mezzanine Front-end.....	29
3.2.4. Bandwidth and Rise Time	29
3.2.5. Input Voltage and Offset.....	29
3.2.6. Vertical Resolution.....	30
3.2.7. DC Accuracy	30
3.3. Trigger.....	30
3.3.1. Trigger Source.....	30
3.3.2. Trigger Coupling.....	30
3.3.3. Trigger Level.....	30
3.3.4. Trigger Slope.....	31
3.3.5. External Trigger Output	31
3.3.6. Trigger Status	31
3.3.7. Trigger Veto and Timeout (AP101/AP201 ONLY)	31
3.4. Data Acquisition – Common.....	32
3.4.1. Timing.....	32
3.4.2. Sampling Rate	32
3.5. Data Acquisition - Digitizer Mode.....	32
3.5.1. Acquisition Memory and Time Base Range - Digitizer Mode.....	32

3.5.2.	Pre- and Post-Trigger Delay - Digitizer Mode	32
3.5.3.	Single and Sequence Acquisitions - Digitizer Mode	33
3.6.	Data Acquisition – Analyzer mode	33
3.6.1.	Random and Synchronized Triggers	35
3.6.2.	Time Base Range	35
3.6.3.	Dual-Bank Memory & Buffered Sequence Acquisition	35
3.6.4.	Buffered Acquisition Mode with Data Readout	36
3.7.	AP235/AP240 Data processing	36
3.7.1.	User Gates	36
3.7.2.	Threshold gates	36
3.7.3.	Time stamp information	37
3.7.4.	Peak detection	37
3.7.5.	Interpolation of Peak Position and Amplitude	37
3.7.6.	Histogramming Peak Data	37
3.8.	AP101/AP201 Data processing	38
3.8.1.	User Gates	38
3.8.2.	Simple peak detection	38
3.8.3.	Peak detection with Hysteresis	38
3.8.4.	Peak Interpolation	38
3.9.	External Clock and Reference	39
3.10.	Front Panel Inputs and Controls	39
3.11.	Internal Calibration	39
3.12.	Electrical, Environmental and Physical Specifications	40
3.12.1.	Electrical	40
3.12.2.	Environmental and Physical	40
4.	RUNNING THE AP_SSRDEMO APPLICATION	41
4.1.	Getting Started with AP_SSRDemo	41
4.2.	Editing Fields	42
4.3.	Digitizer Control Panel and Functions	42
4.3.1.	Input Voltage Scale, Offset, Coupling and Bandwidth	43
4.3.2.	Digitizer Timebase and more	43
4.3.3.	Operation Mode	43
4.3.4.	Trigger	43
4.3.5.	External Clock	44
4.3.6.	Control IO Settings	44
4.3.7.	Acquisition Mode	44
4.3.8.	Display options	44
4.4.	Analyzer Control Panel and Functions	44
4.4.1.	Analyzer Timebase and Memory	45
4.4.2.	User Gate definition	45
4.4.3.	Threshold Gate definition	45
4.4.4.	Peak ^{TDC} mode - Peak and Histogram definition	45
4.4.5.	Additional Controls & Time Stamp data	46
4.4.6.	Waveform storage	46
4.4.7.	Control panel reactivation	46
4.4.8.	Temperature display	47
4.4.9.	Zoom display control	47
5.	RUNNING THE APX01DEMO APPLICATION	48
5.1.	Getting Started with APx01Demo	48
5.2.	Editing Fields	49
5.3.	Digitizer Control Panel and Functions	49
5.3.1.	Input Voltage Scale, Offset, Coupling and Bandwidth	50
5.3.2.	Digitizer Timebase and more	50
5.3.3.	Trigger	51
5.3.4.	Acquisition Mode	51
5.3.5.	Persistence display	52
5.4.	Analyzer Control Panel and Functions	52

5.4.1.	Analyzer Timebase and Memory	52
5.4.2.	Gate definition.....	52
5.4.3.	Trigger Veto and Timeout.....	52
5.4.4.	Hysteresis Mode Peak Parameters	53
5.4.5.	Display of Peaks.....	53
5.5.	Display Window Functions.....	53
5.5.1.	Waveform storage	53
5.5.2.	Control panel reactivation.....	53
5.5.3.	Temperature display.....	54
5.5.4.	Zoom display control	54
6.	APPENDIX A: XA100 BNC INPUT OVERVOLTAGE PROTECTION	55
7.	APPENDIX B: XP102 FAN UNIT FOR THE AP FAMILY OF MODULES.....	56

1. Out of the Box

1.1. Message to the User

Congratulations on having purchased an Agilent Technologies Acqiris data conversion product. Acqiris Digitizers/Analyzers are high-speed data acquisition modules designed for capturing high frequency electronic signals. To get the most out of the products we recommend that you read this manual carefully. We trust the product you have purchased will meet with your expectations and provide you with a high quality solution to your data conversion applications.

1.2. Using this Manual

This guide assumes you are familiar with the operation of a personal computer (PC) running a Windows 2000/XP or other supported operating system. It also assumes you have a basic understanding of the principles of data acquisition using either a waveform digitizer or a digital oscilloscope.

The manual is divided into 5 separate sections. To understand the elements of operation for the module it is essential that you read them as appropriate.

- Chapter 1 ***OUT OF THE BOX***, describes what to do when you first receive your new Acqiris product. Special attention should be paid to sections on safety, packaging and product handling. Before installing your product please ensure that your system configuration matches or exceeds the requirements specified.
- Chapter 2 ***INSTALLATION***, covers all elements of installation and performance verification. Before attempting to use your Acqiris product for actual measurements we strongly recommend that you read all sections of this chapter.
- Chapter 3 ***PRODUCT DESCRIPTION***, provides a full description of all the functional elements of the Analyzer.
- Chapter 4 ***RUNNING THE AP_SSRDEMO APPLICATION***, describes the operation of an application that enables basic operation of Acqiris AP235/AP240 digitizers and analyzers in a Windows 2000/XP environment.
- Chapter 5 ***RUNNING THE APX01DEMO APPLICATION***, describes the operation of an application that enables basic operation of Acqiris AP101/AP201 digitizers and analyzers in a Windows 2000/XP environment

For information necessary for writing you own software to control Acqiris products you should refer to the **Programmer's Guide** and the **Programmer's Reference Manual**.

1.3. Conventions Used in This Manual

The following conventions are used in this manual:

-  This icon to the left of text warns that an important point must be observed.
- WARNING** Denotes a warning, which advises you of precautions to take to avoid being electrically shocked.
- CAUTION** Denotes a caution, which advises you of precautions to take to avoid electrical, mechanical, or operational damages.
- NOTE** Denotes a note, which alerts you to important information.
- Italic* text denotes a warning, caution, or note.
- Bold Italic** text is used to emphasize an important point in the text or a note
- `mono` text is used for sections of code, programming examples and operating system commands.

1.4. Disclaimer and Safety

The AP Series PCI Analyzer cards have been designed to operate in a standard PCI slot found inside most personal computers. Power for the modules is provided by plugging them into a free slot (refer to the installation procedure).



CAUTION: *Do not exceed the maximum input voltage rating! The maximum input voltage for 50 Ω input impedance is ± 5 V.*

1.5. Warning Regarding Medical Use

The AP Series PCI Analyzer cards are not designed with components and testing intended to ensure a level of reliability suitable for use in treatment and diagnosis of humans. Applications of these cards involving medical or clinical treatment can create a potential for accidental injury caused by product failure, or by errors on the part of the user. These cards are *not* intended to be a substitute for any form of established process or equipment used to monitor or safeguard human health and safety in medical treatment.



WARNING: *The modules discussed in this manual have not been designed for making direct measurements on the human body. Users who connect an Acqiris module to a human body do so at their own risk.*

1.6. Packaging and Handling

Your Analyzer is shipped with the following components:

- A small 3.5” compact disc in an Acqiris paper CD envelope that includes
 - 10 product user manuals in electronic form (Signal Analyzers, Streamer Analyzers, Family of 8-bit Digitizers, Family of 10-bit Digitizers, Family of 12-bit Digitizers, Family of Averagers, Family of Analyzers, Family of Time Counters, CC10X Family of CompactPCI Crates and CC121 CompactPCI Crate),
 - 1 Programmer’s Guide and 1 Programmer’s Reference Manual,
 - device drivers with sample software for different operating systems, environments and languages,
 - the AcqirisAnalyzers application, a demonstration program for the AC/SC Analyzer products,
 - the AcqirisLive application, a demonstration program for our digitizer and averager products,
 - the DemoSSR application, a demonstration program for the Acqiris AP235/AP240 Analyzers,
 - the DemoAPX01 application, a demonstration program for the Acqiris AP101/AP201 Analyzers,
 - the DemoTC application, a demonstration program for the Acqiris TC840/TC890 Time Counters,
 - product data sheets,
 - full installation procedures for use with Microsoft Windows, National Instruments LabVIEW RT, Wind River VxWorks, IVI-COM/C, and Linux software.
- A declaration of conformity
- Optional documentation such as a model-dependent document giving Specifications & Characteristics, a Calibration Certificate, or a Performance Verification

After carefully unpacking all items, inspect each to ensure there are no signs of visible damage. Also check that all the components received match those listed on the enclosed packing list. Agilent Technologies cannot accept responsibility for missing items unless we are notified promptly of any discrepancies. If any items are found to be missing or are received in a damaged condition please contact the Customer Support Center or your local supplier immediately. Retain the box and packing materials for possible inspection and/or reshipment.

1.7. Warranty

All Agilent Technologies Acqiris Digitizer products are warranted to operate within specification, assuming normal operation, for a period of three years from the date of shipment. It is recommended that yearly calibration be made in order to verify product performance. All repairs, replacement and spare parts are warranted for a period of 3 months. A 5-year repair warranty is available as an option.

Agilent Technologies endeavors to provide leading edge technology that includes the latest concepts in hardware and software design. As such software and firmware used with the products is under continual refinement and improvement. All software and instrument firmware is supplied “as is” with no warranty of any kind. Software and firmware is thoroughly tested and thought to be functional at the time of shipment. At Agilent Technologies’ discretion software and firmware may be revised if a significant operational malfunction is detected.

Products supplied but not manufactured by Agilent Technologies are covered solely by the warranty of the original equipment manufacturer.

In exercising this warranty, Agilent Technologies will repair or replace any product returned to the Agilent service center, within the warranty period. The warranty covers all defects that are a result of workmanship or materials. This excludes defects that are caused by accident, misuse, neglect, or abnormal operation.

The purchaser is responsible for returning the goods to the nearest Agilent service center. This includes transportation costs and insurance. Agilent Technologies will return all warranty repairs with transportation prepaid.

1.8. Warranty and Repair Return Procedure, Assistance and Support

Agilent Technologies Inc. acquired Acqiris SA and its product lines in December 2006. Before returning any Agilent Technologies Acqiris product for repair please contact your nearest Agilent Technologies Service Center.

You can find information about technical and professional services, product support, and equipment repair and service on the Web, see <http://www.agilent.com/find/service> (or <http://www.agilent.com/> and after selecting your country click on **Contact Us**). The service center will ask for your name, company, phone number and address, the model and serial numbers of the unit to be repaired, and a brief description of the problem.

Before issuing a Service Order they will ask you to communicate with us by phone or eMail so that we can learn as much as needed about the problems observed. If a unit returned under guarantee is found to be working normally and this procedure was not followed we reserve the right to charge you for the work done.

For your nearest customer support center contact Acqiris at 1-877-ACQIRIS in the USA, +41 22 884 32 90 in Europe or +61 3 9210 2890 in the Asia-Pacific region. Alternatively, please contact Acqiris Technical Support (support@acqiris.com) or come visit our web site at <http://www.acqiris.com>. The Agilent Technologies Support Centers can also help redirect you for any questions concerning the installation and operation of your equipment.

1.9. System Requirements

In order to obtain reasonable performance from your digitizer, Acqiris products need the following minimum PC System Requirements:

Processor: 150 MHz Pentium (higher recommended)

Memory: 64 Mbytes RAM. The previous number is a very rough estimate. Assuming that you are using AcqirisLive or an application of your own that operates on the acquired data it seems reasonable to ask for 10 times the total acquisition memory that you will be using at the same time in the application. Performance is likely to be degraded if less memory is available.

Display resolution: At least 800 x 600 pixels and 256 colors for use of AcqirisDemo programs.

Operating System: Microsoft Windows 2000/XP including 2003 Server, Wind River VxWorks 5.5.1, and Linux. The following Linux versions have compiled loadable kernel modules:

- RedHat Enterprise Linux Version 3 - Kernel Version 2.4.21-4.EL
- RedHat Enterprise Linux Version 3 - Kernel Version 2.4.21-4.ELsmp,
- Suse Linux 9.2 - Kernel version 2.6.8-24-smp,
- Debian Sarge 3.1 2.6.8-3-686,
- Debian 4.0 etch Kernel version 2.6.18-4-686,
- Scientific Linux 4.4 Kernel version 2.6.9-11.EL).

Linux kernel driver source code is available for compilation. Support for Windows 95/98/NT4 is included “as is” since these operating systems are no longer supported by Microsoft.

Hard Drive Space: 20 Mbytes Minimum

CD Drive (or any method to copy Acqiris Software installation files from CD to the hard drive such as LAN, floppy drive, etc.)

LabVIEW: The Acqiris LabVIEW driver is available for National Instruments LabVIEW versions 7.1 or 8.0

LabVIEW RT: The Acqiris LabVIEW RT driver is available for National Instruments LabVIEW RT version 7.1 or higher. The VISA driver must be version 3.0 or higher.

MATLAB: The Acqiris MEX interface can be used with MathWorks MATLAB 6.5 or a newer version.

Visual BASIC: The interface files and examples are available for Microsoft Visual Basic versions 5 or 6 and the interface files only for .NET.

Tornado: The example files are useable with Wind River Tornado 2.2.1

1.10. Transport & Shipping



CAUTION: Cards can be safely transported in their original shipping packages. DC cards can be transported when properly mounted in a CompactPCI crate. The transport of DP cards mounted in a PC is a more delicate issue. Because of their mass the cards can vibrate loose unless they are properly secured and braced. DP cards held only in the front and on the bottom should not be shipped in their PC. However, properly mounted DP cards with XP103 or XP105 fans can be sufficiently well held; the Adjustable retainer must be used so as to hold the card in place.

To package the instrument for shipping:

Step

Notes

1. Place the instrument in its original packaging materials.
 - If the original packaging materials are not available, use a professional packaging service. Contact your Agilent Service Center for more information.
2. Surround the instrument with at least 3 to 4 inches (8 to 10 cm) of its original packing material or bubble-pack to prevent the instrument from moving in its shipping container.
3. After wrapping it with packing material, place the instrument in its original shipping container or a strong shipping container that is made of double-walled corrugated cardboard with 159 kg (350 lb) bursting strength.
 - The shipping container must be large and strong enough to accommodate your instrument and allow at least 3 to 4 inches (8 to 10 cm) on all sides for packing material.
4. Seal the shipping container securely with strong nylon adhesive tape.
5. Mark the shipping container “FRAGILE, HANDLE WITH CARE” to help ensure careful handling.
6. Use the address obtained from your Agilent Technologies Service Center.
7. Retain copies of all shipping papers.



CAUTION: Damage can result if the original packaging materials are not used. Packaging materials should be anti-static and cushion the instrument on all sides. NEVER USE STYRENE PELLETS IN ANY SHAPE AS PACKAGING MATERIALS. They do not adequately cushion the instrument or prevent it from moving in the shipping container. Styrene pellets can also cause equipment damage by generating static electricity or by lodging in fan motors.

1.11. Maintenance

The cards do not require any maintenance. There are no user serviceable parts inside. A periodic calibration can be obtained on request.

1.12. Cleaning

Cleaning procedures consist only of exterior cleaning.

Clean the exterior surfaces of the module with a dry lint-free cloth or a soft-bristle brush. If any dirt remains, wipe with a cloth moistened in a mild soap solution. Remove any soap residue by wiping with a cloth moistened with clear water. Do not use abrasive compounds on any parts.

1.13. Disposal and Recycling

Electronic equipment should be properly disposed of. Acqiris Analysers/Digitizers and their accessories must not be thrown out as normal waste. Separate collection is appropriate and may be required by law.



2. Installation

This chapter describes how to install the Acqiris hardware and software for Windows 2000/XP, National Instruments LabVIEW RT, Linux, or Wind River VxWorks.



NOTE: For a first time installation we strongly recommend installing the software **before** inserting the hardware into the PC.

2.1. Installing the Software under Windows

2.1.1. Warnings

If Setup detects a previous installation of Acqiris software on your system, a warning screen will be displayed. It is recommended to exit Setup and uninstall older versions.

The installer from software releases prior to **Acqiris Software 2.0** installed the Digitizer Driver DLL files into the System directory. These will be removed by Setup. If you wish to keep the old installation on your system, you should exit Setup, and move all Acqiris driver files (acqiris*, acqrs* and acqir*) to some archive directory.

The DLL files will be installed into the bin subdirectory of the Acqiris software root, and the corresponding path will be added to the PATH environment variable.

2.1.2. Multiple Versions

With the software installation from **Acqiris Software 2.0** (or above), it is possible to keep multiple versions on the same system, but you must specify a different root directory (i.e. Install Folder). If you keep the same directory, Setup will overwrite your previous installation.

To go back to a previous version, you must change the PATH environment variable and reinstall the Kernel driver. Under Windows 2000/XP:

1. Copy the SYS file from <old_AcqirisSoftware_root>\bin\kernel to the Windows\System32\drivers directory.
2. Change the AcqirisDxRoot, AcqirisDxDir and PATH environment variables to the old root.
3. Reboot the computer.


2.1.3. Installation

Before installing the Acqiris hardware, you should complete the following steps to install the software for Windows 2000/XP.




NOTE: You will need administrator privileges to complete the software installation under Windows 2000/XP.


1. Insert the *Acqiris Software CD* into the CD-ROM drive of your computer. If the Autorun program does not start automatically (Autoplay disabled), you can start it manually, or navigate to the *AcqirisSoftware\Windows* folder in order to display the files included.
2. Choose **Install AcqirisSoftware for Windows 2000/XP** (or run Setup.exe from the *AcqirisSoftware\Windows* folder). After several seconds for initialization the following screen will appear.

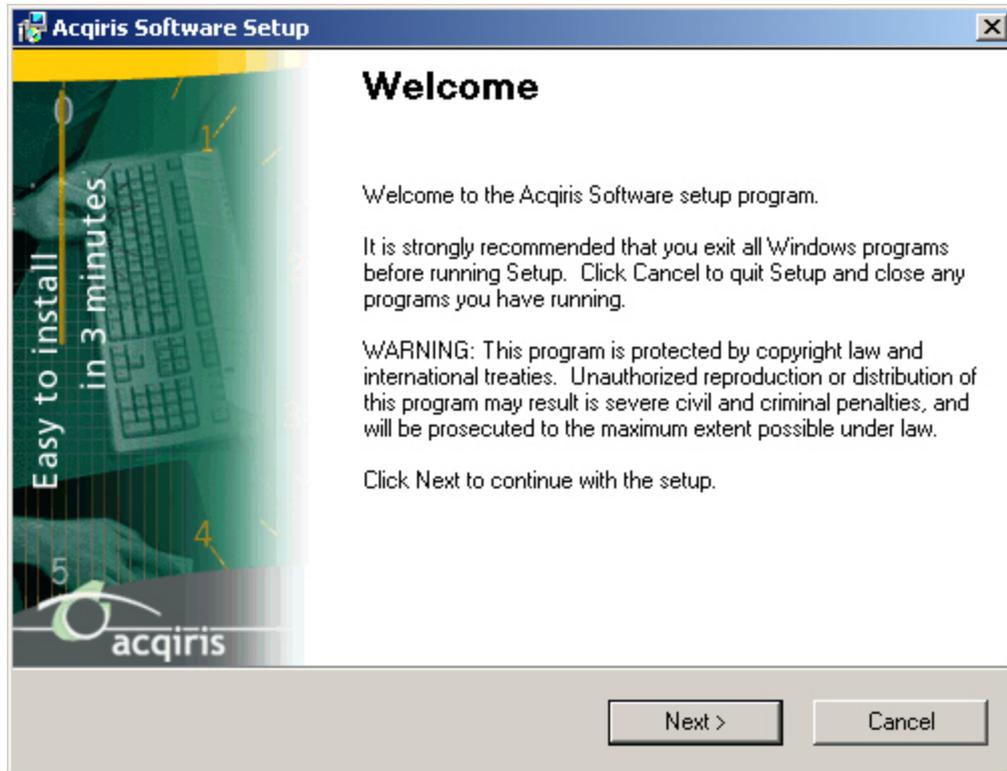


3. To proceed with the installation click "Next". The following screen will be shown.

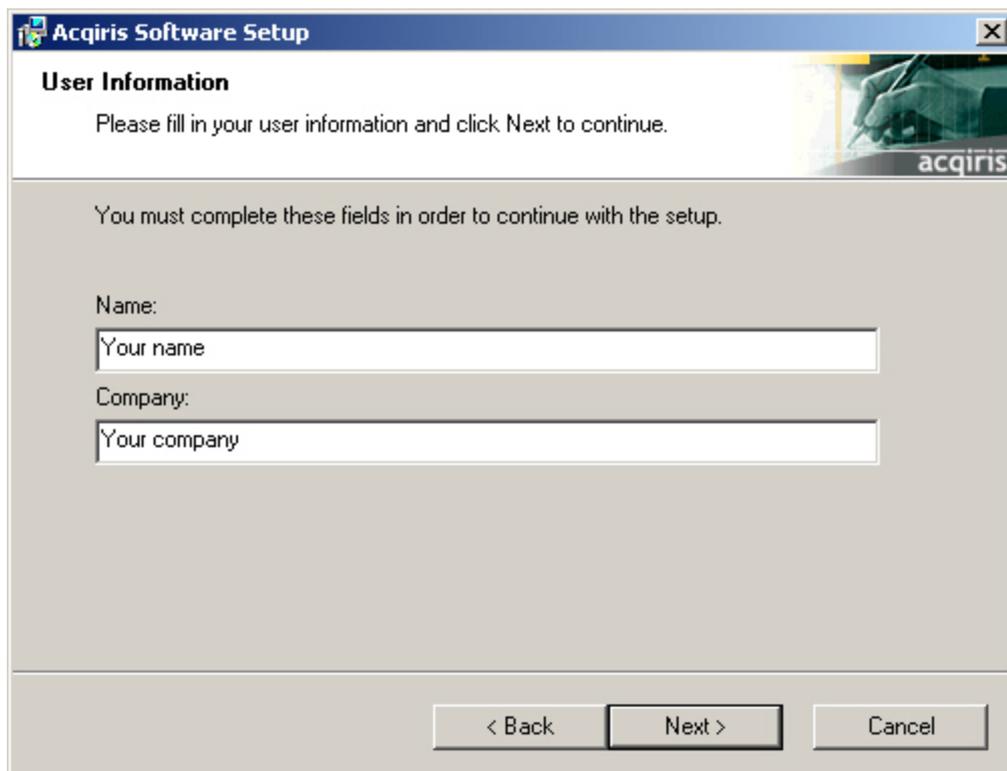


4. Please enter your user information and click "Next" to continue. If the program finds that there is still Acqiris software installed on your machine a warning panel (not shown) will appear. Otherwise, the next screen shows:

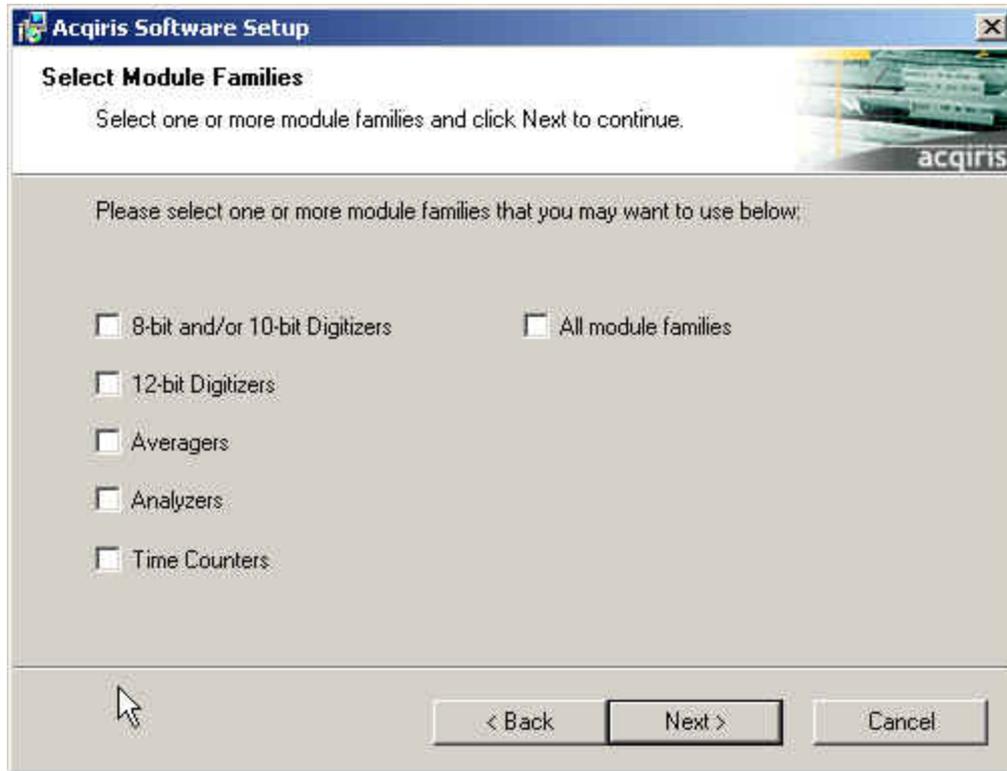


5. Pick the desired combination of module families and click "Next" to continue.

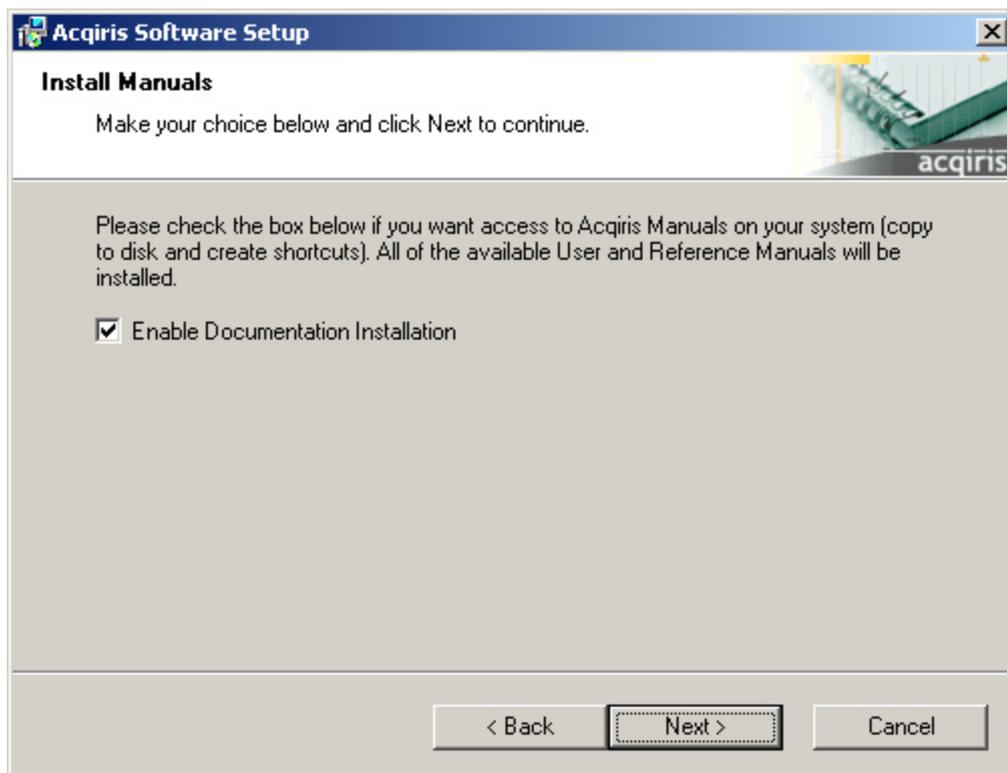


5. The screen above will normally allow the documentation to be installed. Remove the check from the box if you do not want online access to the manuals.
6. The next screen allows you to enable LabVIEW RT and/or Wind River VxWorks support. By default there will be none but if desired you can install any of them together with Windows support or without Windows support. Click "Next" to continue.

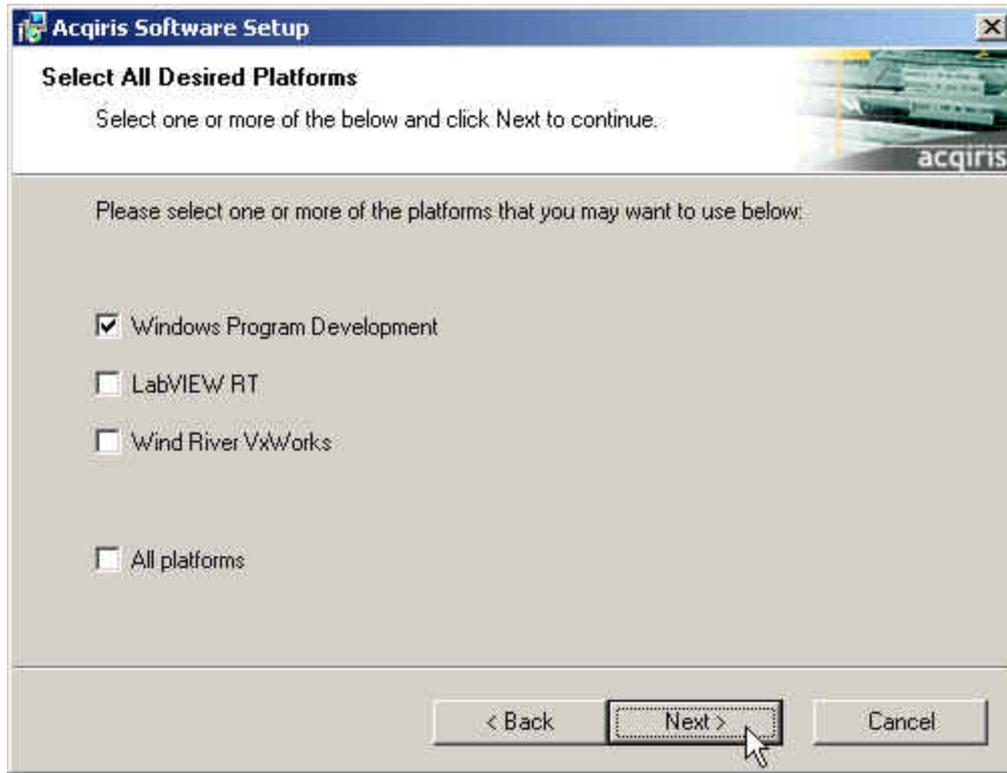


7. The screen below asks for the desired installation type. After having made your choice Click “Next” to continue.

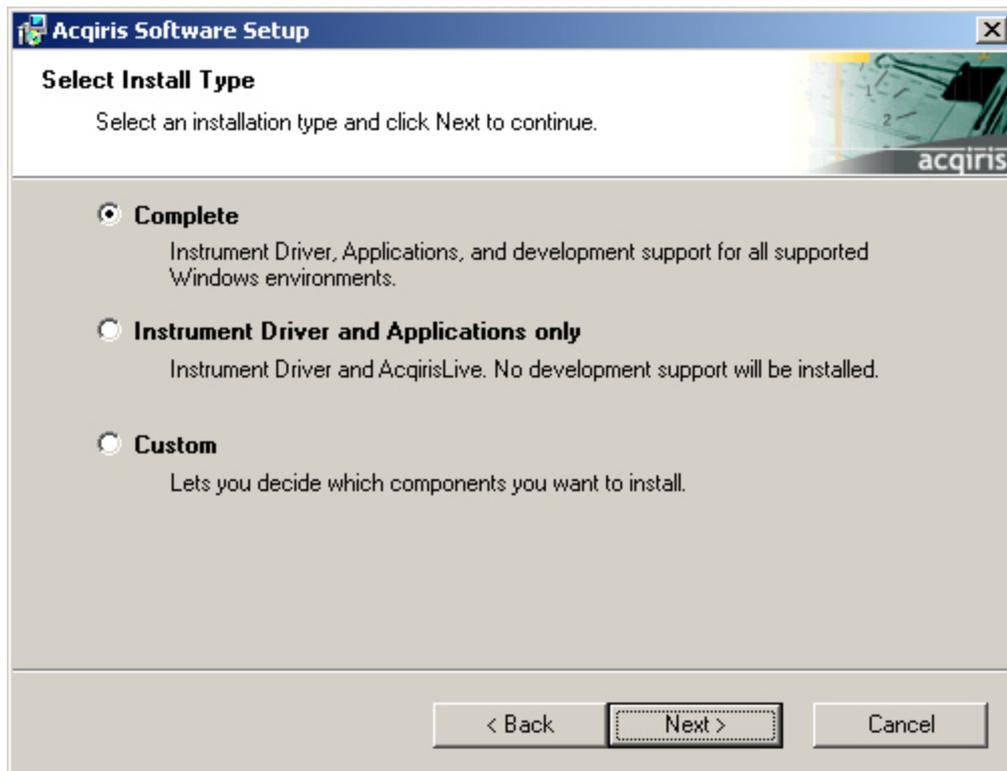


8. If you chose the Custom installation, the following screen will let you select each package individually. Note that the space indicated for LabVIEW, Firmware and UserManual packages is incorrect. The correct values are 4.1 MB, 21 MB, and 27 MB respectively.

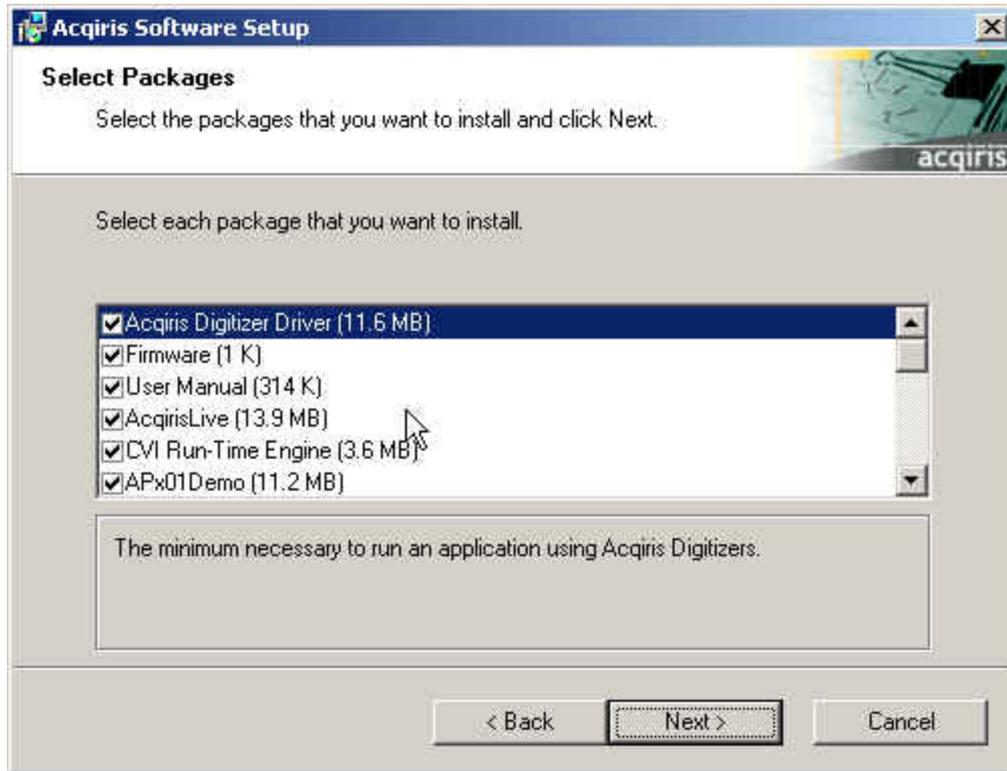


9. If MATLAB is installed on your machine, you will be asked to point the installer to the MATLAB root directory. You should do this if you want the installer to modify the standard startup.m file to incorporate the Acqiris adaptor.
10. Now that the elements of the installation have all been decided you will be prompted for the installation folder. This will be the root directory of the Acqiris software installation. If User Manuals (27 MB) and Firmware (21 MB) are loaded more space than indicated here will be required on the drive. For the case of a Tornado 2.2 installation the folder name should not contain any spaces.

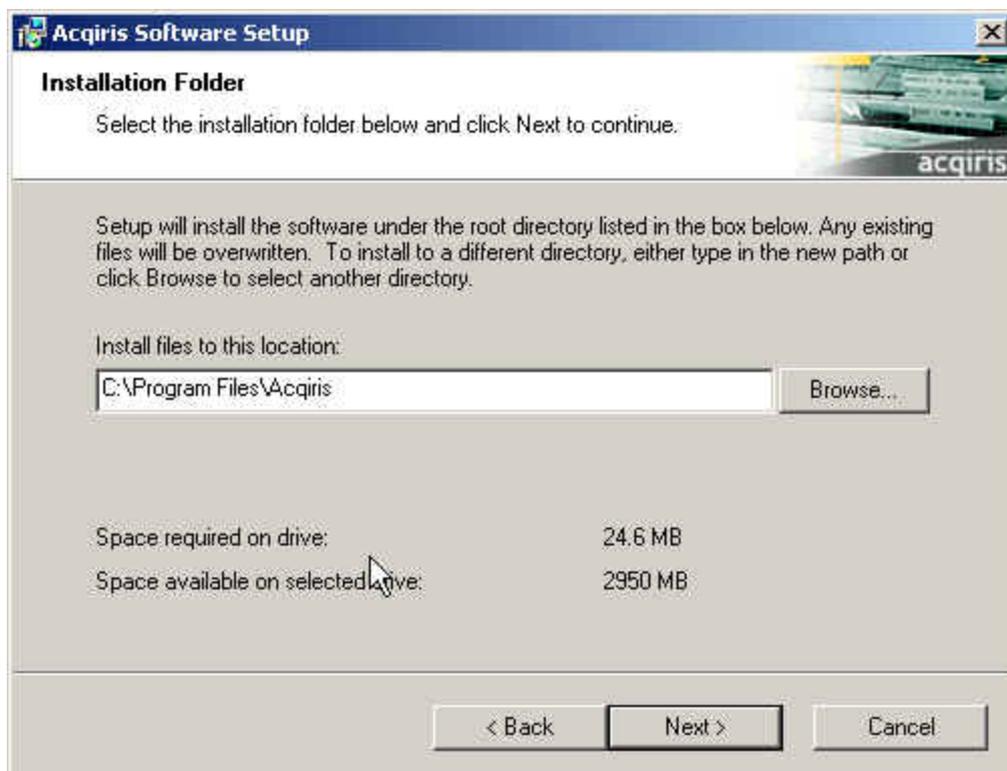


11. Furthermore, you should give a name to the shortcut folder. This is the menu entry under **Start** → Programs where you will find the shortcuts for AcqirisLive, manual(s), etc.

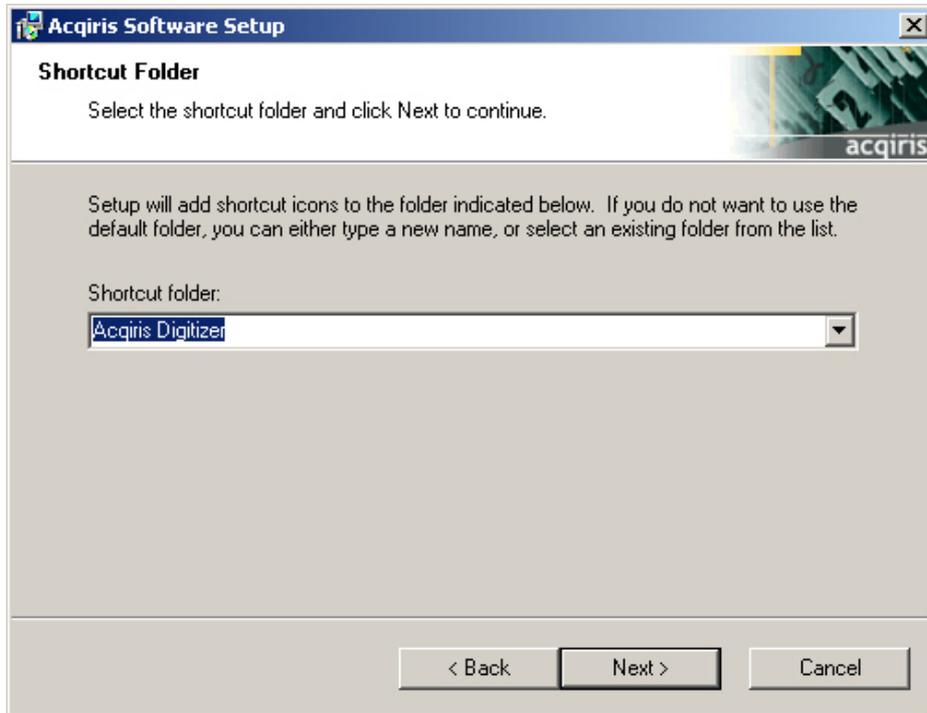


12. If you have enabled the installation of 12-bit Digitizers, Averagers, or Analyzers, the next screen lets you change the name of the directory where the FPGA firmware will be installed.

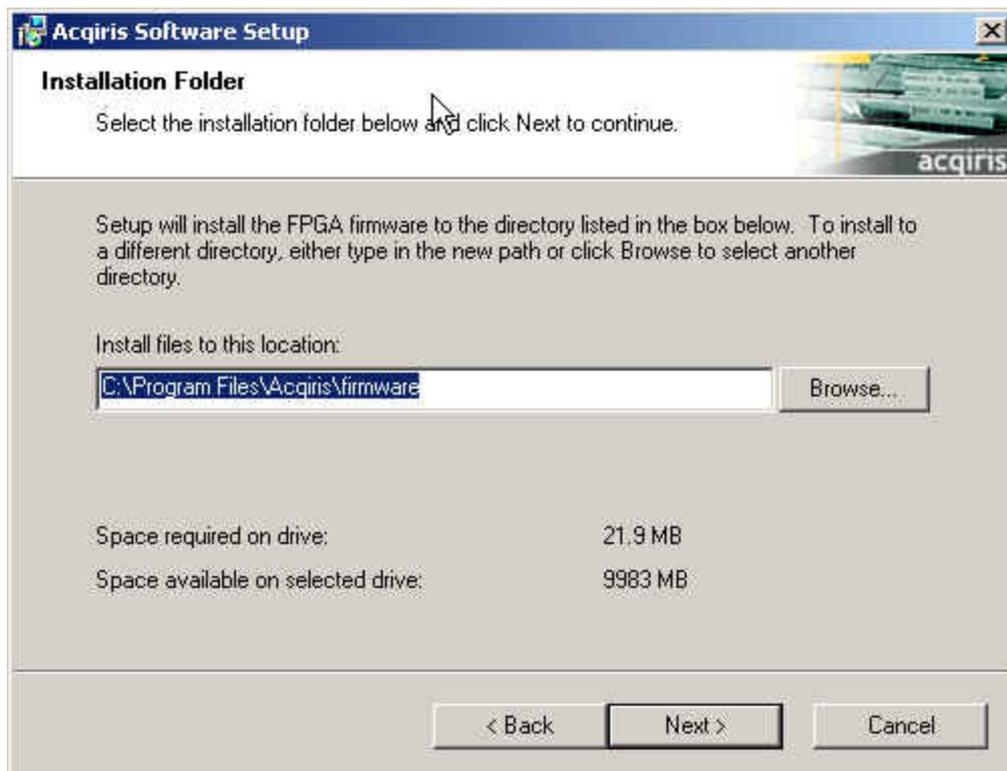


13. AcqirisLive needs the LabWindows/CVI 7.0 Run-Time Engine to run. If Setup has detected that a LabWindows/CVI Run-Time Engine is already installed on your system, it will ask you if you would like to install it locally for AcqirisLive anyway. If you are not sure about the version of the CVI Run-Time Engine on your system, it is recommended to install it locally. Click "Next" to continue.

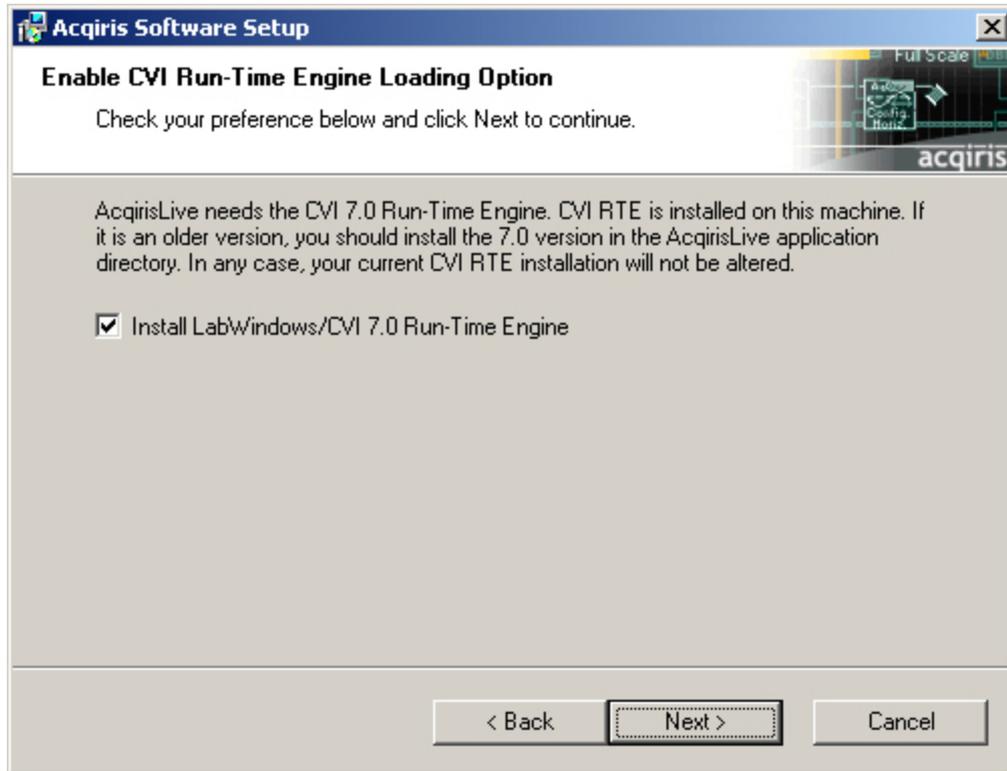


14. Depending on the install type, you may be asked which LabVIEW version format you want for the LabVIEW files. Select the appropriate format and click “Next” to continue.

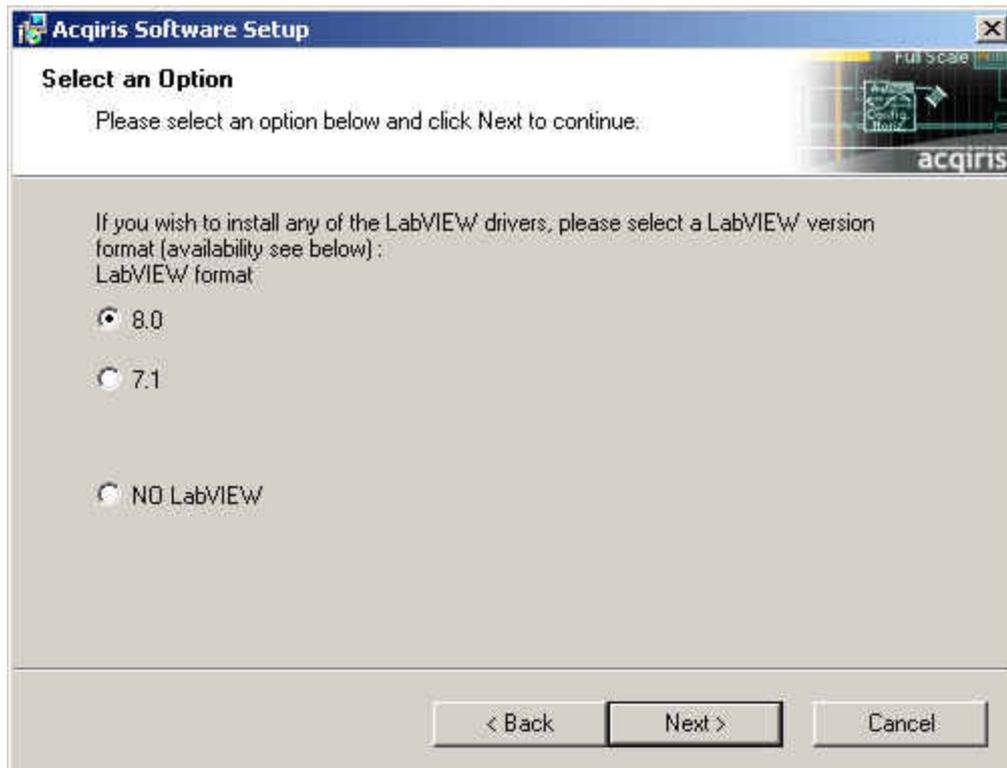

15. A summary will be shown to allow you to check what you have asked for

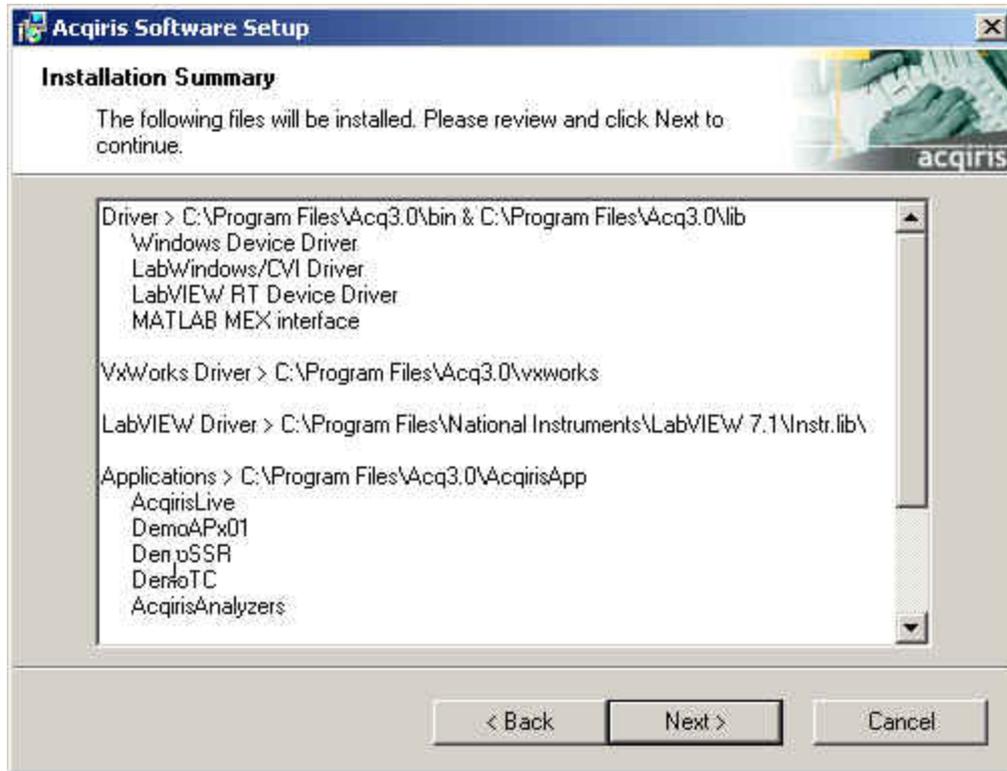


16. You are now ready to install. You may still go back to any previous screen to modify your selection. Click "Install" when ready.

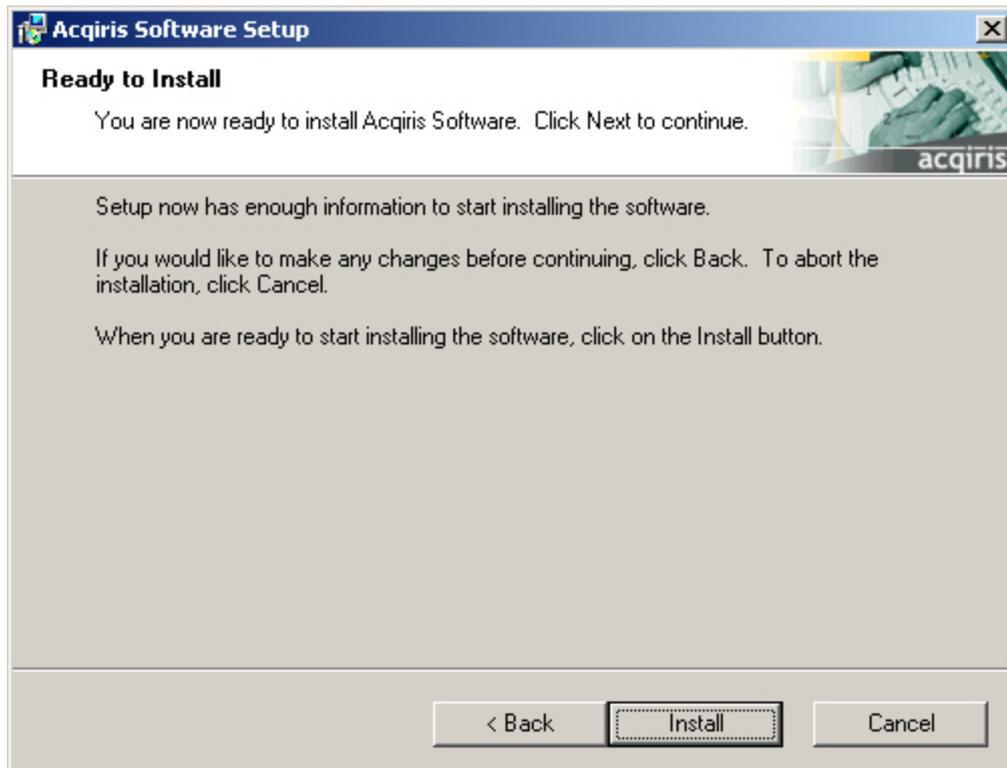

17. Setup will now copy the files and make the necessary changes to your system. When done, an information screen will be displayed. Please read this carefully.



18. Registration of your installation will help us provide you with better support. You will also be notified of updates and upgrades. All information submitted to us will be treated confidentially and never be disclosed outside the company.

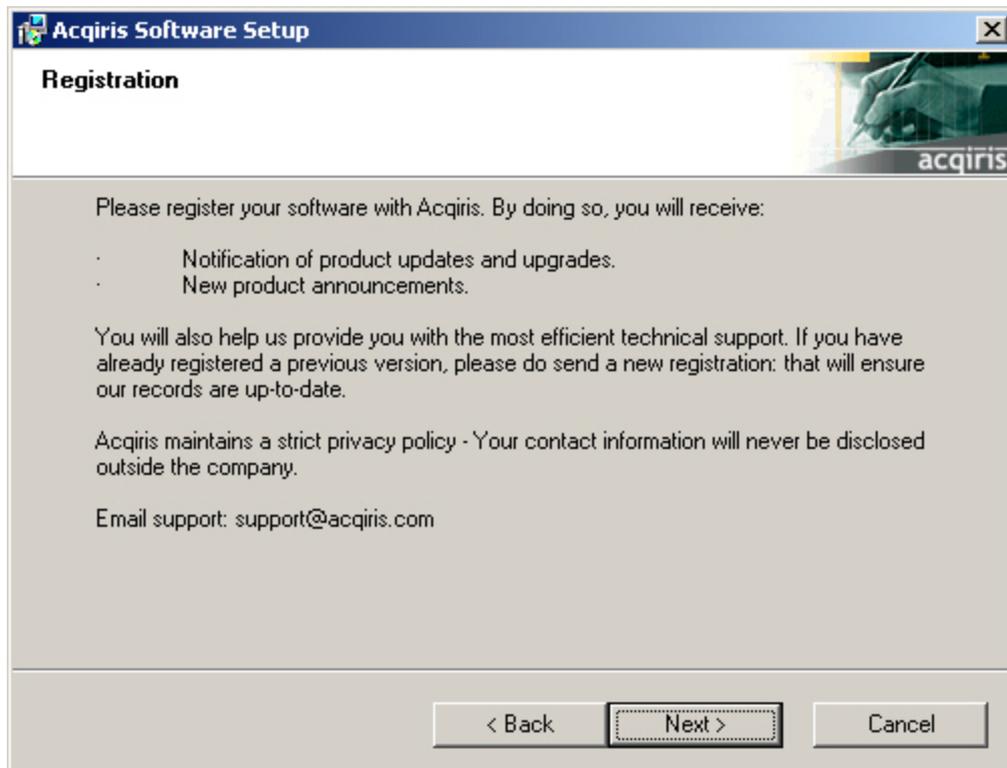




19. Setup will prepare a registration e-mail in your e-mail client application upon termination of the setup procedure. You can then decide whether or not you wish to send it. You may also add comments. Uncheck the box if you do not want to register your installation.



20. Click "Finish". The software installation is now complete.
21. You can now either accept the suggestion to restart the computer or you should shutdown your computer and proceed with the hardware installation.

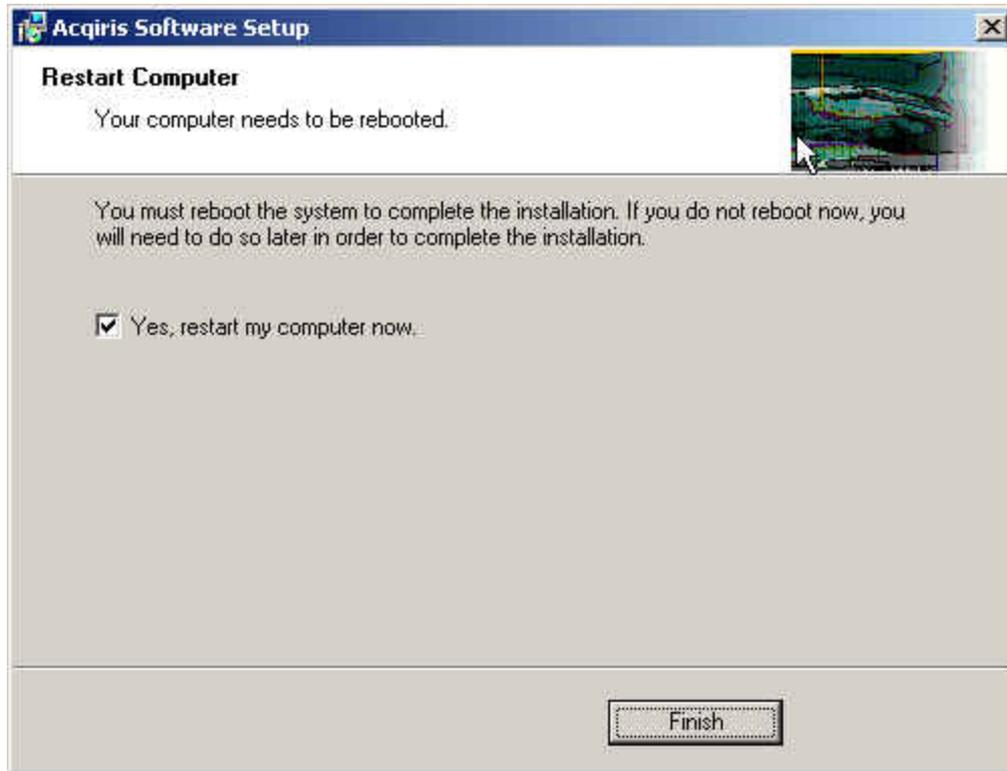

2.2. Installing the Software for Linux

The Acqiris Software is ready to install and run on Linux systems with any of the following:

- RedHat Enterprise Linux Version 3 - Kernel Version 2.4.21-4.EL
- RedHat Enterprise Linux Version 3 - Kernel Version 2.4.21-4.ELsmp,
- Suse Linux 9.2 - Kernel version 2.6.8-24-smp,
- Debian Sarge 3.1 2.6.8-3-686,
- Debian 4.0 etch Kernel version 2.6.18-4-686,
- Scientific Linux 4.4 Kernel version 2.6.9-11.EL)

The following tar files exist with the driver and library compiled with the **GNU gcc** version shown:

- AcqirisLinux-3.1a-gcc-3.3.tar.gz - compiled under Suse 9.2 with GNU gcc 3.3.
- AcqirisLinux-3.1a-gcc-3.4.tar.gz - compiled under SL 4.4 with GNU gcc 3.4.
- AcqirisLinux-3.1a-gcc-4.1.tar.gz - compiled under Debian etch with GNU gcc 4.1.

The Kernel Driver source code is available for recompilation for users with other Linux systems.

The appropriate tar file should be copied to a local directory (e.g. your home directory) and then unpacked by using the following command

```
tar xzf AcqirisLinux-3.1a-gcc-X.x.tar.gz
```

The resulting directory **AcqirisLinux** contains an install script **drv-install** and a graphical Demo program **demo/AcqirisDemo**.

Before installing the driver you have to put the correct version of the Kernel mode driver into acqirisPCI.o. Type

```
cd AcqirisLinux/lib/modules
```

and either

```
cp acqrsPCI.o[.n] acqrsPCI.o (Linux kernel 2.4)
```

or

```
cp acqrsPCI.ko[.n] acqrsPCI.ko (Linux kernel 2.6)
```

where [n] is the running version of the Linux kernel. This can be obtained by running the command

```
uname -r
```

The following versions exist:

```
acqrsPCI.o.2.4.21-4.EL for Red Hat Enterprise Linux version 3
```

acqrsPCI.ko.2.6.8-3-686	for Debian 3.1 'sarge'
acqrsPCI.ko.2.6.8-24-smp	for Suse 9.2
acqrsPCI.ko.2.6.9-11.EL	for Scientific Linux 4.4
acqrsPCI.ko.2.6.18-4-686	for Debian 4.0 'etch'

To install the driver and the load script to the system, you have to get super user privileges and execute the driver install script by typing,

```
cd AcqirisLinux
./drv-install add
```

You can check that the driver is loaded properly with **lsmod** or **dmesg**. The **drv-install** script has the following additional functionality:

drv-install rem to remove the driver and load script

drv-install res to restart the driver

If you intend to use an Averager, an Analyzer, or a 12-Bit Digitizer you must install the Firmware .bit files. You can either copy them from the Firmware directory of the CD-ROM or download them from our WEB site. They should be placed in the Firmware subdirectory of AcqirisLinux. In the case of the CD-ROM, this can be done by continuing the above command sequence with the following:

```
cp /mnt/CDRom/Firmware/*.bit Firmware/
```

2.2.1. Kernel Mode Driver Compilation

For some special kernels, for example SMP, the kernel mode driver will have to be recompiled on the target system. All files needed to compile a new kernel mode driver are in AcqirisLinux/linuxdriverpci for Kernel 2.4 or in AcqirisLinux/linux2.6driverpci for Kernel 2.6

The makefile can be used for compilation under a standard compatible distribution. To compile the driver, the kernel header files need to be installed. Furthermore, depending on the Linux distribution, the kernel source files may need to be installed. The makefile will try to guess where they are located. If it fails, see instructions inside the makefile itself to specify manually where these kernel header files are located.

To compile a new kernel under **linux-2.4**, issue the following commands:

- **cd linuxdriverpci**
- **make clean** to remove all *.o files in AcqirisLinux/linuxdriverpci.
- **make** to generate a new kernel mode driver acqrsPCI.o and also copy it to the directory AcqirisLinux/lib/modules where the install script (drv_install) can access it.
- **drv_install rem** to remove the previous installed kernel mode driver.
- **drv_install add** to install the new driver.

To compile a new kernel under **linux-2.6**, issue the following commands:

- **cd linux2.6driverpci**
- **make clean all** to generate a new kernel mode driver acqrsPCI.ko from scratch.
- **make install** to copy the kernel module where it should reside.

Note that this version of the loadable kernel module had been tested on Linux kernel versions up to 2.6.18. However, starting with kernel 2.6.17, a few kernel functions relating to the device classes are no longer available to proprietary modules. Thus automatic creation of the /dev/acqrsPCI node has to be done manually in the /etc/rc.local file.

For Debian users, in order to compile the kernel module on a system, a few packages that match the current kernel-image must be installed. For instance, on a Debian 3.1 sarge system 686 the required packages are:

- kernel-image-2.6-686
- kernel-image-2.6.8-3-686
- kernel-headers-2.6-686
- kernel-headers-2.6.8-3
- kernel-headers-2.6.8-3-686
- kernel-kbuild-2.6-3

2.2.2. Special cases

If you are running a Linux distribution that doesn't use the standard paths for the load scripts, you can load the driver with the command

```
cd AcqirisLinux
./drv-install load
```

Copy the driver to the module directory as follows:

```
cp lib/modules/acqrsPCI.o[.n] /lib/modules/{kv}/ACQIRIS/acqrsPCI.o
```

where [.n] is the optional designation of the desired version of the acqrsPCI driver as given above and {kv} is the appropriate system kernel version (obtained with **uname -r**, i.e. 2.4.20).

Then add the following command to **/etc/rc.local**.

```
/sbin/insmod -f /lib/modules//{kv}/ACQIRIS/acqrsPCI.o
```

For systems using **udev** (the dynamic **/dev/** hierarchy) it may be convenient to automatically create the **/dev/acqrsPCI** node at boot. To do this, add the following command to the **/etc/rc.local** file:

```
/sbin/mknod -m 666 /dev/acqrsPCI c 124 0
```

2.2.3. Environment variables for the Firmware

Automatic loading of the firmware needed by 12-bit, analyzer, and averager modules relies on the environment variable **AcqirisDxDxDir** pointing to the directory containing the file **AqDrv4.ini** which in turn points to the directory containing the Firmware **.bit** files. Therefore, assuming that your Firmware is in **/usr/local/AcqirisLinux/Firmware** and that **AqDrv4.ini** is in **/usr/local/AcqirisLinux/demo** then you must edit **AqDrv4.ini** so that it contains the line

```
fpgaPath=/usr/local/AcqirisLinux/Firmware
```

Then, if your shell is **csh** or **tcsh** modify the **/etc/csh.login** file to contain the line

```
setenv AcqirisDxDxDir /usr/local/AcqirisLinux/demo
```

or, if your shell is **bash**, **ksh**, **zsh** or **sh**, modify the **/etc/profile** file to contain the lines

```
AcqirisDxDxDir=/usr/local/AcqirisLinux/demo
```

```
export AcqirisDxDxDir
```

2.3. Installing the Hardware

1. Turn off the power of the PC and the crate in the case of a CompactPCI module.



CAUTION: For PCI modules the PC may have to be unplugged to ensure that the PCI bus has no power available. However, CompactPCI crates can be left plugged in since this ensures proper grounding.



CAUTION: Touch the antistatic package to a grounded object before removing the card from the package. Electrostatic discharge can damage the card.

2. Open the PC, identify a free PCI slot and carefully insert the DP Series card into it. Be sure to ground yourself by touching the grounded PC frame and avoid touching any components on the DP Series card. Make sure that the grounding of the card's mounting bracket to the back panel rail of the computer is done correctly. If present make sure that the fan's adjustable retainer is correctly positioned and tightened for mechanical support. Close the PC.
3. Turn on the power of the PC and start the operating system.



NOTE: Acqiris Digitizers are equipped with a LED. If this LED is not glowing orange or red when the power is applied there is a severe problem. Either the module is broken or the necessary voltages for its use are not available.



Devices that were installed using a previous version of **Acqiris Software** the instruments in these logical positions will still appear as *Unknown Devices*. This can be changed to the new Acqiris type category with the Grey Diamond icon by **Uninstalling** the device and then **Installing** again. Instructions on this procedure can be found in the **ReadMe.txt** file in the **manuals** folder of your Acqiris Software installation.

2.4. After Restarting

2.4.1. Windows 2000

Under Windows 2000, you *must* login with administrator privileges after the first boot following the hardware installation; the Plug&Play system must have the appropriate privileges to be able to complete your hardware installation successfully. After a successful hardware installation, you will be able to use your Acqiris Digitizer(s) with normal privileges.

At the first boot following the hardware installation, Windows will detect the new hardware and will install the devices automatically. The following image will appear.



 **NOTE:** In some systems an application program (such as AcqirisLive) will not yet work correctly at this point. One additional boot cycle may be needed if this is the first time that a hardware board is being installed.


2.4.2. Windows XP

Under Windows XP, you *must* login with administrator privileges after the first boot following the hardware installation; the Plug&Play system must have the appropriate privileges to be able to complete your hardware installation successfully. After a successful hardware installation, you will be able to use your Acqiris Digitizer(s) with normal privileges.


If you login with administrator privileges after the first boot following the hardware installation, Windows will detect the new hardware and start the “Found New Hardware Wizard” after a few seconds.



To proceed with the installation click "No, not this time" and then “Next”. The following window will appear:

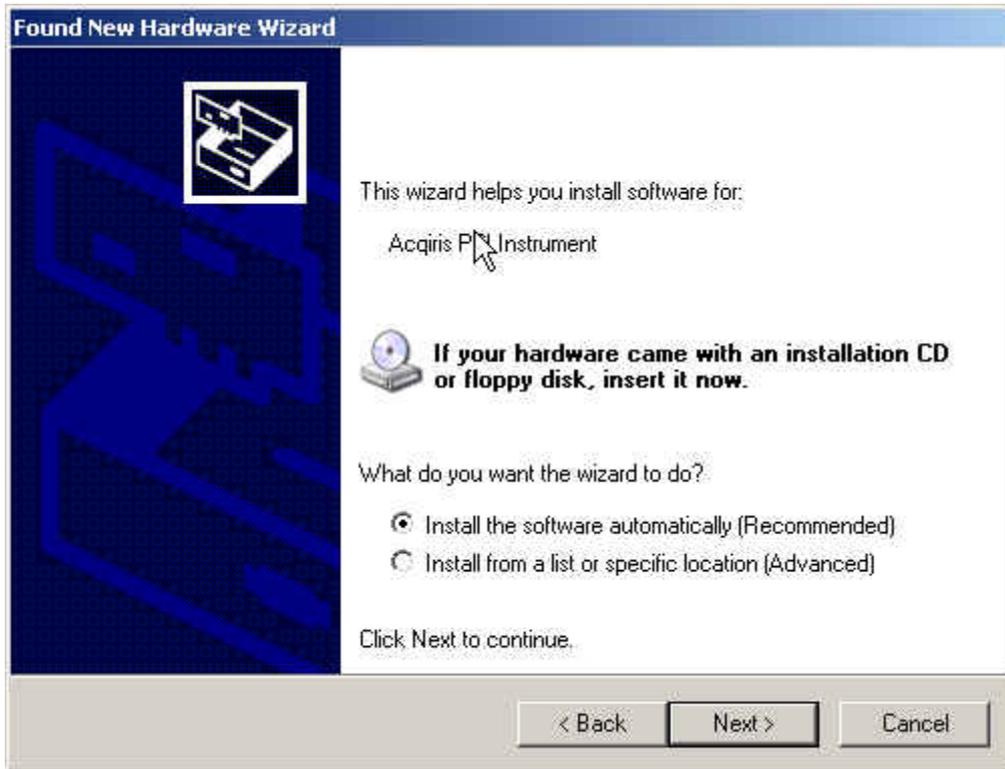

To proceed with the installation click “Next”. You should then see

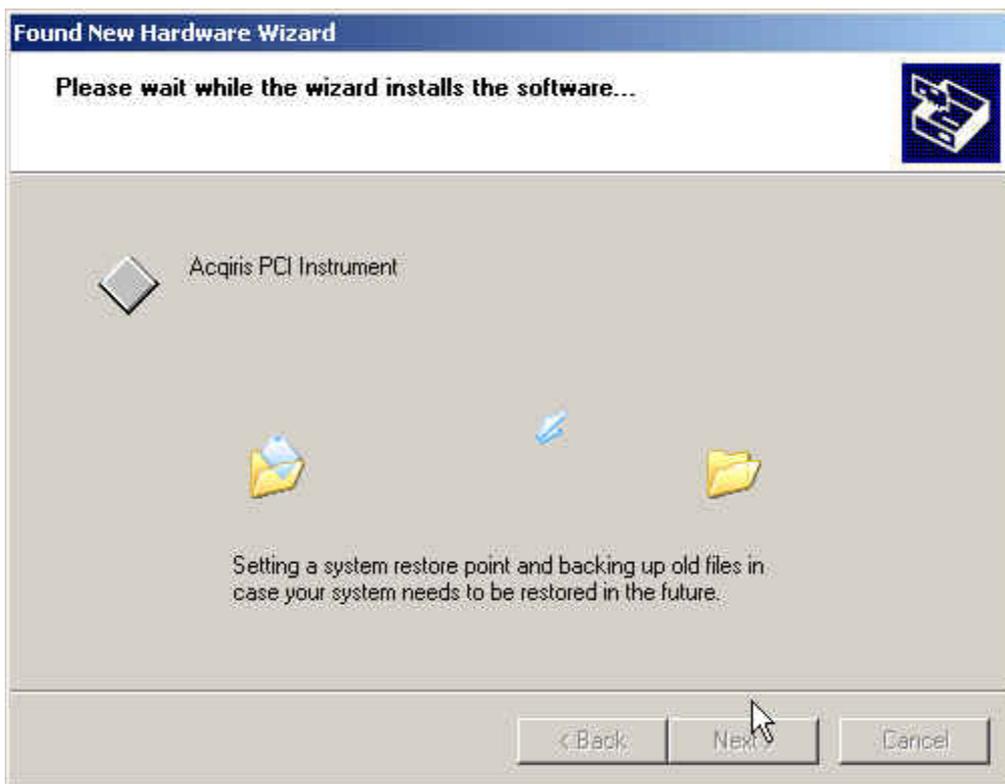

and then



 **NOTE:** In some systems an application program (such as AcqirisLive) will not yet work correctly at this point. One additional boot cycle may be needed if this is the first time that a hardware board is being installed.

2.5. LabVIEW RT

During program development you can choose whether you use LabVIEW or LabVIEW RT compatible libraries by switching the version present in **National Instruments\LabVIEW m.n\instr.lib**. This swap can be facilitated by using the **Install VI library for LabVIEW or LabVIEW Real-Time** shortcut available in the Shortcut folder under **Start** → **Programs**.

There is only one Acqiris Driver. It supports all Acqiris Instruments. The instructions below concern LabVIEW RT as used in NI PXI processors.

The Aq_RT.inf and AqRT_4.ini files must be uploaded to the target. To do this,

- start the MAX application,
- right click on the target
- select file transfer
- select the Aq_RT.inf file on your host machine and upload ('To Remote') to the LabVIEW RT working directory (/NI-RT/system) on the target
- select the Aq_RT.ini file on your host machine and upload ('To Remote') to the LabVIEW RT working directory (/NI-RT/system) giving it the name AqDrv4.ini

For Acqiris modules which need FPGA files you should,

- create the folder \firmware in the /NI-RT/system directory using the file transfer application
- select the FPGA files (from <AcqirisDxRoot>\Firmware) you want to copy to the target and upload them into the firmware directory

Restart the target after finishing the file transfers.

Restart the MAX Explorer and you should have Acqiris digitizers detected in your PXI system.

2.6. Installing the IVI-COM/C Driver

Please install Acqiris Software for Windows first. Then the *Acqiris Software CD* Autorun program gives access to two installers:

- IVI Shared Components 1.4

- IVI-COM/C IviAqD1 driver

These must be installed in the order shown above. For more information you can then consult the Readme.txt file in the IVI\Drivers\IviAqD1 folder or the documentation through the program shortcut present under Ivi/IviAqD1.

2.7. Distribution for Windows 2000/XP and Linux

The manuals/ReadMe.txt file contains a list of files to be found after a complete installation of Acqiris Software on Windows systems. Similarly the ReadMeLinux file gives the list of files corresponding to that installation.


3. Product Description

3.1. Overview


Acqiris Analyzers are designed to provide superior measurement precision and accuracy. Key acquisition specifications (such as DC accuracy, integral and differential non-linearity) have been optimized to deliver maximum measurement fidelity. Careful circuit layout, custom IC's and special packaging techniques have all been employed to reduce the overall random and correlated system noise. The use of custom IC's also dramatically reduces the total number of discrete components required. This has tremendous benefits on reliability and also allows the modules to use a minimal amount of power.

The AP series Analyzers are PCI modules that plug directly into any vacant PCI slot in a PC. They are fully programmable over the PCI bus and deliver 2 GS/s (AP240/AP201) and 1 GS/s (AP235/AP101) buffered acquisition performance. An acquired data sequence can be processed and transferred to a host processor over the PCI bus at speeds up to 100 Mbytes/s, *while simultaneously acquiring the next data sequence*. This explains the use of the term Sustained Sequential Recording (SSR). When appropriately configured, an Analyzer can reliably acquire all incoming triggers up to a rate considerably higher than that achievable by a conventional digitizer. This is particularly useful for the case of constant trigger rates in the kHz range.

For complete technical specifications concerning your particular analyzer please refer to the product's Specifications & Characteristics. In addition, we maintain up-to-date versions of all product data sheets on our web site (www.acqiris.com). The data sheets are available in pdf format and are best viewed using Adobe Acrobat software. If you have trouble accessing our web site, or viewing the data sheets, please contact your nearest sales office.

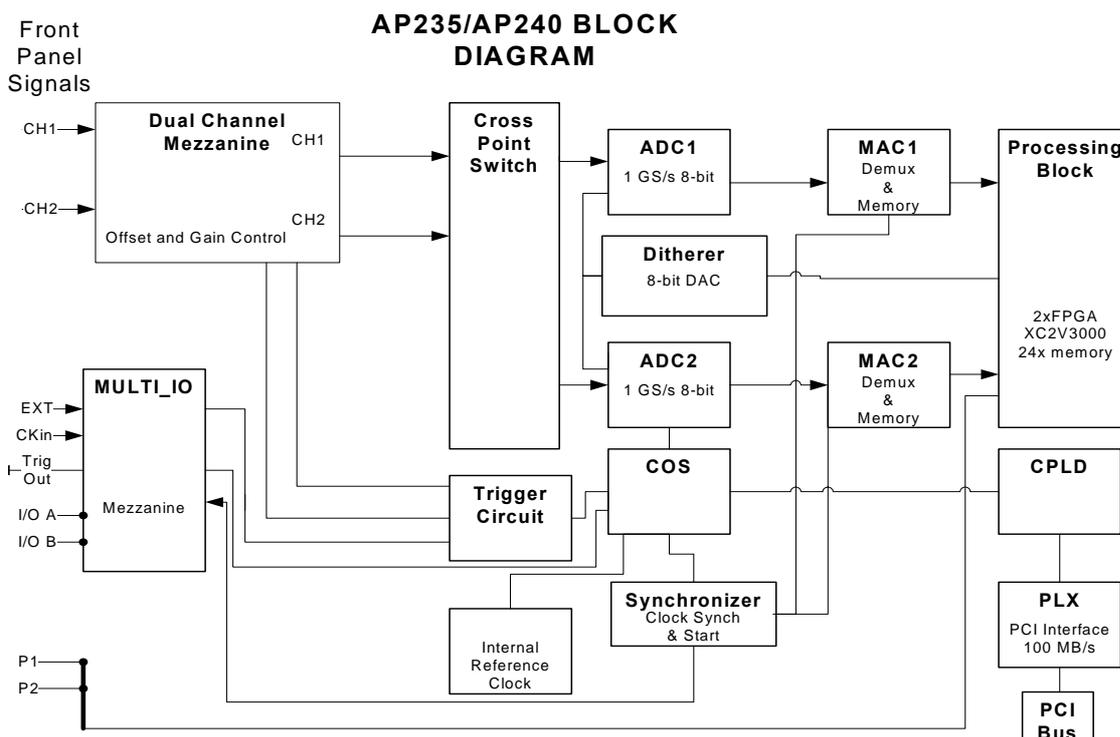

AP201 Block Diagram

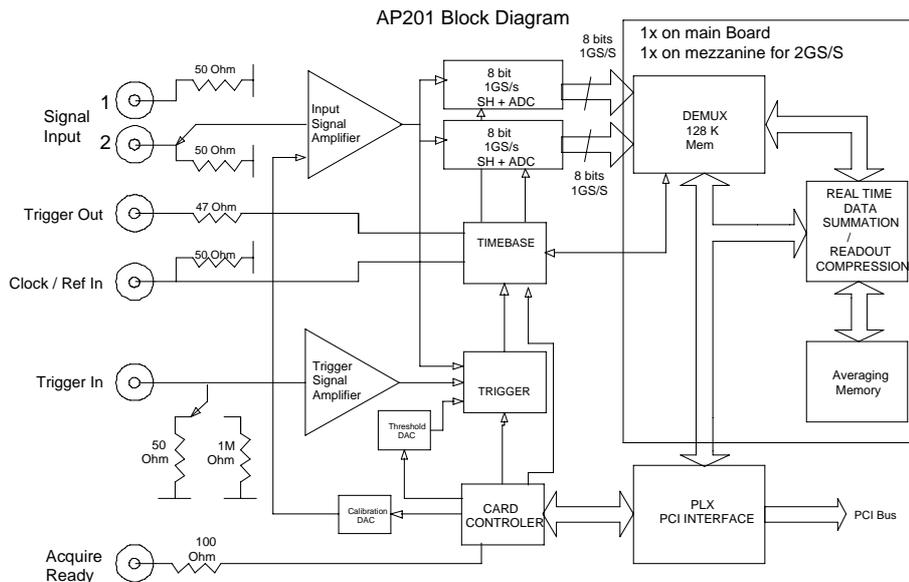

Key features

- High Speed and Precision** – The Acqiris Analyzers provide both high-speed acquisition capability and exceptionally high **simultaneous** data transfer rates to the host processor. The AP240 delivers real-time sampling rates of up to 2 GS/s with 1 GHz of analog bandwidth, while the AP201 delivers the same real-time sampling rates with 500 MHz of analog bandwidth, and the AP235 and AP101 offer sampling to 1 GS/s with 500 MHz of analog bandwidth.
- Signal conditioning** – The AP235/AP240 models include all of the required input signal conditioning such as gain, from 50 mV to 5 V Full Scale, offset adjustment, bandwidth limiter, and 50 Ω DC coupling. The models AP101/AP201 Analyzer offer similar functionality with gains ranging from 50 mV to 500 mV Full Scale. The inputs are fully protected, and calibration is a built-in feature. A crystal controlled time base and advanced trigger circuits are used to allow fast and stable acquisitions.
- Buffered Operation** – All of the Acqiris Analyzers offer a buffered acquisition mode. A sequence of waveforms can be acquired while simultaneously processing and reading out a previously acquired waveform sequence. Data reduction through selective readout of pre-defined *gates* in each waveform segment is also possible. Furthermore the AP235/AP240 models implement a threshold gate readout mode that allows data driven readout of interesting data only. AP235/AP240 models with the Peak^{TDC} option and the AP101/AP201 models support data processing to find peaks in the gates.
- Front panel control signals** – In order to ease the integration and synchronization of the Analyzer in a large variety of systems, a number of input and output signals are made available on the instruments front panel. Two MMCX connectors are dedicated to the Trigger Output and the Clock/Reference input, while two other MMCX connectors can be programmed for a variety of different usages, such as Trigger Ready, 10 MHz Reference Out or Sequence Stop. The AP235/AP240 offer the possibility of hardware reset of the segment timestamps through the P1 or P2 MMCX connectors.
- High Data Throughput** – All of the Acqiris Analyzers can easily be integrated into any standard computer with 1 free long PCI slot that provides the 3.3 V supply and adequate power required by the cards. Acquired data can be transferred, in DMA mode, directly to the processor over the PCI bus at rates of up to 100 Mbytes per second (4 ms for a 100 kpoints record, 10 ms for a 250 kpoints record, including DMA set up time overhead).

3.2. Channel Input

The principle characteristics of the input channels are given in the table below:

Model Agilent #	Bandwidth into 50 Ω	Minimum Rise Time	BW Limiter selections	Maximum Full Scale	Maximum Offset	Other particularities
AP240 U1082A	1 GHz	0.35 ns	20, 200, 700 MHz	5 V	5 V	2 channels
AP235 U1082A	0.5 GHz	0.7 ns	20, 200 MHz	5 V	5 V	2 channels
AP101 AP201 U1081A	0.5 GHz for > 50 mV FS 200 MHz @ 50 mV FS	0.7 ns 1.8 ns	25 MHz	0.5 V	2 V	1 channel with 2 selectable inputs

3.2.1. Coupling & Impedance

Two 50 Ω signal input BNC connectors are provided. The AP240 and AP235 are two channel modules. In the case of the AP101 and AP201 these two connectors are multiplexed to the module's input. The active input may be selected by software to allow connection of one of the two possible signal sources. The 50 Ω coupling offers high quality termination with better than $\pm 1\%$ precision. It is ideally suited for use with 50 Ω transmission lines.

All modules offer DC coupling. The AP240 and AP235 also provide AC coupling with a 32 Hz cutoff frequency.

3.2.2. Input Protection

The input amplifiers are fully protected against over-voltage signals. The device can tolerate an input current of up to 40A for 1 μ s without damage to the front-end electronics. For extreme cases, such as detector sparking, an external coaxial spark gap and attenuator can be supplied on request (See **Appendix A: XA100 BNC Input Overvoltage Protection**).

3.2.3. Mezzanine Front-end

The front-end electronics are all mounted on a removable mezzanine card. In the event of accidental damage, or as components fatigue over time (e.g. relays in high duty cycle automated testing applications), the mezzanine card allows for fast and efficient replacement.

3.2.4. Bandwidth and Rise Time

The bandwidth specification indicates the frequency at which an input signal will be attenuated by 3 dB (approximately 30% loss of amplitude). The bandwidth also affects the minimum rise and fall times that can be passed through the front-end electronics. A pulse with a very sharp edge will be observed to have a minimum rise time (τ_{\min}) determined by the front-end electronics. In general a pulse with a given 10-90% rise time ($\tau_{10-90\text{real}}$) will be observed with a slower value given by:

$$\tau_{10-90}^2 = \tau_{10-90\text{real}}^2 + \tau_{\min}^2 \quad \text{where} \quad \tau_{\min} (\text{ns}) \approx 0.35 (\text{GHz}\cdot\text{ns}) / \text{BW} (\text{GHz})$$

If desired, hardware bandwidth limiters can be selected.

3.2.5. Input Voltage and Offset

The input channel provides a fully programmable amplifier with variable input voltage and offset. Full Scale (FS) input voltages are selectable from 50 mV to the maximum indicated in the table above in a 1, 2, 5 sequence. Care should be taken to select an input voltage range that will allow the signal to be recorded using as much dynamic range of the digitizer as possible. Signals going outside of the FS Range will be clipped and data values for the clipped portion of a signal should be regarded as erroneous. The Variable Offset is programmable in the range of ± 2 V when using an FS Input Voltage setting of 500 mV or below, increasing to ± 5 V for FS settings above 500 mV. The raw 8 bit ADC data values are in the range [-128,+127] with the first and last values reserved for underflow and overflow respectively. The midpoint value, 0, of the range corresponds to the negative of the offset voltage. Thus the Full Scale Range (FSR) goes from

–Offset Voltage – (FS/2) to –Offset Voltage + (FS/2)

Signals going outside of the FSR will be clipped and data values for the clipped portion of a signal should be regarded as erroneous.

3.2.6. Vertical Resolution

The AP Series Analyzers use an ADC system with 8 bits of vertical resolution (256 levels). The dynamic range of the ADC covers the Full Scale (FS) of the Input Voltage setting. For example, if the Input Voltage is set to 500 mV then the ADC resolution is equivalent to 1.95 mV. To obtain the best dynamic range from the ADC care should be taken to ensure that the input signal varies over more than 50% of the Input Voltage Full Scale (FS) setting. The highest and lowest levels of the ADC correspond to underflow and overflow conditions.

3.2.7. DC Accuracy

The AP Series Analyzers use low noise front-end electronics in order to ensure voltage measurement is made with accuracy and precision. DC voltage accuracy is better than $\pm 2\%$ ($\pm 1\%$ typical) of the input voltage full-scale. The differential linearity is better than ± 0.9 LSB for the AP235/AP240 and ± 0.7 LSB for the AP101/AP201.

3.3. Trigger

3.3.1. Trigger Source

The trigger source can be a signal applied to either of the Input Channels (for internal triggering) or the External Trigger Input.

The modules provide a front panel BNC External Trigger Input. The External Input provides a fully functional trigger circuit with selectable coupling, level and slope. The AP240 and AP235 have 50 Ω termination impedance with diode protection against overload. A ± 5 V limit on trigger signals should be respected, although somewhat higher voltages for short time periods will not damage the unit.

The AP101 and AP201 have both 1 M Ω and 50 Ω termination impedance. Overload protection will automatically switch the coupling from 50 Ω to 1 M Ω if the signal is greater than ± 5 V DC. If 50 Ω termination is selected a ± 5 V limit on trigger signals should be respected, although somewhat higher voltages for short time periods will not damage the unit. For 1 M Ω input, signals up to ± 100 V (DC + peak AC < 10 KHz) are allowed.

3.3.2. Trigger Coupling

Trigger coupling is used to select the coupling mode applied to the input of the trigger circuitry. The AC LF Reject mode couples signals capacitively and removes the input signal's DC component and frequencies below 50 Hz for the AP240 and AP235 models (and 50 KHz for the AP101 and AP201). DC mode allows all signal components to be passed through to the trigger circuit. The AP240 and AP235 models have an HF Reject mode that removes signal components above 50 KHz. They also implement an HF trigger that allows triggers to be reliably accepted at rates above ~ 1 GHz. In this mode, triggers occur on every fourth positive edge. In the HF mode negative slope and window triggers are not available.

3.3.3. Trigger Level

The trigger level specifies the voltage at which the selected trigger source will produce a valid trigger. The trigger level is defined as a set voltage. Using the internal trigger with DC coupling, the level is set with respect to the midpoint voltage ($V_m = -$ Offset voltage) of the digitizer's vertical scale. All trigger circuits have sensitivity levels that must be exceeded in order for reliable triggering to occur.

The AP240 and AP235 offer level control for all trigger coupling modes. Internal trigger level settings (expressed in %) must be within $V_m \pm 0.5$ FS, where FS is the channel Full Scale. In addition, they implement a Window trigger. Two trigger level thresholds are used to define the desired range. The trigger can then be chosen to occur either when the signal exits or enters the window range. This mode can be thought of as the appropriate OR of two edge triggers of opposite slope.

For the AP101 and AP201 the AC coupled mode is implemented with an auto-level trigger. Internal trigger level settings for DC coupling must be within $V_m \pm 0.6$ FS, where FS is the channel Full Scale.

The AP240 and AP235 models allow the user to choose the external trigger Full Scale from the set of values 0.5, 1.0, 2.0 or 5.0 V. The external trigger level can then be set to values in the range ± 0.5 FS. The AP101 and AP201 models have an external trigger range of ± 3 V.

The AP series Analyzers will trigger on signals with a peak-peak amplitude > 15% FS from DC to their bandwidth limit.

3.3.4. Trigger Slope

The trigger slope defines the direction of the signal that will be used to initiate the acquisition when it passes through the specified trigger level. Positive slope indicates that the signal is transitioning from a lower voltage to a higher voltage. Negative slope indicates the signal is transitioning from a higher voltage to a lower voltage.


3.3.5. External Trigger Output

When the module is ready to be triggered and a valid trigger signal occurs, a trigger output is generated for external use. It is always available on the Front Panel Trigger Out MMCX connector.



NOTE: The External Trigger Output functionality is implemented in the hardware. No Trigger Out signal occurs for software-generated triggers such as those of the AUTO mode of APx01Demo or through the use of the function `AcqrsD1_forceTrigger`.

Trigger Output Block diagram:

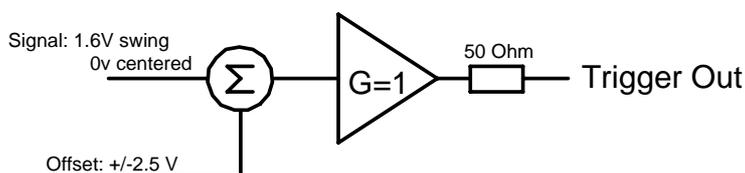


The output swing is 1.6 V (± 0.8 V) when unloaded and 0.8 V when terminated on 50 Ω . The rise and fall times are 2.5 ns typical. The offset can be adjusted, by software control in the range [-2.5 V, +2.5 V] unloaded, or [-1.25 V, +1.25 V] into 50 Ω . The maximum output current capability is ± 15 mA. As the output is retro-terminated, it is possible to drive a 50 Ω line unterminated (HiZ) without loss of performance.

For a TTL compatible signal, set the offset to 1.0 V and the swing at destination will be +0.2 to +1.8 V.

For an ECL compatible signal, terminated on 50 Ω to -1.2 V, set the offset to -1.2 V and the output will be in the range [-0.8 V, -1.6 V]).

Alternatively, to reduce the current drawn from the digitizer, the terminations below can be used:

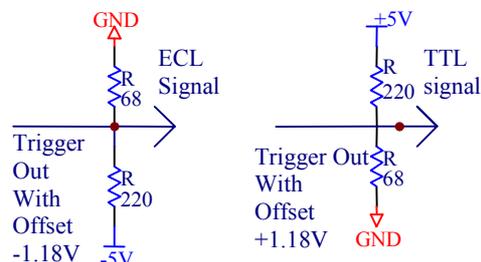

3.3.6. Trigger Status

The front panel includes a tri-color LED indicator to show the status of the trigger. When the LED is green it indicates the trigger is armed and waiting for a valid trigger to occur. Red indicates that the trigger has occurred, the acquisition is complete and the data is waiting to be readout. The user can override the default functions and program the LED color in an application-specific manner.

3.3.7. Trigger Veto and Timeout (AP101/AP201 ONLY)

Certain applications can generate a *Prepare for Trigger* signal that can be generated just before a window in time for an acceptable trigger. An Analyzer in the buffered data acquisition mode (see `AcqrsD1_configMode`) can be configured to recognize the *Prepare for Trigger* signal on the MMCX IO A or B connector (see `AcqrsD1_configControlIO`) and wait a software selectable time before accepting a trigger. This veto time can be set in the range [0 s, 1 s] in steps of 30 ns. The Trigger Veto feature is useful in applications where an early echo (response signal) should be ignored.

Furthermore, if desired, the Analyzer can be asked to generate an artificial trigger if no real trigger occurs after the window has been open for a time that can also be set in that same range. This is a way to ensure that there will be data associated with every *Prepare for Trigger* signal even if no real trigger occurs. A *Prepare for Trigger* Signal can only be accepted after the ACQRDY has been generated. The data values for an artificial trigger will all be -128 = 0x80 and the time taken to generate this special acquisition will be the same as for a real trigger. The Timeout

feature is useful to ensure the correct pairing of stimulus and response and to avoid unwanted “hung” states where the Analyzer is waiting for a trigger that will never occur.

3.4. Data Acquisition – Common

The AP series Analyzers can be used either for buffered acquisition with analysis or normal digitizer data acquisition. These two modes are quite different in their functioning. After describing some common characteristics they will be treated individually in the sections below. The digitizer mode is particularly useful when trying to configure an AP module for later use as an analyzer.

3.4.1. Timing

A crystal controlled time base is used to generate the clock of the digitizer. Clock accuracy is better than ± 2 ppm, and individual sample points are recorded with better than ± 10 ps (< 2 ps rms) of sampling jitter.

3.4.2. Sampling Rate

The AP series Analyzers include an analog-to-digital converter (ADC) that can sample waveforms in real time, at rates from 1 or 2 GS/s (0.5 ns per point) down to 100 S/s (10 ms per point). The sampling rate can be programmed and is selectable in a 1, 2, 2.5, 4, 5 sequence (i.e. 100 MS/s, 200 MS/s, 250 MS/s, 400 MS/s, 500 MS/s, 1GS/s, 2GS/s). The data of all of the active channels is acquired synchronously; all of the ADC's are acquiring data at the same time, to within a small fraction of the maximum sampling rate.

3.5. Data Acquisition - Digitizer Mode

3.5.1. Acquisition Memory and Time Base Range - Digitizer Mode

Data from the ADC are stored in on-board acquisition memory. The amount of memory in use for acquisition can be programmed and is selectable from 2 points to 261350 points, the full amount of acquisition memory available. Note that the analyzer mode memory is much greater.

For technical reasons, a certain memory “overhead” is required for each waveform, reducing the available memory by a small amount. In order to simplify programming, an interface function recommends the best sampling rate and the maximum possible number of data points, taking into account the available memory, the requested time window, the number of segments (in Sequence mode), as well as the required memory overhead.

The Time Base Range defines the time period over which data is being acquired. For example, the AP201 in the digitizer mode has an acquisition memory of just under 256 Kpoints and maximum sampling rate of 2 GS/s. Thus, at the maximum sampling rate, the digitizer can record a signal over a time period of up to 130 μ s (256 Kpoints * 0.5 ns/point). The time base range can be adjusted by varying the amount of acquisition memory or the sampling rate of the digitizer.

3.5.2. Pre- and Post-Trigger Delay - Digitizer Mode

To increase trigger flexibility a pre- or post-trigger delay can be applied to the trigger position. However, pre-trigger is not available in the analyzer mode.

The amount of pre-trigger delay can be adjusted between 0 and 100% of the acquisition time window (i.e. sampling interval x number of samples), whereas the post-trigger delay can be adjusted between 0 and 200 million samples.

Pre- or post-trigger delays are just different aspects of the same trigger positioning parameter:

- The condition of 100% pre-trigger indicates that all data points are acquired prior to the trigger, i.e. the trigger point is at the **end** of the acquired waveform.
- The condition of 0% pre-trigger (which is identical to a post-trigger of 0) indicates that all data points are acquired immediately after the trigger, i.e. the trigger point is at the **beginning** of the acquired waveform.
- The condition of a non-zero post-trigger delay indicates that the data points are acquired after the trigger occurs, at a time that corresponds to the post-trigger delay, i.e. the trigger point is **before** the acquired waveform.

The digitizer hardware accepts pre- and post-trigger adjustments in increments of 16 samples. By definition post-trigger settings are a positive number and pre-trigger settings are a negative number.

Thus it is only natural that the software drivers provided treat pre- and post-trigger delays as a single parameter in seconds that can vary between $-nbrSamples * samplingInterval$ (100% pre-trigger) and $+maxPostTrigSamples * samplingInterval$ (max post-trigger). Since the Acqiris software drivers provide very accurate trigger position

information upon waveform readout, the accepted resolution of the user-requested pre-/post-trigger delay is much better than 16 samples. For more details, refer to the **Programmer's Reference Manual**.

3.5.3. Single and Sequence Acquisitions - Digitizer Mode

Digitizers acquire waveforms in association with triggers. Each waveform is made of a series of measured voltage values (sample points) that are made by the ADC at a uniform clock rate. To maximize sampling rates and utilize memory as efficiently as possible, the digitizers include both Single and Sequential storage modes.


The Single Acquisition mode is the normal operation of most digitizer products. In this mode an acquisition consists of a waveform recorded with a single trigger. The user selects the sampling rate and acquisition memory size and sets the number of segments to 1 (default value).

The modules also include a built-in Trigger Time Interpolator (TTI) that measures the time from the trigger point to the first sample point. This information is essential for determining the precise relation between the trigger or other event of interest and the digitized samples of the signal. The TTI resolution is 80 ps.

The Sequence Acquisition mode allows the capture and storage of consecutive “single” waveforms. Sequence Acquisition mode is useful as it can optimize the digitizer’s sampling rate and memory requirements for applications where only portions of the signal being analyzed are important. The mode is extremely useful in almost all impulse-response type applications (RADAR, SONAR, LIDAR, Time-of-Flight, Ultrasonics, Medical and Biomedical Research, etc.).

In Sequence Acquisition mode the acquisition memory is divided into a pre-selected number of segments. Waveforms are stored in successive memory segments as they arrive. Each waveform requires its own individual trigger. The memory can be divided into any number of segments between 2 and 200. In Sequence Acquisition mode the user needs to specify the sampling rate, the total acquisition memory and number of segments. Note that the Single Acquisition mode is just a special case of the Sequence Acquisition mode with the number of segments set to 1.

.Sequence acquisition enables successive events, which can occur within a very short time, to be captured and stored without loss. A crucial feature of Sequence Acquisition mode is that it has a very fast trigger rearm time. A fast trigger rearm helps produce very low “dead time” (less than 800 ns for the highest available sampling rates) between the segments of a sequence acquisition. The “dead time” is the period after the end of an event when the card cannot digitize data for a new trigger event. To complement this mode of operation the digitizer can also measure and store the arrival time of each trigger using the information from the on board TTI (Trigger Time Interpolator). Readout of

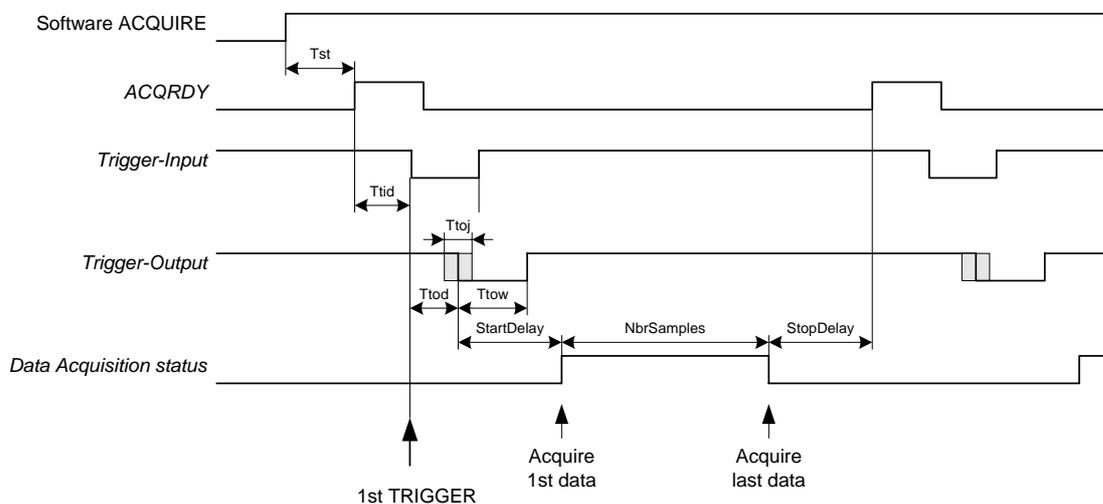

the individual trigger time stamps makes it possible to determine the time from one trigger to any other trigger in the sequential acquisition. The TTI resolution sets the resolution of the trigger time stamps.

3.6. Data Acquisition – Analyzer mode

Before the Analyzer starts digitizing signals it must be armed and in a ready state. Arming the module is done by software (ACQUIRE command). If the Analyzer is in *buffered mode* (also called *dual-memory mode*), it will wait for a valid trigger before recording any data (see the timing diagram below).

When the module is armed (via the **AcqrsD1_acquire** command) and ready to acquire data, it generates, after a time T_{st} (which may vary between 50 μ s and 10 ms, depending on the processor and the operating system used), a Ready signal (ACQRDY). This TTL signal is available on the front panel I/O A MMCX connector. For the AP101/AP201 modules it can also be found on the front panel multi-pin connector. The Ready signal is normally used to activate an external event, such as firing a laser.

This external event in turn can produce a signal that can be used to generate the trigger for the module. The trigger-input pulse can have a delay *Ttid* that can vary from zero to as long as required. The Ready signal (ACQRDY) is reset with the arrival of the trigger pulse. When the acquisition is running, or when the module has not been enabled to start the acquisition, ACQRDY is set low.

When the trigger input pulse is received, the Analyzer generates a trigger output which internally starts the acquisition after a predetermined delay shown as *StartDelay* in the diagram above. The trigger output is provided for applications where an external signal source must be synchronized with the internal sampling clock. The Analyzer's trigger output edge is synchronized to the sample clock to within <100 ps. The trigger output accurately preserves the sample timing with respect to the event of interest from shot-to-shot. Ideally, the trigger output is used to synchronously activate an event serving as the zero-time reference for the time-dependent phenomena that are measured in each accumulated waveform. As an example, the trigger output could activate the application of an electric or magnetic field to a sample under test. This method ensures that the application of the field is identically synchronized to the sample clock for each shot contained within the average. It also means that each sample in the accumulated sum is identically positioned shot-to-shot with respect to a fixed zero-time reference, provided that the time between the trigger output and the event waveform is kept constant.

The trigger output has, with respect to the trigger input, a delay of around $Ttod \pm Ttoj/2$, where *Ttod* is about 10 ns and *Ttoj* is 1 clock pulse. The digitizing clock is started at a precise, user configurable time after the trigger output pulse. This is shown as *StartDelay* in the diagram above. The allowed memory lengths, start delay values, and increments are shown in the table below.

Model - Mode	Channels	Maximum Sampling Rate	Maximum length for Acquisition Default Memory	Maximum length for Acquisition Optional Memory	Memory Minimum/ Increment	Start Delay Maximum
AP240 – Single-channel	1	2 GS/s	4M points	16M points	32	16M points
AP240 – Dual-channel	2	1 GS/s	2M	8M	16	16M
AP235 – Single-channel	1	1 GS/s	4M	16M	32	16M
AP235 – Dual-channel	2	0.5 GS/s	2M	8M	16	16M
AP201	1	2 GS/s	4M	NA	32	32M
AP101	1	1 GS/s	2M	NA	16	16M

The minimum *StartDelay*, that achieved when the user configurable time is set to 0, depends on the sampling rate as shown in the table below.

StartDelay /Sampling Rate	AP240 Single-channel	AP235 Single-channel	AP240 Dual-channel	AP235 Dual-channel	AP201	AP101
2GS/s	18.8 ns	NA	NA	NA	18.8 ns	NA
1	22.6	22.6 ns	18.8 ns	NA	22.6	18.8 ns
0.5	38.2	38.2	22.6	22.6	38.2	22.6
0.4			25.8	25.8		25.8
Zero points						
	10	10	5	5	10	5

It should also be noted that when the user configurable *StartDelay* is 0, the first few data points of data, called Zero points in the table above, are automatically set to -128 (bottom of the range) and should be ignored.

The module can be used to limit the system's overall event rate. A rate adjustment delay *StopDelay* can be programmed in steps from zero to the model maximum with the same range and resolution as *StartDelay*. For maximum acquisition rates *StopDelay* should be set to zero.

3.6.1. Random and Synchronized Triggers

The Analyzers have two trigger synchronization modes, **Non-Resync** or **Resync**. Normally the user's trigger input is not synchronized with respect to the internal sampling clock.

In **Non-Resync** mode, the internal trigger signal, which generates the trigger output and starts the digitizing clock after *StartDelay*, is not resynchronized with respect to the internal sampling clock. The trigger output is synchronized with respect to the sampling clock to better than ± 100 ps. The jitter, i.e. the variation of the position of the clock with respect to the input trigger time, as a function of sampling rate is shown in the table below. In order to preserve the timing accuracy of successively summed acquisitions, the AP101/AP201 also offer a **Resync** mode, where the jitter is limited to ± 200 ps (at the highest sampling rate).

Sampling Clock Synchronization with respect to Trigger Input				
Sample Rate	Resynchronized Mode		Non-Resynchronized Mode	
	AP201 or AP240/AP235 Single-channel	AP101 or AP240/AP235 Dual-channel	AP201 or AP240/AP235 Single-channel	AP101 or AP240/AP235 Dual-channel
2 GS/s	± 200 ps	--	± 500 ps	--
1 GS/s	± 325 ps	± 200 ps	± 1 ns	± 500 ps
500 MS/s	± 325 ps	± 325 ps	± 1 ns	± 1 ns
400 MS/s	± 400 ps	± 400 ps	± 1.25 ns	± 1.25 ns
≤ 250 MS/s	± 400 ps	± 400 ps	± 1.25 ns	± 1.25 ns

The **Resync** mode is the preferred method for applications where the acquisition jitter (and the trigger out jitter) must be minimized to avoid the reproducibility problems possibly associated with the higher jitter found when using the **Non-Resync** mode.

3.6.2. Time Base Range

The Time Base Range defines the time period over which data is being acquired. For example, the AP240 can have an acquisition memory of 16 Mpoints and maximum sampling rate of 2 GS/s. Therefore, at the maximum sampling rate, the Analyzer can record a signal over a time period of up to 8 ms (16 Mpoints * 0.5 ns/point). The time base range can be adjusted by varying the amount of acquisition memory or the sampling rate of the Analyzer.

3.6.3. Dual-Bank Memory & Buffered Sequence Acquisition

Digitizers acquire waveforms in association with triggers. Each waveform is made of a series of measured voltage values (sample points) that are made by the ADC at a uniform clock rate. To maximize sampling rates and utilize memory as efficiently as possible, the digitizers include both Single and Sequential storage modes.

The Buffered Sequence Acquisition mode allows the capture and storage of consecutive "single" waveforms. Sequence Acquisition mode is useful as it can optimize the digitizer's sampling rate and memory requirements for applications where only portions of the signal being analyzed are important. The mode is extremely useful in almost all impulse-response type applications (RADAR, SONAR, LIDAR, Time-of-Flight, Ultrasonics, Medical and Biomedical Research, etc.). Buffered Sequence acquisition enables successive events, which can occur within a very short time, to be captured and stored without loss. A crucial feature of Sequence Acquisition mode is that it has a very fast trigger rearm time. A fast trigger rearm helps produce very low "dead time" (less than 800 ns for the highest available sampling rates) between the segments of a sequence acquisition. The "dead time" is the period after the end of an event when the card cannot digitize data for a new trigger event.

In the *buffered* mode, data from the ADC are stored in the on-board 'dual-bank' memory. A typical acquisition consists of a *sequence* of waveforms of equal size. The actual waveform size and the number of waveforms in an acquisition sequence can be programmed as shown earlier in the table. The number of waveforms in a sequence is adjustable from 1 to 8191 for the AP235/AP240 and 16384 for the APx01. Of course, the number of waveforms *times* their length may not exceed the actually available memory. The driver software protects against invalid requests, by adjusting the waveform length appropriately. Note that the Single Acquisition mode is just a special case of the Sequence Acquisition mode with the number of segments set to 1.

Note that when in *digitizer* mode, the Analyzer stores the ADC data into a different (non-buffered) memory. It is then restricted to 128 Kpoints (AP235/AP240 Dual channel mode and AP101) or 256 Kpoints (AP235/AP240 Single channel mode and AP201) and 200 segments.

3.6.4. Buffered Acquisition Mode with Data Readout

In *buffered* mode, multiple waveforms are acquired and stored, with minimal dead time between them. All waveforms in such a sequence are deposited into one of two possible memory banks. The user program can be either written to explicitly control the processing, start of next acquisition, and readout of data, or to automatically start the processing and the next acquisition as soon as the previous acquisition has completed. This latter case will be called the *autoswitch* mode while the former can be called the *explicit* mode. At start-up of the *autoswitch* mode the program will start the first acquisition and immediately give the command to allow the second acquisition to start as soon as possible.

After the pre-programmed number of waveforms is acquired, the Analyzer terminates the acquisition cycle. At this point in the *autoswitch* mode, if the program has already indicated that the acquisition/readout process is to be continued, the next acquisition, as well as the processing of the newly acquired data will be started. The user program can then initialize readout when the processing has completed and, if so desired, issue the command to allow the subsequent acquisition to be started as soon as possible. In the *explicit* mode after the program learns that the acquisition is terminated, typically by interrupt, it responds by switching to the other memory bank, restarting a new acquisition sequence and then reading out the previously acquired waveform sequence. The readout occurs *while* a new acquisition sequence is in progress.

The readout process may be restricted to pre-defined sections of the acquired waveforms (defined by gates). Alternatively, in the AP101 and AP201 the data of interest may be processed to give peak information before readout.

3.7. AP235/AP240 Data processing

Processing is either

- done automatically by the Analyzer after each acquisition is completed and while the next acquisition is under way for the case of AP235/AP240 with the Peak^{TDC} option using histogramming
- done under control of the host computer after the acquisition is complete and before the readout is started.

In the simplest case this processing does nothing to the data and terminates immediately. If real processing is needed, the throughput is 260 MS/s, ignoring start-up effects associated with each gate and overhead for each peak found. The actual amount of data to be processed can be reduced significantly by intelligent use of the gate definition and of the processing parameters for the peak detection.

3.7.1. User Gates

Gates can be used to define regions of interest within an acquisition segment. This can be used to limit the amount of data to be transferred or as a starting point for peak detection processing. A gate is defined as a range of data defined by a first sample and the number of samples in the gate. Both the first sample and the number of samples must be multiples of 4 (except 8 for an AP201 or an AP235/AP240 in single channel mode). An arbitrary number of gates, up to 4095, can be defined. They are used for each segment in the acquisition.

3.7.2. Threshold gates

The data itself can also be used to define the regions of interest. The threshold gate mechanism selects regions of data to be transferred if the data value exceeds a user-defined threshold. This selection algorithm can also be used for data less than the user-defined threshold if the data inversion option is chosen. The number of samples selected as interesting will always be a multiple of 4 and will start and end on a block of 4 samples defined from the acquisition start. The interesting region will start with the first block of 4 which has a data value above threshold and end after the block of four in which the data drops down below threshold. In dual-channel acquisition, if the below threshold uninteresting region between two selected points does not contain a data point interval $[16*n+1, 16*n+32]$ (for some value of n) the data involved will be combined into a common gate and all of the points in between will be selected. In the case of single-channel acquisition this interval becomes $[32*n+1, 32*n+64]$. In other words, for the dual-channel case, if the interval between two points above threshold does not contain an aligned (to a multiple of 16) block of 32 below threshold points the data will be combined into a single gate. Given that a region is selected for transfer the user can select an additional number of samples (0, 4, 8, 12, or 16) both before and after it for the final transfer. Alternatively, the number of samples before and the total number of samples can be chosen. The subtleties of these definitions imply that all gates should be examined carefully for multiple pulses and that data of interest may be lost if the total number of samples demanded is less than those in a combined gate region. There is a limit of 8192 on the number of gates that can be transferred for a single acquisition on a single channel. Thus any threshold gate waveform with 8192 gates may be incomplete. To help avoid this overflow condition a limit can be placed on the number of threshold gates allowed for each segment. If this limit is active with the value m , then the first $m-1$ gates and the last gate of each segment will be transferred.

3.7.3. Time stamp information

The 10 MHz reference clock is used to increment a register that is recorded as a time stamp for each trigger accepted by the Analyzer in gated mode. These time stamps can be used to fix the absolute time of a trigger. If desired, the time stamp register can be reset with a hardware signal on the MMCX DPU P1 or P2 connectors. It can also be read by the controlling computer at any time when data is available to be read..

3.7.4. Peak detection

The data selected for treatment by either of the gate mechanisms described above can be further processed to identify the position and amplitude of multiple peaks in each region. The peaks are defined using two thresholds to provide some Hysteresis. A peak signal must first increase by a given threshold above its starting level and then decrease below its maximum value by the second threshold value. The processing algorithm for positive peak search is as follows:

- P1. Starting at the first data point of the gate ($V[0]$) skip through the data until the data point voltage exceeds the minimum value observed by *StartDeltaPosPeakV*. The search terminates at the end of the gate. Let i_s be the data point index satisfying this test, then

$$V[i_s] \geq \min(V[0] \dots V[i_s]) + \textit{StartDeltaPosPeakV}.$$

- P2. Continue scanning through the data until the data point voltage is less than the maximum value observed less *ValidDeltaPosPeakV*. The search terminates at the end of the gate. Let i_v be the data point index satisfying this test and i_m be the data point index of the first maximum, then

$$V[i_m] = \max(V[i_s] \dots V[i_v]) \quad \text{and} \quad V[i_v] \leq V[i_m] - \textit{ValidDeltaPosPeakV}.$$

- P3. Record the value and position of the accepted maximum data point (if any). There will be a result if condition 1 above is satisfied.
- P4. Go back to step P1 but start at the point $V[i_v+1]$.

The same algorithm can be used for negative peak search if data inversion has been enabled. This is done with the obvious modifications, namely:

- N1. Starting at the first data point of the gate skip through the data until the data point voltage falls below the minimum value observed by *StartDeltaNegPeakV*. The search terminates at the end of the gate. This means that

$$V[i_s] \leq \max(V[0] \dots V[i_s]) - \textit{StartDeltaNegPeakV}.$$

- N2. Continue scanning through the data until the data point voltage is greater than the minimum value observed plus *ValidDeltaPosPeakV*. The search terminates at the end of the gate. This means that

$$V[i_m] = \min(V[i_s] \dots V[i_v]) \quad \text{and} \quad V[i_v] \geq V[i_m] + \textit{ValidDeltaNegPeakV}.$$

- N3. Record the value and position of the accepted minimum data point (if any). There will be a result if condition 4 above is satisfied.
- N4. Go back to step N1 but start at the point $V[i_v+1]$.

The peak processing is done after the acquisition is completed. All of the gated data has to be examined at a rate of 260 MB/s. Furthermore there is an overhead of around 300 ns/acquisition and 250 ns/peak. In some cases this can reduce the maximum useable trigger rate for the histogram mode described in the next section.

3.7.5. Interpolation of Peak Position and Amplitude

If desired the simple peak values described above can be refined using interpolation. The time and amplitude resolutions can each be increased by a factor of 16. The interpolation is done with a parabolic spline calculated using the three data points around the peak.

3.7.6. Histogramming Peak Data

The next logical step in the Peak^{TDC} option is the capability of histogramming the peak data. For each peak found the histogram can be incremented either by 1 or by the ADC value of the peak. This can be done for a user-defined number of acquisitions. If segmented acquisitions are being used the user can choose between one histogram for all of the segments or a histogram for each segment. Since the buffered acquisition mode is being used the minimum time between triggers is about 1 μ s greater than the acquisition time for the data.

If interpolated peaks are being calculated the histogram's horizontal bin width can be selected to be a fraction of a sample interval. The allowed choices are $\frac{1}{2}$, $\frac{1}{4}$, $\frac{1}{8}$, and $\frac{1}{16}$ of the nominal bin width. Similarly the granularity of the contents of each bin can be chosen to be $\frac{1}{2}$, $\frac{1}{4}$, $\frac{1}{8}$, or $\frac{1}{16}$ of an ADC LSB.

A histogram can be accumulated for a desired number of acquisitions. After this the acquisition process will be stopped until the program explicitly starts it again. Histogram readout can only be done when the acquisition is inactive.

3.8. AP101/AP201 Data processing

Processing is done under control of the host computer after the acquisition is complete and before the readout is started. In the simplest case this processing does nothing to the data and terminates immediately. If real processing is needed, the throughput is 130 MS/s, ignoring start-up effects associated with each gate and overhead for each peak found. The actual amount of data to be processed can be reduced significantly by intelligent use of the gate definition and of the processing parameters for the peak detection. The scans for positive and negative peaks are done in parallel.

3.8.1. User Gates

Gates can be used to define regions of interest within an acquisition segment. This can be used to limit the amount of data to be transferred or as a starting point for peak detection processing. A gate is defined as a range of data defined by a first sample and the number of samples in the gate. Both the first sample and the number of samples must be multiples of 4. An arbitrary number of gates, up to 64, can be defined. They are used for each segment in the acquisition.

3.8.2. Simple peak detection

In this mode the data for each user gate region will be examined and the minimum and maximum values and their positions will be defined as the negative and positive peaks, respectively.

3.8.3. Peak detection with Hysteresis

This mode provides a more sophisticated method of analyzing the data in each user gate region. The final result will be at most one positive and one negative peak per gate. The processing algorithm for the positive peak search is as follows:

- p1. Starting at the first data point of the gate ($d[0]$) skip through the data until the data point value exceeds the minimum value observed by *StartDeltaPosPeak*. The search terminates at the end of the gate. Let i_s be the data point index satisfying this test, then

$$d[i_s] \geq \min(d[0] \dots d[i_s]) + \textit{StartDeltaPosPeak}.$$

- p2. Continue scanning through the data until the data point value is less than the maximum value observed less *ValidDeltaPosPeak*. The search terminates at the end of the gate. Let i_v be the data point index satisfying this test and i_M be the data point index of the maximum, then

$$d[i_M] = \max(d[i_s] \dots d[i_v]) \quad \text{and} \quad d[i_v] \leq d[i_M] - \textit{ValidDeltaPosPeak}.$$

- p3. Record the value and position of the accepted maximum data point (if any). There will be a result if condition p1 above is satisfied.

The same algorithm is repeated for the negative peak search with the obvious modifications, namely:

- n1. Starting at the first data point of the gate skip through the data until the data point value falls below the minimum value observed by *StartDeltaNegPeak*. The search terminates at the end of the gate. This means that

$$d[i_s] \leq \max(d[0] \dots d[i_s]) - \textit{StartDeltaNegPeak}.$$

- n2. Continue scanning through the data until the data point value is greater than the minimum value observed plus *ValidDeltaPosPeak*. The search terminates at the end of the gate. This means that

$$d[i_M] = \min(d[i_s] \dots d[i_v]) \quad \text{and} \quad d[i_v] \geq d[i_M] + \textit{ValidDeltaNegPeak}.$$

- n3. Record the value and position of the accepted minimum data point (if any). There will be a result if condition n1 above is satisfied.

3.8.4. Peak Interpolation

Data in the bins around the simple peak positions can be used to give better values of the peak position and amplitude. This additional processing is done in the host computer using data around the identified peaks. It can be invoked independently of how the peaks are detected. A 3-point spline calculation of the peak position and amplitude is used.

3.9. External Clock and Reference

For applications where the user wants to replace the internal clock of the Analyzer, and drive the ADC with an external source, an External Clock or Reference signal input is available. The Clock or Reference signals can be entered into the Analyzer via the MMCX CK connector on the front panel.

When using an External Clock, the user must ensure that the input signal has a frequency between 20 MHz and 2000 MHz for the AP235/AP240 models (or between 10 MHz and 500 MHz for the AP101/AP201 models) and a minimum amplitude of at least 1 V peak-to-peak (2 V for the AP101/AP201 models). The External Clock allows the Analyzer to make a voltage measurement when the clock signal passes through a predefined threshold.

For the AP240/AP235 single-channel mode and the AP101 the external clock sampling rate is the same as the external clock frequency. However for the AP240/AP235 dual-channel mode the sampling rate is ½ of the clock frequency and for the AP201 it is 2 times the clock frequency. The AP240/AP235 single-channel mode cannot be used with a Start/Stop external clock. The threshold range is variable and user selectable between ± 2 V. The signals should not exceed ±5 V amplitude.

For applications that require greater timing precision and stability than is obtainable from the internal clock, a 10 MHz Reference signal can be used. The amplitude and threshold conditions, for an External Reference, are the same as for the External Clock. If phase synchronization between several Analyzers is required, the reference signal should be applied to all of them. The allowed frequency range is [9.0 MHz, 10.2 MHz].

When using an External Clock, the **Resync** trigger mode is not functional.

3.10. Front Panel Inputs and Controls

Name		Function		Connector
		AP101/AP201	AP240/AP235	
INPUT A	Channel 1	Signal input A	Channel 1 input	BNC 1
INPUT B	Channel 2	Signal input B	Channel 2 input	BNC 2
TRIGGER IN		Trigger input		BNC 3
		User configurable (see below).		MMCX I/O A
TRIGOUT		Signal occurs after an accepted TRIGGER. It is synchronous to the acquisition Clock and can be used to trigger events synchronously to the acquisition clock.		MMCX TR
		User configurable		MMCX I/O B
CLKEXT		50 Ω External clock / External reference Input		MMCX CK
		NA	Timestamp reset	MMCX DPU P1
		NA	Timestamp reset	MMCX DPU P2

The I/O A, I/O B, P1 and P2 signals are 3.3 V compatible CMOS. This means that, on input, low is < 0.7 V and high must be in the range [1.7 V, 5.0 V]. An unconnected signal will be high. This definition ensures TTL compatibility. On output, the low level will be in the range [0 V, 0.7 V] and the high level in the range [1.7 V, 3.3 V] for HiZ. The high level output will typically give 0.8 V into 50 Ω. The I/O A or B connectors can be used for the following output signals: *10 MHz reference clock*, *Acquisition is active*, *Trigger is armed*. They can also be used for *Enable trigger input* for the AP101/AP201.

3.11. Internal Calibration

The software drivers supplied include calibration functions for the timing, gain and offset settings, which can be executed upon user request. The Analyzers are never calibrated in an “automatic” way, i.e. as a side effect of another operation. This ensures programmers have full control of all calibrations performed through software in order to maintain proper event synchronization within automated test applications.

The model AP series Analyzers include a high precision voltage source and a 16-bit DAC, used to determine the input voltage and offset calibration.

For accurate time and voltage measurements it is recommended to perform a calibration once the module has attained a stable operating temperature (usually reached with a few minutes after power on). Further calibration should not be necessary unless temperature variations occur.

3.12. Electrical, Environmental and Physical Specifications

3.12.1. Electrical

Power Consumption (W)		CURRENT REQUIREMENTS (A)			
Model	Maximum / Typical	+12V	+5V	+3.3V	-12 V
AP240	56 / 47	1.1	2.8	8.4	0.04
AP235	56 / 47	1.1	2.8	8.4	0.04
AP101	39	1.3	2.6	2.8	0.04
AP201	48	1.3	3.3	4.5	0.04

Except for the high total power consumption all modules conform to Revision 2.2 of the PCI Local Bus Specification. They use the PCI Bus at 33 MHz. All of these modules are capable of DMA transfers at rates ~100 MB/s.

3.12.2. Environmental and Physical

The modules have a Declaration of Conformity according to ISO/IEC Guide 22 and EN45014 and **CE** Marks of Compliance.

The American Standard MIL-PRF-28800F has been used as the source for many of the characteristics specified below. Those indicated with a '*' are compatible with Class 3 as defined in that standard.

Operating Temperature

0° to 40°C

The above values are for the ambient temperature of the room (or equivalent) where the PC is located. The temperature as measured on the board may well be significantly higher. On-board temperatures above 60°C should be avoided

Relative Humidity*

5 to 95% (non-condensing)

Dimensions

All AP modules conform to the PCI standard and require a long card slot. When used with the XP102 Fan the overall length can be adjusted from ~325 mm up to the long card slot maximum of 342.3 mm.

Shock*

30 G, half-sine pulse

Vibration*

5-500 Hz, random

Safety

Complies with EN61010-1

EMC Immunity

Complies with EN61326-1: Industrial Environment

EMC Emissions

Complies with EN61326-1: Class A for radiated emissions

Required Airflow

XP102 Fan unit required. AP Averagers are delivered with this Fan unit.

4. Running the AP_SSRDemo Application

4.1. Getting Started with AP_SSRDemo


Once the software and hardware installation described in Section 2 of this manual is complete, you will find AP_SSRDemo in the start menu of your computer. Start AP_SSRDemo from the menu. An initial panel of controls and a waveform display window will appear.

During startup, AP_SSRDemo searches for all Acqiris Analyzers on the PCI bus. If none are found, AP_SSRDemo will display an error message indicating this fact and stop. Turn off the computer, install and turn ON the hardware on the PCI bus and, then, restart the computer.

AP_SSRDemo is intended to verify the functionality of the card and to serve as an easy to use and convenient program for capturing, displaying, and archiving captured data in ASCII. AP_SSRDemo is provided free with the purchase of any Acqiris Analyzer.


We supply drivers on the CDROM for National Instruments LabWindows/CVI, LabVIEW, Microsoft Visual C++, Visual Basic, and Visual Basic .NET. We also provide some examples (GetStarted applications) as a guide or starting point for programmers developing their own applications to interface with the AP Analyzers or Averagers and the DC or DP series digitizer cards.

When you start AP_SSRDemo a display like that shown below appears. It has two independent windows.

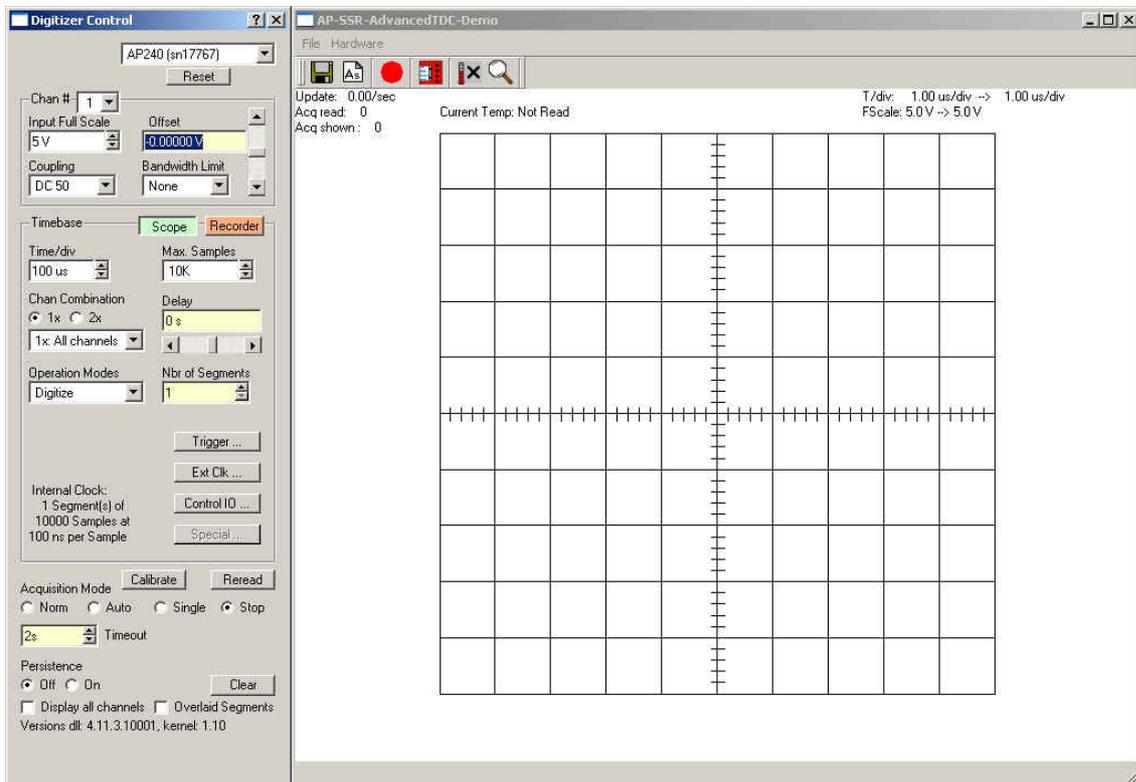

When the operation mode is changed to Peak^{TDC} or SSR, the program displays three independent windows, a Digitizer Control window, an Analyzer Control window, and a display window. The Control windows contain functions that allow you to manipulate the acquisition parameters of the card. In order to display an acquired waveform, use the Operation Mode: Digitizer, connect a signal to the input, and then select the Auto acquisition mode at the bottom of the Control panel. A waveform should appear in the display.

Other controls in the Control panel allow you to select the optimal acquisition settings for your particular signal. Each of the other Control panel functions will be discussed in detail in the following Control Panel Functions section. Here is a sample result.

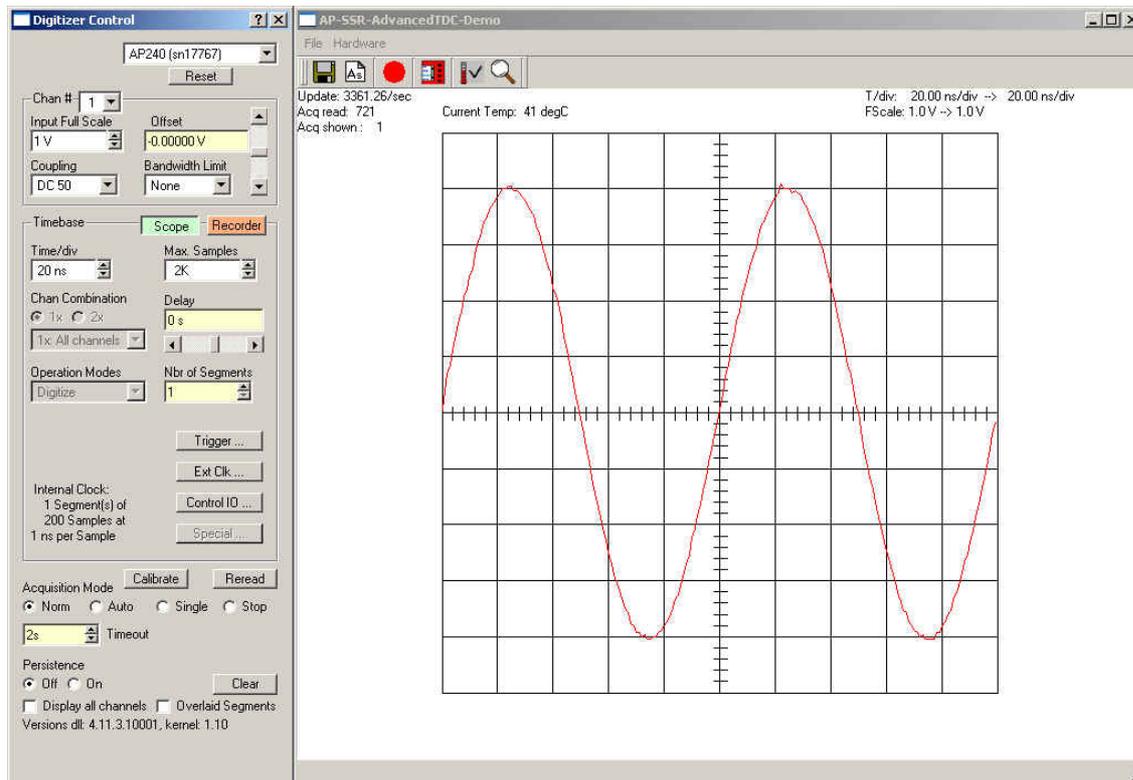

4.2. Editing Fields

Several different possibilities for editing the contents of each field are available. Here is a list:

- Click to select or do: The Scope/Recorder switch and the Acquisition mode are of this type. Select the desired state. When this is used with a text of the form A/B, A is chosen when the Black dot is present and B is chosen when it is not.
- Click to change or do: The Clear and the Calibrate action are of this type.
- Choose from a list: Coupling, and Bandwidth Limit are of this type. Do not attempt to use the keyboard on this kind of field.
- Change value up/down: Input Full Scale, Time/Div, Max Samples. Do not attempt to use the keyboard on this kind of field.
- Change value up/down multiple arrows: Nbr of Segments is of this type. The left pair of arrows modifies the value by the smallest increment possible while the right pair modifies the value by a much larger amount. Do not attempt to use the keyboard on this kind of field.
- Access a lower level menu: Trigger ... is of this type.
- Move a slider: The Offset, the Delay time, and the trigger level in the Trig Control menu are of this type. Use the mouse to slide the arrow along the scale. Numeric entry in the value window is allowed. Scientific notation and common abbreviations (m,u,n,p, S, V) can be used.

Values above or to the left of other values have precedence over the latter.


Selection fields that are gray are not active for the current configuration.

4.3. Digitizer Control Panel and Functions

The Digitizer Control panel provides all of the features needed to use the Analyzer in the digitizer mode. Some of these controls are also active when the analyzer mode is in use. All of the controls are described in the sections below:

4.3.1. Input Voltage Scale, Offset, Coupling and Bandwidth

Settings for selecting the channel, voltage scale, offset, coupling, and bandwidth limit are available in the portion of the Digitizer Control panel shown to the right. The settings are also applicable for the Analyzer.


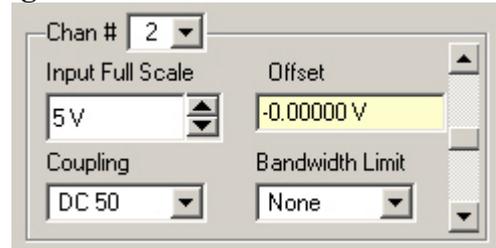


4.3.2. Digitizer Timebase and more

Digitizer timebase characteristics are set using the settings in the Digitizer Control panel shown here. They are not applicable to the Analyzer Operation Mode that will be discussed below. However, the settings in the *Trigger...*, *Ext Clk...*, and *Control IO...* sub-menus are applicable to both modes. The *Ext Clk...* and *Control IO...* sub-menus are for advanced users and will not be presented in detail.

There are two timebase control panel modes available, Scope Mode and Recorder mode. The choice of mode is entirely determined by the preference of the user. Both modes provide the same set of functional capabilities but display the settings slightly differently. In both modes it is possible to set the memory depth and the sample speed independently within the limits of the maximum available memory. However the two modes differ in the following way.

In Scope Mode, shown to the right, the user selects the time window (time per div) and the sampling rate will be adjusted to the fastest possible within the limits of the available memory. The maximum available memory is programmable.

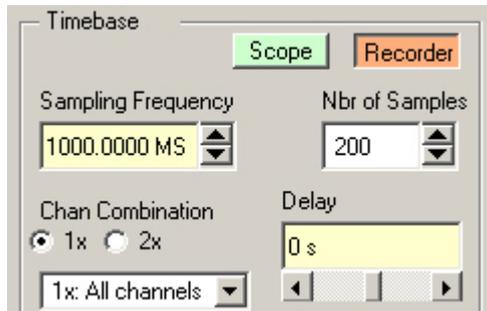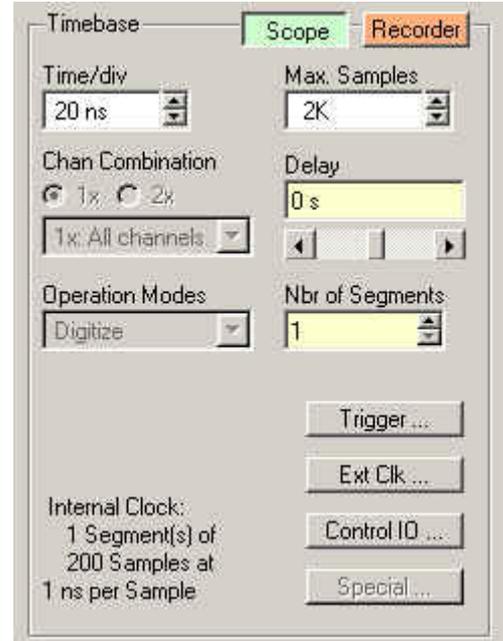

In Transient recorder mode, shown to the left, the user selects the sampling rate and the number of samples. This is closer to the kind of setup used in the analyzer mode.

The trigger *Delay* setting positions the acquisition with respect to the trigger. Pre-trigger values are negative and are selectable in the range of [0 to $-(10 * \text{timebase})$]. The value $(10 * \text{timebase})$ is equivalent to the time window or capture-time. 100% pre-trigger is achieved by setting the Delay to $-\text{capture-time}$. A 100% pre-trigger delay means

the trigger point is at the right edge of the display and all of the data in the acquired waveform were captured just prior to the trigger event.


Maximum Post trigger delay settings depend on the sample interval used. Post trigger delays are adjustable up to 200 Mpoints.

4.3.3. Operation Mode

When the program first appears it is in the Digitize operating mode. Sooner or later you will want to go over to one of the analyzer modes. You can only change operation mode if the acquisition mode is STOP. After this change several fields in the Digitizer Control window will become GREY to show that they are no longer active. The Analyzer Control window should be present after the switch. If it is not, click on the SSR... or TDC... button to make it reappear.

4.3.4. Trigger

Trigger characteristics are set using the settings in the *Trig Control* sub menu. The trigger *source* setting can be either Internal (trigger on any one of the input channels) or External (trigger on an external input). The trigger *class* is always Edge.

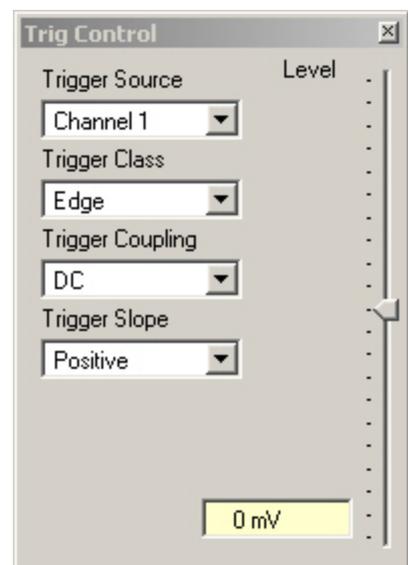

The *Coupling* setting can be AC, DC, or Hfrej coupling. Use the *Level* section to set the desired trigger level. The External Trigger Source menu will also allow selection of the desired Trigger Range and Bandwidth. The relationship between coupling, trigger source, and trigger level was discussed in section 3.3.3, **Trigger Level**.

Use the *Slope* section to choose the appropriate edge for the desired trigger. Alternatively, a window trigger can be defined with two thresholds. Finally, the HF-Divide trigger can also be chosen.

4.3.5. External Clock

The window will allow you to define the desired clock type, threshold, frequency, sampling rate, and delay. You can also ask for a calibration that will take into account a new value of the external clock frequency.

4.3.6. Control IO Settings

The window will allow you to select the desired functionality for the Control I/O A and B MMCX connectors as well as the offset to be used for the Trigger Out MMCX level.

4.3.7. Acquisition Mode

The acquisition mode of the analyzer/digitizer is selected using one of four acquisition mode buttons in the Acquisition section of the Control panel. Available acquisition modes are Auto (digitizer only), Normal, Single and Stop. The use of each of the acquisition modes is described below. In the discussion, a valid trigger indicates a trigger signal meeting the trigger conditions at a time when the device is armed and ready to acquire data.

Stop will stop the acquisition and hold the latest complete acquisition on the display.

Single mode is used in order to capture one event at the first valid trigger or the desired number of segments for that many valid triggers. It freezes the acquisition in the device's memory, and on the display, until the user requests another acquisition. After an acquisition is taken in Single mode, the device will ignore subsequent trigger events until the Single button is pressed again or another acquisition mode is selected. Pressing the single button re-arms the trigger and captures a new acquisition.

Normal mode is used to continuously acquire waveforms into memory for valid trigger events. The display is updated with the new acquisitions. Acquisitions will happen for each valid trigger until either single or Stop is pressed.

If the trigger indicator at the lower left of the waveform display does not blink green, the device is not receiving valid triggers. Check the trigger settings or select Auto mode in order to obtain a display of the waveform.

Auto mode will acquire and display waveforms according to the trigger settings if a valid trigger is present within an adjustable timeout interval. If a valid trigger is not available within this interval, the digitizer generates its own trigger in order to digitize and display whatever signal is at the input at that time. It can only be used when in digitizer mode.

If valid triggers are received at a high enough rate, Auto mode behavior differs little from normal mode behavior on the display. However if the trigger rate is too low or trigger settings are not appropriate for the characteristics of the signal, the input signal will not be stable on the display. This occurs because the internal auto-trigger generated by the digitizer is asynchronous to the input signal.

Auto mode is often used to aid in setup when the input signal must be quickly characterized in order to determine proper trigger settings for Normal or Single mode acquisitions.

4.3.8. Display options

The controls at the very bottom of the Digitizer Control window are used to choose whether the event display of data points is for the last acquisition only (Off) or for all the acquired events (On). A simple Clear button is also present that erases all data points in the display. The display can also be used to show both channels at once. If desired, the individual segments of an acquisition can be shown stacked on top of each other as opposed to spread out along the horizontal axis.


4.4. Analyzer Control Panel and Functions

The analyzer mode differs fundamentally from the digitizer mode since nothing happens before the trigger is received. There is no possibility of acquiring pre-trigger data. The Analyzer Control window can be closed by clicking the **X** in the upper right hand corner of the window.

The functionality of each part of the window will be described in the next few sections.


4.4.1. Analyzer Timebase and Memory

When the analyzer mode is selected the number of data points, the number of segments, and the delays after receiving the trigger and taking the data are set in the Analyzer Control window. Only the sampling rate is taken over from the Digitizer Control window. However, the Digitizer Control window display of the acquisition's horizontal parameters is still valid.

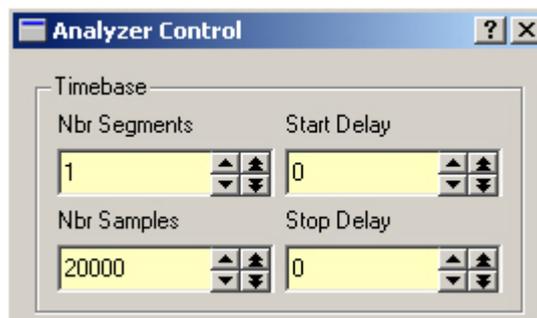

4.4.2. User Gate definition

For further processing or readout data reduction the analyzer mode has the concept of gates, regions of interest after the trigger for each segment of an acquisition. The gates can be defined in the part of the Analyzer Control window shown here. Each channel has its own definition.

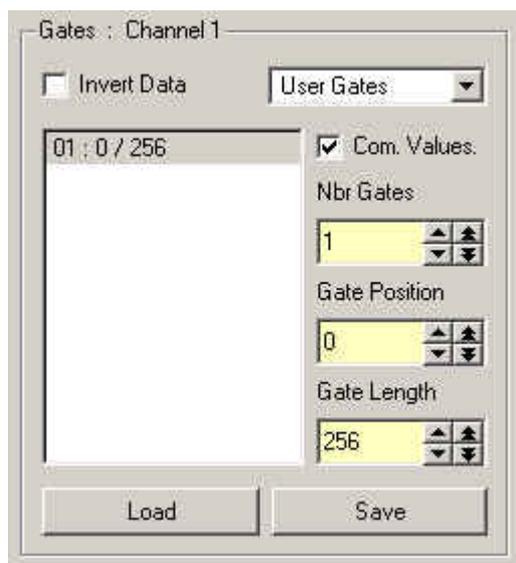

For the User Gate mode shown below the number of gates together with their length and starting position can be modified as needed. As soon as more than one Gate is available the field to the right of the # sign can be used to select a gate for editing. Note that the Gate Length is to the left of the Gate Position, which implies that it has priority. This means that Gate Position will be forced to satisfy the relation

$$\text{Gate Position} \leq \text{Waveform Length} - \text{Gate Length}$$



The Individual gate mode is what will normally be used when there is more than one region of interest. The Common gate mode is useful to judge the readout overhead associated with more data.

The Invert Data option can be enabled to allow the peak definition to select negative-going peaks.

The Load and Save options can be used to store/recall predefined sets of gates.

4.4.3. Threshold Gate definition

For the Threshold Gate mode, the selectors for the Maximum number of gates per segment, the Fixed Number of Samples after the first one satisfying the threshold condition, the number of PreSamples, and the number of PostSamples will be shown. Setting MaxNbrGates = 0 and/or FixedNbrSamples = 0 disables their functions. The Invert Data option can be enabled to allow the Threshold Gates to select data below the Threshold setting and the peak definition to select negative-going peaks.

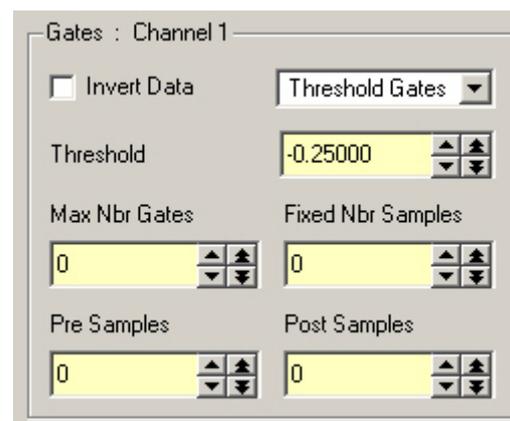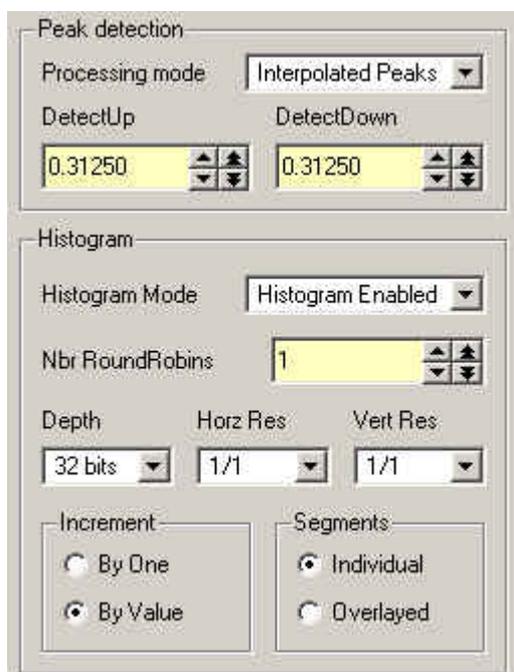

4.4.4. Peak^{TDC} mode - Peak and Histogram definition

The fragment of the Analyzer Control window shown here can be used to select the peak detection parameters and the histogram to be accumulated.

The Processing mode allows a choice between Standard and Interpolated Peaks.

The DetectUp and DetectDown parameters are measured in volts and give the desired StartDeltaPosPeakV and ValidDeltaPosPeakV respectively.

Note that Nbr RoundRobins is the number of acquisitions to be taken in each segment. However, each segment gets only one out of every Nbr Segments triggers.

Depth can be used to trade off between the maximum number of events that can be accumulated and the space required for the histogram. Horz Res and Vert Res define characteristics of the histogram to take advantage of the interpolated peaks better precision. Overlaid segments will accumulate all segments into a common horizontal range in the histogram; individual segments will separated the segments, each in its own histogram.

4.4.5. Additional Controls & Time Stamp data

At the bottom of the Analyzer Control window, several additional options can be set:

the timing of the Trigger Output pulse can be chosen as was discussed in 3.6.1, **Random and Synchronized Triggers**.

when Gates are in use the display can be configured to show the data transferred either


At Acquisition time – the data points will be shown on a horizontal axis covering the entire acquisition or segment

Concatenated full scale – the data points will be shown one after another with no indication of the gap between them. The horizontal scale is the same as in the case above. This is useful to judge the effectiveness of the data reduction.

Concatenated auto scale – the data points will be shown one after another with no indication of the gap between them. The horizontal scale is adjusted for each event to show the data transferred with maximum resolution

The FP I/O P1 and P2 DPU Control can be set to allow for the hardware reset of the trigger time stamp counter.

The time stamp data associated with each trigger can be displayed by clicking on the Show TimeStamps button. Thereafter, the display can be updated on demand.

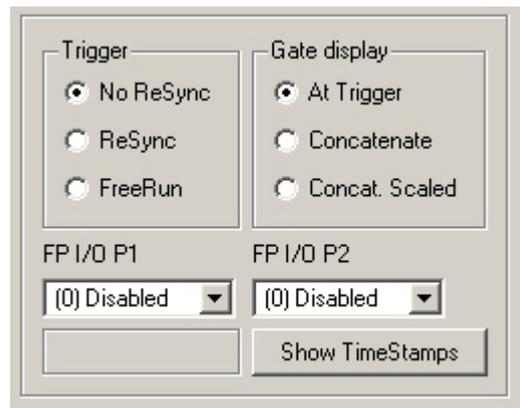

The choices above apply to the selected channel only.

Closing the Display window will terminate the execution of the AP_SSRDemo program.

If desired the Display of waveforms can be turned off. To do this, select the top line menu item **H**ardware to get to click on the option **D**isable display.

4.4.6. Waveform storage

Acquired waveforms can be saved to disk by the AP_SSRDemo program. This feature can be accessed by way of the **F**ile item in the top line menu as **S**ave Waveform or **S**ave Waveform **A**s... Alternatively, the two icons shown here can be used to activate these functions. If this is the first time waveform storage is being used in the current run of the program, the **S**ave Waveform **A**s... variant must be used so that a file name and destination folder can be chosen. Thereafter, the **S**ave Waveform command can also be used. It will save the data by overwriting the file previously created. The data is stored in ASCII format as Raw ADC values. In order to convert the Raw ADC values into Volts the following formula should be used.

$$V = vGain * data - vOffset$$

Where vGain and vOffset can be found in the header of the file.

4.4.7. Control panel reactivation



If the Digitizer Control panel has been closed it can be made to reappear by either, using the top line menu item **H**ardware to get to **D**igitizer Ctrl... or clicking on the icon shown here. If the Analyzer

Control panel needs to be regenerated it can be done by clicking on the **SSR...** button of the Digitizer Control panel when the SSR operation mode is selected.


4.4.8. Temperature display

If desired the current *in-situ* operating temperature of the module can be displayed by either, selecting the top line menu item **Hardware** to get to **Read Temperature**, or clicking on the thermometer icon.

4.4.9. Zoom display control



A waveform zoom feature is available by either, selecting the top line menu item **Hardware** to get to **Zoom...**, or clicking on the magnifying glass icon.

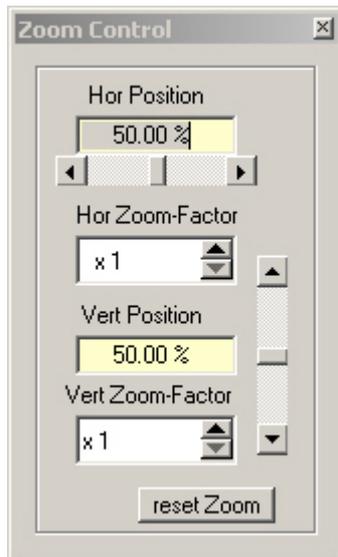

This will activate the Zoom Control menu shown here. The horizontal zoom is controlled by the Hor Position, with 50% in the middle of the screen. The vertical zoom is controlled by the Vert Position, with 0% in the middle of the screen and +50% meaning that the display should be shifted up by $\frac{1}{2}$ the screen height. New zoom values can be used on existing data by pushing the Reread button in the Digitizer Control window.

To clear the display in persistence display mode, press the “Clear” button that appears on the main display window.

5. Running the APx01Demo Application

5.1. Getting Started with APx01Demo


Once the software and hardware installation described in Section 2 of this manual is complete, you will find APx01Demo in the start menu of your computer. Start APx01Demo from the menu. An initial APx01Demo load screen will appear followed by a panel of controls and a waveform display window.

During startup, APx01Demo searches for all Acqiris Analyzers on the PCI bus. If none are found, APx01Demo will display an error message indicating this fact and stop. Turn off the computer, install and turn ON the hardware on the PCI bus and, then, restart the computer.

APx01Demo is intended to verify the functionality of the card and to serve as an easy to use and convenient program for capturing, displaying, and archiving captured data in ASCII format. APx01Demo is provided free with the purchase of any Acqiris Analyzer.


We supply drivers on the CDROM for National Instruments LabWindows/CVI, LabVIEW, Microsoft Visual C++, and Visual Basic. We also provide some examples (GetStarted applications) as a guide or starting point for programmers developing their own applications to interface with the AP Analyzers or Averagers and the DC or DP series digitizer cards.

When you start APx01Demo a display like that shown below appears.

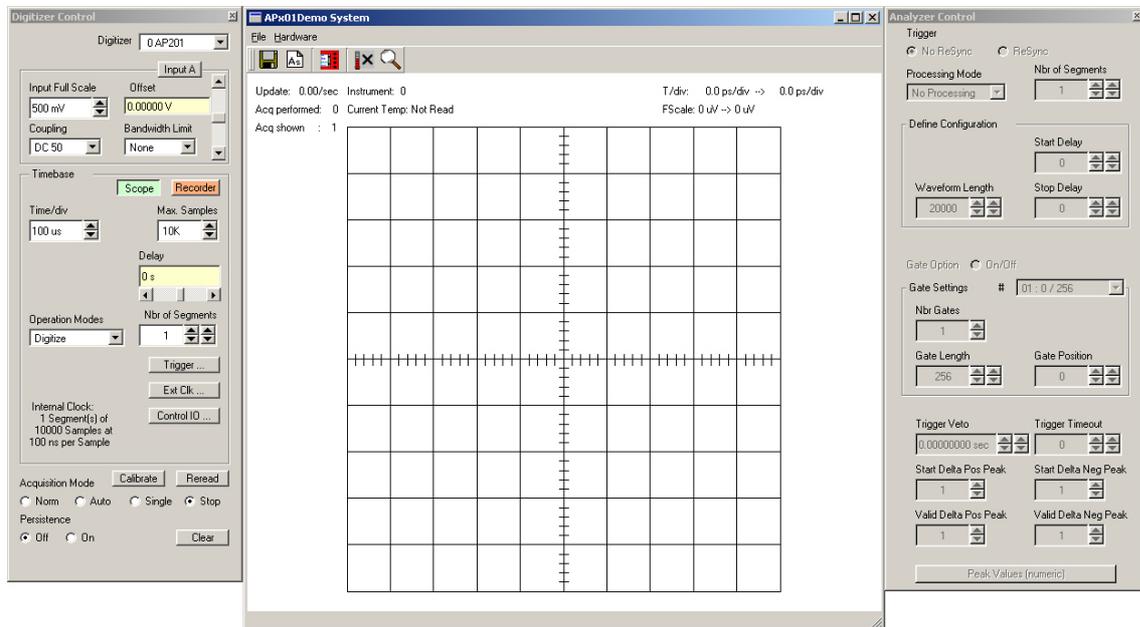

APx01Demo displays three independent windows, a Digitizer Control window, an Analyzer Control window, and a display window. The Control windows contain functions that allow you to manipulate the acquisition parameters of the card. In order to display an acquired waveform, connect a signal to the input, and then select the Auto acquisition mode at the bottom of the Control panel. A waveform should appear in the display. Note that this is best done with the Operation Mode: Digitizer.

Other controls in the Control panel allow you to select the optimal acquisition settings for your particular signal. Each of the other Control panel functions will be discussed in detail in the following Control Panel Functions section. Here is a sample result.

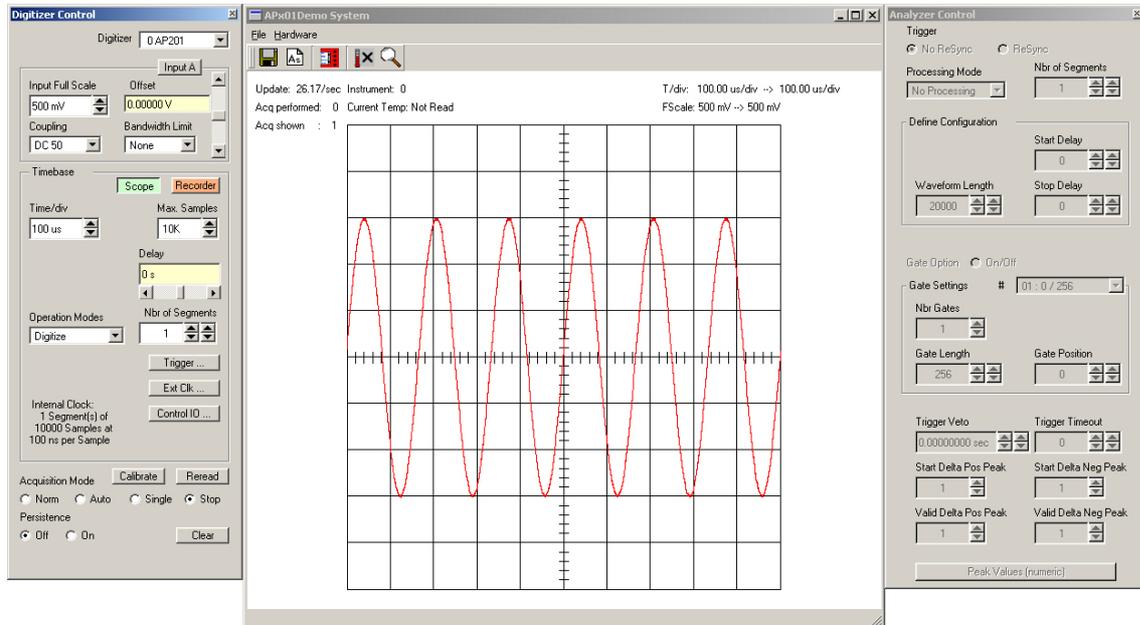

5.2. Editing Fields

Several different possibilities for editing the contents of each field are available. Here is a list:

- Click to select or do: The Scope/Recorder switch and the Acquisition mode are of this type. Select the desired state. When this is used with a text of the form A/B, A is chosen when the Black dot is present and B is chosen when it is not.
- Click to change or do: The Input A/ Input B, and the Calibrate action are of this type.
- Choose from a list: Coupling, and Bandwidth Limit are of this type. Do not attempt to use the keyboard on this kind of field.
- Change value up/down: Input Full Scale, Time/Div, Max Samples. Do not attempt to use the keyboard on this kind of field.
- Change value up/down multiple arrows: Nbr of Segments is of this type. The left pair of arrows modifies the value by the smallest increment possible while the right pair modifies the value by a much larger amount. Do not attempt to use the keyboard on this kind of field.
- Access a lower level menu: Trigger ... is of this type.
- Move a slider: The Offset, the Delay time, and the trigger level in the Trig Control menu are of this type. Use the mouse to slide the arrow along the scale. Numeric entry in the value window is allowed. Scientific notation and common abbreviations (m,u,n,p, S, V) can be used.

Values above or to the left of other values have precedence over the latter.


Selection fields that are gray are not active for the current configuration.

5.3. Digitizer Control Panel and Functions

The Digitizer Control panel provides all of the features needed to use the Analyzer in the digitizer mode. Some of these controls are also active when the analyzer mode is in use. All of the controls are described in the sections below:

5.3.1. Input Voltage Scale, Offset, Coupling and Bandwidth

Settings for selecting the Input A or Input B BNC connector, voltage scale, offset, coupling, impedance, and bandwidth limit are available in the portion of the Digitizer Control panel shown to the right. The settings are also applicable for the Analyzer.


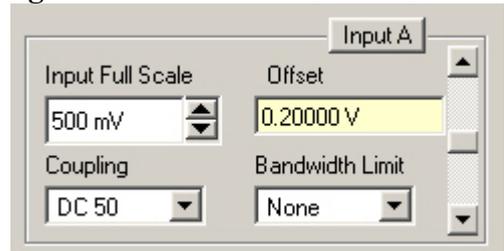


5.3.2. Digitizer Timebase and more

Digitizer timebase characteristics are set using the settings in the Digitizer Control panel shown here. They are not applicable to the Analyzer Operation Mode that will be discussed below. However, the settings in the *Trigger...*, *Ext Clk...*, and *Control IO...* sub-menus are applicable to both modes. The *Ext Clk...* and *Control IO...* sub-menus are for advanced users and will not be presented in detail.

There are two timebase Control panel modes available, Scope Mode and Recorder mode. The choice of mode is entirely determined by the preference of the user. Both modes provide the same set of functional capabilities but display the settings slightly differently. In both modes it is possible to set the memory depth and the sample speed independently within the limits of the maximum available memory. However the two modes differ in the following way.

In Scope Mode, shown to the right, the user selects the time window (time per div) and the sampling rate will be adjusted to the fastest possible within the limits of the available memory. The maximum available memory is programmable.

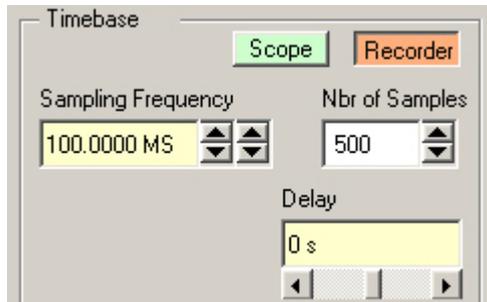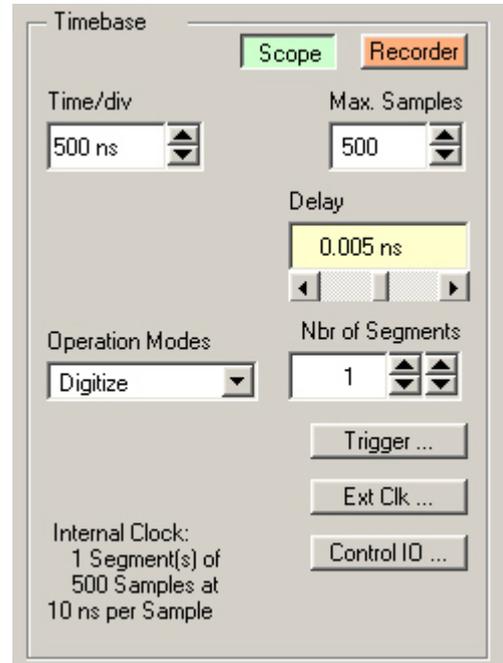

In Transient recorder mode, shown to the left, the user selects the sampling rate and the number of samples. This is closer to the kind of setup used in the analyzer mode.

The trigger *Delay* setting positions the acquisition with respect to the trigger. Pre-trigger values are negative and are selectable in the range of [0 to $-(10 * \text{timebase})$]. The value $(10 * \text{timebase})$ is equivalent to the time window or capture-time. 100% pre-trigger is achieved by setting the Delay to $-\text{capture-time}$. A 100% pre-trigger delay means the trigger point is at the right edge of the display and all of the data in the acquired waveform were captured just prior to the trigger event.

Maximum Post trigger delay settings depend on the sample interval used. Post trigger delays are adjustable up to 200 Mpoints.


5.3.3. Trigger

Trigger characteristics are set using the settings in the *Trig Control* sub menu. The *Trigger Source* setting can be either Internal (trigger on any one of the input channels) or External (trigger on an external input). The *Trigger Class* is always Edge.

The *Coupling* section also contains a selection for AC or DC coupling. AC coupling on the trigger sets the trigger circuit to an AC Low Frequency Reject mode. In this mode the trigger is auto-level and the level value will be ignored.

Use the *Level* section to set the desired trigger level. The relationship between coupling, trigger source, and trigger level was discussed in section 3.3.3, **Trigger Level**.

Use the *Slope* section to choose the appropriate edge for the desired trigger.

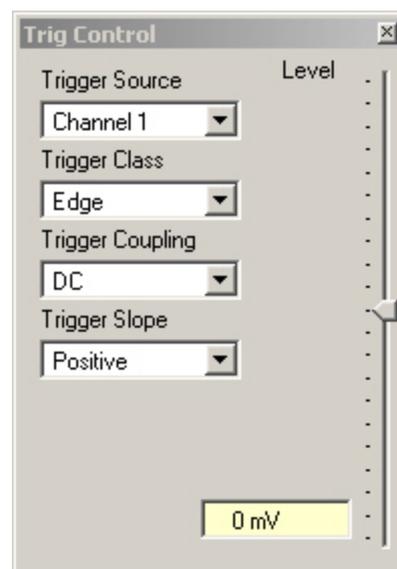

5.3.4. Acquisition Mode

The acquisition mode of the analyzer/digitizer is selected using one of four acquisition mode buttons in the Acquisition section of the Control panel. Available acquisition modes are Auto (digitizer only), Normal, Single and Stop. The use of each of the acquisition modes is described below. In the discussion, a valid trigger indicates a trigger signal meeting the trigger conditions at a time when the digitizer is armed and ready to acquire data.

Stop will stop the acquisition and hold the latest complete acquisition on the display.

Single mode is used in order to capture one event at the first valid trigger or the desired number of segments for that many valid triggers. It freezes the acquisition in the digitizer's memory, and on the display, until the user requests another acquisition. After an acquisition is taken in Single mode, the digitizer will ignore subsequent trigger events until the Single button is pressed again or another acquisition mode is selected. Pressing the single button re-arms the trigger and captures a new acquisition.

Normal mode is used to continuously acquire waveforms into memory for valid trigger events. The display is updated with the new acquisitions. Acquisitions will happen for each valid trigger until either single or Stop is pressed.

If the trigger indicator at the lower left of the waveform display does not blink green, the digitizer is not receiving valid triggers. Check the trigger settings or select Auto mode in order to obtain a display of the waveform.

Auto mode will acquire and display waveforms according to the trigger settings if a valid trigger is present within a timeout interval. If a valid trigger is not available within this interval, the digitizer generates its own trigger in order to digitize and display whatever signal is at the input at that time. It can only be used when in digitizer mode.

If valid triggers are received at a high enough rate, Auto mode behavior differs little from normal mode behavior on the display. However if the trigger rate is too low or trigger settings are not appropriate for the characteristics of the signal, the input signal will not be stable on the display. This occurs because the internal auto-trigger generated by the digitizer is asynchronous to the input signal.

Auto mode is often used to aid in setup when the input signal must be quickly characterized in order to determine proper trigger settings for Normal or Single mode acquisitions.

5.3.5. Persistence display

The controls at the very bottom of the Digitizer Control window are used to choose whether the event display of data points is for the last acquisition only (Off) or for all the acquired events (On). A simple Clear button is also present that erases all data points in the display

5.4. Analyzer Control Panel and Functions


The analyzer mode differs fundamentally from the digitizer mode since nothing happens before the trigger is received. There is no possibility of acquiring pre-trigger data. The Analyzer Control window can be closed by clicking the **X** in the upper right hand corner of the window.

The functionality of each part of the window will be described in the next few sections.

5.4.1. Analyzer Timebase and Memory

When the analyzer mode is selected the number of data points, the number of segments, and the delays after receiving the trigger and taking the data are set in the Analyzer Control window. Only the sampling rate is taken over from the Digitizer Control window. However, the Digitizer Control window display of the acquisition's horizontal parameters is still valid.

In addition, the timing of the Trigger Output pulse can be chosen as was discussed in 3.6.1, **Random and Synchronized Triggers**.

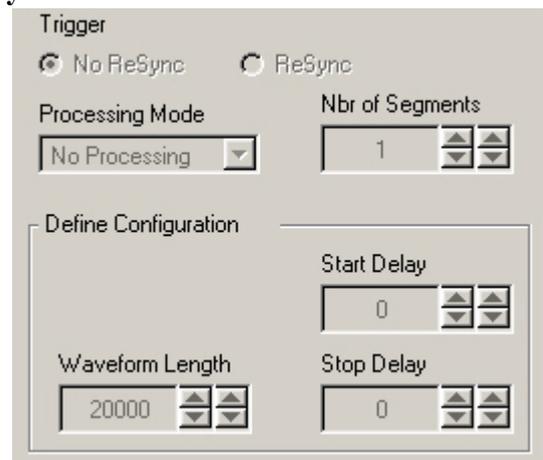


5.4.2. Gate definition

For further processing or readout data reduction the analyzer mode has the concept of gates, regions of interest after the trigger for each segment of an acquisition. The gates can be defined in the part of the Analyzer Control window shown here. The On/Off switch selects whether the readout will be done with/without the gates. This choice is only available for the case of No Processing.

The number of gates together with their length and starting position can be modified as needed. As soon as more than one Gate is available the field to the right of the # sign can be used to select a gate for editing. Note that the Gate Length is to the left of the Gate Position which implies that it has priority. This means that Gate Position will be forced to satisfy the relation

$$\text{Gate Position} \leq \text{Waveform Length} - \text{Gate Length}$$

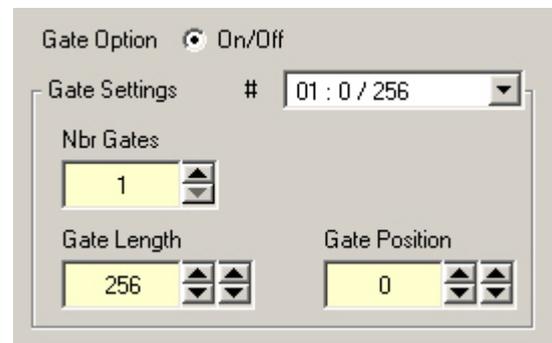


5.4.3. Trigger Veto and Timeout

These parameters were discussed in 3.3.7, **Trigger Veto and Timeout** and can be set in the Analyzer Control window section shown.



5.4.4. Hysteresis Mode Peak Parameters

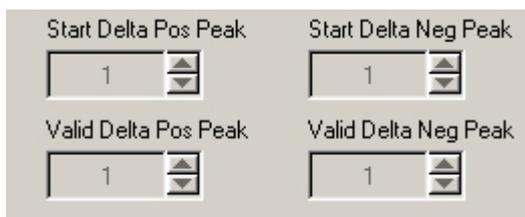


These parameters were discussed in 3.8.3, **Peak detection with Hysteresis** and their values can be set here.

5.4.5. Display of Peaks

A button to activate the display of the position and magnitude of the peaks can be found at the very bottom of the Analyzer Control window. If activated a new window is opened for the peak display. The Update button must be used every time new data is desired. The displayed peaks can be sorted either by Gates within Segments or by Segments for each Gate. Click on the Button to toggle between the two modes.

Note that the positions are given in samples from the start of the gate and the amplitudes are in raw 8-bit ADC counts in the range [-128,+127]. The interpolated values are shown when activated.

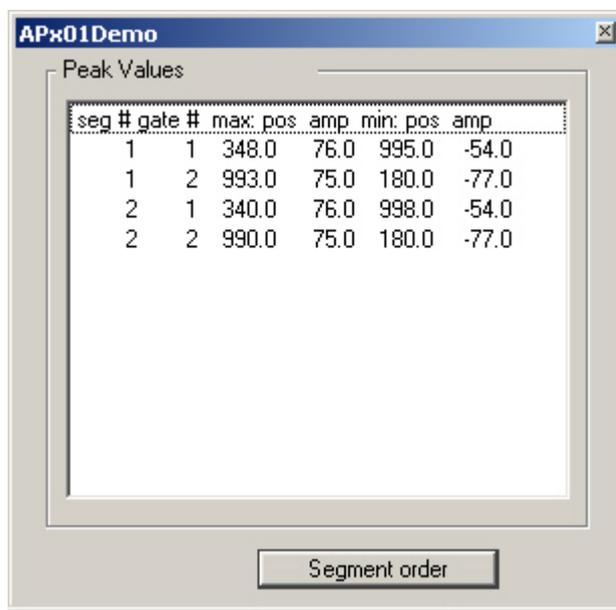

5.5. Display Window Functions

In the digitizer mode, and the analyzer mode with No Processing and Gate Option Off, the display will show the data points as a function of time. If several segments are taken they will be shown one next to the other. For the analyzer mode with No Processing and Gate Option On only the data for the gates will be shown for each segment. If Extrema or Hysteresis peak processing is enabled the position of each of the peaks will be indicated by the appropriate arrow on the display.

Closing the Display window will terminate the execution of the APx01Demo program.

If desired the Display of waveforms can be turned off. To do this, select the top line menu item **H**ardware to get to click on the option **D**isable display.

Other features of the Display window are discussed in the following sections.

5.5.1. Waveform storage

Acquired waveforms can be saved to disk by the APx01Demo program. This feature can be accessed by way of the **F**ile item in the top line menu as **S**ave Waveform or **S**ave Waveform **A**s... Alternatively, the two icons shown here can be used to activate these functions. If this is the first time waveform storage is being used in the current run of the program, the **S**ave Waveform **A**s... variant must be used so that a file name and destination folder can be chosen. Thereafter, the **S**ave Waveform command can also be used. It will save the data by overwriting the file previously created. The data is stored in ASCII format as Raw ADC values. In order to convert the Raw ADC values into Volts the following formula should be used.

$$V = vGain * data - vOffset$$

Where vGain and vOffset can be found in the header of the file.

5.5.2. Control panel reactivation



If either one of the Control panels have been closed they can be made to reappear by either, using the top line menu item **H**ardware to get to **D**igitizer Ctrl... or **A**nalyzer Ctrl..., or clicking on one of the two icons shown here.


5.5.3. Temperature display

If desired the current *in-situ* operating temperature of the module can be displayed by either, selecting the top line menu item **H**ardware to get to **R**ead **T**emperature, or clicking on the thermometer icon.

5.5.4. Zoom display control



A waveform zoom feature is available by either, selecting the top line menu item **H**ardware to get to **Z**oom..., or clicking on the magnifying glass icon.

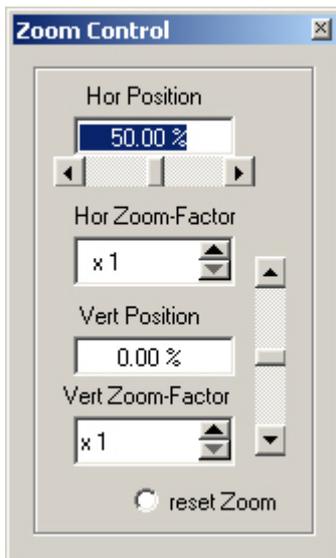

This will activate the Zoom Control menu shown here. The horizontal zoom is controlled by the Hor Position, with 50% in the middle of the screen. The vertical zoom is controlled by the Vert Position, with 0% in the middle of the screen and +50% meaning that the display should be shifted up by 1/2 the screen height. New zoom values can be used on existing data by pushing the Reread button in the Digitizer Control window.

To clear the display in persistence display mode, press the “Clear” button that appears on the main display window.

6. Appendix A: XA100 BNC Input Overvoltage Protection Specifications and User Instructions

The XA100 Acqiris option contain two elements :

- 3dB BNC Attenuator
- 90V BNC Spark Gap


The XA100 overvoltage protection kit allows the digitizers with 50 Ohm only input (without high impedance input) to withstand a high voltage surge of up to 5 kV.

XA100 SPECIFICATIONS :

- Frequency range : DC to 1GHz
- Input impedance : 50 Ohm
- Attenuation : 3 dB
- VSWR, DC to 1GHz : < 1.25
- Impulse discharge current (8/20 us pulse) : 10 kA
- Weight : 60 g
- Length : 85 mm

ATTENTION :

- **BOTH** the 90V spark gap and the 3dB attenuator **MUST** be used in order to absorb very high current.
- In order to limit the input current, the attenuator **MUST** be inserted between the spark gap and the input of the module, as shown in the figure below.

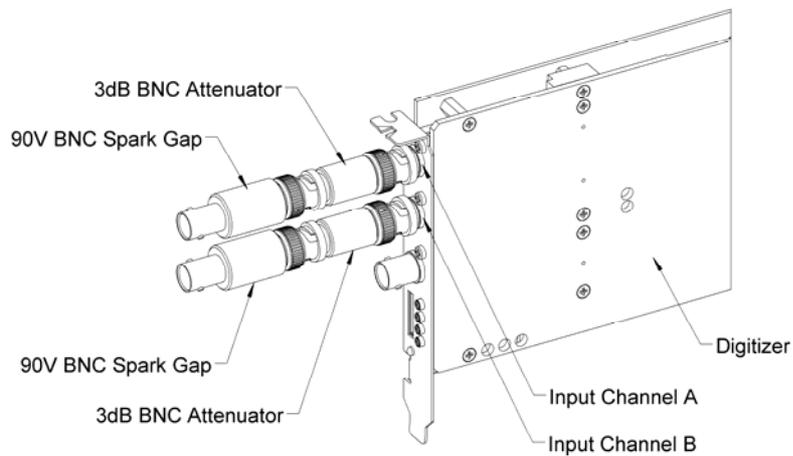

7. Appendix B: XP102 Fan Unit for the AP Family of Modules

Specification and Assembly Instructions

SPECIFICATIONS :

Airflow : 16 m³/h, typical speed 6500 tr/min, maximal noise 30 dB(A).

Nominal Voltage : 5 VDC.

Power input : 2.1 W.

With electronic protection against reverse polarity, impedance-protected against blocking and overloading.

Operating temperature : 0 °C to 70 °C.

Life expectancy (at 25 °C) : 50000 hours.

Safety approvals : CE, VDE, UL, and CSA.

Materials : Box of stainless steel, aluminium and fan of fiberglass reinforced plastic, PBT housing, PA impeller.

HOW TO INSTALL/REMOVE THE FAN ON A DIGITIZER :

1. Remove the digitizer from the PC, prior to installing or removing the fan unit.
2. Twist the fan's cable a few times in order for the 2 wires to stay together.
3. Introduce the cable between the spacers "A" and "B" at the back of the digitizer.
4. Push the fan unit into the back of the digitizer so that the two springs lock into place.
5. Insure that the two springs are safely attached behind the two spacers and that the fan's box is in contact with the digitizer cover.
6. Connect the fan's cable as shown in the figure below. The position of the power connector can vary between different digitizer models.
7. To remove the fan unit, first unplug its power cable and then press on both fan springs and pull it away from the digitizer.

HOW TO INSTALL/REMOVE THE DIGITIZER IN A PCI SLOT OF A PC :

1. Install the fan on the digitizer as described above.
2. Release the three screws "R" from the adjustable retainer and move the retainer against the fans.
3. Introduce the digitizer into the PCI slot and tighten the front panel screw.
4. Push the retainer into the PCI slot guide and tighten the three screws "R". The digitizer is now installed.
5. To remove the digitizer it is necessary to first release the three screws "R" and pull back the adjustable retainer.

