

Using the BYETTA Pen^{1,2}

NEW PEN SETUP – DO THIS ONE TIME ONLY



STEP A Check the Pen



- Wash hands prior to use.
- Check pen label to make sure it is the right dose – 5 µg or 10 µg as prescribed.
- Pull off the blue pen cap.



- Check ^{Pr}BYETTA™ (exenatide) in the cartridge. The liquid should be clear, colourless and free of particles. If it is not, do not use.

Note: A small air bubble in the cartridge is normal.

STEP B Attach the Needle



- Remove paper tab from outer needle shield.
- Push outer needle shield containing the needle straight onto the pen, then screw needle on until secure.



- Pull off outer needle shield. Do not throw away.



- Pull off inner needle shield and throw away. A small drop of liquid may appear. This is normal.

STEP C Dial the Dose



- Check that the **5** is in the dose window. If not, turn dose knob clockwise until it stops and the **5** is in the dose window.



- Pull dose knob out until it stops and the **5** is in the dose window.



- Turn dose knob clockwise until it stops at **5** or **10**. Make sure that the **5** or **10** with the line under it is in the centre of the dose window.

STEP D Prepare the Pen



- Point the needle of the pen up and away from you.



- PUSH and HOLD**
- Use thumb to firmly push the injection button in until it stops, then continue holding the injection button in while slowly counting to 5.
 - If you do not see a stream or several drops come from the needle tip, repeat Steps C and D.



- Pen preparation is complete when the **5** is in the centre of the dose window AND you have seen a stream or several drops come from the needle tip.

STEP E Complete New Pen Setup



- Turn dose knob clockwise until it stops and the **5** is in the dose window.

- New Pen Setup is now done. Do not repeat Section 2 – New Pen Setup for routine use; if you do, you will run out of BYETTA before 30 days of use.

- Go to ROUTINE USE Section, Step 3 for instructions on how to inject your first routine dose.
- inject your first routine dose.

BYETTA is used to improve blood sugar control in adults with type 2 diabetes in combination with metformin and/or a sulfonylurea when metformin and/or a sulfonylurea plus diet and exercise have failed to adequately control blood sugar levels. Continue to follow your diet and exercise plan.

BYETTA helps your body release more insulin when your blood sugar is high. This helps to improve your blood sugar control.

The most common side effects with BYETTA include nausea, vomiting, diarrhea, feeling jittery, dizziness, headache, acid stomach and heartburn.

References: 1. BYETTA 1.2 mL Pre-filled Pen: User Manual. Eli Lilly Canada Inc. January 11, 2011. 2. BYETTA 2.4 mL Pre-filled Pen: User Manual. Eli Lilly Canada Inc. January 11, 2011.

BYETTA is a trademark of Amylin Pharmaceuticals, Inc. and is used under license by Eli Lilly Canada Inc. © 2011 Amylin Pharmaceuticals, Inc. and Eli Lilly and Company. All rights reserved.

ROUTINE USE

STEP 1 Check the Pen




- Wash hands prior to use.
- Check pen label to make sure it is the right dose pen.
- Pull off the blue pen cap.



- Check BYETTA in the cartridge.
- The liquid should be clear, colourless and free of particles. If it is not, do not use.

Note: A small air bubble will not harm you or affect your dose.

STEP 2 Attach the Needle

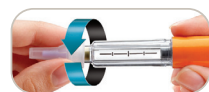

- Remove paper tab from outer needle shield.
- Push outer needle shield containing the needle straight onto the pen, then screw needle on until secure.



- Pull off outer needle shield. Do not throw away.



- Pull off inner needle shield and throw away. A small drop of liquid may appear. This is normal.

STEP 3 Dial the Dose



- Check that the **5** is in the dose window. If not, turn dose knob clockwise until it stops and the **5** is in the dose window.



- Pull dose knob out until it stops and the **5** is in the dose window.



- Turn dose knob clockwise until it stops at **5** or **10**. Make sure that the **5** or **10** with the line under it is in the centre of the dose window.

STEP 4 Inject the Dose



- Grip pen firmly.
- Insert needle into skin using hygienic injection technique recommended by your healthcare professional.



- PUSH and HOLD**
- Use thumb to firmly push injection button in until it stops, then continue holding the injection button in while slowly counting to 5 in order to get a full dose.
 - Remove needle from skin.



- Injection is complete when the **5** is in the centre of the dose window.
- The pen is now ready to be reset.

STEP 5 Reset the Pen



- Turn dose knob clockwise until it stops and the **5** is in the dose window.

STEP 6 Remove and Dispose of the Needle



- Remove the needle after each injection.
- Carefully put the outer needle shield back over the needle.



- Unscrew the needle.
- Replace blue cap on pen before storage.



- Throw away needles in a puncture-resistant container or as recommended by your healthcare professional.

STEP 7 Store Pen for Next Dose

- Store your BYETTA pen properly. For more information, see About the BYETTA pen section.
- When it is time for your next routine dose, go to ROUTINE USE Section, Step 1, and repeat Steps 1-7.

NOTE: The instructions to use both the 5 µg and 10 µg pens are the same.