

GraFit User's Guide

Version 7

The screenshot displays the GraFit software interface. The main window is titled "GraFit - Cover" and features a menu bar with "Home", "Text", "View", and "Graph". Below the menu bar is a toolbar with various icons for file operations (Import Data, Add, Store Stock Item), editing (Paste, Undo, Redo, Cut, Copy, Clear), and graphing (Select All, Edit Text, Object). The interface is divided into several panes:

- Items:** A vertical sidebar on the left containing icons for "Read Me", "Notes", "Data", "Graph", "Non-linear Fit", "Enzyme Kinetics", and "Lineweaver Burk".
- Data:** A central data table with columns for "Type", "Column Name", and two data columns labeled "1" and "2".
- Graph:** A window titled "Enzyme Kinetics Data" showing a plot of "Rate" versus "[Substrate]". The plot includes a red curve representing the Michaelis-Menten equation and an inset Lineweaver-Burk plot showing a linear relationship between $1/\text{Rate}$ and $1/[\text{Substrate}]$.
- Notes:** A window at the bottom right containing "Instructions" for data entry.

Type	1	2
Column Name	[Substrate]	Rate
1	1.0000	3.0000
2	2.0000	4.0000
3	3.0000	5.0000
4	4.0000	6.0000
5	5.0000	6.3000
6	6.0000	6.8000
7	7.0000	7.0000
8		
9		
10		
11		
12		
13		
14		
15		
16		
17		
18		
19		

Parameter	Value	Std. Error
Vmax	9.5250	0.4458
Km	2.5114	0.3092

Instructions

1. Enter substrate, rate data into the data table
2. Watch...

**Erithacus
Software**



**Erithacus
Software**



GraFit User's Guide

Version 7

© Erithacus Software Limited, 1989–2011

Program and Manual Designed and Written by

Robin J. Leatherbarrow

Copyright

Copyright © Erithacus Software Limited, 1989–2011. All rights reserved.

Portions of the software described in this document copyright © Microsoft Corporation 1985–2011. All rights reserved.

Information in this document is subject to change without notice and does not represent a commitment on the part of Erithacus Software Limited. The software is furnished under a license agreement, and the software may be used or copied only in accordance with the terms of the agreement. The purchaser may make one copy of the software for backup purposes. No part of this publication may be reproduced, transmitted, stored in a retrieval system, or translated into any language in any form or by any means, including photocopying and recording for any purpose other than the purchaser's personal use without the written permission of Erithacus Software Limited.

Erithacus Software Limited
PO Box 274, Horley
Surrey, RH6 9YJ, UK

VERSION 7

Last revision January 2011

Trademarks

GraFit is a trademark of Erithacus Software Limited. Microsoft, Windows, NT are trademarks or registered trademarks of Microsoft Corporation.

Other brands and products are trademarks or registered trademarks of their respective holders and should be noted as such.

Citing GraFit in Publications

When citing GraFit, simply treat this software publication as you would a book. For example:

Leatherbarrow, R.J. (2010) *GraFit Version 7*, Erithacus Software Ltd., Horley, U.K.

Contents

Introduction.....	13
Welcome.....	13
Key Features of GraFit	13
About this Guide.....	14
Notational Conventions	14
Contacting Erithacus Software	15
Technical Support	15
Internet	15
Installation and Activation	17
System Requirements.....	17
The GraFit Package	17
Installing GraFit.....	17
Online Installation	17
CD-ROM Installation.....	17
Installation Process	18
Activating the GraFit Program	18
Activation Dialog Box	18
Running in Evaluation Mode.....	19
Purchasing a Full Activation Code	20
Deactivating GraFit	21
Deactivation	21
Changing Computers	21
Checking for Software Updates.....	21
Uninstalling GraFit	21
Different GraFit Editions	22
Academic Licenses	22
Comparison of GraFit Editions	23
Single User versus Multiple User Licensing of GraFit	24
Activating GraFit for Single User or Multiple User Licensing	24
Multiple User Licensing Operation	24
Initial Activation	24
Setting up the License Control Panel	25
Adding a Multi-user License to your Administration Account	26
Managing a Multi-user License	27
Adding Extra Users	28
Reusing Multi-user Seats.....	28
Expiry of Multi-user Licenses	29
Frequently asked Questions About Activation and Licensing.....	31
Single User Licenses	31

Multiple User Licenses	32
GraFit Basics	35
Running GraFit	35
Ribbon or Menu and Toolbars	35
Office 2007-style Ribbon Bar	35
Legacy-style Menu and Toolbars	36
Swapping Between Ribbons and Menu/Toolbars	37
The GraFit Workspace	38
GraFit Items	39
Managing Items	39
Item Representation	39
Editing the Item Name	40
Changing the Item Icon	40
Item-specific Commands.....	40
Printing the Item Contents	40
Page Settings.....	41
Stock Items	42
GraFit Templates	42
Accessing Help	42
Window Display Styles.....	43
Working with Files.....	45
GraFit Workspace Files	45
Opening an Existing File	45
Opening a Recently-Used File	47
Saving a File.....	48
Loading and Saving Individual GraFit Items	49
Creating a New File	49
Importing Data Files	49
GraFit Template Files.....	49
Template Directories.....	50
Data Tables.....	53
The Data Table Item	53
Data Types	55
Ignoring Individual Data Values	56
Moving the Data Entry Position.....	57
Selecting Regions of the Data Table	58
Selecting Part of the Data Table.....	58
Selecting Entire Columns	59
Copying and Pasting Data	59
Copying and Pasting Regions of the Data Table.....	59
Copying and Pasting Columns	59
Removing Data	61

Clearing a Region of Data	61
Deleting a Region of Data	62
Inserting Values into the Table	62
Sorting Data	63
Data Table Appearance	63
Colors and Fonts	63
Number Format.....	63
Magnification	65
Importing and Exporting Data	66
Importing Data	66
Importing Data as a New Data Table	66
Exporting Data	66
Column Types	67
(General) Column Type	69
Data Types.....	69
Column Combinations.....	71
Statistics	71
Data Fitting.....	75
Calculated Values	78
Data Processing.....	82
Locking a Data Column	87
Non-linear Fit Item	89
Overview of Non-linear Fitting	89
Status Bar Indicators	89
Settings Pane	89
Equation	90
Weighting	90
Data Table	90
Data Layout	90
X Data Column	92
Y Data Column.....	93
Error Column	93
Constants Pane	93
Linking Values to Data Table Entries	93
Estimatrix Pane	94
Linking Values to Data Table Entries	95
Interactive Estimator.....	95
Fitting Criteria Pane	96
Fitting Criteria	96
Results Pane	97
Number Format Pane	99
Fonts and Indents Panes.....	100
Formatting Item Windows	101

Formatting Options	101
Fonts Pane	101
Indents Pane	102
Linear Fit Item	103
Overview of Linear Fitting	103
Status Bar Indicators	103
Settings Pane	104
Equation	104
Weighting	104
Data Table	104
Data Layout	104
X Data Column	105
Y Data Column	105
Error Column	105
Results Pane	105
Fonts and Indents Panes	107
Equation Item	109
Overview of Equation Definitions	109
Status Bar Indicators	109
General Pane	110
Changing Font Styles	110
Inserting OLE Objects	110
Editing OLE Objects	111
X Variables Pane	111
Naming Rules	112
Number of X Variables	113
Parameters Pane	113
Constants Pane	114
Definition Pane	115
Equation Definition Syntax	116
Multi-line Equation Definitions	119
Use of Temporary Variables	119
Conditional Statements	120
Comments	120
Looping Statements	121
Goto Statement	121
Estimates Pane	121
The Need for Starting Estimates	121
Method Used to Obtain Initial Estimates	121
The Estimates Pane	123
Compiling the Equation	124
Transformation Item	125

Overview of Transformation Definitions	125
Status Bar Indicators	126
General Pane	126
Data Points Pane.....	126
Data Transformation Definitions	127
Axis Name Settings.....	128
Curve Pane.....	128
Transformed Curve Definition.....	129
Fonts and Indents Panes.....	129
Compiling the Transformation.....	129
Notepad Item.....	131
Uses of Notepad	131
The Notepad Item Window	131
Alternatives to the use of Notepad Items	131
Creating Graphs	133
How to Create a New Graph.....	133
Available Graph Types	133
X/Y Scatter Graph.....	134
Column Chart	135
Bar Chart	136
Histogram.....	136
Graph Item.....	139
Overview of Graph Items.....	139
Controlling the Page View.....	139
Adding Drawing Objects.....	140
Magnification	141
Manipulating Drawing Objects	142
Selecting Objects.....	142
Moving Objects	142
Sizing Objects	143
Common Actions.....	144
Grouping Objects	145
Aligning Objects	145
Setting the Z-Order	146
Setting the Object Color.....	146
Setting Text Properties.....	147
General Object Properties	150
Fill Properties	150
Line Properties	151
Rotation Properties.....	151
Location Properties	152
Arrows Properties	152

Font Properties	153
Paragraph Properties	154
Exporting Graphs	154
X/Y Scatter Graphs	157
The XY Graph Wizard	157
Select Data Source Page	157
Set Initial Scaling	159
Select Tick Styles	160
Select Data Point Style	161
Enter Curve Parameters	162
Enter Curve Constants	162
Select Fitting Curve Style	163
Select Join Style	164
Modifying Graphs	165
Selectable Regions of a Graph	166
Main Graph Region	166
General Page	167
Fill Page	167
Line Page	167
X Ticks Page	168
X Scale Page	170
Y Ticks Page	171
Y Scale Page	172
Location Page	172
Axis Scale Regions	172
Scale Page	172
Axis Text Regions	174
Legend Box Region	174
Legend Page	175
Fill Page	176
Line Page	176
Location Page	176
Data/Curve Properties	176
Data Points Page	176
Join Points Page	177
Fitted Curve Page	178
Data Source Page	179
Data Source Page	180
Error Bars Page	181
Axis Association Page	183
Equation Page	184
Parameters Page	185
Constants Page	185
Legend Text Regions	186

Adding Extra Data to the Graph.....	186
Set Axis Associations Page	186
Using Multiple Graph Axes.....	189
When to Use Multiple Axis Scaling	189
Overlaying Data that have Different Scales	189
Split Axes	189
Creating and Managing Multiple Axes.....	190
Running the Axis Manager	190
Adding a New Axis.....	191
Splitting an Existing Axis.....	192
Removing an Axis	193
Working with Multiple Axes	194
Axis Text and Axis Scale Text.....	194
Data Association.....	194
Axis Scaling and Sizing.....	195
Column and Bar Charts	199
The Chart Wizard	199
Select Data Source Page.....	199
Set Initial Scaling	201
Select Tick Styles	202
Select Bar Style.....	203
Select Data Point Style	204
Select Join Style.....	205
Histograms.....	207
The Histogram Wizard	207
Select Data Source Page.....	207
Set Initial Scaling	208
Select Tick Styles	209
Select Bar Style.....	210
Embedding and Linking	211
What is Embedding and Linking?.....	211
Original Cut and Paste Model	211
Embedded Objects	211
Linked Objects.....	212
Using Embedding and Linking.....	212
Comparison of Data Transfer Types.....	213
Embedding and Linking without using the Clipboard	214
Embedding and Linking with GraFit.....	214
GraFit as an OLE container	214
GraFit as an OLE source.....	215

Data Fitting Basics	217
When Is Data Fitting Appropriate?	217
Types of Data	217
Regression Analysis.....	219
Practical Considerations.....	219
Types of Regression Analysis	220
Linear Regression.....	220
Polynomial Regression	220
Non-linear Regression	220
Additional Considerations for Non-linear Regression	220
Weighting Scheme	221
Linear Regression of Transformed Data	221
Information Obtained from Data Fitting	222
Multiple Regression	222
Uses for Multiple Regression	223
Advantages of Multiple Regression	224
Disadvantages of Multiple Regression	224
Defining Equations	224
Fitting Data by Multiple Non-linear Regression	224
Data Fitting.....	225
Plotting Results from Multiple Regression.....	226
 Theory	 227
Regression Analysis.....	227
Assumptions in Regression Analysis	228
Linear Regression.....	230
Weighting.....	231
Weighting Types	231
Non-linear Regression	232
Standard Errors of the Calculated Parameters	233
Robust Weighting	233
Determining the Most Appropriate Weighting	234
Plot of Residual Errors.....	235
Testing Goodness of Fit — F Test	235
Nested Models	236
Comparison of Equations with Equal Numbers of Parameters.....	236
Mean and Standard Deviation	236
 Worked Data Fitting Examples.....	 237
Using One of the Equations Provided	237
Fitting the Data	238
Checking the Weighting.....	242
Drawing Residual Plots.....	242
Using Robust Weighting	245
Defining a New Equation (1).....	246

Defining the Oxygen Binding Equation	247
Defining a New Equation (2)	249
Defining the Equation	250
Fitting the Equation.....	250
Multiple Regression	251
Displaying Multidimensional Data	252
Reading Values from a Calibration Curve	254
Choosing an Equation to Fit the Data	254
Reading Values from the Fitted Curve	256
Fitting Several Data Sets with Some Parameters Common and Some Specific	256
Method Used to Share Variables between Data Sets	257
Example — Paired Denaturation Curves	259
Enzyme Inhibition Samples	260
Competitive Inhibition Fitting	260
Production of a Dixon Plot	261
Comprehensive Inhibition Modeling.....	262
Equations Provided	265
2 Substrates	266
Enzyme Inhibition	267
Enzyme Kinetics	270
Exponential	271
IC 50	273
Ligand Binding	274
pH	278
Rate.....	281
Standard Templates	285
About the Standard Templates.....	285
Accessing Templates	285
Template Layout	286
Using Your Own Templates	286
Description of the Standard Templates	287
Enzyme Kinetics	287
IC 50	288
IC 50 Back Corrected	289
Ligand Binding.....	290
Linear Regression	291
Exponential	292
Exponential Offset.....	293
Exponential Decay Various Models.....	294
Ligand Binding Various Models	295
Passing Data between Applications	297
Beyond Embedding and Linking	297

Security Warnings When Using OLE Automation	297
Example: Passing Data to and from Microsoft Excel.....	298
Transferring Data using the "GraFit Fitter" Excel Template	300
Using "GraFit Fitter" to Fit the Example File	300
"GraFit Fitter" Reference	302
Transferring Data using Visual Basic (Advanced!)	309
How is Data Transfer Done?	312
Sample Data Files Installed	313
GraFit Automation Reference (Advanced!)	315
Using GraFit Automation with Microsoft Excel	315
GraFit Automation Routines	317
GraFit Helper Routines.....	320
Hard-coded Equations	325
Overview of Hard-coded Equations.....	325
What are Hard-coded Equations?.....	325
Enabling Hard-coded Equations.....	325
Inserting a Hard-coded Equation Item.....	326
General Pane	327
Locating the Equation Library File.....	327
X Variables Pane	327
Parameters Pane	328
Constants Pane	329
GraFit Symbols Font	331
GraFit Symbols Font.....	331
Index	333

Introduction

Welcome

GraFit version 7 for Microsoft® Windows™ combines presentation-quality scientific graphs with comprehensive data fitting. The result is a powerful package that is easy to use, produces stunning output, and can also fulfill the needs for quantitative data analysis. Data can be fitted using various types of theoretical models, and your own equations can easily be added to those supplied with the program. Comprehensive graph plotting facilities allow you to display your data on-screen and edit the graphs interactively.

Key Features of GraFit

GraFit provides the experimental scientist with a tool for visualization and analysis of their data. For the occasional user, GraFit can be configured to perform routine analyses with minimal intervention—in many cases it is possible just to enter the data and allow GraFit to calculate the results and plot the graphs in the background. Advanced users can utilize the full power of GraFit for analyzing complex equations, and may create custom fitting templates for routine data processing.

About this Guide

This manual describes how you can use GraFit to examine your data, and to plot graphs. We recommend that you read the manual fully, and try out some of the examples described.

Notational Conventions

The following notational conventions are used throughout this manual:

Convention	Meaning
SMALL CAPITALS	Small capitals indicate keys to be pressed. For instance, A indicates "press the 'A' key". If more than one key needs to be pressed this is written as (for instance) SHIFT END which would mean "hold down the 'Shift' key and then press the 'End' key".
Monospace	Monospace type indicates text that should be entered on the computer keyboard. For example, to obtain a directory of drive A: when running a Command Prompt window, you would enter DIR A: RETURN
Bold	Bold text is used for menu item commands, e.g. the Edit menu. Where cascading menus are present, an arrowhead (▶) is used to separate the sub-menus. For example, the Windows Control Panel application is started from the desktop Start ▶ Control Panel command.

Note *Important information that should be noted is written in italics, and separated from other text by horizontal rules, as here.*

Contacting Erithacus Software

Technical Support

At Erithacus Software we are working hard to provide you with intuitive technical products. However, if you find yourself with a technical question that you cannot resolve, please contact our technical support department. Registered users can obtain technical support from one of the following sources:

- Fax: [+44] (0)1342 841 939
- Email: support@erithacus.com
- Internet: <http://www.erithacus.com/grafit>

Note Always quote your serial number and the product version, which are shown on the About screen (**Help ▶ About GraFit** command).

Internet

Our web site can be found at <http://www.erithacus.com> and here you will find latest information about our software products. In addition, we post technical notes, tips and answers to common problems.

Tip Be sure to check out <http://www.erithacus.com/grafit> to download the latest fixes and updates. GraFit will check automatically for updates when the program starts. You can also use the **Help ▶ Check for Updates** command to ensure you have the latest release.

Installation and Activation

System Requirements

Software

GraFit version 7 requires Windows XP or later and is fully compatible with Windows Vista and Windows 7.

Hardware

Any computer that is capable of running the Windows operating system in use.

The GraFit Package

The GraFit package comprises the following:

- GraFit program on CD-ROM or via electronic download.
- The GraFit User's Guide (this manual). Depending upon the GraFit package that was purchased, this guide may only be available in electronic form.

Note When GraFit is first run you will be able to register the product. Registration is optional, but will ensure that you are informed of the latest updates.

Installing GraFit

Online Installation

Download and run the installation program from the GraFit website (www.erithacus.com/grafit).

CD-ROM Installation

Insert the GraFit CD-ROM into your CD or DVD drive. The Welcome screen will load automatically.

Note Some systems are configured so that automatic running of CD-ROM software is disabled. If your system is set up in this way, the Welcome screen can be displayed by selecting the **Run** command from the Windows **Start** menu, and entering the command `D: \CDRUN` where `D` is the drive letter of your CD-ROM.

Installation Process

GraFit

Select the **Install GraFit** option to install GraFit onto your computer. The installation program will copy the files to your hard disk, and place an entry for the GraFit program on the Windows **Start** menu, under the **Programs** entry.

Adobe Acrobat Reader


The GraFit installation includes this user manual saved in Adobe Portable Document Format (PDF). To read this file, it is necessary to have the Adobe Acrobat Reader (version 6 or later) installed. If you do not already have Acrobat installed, you can do so by selecting **Install Acrobat Reader**. Installation of Acrobat Reader requires your computer to be connected to Internet.

Activating the GraFit Program

After installing GraFit it is possible to run the software in evaluation mode for 14 days. During this time the program is fully functional and is equivalent to a single-user activation. To continue running the program after the evaluation period expires, it is necessary to activate the program online by entering a valid activation code (serial number). Activation codes are purchased; after the program has been activated it will run without showing the activation dialog box again.


Activation Dialog Box

Before GraFit has been activated, running the program will produce the following Software Activation dialog box:

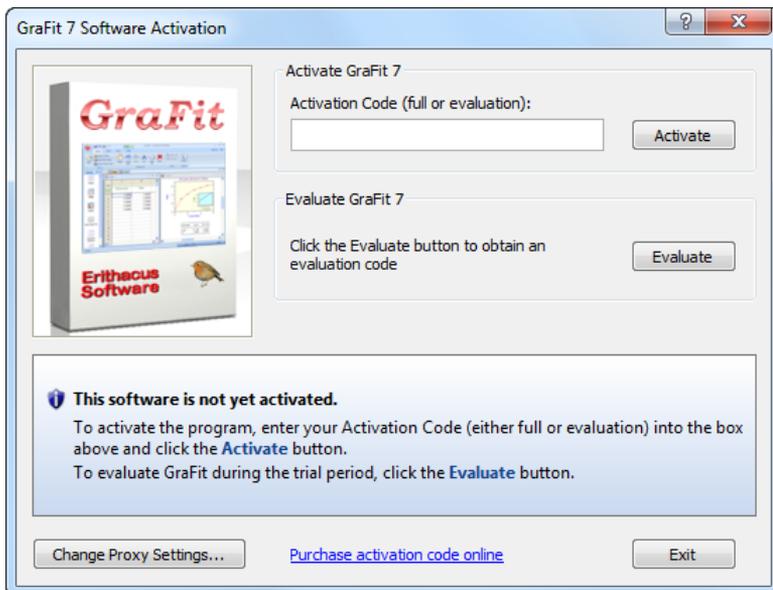

Software Activation dialog box

While the program is being evaluated this dialog box shows the remaining evaluation period; if you are still within the 14 day evaluation limit then you can run the program by clicking the **Evaluate** button. Once the evaluation period has expired you will not be able to run the software until you have entered a full activation code.

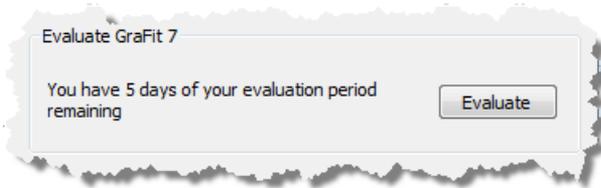


*During the evaluation period the number of days remaining is shown. Click **Evaluate** to run the program.*

Setting Proxy Settings

For some users, access to the internet needs to be via a Proxy Server. Providing you can access the internet using a browser such as Microsoft Internet Explorer, GraFit will automatically make use of any proxy settings that might be needed. In some rare circumstances it may be necessary to use custom proxy settings in order to access the internet; these can be set by the **Change Proxy Settings** button.

Running in Evaluation Mode

The software can be evaluated for 14 days prior to purchase. To start your evaluation period, click the **Evaluate** button on the Software Activation dialog box. In order to start the evaluation your computer must be able to access the internet. The first time you evaluate the software, you will need to supply information about where the Evaluation code should be delivered. This code will be emailed to you at the address entered into the following dialog box:

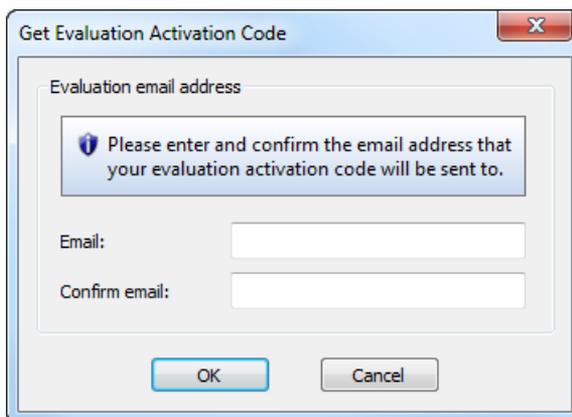

When you first start the evaluation period you must enter the email address to which the evaluation activation code will be sent

Click **OK** to request an evaluation activation code. This code will be sent immediately, but you should allow a few minutes for your email system to process the message.


Note Most email systems are set up to move unsolicited email to a "Junk" email folder. In the event that you cannot find your evaluation code email, try looking in your Junk folder in case it has been moved there.

Using the Evaluation Activation Code

When you receive the email with your code, paste this into the **Activation Code** box and click the **Activate** button.

During activation you will be asked to provide a password, which will be needed in future if you reactivate or deactivate the software and an email address that is used to allow automatic reminders to be sent if you forget your password.

During the evaluation period the software is fully functional. After the evaluation period has ended it will no longer be possible to run the software unless you purchase and enter an evaluation code.

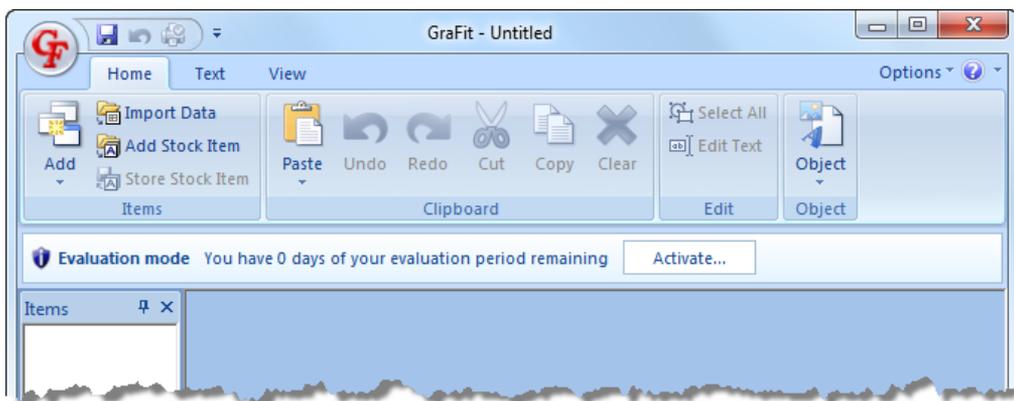

GraFit running in evaluation mode

Purchasing a Full Activation Code

Full Activation codes can be purchased online by following the link on the Activation dialog box. There are two distinct types of full activation code: single-user activation codes and multi-user activation codes. More information about multi-user codes can be found later in this Chapter.

Activating the Software

Once you have purchased a full activation code, paste this into the **Activation Code** box and click the **Activate** button. If you purchased GraFit on CD-ROM the activation code (serial number) will be present on the software case; if you have downloaded the software you may purchase an activation code online.

During activation you will be asked to provide a password, which will be needed in future if you reactivate or deactivate the software. You will also be asked to provide an email address. This address is used to allow automatic reminders to be sent if you forget your password.

Activation is done via the internet; if your computer is not connected to the internet it is possible to carry out offline activation. See the online help for a full description of how offline activation is carried out.

Deactivating GraFit

Your GraFit software needs to be activated on-line, during which time a count is made of the number of times an individual activation code is used. You will not be allowed to activate software on more computers that is allowed by your license agreement. This means that if you wish to move your software to a different computer you should first deactivate it on the original computer.

Note Once you have deactivated your software it will no longer run until it is once more reactivated.

Deactivation

To deactivate your previously-activated software, select the **Help About** dialog and select the **Deactivate** button. You will need to enter the password that you gave when you originally activated the product. If you have forgotten your password, clicking the **Forgotten password?** button will automatically email this to you

Multi-user Installations

If you have a multi-user installation then there is no **Forgotten password?** button as passwords are administered by the designated License Administrator (see later). The License Administrator is also able to deactivate any multi-user installation remotely.

Changing Computers

If you change computer, it is advisable to deactivate your GraFit installation from the old computer before reinstalling the program on the new computer.

Checking for Software Updates

You can check to see whether any program updates are available by running the **Help ► Check for Updates** command. If updates are available online you have the opportunity of downloading and installing them.

Uninstalling GraFit

GraFit can be removed from your system by the following procedure.

1. From the Windows **Start** menu, launch the Control Panel application.
2. Double-click the **Add or Remove Programs** icon.
3. Choose GraFit 7 to uninstall, and click the **Change/Remove** button.
4. Follow the instructions


Note We strongly recommend that you deactivate GraFit, as described above, prior to uninstalling the program.

Different GraFit Editions

GraFit comes in a variety of editions that share key functionality but have different licensing and activation models. The different editions are listed below.

Edition	Features	License Validity
GraFit Pro	Multi-user version with web-based management of activation codes.	Rolling 1-year subscription, extendable
GraFit Pro (Academic)	As above for use in non-profit organizations only.	Rolling 1-year subscription, extendable
GraFit	Single-user stand-alone version of GraFit	Unlimited
GraFit (Academic)	As above for use in non-profit organizations only	Unlimited
GraFit (Student)	Entry-level single-user student edition, for use in non-profit organizations only.	1-year subscription, single-user license

You can tell which edition of GraFit you have from the appearance of the program title bar.

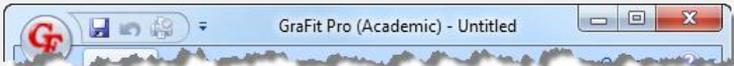Edition	Appearance
GraFit Pro	 The title bar for GraFit Pro shows the GraFit logo on the left, followed by a menu icon, a help icon, and a dropdown arrow. The text in the title bar is "GraFit Pro - Untitled".
GraFit Pro (Academic)	 The title bar for GraFit Pro (Academic) shows the GraFit logo on the left, followed by a menu icon, a help icon, and a dropdown arrow. The text in the title bar is "GraFit Pro (Academic) - Untitled".
GraFit	 The title bar for GraFit shows the GraFit logo on the left, followed by a menu icon, a help icon, and a dropdown arrow. The text in the title bar is "GraFit - Untitled".
GraFit (Academic)	 The title bar for GraFit (Academic) shows the GraFit logo on the left, followed by a menu icon, a help icon, and a dropdown arrow. The text in the title bar is "GraFit (Academic) - Untitled".
GraFit (Student)	 The title bar for GraFit (Student) shows the GraFit logo on the left, followed by a menu icon, a help icon, and a dropdown arrow. The text in the title bar is "GraFit (Student) - Untitled".

Academic Licenses

If you work in a not-for-profit organization then you can purchase an “Academic” GraFit license at a reduced price. Academic licenses are restricted in that they cannot be used in a commercial organization; they are available for both single-user and multi-user installations.

Student Licenses

The GraFit (Student) edition is a low-cost entry version that is specifically designed to support use by students. This edition is restricted to use only in not-for-profit organizations. The Student edition has a number of in-built licensing restrictions: it can be activated only on a single computer at any one time and the license expires after 1 year. Continued use beyond that time will require purchase of a new license code.

Comparison of GraFit Editions

A comparison of the different GraFit editions is given below.

	Single User Licensing	Student Licensing	Multiple User Licensing
Editions	GraFit, GraFit (Academic)	GraFit (Student)	GraFit Pro, GraFit Pro (Academic)
Number of users per license purchased	1 user per license	1 user per license	Depends on license size
Installations allowed	A single user may install the program on a maximum of two separate computers ¹	Installation on a single computer only	The number of installations allowed is equal to the number of user licenses purchased
Cost per seat		Lowest entry price	Lower cost per seat than single user licensing
Remote license management	No	No	Yes
Remote deletion of licenses	No	No	Yes
Installation password	A separate password is used for each user	A separate password is used for each user	A single password is used to install all users ²
Duration of license	Unlimited	Expires after 1 year	Yearly renewal of subscription required
License can be renewed	(Not applicable)	No	Yes
Free upgrade to later versions	No	No	Yes
Internet connection required for activation	Yes	Yes	Yes
Internet connection needed to run program	No	No	Yes ³
Serial numbers	70xxxxxxxxxxxxxxxxx, 7Axxxxxxxxxxxxxxxxx	7Sxxxxxxxxxxxxxxxxx	MUxxxxxxxxxxxxxxxxx, MAxxxxxxxxxxxxxxxxx

¹ These installations must only be to support the designated user.

² The first person to install the program provides an email address for someone who becomes the License Administrator. This person is responsible for managing passwords etc.

³ An internet connection is required periodically to verify the license validity and test for expiry and remote deletion.

Single User versus Multiple User Licensing of GraFit

There are two different licensing models for GraFit: single-user and multi-user. The licensing model does not affect the functionality of the program, but there are some differences regarding purchase, installation and software maintenance.

Single User Licensing

Single user licensing is the traditional way of purchasing software: a separate software license is purchased for each user of the program. Activation of the program requires a separate activation code (serial number) for each installation.

Multiple User Licensing

Multiple user licenses are advantageous in larger workgroups or organizations. A single software license is purchased which allows a far larger number of installations to be performed using the same activation code. In addition, there is a centralized mechanism to monitor and regulate the software use.

Activating GraFit for Single User or Multiple User Licensing

GraFit will automatically adopt the appropriate licensing model depending on the activation code that is entered. GraFit activation codes are always 18 characters long, with the first two letters indicating the type of license code. If the first character is 'M' then this is a multi-user code; if the first character is a number then this is a single user code that is specific for a particular software version. For example, 7xxxxxxxxxxxxxxxxx is a single-user code for version 7 of GraFit whereas Mxxxxxxxxxxxxxxxxx is a multi-user code that will activate any release of GraFit from version 7 onwards.

Purchasing Single User or Multiple User Licenses

Our online software shop allows you to purchase either single-user licenses or multi-user licenses. Please visit www.erithacus.com/grafit for more information.

Converting between Single User and Multiple User Modes

Although it is unlikely that you will need to do this frequently, it is possible to switch modes (though you will need to have purchased appropriate licenses in order to allow this). To change mode:

1. Deactivate your current license by clicking on the **Deactivate** button on the **Help ▶ About** dialog box.
2. Exit the program
3. Restart GraFit
4. Enter the new activation code and click the **Activate** button

Multiple User Licensing Operation

Initial Activation

When you first activate a copy of GraFit using a multi-user activation code, you will see the following Activation Information dialog box.

Activation Information dialog box for multi-user activation codes: first activation

The activation password that is entered will be used to activate all further seats associated with this license code and the email address entered will be that of the License Administrator who will be able to monitor license use via the License Control Panel (LCP, see below).

Note Activation of the first copy of GraFit will use up one of the available license seats. If the License Administrator is not a regular GraFit user (for example, in a large corporation this could be an IT professional) then it is possible to regain use of the seat by deactivating this user installation via the License Control Panel. Deactivating an individual user in the License Control Panel will not remove their registration as the License Administrator.

Setting up the License Control Panel


After the initial activation, the License Administrator will receive an email (with the subject “License Control Panel Information”) telling them how to create an online account that allows them to monitor license use. Our multi-user licensing is implemented using software produced by a company called softWORKZ, who provide the online licensing servers. This email will remind you of the multi-user activation code that has been used and provide a Confirmation Key that will let you add this code to your personal License Control Panel

account. A single account can monitor any number of difference activation codes, allowing one or more GraFit licenses to be regulated by a common administrator if required.

Note Most email systems are set up to move unsolicited email to a "Junk" email folder. In the event that you cannot find your License Control Panel Information email, try looking in your Junk folder in case it has been moved there.

Registering for a New License Control Panel Account

Follow the link in the email and click on the **New Account** button.

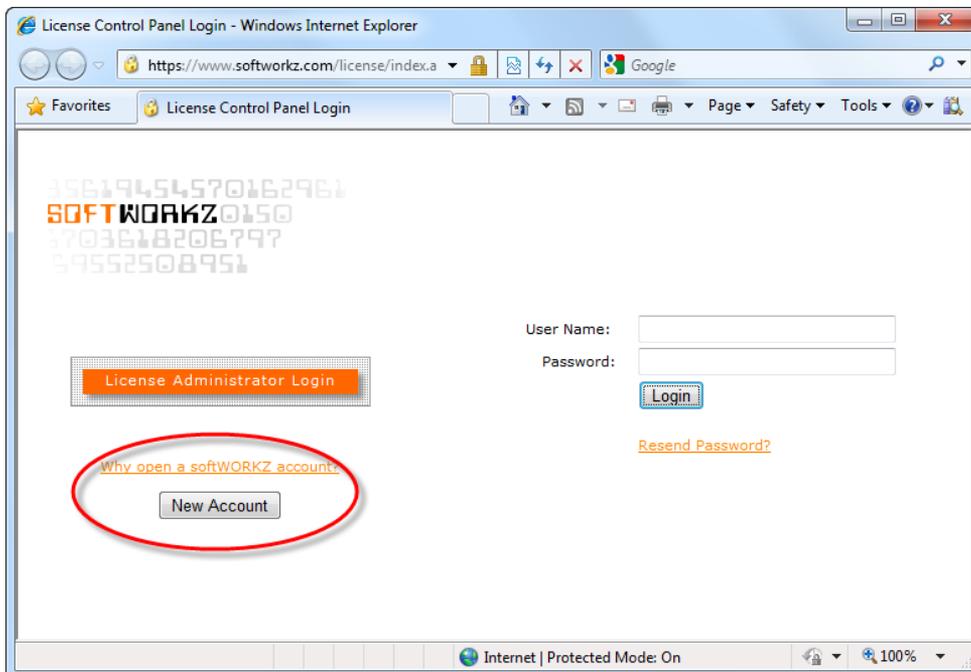

The License Control Panel login screen

Follow the instructions to create your account. Once created you can then log in the your personal account

Adding a Multi-user License to your Administration Account

When logged into you LCP account, click **Add New Software** to add a multi-user license to your account.

Add New software to your Account

To add new software to your account, you will need the **Confirmation Key** that was sent to you by email. This email is sent to you after you activate your software for the first time.

Please 'Copy and Paste' the Confirmation Key from the email here:

Enter Confirmation Key

Adding new software to your License Control Panel account

You will have received the Confirmation Key in the initial email; enter this and click **Submit**. Once added you will be able to administer each multi-user license that you register by clicking on the entry in the list.

Your Software Licenses

single user software - click on code to view details and manage

Software Description	Company	Activation Code	Activation Date	Activation email	Activated Computers	Maximum Computers
multiple seat software - click on code to view details and manage						
Software Description	Company	Activation Code	Activation Date	Activation email	Activated Seats	Maximum Seats
GraFit Pro Academic	Erithacus Software	MA[REDACTED]	12/28/2009	[REDACTED]@erithacus.com	0	10
GraFit Pro Academic	Erithacus Software	MA[REDACTED]	12/28/2009	[REDACTED]@erithacus.com	1	10

[Help](#)

List of all software licenses administered via an individual License Control Panel account

Managing a Multi-user License

To manage an individual multi-user license, clicking on the entry in the list above. This will produce a screen that shows information about each activated user.

Software details for Activation Code MA[REDACTED]

Software	From Company	Activation Date	Activation Email	License Type	Expiry Date
GraFit Pro Academic	Erithacus Software	12/28/2009	[REDACTED]@erithacus.com	Multi User - 10 Seats	12/28/2010

Computer ID	Seat Activation Date	De-activate ?	Computer ID	Seat Activation Date	De-activate ?
ROBIN-ADVENT	12/28/2009	Yes			

Management of an individual multi-user license

All activated seats are listed by Computer ID, which is the name of the individual computer that has activated the software.

Changing the Activation Password and Administrator Email Address


To change the Administrator Email address, click the **Change Email** button and to update the password click **Change Password**. If you change the password then all new activations will need the updated password. Changing the password requires knowledge of the current password—if you are unsure of this click **Resend Password** to email this to the current Administrator Email address.

Deleting a User from the License

The License Administrator can deactivate the software on any of the installed computers by clicking on the **Yes** link under the “De-activate” entry. Once deactivated the user of that computer will no longer be able to run GraFit unless they re-enter a valid activation code and password. Deactivation will typically take 24 hours to take effect.

Adding Extra Users

Additional users, up to the maximum specified by the license, are added in the same way as the initial user. The GraFit software must be installed on the user's computer and then when the software is run, the multi-user license code is entered into the activation information dialog box. For subsequent users this will produce the following dialog box.

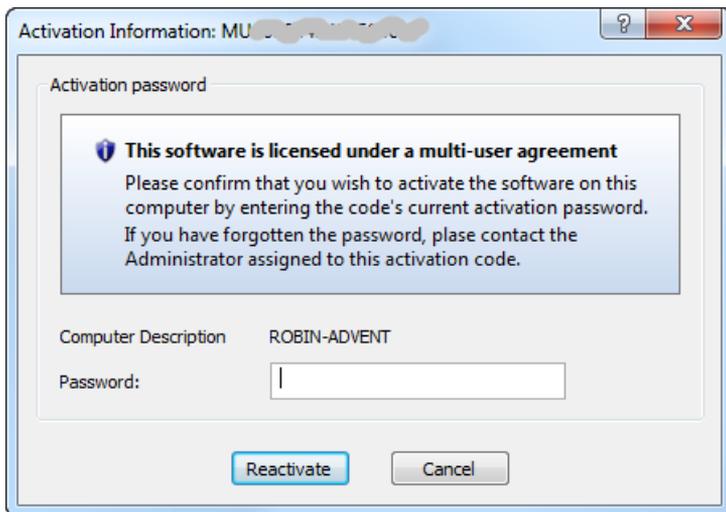

Activation Information dialog box for multi-user activation codes: second and subsequent activations

Activation of the software will need the password that was specified during the first activation and which is available from the License Administrator who first activated the software.

Reusing Multi-user Seats

The seats in a multi-user license can be re-allocated to other users by deactivating one or more of the installed seats. Deactivation can either be performed by the License

Administrator from the LCP, or by the using the **Deactivate** button that is found on the **Help ► About** dialog box.

Note Deactivation of a seat will typically take 24 hours to process. You will therefore not be able to add a new user seat until this time has elapsed.

Expiry of Multi-user Licenses

Multi-user licenses are sold on a time-limited basis and expire one year after the initial activation. Extension of the license needs to be performed before the expiry date in order to ensure all users are able to continue using the software.

When does a License Expire?

The length of time remaining on the license can be found on the **Help ► About** dialog box.



The expiry of a multi-user license is shown on the Help ► About dialog box


Three weeks before your multi-user license is due to expire, all users of the license will see a warning banner whenever GraFit is run.



Multi-user license expiry warning

What Happens if my Multi-user License Expires?

If a multi-user license expires then the next time a user of code tries to run GraFit they will see the following dialog box.

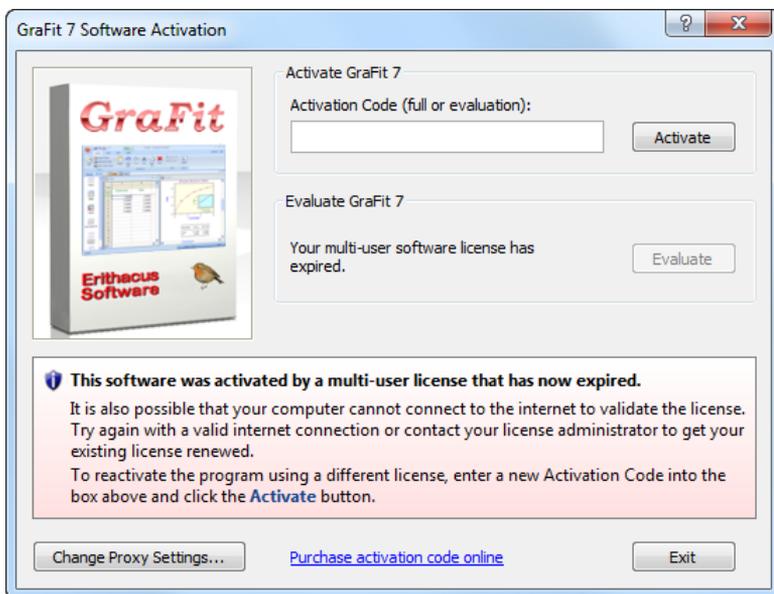

Software Activation dialog box displayed when a multi-user license has expired

This dialog box is displayed either if the multi-user license has expired. To continue working with GraFit either renew the existing multi-user license, at which point all users of this license will then be able to use the software again, or purchase a new single-user or multi-user license and reactivate the program using the new license.

Note If your computer is not connected to the internet for an extended period then this will prevent license verification, which will also result in apparent expiry of the license. Reconnect to the internet and try again in order to rectify this.

How are Licenses Renewed?

Multi-user licenses can be “topped-up” at any time from the License Control Panel. This involves the following process:

1. A new multi-user license is purchased, typically from our online store (see www.erithacus.com/grafit).
2. The activation code will be then be received as normal by email.
3. Instead of using this license to activate GraFit as shown above, it is used to add more time to an existing multi-user installation, using the License Control Panel.

From the License Control Panel, choose the existing activation code that you wish to top-up then click the **Extend License** button.



Extending the expiry date for an individual multi-user license

Enter the code for the new multi-user license that you have just received. The time available on the new license will be added to that of the existing one.

*Note After extending your license it will typically take 24 hours before the new expiry date is shown correctly on the **Help ► About** dialog box. However, the License Control Panel information will be updated immediately.*

License extension can be performed at any time. An activation code can only be used to top-up another code if the top-up code has not already been activated, and can only be used once for this purpose.

Note It is only possible to top-up a multi-user license using one of the same type. For example, a 10-user license can only be topped up with another 10-user license. Please ensure that you purchase the correct license type when you wish to extend your license.

Frequently asked Questions About Activation and Licensing


Single User Licenses

What happens to my installation if my computer is damaged or replaced?

Where possible you should explicitly deactivate your GraFit installation on any old computer, then reinstall on the new one. However, there are circumstances where this is not possible, typically if your computer breaks or is stolen. Under these circumstances you simply need to reinstall GraFit on the new computer using your activation code and password. If you have forgotten your password then this can be resent to you by clicking the **Forgotten Password?** button on the installation dialog box. In the event that the new activation exceeds the permitted activation count you will still be able to proceed though you will also have to reactivate any GraFit installations on other computers (see “*What happens if I try to activate more copies than allowed by my license?*”).

What happens if I try to activate more copies than allowed by my license?

If you try to activate GraFit on an additional computer that results in the installation count for this activation code being greater than permitted by the license then the latest activation will be allowed, but all previous installations will automatically be deactivated.

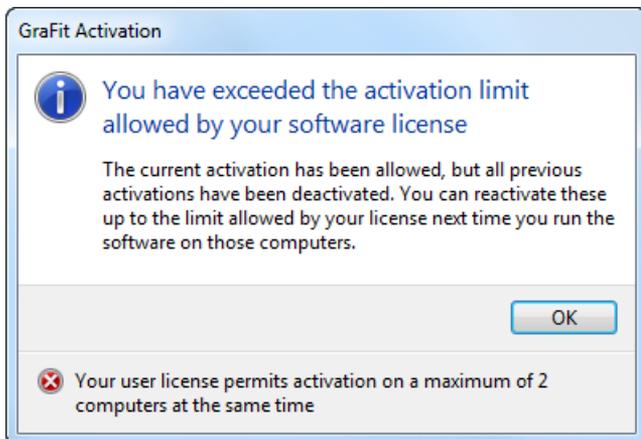

Message shown when activation of GraFit exceeds the license limit for the code being used

You can reactivate the previous installations, up to the permitted license limit, next time you run GraFit on those computers.

What happens if I have lost my activation code?

It is your responsibility to keep your activation code safe. In the event that you lose this code contact Erithacus Sales via our web site and we will do our best to assist you.

Multiple User Licenses

Note You should also consult your License Administrator when dealing with multi-user installations

What happens if I try to activate more copies than allowed by my license?

You will not be able to activate more computers than allowed by the license for the activation code you are using. Please see your License Administrator if you need to add more seats.

What happens to my installation if my computer is damaged or replaced?

Each activation is associated with an individual computer. In the event that computer is no longer running a copy of GraFit then the License Administrator can remotely delete that seat from the license, allowing you to reinstall GraFit on your new computer without affecting the number of available seats remaining on the license. Please see your License Administrator if you need to remove an unwanted installation remotely.

What happens if I have lost my activation code?

Activation codes and passwords are managed by your designated License Administrator. Please see them for assistance.

I am the License Administrator and have lost my administration details!

Please contact Erithacus Sales via our web site and we will do our best to assist you.

My license has expired or is about to expire, what can I do to extend it?

Management of activation codes is done by your License Administrator. They can extend the expiry time by purchasing a new license code and using this to top-up the time remaining on your original multi-user code.

Why is the license time-limited? I don't want my software to expire!

Our multi-user licensing model is based upon defined expiry periods coupled with web-based management of installed seats. It is not restricted to a particular version of GraFit and if a new version is released then your multi-user license will also activate the new release. In addition, the pricing of multi-user licenses factors in the reduced lifetime, making this an affordable and managed method of software use. However, this licensing model is not suitable for all users and if you would prefer not to operate with time-limited software then you should instead purchase one or more a single user licenses.

GraFit Basics


Running GraFit

The GraFit installation program will add a GraFit entry to the Programs list on your Windows **Start** menu. When you run the GraFit program, the main GraFit window will open.

Ribbon or Menu and Toolbars

GraFit 7 uses by default an Office 2007-style ribbon bar to access the various commands. The ribbon bar combines the functionality of both toolbars and menu bars, which are used in more traditional Windows programs.

Office 2007-style Ribbon Bar

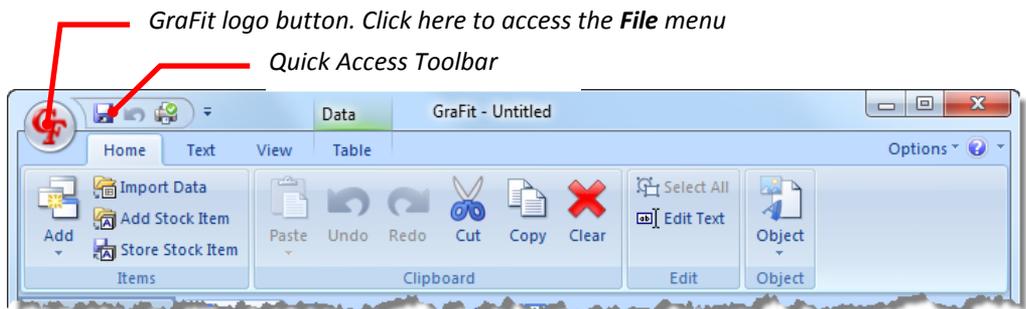


GraFit 7 using the Office 2007-style ribbon bar

The ribbon bar displays commands grouped on a series of tabs, giving a less cluttered layout than the traditional display, together with a “Quick Access Toolbar” that can be customized to hold commonly used commands in a consistent manner. The **File** menu is accessed via the GraFit logo button, found in the top left of the GraFit window.

Using the Ribbon Bar

The ribbon bar comprises a set of tabs that hold related commands. Each tab has one or more ribbon groups that collect toolbar buttons. Click on the button to activate the command.

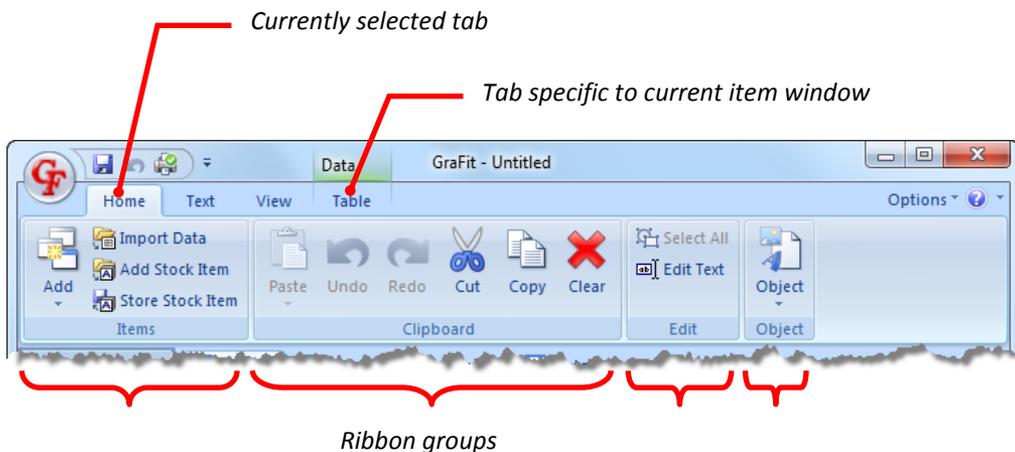

The **Home** tab holds frequently used commands; the **Text** tab has commands that are used for text formatting; the **View** tab has commands that relate to the way information is displayed. These three tabs are always present, but in addition there are some tabs that hold commands specific to particular item windows. When a data table windows is selected (as above), a **Table** tab is present on the ribbon bar; when a graph item window is selected a **Graph** tab will be found.

Using the Quick Access Toolbar


The Quick Access Toolbar can be used to hold any commands that you frequently access. In the example below, various formatting commands have been added.



The **Customize** button allows various aspects of the toolbar to be specified, including whether the toolbar is drawn above (as shown) or below the ribbon bar. To add commands to this toolbar, either use the **Customize** button or right-click on any command button that is found on the ribbon bar and select **Add to Quick Access Toolbar**.

Legacy-style Menu and Toolbars

It is also possible to run GraFit using a legacy-style combination of menu and toolbars, giving an appearance that is similar to GraFit version 5 and earlier. The legacy style may be preferred by existing users of the program.

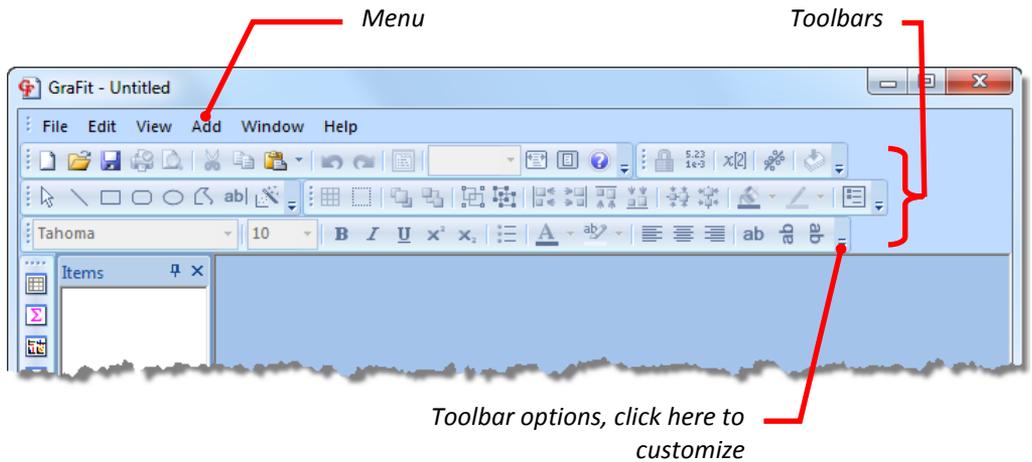


GraFit 7 using legacy-style menu and toolbars

Customizing Menus and Toolbars


To customize menus and toolbars, click the **Toolbar Options** button that is found on the right of each toolbar. You can completely customize the commands that are found on each toolbar and menu by dragging and dropping the various command button icons.

Swapping Between Ribbons and Menu/Toolbars

The **Appearance** command is used to change the layout of the program. The location of this command is found on the **Options** menu of the ribbon bar, or the **View** menu when in menu and toolbar mode.

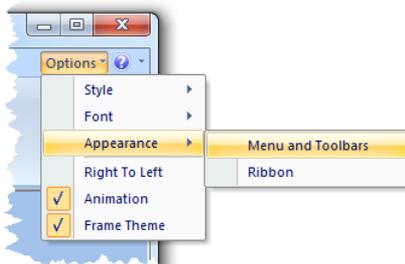


Changing from ribbon to menu and toolbars

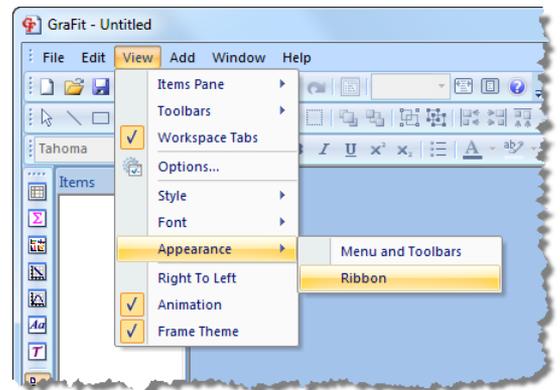

Changing from menu and toolbars to ribbon

The GraFit Workspace

When initially opened, the GraFit workspace is empty. A GraFit file can comprise one or more data tables, graphs, data analyses, equation definitions, transformation definitions, notepads and embedded objects. Collectively, these items make up the GraFit workspace, which are stored to disk as a single file. After adding a number of items to the workspace, the appearance of the GraFit window is shown below.

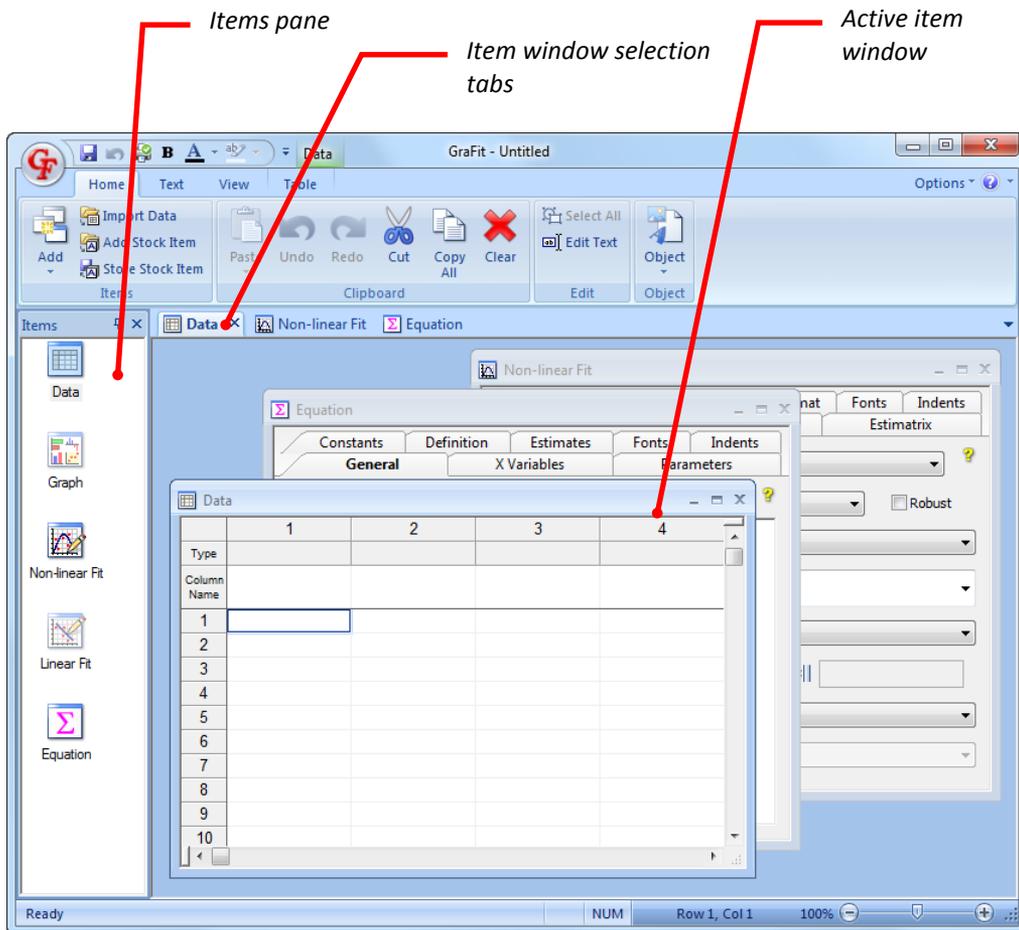

A GraFit workspace with a number of items added

In this example, one each of the following items has been added: a data table; a graph window; a non-linear fitting item; a linear fitting item; an equation definition. The linear fit and graph items have been closed so that their windows are not visible in the workspace; their icons within the Items pane are drawn lighter to reflect this. All the added items are stored inside the same GraFit workspace, and can be saved and loaded from disk using the **File** menu commands. If you run a second copy of GraFit, you can copy items between

workspaces by drag and drop or by using the commands in the **Clipboard** group⁴ to copy and paste items via the Windows Clipboard. Each item in the workspace is represented by an icon in the Items pane.

GraFit Items

GraFit items are initially created in a new window that resides within the main GraFit window, and are identified by the presence of an icon in the Items pane. Closing an item window does not delete the item — it simply makes the window invisible. This allows you to keep open only those items that you are currently working on.

Managing Items

Adding New Items

To add a new item to the workspace:

- Use the **Add** command on the **Home** tab⁵ and select the type of item that you require.

Switching Between Items

To switch between items:

- Click on the window to activate, or click on the icon that represents this item in the project window, or use the item window selection tabs.

Removing Items

To delete an item from the workspace:

- Select the icon that represents the item in the project window, and choose the **Clear** command⁶.

Item Representation

When an item window is closed, the contents of the item still remain present in the GraFit workspace file. The item window can be reopened by double-clicking on the icon, which is found in the Items pane. The icon that represents the item is shaded so that it indicates the current state of the item, as shown below.



Linear Fit, open

GraFit item, open. The icon is drawn normally.



Linear Fit,
closed

GraFit item, closed. The icon is drawn as a ghosted image, which is lighter than normal.

⁴ Found on the **Edit** menu in legacy mode

⁵ The **Add** menu in legacy mode

⁶ The **Edit** ► **Clear** command in legacy mode



Embedded object, open for editing. The object icon is shown overdrawn with hatched lines.



Embedded object, closed. The icon is drawn ghosted.



Linked object, closed. A shortcut overlay (☒) is added to the icon.

Switch between item windows to edit their contents. Information about using each of the GraFit item types is found in Chapters 5-13.

Editing the Item Name

Items are referred to by the name shown in the project window, so it is useful to change this from the default to something more meaningful. To edit the item name, select the icon in the project window, choose the **Edit ▶ Rename** command and type in the new name.

Item names can contain any characters, but must be different from any other item name in the workspace.

Changing the Item Icon

The icon that represents an item may also be changed. This can be useful to allow similar items to be visually distinguished. To change the icon, select it and choose the **Icon** command⁷ (on the **View** tab of the ribbon). For each item type there are a number of alternative images that can be used. In addition, it is possible to use an icon image that resides in any executable or DLL file on your computer.


Item-specific Commands

Although some commands apply to all GraFit items, many are item-specific. For example, the **Fit Data** command applies to the Linear fit and Non-linear fit items, but not to the Data table or Graph items. Commands that are not available are shown grayed. In addition, each item window has one or more specific menus that are only present when this window is selected. For example, only graph items have **Graph** and **Object** menus.

Printing the Item Contents

The contents of any item window can be printed using the **File ▶ Print** command, and the printing previewed using the **File ▶ Print ▶ Preview** command.


⁷ Found on the **Window** menu in legacy mode

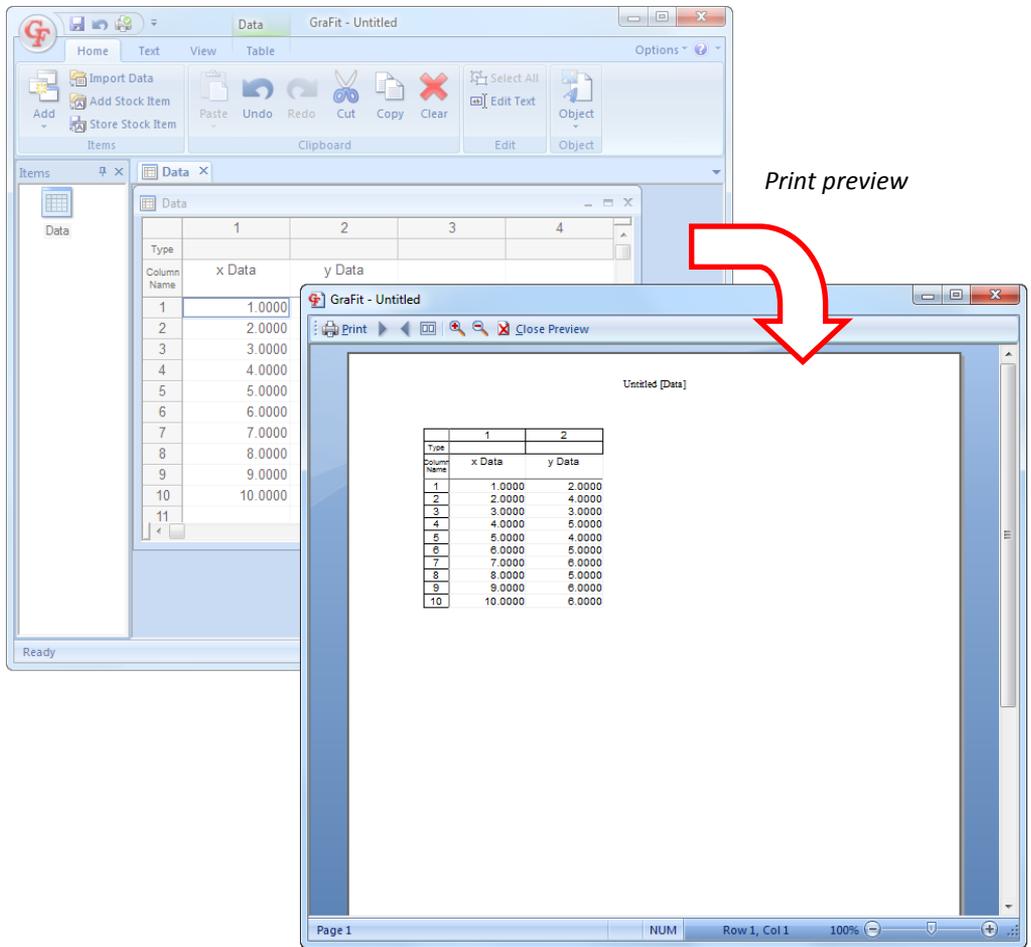

Page Settings

The **File ▶ Print ▶ Page Setup** command⁸ allows the paper size and orientation to be specified, and the margin sizes to be entered. These settings are stored separately for each item in the workspace.

⁸ **File ▶ Page Setup** command in legacy mode

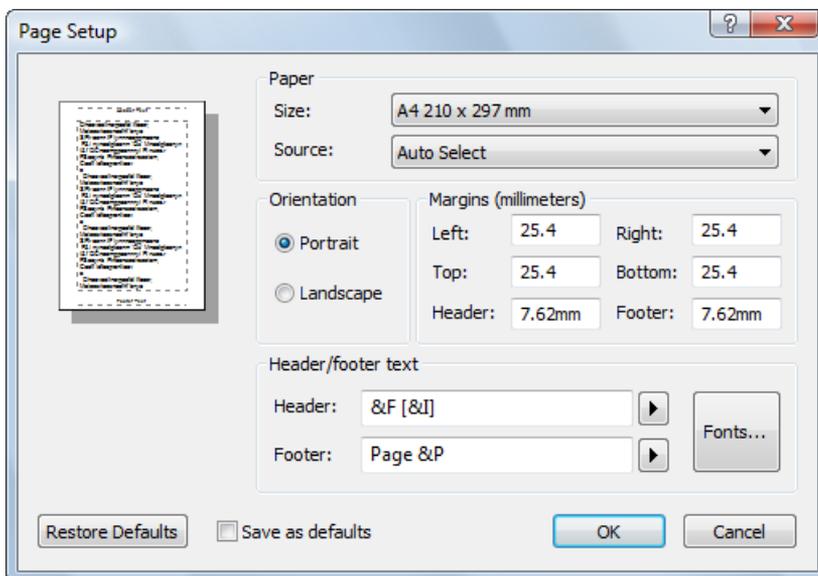

A header and footer are printed at the top and bottom of the page, respectively. The header and footer text can contain special strings such as the file name or the time that the item was printed. These special strings are specified by clicking on the  button. To omit the header and/or footer from the output, leave the entries for the **Header text** and/or **Footer text** blank. The font used to display the header and footer text can be selected by clicking on the **Fonts** button.

Changing the Printer

To change the printer from the default, click on the **Printer** button.

Stock Items

Individual GraFit items can be stored for future reuse by selecting an item and using the **Store Stock Item** command. They are reloaded using the **Add Stock Item** command⁹. Storing and loading of items is especially useful for saving and restoring equation definitions, but can be used for any GraFit item. See Chapter 4 for more information about the location of item files.

GraFit Templates

Workspace files can be created by adding a number of items as described above. Alternatively, it is possible to load a pre-built workspace from a template, using the **File ▶ New** command. See Chapter 4 for more information about using templates.

Accessing Help

GraFit includes a complete on-line help system to guide you through using the program.

⁹ In legacy mode these commands are **File ▶ Store Stock Item** and **File ▶ Add Stock Item** respectively.

To view the help topics

- From the **Help** menu, select the **Contents** command.



Help menu

To obtain specific help for a dialog box


- Click on the Context Help button in the top right of the dialog box to display the help that is available for that dialog box.




Context Help menu

Window Display Styles

By default, GraFit uses a similar window style to that of the Microsoft Office 2007 suite of programs. However, a variety of alternative window styles can be used if desired. To change the window style, use the **Options ► Style** command¹⁰. The various style colors are shown below.

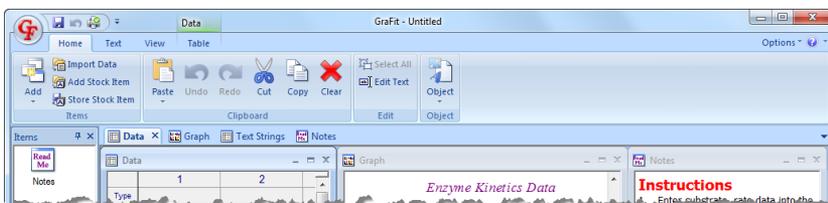

Blue (Office 2007 style)

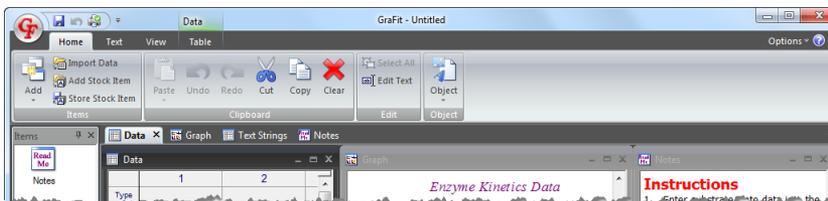


Black




Silver


¹⁰ **View ► Style** command in legacy mode.

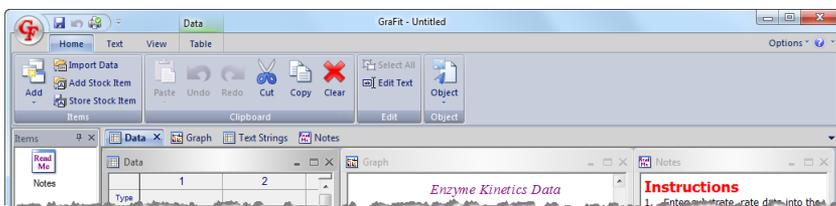

Aqua

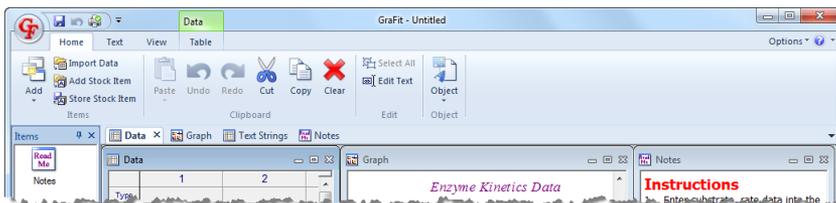

*Scenic
(Windows 7
style)*

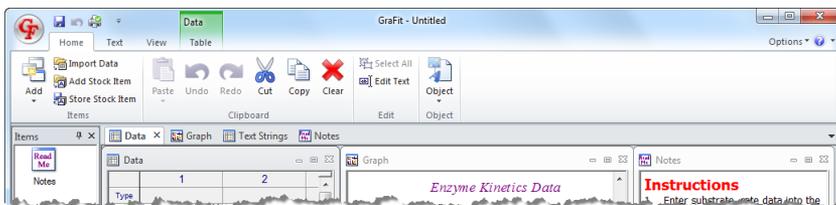

*White (Office
2010 style)*

Working with Files

GraFit Workspace Files


The GraFit workspace consists of all the items displayed in the main GraFit window (data tables, graphs etc.). The contents of the workspace are stored together in a GraFit workspace file, which has the .GFF extension. In addition to the contents of the workspace, the current window sizes and positions are also stored in the workspace file.

Opening an Existing File

Choose the **File ▶ Open** command to open an existing GraFit workspace file. This will bring up the Open dialog box.

Windows Vista


When running GraFit on Windows Vista, the program uses the Vista-style Open dialog box, as shown below.

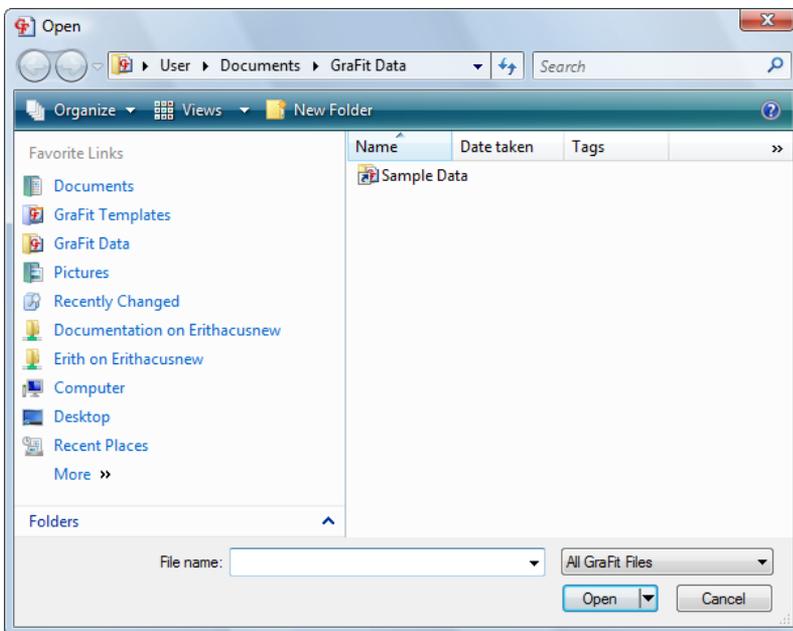

The locations of your **GraFit Data** and **GraFit Templates** directories are added to your *Favorite Links* entries, allowing you to navigate quickly to these. (The locations of these special directories are chosen from the **Options** command, on the **View** tab¹¹.)

Windows 7


With Windows 7 the Open dialog box includes the shortcuts to your **GraFit Data** and **GraFit Templates** directories under a *GraFit 7.0* entry. (The locations of these special directories are chosen from the **Options** command, on the **View** tab¹¹.)

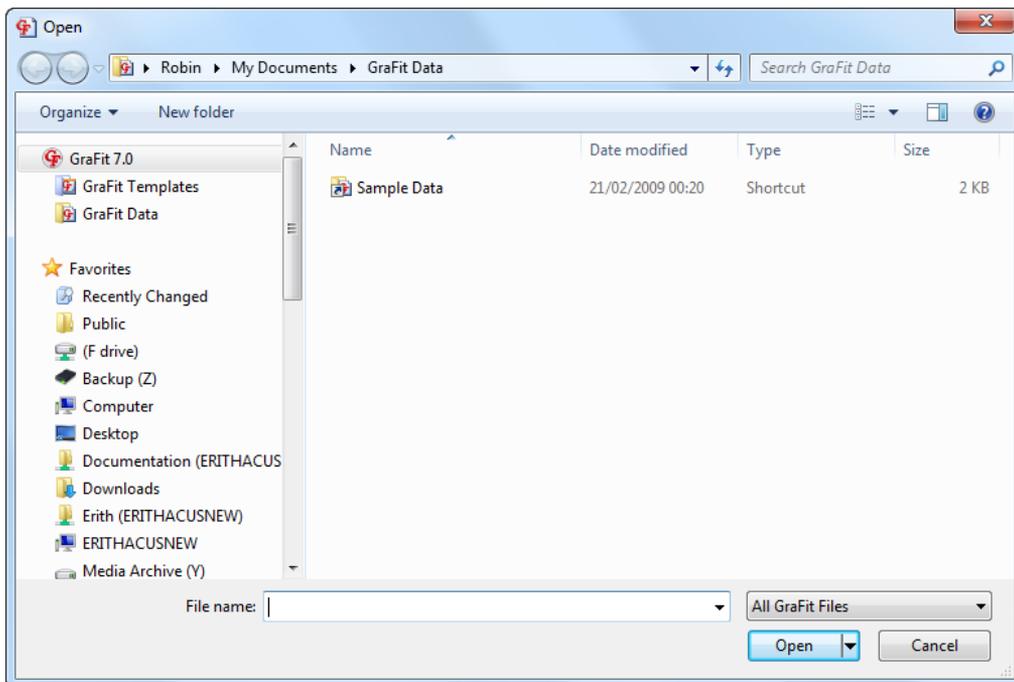

Windows 2000 and Windows XP

Earlier versions of Windows use the older-style file open dialog box.


¹¹ The **View > Options** menu in legacy mode

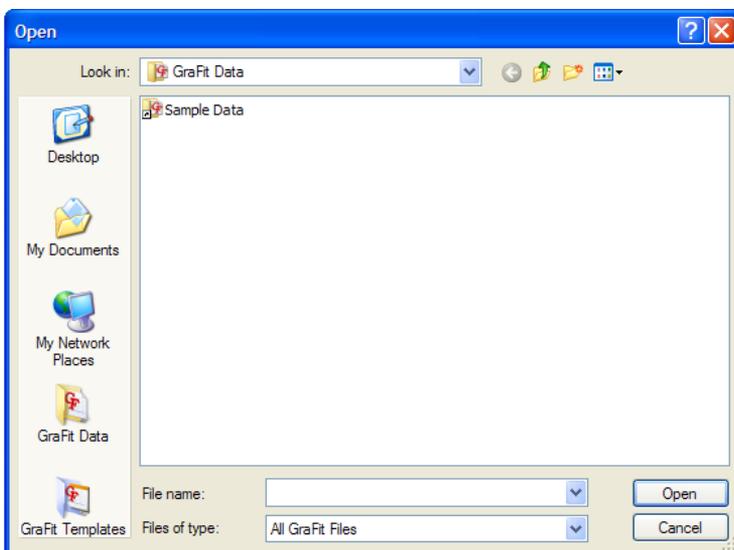

The buttons marked **GraFit Data** and **GraFit Templates** allow you to rapidly to the default data and user template directories respectively. (The locations of these special directories are chosen from the **Options** command, on the **View** tab¹¹.)

Opening a Recently-Used File

GraFit maintains a list of your most recently-used documents. To see this list, open the **File** menu by clicking on the GraFit logo button.

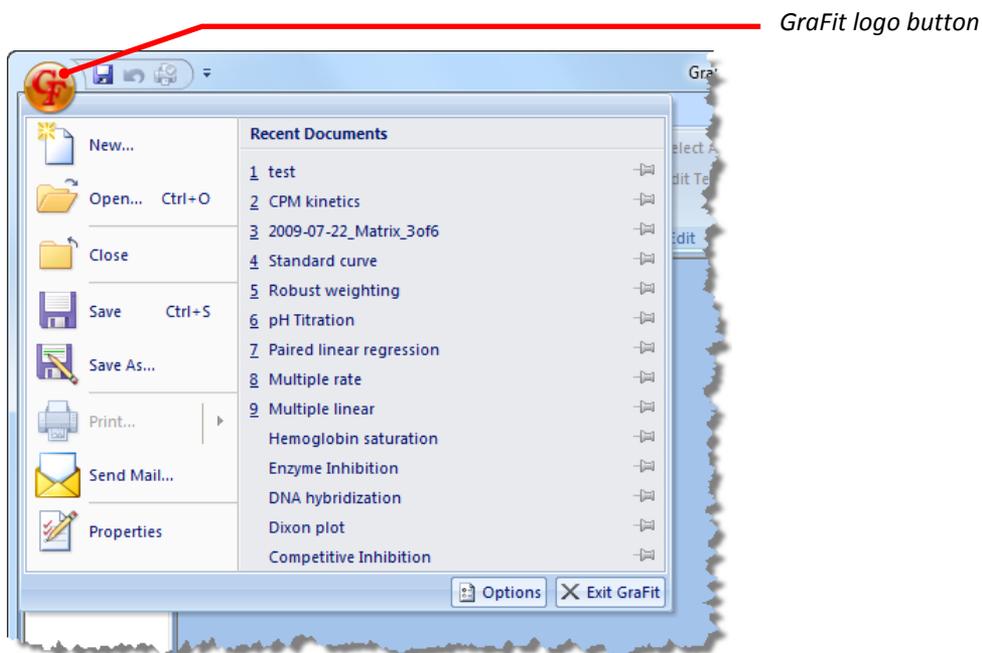

Select a document from the list to re-open it.

Pinning files to the list

The list of recent documents shows files that have been most recently opened by the program. If you wish to ensure that a particular document is kept on the list, so that it will always be available, you can pin any file to the list by clicking on the push-pin icon.

 *Push-pin icon un-pinned*

 *Push-pin icon pinned*


Click on the icon to toggle between un-pinned and pinned.

Saving a File

GraFit workspace files are saved using the **File** ► **Save** or **File** ► **Save As** commands. The contents of the workspace are saved, along with the current window positions and order. When the file is reloaded, these positions and order are restored.

Saving a File to Use With Older Versions of GraFit

Newer versions of GraFit often have features that are lacking in earlier versions. The latest versions of the program will read all earlier data files, but older versions of the program may not be able to open files saved by the latest release. To save a file using a format that can be read by older versions of GraFit (version 5 and 6), use the **File** ► **Save As** command and set the **Save as type** to “*GraFit Workspace File (versions 5-6)*,” as shown below.

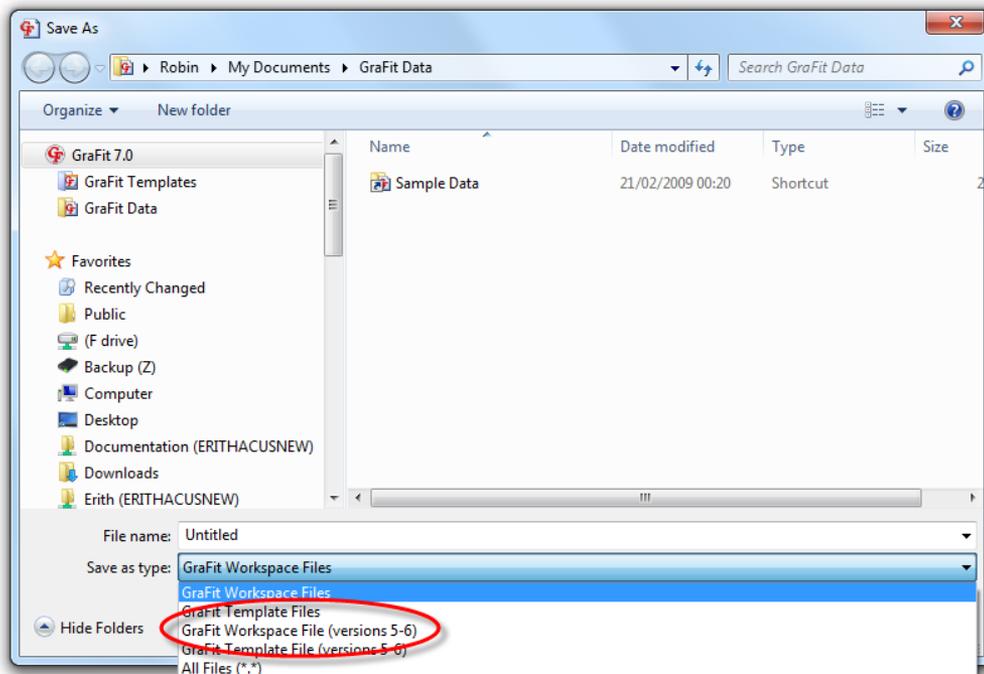

Loading and Saving Individual GraFit Items

Saving Items

The GraFit workspace comprises a number of individual items, such as data tables, graphs and analyses. These individual items are saved in one of the following ways:

- Select the item in the Items pane (or click on the item window), and then select the **Store Stock Item** command¹².

Items are stored as files that have a .GFC file extension. By default, items should be saved in the item store directory, or a subdirectory off this. (The location of the item store directory is specified using the **Options** command, on the **View** tab.)

Loading Items

Stored items can be loaded using the **Add Stock Item** command¹³.

Creating a New File

To reset the GraFit workspace:

- Select the **File ▶ Close** command¹⁴.

Importing Data Files

Data can be imported into GraFit to create a new data table item that holds the imported data. To create this item, select the **Import Data** command¹⁵. Any file for which a filter is available can be imported into GraFit (see the following section on Import/Export Filters).

It is also possible to paste data from a file directly into an open data table. To do this, use the **Paste From** command found on the Table tab¹⁶.

Import/Export Filters

In order to perform conversions between GraFit and other data types, GraFit uses installable import/export filter files. Several standard filter types are supplied by default with GraFit, and other filters may be available from certain instrument manufacturers or from Erithacus Software. These filter files, which have the .GFR file extension, are installed into GraFit from the **Options** command, on the **View** tab, Data Filters page.

Please contact Erithacus Software if you have special requirements for a specific filter file.

GraFit Template Files

Template files are pre-made workspace files that are used to speed up use of common procedures. Template files differ from workspace files in the following ways.

¹² The **File ▶ Store Stock Item** command in legacy mode

¹³ The **File ▶ Add Stock Item** command in legacy mode

¹⁴ In legacy mode, click on the **New** command button, , which is found on the Standard toolbar.

¹⁵ The **File ▶ Import Data** command in legacy mode

¹⁶ The **Edit ▶ Paste From** command in legacy mode

- They are represented by a different icon.
- Double-clicking on a template file opens a copy of the file rather than the file itself.
- Templates stored in the user or group template directories can be rapidly accessed using the **File ▶ New** command.

Creating a Template File


To store a workspace as a template file, select the **File ▶ Save As** command, and save as type “*GraFit Template File*” (to create a template file that can be read by earlier versions of GraFit, save as type “*GraFit Template File (version 5-6)*.” If you wish to access this template using the **File ▶ New** command, ensure that the file is saved into your user or group template directory (see the following section on Template Directories). To help with this, the Save As dialog box contains a button that will jump immediately to your user template directory.

Creating a File Based on a Template

The **File ▶ New** command creates an untitled file based on a selected template. This command is described in the following section. Chapter 23 describes the standard templates that come pre-defined with GraFit.

Template Directories


There are two directories that have special significance when storing template files. These are the user and group template directories, the location of which is specified from the **Options** command, on the **View** tab, Directories page:

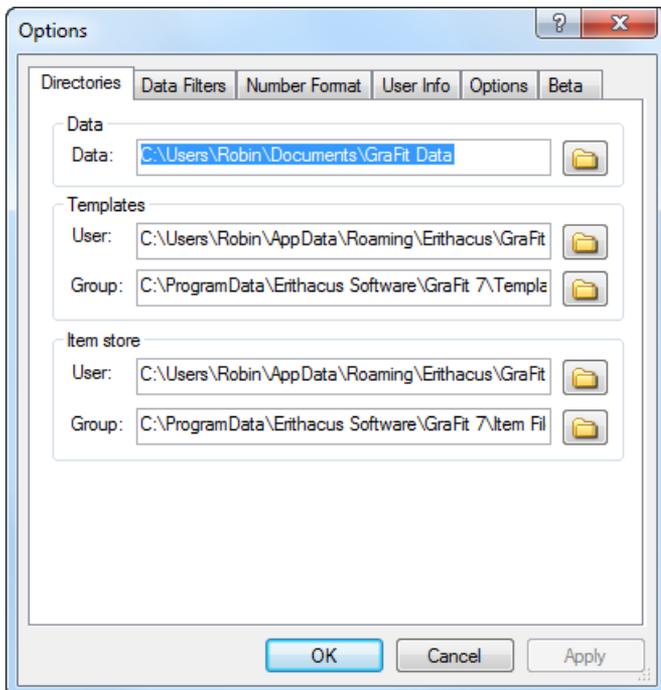

Click on the  button to change a directory setting. By having two separate locations for templates it is possible to keep user-specific files separate from general templates that are used by a group of users. It can be useful, for example, to maintain a group directory on a network server.

Files in a template directory appear on the General tab of the **File ▶ New** dialog box. Any subdirectories have special significance, and result in additional tabs on this dialog; clicking on the tab accesses the files in these subdirectories.

Note Empty subdirectories do not appear in the **File ▶ New** dialog box.

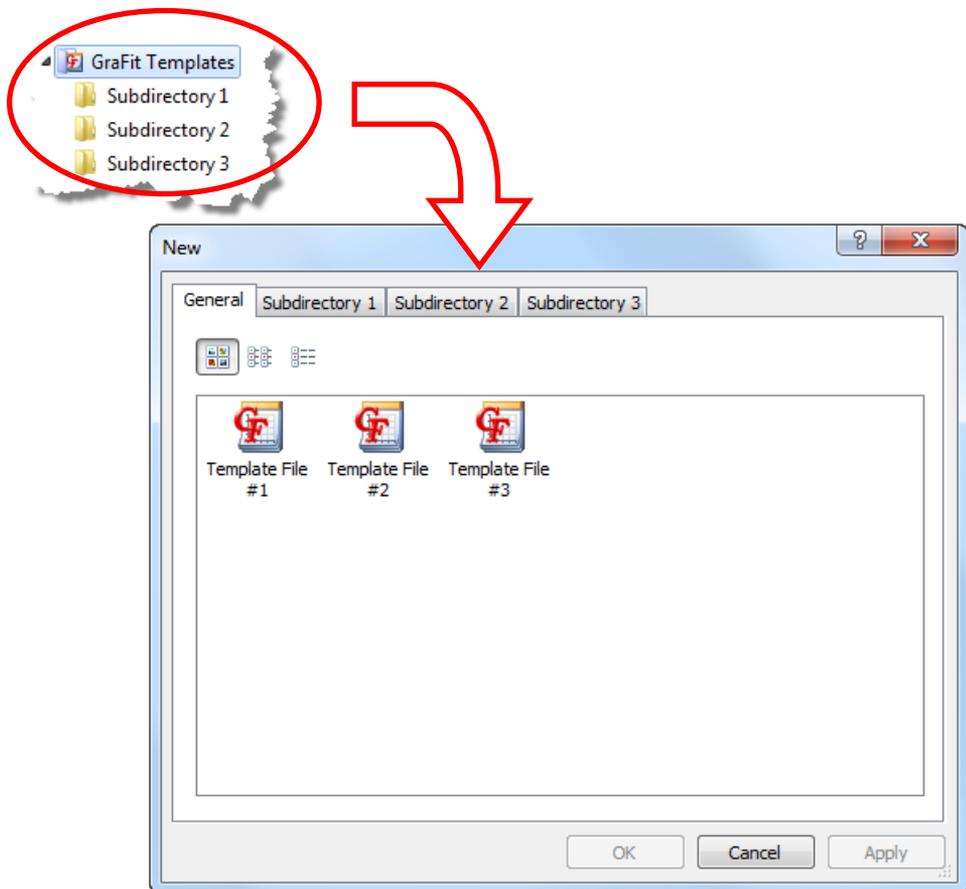

The files displayed on this dialog box are the sum of those present in the user and the group template directories. (The specific location of any individual file can be found by right-clicking on the icon, and choosing the **Properties** command. In the event that a file with the same name is found in both the user and group directories, the user file has precedence.) To create a new workspace file based on one of these templates, select the one to use and click **OK**.


Data Tables

The Data Table Item

To add a new, empty data table to the GraFit workspace:

- Select the **Add ► Data Table** command.

The parts of the data table window are shown in the following illustration.

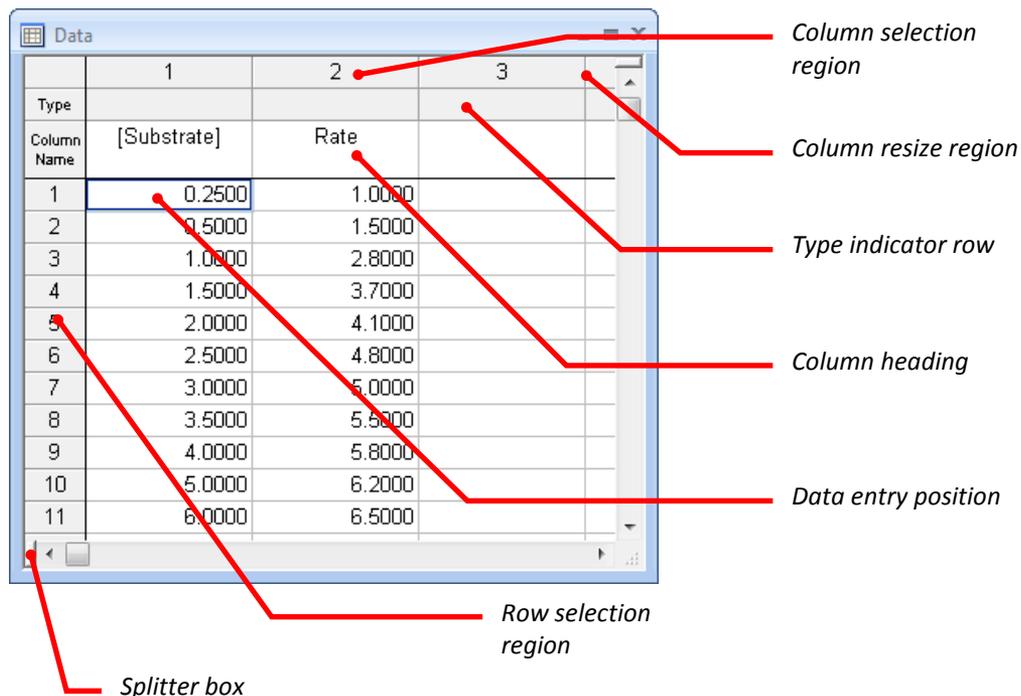

Position row 1 column 1 in the above figure is marked by a dark outline; this is the current data entry position. Data entries are edited *in situ* by typing in the new entry. When the entry is complete, pressing RETURN, TAB or a cursor key, or clicking the mouse on a different position

in the data table, will accept the data, and move the data entry position. To edit the value at the current data entry position, select the **Edit Current Cell** command¹⁷, or press F2.

The following table describes the elements in a data table. Certain sub-regions are distinguished by the mouse cursor changing shape when it is over them.

Element	Description	Cursor
Data table	The data values themselves are listed in the data table. Click the mouse on a value to move the data entry position; drag the mouse to select larger regions of the data table. Newly entered data is placed at the data entry position, which can be moved using the cursor keys or by clicking the mouse on a different part of the data table.	
Column selection region	Click and drag the mouse on the column numbers to select entire columns of data.	
Column resize region	Click and drag the edges of a column to resize it.	
Row selection region	Click and drag the mouse on the row numbers to select entire rows of data within the data table.	
Row resize region	Click and drag the edges of a row to resize it.	
Split table horizontally	Drag the splitter box on the horizontal scroll bar to divide the data table window in two. The two halves can be independently scrolled, allowing distance regions of the table to be viewed simultaneously.	
Split table vertically	Drag the splitter box on the vertical scroll bar to divide the data table window in two. The two halves can be independently scrolled, allowing distance regions of the table to be viewed simultaneously.	

Resizing Rows and Columns

The row height and column width can be set using the **Row Height** and **Column Width** commands on the Table tab¹⁸, respectively, or by using the mouse to drag the edges of the rows or columns within the row or column selection region. Columns can be sized to fit their contents using the **Size Column to Fit** command on the Table tab¹⁹.

¹⁷ The **Edit ► Edit Current Cell** command in legacy mode

¹⁸ The **Table ► Row Height** and **Table ► Column Width** commands in legacy mode

¹⁹ The **Table ► Size Column to Fit** command in legacy mode

Data Types

The number of data points that can be entered is limited only by available memory. The data are grouped into columns, and the contents of one or more columns can be used for data fitting or graph plotting.

Within a column, each individual cell holds a data value that may be one of three different types:

Data Type	Description
Numbers	<p>Numbers are the most common data form that will be entered into GraFit. They are used to plot x/y scatter plots, and for data fitting by regression analysis. Very small or large numbers can be entered in exponential format, for example $1e7$ is equal to 10000000, and $1e-23$ represents 10^{-23}.</p> <div style="border: 1px solid black; padding: 5px; margin-top: 10px;"> <p><i>Note</i> The representation of numbers, particularly the symbol used for the decimal point, differs from country to country. By default, GraFit uses the country-specific settings selected for your Windows installation (using the Control Panel application, Regional Settings). However, you can override these from the Options command²⁰ on the View tab, Number Format page.</p> </div>
Text	<p>Text entries can be used for annotation, and are also used to provide the text for the category axes of charts. Text values are treated as missing data during data fitting.</p>
Missing Data	<p>Experimental data often contain certain values that are not available (perhaps because that particular experimental sample was the one that fell on the floor!). GraFit allows you to insert “missing” values by leaving the data entry blank. All the graphing and data fitting elements of the program will recognize such points, and omit them from the analysis.</p>

The representation of these different data types is illustrated in the following diagram.

²⁰ The **View ► Options** command in legacy mode

	1	2	3	4
Type				
Column Name	This is column 1	This is column 2		
1	1.2345	Text Label		
2		5.4321		
3	1.0000	Another label		
4				
5				

Ignoring Individual Data Values

It is possible to mark one or more data entries to be ignored during data fitting procedures or any calculations performed within the data table. This can be useful if a data value is known to be (or thought to be) erroneous and you do not wish to include it in an analysis. Marking the data value to be ignored is preferable to deleting it as this does not remove the data point completely, allowing it to be reinstated later if required.

To ignore a data value, select it and use the Ignore Data command²¹ on the Table tab.

	1	2	3	4
Type				
Column Name	x Data	y Data		
1	1.0000	4.0000		
2	2.0000	3.0000		
3	3.0000	5.0000		
4	4.0000	21		
5	5.0000	6.0000		
6	6.0000	7.0000		
7	7.0000	8.0000		
8	8.0000	7.0000		
9	9.0000	9.0000		
10				

²¹ The **Table** ► **Ignore Data** command in legacy mode

Numeric values that are flagged to be ignored as treated as Missing Data during any data analysis procedures. To restore the data value, select it and repeat the Ignore Data command.

Moving the Data Entry Position

When you have entered your data value you will need to move to the next data entry position by pressing TAB, RETURN or DOWN after entering the value. The following keys control movement of the data entry position within the data table.

To	Press
Move the entry point to the right	TAB OR RIGHT
Move the entry point to the left	SHIFT TAB OR LEFT
Move the entry point up	UP
Move the entry point down	DOWN OR RETURN
Move to the first column of a row	CTRL LEFT
Move to the last column that contains data entries	CTRL RIGHT
Move to the first row of a column	CTRL UP
Move to the last row that contains data entries	CTRL DOWN

It is often convenient to press RETURN or DOWN between entering data points.

Note If you have a region of the table selected, pressing RETURN will move the data entry position left to right within the selected region. If you prefer to enter $\{x, y\}$ data pairs in sequence " $x_1, y_1, x_2, y_2, \dots$ " then it is convenient first to select the columns (see "Selecting Entire Columns") that will hold the data and then to enter the values, pressing RETURN between each data entry.

Moving the data entry position using the mouse

To move the data entry position using the mouse, move the mouse pointer to the position you require, and click the left mouse button.

Selecting Regions of the Data Table

Selecting Part of the Data Table

To select a portion of the data table, drag the mouse over the required region. Alternatively, hold down the SHIFT key, and use the cursor keys to extend the selection. The selected region is displayed shaded, as shown below.

	1	2	3
Type			
Column Name	[Substrate]	Rate	
1	0.2500	1.0000	
2	0.5000	1.5000	
3	1.0000	2.8000	
4	1.5000	3.7000	
5	2.0000	4.1000	
6	2.5000	4.8000	
7	3.0000	5.0000	
8	3.5000	5.5000	

Clicking the mouse within the row selection area results in the entire row being selected.

To make a more extensive selection:


- Click on the start position.
- Use the scroll bars to move the window as necessary.
- Hold down the SHIFT key, then click on the end position for the selection.

Discontinuous selection

- Hold down the CTRL key and click the required entries.

Selecting Entire Columns

Whole columns of data are selected by clicking on the column selection region. The selected column is displayed inverted.

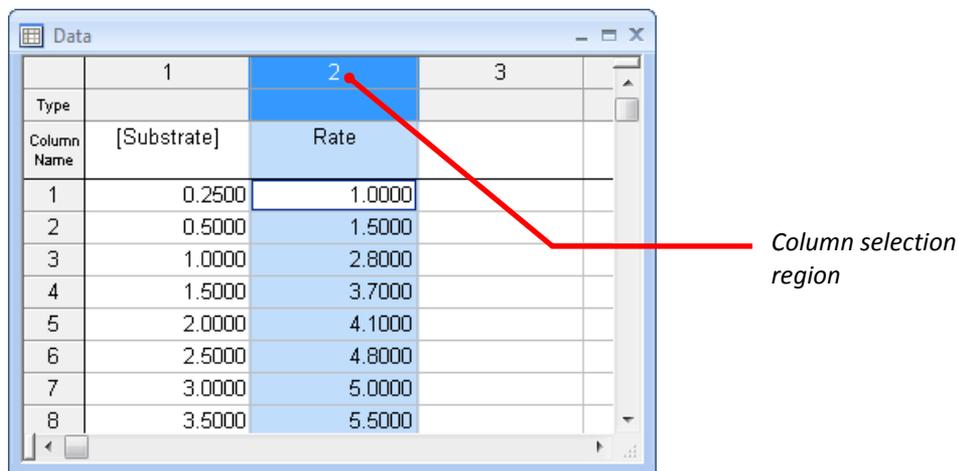

Drag the mouse along the column selection region to extend the selection. In a similar fashion, clicking in the row selection region allows whole rows to be selected.

Copying and Pasting Data

Selected regions of the data table are copied to the Clipboard using the **Copy** command, and pasted using **Paste**²². Data may also be pasted to and from other applications in the same way.

Copying and Pasting Regions of the Data Table

1. Select the region of data that you wish to copy.
2. Select the **Copy** command²³.
3. Move the data entry position to where the data is to be copied.
4. Choose **Paste**. The data will be inserted at the data entry position.


Copying and Pasting Columns

Entire columns of data can also be copied and pasted. Pasting into the data table when a column is selected results in the pasted values being placed in a new column, in front of the selected column.

²² The **Edit ▶ Copy** and **Edit ▶ Paste** commands in legacy mode.

²³ The **Edit ▶ Copy** command in legacy mode


1. Select the column(s) of data that you wish to copy.

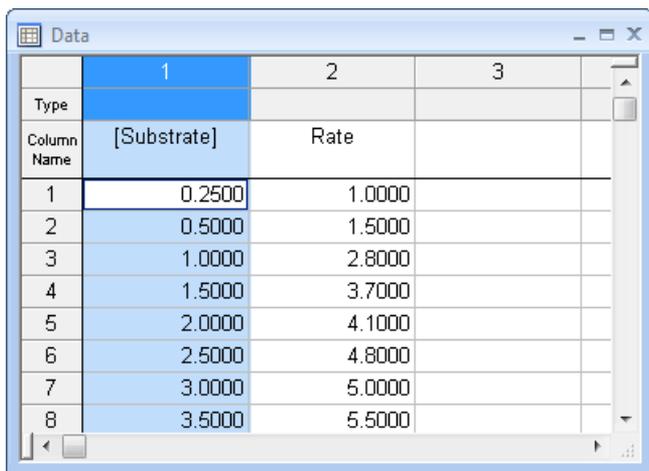

The screenshot shows a window titled "Data" with a table. Column 1 is highlighted in blue. The table has 8 rows and 4 columns. The first two columns are labeled "Substrate" and "Rate".

	1	2	3	
Type				
Column Name	[Substrate]	Rate		
1	0.2500	1.0000		
2	0.5000	1.5000		
3	1.0000	2.8000		
4	1.5000	3.7000		
5	2.0000	4.1000		
6	2.5000	4.8000		
7	3.0000	5.0000		
8	3.5000	5.5000		


2. Select the **Copy** command²³
3. Select the column in front of which the data are to be inserted.

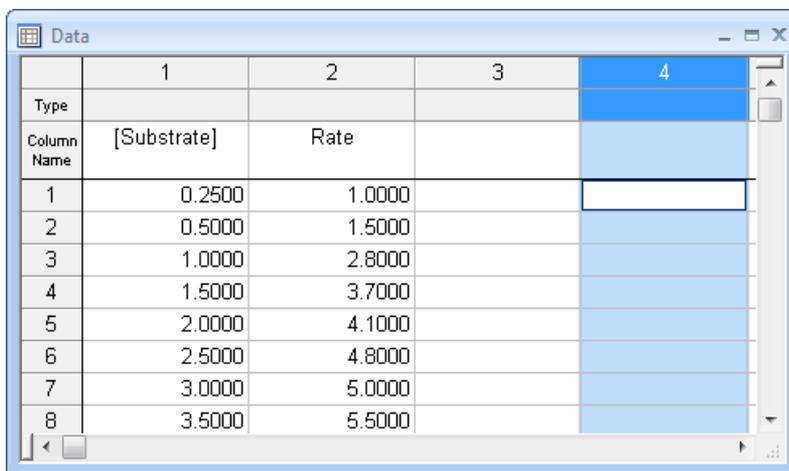

The screenshot shows the same "Data" window. Column 4 is now highlighted in blue. The data in the other columns remains the same.

	1	2	3	4
Type				
Column Name	[Substrate]	Rate		
1	0.2500	1.0000		
2	0.5000	1.5000		
3	1.0000	2.8000		
4	1.5000	3.7000		
5	2.0000	4.1000		
6	2.5000	4.8000		
7	3.0000	5.0000		
8	3.5000	5.5000		

4. Choose the **Paste** command. When pasting into a fully selected column, you have the option to overwrite the selected data or to insert the new data in front of the selection (the overwrite option is only available if the number of columns available to paste is the same size as in the selected region).

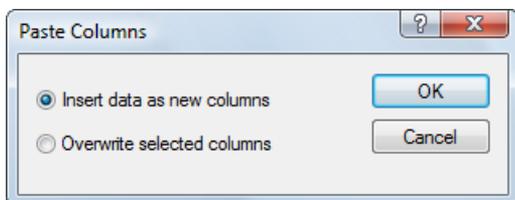

The data are then pasted into the table.

	1	2	3	4
Type				
Column Name	[Substrate]	Rate		[Substrate]
1	0.2500	1.0000		0.2500
2	0.5000	1.5000		0.5000
3	1.0000	2.8000		1.0000
4	1.5000	3.7000		1.5000
5	2.0000	4.1000		2.0000
6	2.5000	4.8000		2.5000
7	3.0000	5.0000		3.0000
8	3.5000	5.5000		3.5000

Removing Data

Clearing a Region of Data

The **Clear** command on the **Home** tab²⁴, or the DEL key, can be used to clear the currently selected data or columns.

1. Select the data value(s) (or columns or rows) to be cleared.

	1	2	3	
Type				
Column Name	[Substrate]	Rate		
1	0.2500	1.0000		
2	0.5000	1.5000		
3	1.0000	2.8000		
4	1.5000	3.7000		
5	2.0000	4.1000		
6	2.5000	4.8000		
7	3.0000	5.0000		
8	3.5000	5.5000		

²⁴ The **Edit** ► **Clear** command in legacy mode

- Choose the **Clear** command, or press DEL.

	1	2	3	
Type				
Column Name	[Substrate]	Rate		
1	0.2500	1.0000		
2	0.5000	1.5000		
3				
4				
5	2.0000	4.1000		
6	2.5000	4.8000		
7	3.0000	5.0000		
8	3.5000	5.5000		

Data in the selected region are cleared, leaving a gap in the table.

Deleting a Region of Data

Select the data value(s) (or columns or rows) to be cleared.

Choose the **Delete** command²⁵ on the **Table** tab.

When deleting a range of cells, the data values are moved up so that no gap appears in the table. When deleting a column, the remaining columns to the right of the selection are moved along to fill the gap.

Inserting Values into the Table

Values are inserted into the data table using one of the **Insert** commands²⁶ on the **Table** tab.


Command	Action
Insert Cells	Inserts blank cells in place of the current selection. Cells are moved down to accommodate the new entries.
Insert Column	Inserts blank columns in place of the current selection. Columns are moved to the right to accommodate the new entries.
Insert Row	Inserts blank rows in place of the current selection. Rows are moved down to accommodate the new entries.

²⁵ The **Table** ► **Delete** command in legacy mode

²⁶ The **Table** ► **Insert** commands in legacy mode

Sorting Data

The data in one or more columns can be sorted into ascending or descending order using the **Sort** command²⁷ on the **Table** tab, after first selecting the columns that you wish to sort.

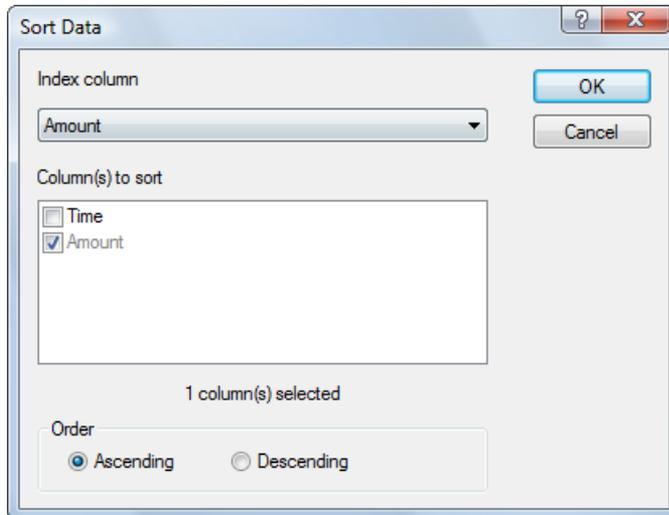

The **Index column** selection specifies the column on which the sort order is based. In addition to this column, any of the other columns can be included in the sort by checking the box to the left of the column name.

Data Table Appearance

Colors and Fonts

The colors used to display the items in the data table, and the font that is used to draw the table, are set using the **Preferences** command²⁸ on the **Table** tab.


Number Format

The manner in which numbers are formatted for display can be set on a column by column basis. The **Column Format** command²⁹ on the **Table** tab is used to control the displayed numerical precision. This command applies to all columns that fall within the current selection.

²⁷ The **Table ▶ Sort Column(s)** command in legacy mode

²⁸ The **Table ▶ Preferences** command in legacy mode

²⁹ The **Table ▶ Column Format** command in legacy mode

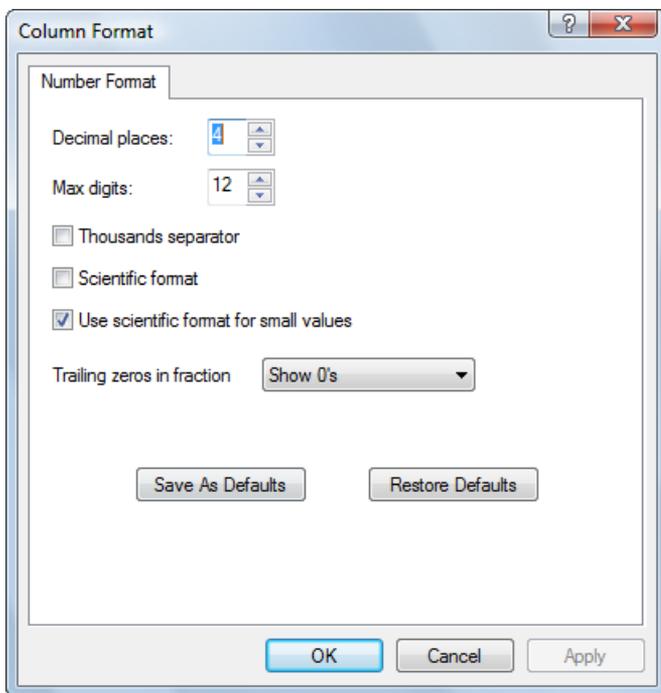

Column Format dialog box

Specification	Meaning
Decimal places	Enter the number of decimal places to be displayed for numbers. <i>Note To ensure that all decimal places are displayed set Trailing zeros in fraction to "Show 0's."</i>
Max digits	Enter the maximum number of digits to use when displaying numbers.
Thousands separator	Check this box to include a grouping separator. The symbol used for this separator, and the number of digits in the group, are set for the Options command ²⁰ on the View tab.
Scientific format	Check this box to display the values using scientific format (e.g. -2.435e003 instead of -2435.000).
Use scientific format for small values	Check this box to ensure that small values, which would otherwise be displayed as 0.0000, are shown in scientific format (e.g. 2.0987e-008).

Specification	Meaning
Trailing zeros in fraction	Chooses the way that trailing zeros are displayed in the formatted numbers. It is possible to show the "0" characters, to omit them, or to replace them with spaces.
Save As Defaults	Stores the current settings as the default for future GraFit sessions.
Restore Defaults	Restores the current GraFit formatting defaults to the selected columns.


Magnification

The magnification of the data table text can be adjusted to allow more or less values to be visible. Magnification is controlled by the **Zoom** command³⁰, found on the **View** tab.

Changing the magnification and altering the font size (**Preferences** command²⁸ on the **Table** tab) both modify the size of text displayed on-screen. However, changing font size using **Preferences** also changes the size of printed text, while **Zoom** only affects the on-screen appearance.


Using the Zoom Slider on the Status Bar

It is also possible to adjust the magnification of the data table text using the zoom slider found on the status bar.

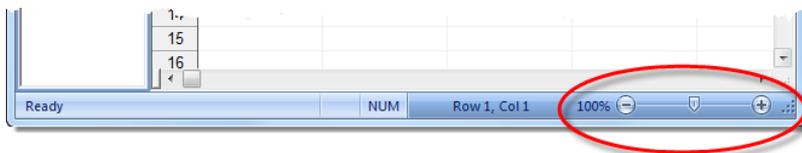

The zoom slider on the status bar

Select the data table then drag the slider to adjust the magnification. The magnification can be set more precisely by clicking on the zoom level text

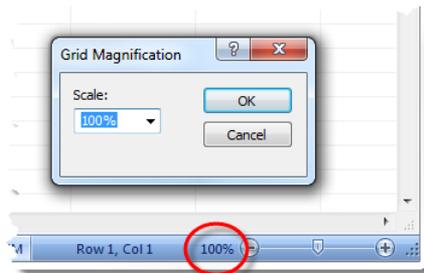

Set the Grid Magnification by clicking on the zoom level text

³⁰ The **View ▶ Zoom** command in legacy mode

Importing and Exporting Data

Importing Data

To import data into a data table, position the data entry position where you want the data values to start and select the **Paste ▶ Paste From** command³¹ on the **Table** tab. Select the file that holds the values to import, and press the **Open** button. The contents of the file will be inserted into the data table.

Importing Data as a New Data Table


It is also possible to import a file into the GraFit workspace as a new data table item. To import a data file in this way, use **Import Data** command³² on the **Home** tab.

Exporting Data

To export data from a data table into a file that can be read by other programs, select the range of values to be exported and use the **Copy To** command³³ on the **Table** tab.

Windows Vista and Windows 7

When running GraFit on Windows Vista or Windows 7, the program uses the new Vista-style Save dialog box, as shown below.

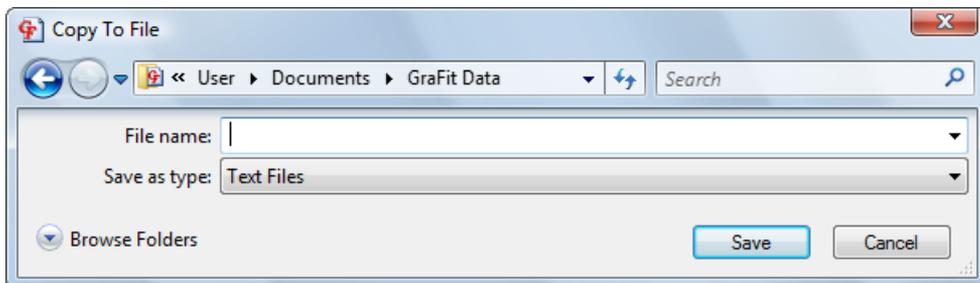

The **Save as type** entry determines the format in which to save the exported data. GraFit supports a number of data export formats, and the available types are listed in this box. GraFit supports export as ASCII (Text), DIF and CSV formats as standard, though additional formats may be available.


Windows 2000 and Windows XP

Earlier versions of Windows use the older-style file save dialog box.

³¹ The **Edit ▶ Paste From** command in legacy mode

³² The **File ▶ Import Data** command in legacy mode

³³ The **Edit ▶ Copy To** command in legacy mode

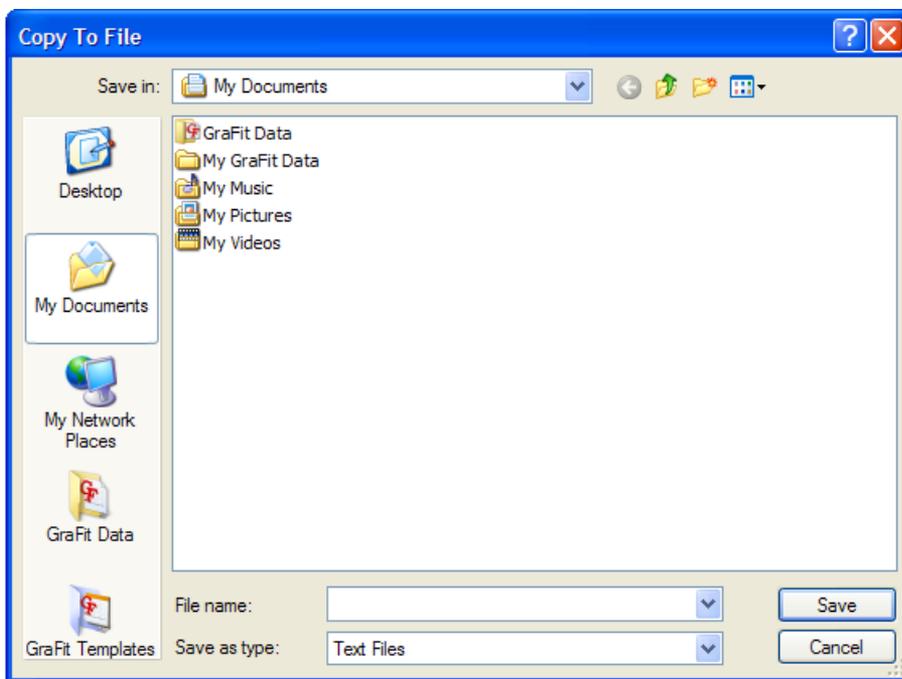


The **Save as type** entry determines the format in which to save the exported data. GraFit supports a number of data export formats, and the available types are listed in this box. GraFit supports export as ASCII (Text), DIF and CSV formats as standard, though additional formats may be available.

Column Types


Columns can hold simple data values, or can be set to hold the results of various calculations or data manipulations. For example, a column can hold the difference between the values in two other data columns. The results of these calculations are automatically updated when the source data change, allowing powerful processing of the underlying data.

Example: setting up a difference calculation

As an illustration, we will show how to set up the following data table in which the third data column holds the difference between the values in the first and second columns.

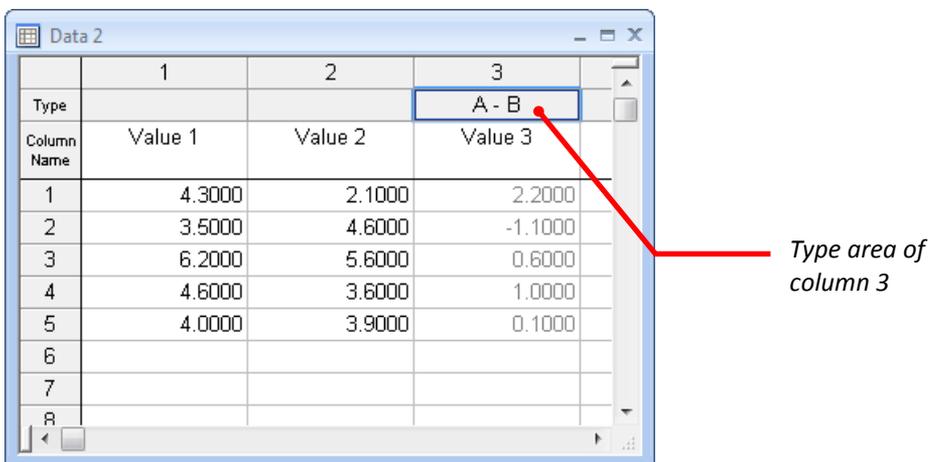

The column type of the third column needs to be set by moving the data entry position to somewhere within this column and choosing the **Column Type ▶ Column Combinations ▶ A - B** command³⁴ from the **Table** tab. Alternatively, the **Column Type** submenu can be accessed by double-clicking in the type area of the column. For this column type setting, the following dialog is displayed.

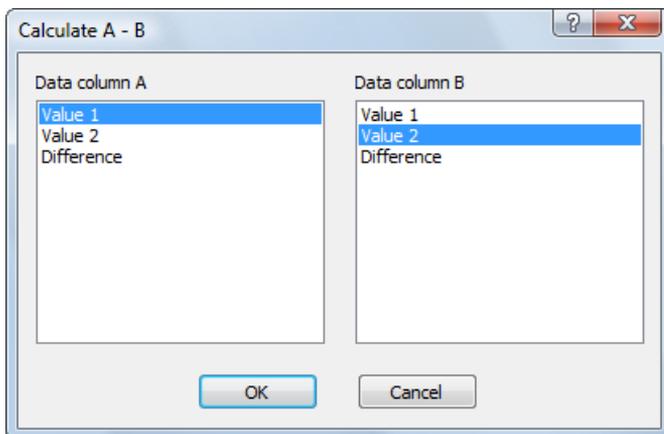

Choosing **Value 1** as the first column and **Value 2** as the second gives the result above. The values that are displayed in this column cannot be overwritten or edited, and are updated whenever the values in the source columns changes. The text color of the calculated values can be adjusted using the **Preferences** command²⁸ on the **Table** tab, but is typically shown a different color than the editable text in order to distinguish it. The dialog settings for the calculated values can be altered by repeating the command, or by choosing the **Properties** command³⁵ on the **Table** tab.

³⁴ The **Table ▶ Column Type ▶ Column Combinations ▶ A - B** command in legacy mode

³⁵ The **Table ▶ Column Properties** command in legacy mode

To reset the column type to the default (non-calculated) settings, choose **Column Type ▶ (General)**.

Note Even though columns that hold calculated values do not allow these values to be edited, it is always possible to modify the column name entry.

(General) Column Type

By default, all columns are set to be of (General) type, and hold general data values that can be edited normally.

Data Types


The column types under this submenu are used to designate that the column holds a certain type of data, which may have associated additional values. The entries in the column are editable, but the data held in the column may have special significance for data fitting.

x Data

You may optionally select this column type if the data in the column are used as the x values for data fitting. Designating a column to contain x data is not essential for data fitting, but can be useful to document the data layout.

x(n) Data

When fitting an equation that has more than one independent (x) variable, one way to hold the data is to have n separate columns for the different x_n values. During the data fitting, GraFit needs to know which column holds which x_n data set. This is done by setting the column type to **x(n) Data**, which gives the following dialog:

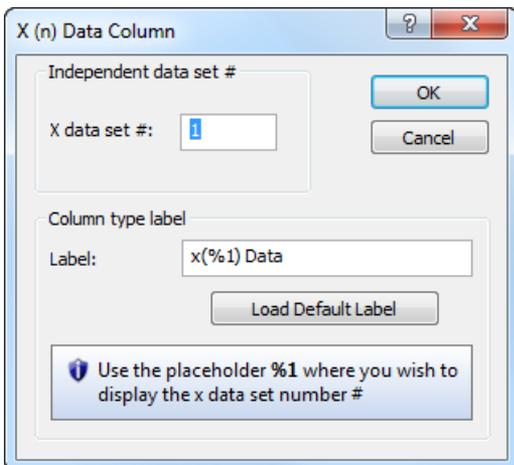


It is necessary to have one **x(n) Data** column for each independent variable in the equation, and these must have their **X data set #** entries set to correspond to the **X variables** number in the equation definition. See Chapter 6 for information about using this option when fitting data.

y Data


You may optionally select this column type if the data in the column are used as the y values for data fitting. Designating a column to contain y data is not essential for data fitting, but can be useful to document the data layout.

y Data with Associated Value

It is common that an equation has two independent (x) variables. Under these circumstances, the most usual way of collecting data is to measure a series of y values at varying x_1 for a fixed x_2 , then to collect a further set of y values at the next x_2 value, etc. For example, the variation of enzyme catalyzed reaction rate might be measured as a function of substrate concentration at a certain fixed inhibitor concentration. The measurements are then repeated at the next inhibitor concentration and so on. Here, the second independent variable (x_2 , in this example is the inhibitor concentration) is the same for an entire set of x_1, y (substrate, rate) data pairs. Under these circumstances, it simplifies the data entry to associate a single x_2 value with a whole column of different y values. This is done by specifying the column type as **y Data with Associated Value**.

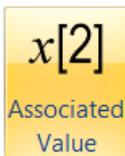

The **y Data with Associated Value** command is also available from the **Table** tab³⁶.

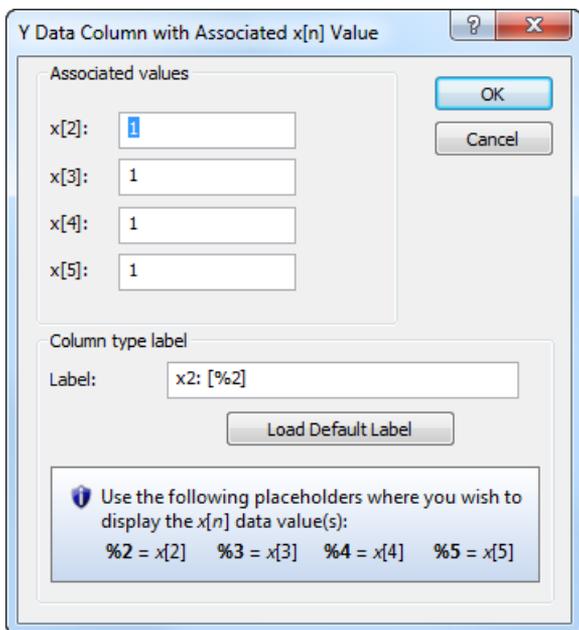

³⁶ It is found on the Data toolbar in legacy mode

The x_2 data value is entered in the above dialog, and is normally displayed in the column type area. Data columns that have associated values are used when fitting to equations with two independent variables. See Chapter 6 for information about using this option when fitting data.

Up to four values can be associated with a data column, allowing x_2 , x_3 , x_4 and x_5 values to be stored, if required, for equations having up to five independent (x) variables. If your equation definition has only 2 independent variables, the values for x_3 , x_4 and x_5 are ignored.

Column Combinations

These column types allow the values within the designated column to be the result of an arithmetic combination of two other columns. In each case, the number of calculated entries depends upon the number of values in two selected columns, and the values are updated whenever one or more data values in the source columns change.

A - B

The column holds the result of subtracting values from the second column from those in the first.

A + B

The column holds the sum of adding values in the first and second columns.

A / B

The column holds the result of dividing values in the first column by those in the second.


A . B

The column holds the product of values in the first and second columns.

Statistics

The column types under this submenu allow the result of various statistical procedures to be displayed in the column.

Mean

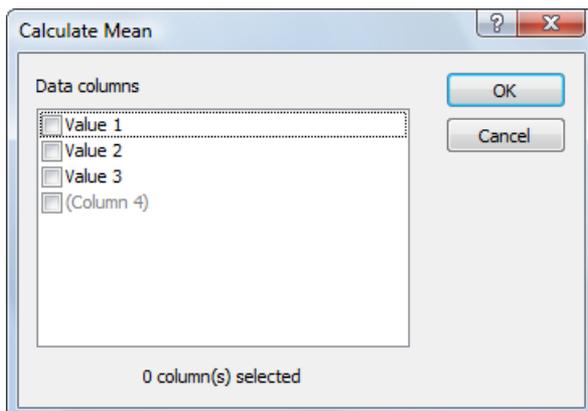

Each row entry in the column is set to be the calculated mean of the corresponding values from the columns selected in the above dialog.

$$\bar{x} = \frac{1}{N} \sum x_i$$

The calculated values are updated whenever one or more values within the selected columns are modified.

No. Data Points

Each row entry in the column is set to be the number of data values (N) present in the corresponding rows from the selected columns. This will be less than the total number of selected columns if missing data values are present.

Standard Deviation

Each row entry in the column is set to be the calculated standard deviation of the corresponding values from the columns selected in the above dialog. Standard deviation, σ , is defined as:

$$\sigma = \sqrt{\frac{1}{N-1} \sum (x_i - \bar{x})^2}$$

Standard Error of the Mean

Each row entry in the column is set to be the calculated standard error of the mean for the corresponding values from the columns selected in the above dialog. SEM, σ_μ , is defined as:

$$\sigma_\mu = \frac{\sigma}{\sqrt{N}}$$

Column Summary

The entries in the column give a summary of the values present within a single selected column.

The screenshot shows a 'Column Summary Statistics' dialog box with 'Value 2' selected in the 'Source data column' list. A red arrow points from this dialog box to the 'Summary' column in the 'Data' table. The 'Data' table has four columns: 1, 2, 3, and 4. The 'Summary' column contains statistical data for 'Value 2'.

Type	1	2	3	4
Column Name	Value 1	Value 2	Value 3	Summary
1	4.3000	2.1000	2.2000	Source: Value 2
2	3.5000	4.6000	-1.1000	Data points 5.0000
3	6.2000	5.6000	0.6000	Minimum 2.1000
4	4.6000	3.6000	1.0000	Maximum 5.6000
5	4.0000	3.9000	0.1000	Range 3.5000
6				Total 19.8000
7				Mean 3.9600
8				Median 3.9000
9				Variance (n) 1.3384
10				Variance (n-1) 1.6730
11				Std. dev. (n) 1.1569
12				Std. dev. (n-1) 1.2934

This summary gives the number of data points, minimum and maximum data values, the data range, sum, mean, median, variance (calculated for a divisor of n and of $n-1$) and standard deviation (calculated for a divisor of n and of $n-1$).

F Test


The entries in the column show the results of performing an F test comparing the goodness of fit of two fitting analyses that use different equations. A fuller description of this test is found in the theory section (Chapter 20), but briefly the F test evaluates whether one fit is significantly better than another. This test can be used to ascertain whether it is realistic to use an equation with additional parameters over a simpler version, or whether one equation fits a data set significantly better than does a second. For example, the analysis shown below fits an exponential decay to a single exponential equation:

$$y = Ae^{-kt}$$

In addition it fits the data to an equation that has an additional background parameter:

$$y = Ae^{-kt} + \text{background}$$

A sample set of data fitted to each of these equations is shown below.

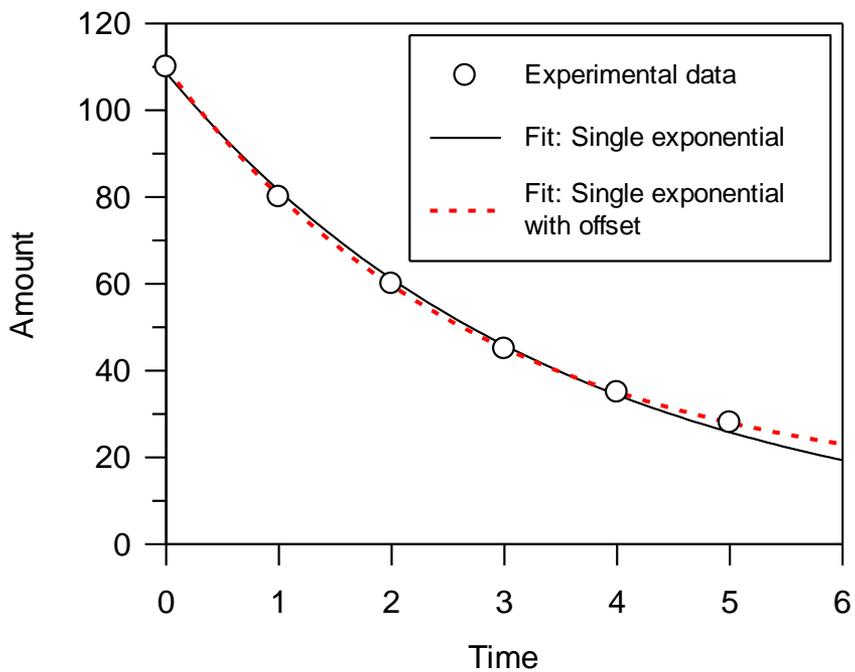

Using the F test, it is possible to test whether the data are fitted significantly better by the presence of the extra parameter. By eye, the fit obtained with the additional background parameter seems slightly better than the one without, but this is subjective and it is better to use a statistical test.

The screenshot shows two windows. The 'F Test' dialog box is in the foreground, with 'Fitting results 1' set to 'Fit Single' and 'Fitting results 2' set to 'Fit Single + Offset'. Below it is the 'Data' table, which contains the results of the F Test. A red arrow points from the 'Fitting results 2' dropdown in the dialog to the 'F Test' column in the data table.

Type	1	2	3	4
Column Name	Time	Amount	Difference	F Test
1	0.0000	110.0000	1.3654	Analysis 1
	1.0000	80.0000	-1.4794	Fit Single
	2.0000	60.0000	-1.1121	Single Exponential Decay
	3.0000	45.0000	-0.8360	Number of parameters: 2.0000
	4.0000	35.0000	0.6216	
6	5.0000	28.0000	2.2151	Analysis 2
7				Fit Single + Offset
8				Single Exponential Decay + Offset
9				Number of parameters: 3.0000
10				
11				Number of data points: 6.0000
12				Chi ² value (1): 11.2814
13				Chi ² value (2): 0.2757
14				F statistic: 119.7674
15				Probability: 0.0016

The results show that the probability that the two fits are equally appropriate is 0.0016 or 0.16%, which is low. This indicates that we can have reasonable confidence in the merit of including the extra parameter (generally we require this value to be <5%, but may prefer a value <1% to be sure).

Note Although the F test results have been shown incorporated into the same data table as the raw data, it is also possible to create the analysis in a separate data table. Using separate data tables can be preferable if you wish to keep the raw data separate from any analysis.

Data Fitting

The column types under this submenu allow the results from data fitting to be displayed in a column, or entry of parameter/constants values associated with an equation. This lets the results be viewed, but more importantly allows the values to be used for subsequent analyses. For example, the results of one set of data fitting might provide the source data for subsequent fits.

Results

This column type lists the results from data fitting.

The image shows two overlapping windows from the GraFit software. The background window is titled "Data Fitting Results" and has a "Fit" dropdown menu set to "Fit Single". Below it, under "Parameters", there are two checked checkboxes: "Initial value" and "Rate constant". At the bottom, under "Display", there are three options: "Parameter value" (selected with a radio button), "Std. error value" (unselected), and "Value name" (checked with a checkbox). The foreground window is titled "Data" and displays a table with 10 rows and 4 columns. The columns are labeled "1", "2", and "3". The first column is labeled "Type" and "Column Name". The second and third columns are labeled "Initial" and "Rate" respectively. The fourth column is labeled "Results". The data in the table is as follows:

Type	1	2	3
Column Name	Initial	Rate	
1	0.0000	110.0000	Initial value 108.6349
2	1.0000	80.0000	Rate constant 0.2876
3	2.0000	60.0000	
4	3.0000	45.0000	
5	4.0000	35.0000	
6	5.0000	28.0000	
7			
8			
9			
10			

A red arrow points from the "Value name" checkbox in the "Data Fitting Results" dialog to the "Initial value" and "Rate constant" labels in the "Data" table.

After the fit to be displayed is selected, the parameters available are listed. Some or all of these can be included by clicking in the check box to the left of the parameter name. If more than one set of data were analyzed in the selected fit, the results for each fit are given. It is possible to display either the fitted parameter value, or the fitted standard error value, as chosen in the **Display** settings. In addition, the name of the parameter can optionally be displayed next to the data value, by checking the **Value name** option.

Information

This column type lists information about a data fitting item.

The 'Data Fitting Information' dialog box has the following options checked:

- X data column name
- Y data column name
- Error data column name
- Number of data points
- Chi² value
- Reduced Chi² value

The 'Display value name' checkbox is also checked.


The 'Data' table shows the following data:

Type	1	2	3
Column Name	Initial	Rate	Fit Info
1	0.0000	110.0000	X column Initial
2	1.0000	80.0000	Y column Rate
3	2.0000	60.0000	# data 6.0000
4	3.0000	45.0000	Chi ² 11.2814
5	4.0000	35.0000	Reduced Chi ² 2.8204
6	5.0000	28.0000	
7			
8			
9			
10			

After the fit to be displayed is selected, the available information items are listed. Some or all of these can be included by clicking in the check box to the left of the parameter name. If more than one set of data were analyzed in the selected fit, the information for each fit is given. The name of the information type can optionally be displayed next to the data value, by checking the **Display value name** option.

Parameter Entry

The Parameter Entry column type gives a convenient way of entering parameter values for a particular equation into a column. Certain graph plotting options allow these values to be used when displaying graphs. The equation to use is selected, after which it is possible to enter values into the column.

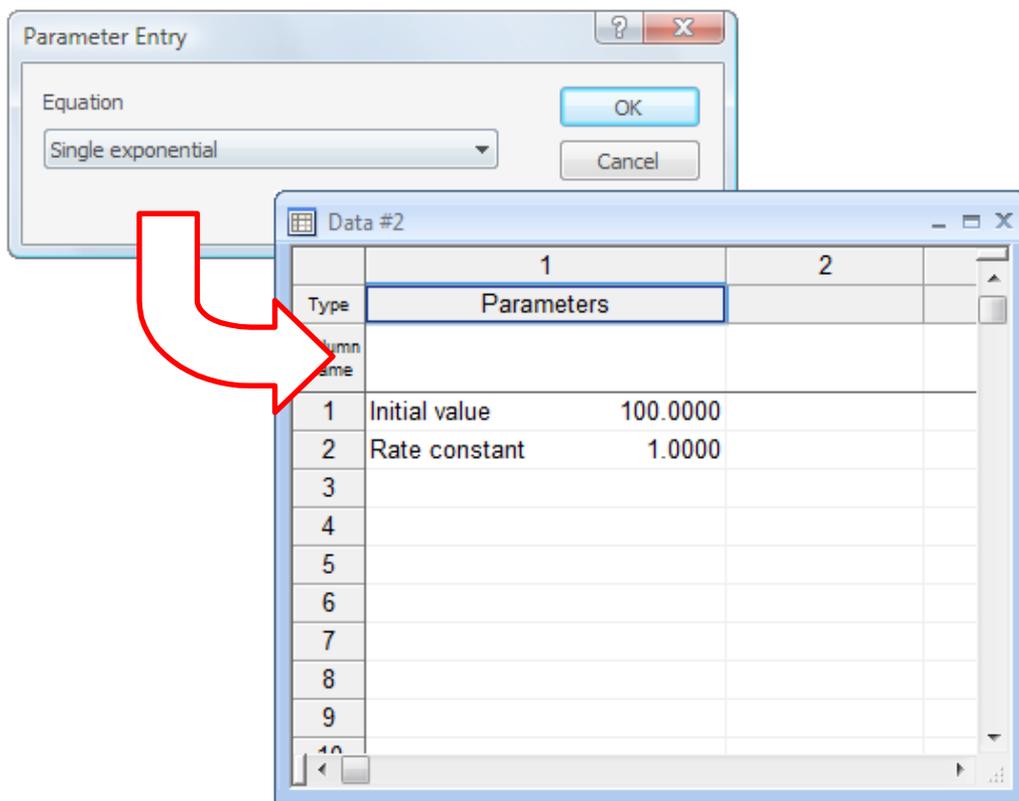

The column display differs from normal in that the name of the parameter remains displayed on the left, with the value on the right. When a new parameter value is entered from the keyboard, the parameter name remains but the revised value is displayed in the table.

Constants Entry

The Constants Entry column type allows the values of any constants associated with an equation to be entered in the data table. This option works in the same way as the Parameter Entry option described above, but shows the constants associated with the selected equation.

Calculated Values

The column types under this submenu allow the x or y values to be calculated using the results from a data fitting item.

x Values

Use this column type to calculate theoretical x values corresponding to user-entered y values, according to the fitting model described in a selected workspace item.

The 'Calculated X Values' dialog box is shown with the following settings:


- Fitting results: Fit Single
- Use fitted values (selected), Use initial estimates
- Fit: 1
- Y data column (y): Rate (Column 3)
- Interpolate between X limits: From: -1, To: 10

The 'Data' table shows the following data:

Type	1	2	3	4
Column Name	Initial	Rate	x = ?	
1	0.0000	110.0000	-0.0434	
2	1.0000	80.0000	1.0637	
3	2.0000	60.0000	2.0638	
4	3.0000	45.0000	3.0640	
5	4.0000	35.0000	3.9377	
6	5.0000	28.0000	4.7134	
7				
8				
9				
10				

The fitting item that holds the results is chosen in the dialog box. It is possible to use either the results from the data fitting or the parameter values used as the initial estimate values (starting values only) when carrying out the calculations. Values are calculated for x values that correspond to the y values present in the selected **Y data column** entry.

A potential problem with calculating x values for a given y value is that some equations will have multiple solutions. For example, consider the curve shown below:

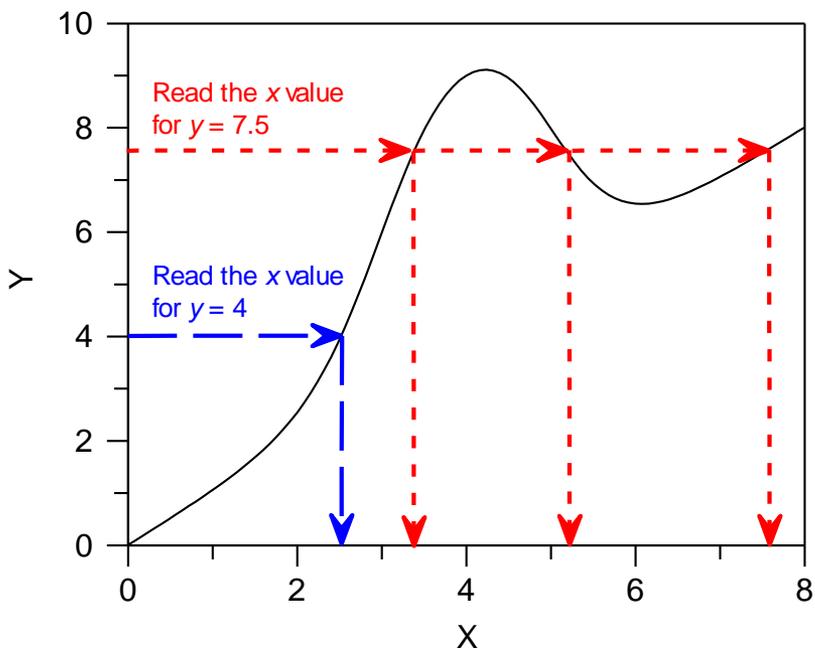


At $y = 4$ the value of x read from the curve is 2.5; but at $y = 7.5$ there are three x values, 3.4, 5.2 and 7.5. Clearly, whether the curve is single or multiple-valued in x depends upon the equation in use and the value of y that is used.

To overcome this complication it is necessary to specify limits to the x values between which the data are to be calculated. The **Interpolate between X limits** entries in the above dialog control the range between which the values are tested. For example, in the plot above the solution $x = 3.4$ at $y = 7.5$ would be correctly found by setting interpolation limits between $x = 2$ and $x = 4$.

y Values

Use this column type to calculate theoretical y values corresponding to user-entered x values, using the fitting model described in a selected workspace item.

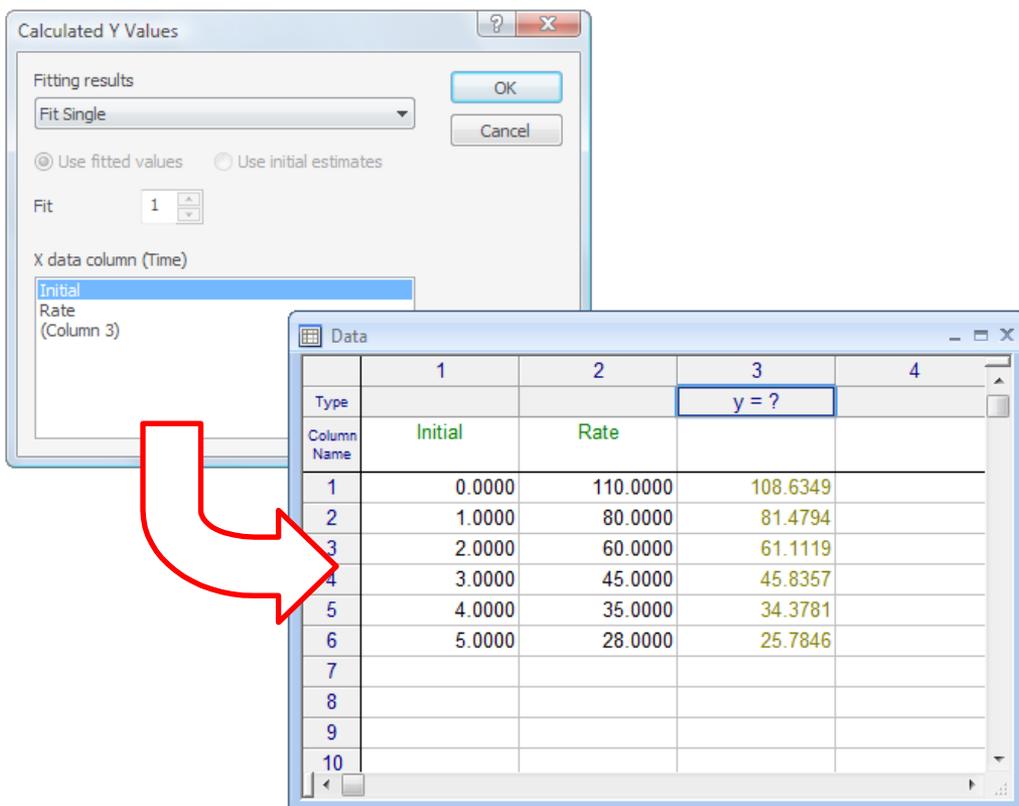


All fitted GraFit equations have only a single y value that corresponds to each x value, and so there is no need to specify interpolation limits as in the previous column type.

When calculating y values, it is possible to choose whether to use the fitted values or the explicit initial estimates that were used in the Estimatrix (see Chapter 6). This allows the generation of a hypothetical data set, which would result from the initial estimate values. When no initial estimates were used, these options are unavailable (as in the figure above, as the estimates were generated automatically for the equation used).

Difference

Use this column type to calculate the difference between the experimental and the calculated y values, using the fitting model described in a selected workspace item.

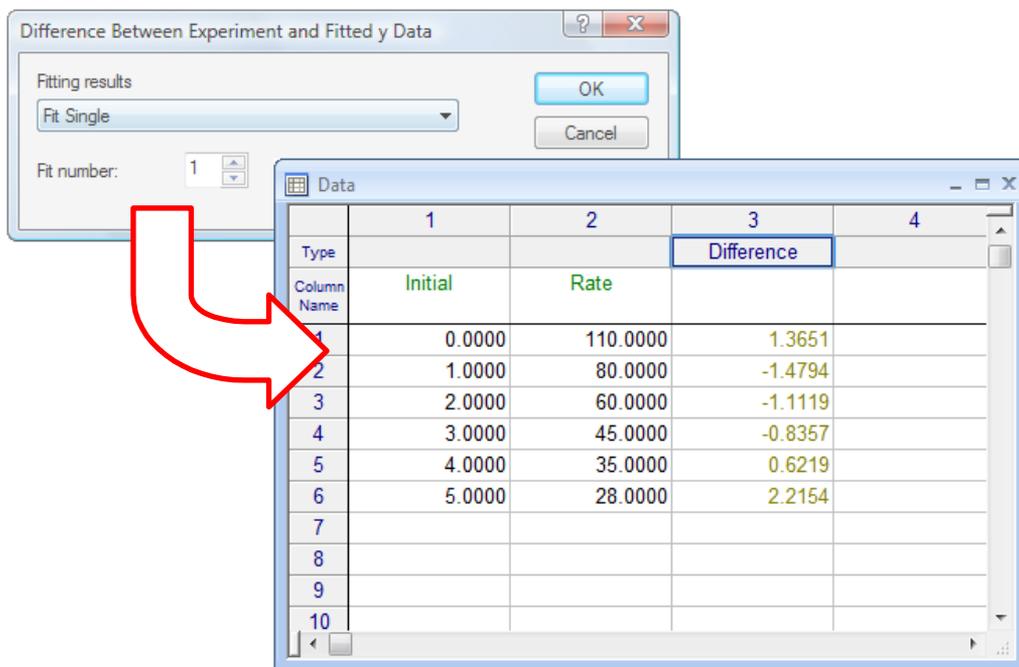


Plots of difference versus x or y values are very useful to indicate whether the calculation accurately reflects the variation in the data. A good fit is characterized by a random distribution of (small) differences between the raw and fitted data.

Data Processing

The column types under this submenu allow various types of processing to be performed to generate data values for the column.

Calculated Function

The calculated function column type fills the selected column with the results of a user-defined function. This column type can be used to perform sophisticated processing of data, including background subtraction, calculation of logarithmic values and complex rescaling. The calculated function dialog box is shown below.

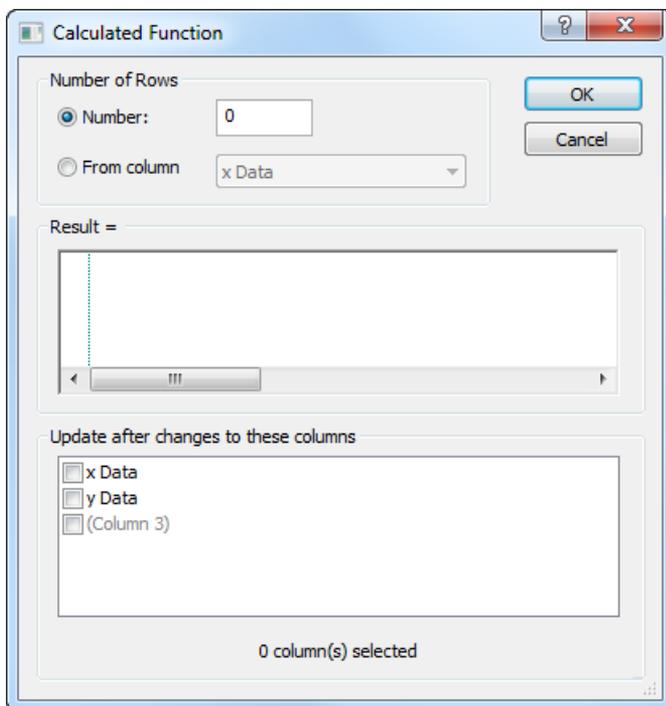

The **Number of Rows** of data to be calculated can be a defined number, or can be set to be the same as the number of entries in one of the other columns. The function definition is entered into the **Result =** field, and uses the same syntax as is used for equation definitions (see Chapter 9). By default, the text entered into this field is used to calculate a result without the need for an explicit assignment statement. However, the name “Result” is predefined so that it is possible to assign values using a “Result =” statement. In addition to the standard arithmetic functions (sin, cos etc.), the following can be included within this definition

Function	Description
Data (x, y)	Returns the data value present within the data table at column x , row y .
ValidData (x, y)	Returns TRUE if a valid number is present in the data table at column x , row y , and FALSE otherwise.
SetBlank	If this statement is encountered, the current cell is set to be blank and the calculation skips to the next row.
Column	Returns the number of the current column. This value is most useful in conjunction with the Data () function.

Function	Description
Row	Returns the number of the current row. This value is most useful in conjunction with the <code>Data()</code> function.
<code>_rnd_</code>	Returns a random number between 0 and 1.
<code>_grnd_</code>	Returns a random number that falls within a Gaussian distribution that has a zero mean and unit variance.


The following examples show how calculated functions might be specified and used.

Process	"Result =" field	Description
Background subtraction	<code>Data(1, Row) - 0.234</code>	Subtracts the value 0.234 from the values in column 1.
Log calculation	<code>log(Data(2, Row))</code>	Takes the logarithm (to base 10) of the values in column 2.
Complex rescaling	<code>(Data(1,Row) * 1.34) - Data(2,Row) / Data(3,2)</code>	For each row in column 1, multiplies this value by 1.34, then subtracts the data in column 2 which has been divided by the value in column 3 row 2.
Multiple line statement	<code>If (ValidData(1, Row) Then Result = Data(1, Row) + 23.3 Else SetBlank EndIf</code>	For each valid data entry in column 1, sets the corresponding entry in the current column to be 23.3 greater.

Note It is essential to inform the program, by checking the appropriate boxes, which columns (if any) the calculation depends upon. Failure to do this will result in the calculation not being updated correctly when these values change.

Derivative


The derivative column type fills the column with the calculated derivative (order 0 to 5) of a selected source data column. The contents of the source column are not affected by this operation.

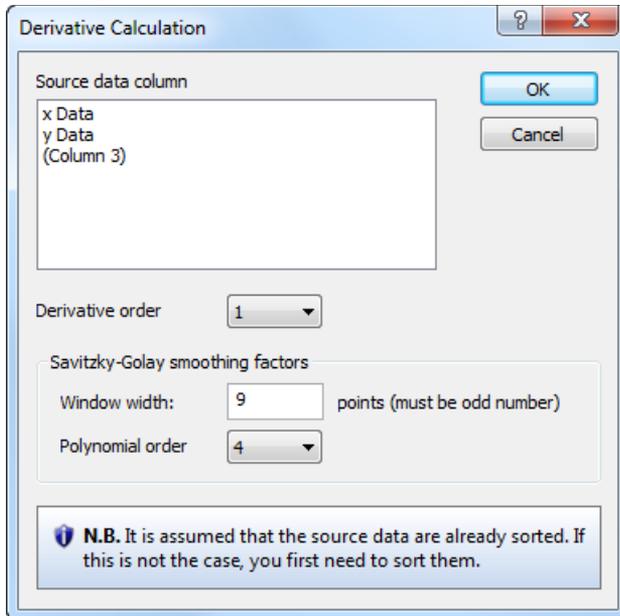

Calculation of derivatives for a set of data is performed after application of Savitzky-Golay smoothing to the source data, and so it is necessary to set the Savitzky-Golay smoothing factors. A full description of Savitzky-Golay smoothing is found under the **Smooth** data processing type.


Series

The series column type provides a convenient way of filling a column with a series of equally spaced values, for example, time-series values.

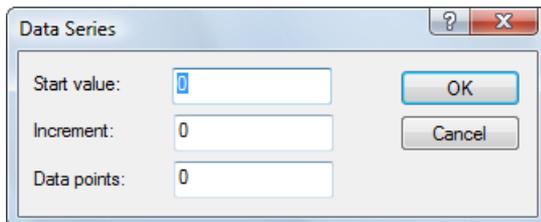

Smooth

This fills the selected column with the result of applying a Savitzky-Golay smoothing function to a source data column. The contents of the source column are not affected by this operation.

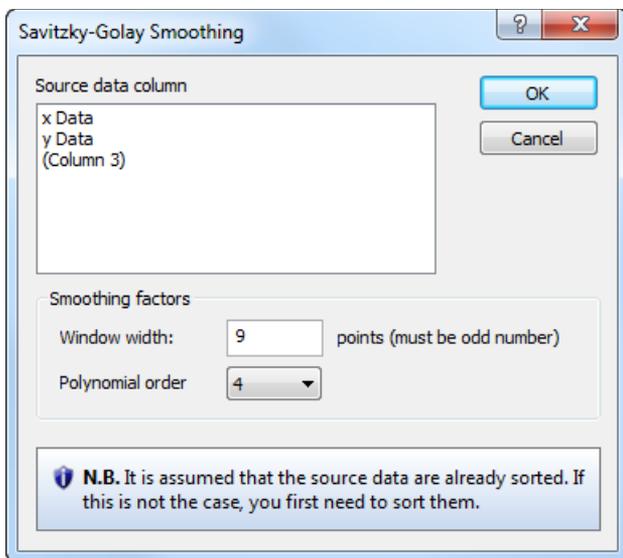


Savitzky-Golay smoothing is probably the most effective means of smoothing a set of data in such a way that noise is removed without too much loss of resolution (as happens with simple moving window averaging). In this smoothing procedure, the underlying function is approximated by a polynomial fit applied to a window of points around each data value. The degree of smoothing is therefore controlled by the window width, and by the order of the polynomial applied. The lower the order of the polynomial or the greater the window width, the smoother the result but at the expense of losing resolution. For a polynomial order of 4, best results are obtained when the window width is between 1 and 2 times the full width at half of maximum of the desired features in the data. An excellent discussion of the effects of various smoothing parameters can be found in Press, W.H., Teukolsky, S.A., Vetterling, W.T. & Flannery, B.P. (1992) *Numerical Recipes in C*. Cambridge University Press.

Note Smoothing of data should always be undertaken with care. Although smoothing can improve the appearance of a set of data, this improvement is cosmetic only. If you are fitting data, it is almost always preferable to fit to the raw data values rather than the smoothed ones.

Subset

Data that are imported from other applications, or obtained from an instrumental source such as a spectrometer, often have more data values than are necessary to define the experimental curve. For example, the data obtained from a spectrometer may contain more than a thousand points; fitting so many points would be extremely time consuming, and would not necessarily give a “better” fitted value than if (say) thirty points spread over the

data range were analyzed. The **Subset** type fills the column with a subset of the data values present in one of the other columns, and provides a convenient means of reducing the number of data values for analysis.

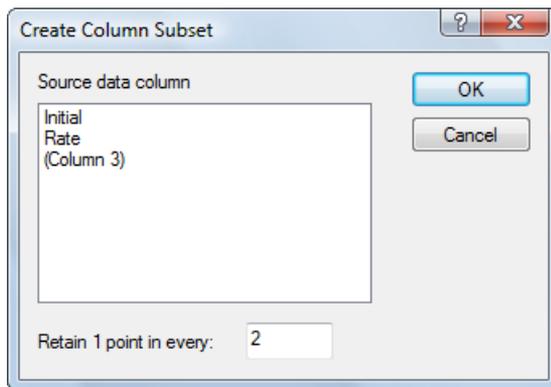

Select the source data column, and enter a value in the **Retain 1 point in every** field. For example, choosing to retain one data point in every 10 of a data set that has 200 data values would fill the column with point 1, 11, 21... giving 20 data points in this column. The contents of the source data column are not affected by this operation.

Locking a Data Column

It is possible to lock a data column to prevent any further modifications being made to the values that are present. To lock one or more columns:

- Move the data entry point to the column that you wish to lock, or select a range of values encompassing the columns that you wish to lock.
- Choose the **Lock Column** command³⁷ on the **Table** tab.

Repeat these actions to unlock the column(s).

³⁷ The **Table** ► **Lock Column** command in legacy mode

Non-linear Fit Item

Overview of Non-linear Fitting

A non-linear fit item is used to control data fitting by non-linear regression. This item specifies the data source and equation to use, and also presents the results of the data fitting. The item is created from the main GraFit window using the **Add ▶ Non-linear Fit** command. The non-linear fit window itself contains a tabbed dialog box with several panes that are used to specify how the fit is calculated. Each of these panes is discussed below. Data fitting involves making appropriate selections in the Settings pane (and possibly one or more of the other panes as described below), then selecting the **Fit Data** command³⁸ on the **Home** tab. After fitting the data the results can be viewed from the Results pane.

Differences from Linear Fitting

Non-linear regression allows much more general equations to be used for analyzing data. The biggest practical difference, however, is that non-linear regression requires starting estimates of the unknown parameter values; linear regression needs no such estimates.

Status Bar Indicators

When a fitting item is selected, the status bar shows the current fitting status, indicating whether the item is **Fitted** or **Not Fitted**. In addition, the status bar shows when a calculation is in progress. Whenever any calculation is being processed, the calculation pane shows an animated abacus.




Fitting status pane

Calculation pane

Settings Pane

The Settings pane is used to specify the main details of the fit.

³⁸ The **Calculate ▶ Fit Data** command in legacy mode

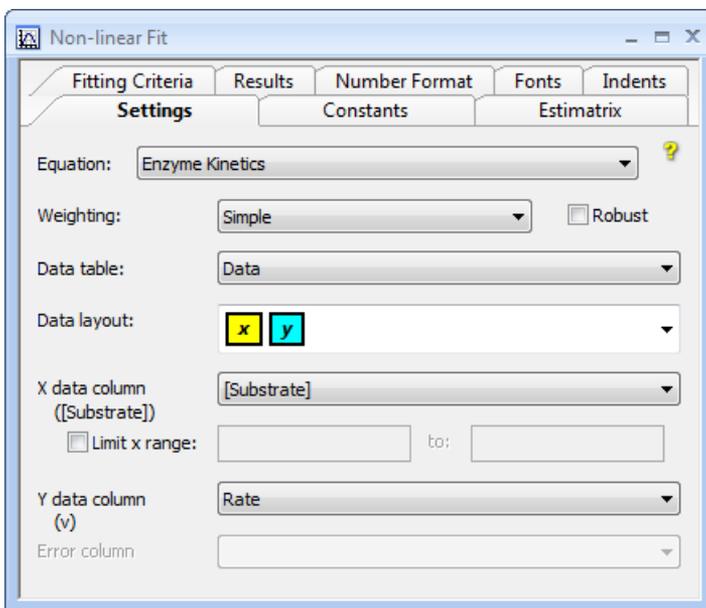

The Settings pane of the Non-linear Fit item

Equation

Choose the equation to use for the data analysis. The equation must exist as a separate item within the main GraFit window. See Chapter 9 for information about equation items.

Weighting

All types of regression analysis require information about the error distribution in your data. The reason for this is explained fully in Chapter 20, but briefly it is necessary to specify whether the inaccuracies are likely to be the same for each data point (“Simple” weighting), are proportional to the value of the y data points (“Proportional” weighting) or have been individually estimated (“Explicit” weighting). The default Simple weighting is often the most appropriate setting.



Robust weighting is additional weighting that provides a way of removing outliers (see Chapter 20). Check the **Robust** box to add robust weighting to the fit.

Data Table

The data to be fitted must reside in a separate data table item within the GraFit workspace. Select the table that holds the data from this drop-down list.


Data Layout



The data layout specifies how the data are arranged for data fitting. The default layout, as shown in the figure above, is that the x and y data columns are to be selected individually. This is used to fit a single {x, y} data set. It is also possible to specify that multiple data sets are present, which will result in several sets of results being generated. The available options for a standard equation having a single independent variable are given below.

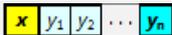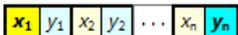Layout	Meaning
	A single data set is present and the columns holding the x and y data values are specified individually from the X data column and Y data column entries. The x and y data columns do not need to be in any particular order within the data table.
	Multiple data sets are present, although all data sets have common x data values. The x and y data columns are consecutive in the data table. The location of the x column is set from the X data column entry; the location of the final y column is set from the Last Y data column entry. The number of data sets present is determined automatically.
	Multiple data sets are present, and the data sets have different x data values. Pairs of x, y data columns are present consecutively in the data table. The location of the first x column is set from the First X data column entry; the location of the final y column is set from the Last Y data column entry. The number of data sets present is determined automatically.

Equations with Multiple Independent Variables

Equations that have multiple independent variables (x axes) need to access multiple x data columns. They therefore have a different set of data layout options.

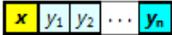Layout	Meaning
$x[n]$ from y header 	This option is only available if there are two independent variables. In these circumstances, it is often the case that a series of y data values are collected at a fixed value of the second x variable. For example, the progress of a reaction might be collected with time at a number of different reagent concentrations. The first x variable would be time, the second would be reagent concentration. Within the data table it is possible to associate a value with an entire column; when this layout is selected the associated value is taken as the second x variable. The associated value is set by setting the column type to “y Data with Associated Value” and entering a numeric x_2 value as appropriate. See Chapter 5 for further details. The x_1 and y data need to be arranged contiguously in the data table, and multiple y data sets can be present. The location of the X data column and the Last Y data column are selected,

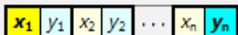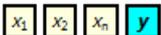Layout	Meaning
<p>$x[n]$ from y header </p>	<p>and from the program can work out how many data sets are present. All the y data sets need to be in columns that are of the type “y Data with Associated Value.”</p> <p>This option is only available if there are two independent variables. As with the previous layout, the second x variable is obtained from the column type header. The difference is in the arrangement of data columns. For this layout, pairs of x, y data columns are arranged contiguously in the data table, and a separate set of x data values are used for each y data set. All the x data columns hold values corresponding to the first independent variable as defined in the equation; values for the second independent variable are obtained from the type header. The First X data column and the Last Y data column are selected, and from this the program can work out how many data sets are present.</p>
<p>$x[n]$ from table </p>	<p>This layout is available whenever there is more one independent variable, and is the only layout type when three or more independent variables are present. When using this layout, the program determines which column is associated with which independent variable by examining information present in the column type headers of the data table. The column types of the various independent data sets (x axes) need to be set to “x(n) Data” with the “n” value representing the order of the x variable found in the equation definition. See Chapter 5 for further details on column types. The position of the y data column is selected from the Y data column entry.</p>

X Data Column

This entry specifies the column that holds the x data values. See the **Data Layout** section (above) for further information.

Limiting X Data Range

It can sometimes be useful to restrict data fitting to a subset of the values in the columns. For example, it can be useful to fit just the initial portion of a data set. Checking the **Limit X range** box allows a lower and upper limit to be set on data that will be included in the analyses.

Y Data Column


This entry specifies the column that holds the y data values. See the **Data Layout** section (above) for further information.

Error Column

If explicit weighting is chosen, the error values are read from a separate column, which is selected here.

Constants Pane

Some equations are defined to include Constants, the values of which are passed to the equation when fitting data. If the selected equation includes Constants values, these are entered into the grid found on this pane. For equations that lack Constants, this pane is empty.

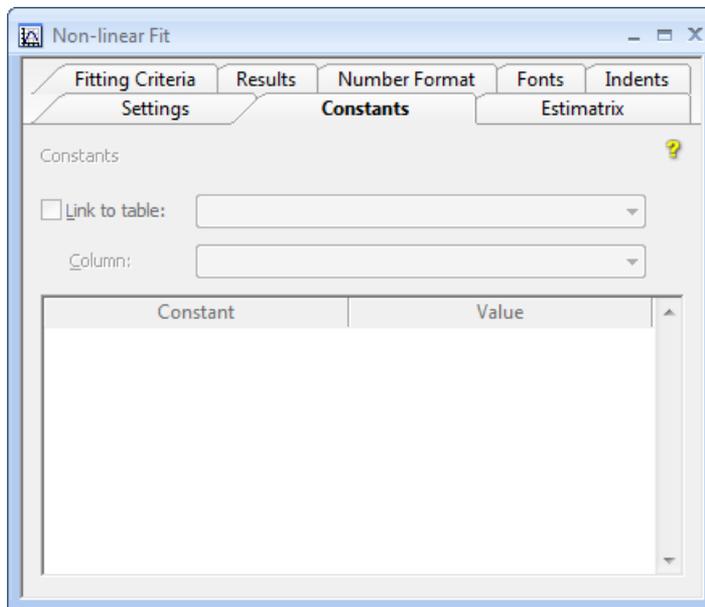


The Constants pane of the Non-linear Fit item

Linking Values to Data Table Entries

It is possible to obtain the constants values from entries within a data table. Data table entries can be updated as a result of calculations and can even further linked to values in external data files. To link values in this way, click the **Link to table** checkbox and select the data table that holds the data plus the columns that are the source of the data values.

Estimatrix Pane

Non-linear regression analysis requires approximate starting estimates for the parameter values. These starting values are then improved upon by iteration, until the best-fit curve is calculated. Some equations provide a means of generating estimates automatically; for other equations it is always necessary to enter explicit estimates (see Chapter 9). For the former case, this pane allows supplementary estimates to be entered; for the latter it is essential to enter initial estimates into the grid found on this pane.

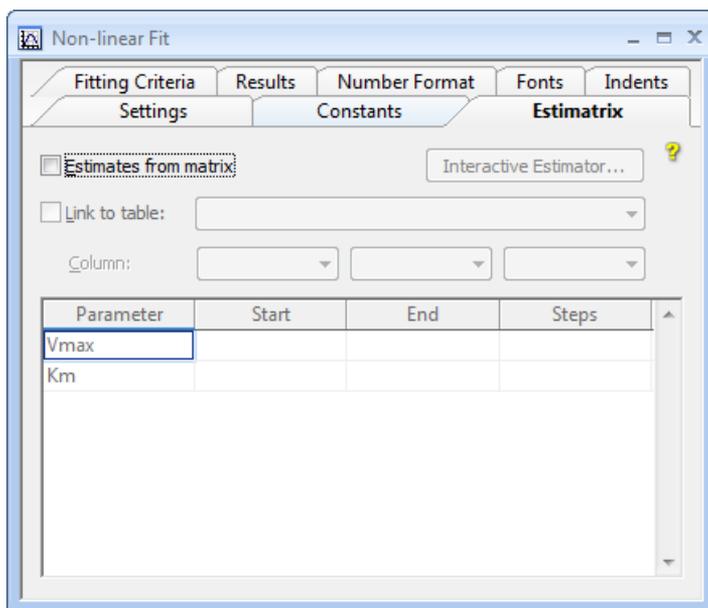

The Estimatrix pane of the Non-linear Fit item

For equations that automatically provide estimates, the **Estimates from matrix** box needs to be checked in order to allow supplementary estimates to be entered into the grid. If automatic estimation is not available then this box is always checked. The grid is used for entry of estimates, which can be done in one of two ways.

1. Just the **Start** value is specified. In this case the value entered is taken as the initial estimate for this parameter.
2. Values are entered under the **Start**, **End** and **Steps** columns. In this case a range of estimates, total number Steps, is generated between the Start and End values inclusive.


It is possible to enter just a Start value for some parameters and a range for others. Whenever a range of values is involved, a matrix of permutations is tested. It is important to limit this grid search to a sensible number of combinations, for example if 10 parameters were present and 10 steps were selected for each, this would involve 10^{10} values to be tested – clearly an unrealistic number. When data fitting starts, each permutation is tested and the combination that gives the best fit (lowest χ^2 value) is used as the starting estimates for non-linear regression. If both automatic and manual estimates are present, each is tested and the best of the two is used as the starting estimates for regression.

Linking Values to Data Table Entries

It is possible to obtain the estimate values from entries within a data table. Data table entries can be updated as a result of calculations and can even further linked to values in external data files. To link values in this way, click the **Link to table** checkbox and select the data table that holds the data plus the columns that are the source of the data values.

Interactive Estimator

Finding suitable initial estimates can be a laborious task, and to simplify things GraFit provides an interactive visual tool that is accessed by pressing the **Interactive Estimator** button. The Interactive Estimator dialog is shown below:

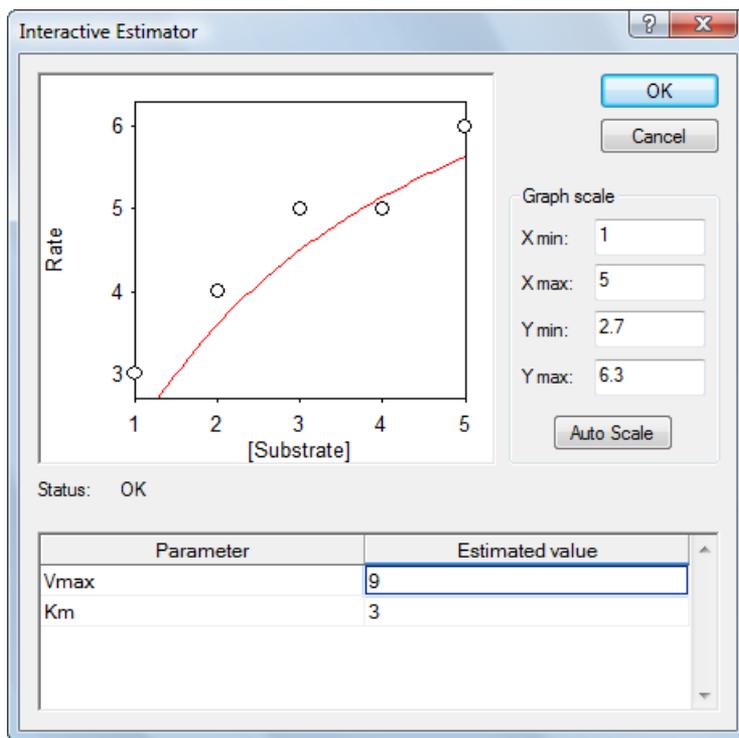


The Interactive Estimator dialog

A simple graph of the data, selected on the Settings pane, is shown. The graph is automatically scaled by default, although the scale can be altered if desired. To see the curve generated by a particular set of parameter values, type in values into the **Estimates value** column of the grid. When you move the data entry position, the graph will be updated to show the curve generated by the current set of values. These estimates should be such that the calculated curve falls close to the data values. Just how accurate these parameter estimates need to be depends upon the equation in use and the data present. However, it is usually only necessary to produce a curve that roughly follows the data (as shown above) for the subsequent data fitting to find the best-fit values.

When acceptable parameters have been found, pressing **OK** will return these values into the **Start** column of the Estimatrix pane.

Fitting Criteria Pane

As described in Chapter 20, non-linear regression calculates best-fit parameters by a series of iterations that produce successively more accurate results. It is up to the program to decide at what point the values are deemed acceptable, and at which point the calculations stop. The criteria used to decide this are set from the Fitting Criteria pane. It is likely that for most situations you will not need to change the default values on this page.

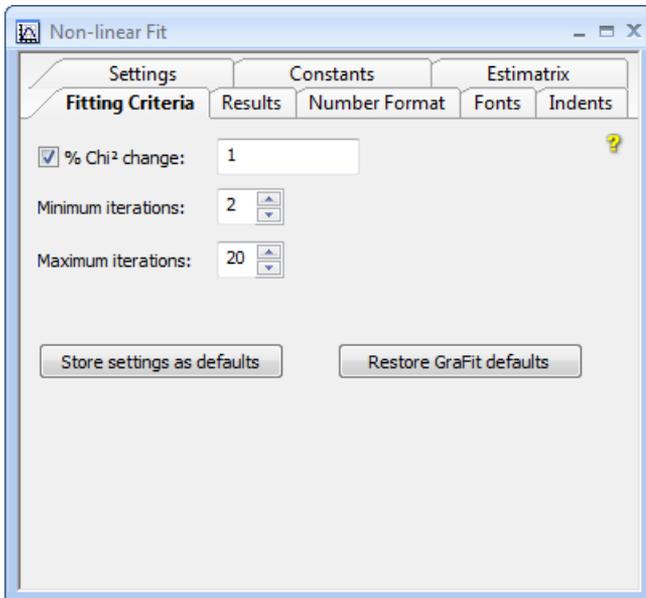

The Fitting Criteria pane of the Non-linear Fit item

Fitting Criteria

GraFit uses three separate criteria to decide when the calculations are acceptable.


Criterion	Meaning
% Chi² change	This setting is optional, and is activated by checking the % Chi ² change box. When selected, fitting will terminate when the change in the χ^2_v value is less than the percentage value entered (providing the minimum number of iterations has been reached).
Minimum iterations	GraFit will always perform the specified number of iterations.

Criterion	Meaning
Maximum iterations	Fitting will terminate when the specified number of iterations has been reached, irrespective of whether the designated change in χ^2_v has taken place. The maximum iterations limit serves to guard against the possibility that the calculations will never reach a minimum, and so forces control to return to the program eventually.

The settings can be stored as the defaults by clicking the **Store settings as defaults** button; GraFit will remember these defaults next time the program starts.


Results Pane

The Results pane holds a summary of the fitting results.

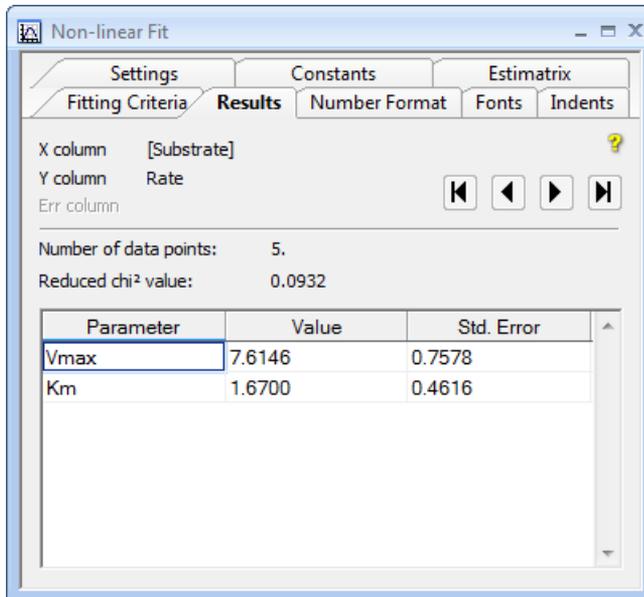

The Results pane of the Non-linear Fit item

If more than one set of data is specified in the Settings pane, each fit can be viewed by clicking the arrow buttons. The **File ▶ Print** and **File ▶ Print ▶ Print Preview** commands display the results as illustrated below:

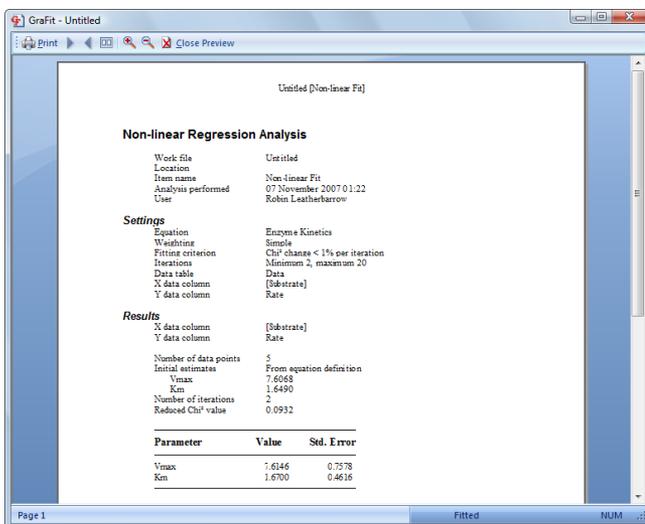

Using Print Preview from the Non-linear Fit item

The formatting of these results is controlled from the Number Format, Fonts and Indents panes.


Showing the Fitted Results on a Graph

A listing of the fitted parameter values can be displayed in a graph item window by selecting the **Results Listing** command³⁹ on the **Graph** tab. See Chapter 12 for further information.

³⁹ The **Draw ► Results Listing** command in legacy mode

Number Format Pane

The numeric result values are displayed formatted as is specified on the Number Format pane.

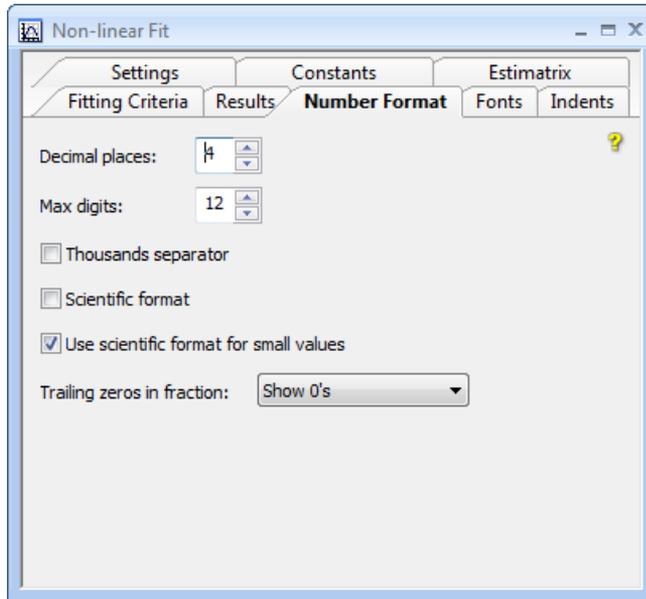

The Number Format pane of the Non-linear Fit item

Specification	Meaning
Decimal places	Enter the number of decimal places to be displayed for numbers in the results listing. <i>Note To ensure that all decimal places are displayed set Trailing zeros in fraction to "Show 0's."</i>
Max digits	Enter the maximum number of digits to use when displaying numbers in the results listing.
Thousands separator	Check this box to include a grouping separator. The symbol used for this separator, and the number of digits in the group, are set for the View ► Options command on the Home tab.
Scientific format	Check this box to display the values using scientific format (e.g. -2.435e003 instead of -2435.000).
Use scientific format for small values	Check this box to ensure that small values, which would otherwise be displayed as 0.0000, are shown in scientific format (e.g. 2.0987e-008).

Specification	Meaning
Trailing zeros in fraction	Chooses the way that trailing zeros are displayed in the formatted numbers. It is possible to show the "0" characters, to omit them, or to replace them with spaces.

Fonts and Indents Panes

These panes control the way that the results are formatted. The contents of these panes are common to several different item windows, and are described in Chapter 7.


Formatting Item Windows

Formatting Options

Many GraFit Item windows are built from a tabbed dialog box that allows various settings to be entered. The Linear Fit, Non-linear Fit, Equation, Hard-coded Equation and Transformation definition items are all built in this way. In each case, the settings in the window can be displayed using the **File ▶ Print ▶ Print Preview** command⁴⁰, and printed by the **File ▶ Print** command. Each of these item windows includes a Fonts and an Indents pane that control the way the information is displayed when printing.

Fonts Pane

The Fonts pane selects the fonts to use when displaying the contents of the item window.

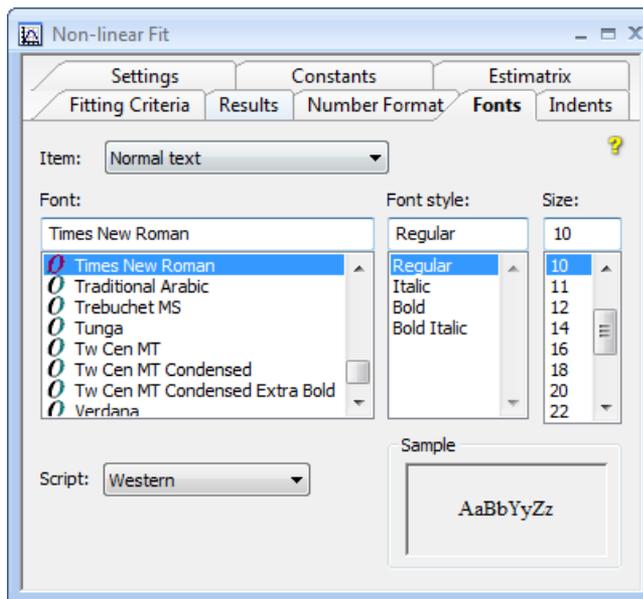


The Fonts pane of the Non-linear Fit item

⁴⁰ The **File ▶ Print Preview** command in legacy mode

This pane behaves in the same way as a standard Windows font selection dialog. The fonts to use can be set for a variety of different text elements, specified by the **Item** selection. To see the effect of changing the font settings, use the **File ▶ Print ▶ Print Preview** command⁴⁰.

Indents Pane

The Indents pane controls the organization of information for printing.

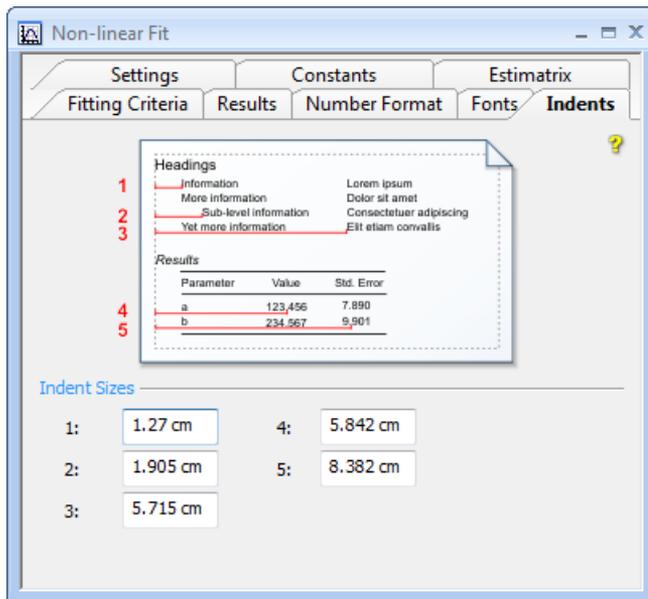

The Indents pane of the Non-linear Fit item

Each of the different indent levels can be individually set. To see the effect of changing the indent settings, use the **File ▶ Print ▶ Print Preview** command⁴⁰.

Linear Fit Item

Overview of Linear Fitting

A linear fit item is used to control data fitting by linear regression. This item specifies the data source and equation to use, and also presents the results of the data fitting. The item is created from the main GraFit window using the **Add ▶ Linear Fit** command. The linear fit window itself contains a tabbed dialog box with several panes that are used to specify how the fit is calculated. Each of these panes is discussed below. Data fitting involves making appropriate selections in the Settings pane (and possibly one or more of the other panes as described below), then selecting the **Calculate ▶ Fit Data** command⁴¹ on the **Home** tab. After fitting the data the results can be viewed from the Results pane.

Status Bar Indicators

When a fitting item is selected, the status bar shows the current fitting status, indicating whether the item is currently **Fitted** or **Not Fitted**. In addition, the status bar shows when a calculation is in progress. Whenever any calculation is being processed, the calculation pane shows an animated abacus.




Fitting status pane

Calculation pane

⁴¹ The **Calculate ▶ Fit Data** command in legacy mode

Settings Pane

The Settings pane is used to specify the main details of the fit.

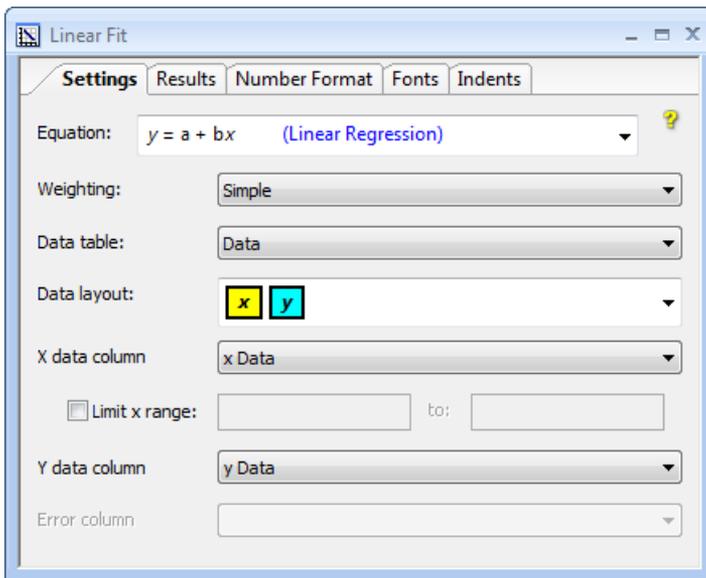

The Settings pane of the Linear Fit item

Equation

Choose the equation to use for the data analysis. All the available equations are linear, and are of the form given below:

$$y = a + bx + cx^2 + dx^3 + \dots$$

The default linear regression equation is $y = a + bx$.

Weighting

All types of regression analysis require information about the error distribution in your data. The reason for this is explained fully in Chapter 20, but briefly it is necessary to specify whether the inaccuracies are likely to be the same for each data point (“Simple” weighting), are proportional to the value of the y data points (“Proportional” weighting) or have been individually estimated (“Explicit” weighting). The default Simple weighting is often the most appropriate setting.



Data Table

The data to be fitted must be in a separate data table item within the GraFit workspace. Select the table that holds the data from this drop-down list.

Data Layout

The data layout specifies how the data are arranged for data fitting. The default layout, as shown in the figure above, is that the x and y data columns are to be selected individually.

This is used to fit a single $\{x, y\}$ data set. It is also possible to specify that multiple data sets are present, which will result in several sets of results being generated. The available options for a standard equation having a single independent variable are given below.

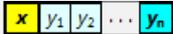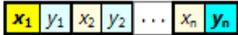Layout	Meaning
	A single data set is present and the columns holding the x and y data values are specified individually from the X data column and Y data column entries. The x and y data columns do not need to be in any particular order within the data table.
	Multiple data sets are present, although all data sets have common x data values. The x and y data columns are consecutive in the data table. The location of the x column is set from the X data column entry; the location of the final y column is set from the Last Y data column entry. The number of data sets present is determined automatically.
	Multiple data sets are present, and the data sets have different x data values. Pairs of x, y data columns are present consecutively in the data table. The location of the first x column is set from the First X data column entry; the location of the final y column is set from the Last Y data column entry. The number of data sets present is determined automatically.

X Data Column

This entry specifies the column that holds the x data values. See the Data Layout section (above) for further information.

Limiting x Data Range

It can sometimes be useful to restrict data fitting to a subset of the values in the columns. For example, it can be useful to fit just the initial portion of a data set. Checking the **Limit x range** box allows a lower and upper limit to be set on data that will be included in the analyses.

Y Data Column


This entry specifies the column that holds the y data values. See the Data Layout section (above) for further information.

Error Column

If explicit weighting is chosen, the error values are read from a separate column, which is selected here.


Results Pane

The Results pane holds a summary of the fitting results.

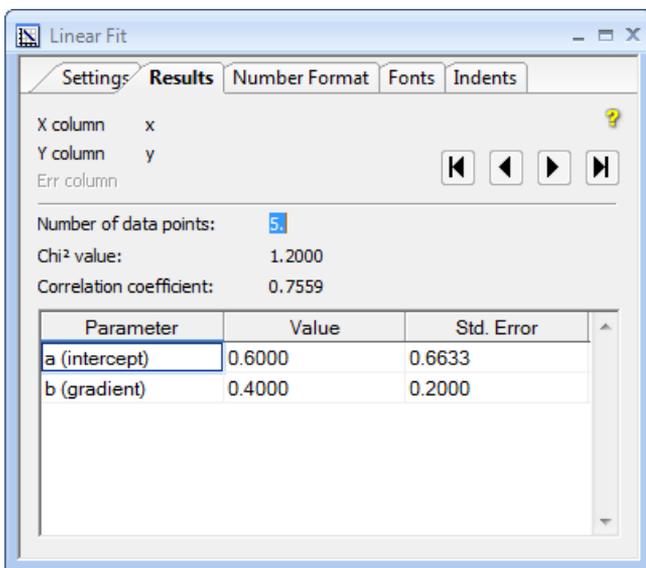

The Results pane of the Linear Fit item

If more than one set of data is specified in the Settings pane, each fit can be viewed by clicking the arrow buttons. The **File ▶ Print** and **File ▶ Print ▶ Print Preview** commands show the results in a more stylized manner, as illustrated below:

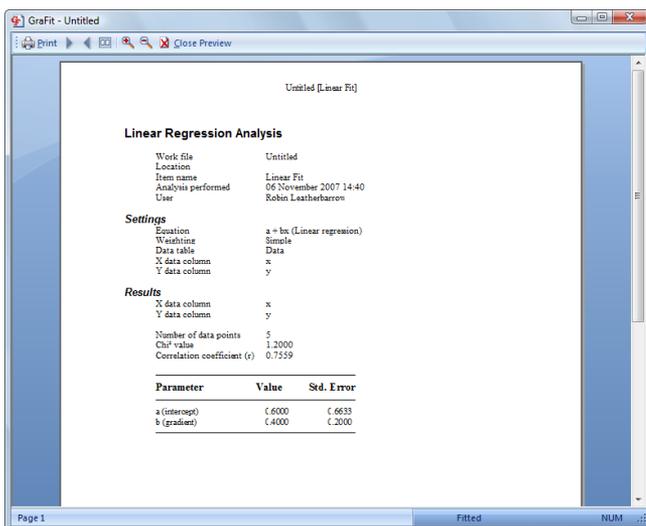

Using Print Preview from the Linear Fit item

The formatting of these results is controlled from the Number Format, Fonts and Indents panes.

Showing the Fitted Results on a Graph

A listing of the fitted parameter values can be displayed in a graph item window by selecting the **Results Listing** command⁴² on the **Graph** tab. See Chapter 21 for further information.

Fonts and Indents Panes

These panes control the way that the results are formatted. The contents of these panes are common to several different item windows, and are described in Chapter 7.

⁴² The **Draw** ► **Results Listing** command in legacy mode

Equation Item

Overview of Equation Definitions

Equation definitions are used for non-linear fitting (see Chapter 6 for more information about non-linear fitting items). The item is created from the main GraFit window using the **Add ▶ Equation** command. The equation window itself contains a tabbed dialog box with several panes that are used to define the equation. After defining the equation, it must be compiled using the **Compile** command⁴³, found in the **Equation** group of the **Home** tab, before it can be used. Each of the panes in this item window is discussed below.

Status Bar Indicators

When an equation definition item is selected, the status bar shows whether the equation is **Compiled** or **Not Compiled**.




Compilation status pane

⁴³ The **Equation ▶ Compile** command in legacy mode

General Pane

The general pane allows a description of the equation to be entered.

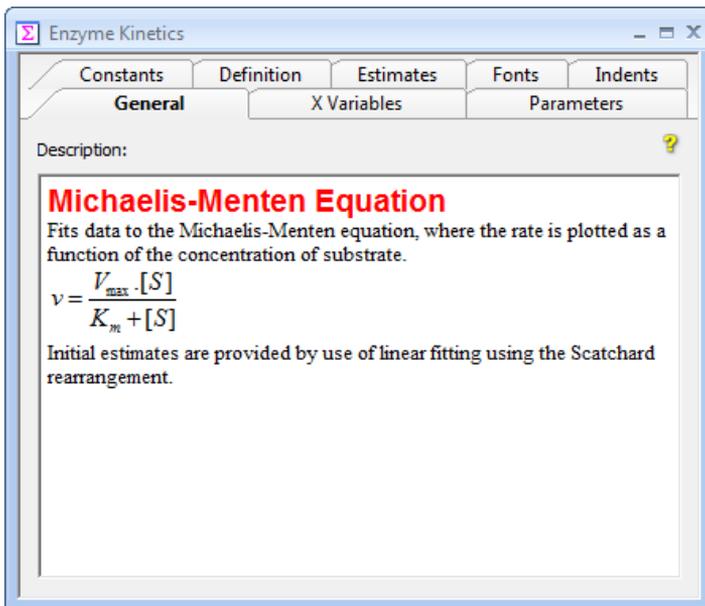

The General pane of the Equation item

The description is optional, but adding a text description of your equation definitions is strongly recommended. As illustrated in the figure above, the descriptive text allows rich formatting options, with multiple font types, styles, sizes and colors allowed. It is also possible, as shown above, to include an OLE object embedded in the text. In this particular example, a Microsoft Equation object has been embedded within the text.

Note The Microsoft Equation editor, which comes with Microsoft Word, is not provided with GraFit.

Changing Font Styles


To alter the style of a selected region of text, use the commands available on the **Font** tab⁴⁴.

Inserting OLE Objects

To insert an OLE object into the description text, position the cursor where you wish the object to be placed, and select the **Object ► Insert New Object** command⁴⁵ from the **Home** tab.

⁴⁴ The **Edit ► Font** command in legacy mode

⁴⁵ The **Edit ► Insert New Object** command in legacy mode

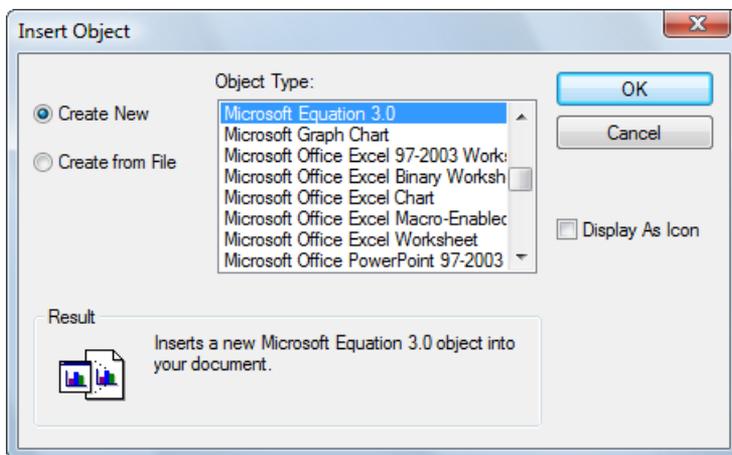

Choose the object type to create, and then click **OK**.


Editing OLE Objects

An inserted OLE object can be edited by double-clicking on the object, or by choosing the appropriate menu command from the **Object** menu⁴⁶ on the **Home** tab.

X Variables Pane

The X Variables pane is used to define the names of the independent (x) variables that are used in the equation. All equations must have at least one independent variable, and can include more than one.

⁴⁶ The **Edit** menu in legacy mode

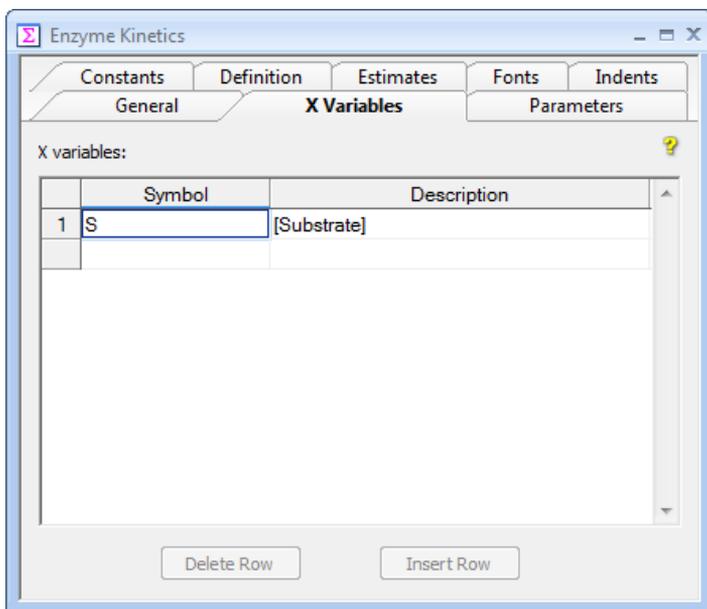

The X Variables pane of the Equation item

Naming Rules

All name definitions for use in an equation definition are entered into a grid that has space for a Symbol and a Description entry. All names must have a Symbol entry, but the Description field is optional.

Description Entry

If a Description entry is made, GraFit will use this long name whenever it refers to the name, for example when listing results. If no Description entry is present, GraFit will use the Symbol entry. Description entries can incorporate any characters, including spaces and punctuation characters that are illegal in Symbol entries.

Symbol Entry

The symbol entry is the name used in the equation Definition page. Symbol names must obey a number of rules, in a similar way to variable names within a computer program.

- The name must begin with a letter, and contain only letters and numbers. In particular, a symbol name cannot contain the space character.

The following names are legal symbols:

```
Fred
k_on
x23
```

The following names, however, are invalid:

```
2nd          (starts with a number)
k off       (contains a space)
```

test-2	(contains a minus sign)
ten:	(contains a punctuation character)

- No distinction is made between upper and lower case characters within a symbol definition. Therefore the symbol names `grafit`, `GraFit`, `GRAFIT` and `gRAfIT` are all considered identical.

We recommend using meaningful names wherever possible when defining symbol names. Using `Final` to represent the final concentration of a substance rather than `z` will make understanding your definitions far easier for others as well as yourself. In addition, it is often useful to add a full descriptive name.

Number of X Variables


After entering an x variable definition, the data entry grid will show a blank entry position immediately below the defined variable. If you enter a second definition in this location, a further blank entry position will appear. The number of rows of information that are entered into the grid therefore defines the number of x variables in the equation.

When to use Multiple X Variables

Most equations are of the form $y = f(x)$, and so have just a single independent (x) variable. In equations with multiple independent variables, the observed value (y) varies as a function of more than one variable condition. For example, the rate of an enzyme reaction may vary as a function of both substrate and inhibitor concentrations. In this case, there are two independent (x) variables. In principle it is possible to have any number of independent variables, although in practice it is rare to fit data where more than two independent variables are present.

Parameters Pane

The Parameters pane is used to define the parameter values in the equation definition. The parameter values are the “unknowns” in the equation, and non-linear regression allows data to be fitted in order to find optimal values for these parameters.

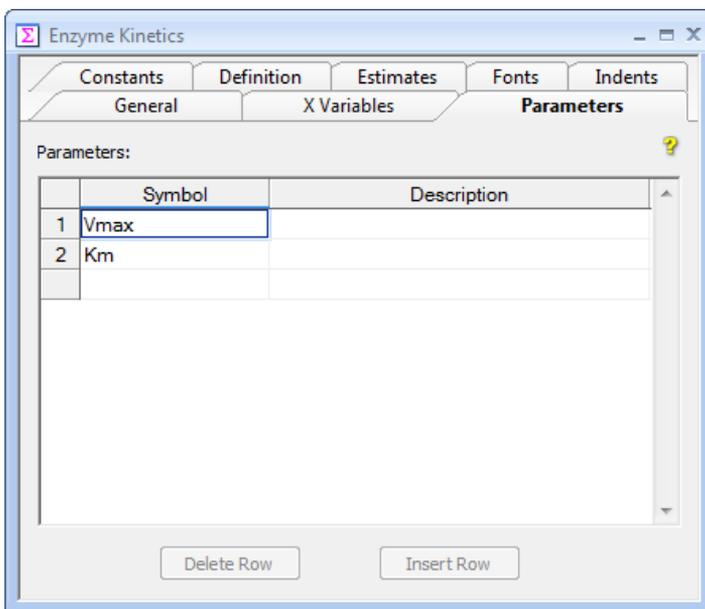


The Parameters pane of the Equation item

The organization of the Parameters pane is very similar to that of the X Variables pane, and the definition of Symbol and Description names is performed in an analogous manner (see the X Variables Pane section above). The number of rows of information that are entered into the grid defines the number of parameters in the equation. At least one parameter must be defined. There is no formal limit to the number of parameters allowed in a GraFit equation definition, but it is rare to find equations that involve more than four or five parameters.

Constants Pane

The Constants pane is used to define constants in the equation definition. Constants are used to allow the entry of a value into an equation definition that is fixed in any single experiment, but varies from experiment to experiment. The presence of constants is optional, but when present an equation definition can include any number of constants. When constants are included in an equation used for non-linear regression, their values are entered into the non-linear fit item window as described in Chapter 6.

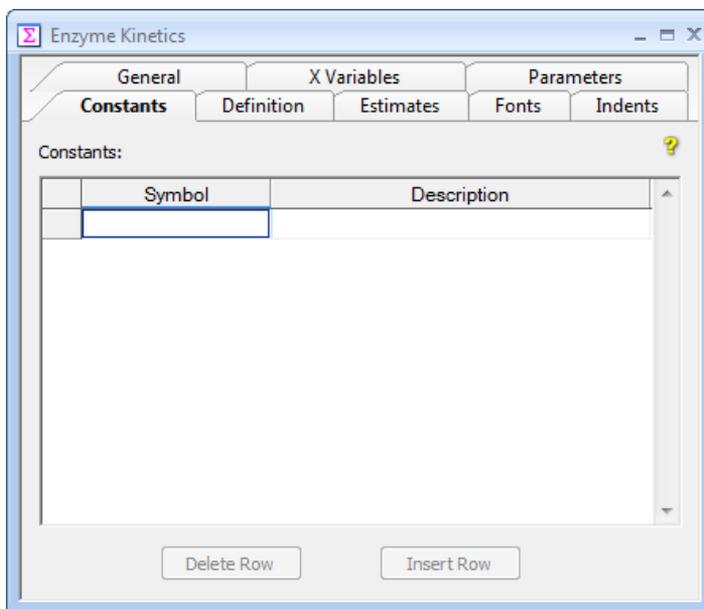


The Constants pane of the Equation item

The organization of the Constants pane is very similar to that of the X Variables pane, and the definition of Symbol and Description names is performed in an analogous manner (see the X Variables Pane section above). The number of rows of information that are entered into the grid defines the number of parameters in the equation. If no constants are required, as will be the case for many equations, this page should be left blank.

Definition Pane

The definition pane describes the functional form of the equation using the Symbol names of the x variables, the parameters and constants (if any) entered as on the earlier panes.

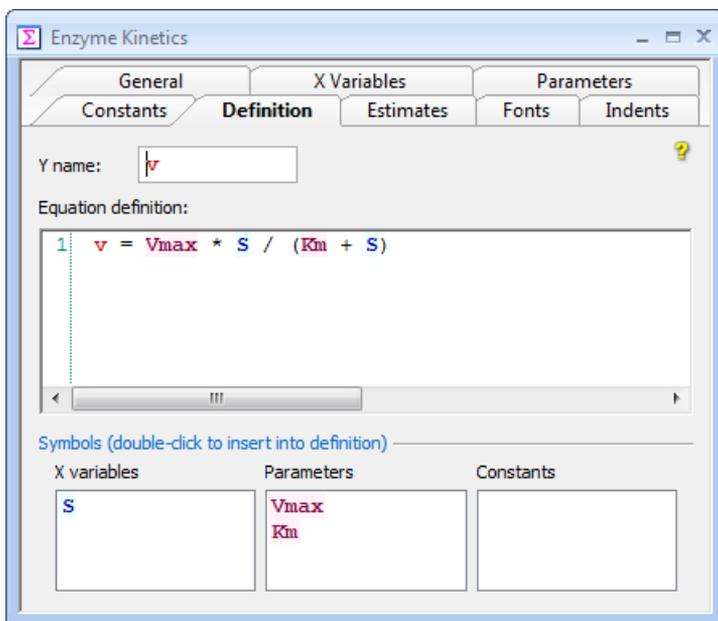

The Definition pane of the Equation item

All equation definitions that are used within GraFit must have the functional form $y = f(x)$, i.e. it must be possible to describe the variation of y as some defined function of x . The information on this pane describes the required equation in symbolic form. At the top of the pane is an entry for **Y name**, which allows the y axis to be given a Symbol name (using the rules described earlier).

Equation Definition Syntax

The equation definition is built up from the symbols representing parameters and constants (if any) with whatever arithmetic operators and mathematical functions are required. For the example shown in this Chapter, we have used the following equation:

$$v = \frac{k_{cat} \cdot [E_0] \cdot [S]}{K_m + [S]}$$

In the **Equation definition** entry, this becomes:

$$v = kcat * Eo * S / (Km + S)$$

The syntax used is similar to that of a programming language, and uses the rules described below.

Arithmetic Operators

Symbol	Operation	Example	Explanation
+	Addition	$x + y$	Add x and y
-	Subtraction	$lim - off$	Subtract off from lim
*	Multiplication	$m * x$	Multiply m by x
/	Division	$a/b52$	Divide a by $b52$
^	Exponentiation	23^N	Calculate 23^N
=	Assignment	$y = x*2$	Assign y to be twice the value of x

Logical Operators


Symbol	Operation	Example
>	Is greater than	$k > 7$
<	Is less than	$x_value < con1$
>=	Is greater than or equal to	$test \geq 23.87$
<=	Is less than or equal to	$rate \leq maxrate$
=	Is equal to	$i76 = 0$

Note The "=" symbol is used for both assignment and as a logical operator. Which action is taken depends on the context.

Logical operators allow a comparison to be made. The result of the comparison is either 0 if FALSE, or -1 if TRUE. To see how these operators work, consider a graph of the form:

$$y = (x < test) * -1 + 3$$

Here, x is the x axis, and $test$ is a constant that has the value 2. This function is plotted below.

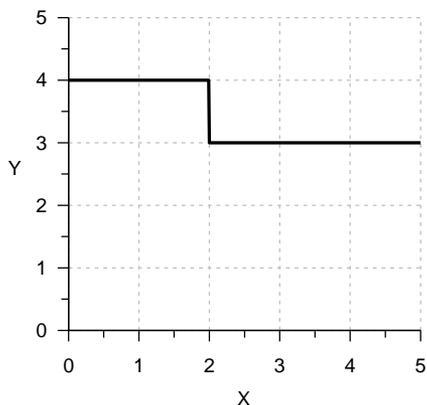

These comparisons can be used to introduce a discontinuity into a curve, and though they will seldom be used in equations, are extremely useful for displaying limits etc. to annotate a graph. An alternative method to obtain discontinuities is to use the `If Then Endif` constructs described below.

Note To obtain a sharp discontinuity it is necessary to create the curve with **Resolution** set to **High** (see page 163).

Numbers can be entered directly into the equation definitions, and can include exponents, e.g. $27.87e^{-12}$.

Functions

A large number of functions are provided for use in equation definitions.

Function	Returns
<code>abs(x)</code>	Absolute value of x , i.e. if $x < 0$ returns $-x$, otherwise returns x .
<code>acos(x)</code>	Cos^{-1} of x .
<code>alog(x)</code>	Antilogarithm to base 10 of x .
<code>asin(x)</code>	Sin^{-1} of x .
<code>atan(x)</code>	Tan^{-1} of x .
<code>cos(x)</code>	Cosine of x .
<code>exp(x)</code>	e^x .
<code>frac(x)</code>	The fractional part of x . The fractional part of 13.387 is 0.387.
<code>int(x)</code>	The integer part of x . The integer part of 13.387 is 13.
<code>log(x)</code>	The logarithm to base 10 of x .

Function	Returns
<code>ln(x)</code>	The natural logarithm of x .
<code>sin(x)</code>	Sine of x .
<code>sqr(x)</code>	The square of x , x^2 .
<code>sqrt(x)</code>	The square root of x , \sqrt{x} .
<code>tan(x)</code>	Tangent of x .

The expression within the brackets may be of arbitrary complexity, for example

```
exp(-k * sqrt(j - 1))
```

All the geometric functions (sin, cos etc.) use radians as their units.

The following predefined constants are also available:

Constant	Value
<code>_pi_</code>	$\pi = 3.1415926536$
<code>_e_</code>	$e = 2.7182818285$
<code>_rnd_</code>	A random number between 0.0 and 1.0.
<code>_grnd_</code>	A random number with Gaussian error distribution. The number has a zero mean and unit variance.

Multi-line Equation Definitions

In its simplest form as described earlier, the equation definition is a single line that calculates the y value. For example, if the equation to be used was $y = ax + b$, then the equation definition might be written as:

```
y = a * x + b
```

Here, a and b have been declared to be the symbols that represent the parameters, x has been declared to be the symbol that represents the x data and y the symbol that represents the y data. Note that the multiplication between a and x needs the explicit inclusion of the multiplication operator, $*$. The equation definition can also be defined over several lines in order to simplify a complicated definition. This allows the calculation to be broken up into several steps.

Use of Temporary Variables

Temporary variables can be used to hold intermediate calculation results and simplify definitions. For example, the following definition uses three temporary variables, a , b and c .

```
a = 1
b = Kd + total + Capacity
c = total * Capacity
y = -( -b + sqrt( sqr(b) - 4 * a * c ) ) / ( 2*a )
```

The remaining symbols (**Kd**, **total**, **Capacity**) are defined in the equation editor as parameters or constants. Values are assigned to the temporary variable using the assignment operator, =, after which they may be used in expressions as normal. Use of temporary variables makes the above definition more easily read, and speeds up the calculation as the value for **b** need only be calculated once. Each line of the definition is separated from the next by pressing the RETURN key.

Conditional Statements

The values that are evaluated can be controlled using the **If** statement. This takes the form

```
If condition Then
.....statements.....
Endif
```

or alternatively,

```
If condition Then
.....statements.....
Else
.....statements.....
Endif
```

For example, it is possible to return different values to the program depending upon the value of a parameter, constant or independent variable:

```
If extra=1 Then
  Return( a1 + b*x )
Endif
If extra=2 Then
  Return( a2 + b*x )
Endif
If extra=3 Then
  Return( a3+ b*x )
Endif
Return( a4+ b*x )
```

Comments

Comments can be added to annotate the equation definitions.

```
' Anything after an apostrophe is a comment
; as is anything after a semi-colon
```

Looping Statements

The `For` statement has the following syntax.

```
For varname = startValue To endValue
.....statements.....
Next
```

or alternatively

```
For varname = startValue To endValue Step stepValue
.....statements.....
Next
```

Goto Statement

Another statement that will be familiar to programmers is the `Goto` statement. In GraFit equations this is implemented in the following manner:

```
:labelName
.....statements.....
Goto labelName
```

`labelName` can be any valid name, using the same rules as apply to the names for symbol names. The label may come before or after the `Goto` statement.

Note Indiscriminate use of `Goto` could result in an infinite loop.

Estimates Pane

The Need for Starting Estimates

Non-linear regression works by taking initial estimates for the best-fit parameter values, and by an iterative process calculates successively better values until the best-fit values are reached. It is therefore necessary to provide the program with suitable starting estimates, which should generally be of the same order of magnitude as the actual ones. If an equation definition alone is provided, the program cannot know what initial estimates are appropriate. In this event, when the equation is used the program will ask you to enter initial estimates before fitting the data.

To relieve users the inconvenience of estimating parameter values, GraFit provides a means of generating these estimates automatically. It should be noted that the method used is not foolproof, and requires familiarity with manipulating equations. It does, however, make data fitting simpler, and can be over-ridden if the estimates provided are inappropriate. In addition, equations can be defined without provision for automatic initial estimation.

Method Used to Obtain Initial Estimates

The method used is based on performing linear regression on rearranged data. It is often possible to rearrange data in such a way that a linear plot is obtained, with the parameter

values being related to the gradient and intercept. As linear regression needs no initial estimates, this can be used to estimate these parameter values providing the program knows the following.

1. The rearrangement that is needed.
2. How the parameter values relate to the gradient and intercept.

This process is illustrated using the previous example, where the equation is:


$$v = \frac{k_{cat}[S][E_0]}{K_m + [S]}$$

We can rearrange this to a linear form by plotting $v/[S]$ against v :

$$\frac{v}{[S]} = \frac{k_{cat}[E_0]}{K_m} - \frac{v}{K_m}$$

This is the so-called Eadie-Hofstee equation. Plotting $v/[S]$ against v allows the slope and intercept to be found, from which it is possible to calculate k_{cat} as $(-\text{intercept}/\text{gradient})/[E_0]$, and K_m as $-1/\text{gradient}$. It must be stressed that this approach does *not* provide strictly correct values for k_{cat} and K_m —this is why we wish to use non-linear regression in the first place! (Further information can be found in Chapter 19.) The values provided are suitable, however, to use as initial estimates for the non-linear regression analysis. (It is instructive to compare the initial estimates provided by this method with the final values produced by the non-linear regression analysis. As they are demonstrably different, this illustrates the danger of using linear regression of a linearized function for analyzing data.)

The Estimates Pane

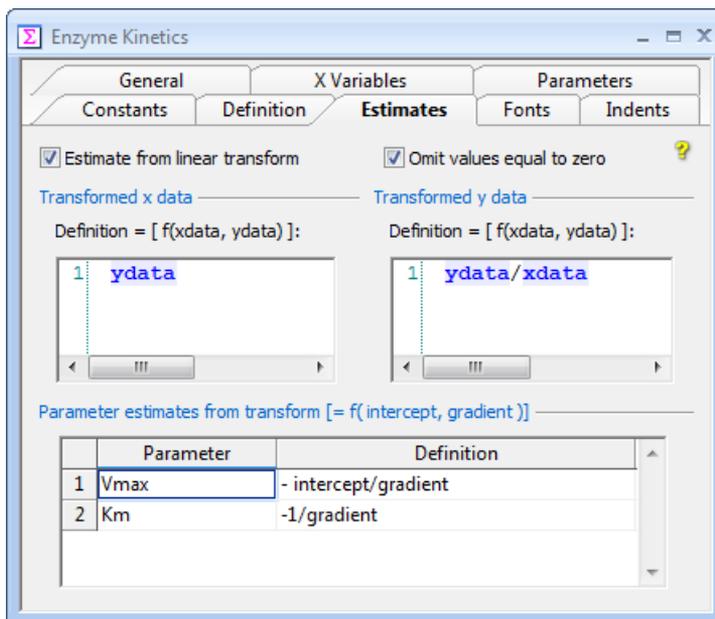

The Estimates pane of the Equation item

To specify that initial estimates are provided for an equation, the **Estimate from linear transform** box should be checked. If you do not want the program to calculate initial estimates, simply leave the box unchecked.

Defining How to Rearrange the Data Points

For our current equation, the rearranged plot that we require is of $v/[S]$ against v , i.e., y/x versus y . The **X Data** and **Y Data** entries are used to define the data points for such a plot (note that the stored data points in the data table are, of course, unaffected). For this equation the definitions are:

$ydata$

and

$ydata/xdata$

respectively. The symbols “ $xdata$ ” and “ $ydata$ ” are predefined by the program to represent the x and y data points.


These definitions allow the program to perform the rearranged plot, and any valid single line expression can be used. The grid at the bottom of this page is used to define how we obtain our initial estimates from linear regression on the rearranged plot. Shown in the previous figure is the definition of k_{cat} as $(-intercept/gradient)/E_0$, and K_m as $-1/gradient$. Use the scroll bar to bring any further parameters into view. The symbols “ $intercept$ ” and “ $gradient$ ” are predefined, and represent the intercept and gradient of the rearranged linear plot.

Other predefined symbol names may also be used when defining initial estimates of your variables. The symbol `min` will give the value of the minimum y data point, and `max` the value of the maximum y data point. The corresponding minimum and maximum x data points are `minx` and `maxx`, respectively. These values can be of benefit if you need to obtain rough estimates of limits, where the limiting value can be approximated by the maximum recorded value. It is also possible to enter numeric values directly as estimates. In some situations, rough numerical values provide suitable estimates without recourse to the more complex methods described above.

If the **Omit values equal to zero** box is checked then any data values (x or y) that are equal to zero will not be used during the data transformation. This is useful with several transformations that involve taking logarithms or dividing by one of the data values (as here). If this box is not checked then a data point which happens to have a value of zero might cause an arithmetic error, resulting in the estimates routine failing.

Compiling the Equation

To compile the equation, choose the **Compile** command⁴³ from the **Equation** group on the **Home** tab. Compilation translates your definition into an internal form that can be executed rapidly. During compilation the definition is comprehensively checked to ensure that the correct syntax has been used. If an error is detected a message box will be displayed informing you of the mistake; after pressing the **OK** button to continue the cursor will be placed at the most likely place of error.

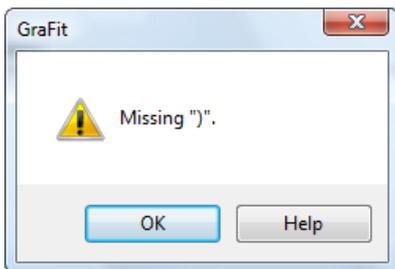

An Error Message Box


To obtain further information on the specific error, click the **Help** button, or press the **F1** key.

Once the equation has been correctly compiled, it can be used for data fitting and/or graph plotting.

Transformation Item

Overview of Transformation Definitions

Transformation definitions are used in order to allow data and curves to be displayed on a graph after application of a user-defined transformation. For example, one use of a transformation definition would be to allow a double-reciprocal plot to be displayed. Transformations describe methods for rearranging data and equations that can be applied during graph drawing without the need to create a separate set of rearranged data.

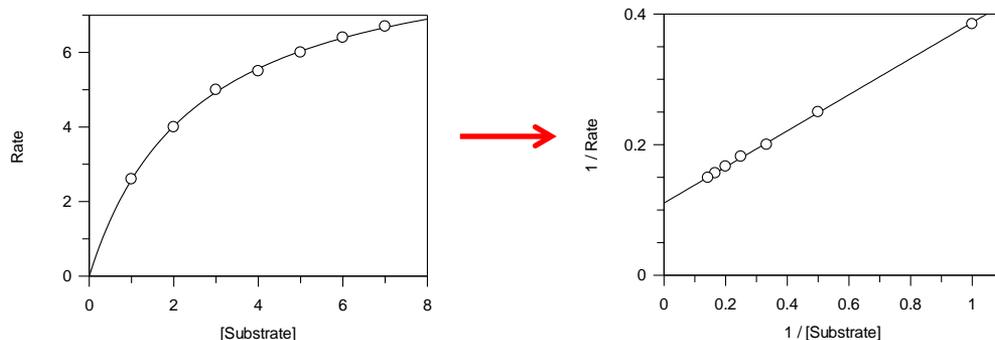

An example of a transformed double-reciprocal plot. The transformation definition defines the way that the data points and the fitted curve are rearranged

The transformation item is created from the main GraFit window using the **Add ▶ Transformation** command. The transformation window itself contains a tabbed dialog box with several panes that are used to define the transformation. After defining the transformation, it must be compiled using the **Compile** command⁴⁷, found in the **Transform** group of the **Home** tab, before it can be used. Each of the panes in this item window is discussed below.

⁴⁷ The **Transformation ▶ Compile** command in legacy mode

Status Bar Indicators


When an equation definition item is selected, the status bar shows whether the equation is **Compiled** or **Not Compiled**.



Compilation status pane

General Pane

All transformations are associated with an equation definition item, which must also be present in the GraFit workspace. See Chapter 9 for more information about defining equations. The general pane also allows a description of the transformation to be entered.

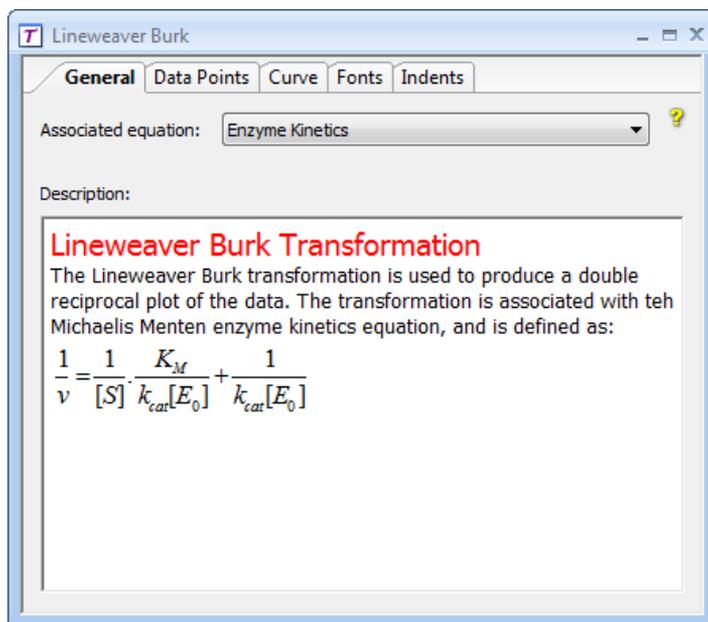


The General pane of the Transformation item

The description is optional, but adding a text description of your transformation definitions is strongly recommended. As illustrated in the figure above, the descriptive text allows rich formatting options, with multiple font types, styles, sizes and colors allowed. It is also possible, as shown above, to include any OLE object in the text. In this particular example, a Microsoft Equation 3.0 object has been embedded within the text.

Data Points Pane

The Data Points pane describes how data points are transformed for plotting. If the **Transform data points** box is unchecked, no data point transformation is used. Otherwise,

the remainder of the window describes how the data are to be rearranged when this transformation is applied.

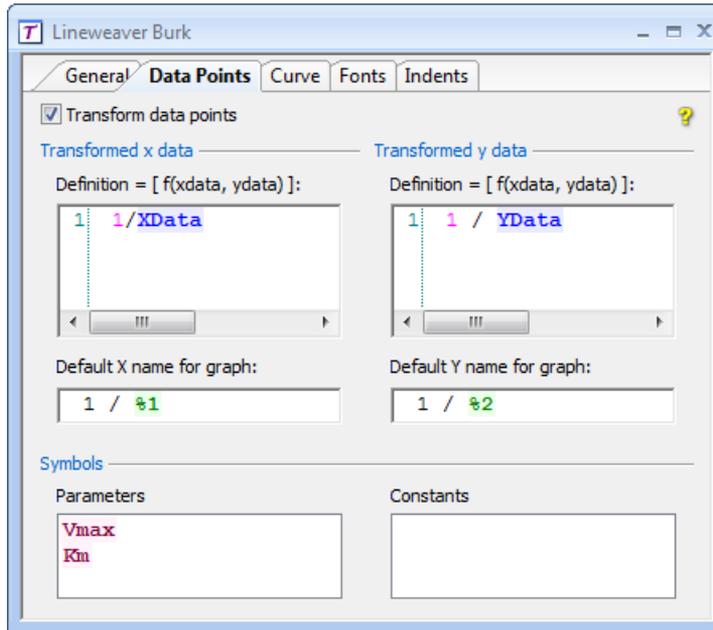

The Data Points pane of the Transformation item

In the sample shown above, the data transformation involves plotting the reciprocals of the x and y data.

Data Transformation Definitions

Two transformation definitions need to be entered, one each for the transformed x and y data. Use the edit fields to enter the transformation definitions. The rules for these transformation definitions are the same as described in Chapter 9 for equation definitions, although the definition is limited to a single line and only the right hand side of the “Definition =...” equation is entered. The definition can be of arbitrary complexity, and can include mathematical functions, parameter names etc. The pre-defined symbols $xdata$ and $ydata$ represent the original x and y data values respectively. The symbol $xdata2$ can be used to access any data values associated with the y data column (“ y Data with Associated Value” column type—see the discussion of column types in Chapter 5). In the example illustrated above we wish to take reciprocals of these data; the definitions are therefore $1/xdata$ and $1/ydata$. Any algebraic expression may be entered, for example

$$xdata / (ydata + \text{sqrt}(xdata))$$

is a perfectly valid transformation definition. It is permissible to access the values of parameters and constants in the equation, if desired. The symbol names of these are listed at the bottom of the Data Points pane.

Generating Residuals Plots

A very useful symbol that is specific to the data transformation definition is `ycalc`. This returns the y value that is calculated from the fitted parameters. Therefore, to generate a residuals plot, define the transformed y data as:

$$ydata - ycalc$$


Axis Name Settings

When first creating a transformed plot, the text used to name the x and y axes is taken from the **Default X name for graph** and **Default Y name for graph** entries. The following special entries can be used:

Entry	Effect
%1	All occurrences of this are replaced by the name of the untransformed x data column.
%2	All occurrences of this are replaced by the name of the untransformed y data column.

Curve Pane

The Curve pane describes how fitted curves are transformed for plotting. If the **Transform curve** box is unchecked, no curve transformation is used. Otherwise, the remainder of the window describes the transformed curve definition.

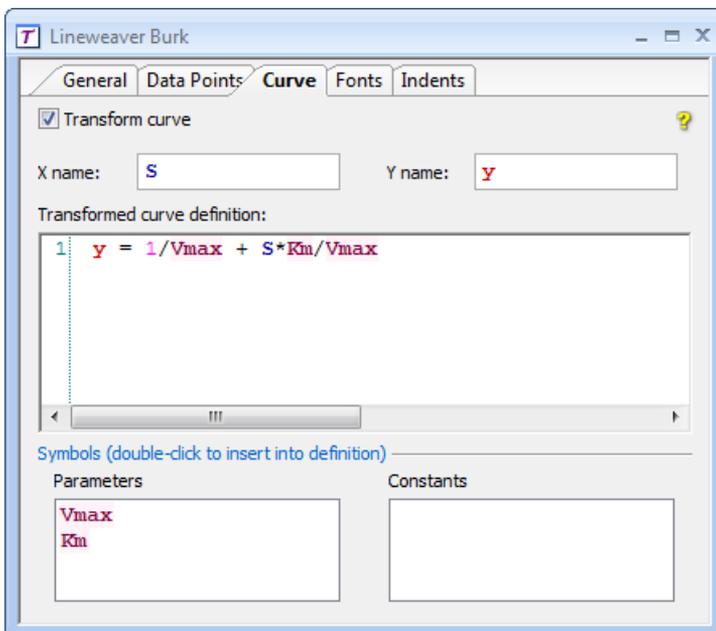

The Curve pane of the Transformation item

In the sample shown above, the curve transformation displays a line through the transformed data points making use of the fitted parameter values

Transformed Curve Definition

To draw the transformed curve (which in the example above gives the linear Lineweaver Burk plot) it is necessary to have an appropriate equation definition. For the example that we are using, the transformed equation is:

$$\frac{1}{v} = \frac{1}{[S]} \cdot \frac{K_m}{k_{cat}[E_0]} + \frac{1}{k_{cat}[E_0]}$$

What is being plotted is $1/v$ (the y axis) versus $1/[S]$ (the x axis). To show this line, we therefore need an equation that is defined as:

$$y' = x' \cdot \frac{K_m}{k_{cat}[E_0]} + \frac{1}{k_{cat}[E_0]}$$


Here, x' and y' are the transformed data. The transformed curve definition uses the **X name** and **Y name** entries, which correspond to x' and y' in the above equation. The syntax used for this definition is as described for equations (see Chapter 9), and can occupy more than one line if required.

Fonts and Indents Panes

These panes control the way that the results are formatted. The contents of these panes are common to several different item windows, and are described in Chapter 7.

Compiling the Transformation

To compile the equation, choose the **Compile** command⁴⁷ from the **Transform** group on the **Home** tab. Compilation translates your definition into an internal form that can be executed rapidly. During compilation the definition is comprehensively checked to ensure that the correct syntax has been used. If an error is detected a message box will be displayed informing you of the mistake; after pressing the **OK** button to continue the cursor will be placed at the most likely place of error.

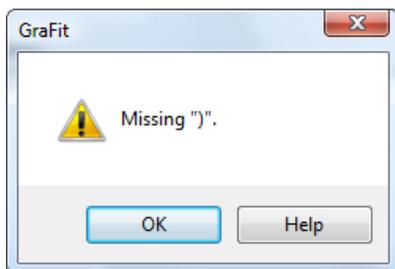

An Error Message Box

To obtain further information on the specific error, click the **Help** button, or press the F1 key.


Once the transformation has been correctly compiled, it can be used for graph plotting.

Notepad Item

Uses of Notepad

The notepad provides a text editor window where information can be stored in the GraFit workspace. Uses include making notes about the data, the analysis of the results, and writing brief documentation on how to use a particular data file or template. The notepad item is created from the main GraFit window using the **Add ▶ Notepad** command.

The Notepad Item Window

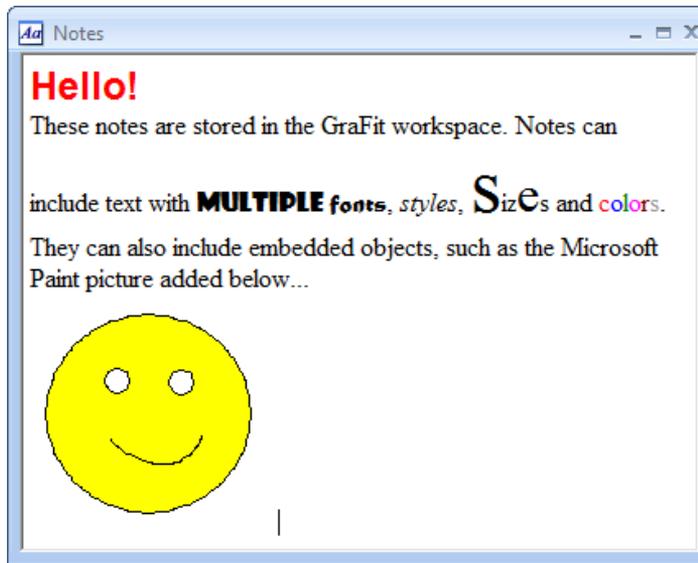

The Notepad item window

The font styles are set from the **Font** tab⁴⁸. Notepad windows automatically wrap their contents so that the line width is equal to the width of the window.

Alternatives to the use of Notepad Items

The notepad item provides a simple text editor that has many of the capabilities of the Windows WordPad application. However, it does not provide the wealth of features present in a dedicated word-processing application such as Microsoft Word. It is, however, possible to

⁴⁸ The **Font** menu, or the formatting toolbar, in legacy mode

harness the power of applications such as Word by embedding documents into the GraFit workspace. For example, a Word document can be embedded within the main GraFit window, as described in Chapter 18. Incorporating notes in this manner can be a more suitable alternative to using a notepad item.

Creating Graphs

How to Create a New Graph

To create a new graph, first add a graph item window to the workspace using the **Add ▶ Graph** command. You will then be given a choice of graph type to create from the **Select Graph Type** dialog (below). The resulting graph will be created using the default size. If you want a graph of a specific size, press the **Cancel** button, which will result in a blank graph page. Then use the New Graph tool to add a graph at a specific location and size to an existing graph item window.



New Graph tool


Adding Graphs to an Existing Graph Item

It is possible to create as many graphs as required on each graph item page. To add an additional graph to an existing graph item window, either select the **New Graph** command⁴⁹ from the **Graph** tab, which will create a graph of a default size, or click the New Graph tool and drag to position the new graph where you wish it to be placed.

Available Graph Types

Although GraFit can create a large variety of different graphs, the initial choice when making a new graph is between four fundamental graph types. These differ in the data that they use; for example, X/Y Scatter Graphs hold data with separate numeric x and y values whereas a Column Chart has data values that associated with named categories.


⁴⁹ The **Graph ▶ New Graph** command in legacy mode

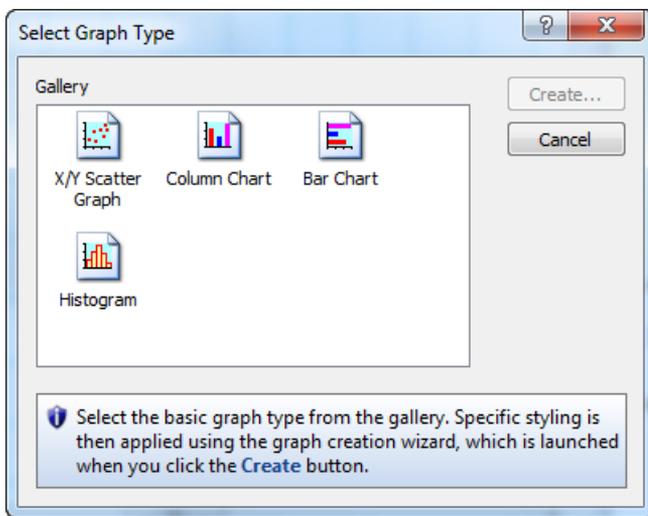


After selecting the type of graph to create, clicking the **Create** button will result in a specific graph creation wizard that allows the graph to be customized. A description of the kind of graphs that can be created for each type is given below, and more detailed information is found in the following chapters.

X/Y Scatter Graph


X/Y Scatter Graphs have numeric ranges on both x and y axes. Some of the types of graph that can be created are illustrated below.

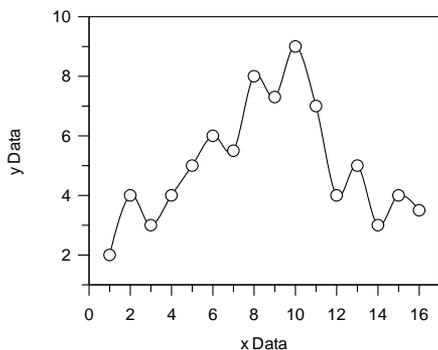


This graph shows data points joined by a smooth Bezier curve

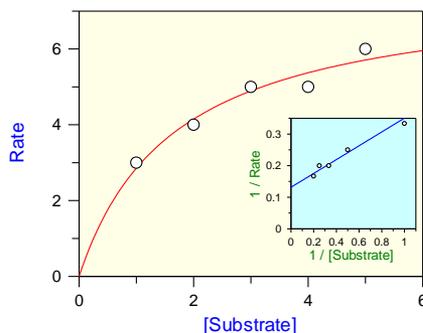


This graph shows data and a fitted curve; inset is a linear transformation of the data

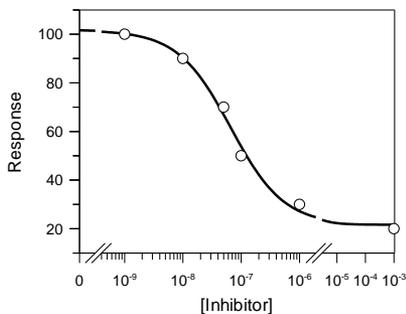


This graph shows a semi-log plot with a doubly split axis scale

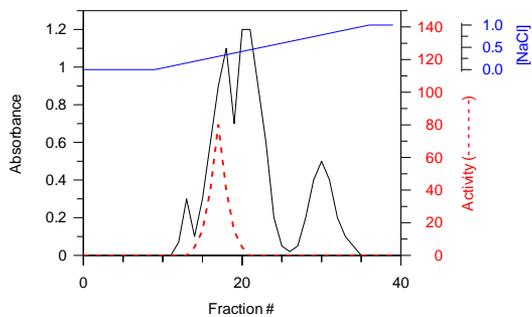

This graph shows a single graph having multiple axis scales

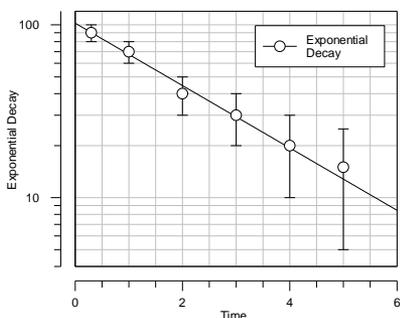


This graph shows offset axes and error bars

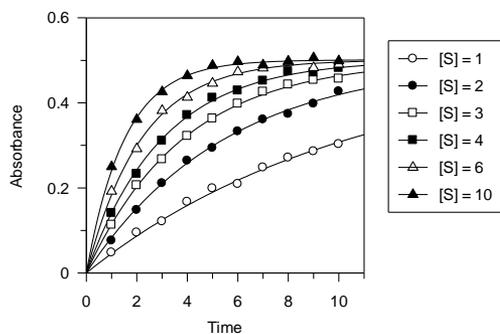


This graph shows multiple data sets fitted simultaneously to an equation with two independent variables

Column Chart


Column Charts are categorized by having a numeric range on the y axis and named categories on the x axis. Some of the types of chart that can be created are illustrated below.

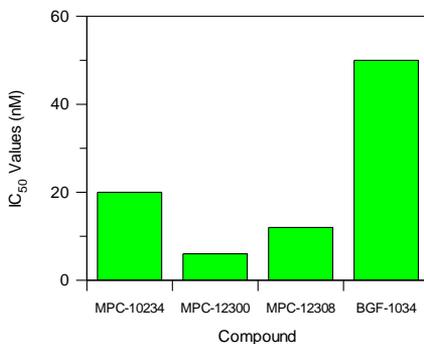


Simple comparison of a set of parameter

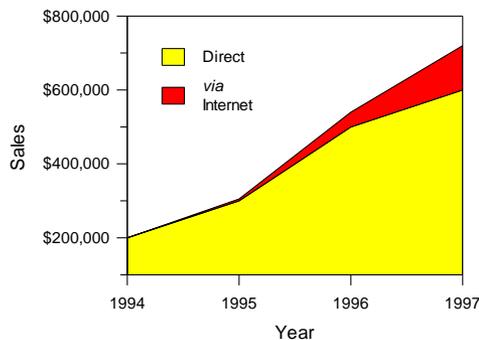

Summed area plot

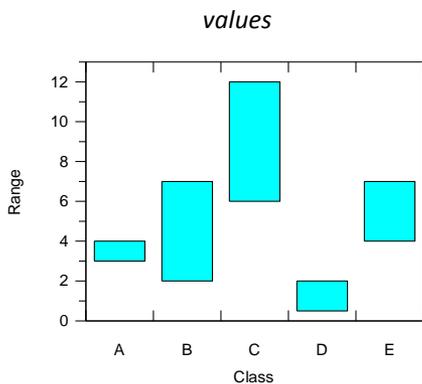


A Low/High column chart

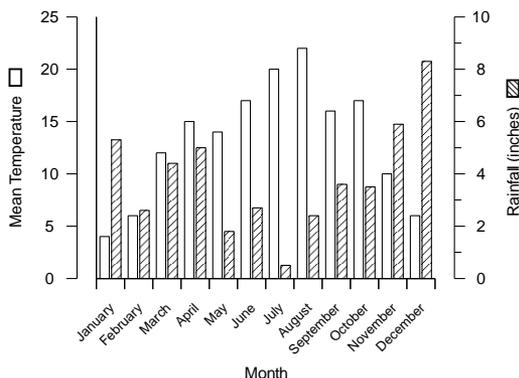


A chart with different y axis scales for the two data sets plotted.

Bar Chart

Bar Charts are categorized by having a numeric range on the x axis and named categories on the y axis, and have the same range of styles available as column charts. Some of the types of chart that can be created are illustrated below.

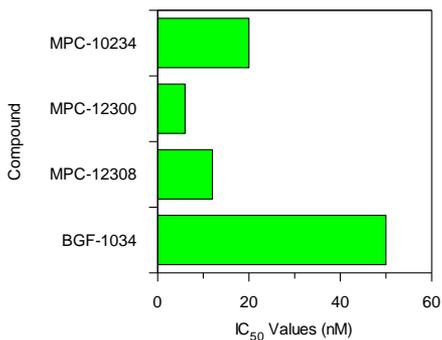

Comparison of values using a bar chart

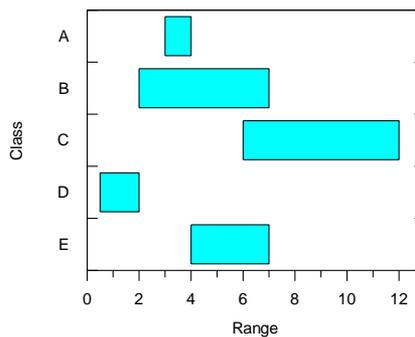


A Low/High bar chart


Histogram

Histograms have numeric ranges on both x and y axes; the x axis represents value ranges and the y axis shows the number of occurrences within the range. Some histogram styles are shown below.

Note Despite the similar appearance, a histogram is fundamentally different to a column chart.

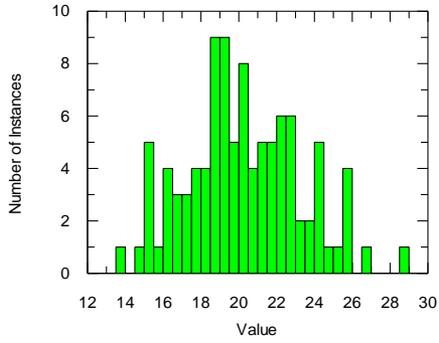

A distribution of data

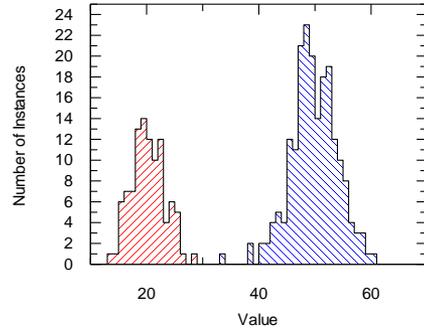


Comparison of two distributions

Graph Item

Overview of Graph Items

Graph items display a page of information that can hold one or more graphs, and may also display results or contain OLE objects. This item window is created from the main GraFit window using the **Add ▶ Graph** command. The appearance of a blank graph item window (at 50% magnification, in landscape orientation) is shown below.

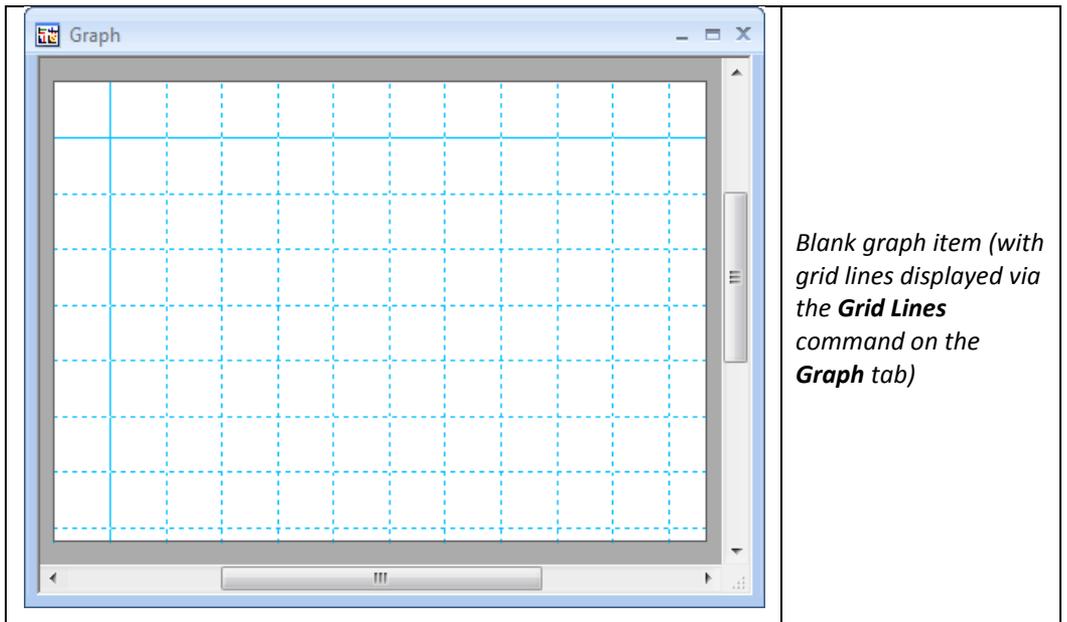

Controlling the Page View

The orientation of the printed page is controlled using the **File ▶ Print ▶ Page Setup** command⁵⁰. The magnification is set from the **View** tab⁵¹.

⁵⁰ The **File ▶ Page Setup** command in legacy mode

⁵¹ **View** menu, or using the Zoom selection on the toolbar, when in legacy mode

Grid Lines

By default, the graph is displayed with a blank background, but if required it is possible to add a background grid to help position the various items that form part of the graph. This grid does not show during print preview, and is not printed. It can be turned on and off by the **Grid Lines** command⁵², on the **Graph** tab, which toggles their display. The color of the grid lines, and whether GraFit graphs show gridlines by default, can be controlled from the **Preferences** command⁵³ on the **Graph** tab. The grid line spacing (inches or centimeters) depends on whether you are using the metric or the US measurement system; this is set from the **Options** command⁵⁴ (Number Format page) on the **View** tab. The top and left margins are shown on the graph as solid rather than dashed lines; the margin is specified using the **File ▶ Print ▶ Page Setup** command⁵⁰.

Element Outlines

By default, text elements and embedded objects are displayed with dotted and solid outlines, respectively. These do not show during print preview or printing, and can be hidden by the **Outlines** command⁵⁵, on the **Graph** tab, which toggles their display. Whether GraFit graphs show outlines by default can be controlled from the **Preferences** command⁵³ on the **Graph** tab.

Adding Drawing Objects

Various drawing objects can be added to the graph, as controlled by the **Draw** command group⁵⁶ on the **Graph** tab.


⁵² The **View ▶ Grid Lines** command in legacy mode

⁵³ The **Graph ▶ Preferences** command in legacy mode

⁵⁴ The **View ▶ Options** command in legacy mode

⁵⁵ The **View ▶ Show Outlines** command in legacy mode

⁵⁶ **Draw** menu or Drawing toolbar in legacy mode

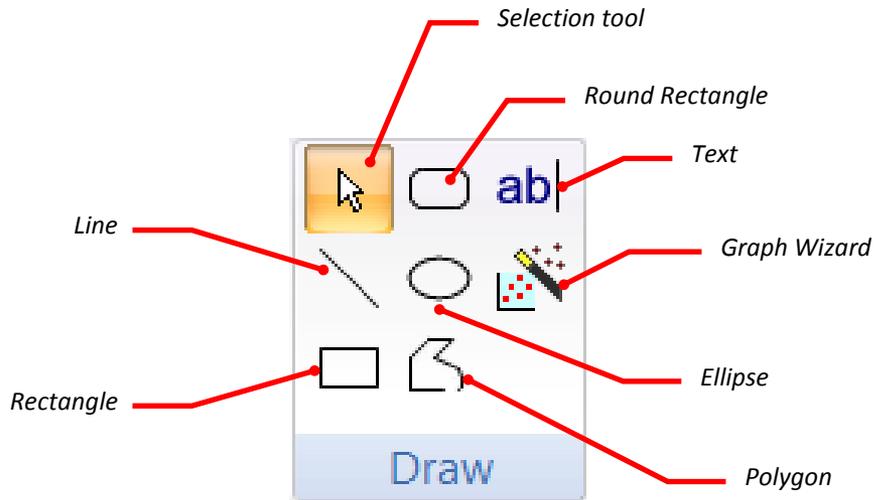

Selection Tool

When the selection tool is selected, the mouse pointer displays as an arrow. Clicking on a drawing object will select it, allowing moving or resizing. Clicking on the graph background and dragging will allow selection of a number of items, which can also be achieved by clicking on items while holding down the **SHIFT** key.

All Other Tools

Use of all other drawing tools follows a simple procedure. The tool is selected, which results in the mouse cursor being shown as a cross. Click the left mouse button where you wish the element to be drawn and drag until it is the correct size. When the button is released, the element is drawn using the current fill, line and font settings (these can be specified, see later). Creating graphs using the Graph Wizard will be discussed in the next chapter.


Magnification

The magnification of the graph can be adjusted to allow more or less of the page to be visible. Magnification is controlled by the **Zoom** command⁵⁷, found on the **View** tab. Note that **Zoom** only affects the on-screen appearance.

Using the Zoom Slider on the Status Bar


It is also possible to adjust the magnification of the data table text using the zoom slider found on the status bar.

⁵⁷ The **View** ► **Zoom** command in legacy mode

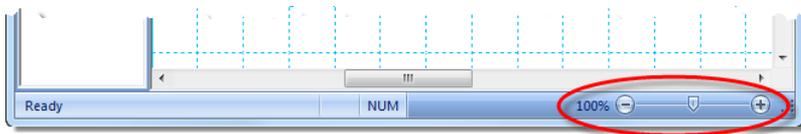

The zoom slider on the status bar


Select the data table then drag the slider to adjust the magnification. The magnification can be set more precisely by clicking on the zoom level text

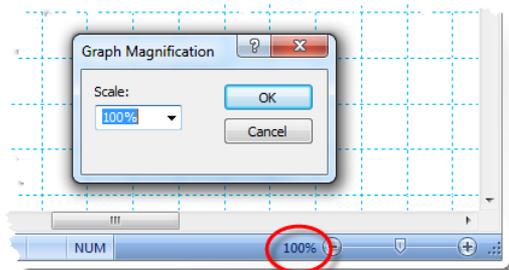

Set the Graph Magnification by clicking on the zoom level text

Manipulating Drawing Objects

Selecting Objects

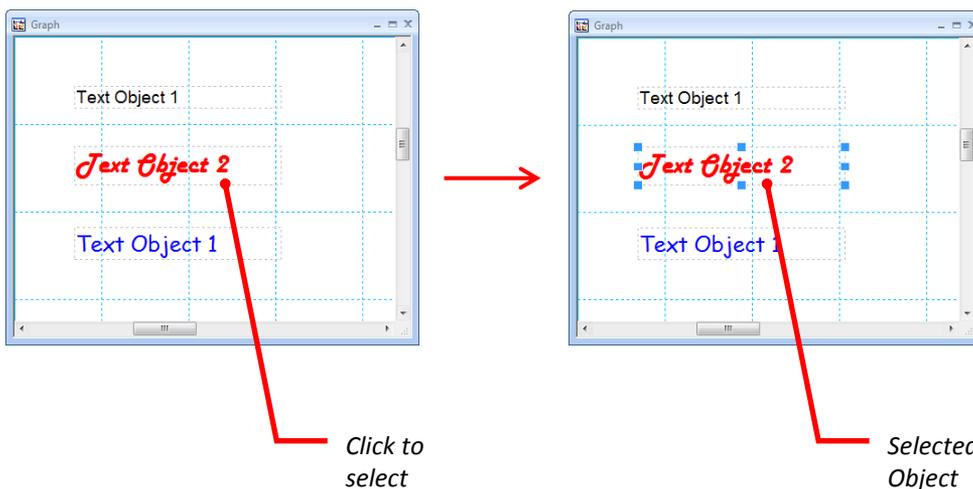

To select an object, move the mouse pointer over it and left-click.

Moving Objects

To move an object or group of objects, select the objects required and then either drag them to where you wish them to be placed, or use the cursor keys to move them around. Holding down the CTRL key while pressing the cursor keys produces smaller step sizes; hold down SHIFT gives larger ones.

Duplicating Objects


If the CTRL key is held down while dragging an object, a duplicate of the object will be moved to the new location.

Sizing Objects


To resize an object, first select it and then drag one of the resizing handles that are drawn in the corners of the object. Solid resizing handles are drawn on objects that can be resized; hollow handles are found on objects that do not permit resizing.



This selected object has solid resizing handles that permit resizing.

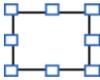

This selected object has hollow resizing handles, and cannot be resized.

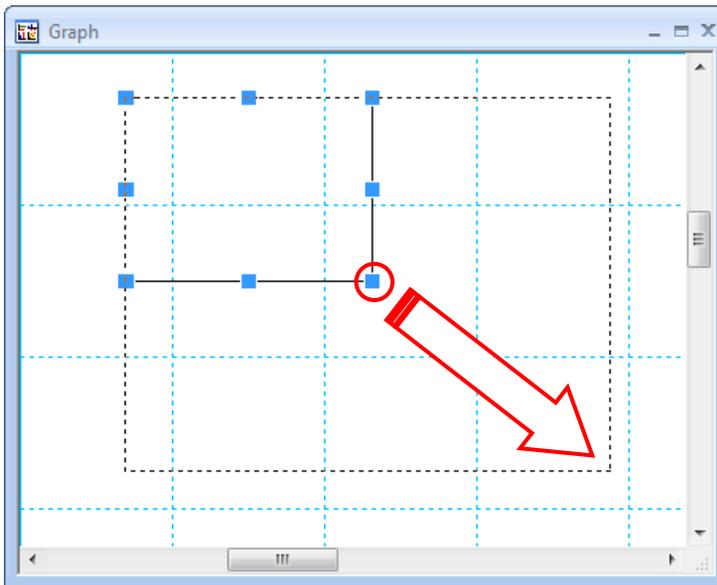

Drag one of the resizing handles to alter the size of a graph object.

The resizing action can be modified by holding down the SHIFT or CTRL keys while dragging the handle. It is possible to combine these effects by holding down both keys together.

Key	Effect on dragging
SHIFT	The aspect ratio of the object is maintained during resizing.
CTRL	The object remains maintains its current center position during resizing.

Resizing Text Objects

Text objects have slightly different resizing behavior to other graph objects. Text objects always wrap their text to fit, so the height of a text object is determined by its width, the text

it contains and the current formatting. Resizing a text object will therefore change the text wrapping.


Common Actions

The following actions are common to many of the objects that are found on the graph.

Action	Effect
Left click	The object under the cursor is selected, and any selected objects are deselected. If the cursor is not over a drawing object, then any selection is removed.
SHIFT + left click	The object under the cursor is added to the current selection.
Right click	An object-specific context menu is displayed.
Left click + drag	The selected object(s) are moved. Dragging can be stopped by right-clicking or by pressing the ESCAPE key.
CTRL + left click + drag	The selected object(s) are duplicated, and a copy of these of these is moved to the new location. This procedure can be stopped by right-clicking or by pressing the ESCAPE key.
SHIFT + left click + drag	The selected object(s) are moved, but the movement is constrained to a purely horizontal or vertical direction, depending upon the initial direction of the drag.
ALT + left click + drag	Certain objects are normally constrained in their locations. For example, graph scale text is normally only allowed to be moved normal to the axis orientation. Holding down the ALT key removes such constraints and allows arbitrary positioning.
Double left click	For most objects, this displays the object Properties dialog. The exceptions are text objects, for which this results in editing of the text, and embedded objects, for which this results in executing the default object action (normally opening the item for editing).

The shape of the mouse cursor gives information about the current action, or about which actions are available.

Cursor shape	Meaning
	Standard cursor. This cursor is shown when the selection tool is chosen, and the cursor is not over a selected object.

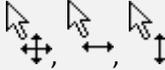Cursor shape	Meaning
	A drawing tool other than the selection cursor has been chosen. Left click where you wish the object to be positioned, then drag until it is the desired size.
	The cursor is over a resizing handle, which can be dragged to alter the shape of the element. There are several different resizing cursors, depending upon which edge of the object is to be resized.
	The cursor is over a selected object that can be moved in any direction, only in a horizontal direction, or only in a vertical direction, respectively.
	The cursor is over a selected object that can be moved in any direction, and the CTRL key is down. This means that dragging the element(s) will result in creation of a copy of the selection.
	(Polygon objects only.) The cursor is over one of the lines of the polygon, and the CTRL key is down. Under these circumstances, left clicking will add an additional control point to the polygon.
	(Polygon objects only.) The cursor is over one of the control points of the polygon, and the CTRL key is down. Under these circumstances, left clicking will delete this control point from the polygon.
	(Text objects only.) The cursor is over the text object, which is currently open for editing.

Grouping Objects

Drawing objects can be grouped in order to ensure that they can be moved as one. To group objects, select those required and choose the **Group** command⁵⁸, found in the **Objects** group of the **Graph** tab. Objects that have been grouped can have their grouping removed using the **Ungroup** command⁵⁹.

Aligning Objects

Selected objects can be aligned by selecting those required, and choosing the appropriate **Align** command⁶⁰, found on the **Graph** tab. Objects are aligned relative to all currently selected.

⁵⁸ The **Object ▶ Group** command in legacy mode

⁵⁹ The **Object ▶ Ungroup** command in legacy mode

⁶⁰ The **Object ▶ Align** command in legacy mode


Setting the Z-Order

Objects are drawn on the graph page in a defined order, with the uppermost item overlaying ones beneath. This z-order can be changed for a selected object by using commands found in the **Objects** group⁶¹ of the **Graph** tab.

Command	Meaning
 Bring to Front	The selected object is brought to the front of all objects on the graph.
 Send to Back	The selected object is placed underneath all objects on the graph. This may result in the object becoming obscured by others on the page.
 Move Forward	The selected object is moved one place forward in the z-order.
 Move Back	The selected object is moved one place back in the z-order.

Setting the Object Color

The fill, pen and text colors of an object are specified separately. Text color, as well as the text formatting, is set from the **Text** tab⁶¹. Pen and fill colors are set from the **Graph** tab⁶¹.

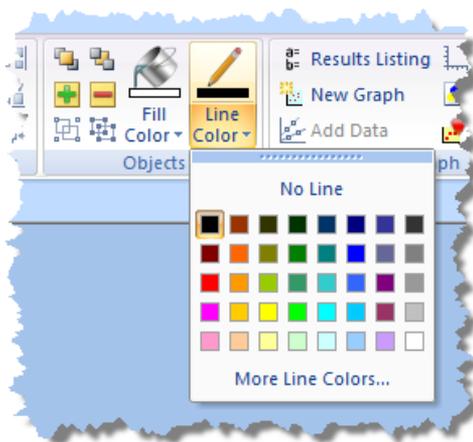

Setting the line color

Note Some objects (notably graph objects) have more complex color settings than can be made using the Color selectors. For these it is possible to select other color options from the object **Properties** dialog⁶², this command found on the **Graph** tab or by right-clicking on the object.

⁶¹ The **Object** menu in legacy mode


⁶² The **Edit ▶ Properties** command in legacy mode

Setting Default Colors

The default fill, line and text colors for new graph objects can be specified by altering the color when no object is selected. The new default colors will be used next time an object is added. Alternatively, the **Preferences** command⁶³ on the **Graph** tab can be used to set these defaults. This command also lets the default settings be stored for subsequent GraFit sessions.

Setting Text Properties

Text properties are specified by commands found on the **Text** tab⁶⁴.

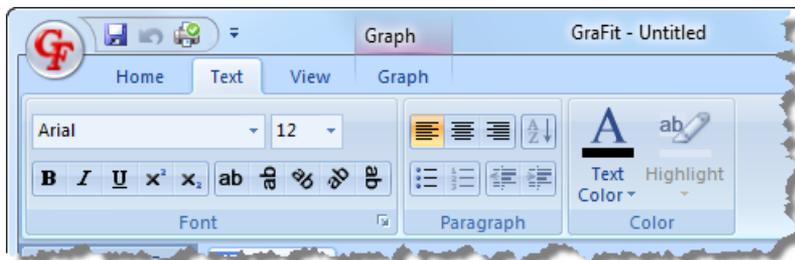

The font, size, style, alignment and orientation of the text can be specified from this toolbar. Further text properties that are specific to particular objects may also be available from the **Properties** dialog⁶² for the selected object.

If a text object is selected by clicking on it, altering text properties such as font face, size and color, will alter the appearance of the entire object. Finer control of text settings is possible for editable text by double-clicking on the text item, then selecting just a portion of the text. Changes to text formatting will then apply only to the selected text.

Setting Default Text Properties

The default font properties for new text objects can be specified by altering a property (for example, changing the font size) when no object is selected. The new defaults will be used next time a text object is added.

Using Field Codes

Text within graphs can contain field codes that are substituted with values such as entries from a data table. This allows graph text to incorporate information from other item windows; such information is updated dynamically. Field codes are entered between curly braces (“{}”) and have the following syntax:

```
{“item name” command parameters}
```

The first entry is that of the item window from which the information is to be obtained. This name must be surrounded by quotes if the item name contains space characters. The next entry is a command name, which is specific to the item window type. Finally, there are one or more optional parameters, which are specific to the command and separated by spaces.

⁶³ The **Graph** ► **Preferences** command in legacy mode

⁶⁴ The **Text** menu, or on the Formatting Style toolbar, in legacy mode

Data Table Command	Parameter(s)	Description
AssValue	col index	Returns the associated value "index" for column "col" (the lowest allowed index number is 2, which is the x_2 value; index number 3 gives the x_3 value, etc.)
ColumnName	col	The parameter col is the number of the column. The value returned is the name of the column within the data table
Data	col row	The data entry from column "col" and row "row" is returned. The number formatting of the data value is as used by the source column.
IfTrue	col row text	If the data value at column "col" and row "row" is not zero, then returns the text "text".
IfFalse	col row text	If the data value at column "col" and row "row" is zero, then returns the text "text".


Result Item Command	Parameter(s)	Description
xColumn	no	Returns the name of the x data column. For data with a single dependent variable, "no" is the fit number; for data with two independent variables, "no" is the data set number; for data with multiple independent variables then "no" is the x data number.
yColumn	no	Returns the name of the y data column. For data with a single dependent variable, "no" is the fit number; for data with two independent variables, "no" is the data set number.
errColumn	no	Returns the name of the error data column. For data with a single dependent variable, "no" is the fit number; for data with two independent variables, "no" is the data set number.

Result Item Command	Parameter(s)	Description
<code>parameter</code>	<code>fitno</code> <code>paramno</code>	Returns the value of a fitted parameter from the data fitting. The number formatting is taken from the result item format specification. “ <code>fitno</code> ” is the fitting number (1 if there is a single fit); “ <code>paramno</code> ” is the number of the parameter (1 for the first parameter, 2 for the second, etc.)
<code>paramname</code>	<code>fitno</code> <code>paramno</code>	Returns the parameter name from the data fitting. See the <code>parameter</code> description for more information.
<code>error</code>	<code>fitno</code> <code>paramno</code>	Returns the error value from the data fitting. See the <code>parameter</code> description for more information.

Example


The following text would allow a data value present in the first row of the second column of table “My Data” to be included within the graph. The graph text is dynamically updated whenever the underlying data value alters.

The answer to life the universe and everything is {"My Data" Data 2 1}

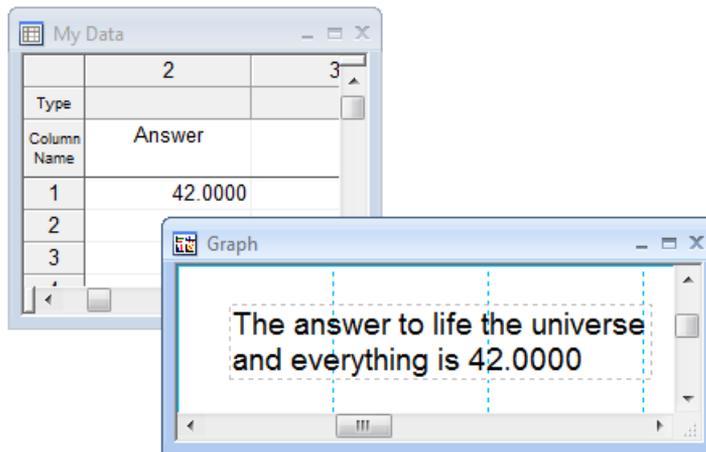**Entering Field Codes into Graph Text the Easy Way**

An easier way of entering field codes is to use the **Field Code** command⁶⁵ on the **Graph** tab. First create a text item on the graph, then use the **Field Code** command where you wish the field code to be placed. The following Text Field dialog box controls the entries.

⁶⁵ The **Edit ▶ Field Code** command in legacy mode

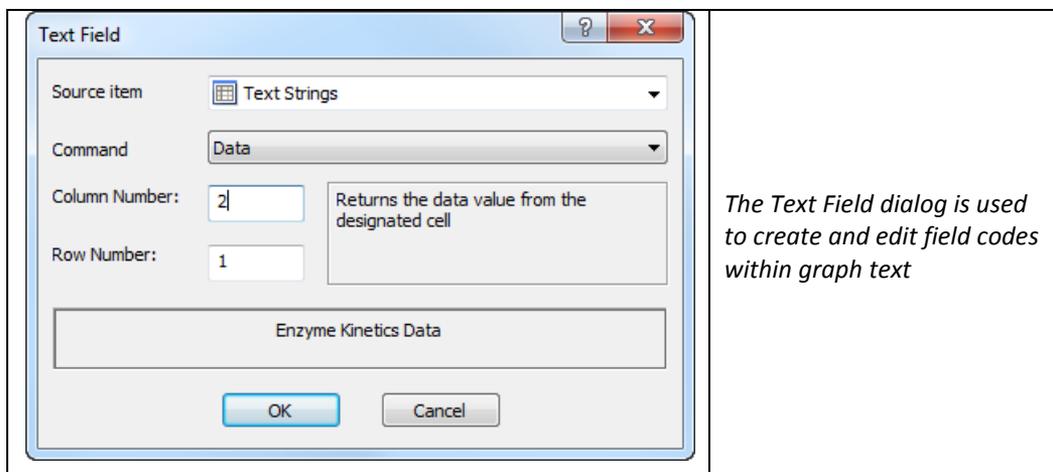


Chose the Source item, set the Command as required, then enter suitable values for the parameters as prompted. To edit an existing text field code, position the cursor between the two curly brackets and enter the **Field Code** command—the Text Field dialog box will be populated with the currently selected field code settings.

General Object Properties


Selecting the **Properties** command⁶² on the **Graph** tab, or double-clicking on a graph object, shows the object properties dialog box. This dialog box contains a number of pages, some which are common to many objects and some which are specific to the object selected. The common property pages are shown below; specific property pages discussed later.

Fill Properties

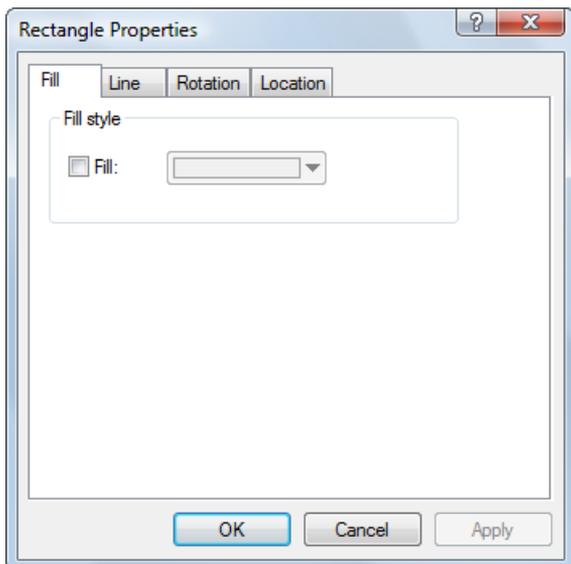


Check the **Fill** box to fill the object with the selected color. If this box is unchecked, the object will be transparent.

Line Properties

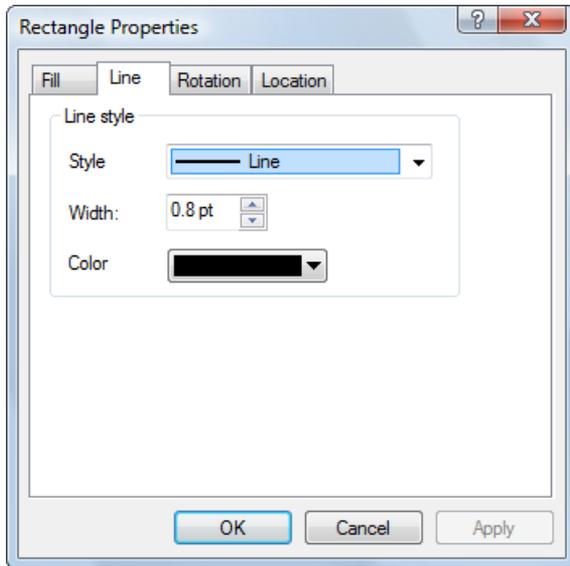


The Line page sets the width, style and color of the line. It is possible to set the Style to "No line" if you do not wish a line to be drawn.

Rotation Properties

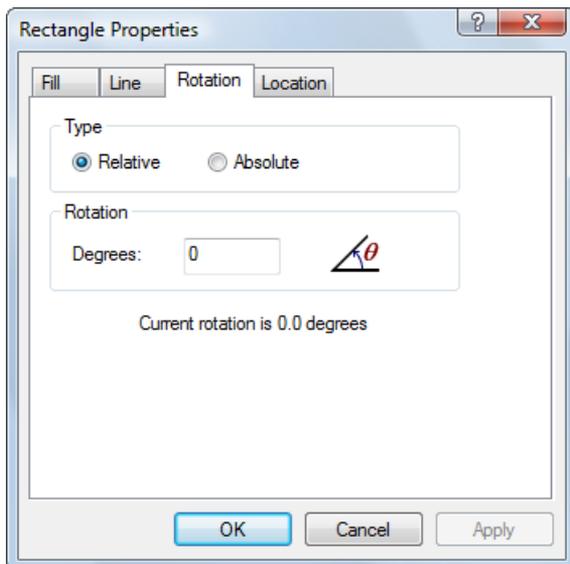


Objects can be rotated, and their rotation is specified from the Rotation page. Rotation is entered in degrees, which is measured in a counter-clockwise direction. The rotation can be set as an absolute figure, or relative to the current rotation.

Location Properties

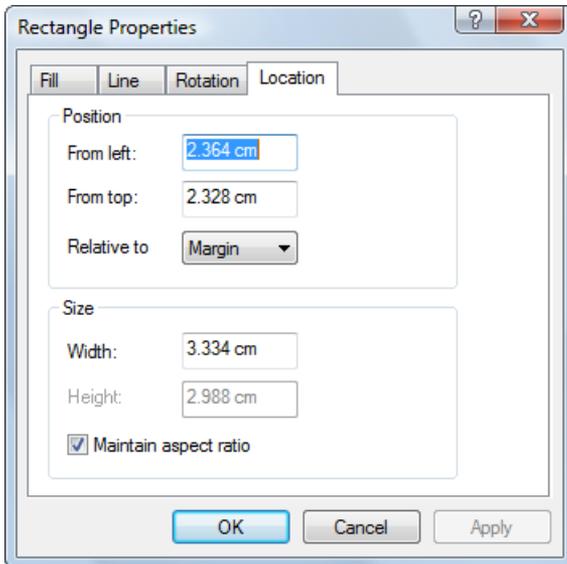


The location of the selected object is set from the Location page. This allows the object position to be specified exactly, relative to the current margin or to the edge of the page.

Arrows Properties

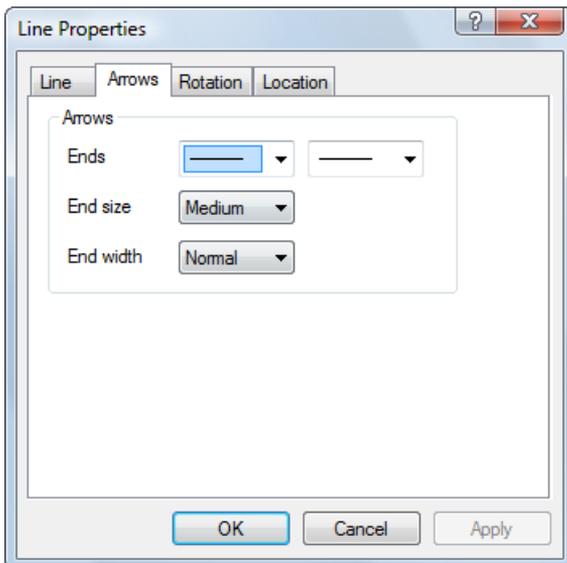


Objects that can have lines with terminating ends allow the end type to be specified via the Arrows page. Choose from a variety of arrow styles for each end of the line, together with the end size and width. The arrowhead color is set to be the same as the line color.

Font Properties

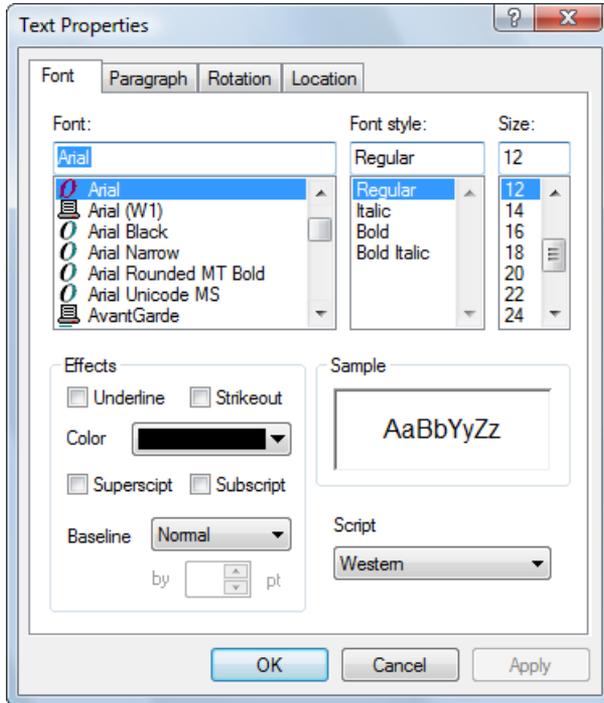

The Font page sets the font to use, together with the style, size and various effects including superscript and subscript options.

Paragraph Properties

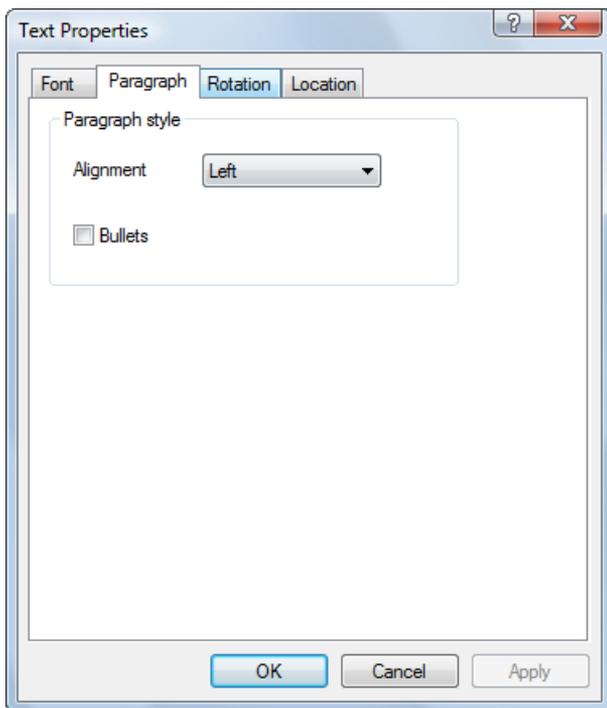

This page sets the formatting for the selected paragraph(s) or text.

Exporting Graphs

While it is possible to insert graphs into other applications via the Clipboard, for publishing purposes it is often important to create image files in a portable format. The **Export Graph** command⁶⁶ on the **Graph** tab allows graph images to be saved using a file format that can be used by other applications. The following graph file types can be generated:


- Windows Metafile (WMF)
- Enhanced Metafile (EMF)
- Bitmap file (BMP)
- Tagged Image Format File (TIFF) (compressed or uncompressed)

For most purposes, the compressed TIFF format is the most suitable.

Windows Vista or Windows 7


When running GraFit on Windows Vista or Windows 7 (or later), the program uses the new - style Save dialog box, as shown below.

⁶⁶ The **File** ► **Export Graph** command in legacy mode

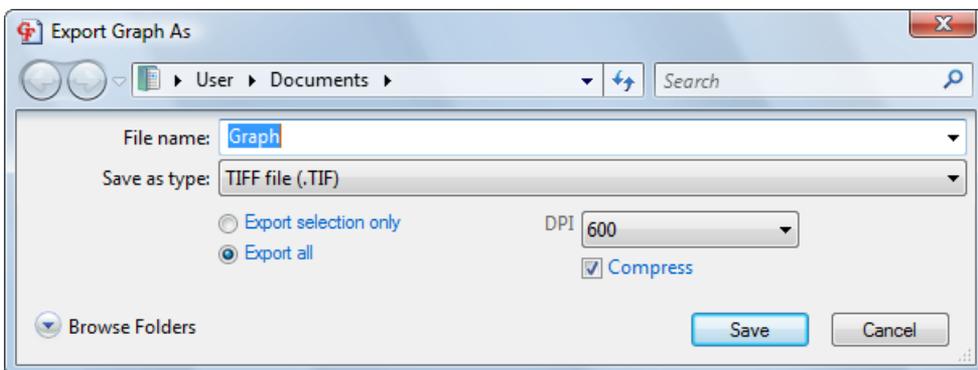

Windows 2000 or Windows XP

Earlier versions of Windows use the older-style File Save As dialog box.

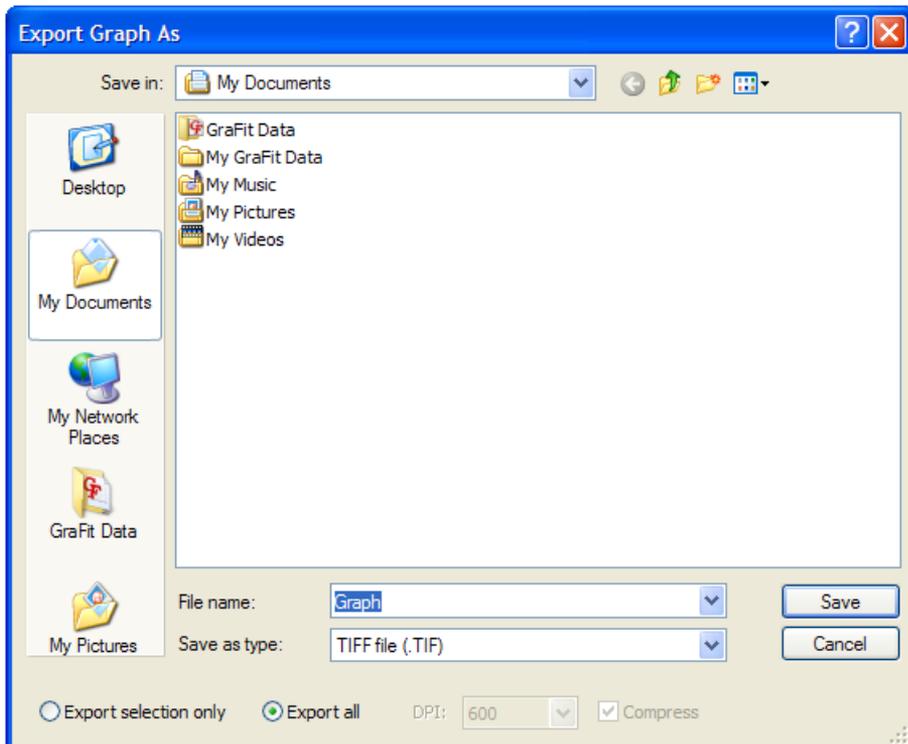


X/Y Scatter Graphs

The XY Graph Wizard

Chapter 12 described the basics of creating graphs. Selection of an X/Y Scatter Graph produces the “XY Graph Wizard” to assist with the procedure (all settings can be adjusted later, if required).

Select Data Source Page

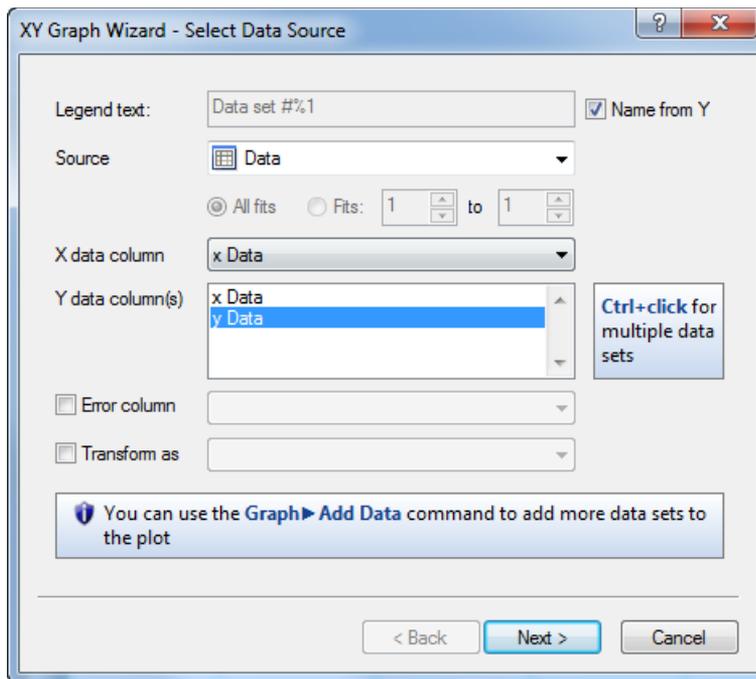


The Select Data Source page allows you to choose the source and type of the data to plot.

Legend text

The data set plotted is associated with an entry in the data legend. The entry shows the symbol used to represent the data, and a descriptive text. By default, the text used is taken from the y data column. However, if **Name from Y** is unchecked then you can enter your own

text into this field. This text field can also contain two special entries that allow the data set number or the fit number to be inserted into the data legend:

Entry	Source Item Type		
	Linear/non-linear Fit	Equation	Data table
%1	Replaced by the fit number	Replaced by the data number	Replaced by the data number
%2	Replaced by the data number	(not applicable)	(not applicable)

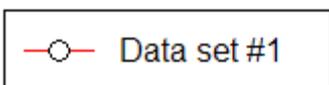

A data legend showing the symbol and descriptive text.

Source

Graphs are plotted using data that are present in one of the items found in the GraFit workspace. The **Source** selection lists all items that contain suitable data for plotting, and one of these must be selected. GraFit graphs can plot three different types of plot on an X/Y graph, depending on the source selected.

Item Type	Effect
 Data	When raw data are plotted, the data table item is selected and the x and y columns are selected for display. The data points can optionally be joined with lines or curves. If the data have been fitted to an equation, it is usually better to select the data fitting item as this will also show the fitted curve.
 Equation	To plot a theoretical curve, select the equation item for the curve you want.
 Linear Fit	Selecting a data fitting item (linear or non-linear fit) results in the data being plotted together with the fitted curve. This is usually the way of plotting data fitting results. Data fitting items may contain multiple fits, so when a data fitting item is selected, you have the option of displaying a single fit (and selecting which you require) or of plotting all fits.
 Non-linear Fit	

X data column / Y data column

When the **Source** is a data table item, it is necessary to select which columns hold the x and y data values to be plotted. For all other source types, these selections are not available.

Error column

Checking the **Error column** box allows a quick way of producing an error plot when the **Source** is a data table item. Select the column that holds the error values from the list. More extensive error plotting options are available once the graph has been produced (see page 181).

Transform as


Check the **Transform as** box to plot the data after transformation using the selected transformation item. All available transformations (which must be present in the workspace) are listed. Transformation of data is only available if the **Source** is a data table or a non-linear fit item. See Chapter 10 for more information about plotting transformed data.

Set Initial Scaling

The screenshot shows the 'XY Graph Wizard - Set Initial Scaling' dialog box. It is divided into two main sections: 'X axis' and 'Y axis'. Each section has a 'Type' dropdown menu set to 'Linear', a radio button for 'Autoscale' (which is selected), and a radio button for 'Manual scaling'. Under 'Manual scaling', there are 'From' and 'To' input fields. For the X axis, 'From' is 0 and 'To' is 11. For the Y axis, 'From' is 1 and 'To' is 7. Below these sections are 'X axis text' and 'Y axis text' fields, both containing 'x Data' and 'y Data' respectively. At the bottom, there are three buttons: '< Back', 'Next >', and 'Cancel'.

The Set Initial Scaling page allows the scale of the x and y axes to be specified. By default the scale is set automatically, and linear axes are chosen. These can be altered as required (it is also possible to adjust the scale settings later). The text used for the x and y axis can also be entered.


Select Tick Styles

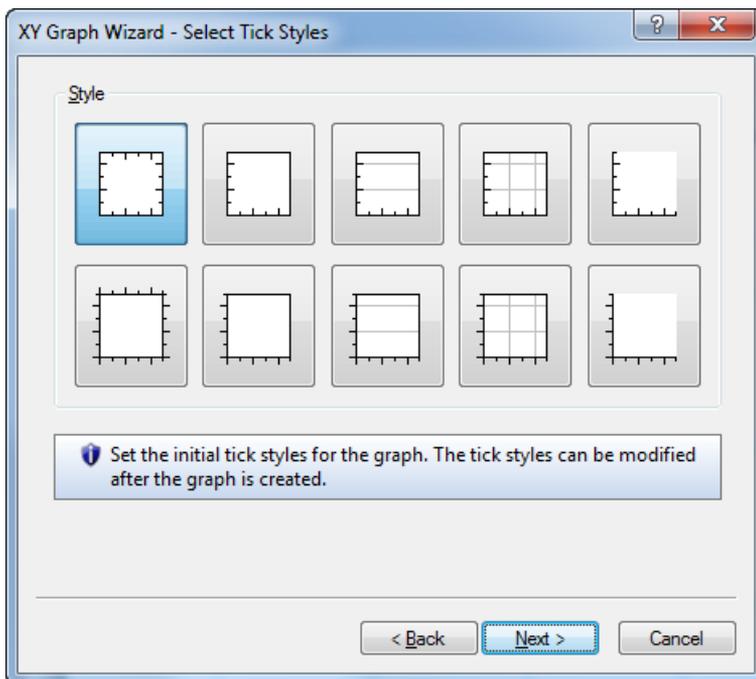

The Select Tick Styles page allows rapid selection of some of the more common axis tick styles. More extensive tick style choices can be made once the graph has been created.

Select Data Point Style

(not equation items)

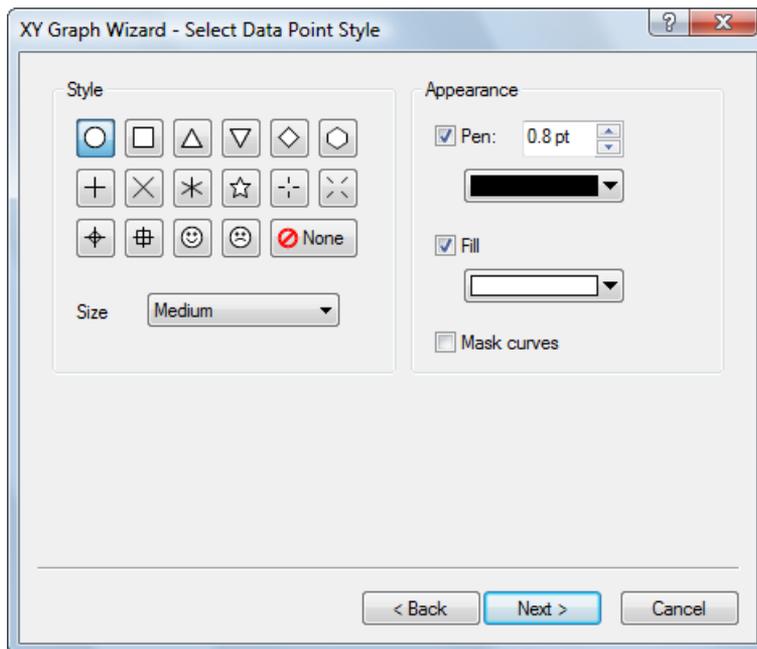

The Select Data Point Style page allows the data point representation to be chosen.

Size

Data points are scaled relative to the overall graph width, and so the absolute size will depend upon the graph dimensions. Choose a size that represents your data appropriately; with more data you should use smaller data points. If very many points are present it can be more effective to set the **Style** as None, and to join the data points with lines (next page).

Pen

Check this box to draw an outline around the data point. The size and color of the outline can be specified.

Fill

Check this box to fill the point with the solid color selected. If **Fill** is not checked, the data points will be transparent.

Mask curves

Check this option to ensure that any curves on the graph are not drawn immediately adjacent to the data points.

Enter Curve Parameters

(Equation items only)

XY Graph Wizard - Enter Curve Parameters

Parameters

Read from table

Column

Direct entry

Parameter	Value
Vmax	
Km	

< Back Next > Cancel

When plotting theoretical curves, this page allows the parameter values that describe the curve to be entered.

Read from table

When this option is selected, the parameter values are read from the data table item and column that are selected. Any column in the data table can be selected, but it is most convenient to use a column that has been specifically created to facilitate parameter entry (see page 77).

Direct Entry

When this option is selected, the parameter values are entered directly into the grid.


Enter Curve Constants

(Equation items only)

If the equation definition uses constants, their values are entered into this page. The operation of this page is analogous to that of the **Enter Curve Parameters** page described above.

Select Fitting Curve Style

(Fitting and equation items only)

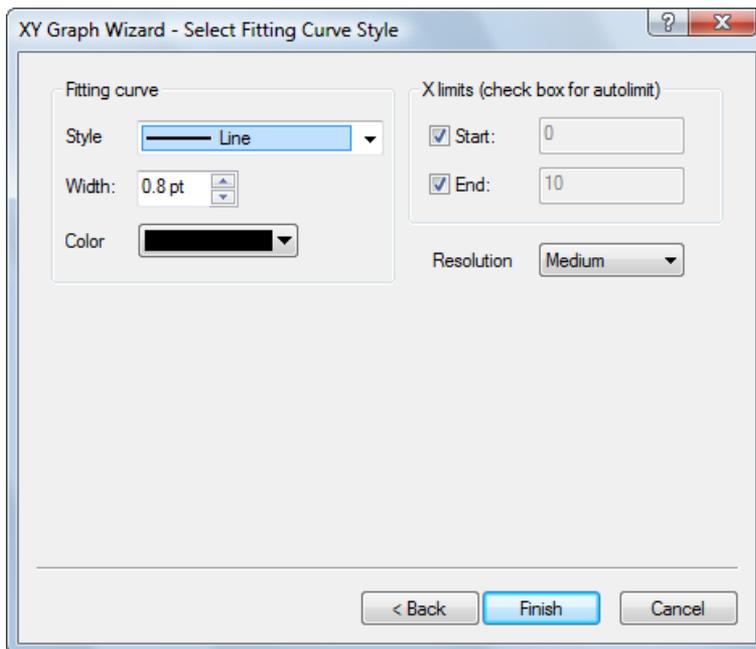

Fitting curve

Choose the width, line style and color of the curve that represents the fit.

X limits


By default, the fitted curve is drawn over the entire range of the x axis. Some equations, however, are undefined over certain ranges (for example, at $x < 0$). When the **Start** or **End** box is checked, GraFit will set the limit automatically, but by unchecking these boxes it is possible to define the range over which the curve should be displayed.

Resolution

Fitted curves are calculated for a range of values between the **Start** and **End** limits. The **Resolution** entry controls how many points are used in this range. High resolution uses more points, which takes longer and uses more memory. In general, use the medium resolution unless the curve has a sharp discontinuity (e.g. when plotting a step function). Under such circumstances use the high resolution setting.

Select Join Style




(not equation items)

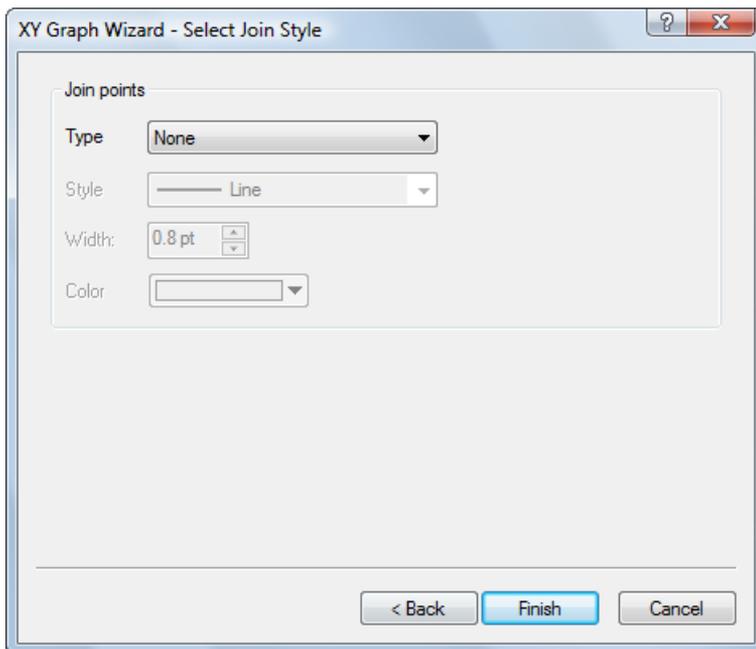

Type

Chooses the type of join to use.

Effect	Join Type
	<p>None</p>

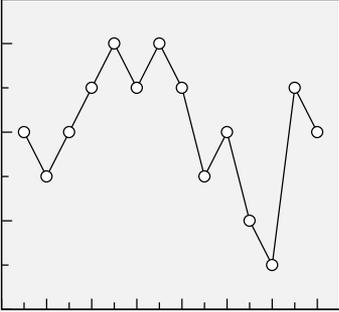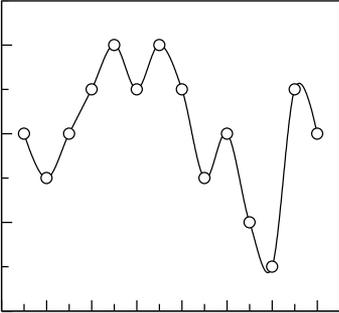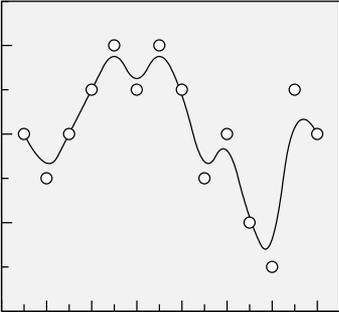Effect	Join Type
	<p>Line</p>
	<p>Bezier</p> <p>Three degrees of curvature are available: slight, medium and strong. The illustration shows a medium curve. Bezier curves give smooth lines between the points, and always pass each data point.</p>
	<p>Spline</p> <p>Four orders of spline curve are available: 3, 4, 5 and 6; the illustration is created using order 3. Higher orders give smoother curves. Spline curves smooth out variation in the data, and do not necessarily pass through the data points.</p>

Width, Style and Color

These options control the appearance of the line that is drawn joining the points.


Modifying Graphs

Once created, graphs can be customized as required. This involves selecting the appropriate part of the graph, and invoking a specific property sheet by double-clicking, right-clicking or choosing the **Properties** command⁶⁷ on the **Graph** tab.

⁶⁷ The **Edit ▶ Properties** command in legacy mode

Selectable Regions of a Graph

The illustration below shows the various parts of a graph that can be selected. For each, a specific property sheet is available that controls the formatting options.

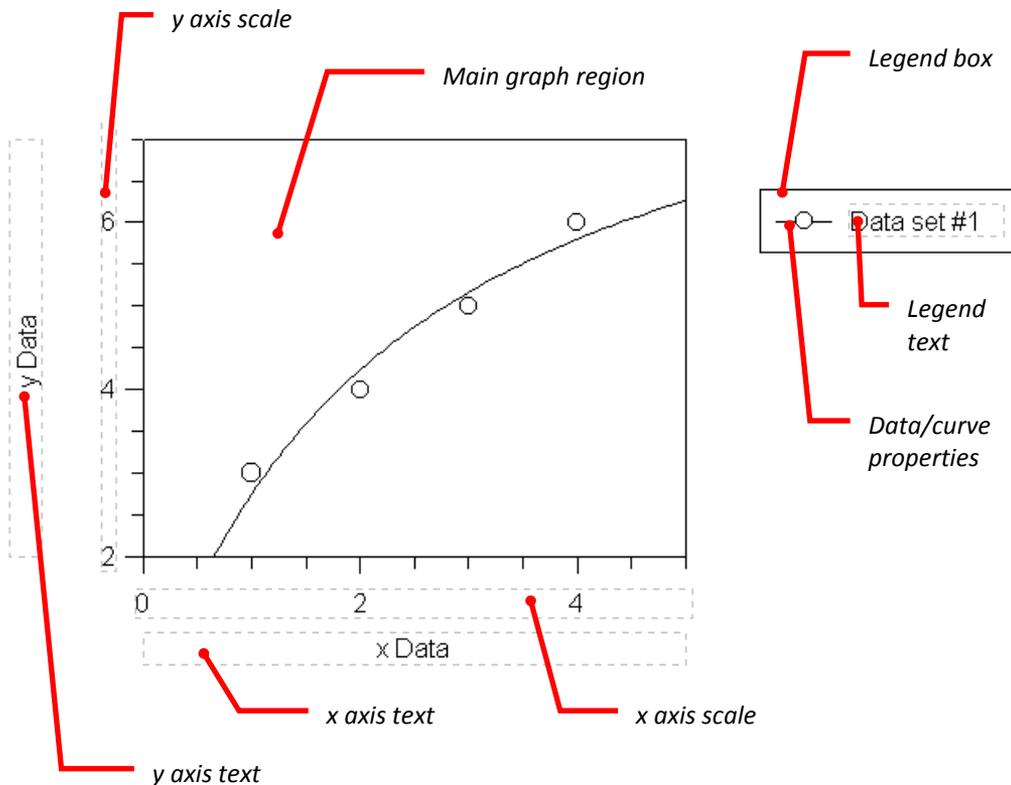


Main Graph Region

The main graph region acts as the “parent” of the other selectable regions and is always visible (the other parts of the graph can optionally be hidden). Moving the main graph region will also move all the other “child” parts of the graph; however, the child regions can be moved individually to position them. This ensures that the whole graph can be moved around the page while maintaining the relative positioning of the constituent parts.

Note *The x and y axis scales have a particular relationship with the main graph, as it is usual to ensure that their numbers line up with the scale marks. For this reason, by default it is only possible to move the x axis scale in a vertical direction, and the y axis scale in a horizontal direction. This effect can be overridden by holding down the ALT key when moving the axis scale texts.*

General Page

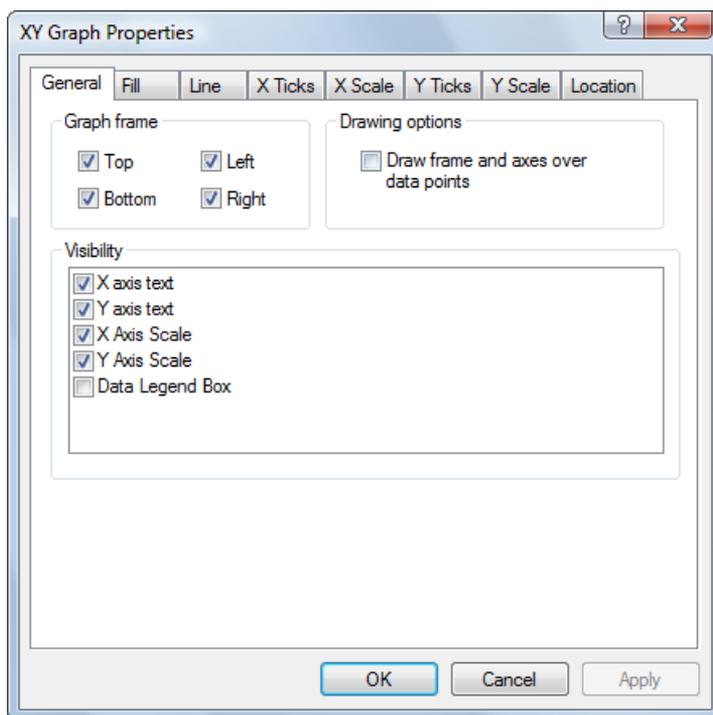

Graph frame

Check the various boxes to draw a line on the various edges of the graph.

Drawing options

Check this box to specify that the graph frame, axes and tick marks should be drawn after the data points, which will result in these items lying on top of the data points. By default, XY graphs have this option unchecked and column graphs have it checked.

Visibility

The components of the graph are listed here, and their visibility can be set by checking the appropriate box. Individual items on the graph can also be hidden selecting the **Hide Object** command⁶⁸ on the **Graph** tab.

Fill Page


The fill style applies to the graph background.

Line Page

The line style applies to the frame drawn around the graph, as set on the **General** page. It does not apply to the color of the tick marks, which are set on the appropriate tick page.

⁶⁸ The **Edit** ► **Hide Object** command in legacy mode

X Ticks Page

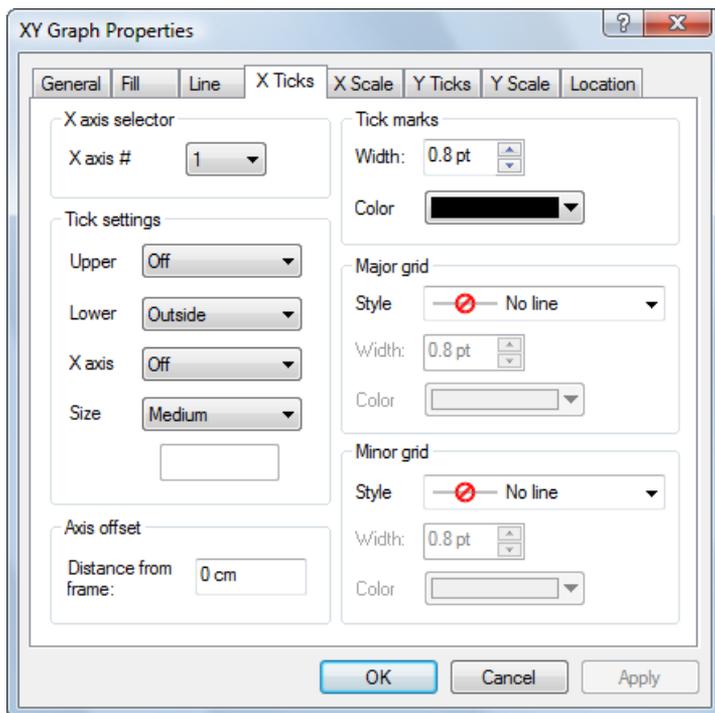


X axis selector


X/Y graphs can have multiple x axes, and the tick settings for each are set individually. Set the x axis number from this box (when first created, only a single x axis will be available).

Tick settings


Ticks can be drawn on the inside, outside, both sides or neither side of the **Upper** graph edge, the **Lower** graph edge, or on the **X axis** itself (assuming the scale settings are such that this lies in the center of the graph). The tick sizes can be set to small, medium or large, which are scaled relative to the overall graph width, or to “specific,” where the exact tick size is entered into the edit field.

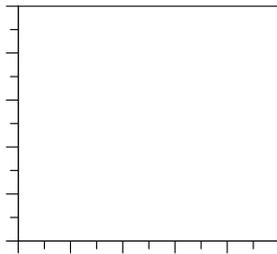

Left/Lower ticks outside

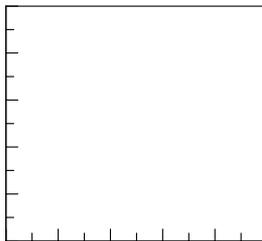


Left/Lower ticks inside

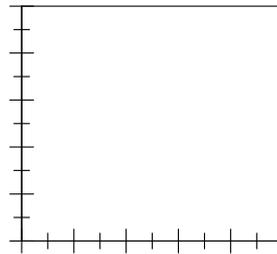


Left/Lower ticks both sides

Axis offset


The axis offset controls the distance that the axis scale is drawn from the frame. If the offset is 0, the scale is drawn immediately adjacent to the central graph frame; increasing offset results in the scale being removed from the frame, as shown below:

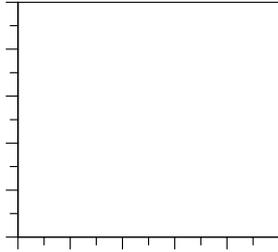

No offset

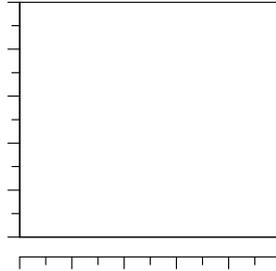

x axis offset

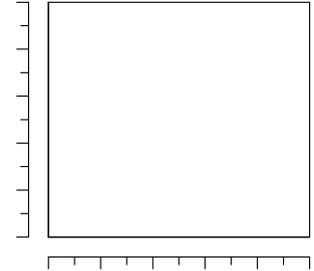

x and y axes both offset


If more than one x axis is present, it can be useful to set each to have a different axis offset. This allows the two scales to be differentiated.

Tick marks


These settings control the width and color of the pen used to draw the tick marks.

Major grid / Minor grid


These settings control the display of grid lines drawn across the axis at intervals determined by the **Major increment** and **Minor increment** settings (Scale Page).

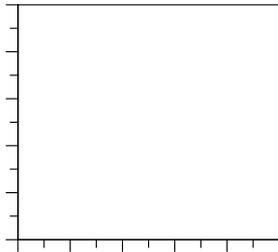

No grid

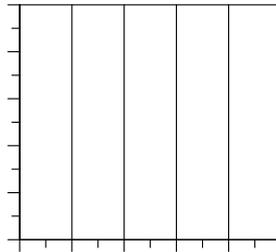


x axis with major grid (solid lines)

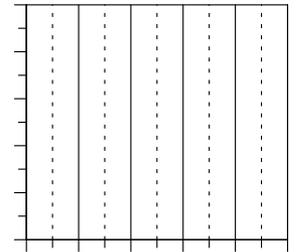


x axis with major grid (solid lines) and minor grid (dotted)

X Scale Page

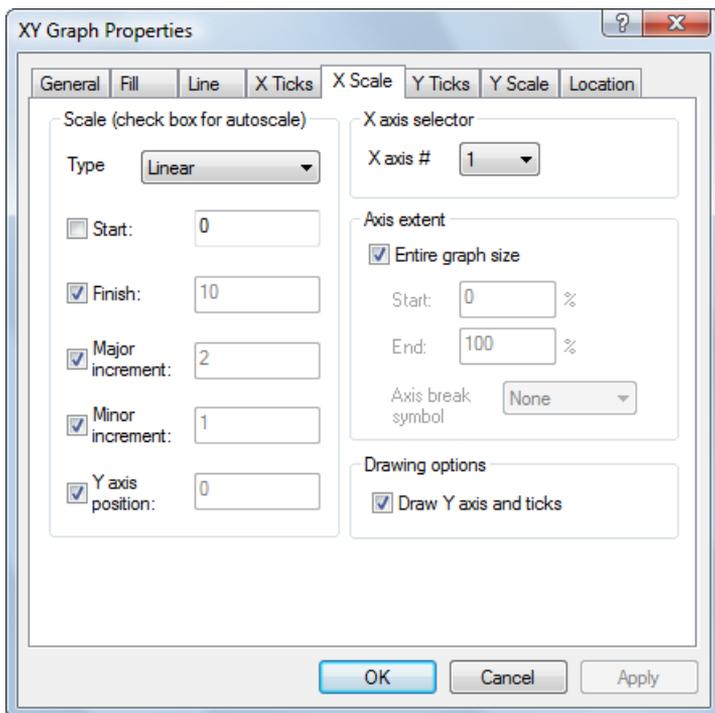


Scale Type

Select linear or logarithmic scaling.

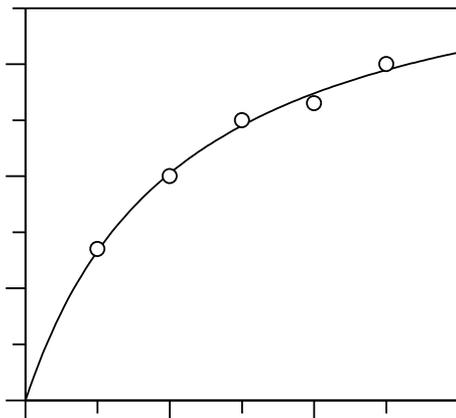

Graph with linear x axis scale

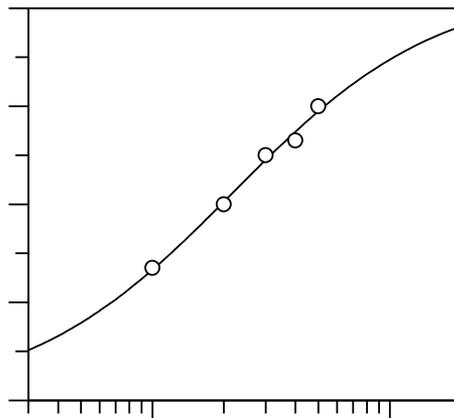

The same graph with a logarithmic x axis scale

Scale Start

If the box is checked, this value is set automatically from the graph data. Uncheck the box to specify the starting value of the x axis (left edge).

Scale Finish

If the box is checked, this value is set automatically from the graph data. Uncheck the box to specify the finishing value of the x axis (right edge).

Scale major increment


If the box is checked, this value is set automatically from the graph data. Uncheck the box to specify the major increment value. The major increment value controls how often scale numbers and major ticks are drawn.

Scale minor increment

If the box is checked, this value is set automatically from the graph data. Uncheck the box to specify the minor increment value. The minor increment value controls how often major ticks are drawn.

Scale y axis position

If the box is checked, this value is set automatically from the graph data. Uncheck the box to specify the axis position, which should lie in-between the **Start** and **Finish** values.

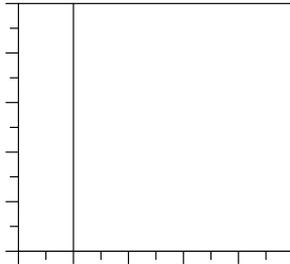

*The y axis position has been set so that it appears between the **Start** and **Finish** values. If required, tick marks can be drawn on the y axis line from the Y Ticks page (Y axis Tick settings).*

X axis selector

X/Y graphs can have multiple x axes, and the scale settings for each are set individually. Set the x axis number from this box (when first created, only a single x axis will be available).

Axis extent

By default, the **Entire graph size** box is checked, which ensures that the axis is drawn over the whole graph width. Unchecking this box allows the axis to be drawn over just a portion of the graph, which can be useful when multiple axes are present. Enter the **Start** and **End** values as percentages of the whole graph width. The **Axis break symbol** setting controls whether a break mark is drawn at the end of the axis. This setting is used when drawing split axis scales (see Chapter 15 for more information about multiple axes).

Y Ticks Page

The Y Ticks page controls the tick settings on the y axis in the same way as the X Ticks page controls the tick settings on the x axis (see page 168).

Y Scale Page

The Y Scale page controls the scale settings for the y axis in the same ways as the X Scale page controls scale settings on the x axis (see page 170).


Location Page

This sets the location of the graph within the graph page.

Axis Scale Regions

The axis scale region controls the display and formatting of the scale numbers.

Scale Page

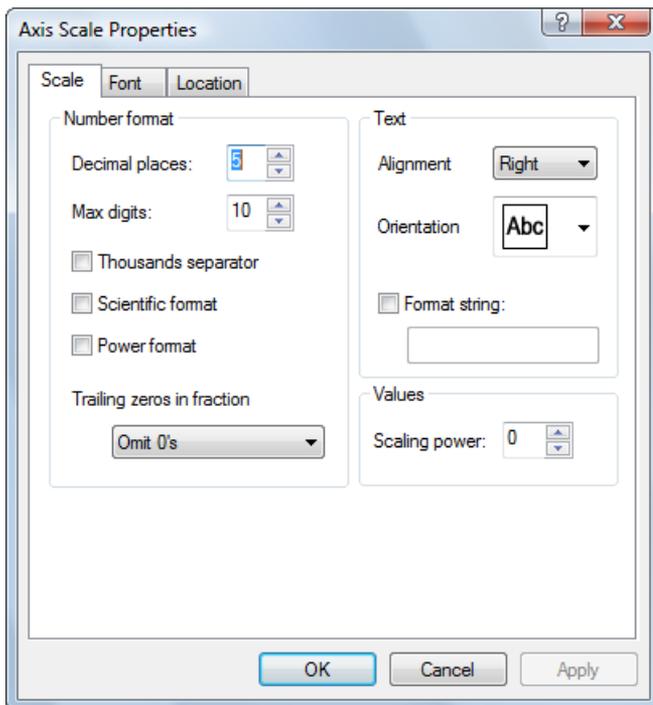




Number format

These settings determine the way that the numbers are displayed.

Setting	Description
Decimal places	Enter the number of decimal places to display. Note that trailing zeros will be removed if the Trailing zeros in fraction setting is "Omit 0's."
Max digits	Enter the maximum number of digits to be used. The displayed values will be rounded to ensure that only the requested precision is displayed.

Setting	Description
Thousands separator	Check this option to separate groups of digits by a separator. Typically, digits are grouped into threes and separated by a comma. However, the grouping and the separator used depend on values set by in the Options command ⁶⁹ , Number Format pane, found on the View tab.
Scientific format	Check this option to display numbers using scientific format. For example, 1.234e+003 rather than 1234.
Power format	Check this option to display numbers using a power representation. For example, 1.234 x 10 ³ rather than 1234.
Trailing zeros in fraction	This controls the way that trailing zeros are displayed in a fraction.

Some examples of different formatting are shown below; all have **Decimal places** set to 4, and **Max digits** set to 10.

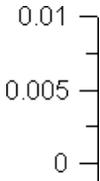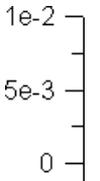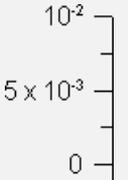Effect	Settings	Effect	Settings
	Omit 0's		Omit 0's, scientific format
	Show 0's		Omit 0's, power format

Text

These settings control the text display.

Setting	Description
Alignment	The scale text can be aligned to the left, right or centrally.
Orientation	Select the orientation of the text. The following orientations are available: 

⁶⁹ The **View ► Options** command in legacy mode

Setting	Description
Format string	Check this option to embed the numeric value within a text string. The number is substituted for any occurrences of “%1” within the entered text.

Values

The **Scaling power** setting allows a scaling factor to be applied to the numbers. Before display, the scale numbers are multiplied by 10 raised to the power of the value entered. The power must be a whole number, and may be positive or negative. For example, an entry of -4 would cause the value of 40,000 to be displayed as 4. An entry of 3 would cause the value 0.002 to be displayed as 2 ($0.002 \times 10^3 = 2$). A scaling power of 0 results in the actual value being displayed ($10^0 = 1$).

Scaling power entries are very useful when dealing with large or small data values. By scaling the displayed numbers it is possible to produce graphs that are correctly scaled, but to enter conventionally ranged values into the data table. For example, time values in nanoseconds could be entered into the data table as whole numbers of nanoseconds then a **Scaling power** of -9 applied.


Axis Text Regions

Axis text properties are controlled in exactly the same way as normal text objects. See Chapter 13 for further information.

Legend Box Region

The legend box lists all the data sets currently present within the graph.

Legend Page

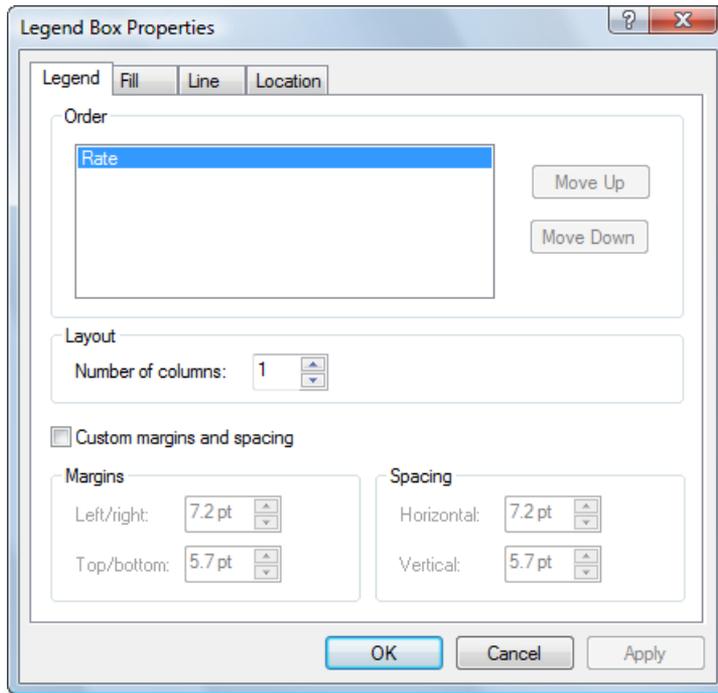


Order


The current arrangement of items within the legend box is displayed. This order can be changed by selecting an entry and clicking the **Move Up** or **Move Down** buttons.

Layout

This section controls the number of columns into which the data entries are arranged.

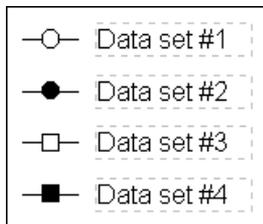

*4 data entries, with **Number of columns** set to 1 (the default setting).*

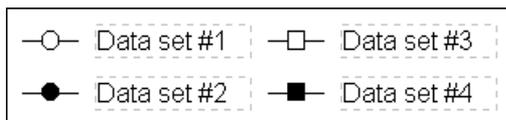

*4 data entries, with **Number of columns** set to 2.*

Custom Margins and Spacing

Checking the **Custom margins and spacing** box allows the spacing of the elements within legend box to be set explicitly.

Fill Page

This page sets the fill settings for the legend box.

Line Page


This page specifies the line used to draw the surround of the legend box.

Location Page

This page controls the position of the legend box within the graph page.


Data/Curve Properties

Data Points Page

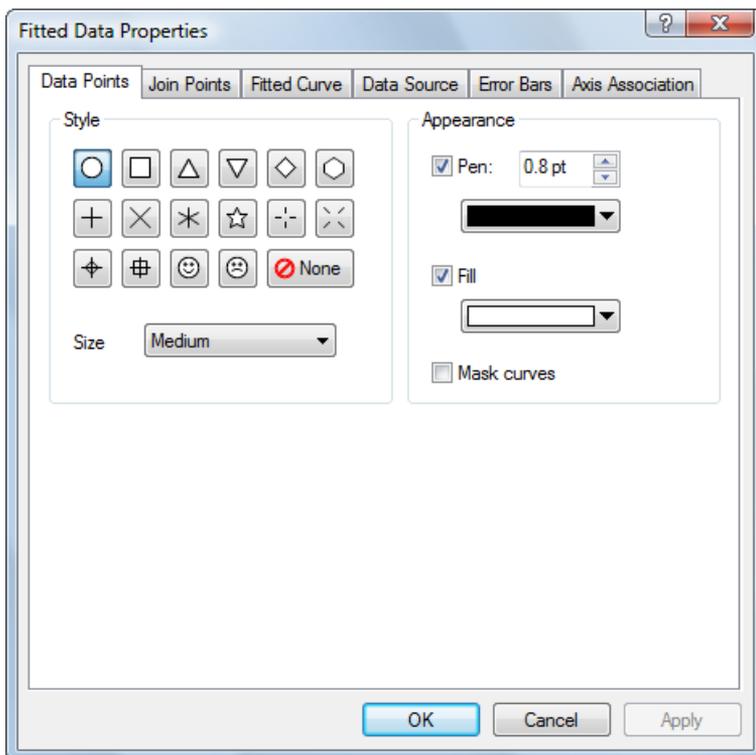

This page is equivalent to the Data Point Style page of the XY Graph Wizard. See page 161 for further information.


Join Points Page

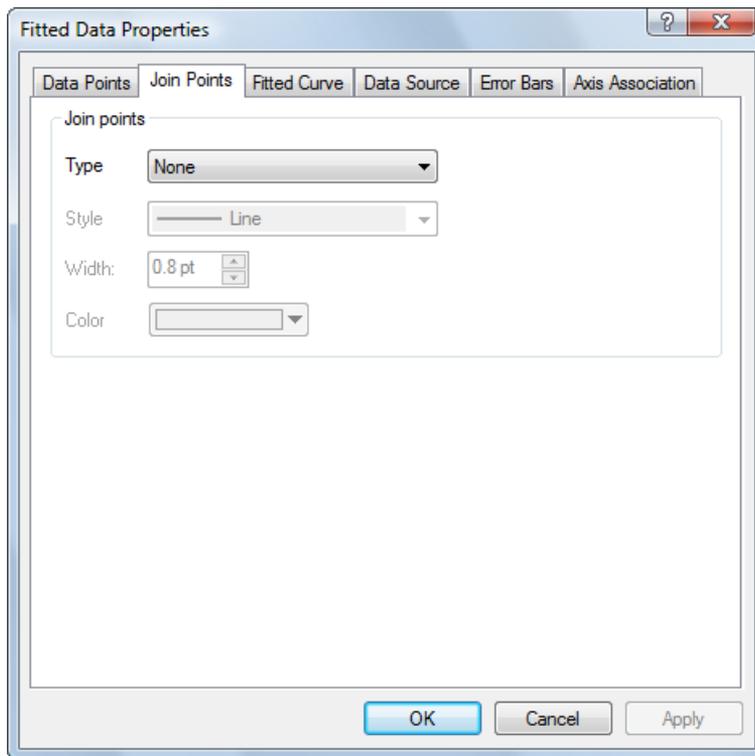

This page is equivalent to the Join Style page of the XY Graph Wizard. See page 164 for further information.

Fitted Curve Page


(fitted data and equations only)

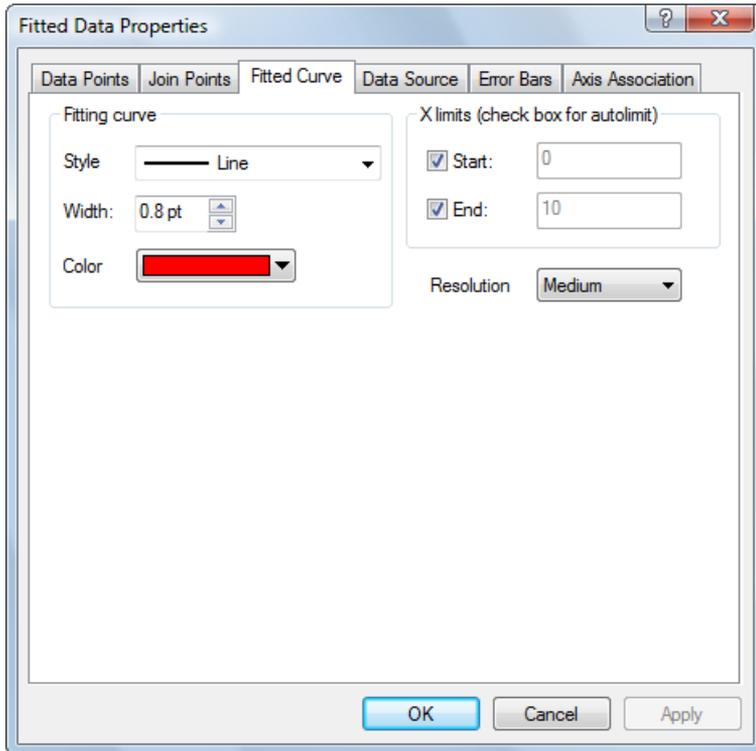

This page is equivalent to the Fitting Curve Style page of the XY Graph Wizard. See page 163 for further information.

Data Source Page

(fitted data only)

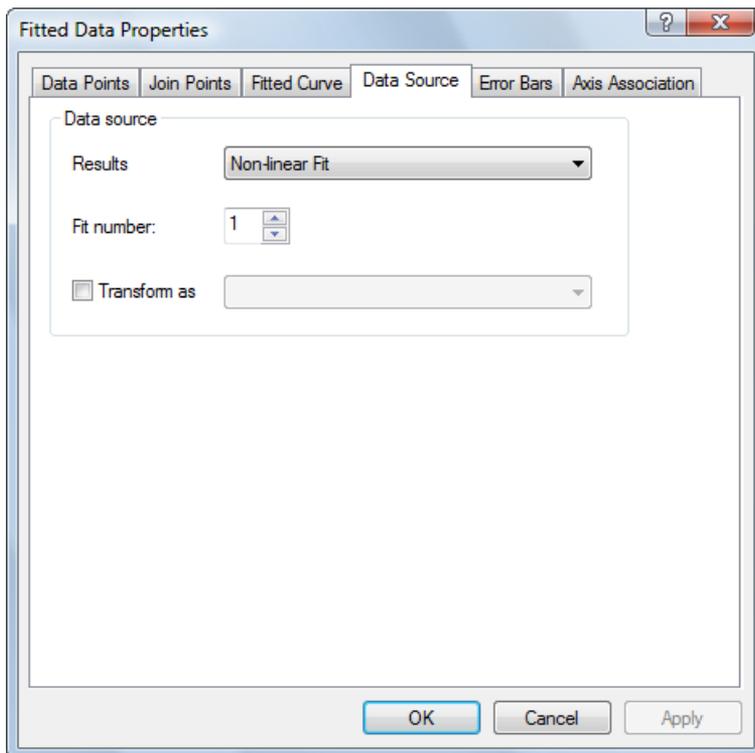


Select the data fitting item to be plotted from the **Results** selection, which lists all the available linear and non-linear fits in the GraFit workspace. Fitting items can hold several fits, and the **Fit number** can be chosen. The first fit is fit 1, which in most cases will be the only fit present.

Check the **Transform as** box if you wish to display the data points and fitted curve transformed as defined in the selected transformation. All the transformations present in the workspace are listed. See Chapter 10 for information about defining transformation items.

Data Source Page


(raw data only)

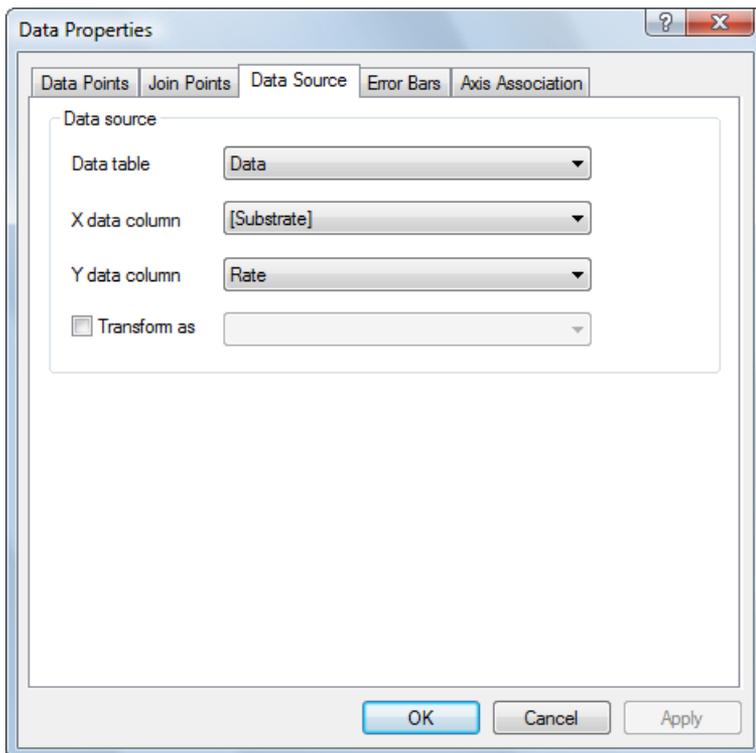

Select the **Data table** item that holds the data, together with the columns that hold the *x* and *y* data.

Check the **Transform as** box if you wish to display the data points transformed as defined in the selected transformation. All the transformations present in the workspace are listed. See Chapter 10 for information about defining transformation items.

Error Bars Page

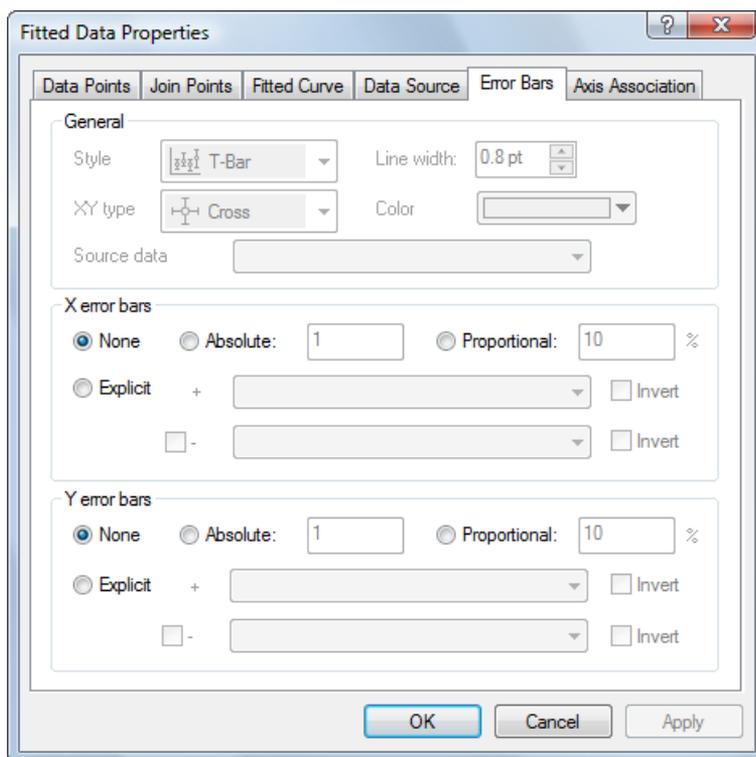

General

This section defines the general types to use for the error bars, including the width and color of the pen used to draw them. The **Style** choice is only available after error bars are selected in one of the **X error bars** or **Y error bars** sections.

Style	Effect
 Line	The error bars are drawn as simple lines.
 T-Bar	The ends of the bars are delimited with a “T” bar.

The **XY type** option is only available if error bars are present on both axes.

XY type	Effect
 Cross	Error bars are drawn for the x and y axes separately.
 Diagonal 1	A single error bar is drawn, which extends diagonally to show the x and y error extents.
 Diagonal 2	A single error bar is drawn, which extends diagonally to show the x and y error extents.


The **Source data** entry must be completed if **Explicit** errors are specified for either x or y errors. Select the data table that is used for the explicit errors from the list given.

X error bars


This section is used to specify the type of errors present on the x axis.

Error type	Effect
None	No error bars are drawn. This is the default setting.
Absolute	The size of the error bar is the same for each data point of the data set; this size is entered into the field to the right of this button.
Proportional	The size of the error bar depends upon the magnitude of the data point; the size is expressed as a percentage of the data point value and is entered into the field to the right of this button.
Explicit	<p>The size of the error for each data point is specified individually. The error values must be stored in a data table, which is specified in the General section, and the size of the positive and negative parts of the error bar is set separately. The data column that holds the positive error values is selected from the list marked "+." If only this option is chosen, the error bar just has a single lobe.</p> <p>To specify additionally the size of the negative lobe, click the "-" button, and select the column that holds the value of this error. If you select the same column as selected for the positive error, this will result in symmetrical error bars.</p> <p>The Invert check boxes allow the direction of the error bars to be reversed, drawing the error bar in the opposite direction.</p>


Some examples of error bar types are shown below.

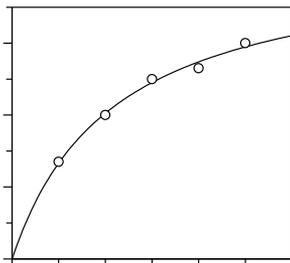


No error bars

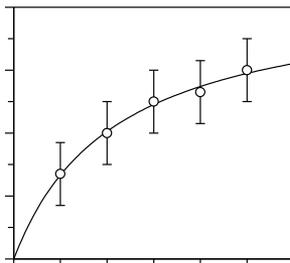


Absolute error bars on the y axis

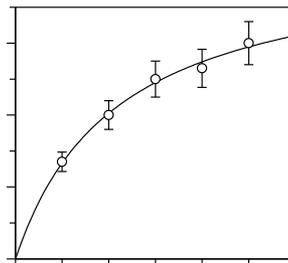


Proportional error bars on the y axis

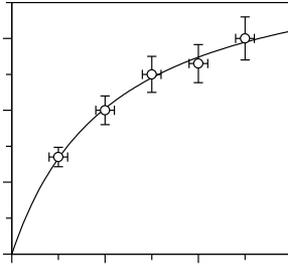

Absolute error bars on the x axis, proportional errors on y

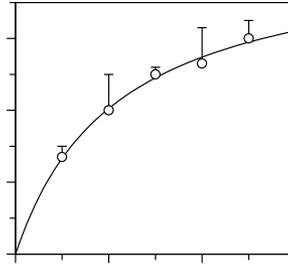

Explicit error bars on the y axis; only the positive lobe is shown

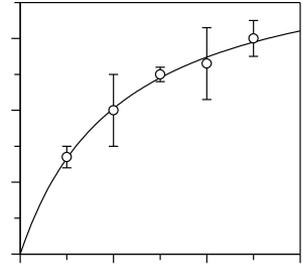


Explicit error bars on the y axis; the positive and negative lobes each use the same error column

Y Error bars

This section is used to specify the type of errors present on the y axis. The settings are analogous to those described above for x error bars.

Axis Association Page

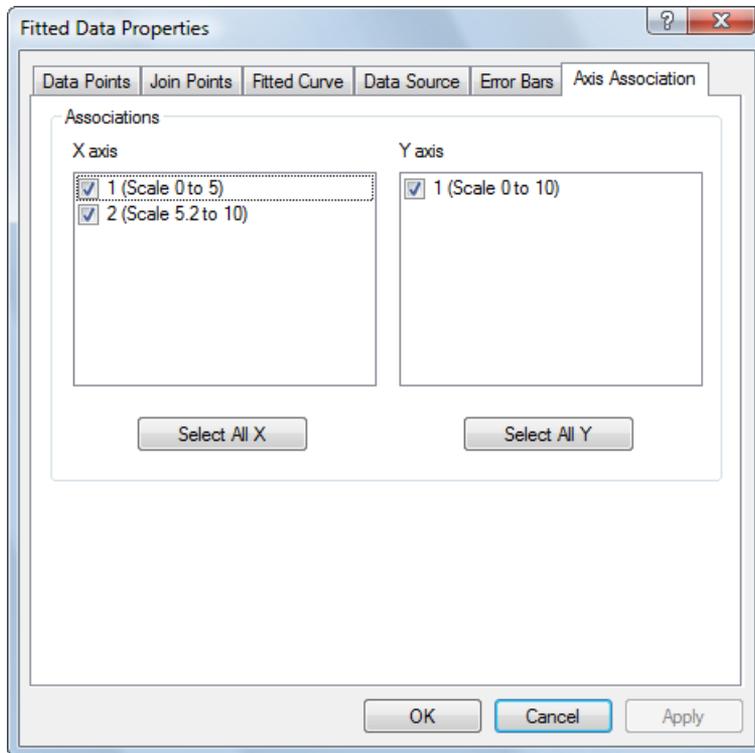


The Axis Association page allows the data to be associated with the various axes that are present (see Chapter 15 for information on multiple axes). XY graphs can contain many

different x and y axes, and the data can be associated with one (the minimum) or many different axis scales. Check the boxes of all axes that you wish to be associated with this data set.

Equation Page


(equation data only)

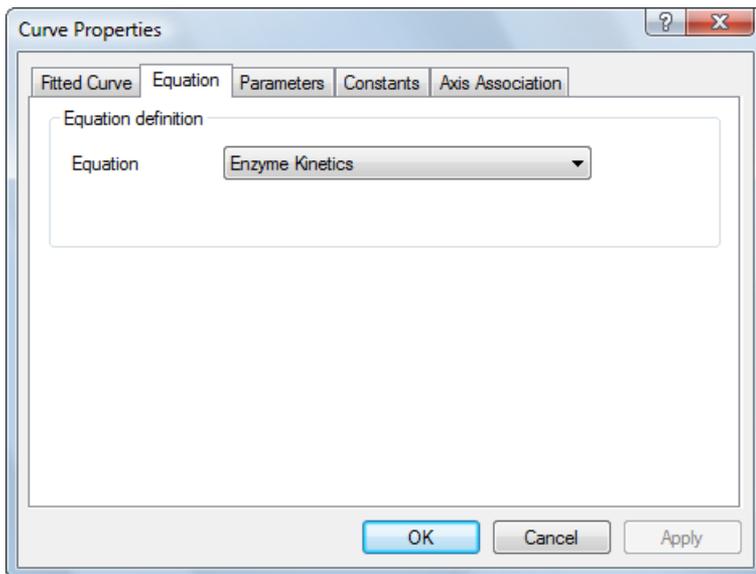

This page selects the equation that describes the curve to be plotted. All equation items present in the GraFit workspace are listed.

Parameters Page


(equation data only)

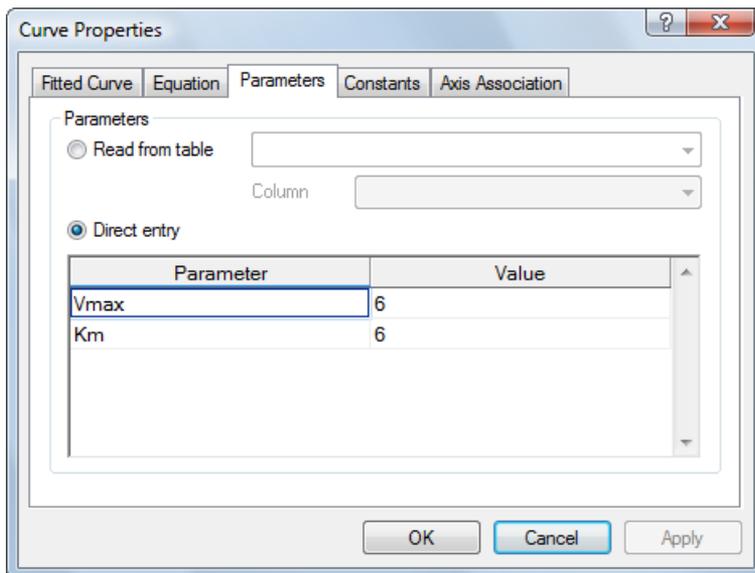

This page is equivalent to the Curve Parameters page of the XY Graph Wizard, see page 162 for further information.

Constants Page

(equation data only)

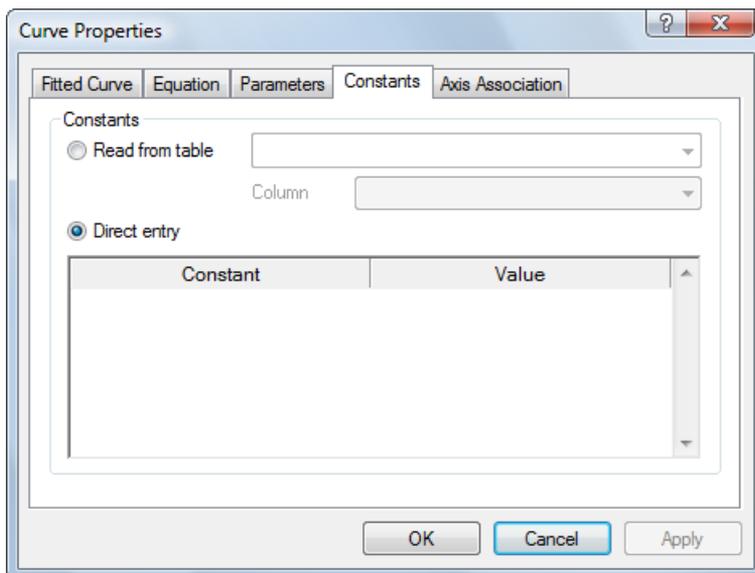

This page is equivalent to the Curve Constants page of the XY Graph Wizard, see page 162 for further information. This page is present even for equations that do not define any constants, as is the case for the figure above.

Legend Text Regions

Legend text properties are controlled in the same way as normal text objects, except that legend text cannot be displayed rotated. See Chapter 13 for further information on text formatting.

Note The overall width of the data legend depends on the size of the individual legend text elements. Each of these can be individually sized using the mouse.


Adding Extra Data to the Graph

Additional data sets, fitted results or equations can be added to the graph by selecting the main graph region, then choosing the **Add Data** command⁷⁰ on the **Graph** tab. This re-invokes the XY Graph Wizard to allow further data to be added. Most of pages are equivalent to those described earlier for creating new graphs, although those pages that control initial scaling are not relevant, and so are not displayed.

Set Axis Associations Page

When graphs are first created, they have a single x and y axis, with which the selected data sets are associated. However, after creation it is possible that a graph may have additional axes, and so the axis association needs to be specified when adding data. This results in an additional Axis Associations page being added to the Wizard.

⁷⁰ The **Graph** ► **Add Data** command in legacy mode

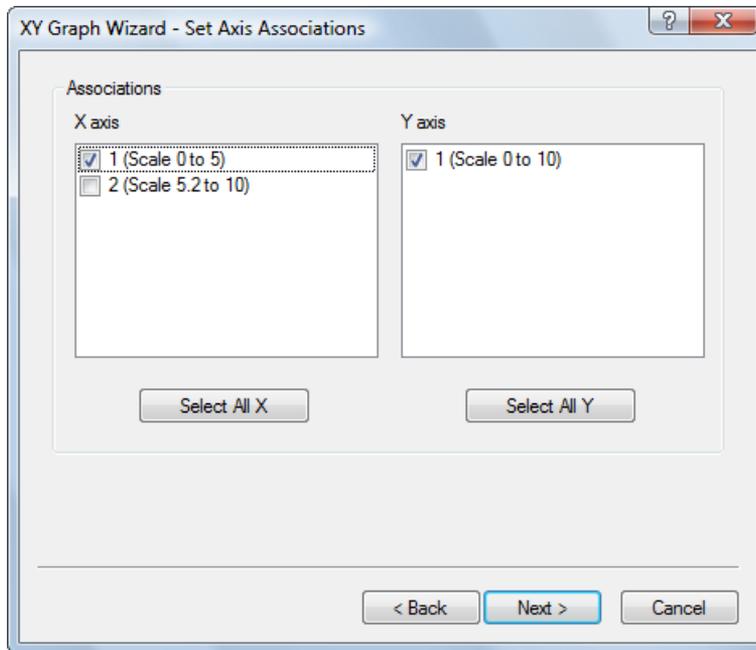

This is equivalent to the Axis Association parameters page described on page 183. Check the boxes of the axes with which your data should be associated. In the figure shown above, the graph has a single x and y axis and so no modification is needed.


Using Multiple Graph Axes

When to Use Multiple Axis Scaling

Graphs created using GraFit can optionally have several different x and/or y axes, which are independently scaled.

Overlaying Data that have Different Scales


Multiple axes are useful when data need to be overlaid, but have quite different numerical ranges. For example, in the plot below there are two separate y axes, but a common x axis. Two different data sets are plotted, each associated with different y scales.

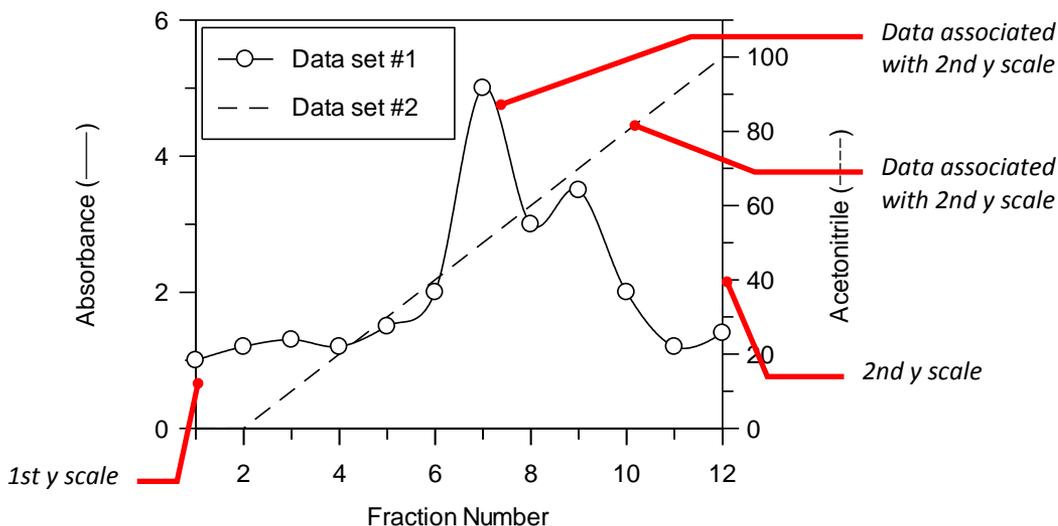

XY scatter graphs can have as many different x or y axes as you require; bar charts and column charts allow unlimited numeric axes (the x and y axes, respectively, for these plots).

Split Axes

Split axes are used when it is necessary to show some region(s) of a data set at a different scale to the remainder. The example drawn below has a single split x axis. GraFit uses multiple axes to construct the splits, and permits any number of splits on both x and y axes.

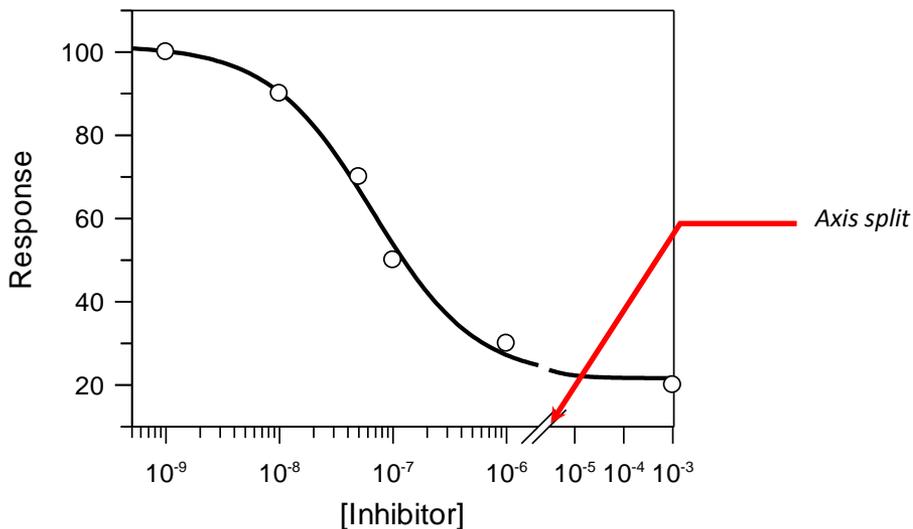


Creating and Managing Multiple Axes

By default, GraFit graphs have just a single x and y axis. Any additional axes must be added to the graph using the Axis Manager, which allows axes to be created or deleted, and simplifies generation of split or multiple axes.

Running the Axis Manager

To activate the Axis Manager, use the mouse to select the graph that you wish to modify, then choose the **Axis Manager** command⁷¹ from the **Graph** tab. This will launch the Axis Manager Wizard.

⁷¹ The **Graph** ► **Axis Manager** command in legacy mode

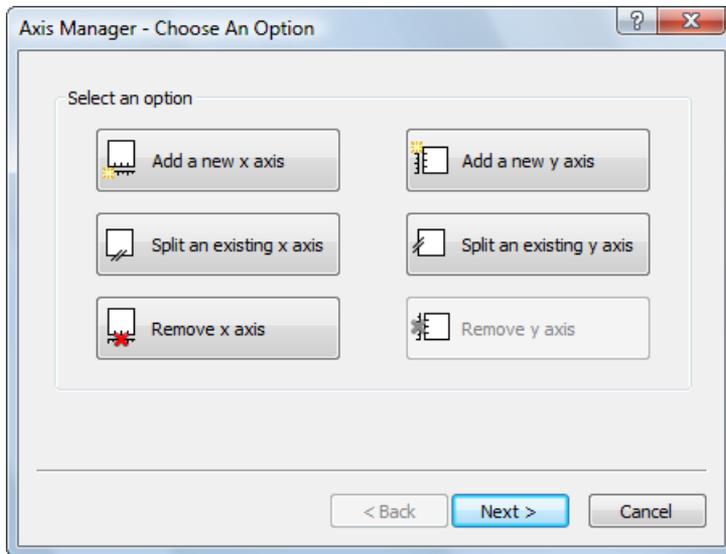


Six options are available on the initial page:

- Add a new x axis
- Add a new y axis
- Split an existing x axis
- Split an existing y axis
- Remove x axis
- Remove y axis

If only a single axis is present, the **Remove axis** options will not be available, as a minimum of one axis must exist. Select one of the options, and click the **Next** button.

Adding a New Axis


Choosing to add a new axis will give the following page.

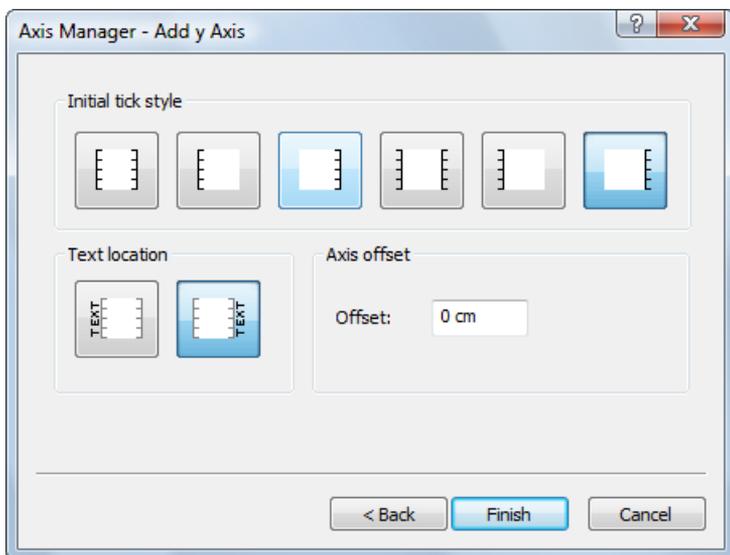


A number of initial tick styles and axis text locations can be chosen, and in addition it is possible to specify that the new axis is created offset from the main graph. All new axes have associated axis and scale text, which needs to be positioned as appropriate. In addition, the new axis ticks are set and scaled independently of any tick settings already present on the graph. When you have made your initial choices, press the **Finish** button to add the new axis.

Splitting an Existing Axis

Choosing to split an existing axis from the Axis Manager will give the following page.

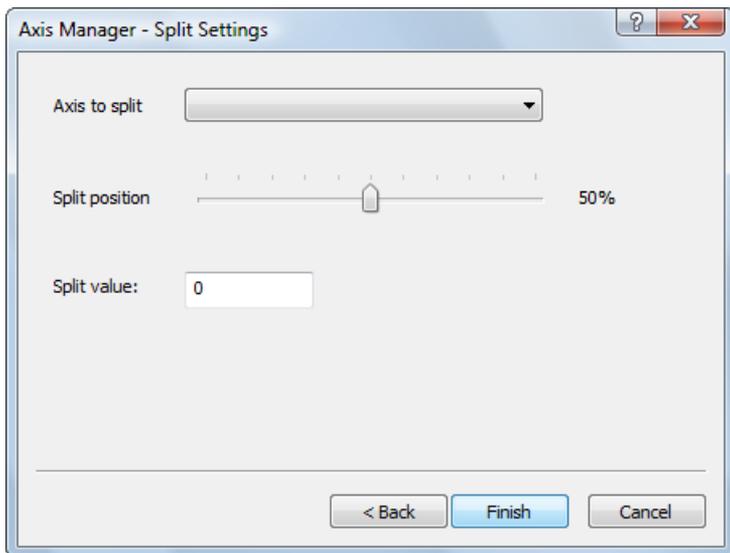


The axis to be split should be selected from the list available. The initial split position can be set using the slider control. By default, the axis is split in the middle, i.e. at 50%. Use the slide to set this to an appropriate position. The axis is split so that the scale for two new axes is set as below, using the **Split value** to set the end of the first scale and the start of the second.

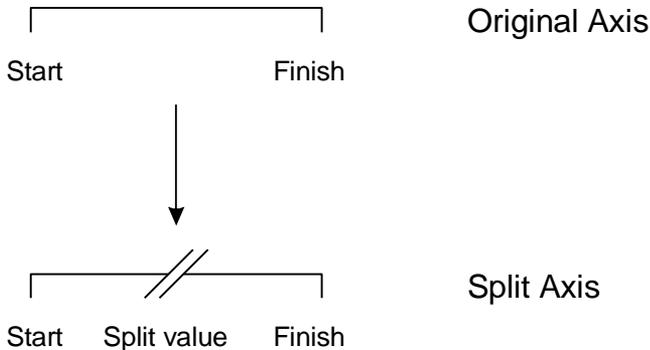

After entering a suitable **Split value**, click the **Finish** button to create the new axis.

Removing an Axis

If you elect to remove an axis using the Axis Manager, you will be presented with a list of all axes that can be deleted.

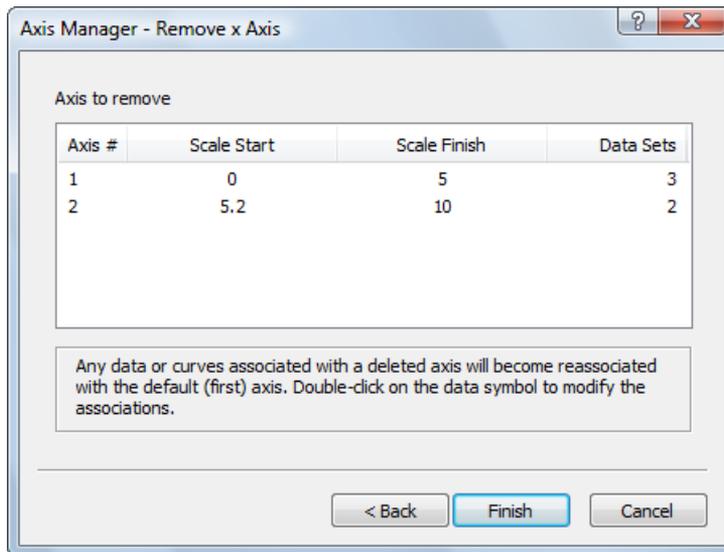


Choose the axis to remove, then click the **Finish** button to delete it from the graph. Any data or curves that are associated with the deleted axis become re-associated with the default (first) axis. You can modify axis associations by editing the properties associated with the data

set. To do so, select the data symbol (found in the data legend), and choose the **Properties** command⁷² from the **Graph** tab (or double-click on the data symbol).

Working with Multiple Axes

Axis Text and Axis Scale Text

Each axis has a separate text legend that can be positioned and edited as required, and a separate scale text that can also be formatted as required.

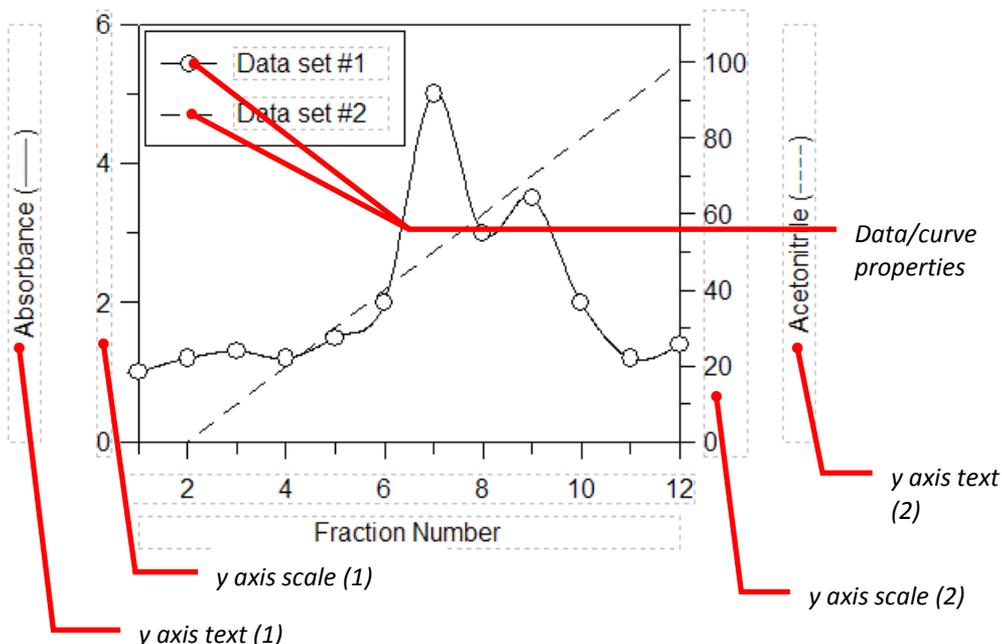

Visibility of Axis Text and Axis Scale Text

There are often cases where it is not necessary to show a separate axis text and/or axis scale text for each axis. The various items can be hidden by selecting the item, then choosing the **Hide Object** command⁷³ on the **Graph** tab. Objects can also be hidden and unhidden from the General page of the graph properties (double-click on the main graph area, or choose the **Properties** command⁷² when the main graph is selected).

Data Association

Any dataset that is added to a graph can be associated with one or all axes. When using multiple axes to display data that have distinct ranges, it is usual to associate each dataset with one specific axis. For split axes, it is usual to associate all datasets with each axis. The data association for an axis is set from the Axis Association Page, found by double-clicking on

⁷² The **Edit ▶ Properties** command in legacy mode


⁷³ The **Edit ▶ Hide Object** command in legacy mode

the appropriate Data/Curve Properties region of the data legend (see figure above). Further information about modifying data properties is found in Chapter 14.

Axis Scaling and Sizing

Each of the multiple axes has separate scaling, size and tick settings. These are modified from the Scale or Ticks pages of the graph properties. To access these, select the graph and choose the **Properties** command⁷², or double-click on the main graph region. For more detailed information about setting these properties, see Chapter 14.


Scaling

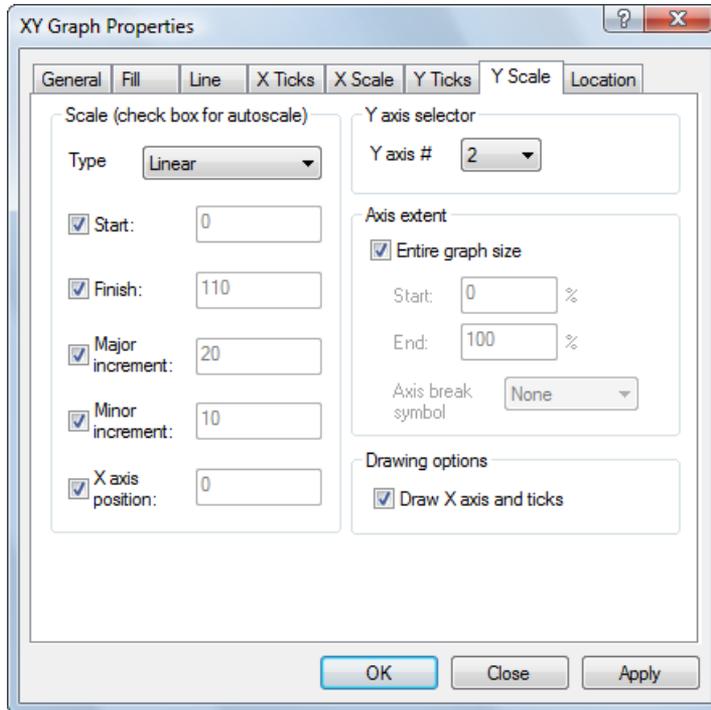

Use the **Axis selector** to choose which axis is to be scaled, then set the scale as appropriate.

Sizing

The **Axis extent** settings on the page shown above control the proportion of the graph over which the current axis extends. When using multiple axes to show data of with different numerical ranges, it is normal to have the axis set to occupy the **Entire graph size**. However, in some circumstances it can be useful to limit the axis to a small region. For example, this can be a good way to overlay residuals, as illustrated below.

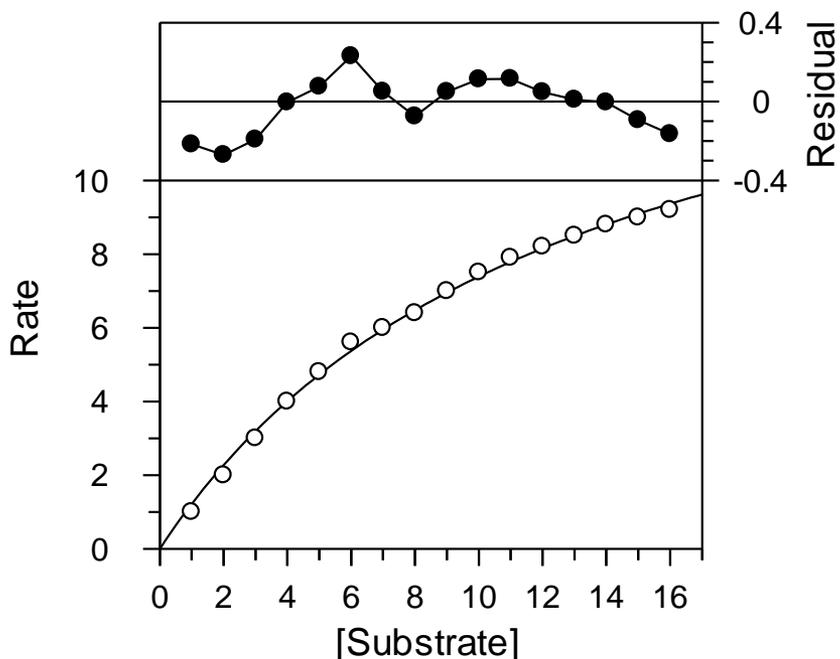


In this graph there are two y axis scales, the first which holds the normal fitted data, and the second that is used to plot the residual error. The first y scale is set to occupy from 0 to 70% of the graph size, the second is set from 70 to 100%.

Setting the Split Position

Split axes are created by using axis extents less than the full graph size, and by including a suitable **Axis break symbol**. In order to modify the split position along the axis, it is necessary to change the **Start** and/or **End** position of the **Axis extent**. This must be done for each axis concerned.

Tick Settings

Each axis has full control over where, and how often, ticks are drawn. For this reason, when a new axis is created it often the case that the initial tick settings overlap existing settings from existing axes. If this is the case, it is necessary to modify the tick settings for one or more axis in order to obtain a suitable appearance.

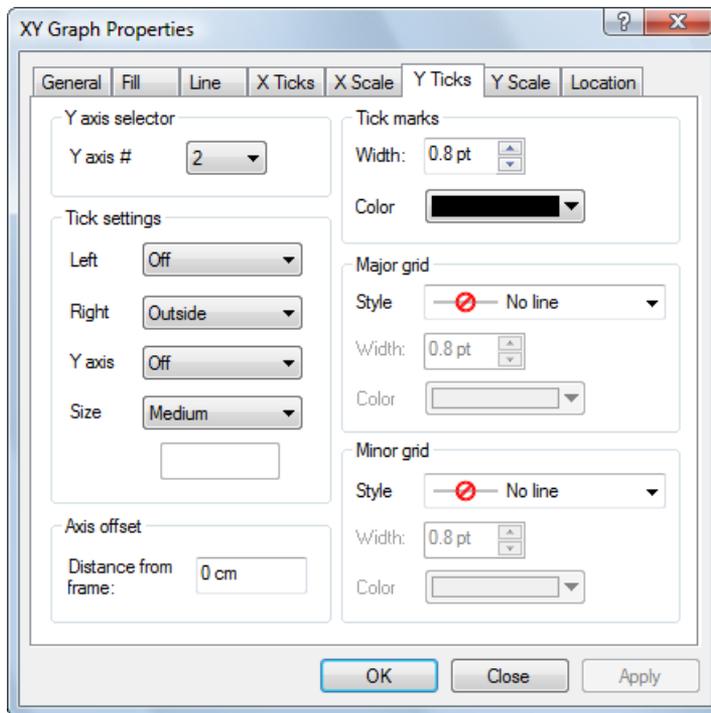


The **Ticks** page also includes an **Axis selector** section, which allows each different axis to be selected. After choosing the axis to modify, alter the tick settings as required.

Column and Bar Charts

The Chart Wizard

Chapter 12 described the basics of creating graphs. Selection of Column Chart or Bar Chart produces the “Chart Wizard” to assist with the setup (all settings can be adjusted later, if required). Column and bar charts differ only in their orientation, and so the same process is carried out for each. The process of editing charts is broadly similar to editing X/Y graphs, and Chapter 14 should be consulted for further information.

Select Data Source Page

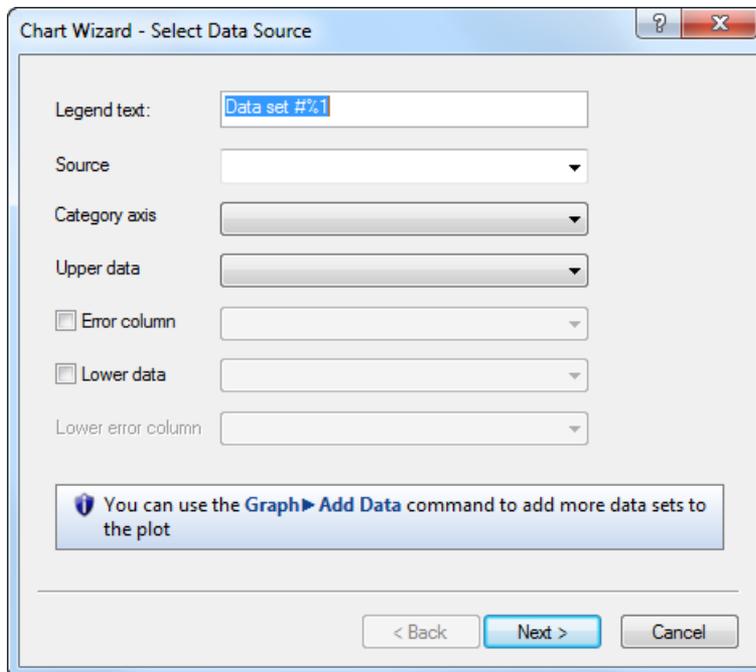

The screenshot shows the "Chart Wizard - Select Data Source" dialog box. It has a title bar with a question mark and a close button. The main area contains several fields and checkboxes:

- Legend text:
- Source:
- Category axis:
- Upper data:
- Error column:
- Lower data:
- Lower error column:

At the bottom, there is a blue box with a shield icon containing the text: "You can use the **Graph ► Add Data** command to add more data sets to the plot". Below this box are three buttons: "< Back", "Next >", and "Cancel".

The Select Data Source page allows you to choose the source and type of data to plot.

Legend text

The data set plotted is associated with an entry in the data legend. The entry shows the symbol used to represent the data, and a descriptive text. The initial text used to describe the entry can be entered into this field.

Source

Graphs are plotted using data that are present in one of the items found in the GraFit workspace. The **Source** selection lists all items that contain suitable data for plotting, and one of these must be selected.

Category axis

Select which column holds the category names that are to be plotted.

Upper data

The **Upper Data** selection is used to choose the data column that holds the values to be associated with the selected categories.

Error column

Checking the **Error column** box allows a quick way of producing an error plot. Select the column that holds the error values from the list. More extensive error plotting options are available once the graph has been produced.

Lower data

Check the **Lower data** box to generate a High-Low chart (where the entries have both an upper and lower value). If this box is unchecked, a standard chart will be produced; if checked then it is necessary to choose the column that holds the lower data values for the High-Low chart.

Lower error column

When the **Error column** and **Lower data** boxes are both checked, the error values for the lower data are selected from this box.

Set Initial Scaling

Chart Wizard - Set Initial Scaling

Data axis

Type: Linear

Autoscale

Manual scaling:

From: 1

To: 7

Data axis text

Text: Values


Category axis text

Text: Categories

< Back Next > Cancel


The Set Initial Scaling page allows the scale of the data axis to be specified. By default the scale is set automatically, and a linear data axis is chosen. These can be altered as required (it is also possible to adjust the scale settings later). The text used for the data and category axes can also be entered.

Select Tick Styles

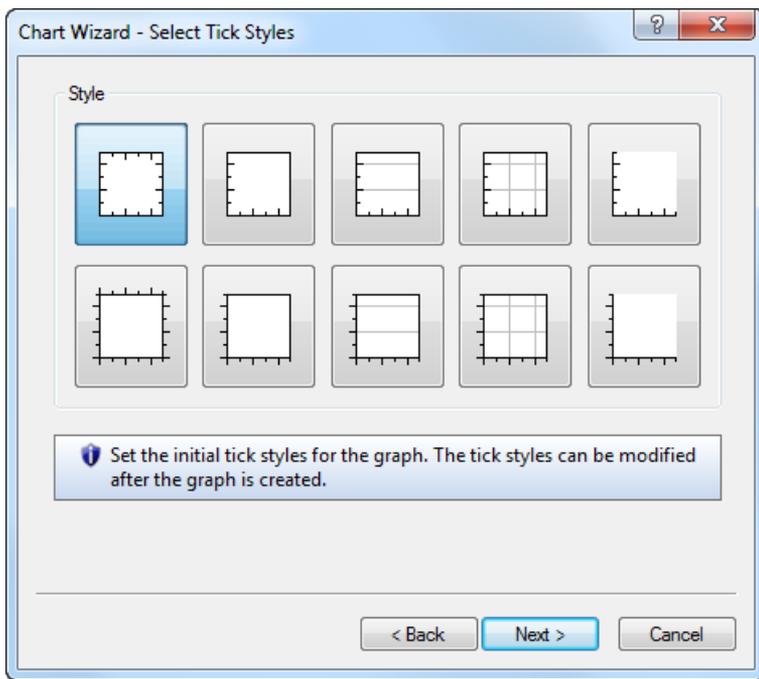


The Select Tick Styles page allows rapid selection of some of the more common axis tick styles. More extensive tick style choices can be made once the graph has been created.

Select Bar Style

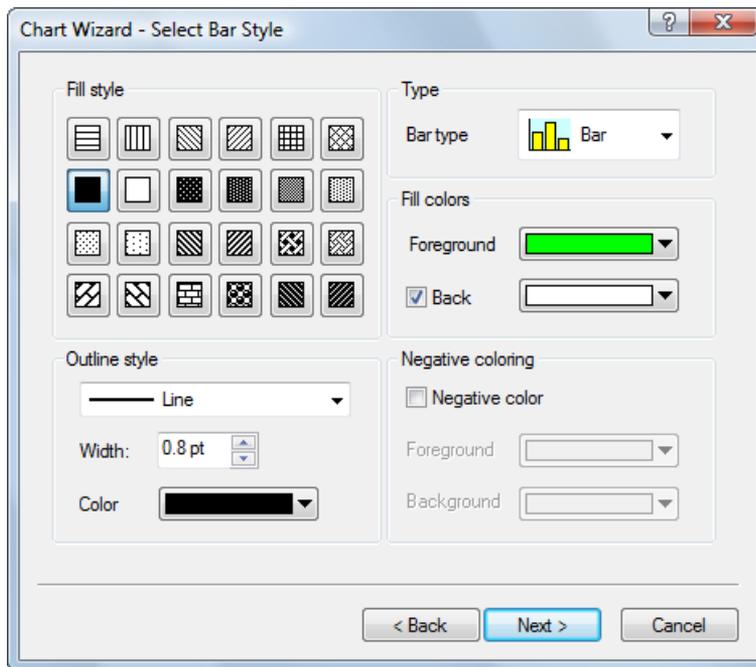


The Select Bar Style page allows the chart bar representation to be chosen.

Select Data Point Style

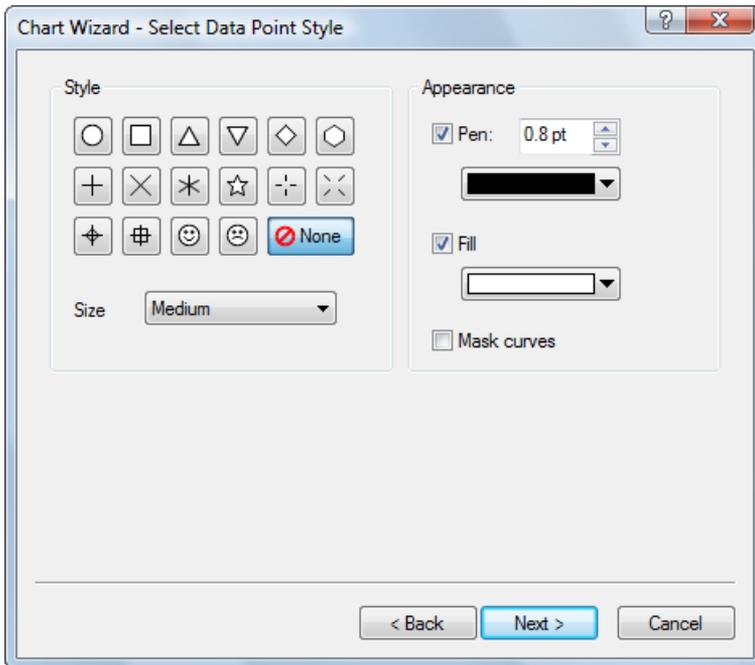

It is possible to use data points to represent the values, if required. By default, no data point style is selected.

Select Join Style

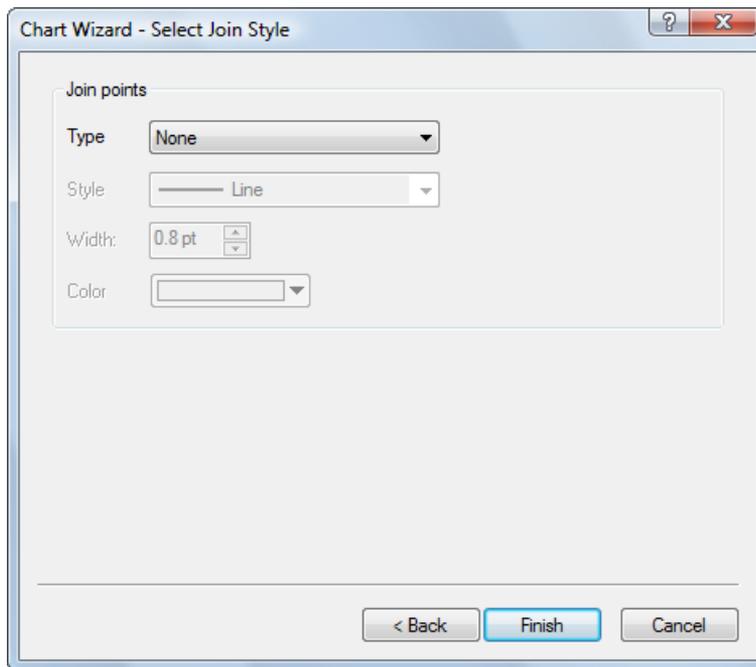


It is possible to draw lines that join the points, if required. By default, no lines are drawn. Click the **Finish** button to create the chart.

Histograms

The Histogram Wizard

Chapter 12 described the basics of creating graphs. Selection of Histogram produces the “Histogram Wizard” to assist with the procedure (all settings can be adjusted later, if required). The process of editing histograms is broadly similar to editing X/Y graphs, and Chapter 14 should be consulted for further information.

Select Data Source Page

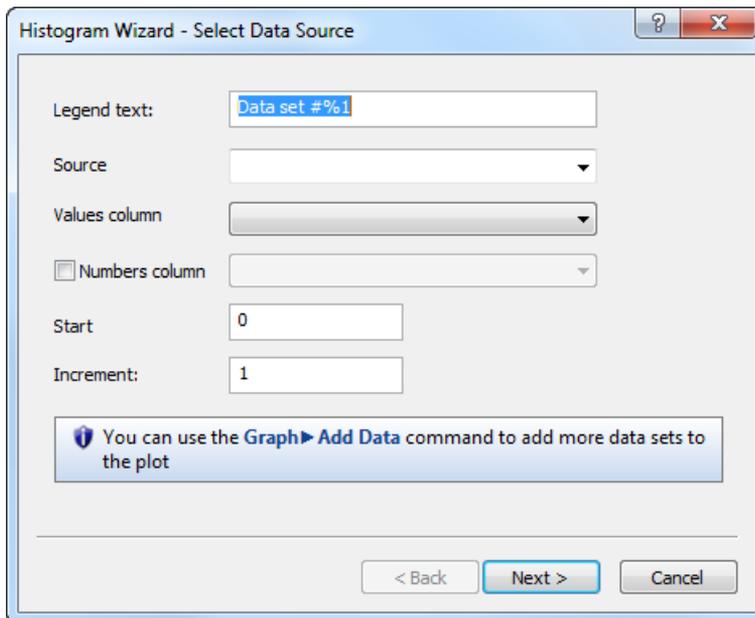

The screenshot shows a dialog box titled "Histogram Wizard - Select Data Source". It contains several input fields and a message box. The "Legend text" field contains "Data set #%%1". The "Source" field is a dropdown menu. The "Values column" field is a dropdown menu. The "Numbers column" field has a checkbox and a dropdown menu. The "Start" field contains "0" and the "Increment" field contains "1". At the bottom, there are three buttons: "< Back", "Next >", and "Cancel". A message box at the bottom states: "You can use the Graph ► Add Data command to add more data sets to the plot".

The Select Data Source page allows you to choose the source and type of data to plot.

Legend text

The data set plotted is associated with an entry in the data legend. The entry shows the symbol used to represent the data, and a descriptive text. The initial text used to describe the entry can be entered into this field.

Source

Graphs are plotted using data that are present in one of the items found in the GraFit workspace. The **Source** selection lists all items that contain suitable data for plotting, and one of these must be selected.

Values column

Select which column holds the data values whose distribution is to be plotted.

Numbers column

Check this box if a second column is present that holds a count of the number of items an item in the **Values column** is present. If unchecked, that it is assumed that each entry in the **Values column** is represented once only; if checked then the column that holds the number of entries must be selected.


Start position

Enter the starting position for the histogram grouping. The histogram records the number of occurrences of data values within a specified data range. This range is defined as having an origin at the **Start position** value, and the group size is set by the **Increment** value.

Increment

Enter the histogram increment value, as defined above.

Set Initial Scaling

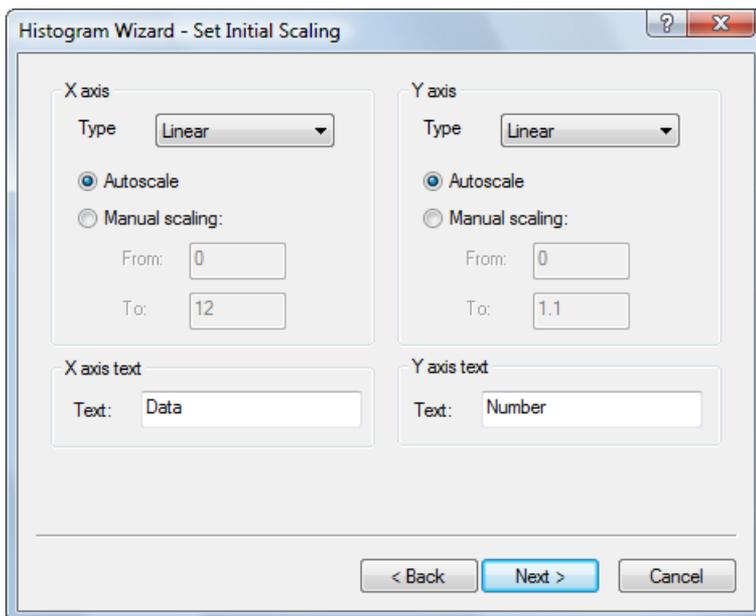


The screenshot shows the 'Histogram Wizard - Set Initial Scaling' dialog box. It is divided into two main sections for the X and Y axes. Each section has a 'Type' dropdown menu set to 'Linear'. Below each dropdown are two radio buttons: 'Autoscale' (selected) and 'Manual scaling:'. Under 'Manual scaling:', there are 'From:' and 'To:' input fields. For the X axis, 'From:' is 0 and 'To:' is 12. For the Y axis, 'From:' is 0 and 'To:' is 1.1. Below these sections are 'X axis text' and 'Y axis text' labels with corresponding text input fields containing 'Data' and 'Number' respectively. At the bottom of the dialog are three buttons: '< Back', 'Next >', and 'Cancel'.

The Set Initial Scaling page allows the scale of the data axis to be specified. By default the scale is set automatically, and linear axes are chosen. These can be altered as required (it is


also possible to adjust the scale settings later). The text used for the x and y axes can also be entered.

Select Tick Styles

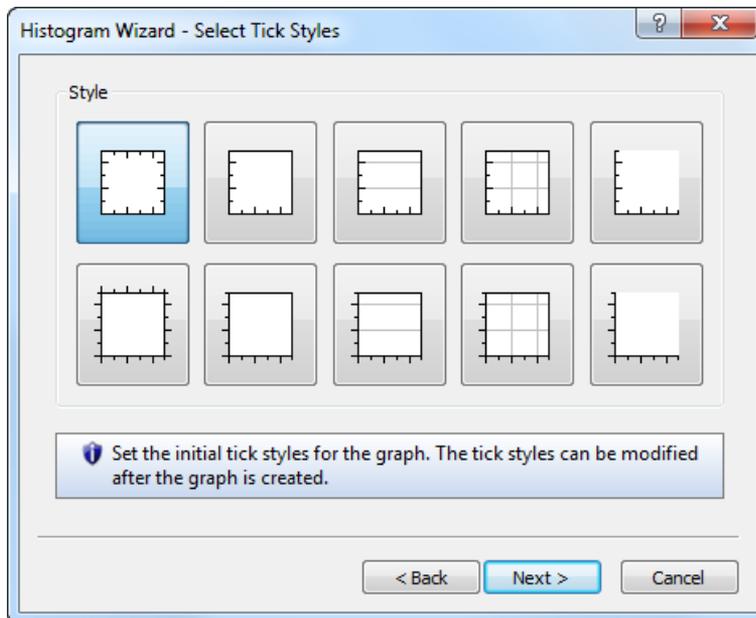

The Select Tick Styles page allows rapid selection of some of the more common axis tick styles. More extensive tick style choices can be made once the graph has been created.

Select Bar Style

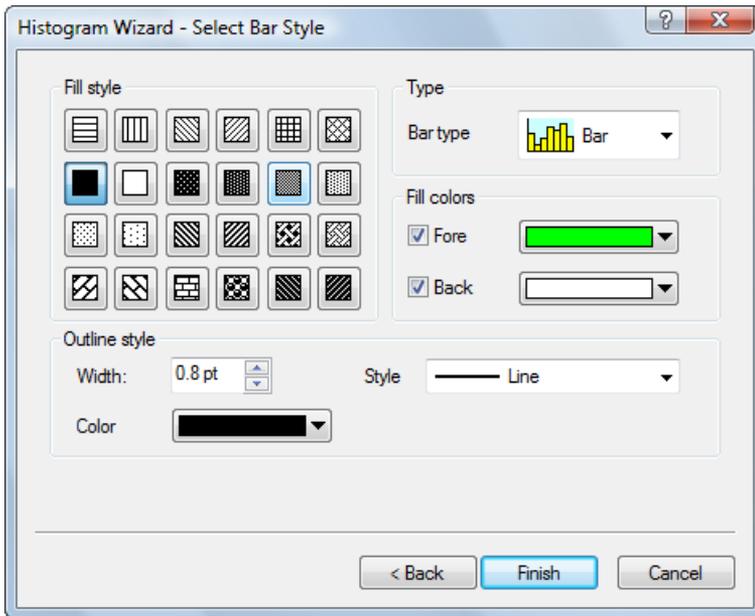

The Select Bar Style page allows the histogram bar representation to be chosen. Click the **Finish** button to create the chart.


Embedding and Linking

What is Embedding and Linking?

Embedding and linking are two related ways that information from one application can be incorporated into a different program. Also known as Object Linking and Embedding (OLE), it provides a uniform mechanism for Windows programs to cooperate and to share information. Most Windows applications allow information to be copied between programs using the Clipboard. For example, a picture could be copied from one application and pasted into a second. The problem with simple cut and paste is that the second application often has no idea what information is being transferred. In non-OLE applications, this often means that once the information is copied it is impossible to make any subsequent changes to the image. OLE gets around this by allowing two cooperating applications to pass information between themselves in such a way that editing is still possible, even after pasting into an application that has no knowledge about the information transferred. Both embedding and linking are extensions of the standard Clipboard cut and paste mechanism.

Original Cut and Paste Model

The original use of the Clipboard was to pass a representation of the object being copied. For example, copying a graph might place an image of the graph on the Clipboard. Pasting would then insert this picture, as distinct from inserting the information required to draw the graph.

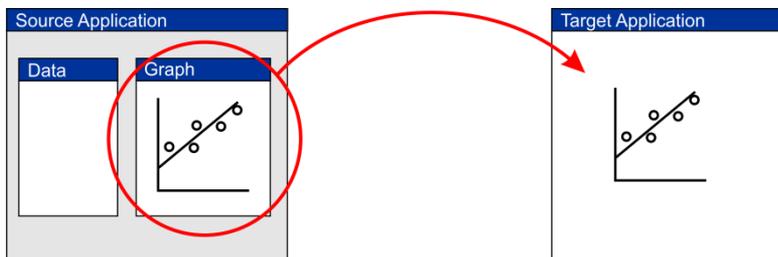


Standard cut and paste inserts an image representing the source object into the target application. This image is simply a picture, and contains no information about the underlying data in the source application.

Embedded Objects

Embedded objects hold both a display representation **and** the underlying data. A copy of the current state of the source object is created and is passed to the target application. This


means that subsequent changes within the source application do not affect the newly pasted object. In turn, although the embedded object can be edited, any changes are not reflected in the original source document. Because a complete copy of the source data is placed in the target application, files containing embedded objects can be large.

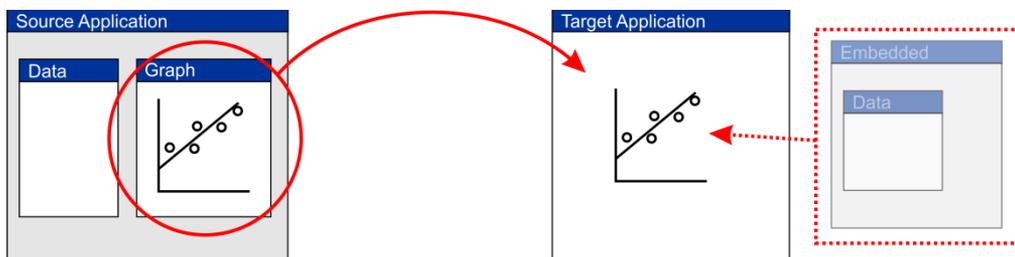

Embedding places a copy of the internal data into the target application, as well as an image representing the source.

Linked Objects

Linked objects also copy a representation of the original data, but do not copy the data itself to the target. Instead, the target application maintains a link to the original source file. Changes to this source data file result in the information copied to the target application also changing. In addition, if the linked object is edited from within the target, this will change the original document. Because linking must always involve files, it is only possible to create a link after the information in the source application has been saved to disk.

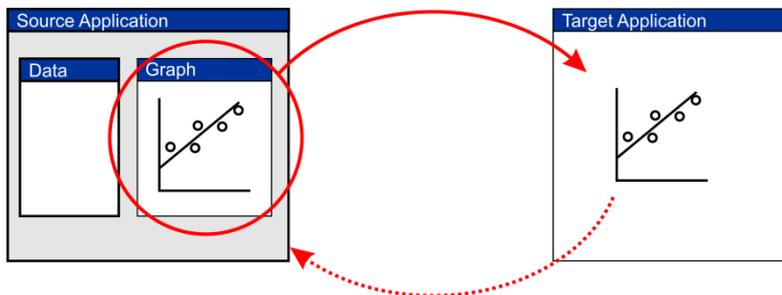

Linking does not copy the original data, but allows the target application to maintain links to the original source file. Changes to the source result in the target application being updated.

Using Embedding and Linking

Windows applications that are able to employ OLE use the convention that **Edit ▶ Copy** from the source program will place the necessary information on the Clipboard, and **Edit ▶ Paste** from the target application will embed an object. Therefore, it is the standard that embedding is used by default. To link to an object, the same **Edit ▶ Copy** command is used in the source program, but, depending upon the target application, the object is inserted with the **Edit ▶ Paste Link** or **Edit ▶ Paste Special** command.

The OLE object can be activated by double-clicking on its representation within the target, or by accessing the object menu that resides on the **Edit** menu (and is often also available by right-clicking on the object).

Comparison of Data Transfer Types


The following table compares the various ways that data can be transferred between applications via the Clipboard.

	Non-OLE Clipboard transfer	Embedding	Linking
How to paste information in this format	For non-OLE applications, Edit ► Paste For OLE applications, Edit ► Paste Special and select a non-OLE format	For OLE applications, Edit ► Paste	For OLE applications, Edit ► Paste Link or Edit ► Paste Special and choose to paste as a link
Advantages	Compact data format	Allows information to be edited after pasting Pasted information is stored within the target document	Allows information to be edited after pasting Since data are kept within the source file, this allows efficient use of space when several documents may wish to link to the same object
Disadvantages	No editing allowed after pasting	No link maintained with the source document Documents containing embedded objects tend to be large	Linked information is kept in a separate file, which may become “lost”

Note The default action for a Windows application that supports OLE is to embed information whenever possible if the **Edit ► Paste** command is used. As described in the table above, embedding is not always the most appropriate way of sharing the information, and linking is frequently more desirable, for which you need to select the **Edit ► Paste Special** command.


Embedding and Linking without using the Clipboard

Data objects can also be inserted into a target application directly, using an **Insert Object** command, which is often found on the **Edit** menu.

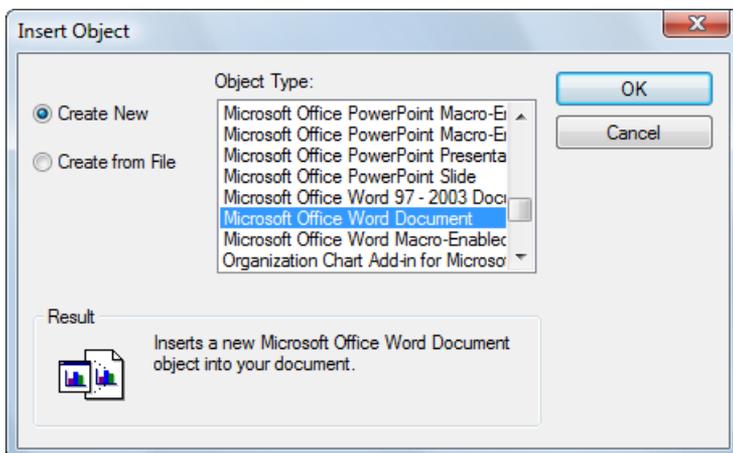

The Insert Object command allows the insertion of OLE objects into a document

When the **Create New** button is selected, an empty object of the selected type is inserted as a new embedded OLE item. If the **Create from File** button is chosen, it is possible to initialize the embedded item using the contents of an existing file, or to create a linked object.

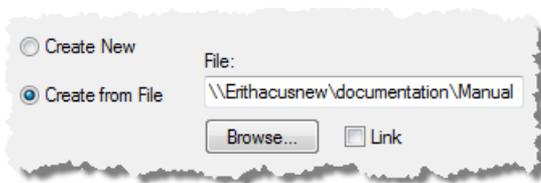

*The Create from File option allows creation of an embedded object based on the contents of an existing file. If in addition the **Link** button is checked, a linked object will be created rather than an embedded one*

Embedding and Linking with GraFit

GraFit can act as both a container of OLE items that come from other source applications, and as the source of OLE items that are inserted into other targets. The situations where GraFit can act as container (target) or source are discussed below.

GraFit as an OLE container

There are several places where it is possible to insert OLE items into GraFit.

GraFit Project Window

The main GraFit project window can hold OLE objects, which are inserted from the Clipboard using the **Paste** or **Paste Special** commands⁷⁴ from the **Home** tab, or created anew using the

⁷⁴ The **Edit** ► **Paste** or **Edit** ► **Paste Special** commands in legacy mode

Object ► Insert New Object command⁷⁵ on the **Home** tab. Objects inserted into the project window are always displayed as icons. Activation of these objects will always launch their source application in a separate window.

Graph Window

All graph windows can hold OLE objects, which can be inserted from the Clipboard as usual or created using the **Object ► Insert New Object** command. Objects inserted in a graph window can be moved and sized as required. In-place editing is available by double-clicking on the item. The size of the object can be reset using the **Object ► Reset Original Size** command⁷⁶ from the **Home** tab.

Equation and Transformation Definition Windows

These GraFit item windows hold textual comments, which can incorporate embedded OLE objects. See Chapters 9 and 10 for further details.

GraFit as an OLE source

The contents of a GraFit workspace can be embedded or linked into other applications. What is displayed is the contents of a particular graph window, although the entire workspace information is copied. To embed a GraFit graph into another application, first activate the graph to use, then choose the **Select All** command⁷⁷ from the **Home** tab to ensure that all graph items are selected.

Note It is extremely important to use the **Select All** command, as no OLE information is copied unless **all** items are selected). When all objects are selected, the **Copy** command is replaced by **Copy All**.

Next choose the **Copy All** command⁷⁸. Insert the information into an OLE-compliant application using **Edit ► Paste** or **Edit ► Paste Special** as appropriate to the application in question.

⁷⁵ The **Edit ► Insert New Object** command in legacy mode

⁷⁶ The **Edit ► Reset Original Size** command in legacy mode

⁷⁷ The **Edit ► Select All** command in legacy mode

⁷⁸ The **Edit ► Copy All** command in legacy mode

Data Fitting Basics

This chapter describes some of the basic principles behind data fitting. A more complete description of the theory is given in Chapter 20, and additional information can be found in several reviews.⁷⁹⁻⁸⁴ Within GraFit, data are fitted using the Non-linear Fit or Linear Fit item windows, which are described in Chapters 6 and 8, respectively.

When Is Data Fitting Appropriate?

Data fitting is the mathematical analysis of a data set in order to analyze trends in the data values. This will usually involve linear or non-linear regression analysis of these data values, in order to define a set of parameter values that best characterize the relationship between the data points and an underlying theoretical model. While many types of data can be usefully presented in graphical form, not all such data are suitable for data fitting.

Types of Data

In general, data can be divided into three groups.

1. No Systematic Trend

There are many examples of such data, and two are displayed below.

⁷⁹Cleland, W.W. (1967) *Adv. Enzymol.* **29**, 1-32.


⁸⁰Bevington, P.R. (1969) *Data Reduction and Error Analysis for the Physical Sciences* McGraw-Hill.

⁸¹Duggleby, R.G. (1981) *Analytical Biochemistry* **110**, 9-18.


⁸²Press, W.H., Teukolsky, S.A., Vetterling, W.T. & Flannery, B.P. (1992) *Numerical Recipes in C*. Cambridge University Press.

⁸³Motulsky, H.J. & Ransnas, L.A. (1987) *FASEB J.* **1**, 365-374.

⁸⁴Leatherbarrow, R.J. (1990) *TIBS* **15**, 455-458.

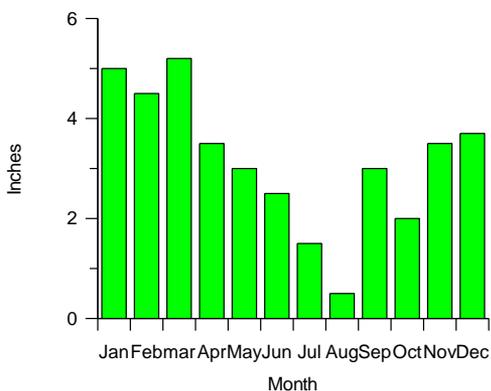

A rainfall chart

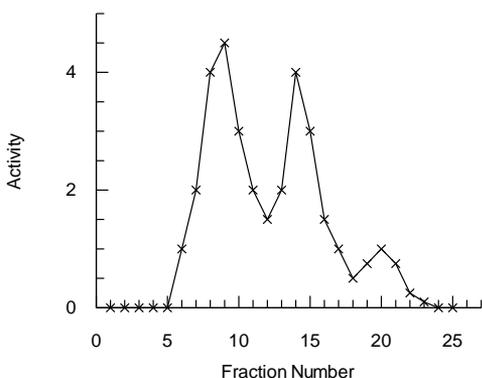

An elution profile from a purification column

Such data are typically presented as bar or column charts, or as scatter plots with the data points joined by lines (or possibly by a smooth curve). As no systematic relationship describes the data, they are not suitable for data fitting.

2. The Trend is described by a Known Equation

This is the situation where data fitting is most useful. In this case, based on knowledge of the experiment, and the laws that govern it, the data are described by a known equation, i.e.

$$y = f(x, p_1, p_2, p_3, \dots)$$

In non-mathematical terms, the observed data, y (also known as the dependent variable), varies as a function of some experimentally controlled value x (the independent variable) and one or more parameters p_1, p_2, p_3 , etc. The exact form of the equation depends on the experiment and the conditions that are used: it might be a simple linear equation,

$$y = ax + b$$

or a more complex function. In either case, the purpose of data fitting is find values for the parameters in this equation that best describe the data. For a simple linear equation, as shown above, this is done by linear regression. In general, however, a more complex equation will apply. Such data are fitted by non-linear regression. GraFit allows non-linear equations of arbitrary complexity to be used for data fitting.

3. A Systematic, But Uncharacterized Trend is Present

Occasionally, data may be collected that show a distinct trend, but no known equation describes this relationship. For example, an empirically determined calibration curve may be collected. If the physical basis for the shape of the data is not well defined, the appropriate equation is unknown. In such situations it can be useful to fit the data to a general-purpose equation, such as a polynomial. This will allow values to be extrapolated from the standard curve, although any parameter values that are obtained will not have physical meaning. The use of data fitting for this purpose is illustrated in Chapter 21.


Regression Analysis

Regression analysis finds the “best fit” line or curve through a series of data. This is considered to be the one that minimizes the sums of the squares deviations of the experimental data points from the theoretical curve. GraFit uses regression to allow the analysis of experimental data.

Practical Considerations

The Errors are in the y Data

It is the convention to assume that the x data values represent the independent variable, and the y data values represent the dependent variable, i.e. the x data are what is varied in the experiment (for example time, concentration, etc.), and the y data are the experimentally observed values.

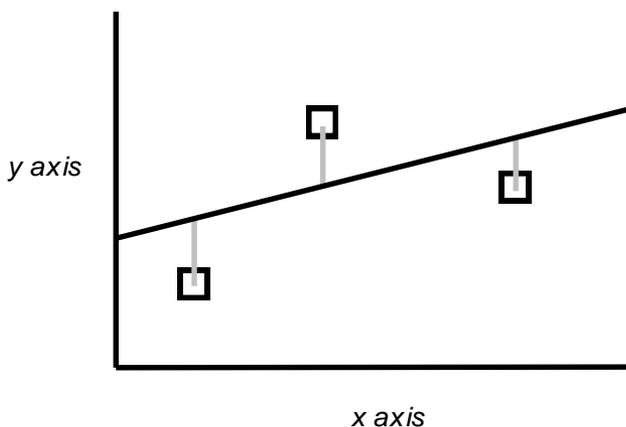

The vertical distance (the gray lines in the above diagram) between the data point and the curve or line is that which is minimized by normal regression analysis, i.e., it is assumed that all errors are present in the y data. An important aspect of experimental design is to ensure that the x data are as accurate as possible so that this assumption is valid.

The Data Values Have a Suitable Numerical Range

A consequence of computer calculation is that rounding errors can occur, particularly if very large or very small numbers are used. This is because computers only store decimal numbers with a limited numerical precision. Although this should not affect most calculations, it is

advisable to scale any set of very large or very small values prior to performing calculations. For example, if a series of time values were entered as 0 to 10^{-15} seconds, the rounding errors could be significant. Entering the data as 0 to 1 femtoseconds may be more appropriate.

Types of Regression Analysis

Regression acts to minimize the sum of the squares deviations of the experimental values from values calculated using some theoretical equation. However, the mechanics of the calculations are different for linear and non-linear equations.

Linear Regression

Linear regression is a technique that should be familiar to most scientists. It is the mathematical equivalent of using a ruler to draw the “best” line through a series of data points to obtain values for the slope and intercept of this line. The calculations needed to perform linear regression are relatively simple, and are even found on many scientific calculators.

Polynomial Regression

Polynomial regression is an extension of linear regression, where the equation is

$$y = a + bx + cx^2 + dx^3 + ex^4 + \dots$$

Non-linear Regression

Non-linear regression allows data to be fitted in the more general case to any equation where the y data value can be described as a function of the x data and a series of parameters, i.e.

$$y = f(x, p_1, p_2, p_3, \dots)$$

As a result, it is applicable to many more situations than linear regression, and is more generally useful for analyzing experimental data. The only drawback to non-linear regression is that the calculations are too complex to perform using a calculator, and require a computer program such as GraFit.

Additional Considerations for Non-linear Regression

The more complex nature of the non-linear regression calculations does require some further considerations compared with linear regression.

Equation to Use

Non-linear regression can employ a vast range of different equations. It is therefore necessary to select an equation that is appropriate to the particular experimental situation before fitting the data.

Initial Estimates

The mathematical method used to perform non-linear regression acts to take starting estimates for the parameters in the equation, and optimize them during the calculations. As described in Chapter 9, GraFit provides a mechanism for providing these initial estimates

automatically for some equations, but in principle it is necessary to give the calculations some rough estimates prior to data fitting.

Convergence Criteria

Unlike linear regression, non-linear regression does not provide an exact solution, but, by an iterative process, calculates successively better parameter values until a suitable tolerance is reached. The criteria for defining this tolerance are determined by the program, and can be specified when using GraFit.

Weighting Scheme

Both linear and non-linear regression analyses incorporate assumptions about the accuracy of the data points. Normally, they assume that all data points are equally accurate, and so carry equal weighting in the analysis. It is also possible to use other weighting schemes that give different weighting to one or more data values. A full description of weighting is found in Chapter 20.

Linear Regression of Transformed Data


A familiar method of analyzing non-linear data graphically is to rearrange the data to a linear form, find the “best line” through the points, and so derive the parameter values. This is illustrated in the following example. Enzyme kinetic data are described by a hyperbolic equation of the form:

$$y = \frac{Ax}{B + x}$$


This can be linearized by plotting $1/y$ versus $1/x$, as

$$\frac{1}{y} = \frac{1}{x} \cdot \frac{B}{A} + \frac{1}{A}$$

The ease of linear regression analysis makes it tempting to use in order to fit data that have been rearranged in this way. However, the results that are obtained are not statistically optimal. This is because the regression calculations make assumptions about the error distribution that are not usually valid after transformation. The distortion of the experimental errors by the above transformation is easily demonstrated in the following graphs.

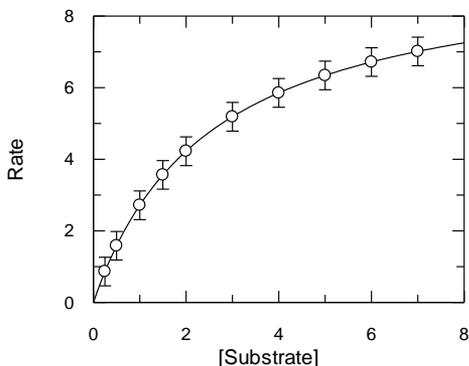

Enzyme kinetics plot. The error is present only in the Rate (y axis) measurement. Errors are simple, i.e. all data points have the same error.

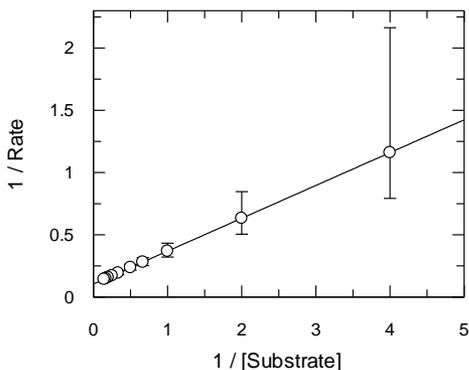

Double reciprocal (Lineweaver Burk) plot of the data above. Note the large distortion of the error bars resulting from this plot.

Here, the error involved in the high data points of the transformed plot is much greater than in the low data points — this is simply a consequence of the reciprocal transformation. Linear regression (unless suitable compensation is made for the distorted errors) is therefore not an appropriate way to analyze the data, and the original data should be fitted to the appropriate non-linear equation.

Information Obtained from Data Fitting

Data fitting provides the following information regarding the experimental data that are analyzed.

1. The values of the parameters in the equation used that best describe the data. These parameters may include (for example) rate constants, inhibition constants etc. whose values need to be determined.
2. The accuracy of these calculated parameters.

From this information it is possible to judge how well the data fit the experimental model (equation) that was used. In most cases, however, it is the numerical values of the parameters that are important, and whose values are required from the analysis.

Multiple Regression

The most common form of non-linear equation that is analyzed by GraFit is of the form

$$y = f(x)$$

where the right hand side of this expression may contain one or several parameters. Such an equation has one dependent (y) and one independent variable (x). In addition, GraFit can be used to analyze equations with more than one independent variable, i.e.

$$y = f(x_1, x_2, x_3, \dots)$$

In this case the analysis is termed multiple non-linear regression. GraFit allows unlimited independent (x) variables, as well as unlimited numbers of parameters. Note, however, in all cases it is assumed that only the y axis data contain errors.

Note Multiple regression is sometimes also referred to as "global analysis."

Uses for Multiple Regression

Multiple regression is essential when the data to be examined are inherently multidimensional. Examples include peak-shape analysis of 2D spectra. However, there are several situations where multidimensional data are collected, but are traditionally analyzed by performing a series of normal two-dimensional (x, y) curve fitting.

Example

The progress of an enzyme-catalyzed reaction at a fixed concentration of substrate is given by the equation:

$$A_t = A_\infty(1 - e^{-kt})$$

Here, A_t is the extent of reaction at time t , A_∞ the maximum extent of reaction, and k the rate constant at this concentration of substrate. The value of k from this equation is related to the kinetic parameters by the equation:

$$k = \frac{V_{max}[S]}{K_m + [S]}$$

In this equation, V_{max} and K_m are the parameters that we need to obtain from data fitting. The usual way of analyzing these data is to collect a series of progress curves (A as a function of t) at varying concentrations of substrate. For each concentration of substrate the data are fitted to the first equation, and a value for k obtained. The resulting k values are then fitted to the second equation to find V_{max} and K_m . A better solution is to recast the equation as a three-dimensional one, where the progress of the reaction varies with t and $[S]$, which are the two independent variables. The equation is obtained by combining the previous ones:

$$A_t = A_\infty \left(1 - \exp \left(- \frac{V_{max}[S]}{K_m + [S]} \cdot t \right) \right)$$

This allows V_{max} and K_m to be determined directly from the raw data, which has two x variables (t and $[S]$).

A sample set of data that has been fitted in this way can be found in the Examples data subdirectory, and is named "Multiple rate."

Advantages of Multiple Regression

Use of multiple regression simplifies the analysis, as the curve fitting need only be carried out once. In the above example, the value of A_∞ is fixed for the entire experiment and only found once for the multidimensional equation. However, if the data are analyzed by several applications of the first equation, the A_∞ value would be determined several times.

Disadvantages of Multiple Regression

The disadvantages of multiple regression are as follows.

- Conceptual. It is often difficult to visualize the problem in more than two dimensions.
- Finding initial estimates. This is sometimes difficult for the same reasons.
- Difficulty of display. The current version of GraFit does not support three-dimensional plotting; the display of four-dimensional data is even more problematic!

Defining Equations

Definition of equations having multiple independent variables is achieved in the same way as the definition of normal two-dimensional x/y equations, but the equations have more than one x variable.

Fitting Data by Multiple Non-linear Regression

For a normal data set that contains x/y data, these data must be present in the data table with the x data and the y data in separate columns. When multiple independent variables are involved, there are special considerations as to how the data should be arranged.

Case 1: The equation has only two independent variables

In many circumstances there will be just two independent variables in an equation. In these circumstances, it is common to collect one set of $\{x_1, y\}$ data at a particular x_2 value, then the next set of $\{x_1, y\}$ data at a different x_2 value etc. In the case of the example described above, this would entail collecting A versus t at one particular $[S]$, then another set of A versus t at the next $[S]$ and so on.

When this arrangement of data is present, it is convenient to enter the $\{x_1, y\}$ data into the data table, and to associate the x_2 value with the y data column. For example, the data can be entered as shown below.

	1	2	3	4
Type		x2: [1.0000]	x2: [2.0000]	x2: [3.0000]
Column Name	x	y 1	y 2	y 3
1	1.0000	0.0480	0.0760	0.1130
2	2.0000	0.0950	0.1480	0.2060
3	3.0000	0.1210	0.2110	0.2670
4	4.0000	0.1670	0.2640	0.3220
5	5.0000	0.1990	0.2940	0.3630
6	6.0000	0.2090	0.3330	0.3980
7	7.0000	0.2480	0.3610	0.4260
8	8.0000	0.2710	0.3740	0.4430
9	9.0000	0.2860	0.3980	0.4540
10	10.0000	0.3030	0.4270	0.4570

Note that the y data each have an associated x_2 value. For this particular data set, all the y data are associated with the same x_1 data, which corresponds to the Non-linear Fit data layout type:

$x[n]$ from y header x y_1 y_2 ... y_n

(See Chapter 6 for more details about the Non-linear Fit window.)

It is also possible to have a separate set of x_1 data for each y data, corresponding to the Non-linear Fit data layout type:

$x[n]$ from y header x_1 y_1 x_2 y_2 ... x_n y_n

Case 2: The equation has several independent variables

When an equation has several independent variables it is necessary to have the x_1 , x_2 , x_n , y data in separate columns. This corresponds to the Non-linear Fit data layout type:

$x[n]$ from table x_1 x_2 x_n y


In order to distinguish the various data types, it is essential that all the x data columns are set to “ $x(n)$ Data” with the “ n ” value representing the order of the x variable found in the equation definition. See Chapter 5 for further details on column types.

Data Fitting

Apart from the considerations above, fitting data that have multiple independent variables is no different from fitting standard equations. A sample set of data that has been fitted in this way can be found in the Examples data subdirectory, and is named “Multiple rate.”

Plotting Results from Multiple Regression

The data that form the basis for multiple regression have three or more dimensions. To represent such data would, strictly speaking, require a multidimensional plot, which is not available. However, the usual representation for such data is to plot multiple curves, representing the raw data as collected, on the same graph. This is particularly simple for the situation noted as "Case 1" above. Choosing to create a graph to display the non-linear fit will automatically generate a suitable plot. In the case of the sample data set, the fitted curve is shown below. The file that generated this curve can be found in the "Samples" sub-directory of your My GraFit Data directory, named "Multiple rate".

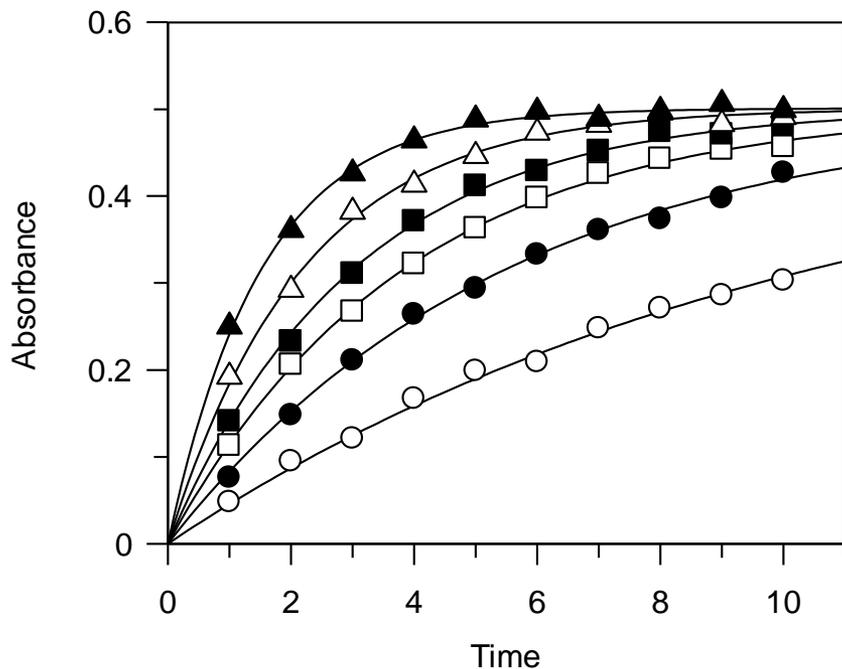

Theory

To make best use of GraFit, it is helpful to know some of the theory behind regression analysis. Chapter 19 provided a general overview; this section describes more of the theoretical background.

Regression Analysis

Regression analysis is a technique used to calculate the “best fit” line or curve through a data set by minimizing the deviation of the data from the curve. Mathematically, it finds values for various unknown parameters in an equation (for example, the slope and intercept of a straight line) that produce the smallest possible sum of the squares deviations of the experimental values from the calculated ones. This quantity that is minimized is denoted by the symbol χ^2 (chi-squared). The value of χ^2 includes the variance of the data point (the variance of data point i is σ_i), so that data points that have greater errors contribute less to the analysis. It is calculated as:

$$\chi^2 \equiv \sum \left(\frac{\Delta y_i}{\sigma_i} \right)^2$$

Δy_i is the difference between the experimental and calculated y data values (the subscript “ i ” is used to indicate the “ i th” data value).

A further useful statistic is χ^2_ν , or the reduced chi-squared. This is defined as below.

$$\chi^2_\nu = \frac{\chi^2}{\nu}$$

Here, ν is the number of degrees of freedom, which is given by $\nu = N - n$; N is the number of data points and n is the number of variables in the equation.

Regression analysis is the mathematical equivalent of finding the best line by eye—it has the advantage, however, of being non-subjective. It is also possible to fit the “best curve” through the data by non-linear regression. This has no manual equivalent (because humans are very poor at judging accurate curvature), but is equally simple mathematically using computational methods.


Assumptions in Regression Analysis

Regression analysis is probably the most powerful method to analyze data. However, such analysis involves several assumptions about the function chosen and the error structure of the data. These assumptions are listed below.

The Correct Equation is used

This assumption is often overlooked, but should be borne in mind when using any data analysis technique. It is important to be aware that deviations from ideality are often caused by use of an inappropriate equation. Graphical methods provide a useful means to see whether the chosen model equation is appropriate. This is because use of an inappropriate equation often results in systematic deviations from the curve, which is easy to see in a graph.

Residual plots of Δy versus x are useful to assess whether the equation is appropriate. Systematic deviations in the residuals can mean that the equation used was inappropriate to describe the data. The differences between experimental and calculated values can be calculated using options available in the data table (see Chapter 5). The data shown below have been fitted to a single exponential decay, and to a double exponential decay, and the residuals from each have been plotted. Clearly, the residuals from the double exponential fit are smaller and more randomly distributed, and show that this fit is superior.

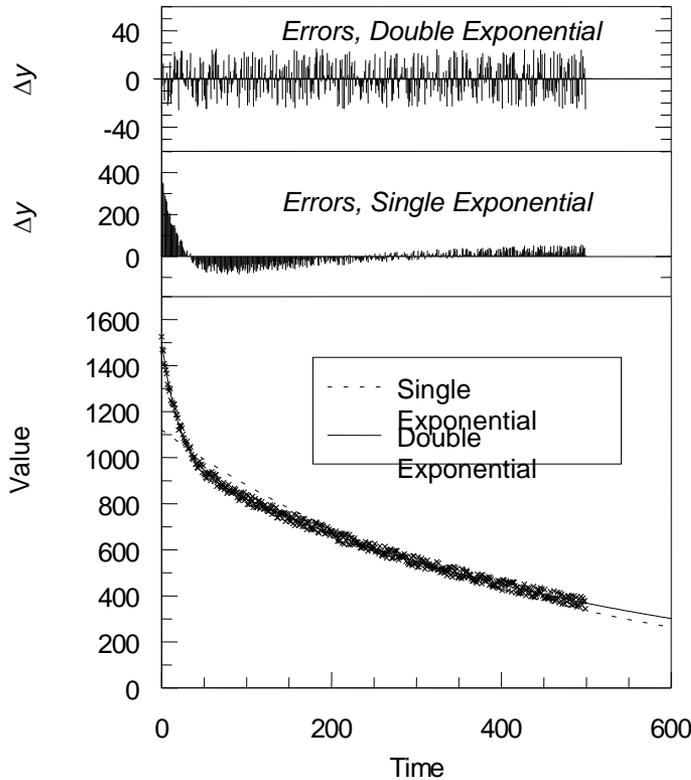

Only the Y Data Contain Errors

Only the dependent variable, which by convention is always the y variable (y_i), is subject to error; the other (x_i) is known exactly. This requirement should be satisfied by suitable experimental design — the x data values often comprise time or concentration values, which should be known with precision.

The Errors are Normally Distributed

The errors in the y data values should follow a normal distribution function, i.e. if sufficient replicates were taken they would have a Gaussian spread about the true y value. Two common factors that invalidate this assumption are the presence of (i) systematic errors, and (ii) outliers (data points that are wildly inaccurate, usually due to operator or machine failure). Systematic errors invalidate any data analysis, and should be removed by correct experimental design. Outliers can be eliminated by use of “robust” techniques, as will be described later in this chapter.

The Correct Weighting is used

An explanation of what “weighting” of data in regression analysis is, why it is incorporated into all such analysis (including the subjective regression performed by drawing the “best line” using a ruler), and how the appropriate weighting can be determined are described later.

Linear Regression

The simplest case for regression analysis is when the equation is linear. For the linear equation:

$$y = a + bx$$

χ^2 is defined as:

$$\chi^2 \equiv \sum \left(\frac{\Delta y_i}{\sigma_i} \right)^2$$

$$\chi^2 = \sum \left(\frac{1}{\sigma_i^2} (y_i - a - bx_i)^2 \right)$$

In this equation, σ_i^2 is the standard deviation of the data point i .

The best fit to this equation will yield parameter values for a and b that produce the minimum possible value of χ^2 . This is equivalent to saying that the first derivative of χ^2 with respect to a and b is equal to zero:

$$\begin{aligned} \frac{\partial}{\partial a} \chi^2 &= \frac{\partial}{\partial a} \sum \left(\frac{1}{\sigma_i^2} (y_i - a - bx_i)^2 \right) \\ &= -2 \sum \left(\frac{1}{\sigma_i^2} (y_i - a - bx_i) \right) = 0 \end{aligned}$$

$$\begin{aligned} \frac{\partial}{\partial b} \chi^2 &= \frac{\partial}{\partial b} \sum \left(\frac{1}{\sigma_i^2} (y_i - a - bx_i)^2 \right) \\ &= -2 \sum \left(\frac{x_i}{\sigma_i^2} (y_i - a - bx_i) \right) = 0 \end{aligned}$$

If we rearrange these equations it is possible to solve for a and b exactly to give values for these parameters that minimize χ^2 . This is shown in the equations below.

$$a = \frac{1}{\Delta} \begin{vmatrix} \sum \frac{y_i}{\sigma_i^2} & \sum \frac{x_i}{\sigma_i^2} \\ \sum \frac{x_i y_i}{\sigma_i^2} & \sum \frac{x_i}{\sigma_i^2} \end{vmatrix}$$

$$b = \frac{1}{\Delta} \begin{vmatrix} \sum \frac{1}{\sigma_i^2} & \sum \frac{y_i}{\sigma_i^2} \\ \sum \frac{x_i}{\sigma_i^2} & \sum \frac{x_i y_i}{\sigma_i^2} \end{vmatrix}$$

$$\Delta = \begin{vmatrix} \sum \frac{1}{\sigma_i^2} & \sum \frac{x_i}{\sigma_i^2} \\ \sum \frac{x_i}{\sigma_i^2} & \sum \frac{x_i y_i}{\sigma_i^2} \end{vmatrix}$$

Weighting

This analysis presupposes that we know the standard deviation, σ_i^2 , for each data point. In general, we will not have such information available. However, only relative uncertainties need be known, and if we assume that all the data points have the same inherent errors, the equations simplify to give:

$$a = \frac{1}{\Delta} \begin{vmatrix} \sum y_i & \sum x_i \\ \sum x_i y_i & \sum x_i^2 \end{vmatrix}$$

$$b = \frac{1}{\Delta} \begin{vmatrix} N & \sum y_i \\ \sum x_i & \sum x_i y_i \end{vmatrix}$$

$$\Delta = \begin{vmatrix} \sum y_i & \sum x_i \\ \sum x_i & \sum x_i y_i \end{vmatrix}$$

This assumption that the errors are constant is often made, for example in the simple linear regression calculations performed by pocket calculators. Although often true, this assumption is not always valid. In many experiments the uncertainty varies with the magnitude of the y value, such that the percentage error is constant. To account for the different types of error distribution we introduce a weighting factor, which is related to $1/\sigma_i^2$ and is used to multiply all elements of the above matrices.

Weighting Types

The following types of error distribution are easily accounted for, and can be selected for use in the GraFit program within the Linear Fit or Non-linear Fit item windows.

Constant error

The standard error, σ_i^2 , is the same for each data point and therefore a weighting factor of 1 is used (“Simple weighting”). This type of error is often given by instrumental variations (for example, a gauge may be read to an accuracy of 0.1 units, etc.).

Statistical error

For statistical variation, where σ_i^2 varies in proportion to y_i , the weighting factor is $1/y_i$. This type of error distribution is uncommon, but is found in systems where the main errors arise from sampling, for example counting low numbers of radioactive disintegrations.

Proportional error

For many experiments the percentage error is constant. In this case, σ is proportional to y , leading to a weighting factor of $1/y_i^2$. A typical example of such errors might be when, for instance, a measuring instrument is accurate to 2% of the delivered volume, and the so error is proportional to the recorded value.

Individually determined error

All the above cases make assumptions about the error of each measurement. However, if sufficient replicates of each data point are available, or a complex error function is applicable, it is possible to calculate σ_i^2 for each point and use these values explicitly. This method can also be used if you know that certain values are less accurate than others (for example if they were measured in a different way) as again only relative σ_i^2 values need be given.

In practice, the use of simple weighting or proportional weighting is most often appropriate. If the errors are small the weighting function chosen has only a small effect on the calculated values. We will see later how it is possible to estimate which type of error function best describes a set of data. If any data values are close to 0, it is not appropriate to use proportional weighting, and simple weighting should be selected.

Non-linear Regression

For an equation that is not a linear function the determination of best-fit parameters becomes slightly more difficult, although the principles discussed above still apply. The difficulty arises because it is not possible to derive an exact analytical expression for minimizing χ^2 . It is possible, however, to take initial estimates for the best-fit parameter values and to improve on them by a series of iterative calculations until values are obtained that approach the minimum χ^2 value.

Apart from this mathematical inconvenience, non-linear regression functions in the same way as linear regression. We still need to know the weighting that has to be applied to the data, and the calculations still aim to minimize χ^2 . The advantage is that non-linear regression may be used to analyze any function that may be written $y = f(x)$, i.e. the y data value is dependent upon some arbitrary single-valued function of the x data value.

GraFit performs non-linear regression using the method of Marquart⁸⁵ using a numerical second order method to calculate partial differentials. In this way it is possible to incorporate

⁸⁵ Marquart, D.W. (1963) *J. Soc. Ind. Appl. Math.*, **11**, 431-41.

The Marquart method is sometimes also referred to as the Marquart-Levenburg method. It is perhaps best explained in reference 2.

equations of arbitrary form. For an excellent explanation of regression analysis the reader is referred to the text by Bevington⁸⁶.

Standard Errors of the Calculated Parameters

As well as giving the best-fit values for the parameters, regression analysis can provide estimates of the standard error of each parameter. GraFit calculates these errors by the matrix inversion method. These estimated standard errors are useful to show how accurate the calculated parameters are. As a rough guide, the standard error should not be greater than about 10% of the parameter value if you are to have confidence in the calculated parameters. High standard errors will be obtained if the data are of poor quality, or if the data do not allow accurate estimation of one or more of the parameters.

Robust Weighting

One of the assumptions of regression analysis is that all the data points have a normal error distribution. However, a common feature of actual experimental data is the presence of one or two observations that are well outside the error range of the other observations. The stray data points are usually due to mistakes by the experimentalist (incorrect addition of reagents, inaccurate keying in of data, etc.) and are termed outliers. If these points are included in the analysis they will disproportionately affect the results, since the square of the residuals for these points will dominate the fit. It is possible to edit the data manually to remove such points, but it is preferable to eliminate them in a less subjective manner. Robust methods achieve this by applying an additional weighting factor varying, for example, from 1.0 for points that have low residuals, to 0 for extreme outliers — effectively eliminating such points.

GraFit employs optional robust weighting using the algorithm of Mosteller & Tukey⁸⁷ as implemented by Duggleby⁸⁸. An additional ‘bisquare’ weight, b , is incorporated into the analysis, and is calculated as:

$$b_i = \begin{cases} (1 - u_i^2)^2 & \text{if } |u_i| \leq 1 \\ 0 & \text{if } |u_i| > 1 \end{cases}$$

Here, with z being the residual weighted by the *a priori* weight and N the number of observations:

$$u_i = \frac{z_i}{6 \sum |z_i| / N}$$


Including robust weighting will result in improved estimates for the parameters in the presence of outliers. This is, however, gained at the expense of more calculation, which will slow down the fitting procedure.

⁸⁶ Bevington, P.R. (1969) *Data Reduction and Error Analysis for the Physical Sciences*. McGraw-Hill Book Company.

⁸⁷ Mosteller, F. & Tukey, J.W. (1977) *Data Analysis and Regression*. Addison-Wesley, Reading MA, USA pp.353-65.

⁸⁸ Duggleby, R.G. (1981) *Analytical Biochemistry* **110**, 9-18.



The effect of using robust weighting is seen in the graph below:

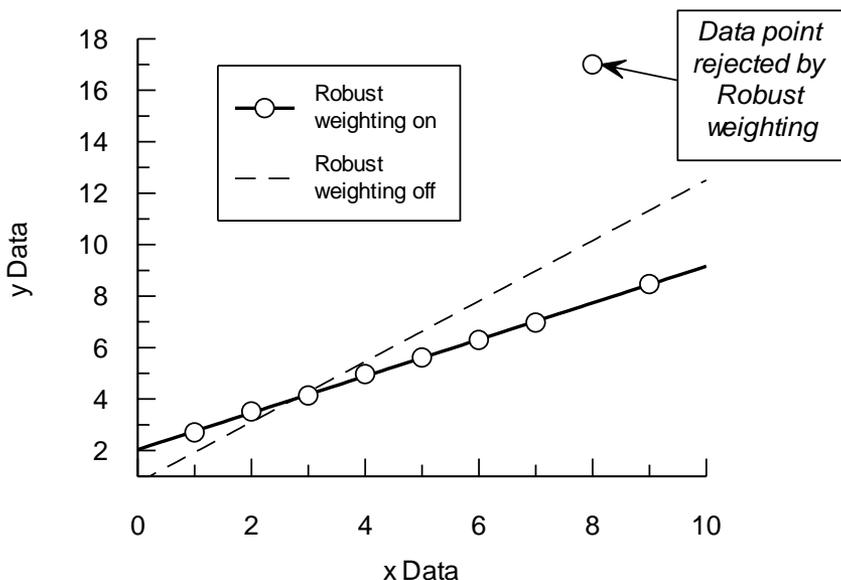

The data, which contain a single outlier, have been fitted to a straight line with and without robust weighting enabled. When robust weighting is omitted, the presence of this outlier markedly affects the calculated results. With robust weighting selected, however, the rogue data point is rejected, and the calculated line now correctly follows the main trend.

Determining the Most Appropriate Weighting

For most experimental data the error is commonly either independent of the value of the observation (simple or constant errors) or is proportional to the y data value (proportional errors). Simple errors may be caused by instrumental uncertainties (for instance, a machine reading may be only accurate to \pm some intrinsic error, independent of the magnitude of the reading). Providing sufficient data points have been sampled, it is possible to determine the correct weighting by examination of the residual error (the experimental minus the calculated value). A plot of the residuals versus the observed value (Δy versus y) is diagnostic of the type of error involved. For simple errors by definition, the residual is independent of y , as is seen in the plot (a). If the errors are of a proportional nature the residual plot will be as plot (c). A plot of the relative residual ($\Delta y/y$ versus y) results in the profiles (b) and (d) for simple and proportional errors respectively, and confirms the error distribution.

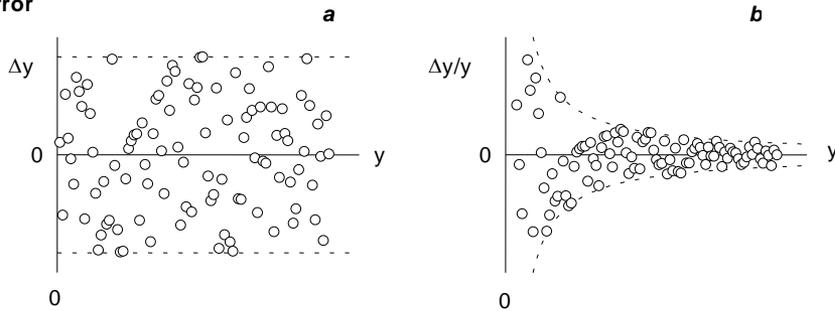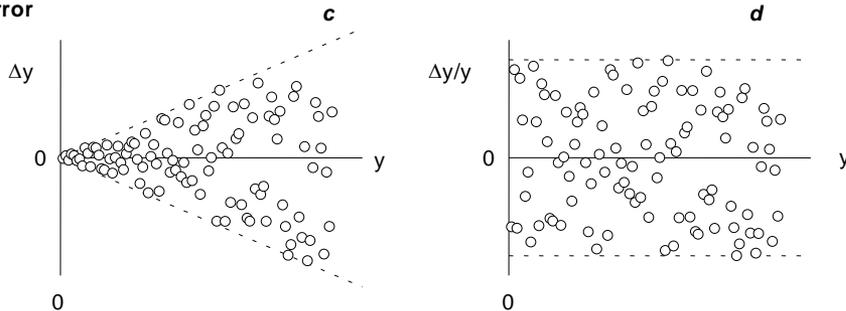Constant Error**Proportional Error****Plot of Residual Errors**

Above is shown the type of residuals plots obtained when the error is constant or proportional in nature. It must be stressed that visualization of clear-cut plots of this kind depends upon having many data values to define the curves. The dotted lines on the plots show the envelope within which all points are expected to be found. Residuals can be calculated using options available in the data table (see Chapter 5).

Testing Goodness of Fit — *F* Test

There are situations where it is necessary to know whether a set of data is fitted better by one particular equation or another. The *F* test is a statistical way of judging the goodness of fit by two different equations. This test compares the χ^2 values obtained from the two fits and the degrees of freedom in each, and calculates the probability that the fits are equally good. A low probability value indicates that one of the two equations (that giving the lower χ^2) fits the data significantly better. To be significant the probability should be lower than 0.1, and preferably lower than 0.05. To perform the *F* test, use the column type options available in the data table (see Chapter 5).

Nested Models

In the case where data have been fitted to two related models, where one is an extension of the other by virtue of having one or more additional parameters, the F statistic is calculated from the equation below.

$$F = \frac{(\chi_1^2 - \chi_2^2)/(v_1 - v_2)}{\chi_2^2/v_2}$$

Comparison of Equations with Equal Numbers of Parameters

Here the F statistic is calculated from the ratio of the χ_v^2 values:

$$F = \frac{\chi_{v(1)}^2}{\chi_{v(2)}^2}$$

If the fits are equally appropriate, this ratio will be 1.0; the F test assesses the significance of the observed differences.

Mean and Standard Deviation

For a full description of means and standard deviations, see any text on statistics (for example, that by Bevington⁸⁶). A brief description follows.

The mean value of a set of data points is defined as:

$$\mu \approx \bar{x} = \frac{1}{N} \sum x_i$$

Here, μ is the mean of the parent distribution, and \bar{x} is the mean of our experimental sample. The uncertainty, or standard deviation, of the data points, σ , is estimated from the data as follows:

$$\sigma \approx s = \sqrt{\frac{1}{N-1} \sum (x_i - \bar{x})^2}$$

The uncertainty of the mean, σ_μ , is given by the equation:

$$\sigma_\mu = \frac{\sigma}{\sqrt{N}} \approx \frac{s}{\sqrt{N}}$$

All the above equations assume that the uncertainties of all the data points are equal. If they are not, for reasons described earlier, we must use equations that incorporate weighting. For example, calculation of the mean becomes:

$$\mu \approx \bar{x} = \frac{\sum(x_i/\sigma_i^2)}{\sum(1/\sigma_i^2)}$$

Worked Data Fitting Examples


This chapter provides some further examples of data fitting to help acquaint you with the program. The examples are presented so that they can be worked through in a step-by-step manner. The data files that are used can be found in the “Samples” sub-directory of your My GraFit Data directory. Don’t worry if the examples do not seem relevant to the type of data that you need to analyze — the same principles are involved in all data fitting.

This chapter assumes that you are familiar with the basic functioning of the program, and if not you may like to consult the earlier chapters.

Using One of the Equations Provided

We will fit the data obtained from an NMR titration curve to obtain the pK_a of an ionizable group.

Background The position of a peak in the NMR spectrum (the chemical shift position) resulting from a magnetic nucleus that is part of a titratable group depends upon the equilibrium between the two ionized forms. For example, in the following equilibrium:

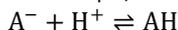

Nuclei close to the ionizable group, A, will have chemical shift values that depend on the ratio of A^- and AH. By definition, the acid dissociation constant, K_a is given by

$$K_a = \frac{[A^-][H^+]}{[AH]}$$

It is usually more convenient to use pK_a values, where

$$pK_a = -\log_{10} K_a$$

The variation of the ratio of A^- and AH is given by the Henderson-Hasselbalch equation:

$$pH = pK_a + \log_{10} \frac{[A^-]}{[AH]}$$

The chemical shift observed for the nucleus varies between that measured at the extreme acid range, δ_{AH} , and that at the extreme alkaline range, δ_{A^-} , depending on the ratio of $[A^-]/[AH]$.

The data obtained are chemical shift values (in units of ppm, or parts per million) at various pH readings. These values need to be fitted to a single ionization curve where the start and end chemical shift values and the pK_a are the unknown parameters.

Fitting the Data

The example data are in the file “pH Titration”, which can be opened using the **File ► Open** command. Two columns are present in the data table: pH and Chemical Shift. No other workspace items are present.

Load the Equation


The first action to take is to load the correct equation. Use the **Add Stock Item** command⁸⁹ on the **Home** tab; select the “pH” tab and choose the “pKa, Single” equation definition item.

Add a Non-linear Fit Item to the Workspace

Use the **Add ► Non-linear Fit** command to add a non-linear fit item to the GraFit workspace.

Fill out the Non-linear Fit Item Settings

Complete the Settings page as shown below:

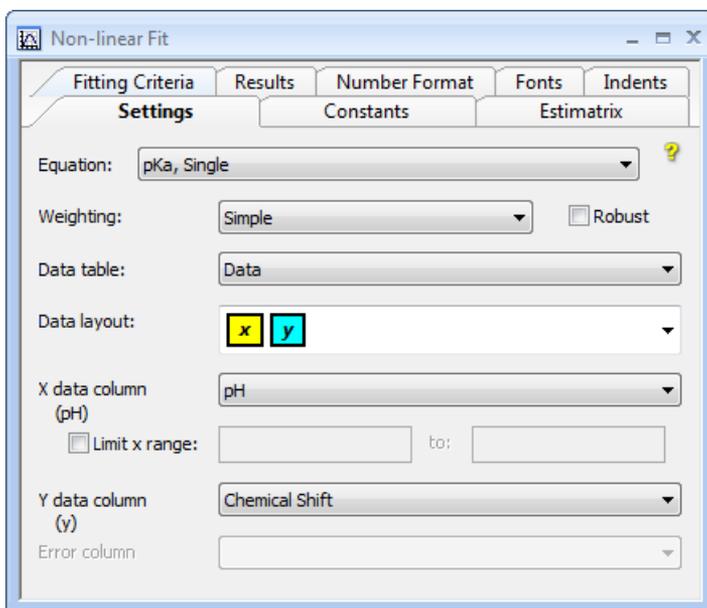

Set up Initial Estimates

The “pKa, Single” equation is one for which initial estimates for the parameter values are difficult to provide automatically. Therefore it is necessary to enter these initial estimates explicitly. For equations that require initial estimates it is essential to set these from the Estimatrix page. Click on the Estimatrix tab, then on the **Interactive Estimator** button.

⁸⁹ The **File ► Add Stock Item** command in legacy mode

The screenshot shows the 'Non-linear Fit' software window with the 'Interactive Estimator' dialog box open. The dialog box contains a plot of Chemical Shift versus pH, a table of estimated values, and various control buttons.

Graph Data:

pH	Chemical Shift
4.0	9.5
4.5	9.5
5.0	9.5
5.5	9.5
6.0	9.5
6.5	9.4
7.0	9.2
7.5	8.8
8.0	8.4
8.5	8.1
9.0	8.0

Estimated Values Table:

Parameter	Estimated value
pKa	7
Limit1	9.5
Limit2	8


Enter values into the **Estimated value** grid; as you move the entry position, the new values will be used to create a theoretical curve (red) on the plot. As described in Chapter 20, the starting estimates for non-linear regression should be reasonably close to the true values, but it is not necessary to be too exact. The program then iterates to find successively better values until a satisfactory fit is found. When you are satisfied with your estimated values, click **OK** to return these as the starting estimates on the Estimatrix page.

Fit the Data

To fit the data, choose the **Fit Data** command⁹⁰ on the **Home** tab. The results are then calculated, and can be displayed using the **File ► Print Preview** command, or viewed from the Results page.

Graph the Fit

To plot a graph of the fitted data, add a new graph to the workspace using the **Add ► Graph** command. Elect to create an X/Y Scatter Graph, then use the newly-created non-linear fit item as the source of the data:

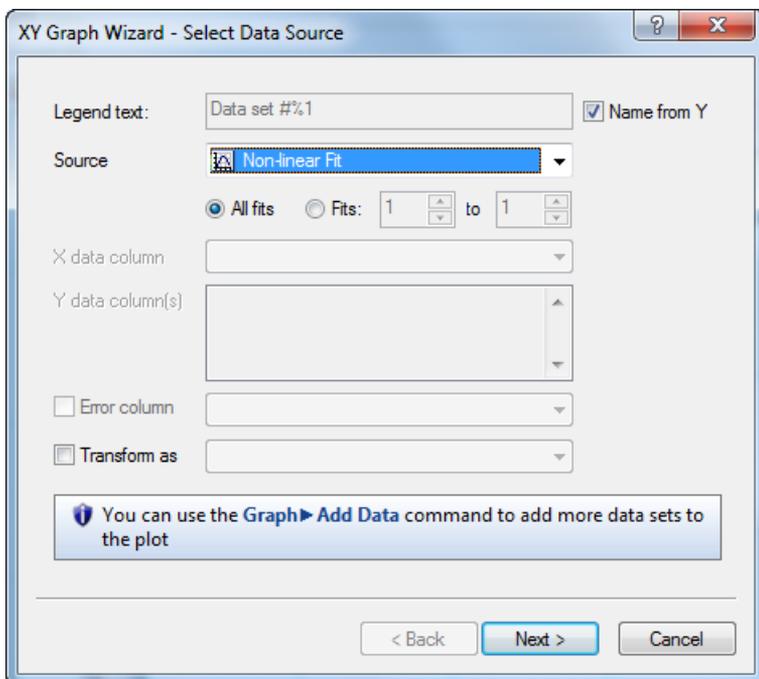


Show the Results on the Graph


It is often a good idea to display the data fitting results on the graph. To do this, select the **Results Listing** command⁹¹ (on the **Graph** tab), and elect to display the results of your non-linear fit item.

⁹⁰ The **Calculate ► Fit Data** command in legacy mode

⁹¹ The **Draw ► Results Listing** command in legacy mode

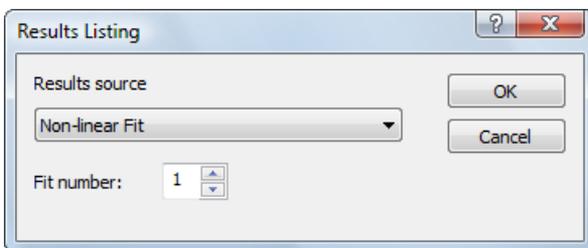


The information displayed in the results listing includes the parameter values for the fitted data, but can also show various other pieces of information about the fit. To customize the display, double-click on the results listing once it has been generated. The following dialog box is displayed:

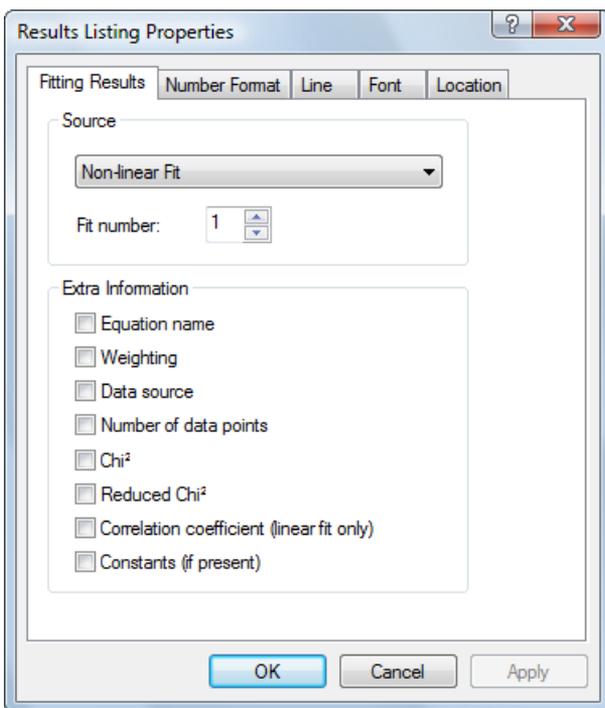

The Final Result

For these data, the final result is shown below. Here, the results listing has been added to the graph, and in addition the pK_a value has been indicated visually using the “pKa, Single, Show pKa” transformation. It is clear in this instance that the fitted curve describes the data well.

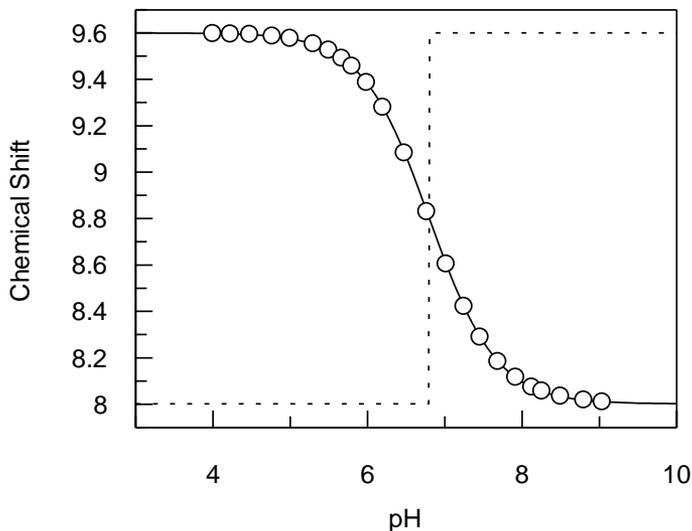

Parameter	Value	Std. Error
pKa	6.8001	0.0006
Limit1	9.6005	0.0003
Limit2	8.0015	0.0003

(To use this transformation, load it using the **Add Stock Item** command⁸⁹, then click on the graph and select the **Add Data** command⁹² from the **Graph** tab. The data source is the Non-linear fit item, and the **Transform As** box should be checked to allow the “pKa, Single, Show pKa” transformation to be selected.)

Checking the Weighting

It is often useful when performing an experiment for the first time to estimate the type of weighting that should be used. It must be stressed that for these checks to be meaningful a large number of data points should be available, and the equation describing the data must be well established.


As was described in Chapter 20, plots of residual error (Δy) versus y or $\Delta y/y$ versus y give characteristic patterns for the common simple and proportional error distributions. We will produce these plots for the data that have been fitted above.

Drawing Residual Plots

We will plot the residual plot Δy versus y . To make this plot we need a column of data containing the Δy values. This can be created after data fitting by selecting the data table,

⁹² The **Graph** ► **Add Data** command in legacy mode


then setting the **Column Type (Table tab⁹³)** of one of the columns to **Calculated Values ▶ Difference**

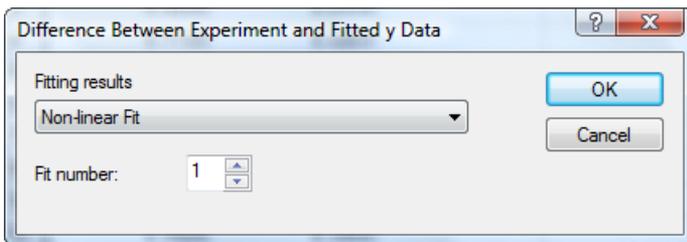

Select the Non-linear fit that was generated above, and click **OK**.


It is advisable to set the column name so that you can identify it; in the current example it has been set to “Differences”

To create a graph of the residuals on a new graph page, select the **Add ▶ Graph** command, and create an X/Y Scatter Graph. Create this graph using the data table as the source, the Chemical Shift for the x data and the Differences column for the y data:

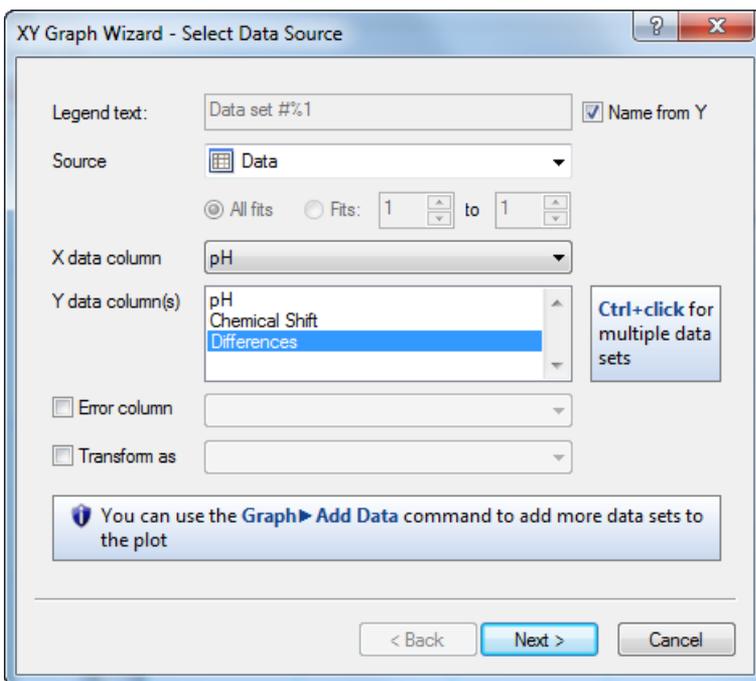


Select to join the points with a straight line, as below.

⁹³ **Table ▶ Column Type** in legacy mode

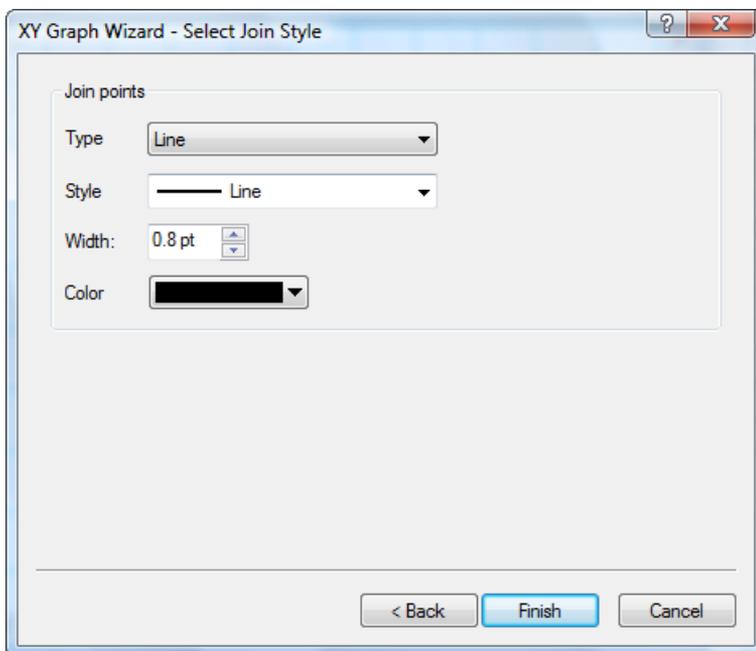


The resulting graph is shown below.

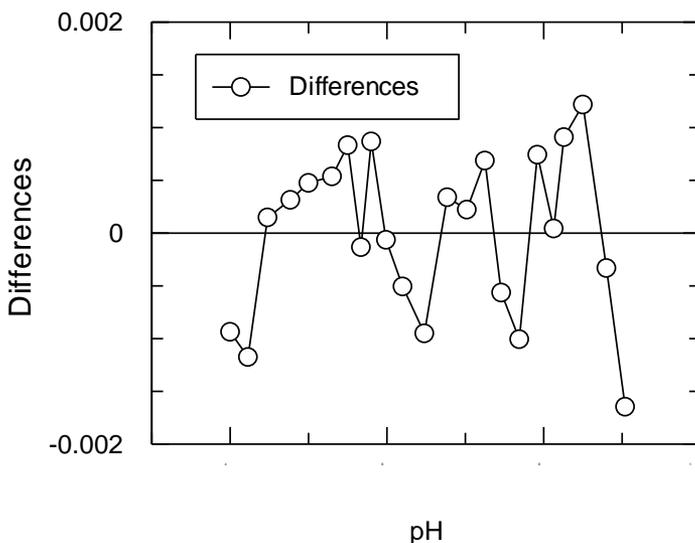

As is clear from the shape of this plot (see the discussions in Chapter 20), it was probably correct to use simple weighting with these data.

Δy versus x plots are easily made by plotting pH against the Differences. To make $\Delta y/y$ versus y plots, it is possible to create a new column to hold the result of dividing the Differences column by the original y data column (Chemical Shift). This can be created after data fitting by

setting the **Column Type** of one of the next available columns to **Column Combinations** ▶ **A/B**.

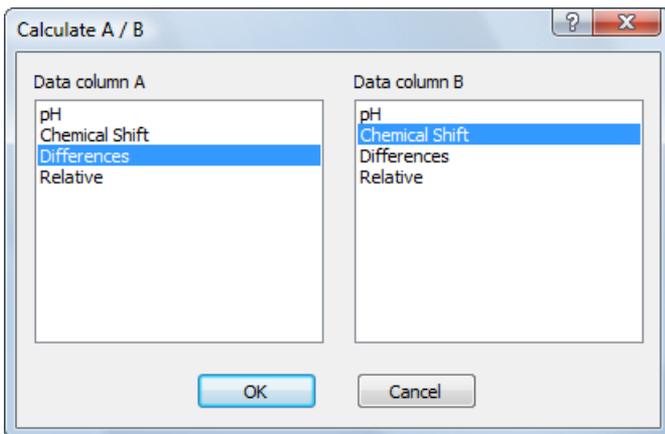


Using Robust Weighting

Robust weighting is additional weighting that is applied to a data set. It acts to reject data points that are very much in error, and so removes the occasional rogue data point or “outlier” that may be present in a data set. A complete description of the background to robust weighting can be found in Chapter 20.

To show the detrimental effects of outliers on curve fitting, and see how robust weighting eliminates these effects, an analysis has been created which can be found in the file “Robust weighting”. The data being have a linear relationship, but are fitted using a linear equation defined as a non-linear equation in order to make use of the robust weighting option that is available for non-linear fitting.

Note It is not possible to use a linear fitting item when using robust weighting.

The graph obtained by plotting the data with the best-fit lines obtained with and without robust weighting is shown below.

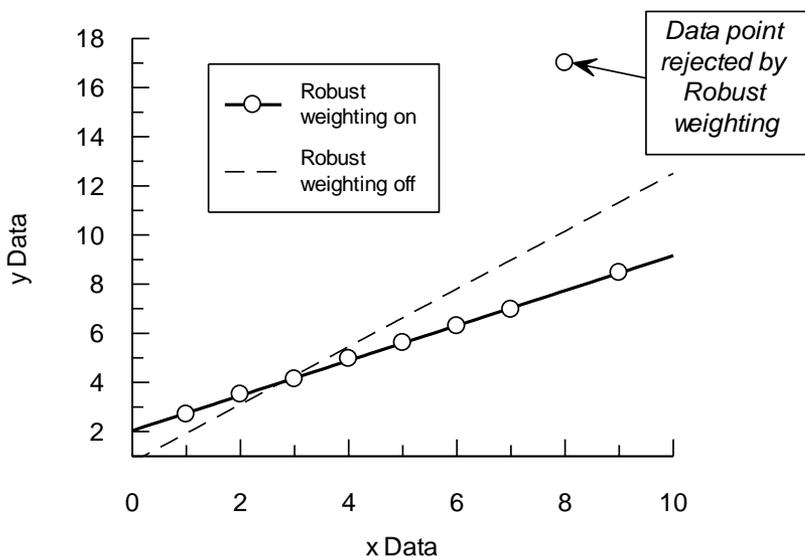

As is apparent, without robust weighting selected the presence of the outlier biases the fit considerably. This is because the large sum-of-squares residual from this point has a disproportionate effect on the regression analysis. However, using robust weighting eliminates the outlier, and the regression line is drawn through the majority of the data points.

Defining a New Equation (1)

There are very few examples found in the literature where both raw data and calculated values are included. However, the paper by Roughton *et al.* (1955)⁹⁴ is notable because the raw data are presented with a set of calculated results. Also the equation to be fitted is somewhat complex. Because of this, these data provide a good test for a data analysis package. The Roughton paper examines the binding of oxygen to sheep hemoglobin, as a function of the pressure of oxygen. Hemoglobin has four non-equivalent sites, and the binding is described by the equation:

$$\text{Saturation} = \frac{K_1 p + 2K_1 K_2 p^2 + 3K_1 K_2 K_3 p^3 + 4K_1 K_2 K_3 K_4 p^4}{4(1 + K_1 p + K_1 K_2 p^2 + K_1 K_2 K_3 p^3 + K_1 K_2 K_3 K_4 p^4)} \times 100\%$$

In this equation, p is the pressure of oxygen, and K_1 to K_4 are the equilibrium constants for the binding of four successive moles of oxygen per mole of hemoglobin.

The data are in the file "Hemoglobin saturation". Read this file in using the **File ▶ Open** command. An interesting feature of these data is that some data points were collected with greater accuracy than others. To account for this it is necessary to use Explicit weighting.

⁹⁴ Roughton, F.J.W., Otis, A.B. & Lyster, R.L.J. (1955) *Proc. Roy. Soc. London Ser B* **144**, 29-54.

Defining the Oxygen Binding Equation

To define this equation we need to create a new definition in our workspace using the **Add ▶ Equation** command. (See Chapter 9 for a complete description of defining equations.)


Definition of the new equation requires that we make appropriate entries in this item window.

X Variables	Symbol	Description
p	p	pO_2
Parameters	Symbol	Description
K_1	$K1$	$K1$ value
K_2	$K2$	$K2$ value
K_3	$K3$	$K3$ value
K_4	$K4$	$K4$ value
Y name:	Sat	
Definition: ⁹⁵	$\text{Sat} = 100 * (K1 * p + 2 * K1 * K2 * p * p + 3 * K1 * K2 * K3 * p * p * p + 4 * K1 * K2 * K3 * K4 * p * p * p * p) / (4 * (1 + K1 * p + K1 * K2 * p * p + K1 * K2 * K3 * p * p * p + K1 * K2 * K3 * K4 * p * p * p * p))$	


The appearance of the equation definition item after entering this definition, and changing the name of the item, is shown below.

⁹⁵ As this definition is for demonstration purposes only, it is defined in this way as this is the nearest form to that of the equation given. Although the definition would function in this form, it is rather inefficient to write it like this, and a better approach would be to enter the mathematically simpler:


$$\text{Sat} = 100 * K1 * p * (1 + K2 * p * (2 + K3 * p * (3 + K4 * p * 4))) / (4 * (1 + K1 * p * (1 + K2 * p * (1 + K3 * p * (1 + K4 * p))))$$

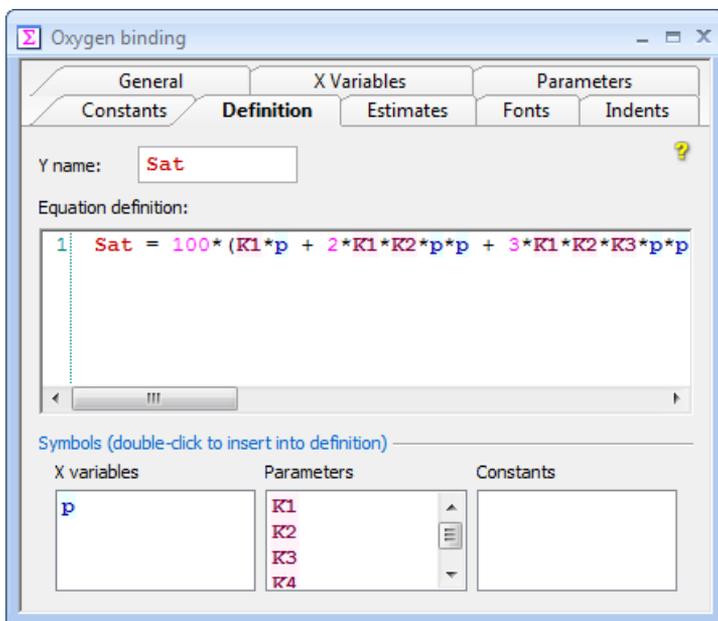

To fit these data, add a new non-linear fitting item using the **Add ▶ Non-linear fit** command, then complete the Settings page as below:

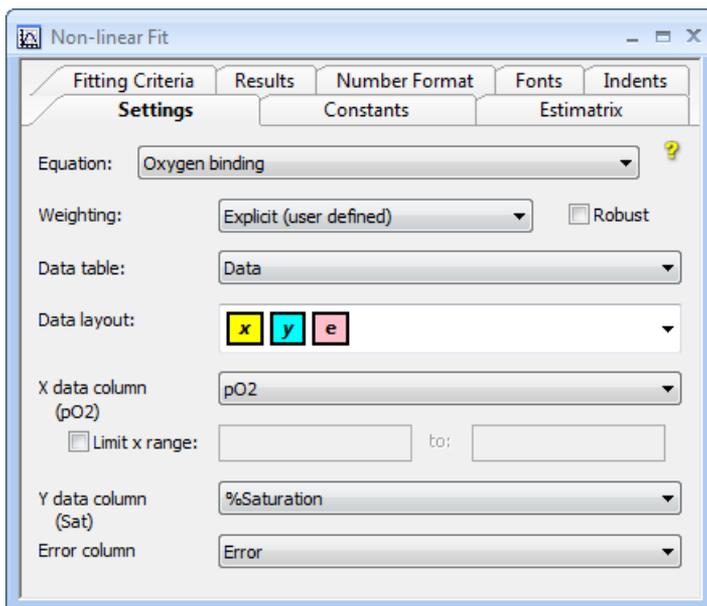

Note the use of Explicit weighting, which requires the selection of the column that holds the error values. The results of the analysis are as follows:

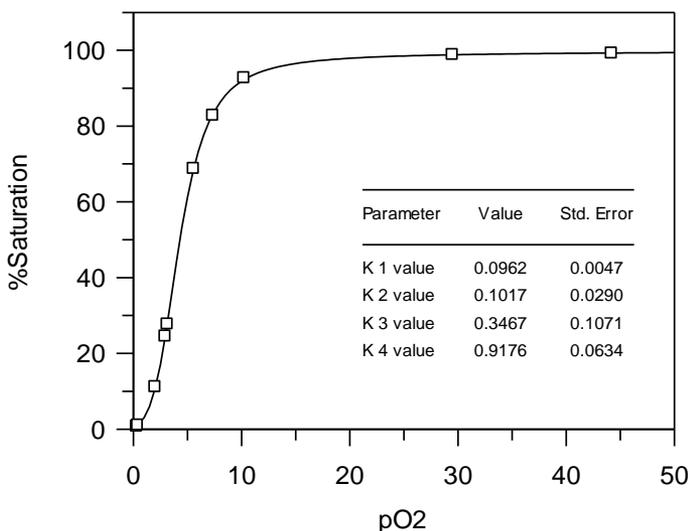

For comparison, the results quoted by Roughton *et al.* are given below.

	Parameter value	Standard Error
K_1	0.0961	0.0049
K_2	0.1040	0.0293
K_3	0.3352	0.1018
K_4	0.9263	0.0648

It should be borne in mind that Roughton *et al.* performed their regression analysis without the aid of personal computers or even pocket calculators!

Defining a New Equation (2)

The hybridization of RNA and cDNA molecules is used to provide information about the genome organization and complexity of an organism. This hybridization can be described by the equation:

$$\frac{d}{D_0} = B + \sum_{i=1}^n P_i \left[1 - e^{-0.693(Rot/Rot_{0.5(i)})} \right]$$

Here, d/D_0 is the fractional hybridization at any given value of Rot (moles second L^{-1} of nucleotides of RNA), B is a background value, P_i is the proportion of the total cDNA hybridized to the i th component, and $Rot_{0.5(i)}$ is the Rot value at which 50% of the cDNA of the i th component is hybridized. In practice, the data obtained may be fitted assuming three binding

components, giving a total of seven independent variables. Green *et al.*⁹⁶ have presented a sample of such data; these can be found in the file "DNA hybridization".

Defining the Equation

The equation can be defined as follows:

X Variables	Symbol	Description
Rot	Rot	Rot value
Parameter	Symbol	Description
B	B	Background
P ₁	P1	Proportion Component 1
Rot _{0.5(1)}	R1	Rot 0.5 (1)
P ₂	P2	Proportion Component 2
Rot _{0.5(2)}	R2	Rot 0.5 (2)
P ₃	P3	Proportion Component 3
Rot _{0.5(3)}	R3	Rot 0.5 (3)
Y name:	Sat	
Definition:	$\text{Sat} = B + P1 * (1 - \exp(-0.693 * \text{Rot} / R1))$ $+ P2 * (1 - \exp(-0.693 * \text{Rot} / R2))$ $+ P3 * (1 - \exp(-0.693 * \text{Rot} / R3))$	

Fitting the Equation


The following initial estimates are appropriate to analyze the data.

B	0.05
P ₁	0.3
Rot _{0.5(1)}	0.2
P ₂	0.3
Rot _{0.5(2)}	40.0
P ₃	0.3
Rot _{0.5(3)}	600.0

It is interesting to note that when Green *et al.* published their results in 1982, using an Apple II computer and employing a "Patternsearch" algorithm they reported that these data were fitted in approximately two hours of calculation. Compare this to the amount of time taken by GraFit using a modern PC, which should fit the data in a fraction of a second.

The calculated result is displayed below.

⁹⁶ Green, S., Field, J.K., Green, C.D. & Beynon, R.J. (1982) *Nucleic Acids Research* **10**, 1411-1421.

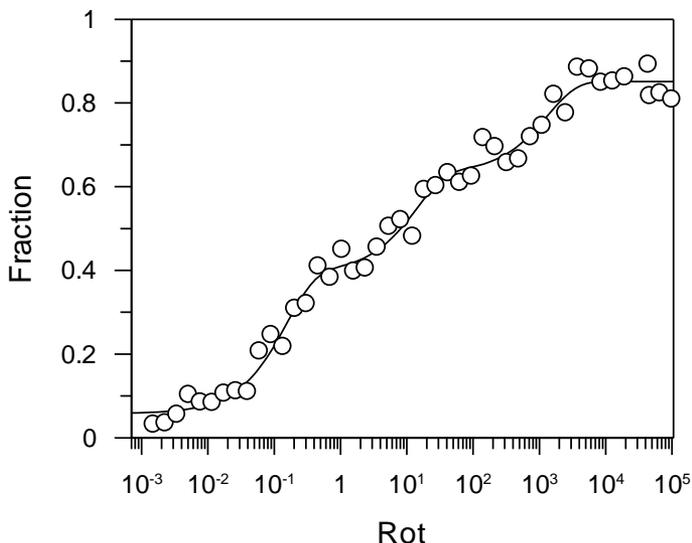

Multiple Regression

Multiple regression is used when the observed data values depend on more than one dependent (x) value. As an example of this, we will examine the data listed by Toney & Kirsch⁹⁷. These data consist of observations of the rate constant of a reaction, k , as a function of the $\text{p}K_a$ of the group involved in the reaction and the molecular volume of this group, Vol. The data are described by the equation:

$$\log k = \beta \cdot K_a + V \cdot \text{Vol} + c$$

Here, β , V and c are unknowns in the equation, which is a multiple linear equation having two independent variables. The data presented in the paper are in the file "Multiple linear." It is possible to fit these data by defining a new specific equation. The required definitions are as follows.

Equation Name:	Multilinear kinetics	
X Variables	Symbol	Description
$\text{p}K_a$	pKa	pKa value
Vol	Vol	Molecular volume
Parameter	Symbol	Description
β	Beta	Beta
V	V	V
c	C	Constant
Definition:	beta*pKa + V*Vol + C	

⁹⁷ Toney, M.D. & Kirsch, J.F. (1989) *Science* **243**, 1485-1488.

The results obtained from fitting these data are:

Variable	Value	Standard Error
β	0.39	0.05
V	-0.055	0.005
c	-0.7	0.5


Displaying Multidimensional Data

Displaying data sets that have more than one x axis can be difficult. One solution is to remove all but one dimension from a plot. For example, in the above data, a plot of $(\log k - V.Vol)$ versus pK_a will result in a linear graph with a slope β and intercept c . We can create the $(\log k - V.Vol)$ data using the data manipulation abilities of GraFit; the result for V is obtained from the curve fitting.

The first step (assuming that the data have already been fitted) is to allow the calculated V value to be accessed. This can be done by creating a data column that displays the fitting results. To do this, select a blank column and choose the **Column Type ▶ Data Fitting ▶ Results** command (**Table** tab⁹⁸), and select the appropriate Non-linear Fit item.

The second step is to select the next blank column and choose the **Column Type ▶ Data Processing ▶ Calculated Function** command to set up a calculation. The following calculated function allows the column to be filled with the $(\log k - V.Vol)$ entries.


⁹⁸ **Table ▶ Column Type ▶ Data Fitting ▶ Results** command in legacy mode

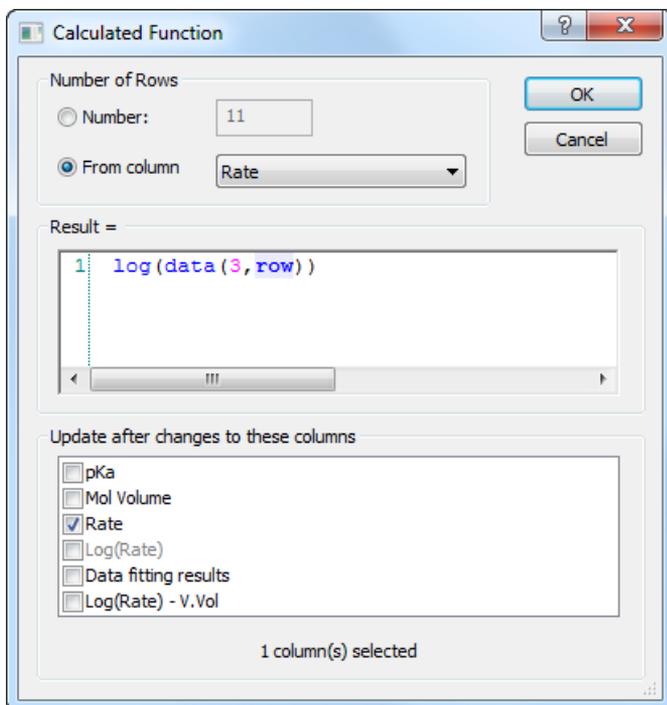

Note the use of the `Data()` function to access data present in the data table. The appearance of the data table that uses this calculated function is shown below.

Type	1	2	3	4	5	6	
Column Name	x(1) Data pKa	x(2) Data Mol Volume	Rate	Calculated Log(Rate)	Results Data fitting results	Calculated Log(Rate) - V.Vol	
1	10.6000	42.1000	40.0000	1.6021	Beta	0.3934	3.9142
2	10.6000	60.9000	1.3000	0.1139	V	-0.0549	3.4587
3	10.5000	79.8000	0.2700	-0.5686	Constant	-0.6977	3.8141
4	10.6000	98.7000	0.0090	-2.0458			3.3750
5	10.0000	74.6000	0.0840	-1.0757			3.0214
6	9.5000	71.5000	0.0780	-1.1079			2.8190
7	9.2000	23.2000	21.0000	1.3222			2.5964
8	9.0000	64.4000	0.2570	-0.5901			2.9469
9	7.7000	70.5000	0.0108	-1.9666			1.9054
10	5.7000	71.6000	0.0073	-2.1367			1.7957
11	5.3000	51.1000	0.0500	-1.3010			1.5055

The result, plotting $(\log k - V.Vol)$ versus pK_a and using a transformation definition to display the appropriate line through the data, is shown below.

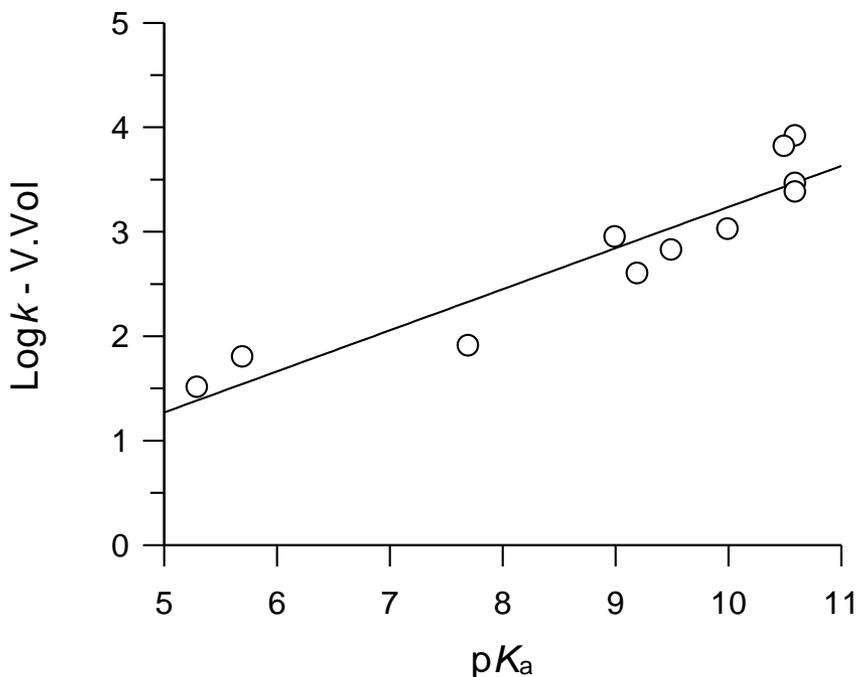

Reading Values from a Calibration Curve

A calibration curve records the variation of an observed value when an experimental condition is altered in a defined manner. For example, an instrumental reading may be measured for various known concentrations of a solution. The standard curve that is produced can then be used to allow, in this example, the concentrations present in an unknown solution to be determined by extrapolating values from the curve.

These processes can be performed using GraFit by the following procedure.

1. Fit the data to an equation that describes the observed trend.
2. Use the **Column Type ▶ Calculated Values ▶ x Values** or **Column Type ▶ Calculated Values ▶ y Values** commands (**Table** tab⁹⁹) to read values from the fitted curve.


The sample data file "Standard curve" holds the analysis described below.

Choosing an Equation to Fit the Data

In many circumstances, the relationship between the x and y data values will be well defined. In these cases, the data will be fitted by linear or non-linear regression in the usual manner. However, some calibration data show a clear trend of how y varies with x, but are not necessarily defined by a well-characterized equation.



⁹⁹ **Table** menu in legacy mode

For example, the data below show the absorbance reading given by varying concentrations of a compound, and a straight line fit to the data.

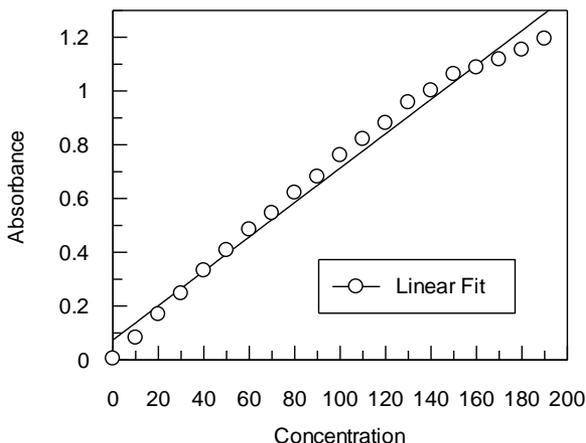


In theory, the plot of absorbance versus concentration should be linear, but the data shown above describe a shallow curve. While it might be of interest to find out why this effect occurs, the only important consideration here is to allow readings to be taken from the calibration curve. Quite clearly, fitting the data to a straight line is not suitable.

Under these circumstances, the best recourse is to fit the data to an arbitrary equation that allows a smooth curve to be drawn through the data. As there is no theoretical basis for the equation, any calculated parameter values will have no physical significance. The choice of equation will depend upon the shape of the data, and it is possible that more than one equation would prove suitable. In the case shown above it might be appropriate to use a polynomial, or any equation that can produce a gentle curve. Below is shown the results of fitting the data to a second order polynomial ($y = a + bx + cx^2$) using a Linear Fit item, and to a first order rate equation by non-linear regression. The latter equation has no possible theoretical justification, but has a suitable mathematical form to allow a smooth curve to run through the data.

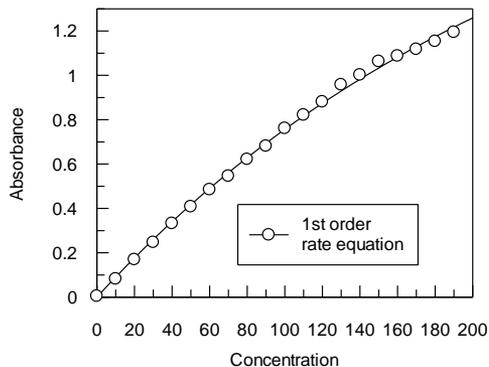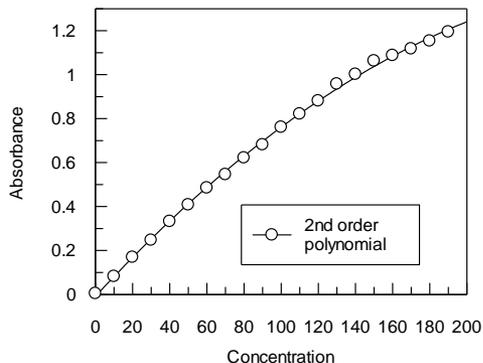

Reading Values from the Fitted Curve


Either of the fitted curves shown above will allow the concentration values that correspond to a series of measured absorbance values to be read from the curve. This is performed by entering the unknown absorbance values into a data column, then generated the calculated values in a separate column using the **Column Type ▶ Calculated Values ▶ x Values** command (**Table tab**⁹⁹).

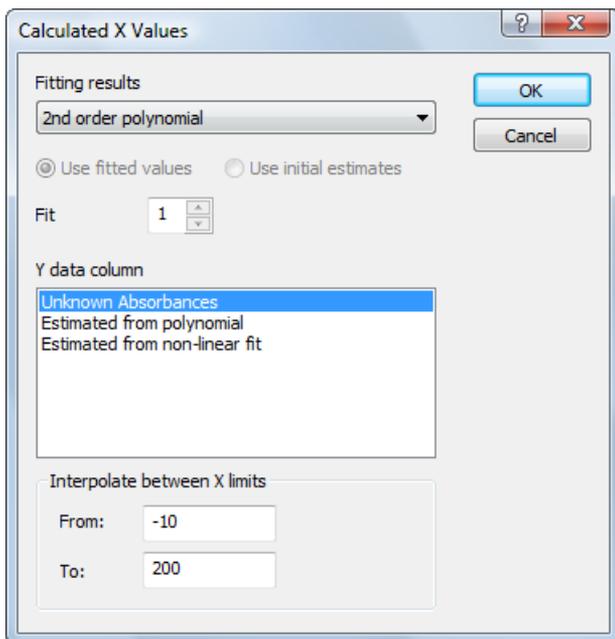

Fitting Several Data Sets with Some Parameters Common and Some Specific

To illustrate, imagine that you have collected two linear data sets, and that you already know that for this experiment the slope is the same for each although the intercepts are different. Using GraFit it is possible to fit these two data sets to find the best slope that is common to each, and the two intercept values that are specific to the two data sets. Indeed, it is possible to fit two or more data sets to any of the normal equations that GraFit can use and to share one or several parameters between the data sets.

The example below shows how this is done for the paired linear data. Two raw data sets are illustrated below (these data are in the example file “Paired linear regression”). Fitting these two data sets individually results in similar, but not identical, values for the slope. This is shown in the figure below.

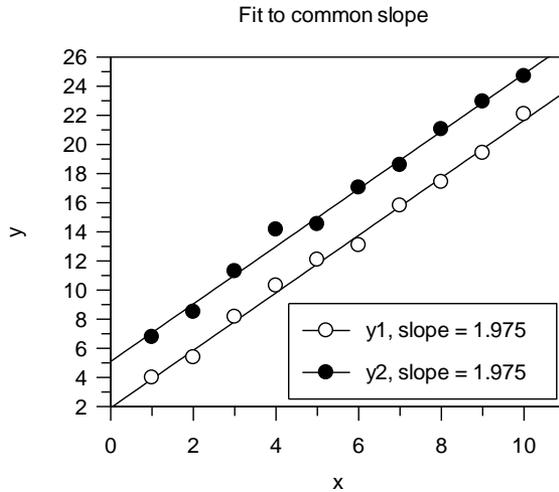

If the data sets are known to have the same slope, this can be incorporated into the data fitting as follows.

Method Used to Share Variables between Data Sets

Sharing variables between data sets involves using an equation with a second independent (x) variable, which is defined in such a way that that this value acts as a switch to decide which parameters apply to the data value. For the simple example described above, the independent variables are x and `dataset`; parameters $a1$ and $a2$ are the two intercepts, and b is the common slope. The equation is then defined as

```
If dataset=1 Then
  y = a1 + b*x
  Return
Endif
y = a2 + b*x
```

How does this work? When data from the first data set are accessed, the value of `dataset` is 1 and the equation returns


$$a1 + b*x$$

When data from the second data set are accessed, the value of `dataset` is 2 and the equation returns

$$a2 + b*x$$

Therefore the value of $a1$ is determined, in effect, from only the first set of data, and $a2$ from the second. However, the slope b is applicable to all the data. The graph below shows these data fitted to a common slope in this manner.

Fit to individual slope

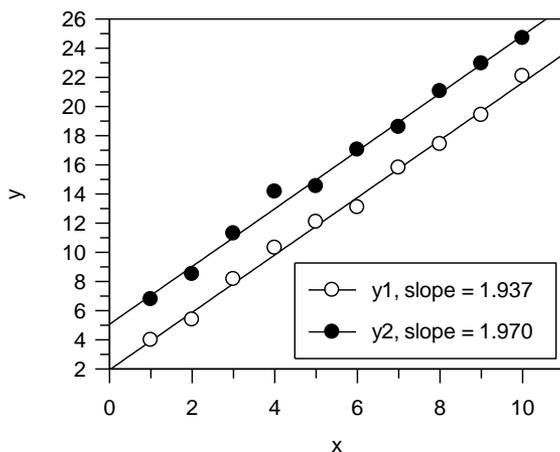

The dataset setting is set using the **Column Type ▶ Data Types ▶ y Data with Associated Value** command (Table tab⁹⁹), resulting in the appearance shown below.

Type	1	2	3
Column Name	x	y1	y2
1	1.0000	3.9814	6.7743
2	2.0000	5.3816	8.5039
3	3.0000	8.1534	11.2915
4	4.0000	10.3019	14.1558
5	5.0000	12.0741	14.5192
6	6.0000	13.0814	17.0320
7	7.0000	15.7954	18.5855
8	8.0000	17.4193	21.0429
9	9.0000	19.4161	22.9451
10	10.0000	22.0809	24.6848

It is simple to extend this approach to many data sets. The equation to find the best slope through four data sets and a separate intercept for each is:

```

If dataset=1 Then
  y = a1 + b*x
  Return
Endif

```

```

If dataset=2 Then
  y = a2 + b*x
  Return
Endif
If dataset=3 Then
  y = a3 + b*x
  Return
Endif
' assumes that dataset must be 4
y = a4 + b*x

```

Example — Paired Denaturation Curves

This technique of sharing parameters between two or more data sets is not limited to linear equations. Any equation may be recast in this manner. For example, the equation below describes the denaturation of a protein as a function of the concentration of a denaturing solvent¹⁰⁰.

$$F = 1 - \frac{\exp\left(\frac{m[D] - \Delta G}{RT}\right)}{1 + \exp\left(\frac{m[D] - \Delta G}{RT}\right)}$$

In this equation, F is the fraction of the protein that is denatured, $[D]$ is the concentration of denaturant (e.g. guanidinium hydrochloride or urea), R is the gas constant, T is the temperature, ΔG is the stability of the protein, and m is a slope factor. Finding the difference in stability between two proteins is often more important than determining the absolute stability. If the proteins are highly related, it can be assumed that the value for m is the same for each, but ΔG differs.

Using the technique described above, it is possible to fit two data sets to a common m value but different ΔG values. In this case, the equation to be used can be defined as

```

If dataset=1 Then
  a = exp((m*D - G1) / (8.314*T))
Else
  a = exp((m*D - G2) / (8.314*T))
Endif
y = (1 - a / (1 + a))


```

The x variable is the concentration of denaturant, D , the slope factor is m and the two values of ΔG are $G1$ and $G2$. The temperature is T , and should be defined as a constant. The

¹⁰⁰ Jandu, S.K., Ray, S.R., Brooks, L. & Leatherbarrow, R.J. (1990) *Biochemistry* **29**, 6265-6269.

dataset column (x_2) must hold the value 1.0 if the x and y columns are from the first data set, and 2.0 if they are from the second. (Clearly, this can be extended to further data sets if required.) The value a in the above definition is a temporary value used to simplify the calculation.

The figure below shows the guanidinium hydrochloride unfolding of a native protein, and of a mutant variant of this protein. The mutation has caused the protein to become less stable, making it unfold at lower concentrations of denaturant. The solid lines indicate the fitted curves, using the above equation. The difference in stability can therefore be determined.




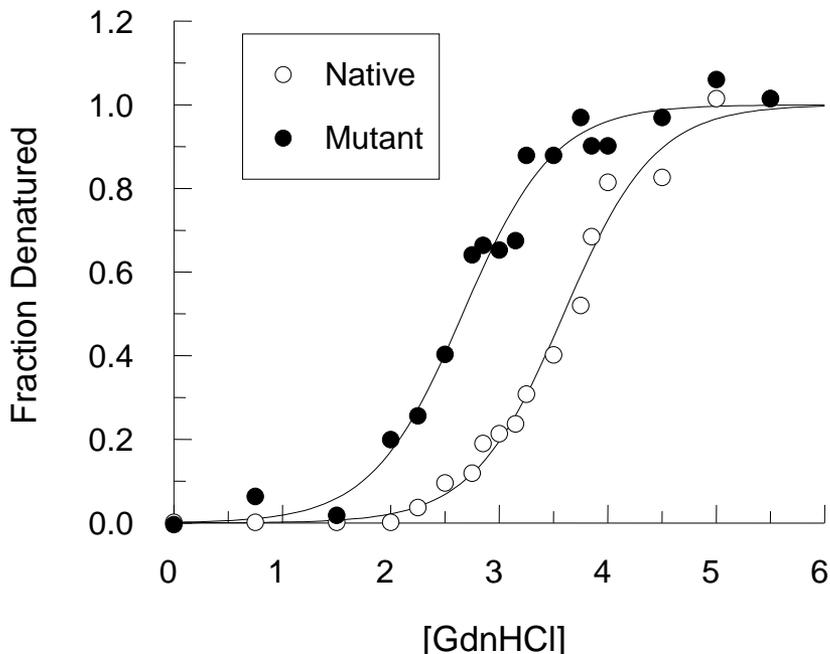


Enzyme Inhibition Samples

The final set of samples concern the analysis of enzyme inhibition data. These sample files are relatively complex and provide the basis for your own enzyme inhibition analyses. If you need to fit inhibition data, we recommend that you examine the files to see how they analyses have been set up.

Competitive Inhibition Fitting

The files "Competitive Inhibition" and "Competitive Inhibition LB" show how it is possible to fit data to a competitive inhibition model. They differ in that the first of these files shows an Eadie Hofstee transformation, while the second uses a Lineweaver Burk transformation.

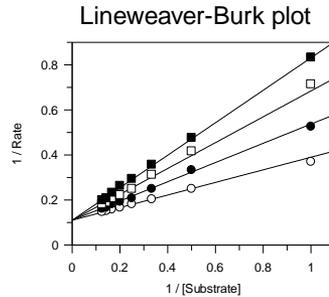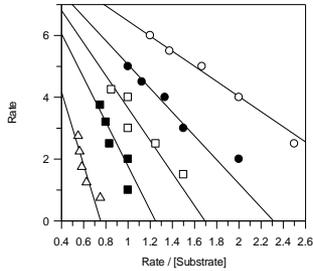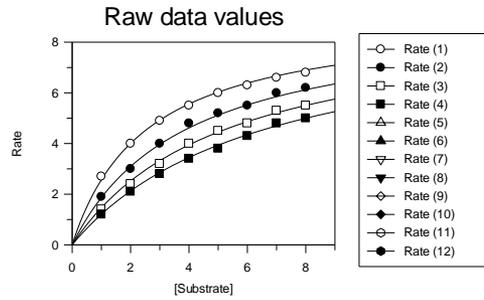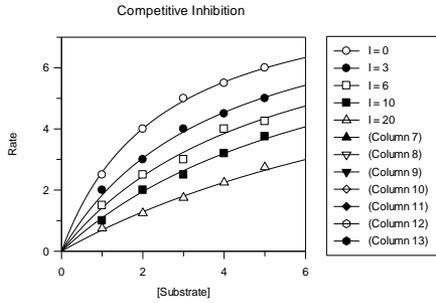

Reduced Chi²: 0.0094

Parameter	Value	Std. Error
Vmax	8.9147	0.2682
Km	2.4453	0.1720
Ki	5.2035	0.2681


Parameter	Value	Std. Error
Vmax	9.0873	0.1213
Km	2.5422	0.0973
Ki	1.9008	0.0710

Contents of "Competitive Inhibition" sample data file

Contents of "Competitive Inhibition LB" sample data file

Production of a Dixon Plot

The sample data file, "Dixon plot", analyses a set of inhibition data to a competitive fit, then displays the data using a Dixon plot. The data file shows how it is possible to generate complicated rearrangements of the underlying data (which is needed to generate the Dixon plot).

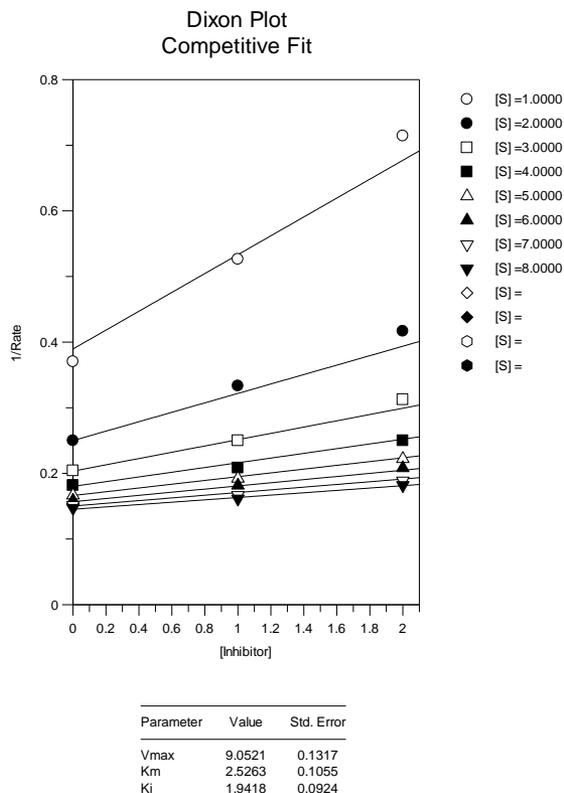


Contents of "Dixon plot" sample data file

Comprehensive Inhibition Modeling

The final sample data file is the most complex that we provide. The "Enzyme Inhibition" file carries out a full analysis of four different inhibition models for the same data set, allowing assessment of whether the data are best described by a Competitive, a Non-competitive, a Mixed or an Uncompetitive inhibition mechanism. The output for the sample data set used is reproduced below. For this particular set of data it can be seen that a competitive inhibition model best fits the data, as this produces the best fit (both by eye and by examination of the reduced χ^2 values). (The mixed type inhibition gives a similar fit quality, but the additional parameter in this model does not significantly improve the fit. The high error value of the additional K_i' value is also a good indication that the extra parameter is not appropriate for these data.)


Competitive model

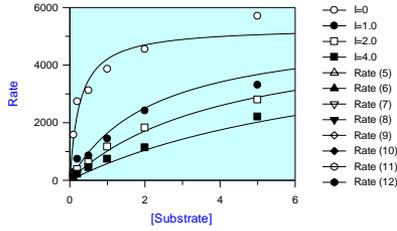

Reduced Chi²: 64387.8893



Parameter	Value	Std. Error
Vmax	5320.8514	211.8898
Km	0.2646	0.0423
Ki	0.1336	0.0203

Non-Competitive model

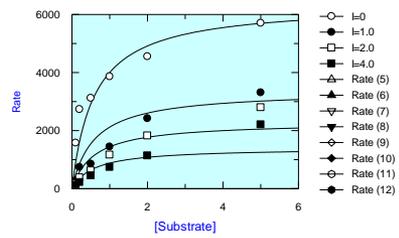


Reduced Chi²: 232468.9108

Parameter	Value	Std. Error
Vmax	6365.0292	519.7938
Km	0.5911	0.1403
Ki	1.1252	0.1826

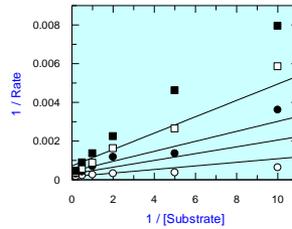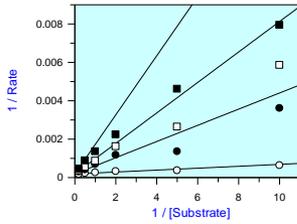


Mixed type model

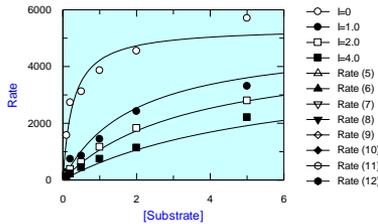

Reduced Chi²: 65843.5622



Parameter	Value	Std. Error
Vmax	5379.6409	226.6466
Km	0.2752	0.0460
Ki	0.1512	0.0330
Ki' value	13.3794	11.8828

Uncompetitive model

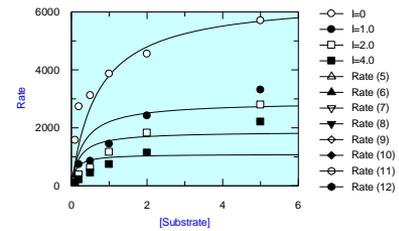

Reduced Chi²: 415505.5294


Parameter	Value	Std. Error
Vmax	6478.1123	819.7069
Km	0.7032	0.2509
Ki	0.8107	0.2008

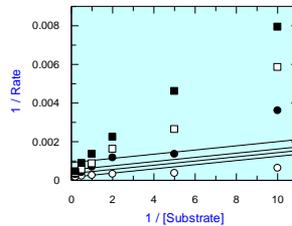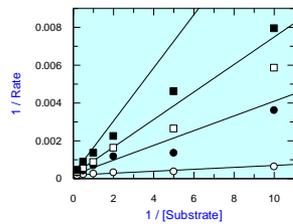

Equations Provided

GraFit comes with a large number of equations provided. These equations may be edited to see how they have been defined, and further equations can be added. This chapter briefly describes the equations that come with the program. It does not attempt to discuss the detailed background of each equation, as it is assumed that the user of a specific equation will be familiar with this.

The equations and transformations that come with the program are installed into the Group Item Store directory, and can be accessed using the **Add Stock Item** command¹⁰¹. Equations are grouped into separate sub-directories, which appear as different pages on the Add Item dialog box.

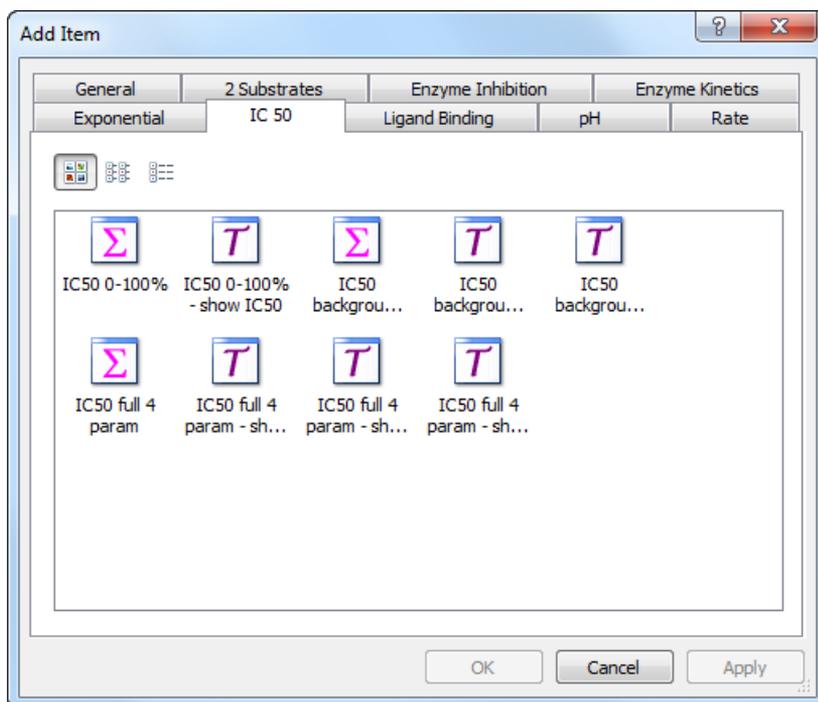

¹⁰¹ The **File** ► **Add Stock Item** command in legacy mode

Equation and transformation definitions are distinguished by the nature of their icons. The following sections describe the contents of these various pages.

2 Substrates

The definitions on this page are designed to analyze enzyme kinetic data involving two substrates.


Σ Ping pong

Use this equation where a ping pong mechanism applies. The equation relating the rate of reaction, v , to the concentration of substrates A and B and the constants K_A and K_B is:

$$v = \frac{V_{max}[A][B]}{K_B[A] + K_A[B] + [A][B]}$$

T Ping pong Lineweaver Burk

This transformation generates a Lineweaver Burk plot for the fitted data.

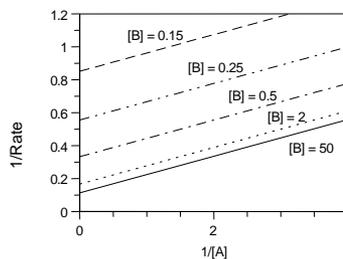



Σ Ternary

Use this equation where a ternary complex is formed. The equation relating the rate of reaction, v , to the concentration of substrates A and B and the constants K'_A , K_A and K_B is:

$$v = \frac{V_{max}[A][B]}{K'_A K_B + K_B[A] + K_A[B] + [A][B]}$$

T Ternary Lineweaver Burk

This transformation generates a Lineweaver Burk plot for the fitted data.

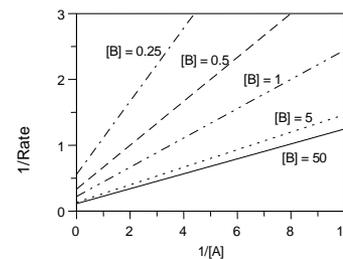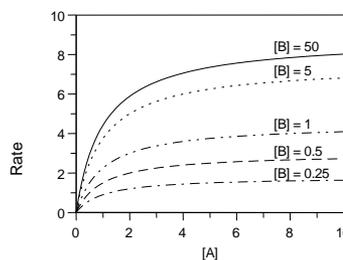


Enzyme Inhibition

Σ Competitive Inhibition


In this multidimensional equation, rate is monitored versus [substrate] and [inhibitor]. The equation for competitive inhibition is:

$$v = \frac{[S][E_0]k_{cat}}{[S] + K_m(1 + [I]/K_i)}$$

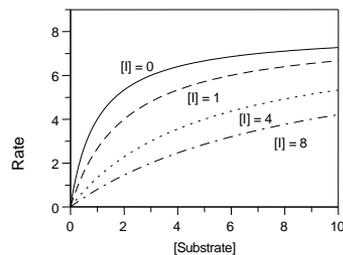

T Competitive Eadie Hofstee

Use this to display the Eadie Hofstee plot.

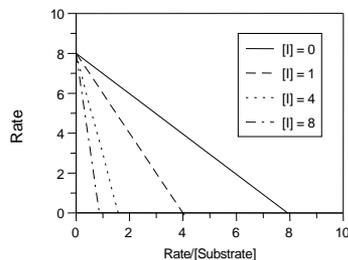

T Competitive Lineweaver Burk

This transformation displays the Lineweaver Burk plot.

Σ Mixed Inhibition

In this most general inhibition equation, rate is monitored versus [substrate] and [inhibitor], using the equation:

$$v = \frac{V_{max} \cdot [S]}{K_m \left(1 + \frac{[I]}{K_i}\right) + \left(1 + \frac{[I]}{K_i'}\right) [S]}$$

T Mixed Lineweaver Burk

This transformation displays the Lineweaver Burk plot.

T Mixed Eadie Hofstee

Use this to display the Eadie Hofstee plot.


Σ No inhibition

This equation is the normal Michaelis Menten enzyme kinetics equation.

Σ Non-competitive Inhibition


In this multidimensional equation, rate is monitored versus [substrate] and [inhibitor]. The equation for non-competitive inhibition is:

$$v = \frac{V_{max} \cdot [S] \cdot \frac{1}{1 + [I]/K_i}}{K_m + [S]}$$

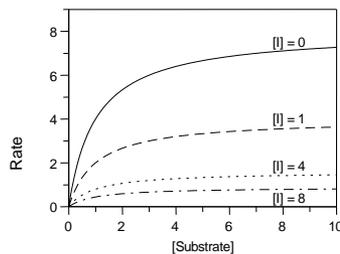

T Non-competitive Eadie Hofstee

Use this to display the Eadie Hofstee plot.

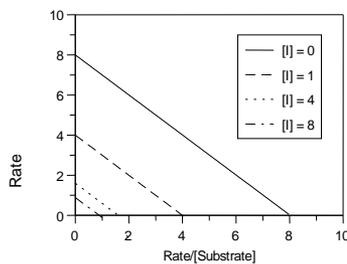

T Non-competitive Lineweaver Burk

This transformation displays the Lineweaver Burk plot.

Σ Non-tight binding

This equation describes standard, or non-tight binding inhibition. Use in contrast to the tight binding inhibition equations (see below).

$$v = v_0 \left(1 - \frac{[I]_0}{K_d + [I]_0} \right)$$

Σ Non-tight binding (offset)

This equation describes standard, or non-tight binding inhibition. Use in contrast to the tight binding inhibition equations.

The equation assumes that inhibition does not proceed to 100%, but that there is a background rate that is not being inhibited. In real-world situations, this can happen (for example, due to residual impurities in the enzyme preparation). If the background value is significant, take care to question whether it is real.

$$v = (v_0 - \text{background}) \cdot \left(1 - \frac{[I]_0}{K_d + [I]_0}\right) + \text{background}$$

Σ Tight binding

This equation describes tight binding inhibition, where [I] is comparable to [E].

Note: The concentration of enzyme is treated as a parameter, rather than a constant. This is because even in circumstances where you “know” the enzyme concentration, it is unlikely that this will be known with sufficient precision. After fitting your data to this equation, always ensure that the fitted [E] value is sensible; if it is not (far higher or lower than the actual [E]) then it is possible that you need to consider an alternative inhibition mechanism.

$$v = v_o \frac{[E]_0 - [I]_0 - K_d + \sqrt{([E]_0 - [I]_0 - K_d)^2 + 4[E]_0 K_d}}{2[E]_0}$$

Morrison, J. F. (1969) Kinetics of the reversible inhibition of enzyme-catalysed reactions by tight-binding inhibitors. *Biochim. Biophys. Acta* **185**, 269-286.

Cha, S. (1975) Tight binding inhibitors - I. Kinetic behavior. *Biochem. Pharmacol.* **24**, 2177-2185.

Σ Tight binding (offset)

This equation describes tight binding inhibition, where [I] is comparable to [E]. This equation assumes that inhibition does not proceed to 100%, but that there is a background rate that is not being inhibited. In real-world situations, this can happen (for example, due to residual impurities in the enzyme preparation). If the background value is significant, take care to question whether it is real.

See the description of the tight binding equation for more information.

$$v = (v_0 - \text{background}) \frac{[E]_0 - [I]_0 - K_d + \sqrt{([E]_0 - [I]_0 - K_d)^2 + 4[E]_0 K_d}}{2[E]_0} + \text{background}$$

Σ Uncompetitive Inhibition

In this multidimensional equation, rate is monitored versus [substrate] and [inhibitor]. The equation for uncompetitive inhibition is:

$$v = \frac{\frac{V_{max}}{\left(1 + \frac{[I]}{K_i}\right)} \cdot [S]}{\frac{K_m}{\left(1 + \frac{[I]}{K_i}\right)} + [S]}$$

T Uncompetitive Lineweaver Burk

This transformation displays the Lineweaver Burk plot.

T Uncompetitive Eadie Hofstee

Use this to display the Eadie Hofstee plot.


Enzyme Kinetics

Σ Allosteric kinetics (Hill)

The Hill equation describes enzyme kinetics for an allosteric enzyme. For equation is:

$$v = \frac{V_{max}[S]^n}{K^n + [S]^n}$$

The Hill coefficient, n , equals 1 if standard kinetic behavior is followed, but differs significantly from 1 if the enzyme kinetics are cooperative.

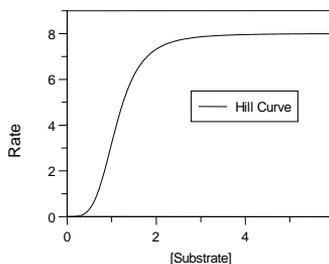

T Allosteric kinetics (Hill), Linear plot

This transformation is used to display the linearized Hill plot where $v/(V_{max} - v)$ is plotted versus $[S]$. The slope is $[S]^n / K^n$.

Σ Enzyme Kinetics


This equation fits the basic Michaelis-Menten enzyme kinetics equation:

$$v = \frac{V_{max}[S]}{K_m + [S]}$$

In this equation, the rate, v , varies with the concentration of substrate, $[S]$, depending upon the maximal rate, V_{max} , and the Michaelis constant, K_m .


T Enzyme Kinetics, Eadie Hofstee

This transformation draws the linear Eadie-Hofstee plot for fitted enzyme kinetic data.

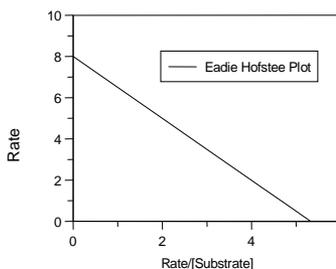


T Enzyme Kinetics, Lineweaver Burk

This transformation draws the linear Lineweaver Burk plot for fitted enzyme kinetic data.

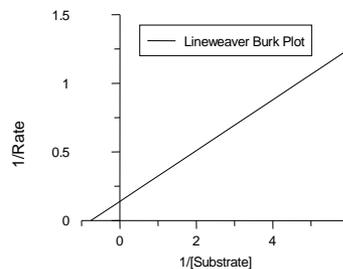

T Enzyme Kinetics, show Km

This transformation displays the Michaelis constant, K_m . For best results set the curve resolution to **High** when using this transformation.

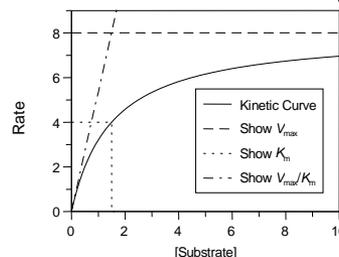

T Enzyme Kinetics, show Vmax

Displays the limiting rate, $V_{max} (= k_{cat}[E_0])$.

T Enzyme Kinetics, show Vmax over Km

This transformation displays the linear portion at the beginning of the graph where $[Substrate] \ll K_m$.

Σ Kinetics, contaminating S

Occasionally the reagents etc. used in an enzyme assay will contain traces of contaminating substrate. Therefore the rate will vary with added substrate, S_{add} as:

$$v = \frac{V_{max}([S_{add}] + [S_{con}])}{K_m + [S_{add}] + [S_{con}]}$$

In this equation, $[S_{con}]$ is the concentration of contaminating substrate.

Exponential


This page contains equations for calculating the rate constants of various exponential decay processes.

Σ Double Exponential

If two decay processes occur simultaneously the observed values decay with time following the double exponential decay equation:

$$y = A_{0(1)} \cdot e^{-k_1 t} + A_{0(2)} \cdot e^{-k_2 t}$$

Here, the two processes have separate A_0 and k values.

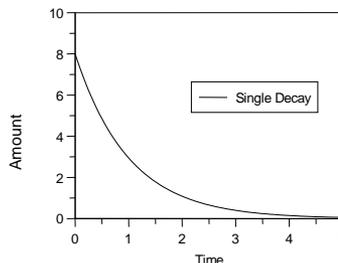

T Double Exponential - phase 1

This displays the first of the two phases in isolation.

T Double Exponential - phase 2

This displays the second of the two phases in isolation.

Σ Double Exponential with offset


This equation is used when a double exponential process decays to a non-zero limit value.

Σ Single Exponential

The equation for a single exponential decay is:

$$y = A_0 \cdot e^{-k_1 t}$$

The observed values decay with time, t , from an initial value, A_0 , at $t = 0$, to 0 at $t = \infty$. The rate constant for the decay is k . (The half-life for the process, $T_{1/2}$, is given by $\log_e 2/k$.) To display a linearized plot of this equation, set the y axis to have a logarithmic scale.

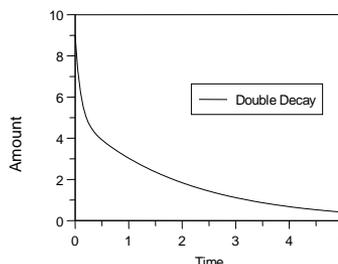

Σ Single Exponential with offset

If the observed values do not decay to 0, but to some defined offset or background value the equation for a single exponential process is

$$y = A_0 \cdot e^{-k_1 t} + \text{offset}$$

In this equation, 'offset' is the limit of the decay.

Σ Triple Exponential

For a triple exponential decay process the equation is:

$$y = A_{0(1)} \cdot e^{-k_1 t} + A_{0(2)} \cdot e^{-k_2 t} + A_{0(3)} \cdot e^{-k_3 t}$$

Σ Triple Exponential with offset

This equation is used when a triple exponential process decays to a non-zero limit value.

IC 50

These equations are used for analyzing radio-immuno assays or Dose-response curves.

Σ IC50 0-100%

This equation fits inhibition data to a 2-parameter equation, where the lower data limit is 0, and the upper data limit is 100, i.e. the data are both background and range corrected.

$$y = \frac{100\%}{1 + \left(\frac{x}{IC_{50}}\right)^s}$$

In this equation, s is a slope factor. The equation assumes that y falls with increasing x .

T IC50 0-100% - show IC50

This transformation displays the IC_{50} value on the plot. For best results set the curve resolution to **High** when using this transformation.

Σ IC50 background corrected

This equation fits inhibition data to a 3-parameter equation, where the lower data limit is 0, i.e. the data are background corrected.

$$y = \frac{\text{Range}}{1 + \left(\frac{x}{IC_{50}}\right)^s}$$

In this equation, s is a slope factor. The equation assumes that y falls with increasing x .

T IC50 background corrected - show IC50

This transformation displays the IC_{50} value on the plot. For best results set the curve resolution to **High** when using this transformation.

T IC50 background corrected - show limit


This transformation displays the limit value on the plot.

Σ IC50 full 4 param

Fits data to a full 4-parameter equation, which is usually the best way to fit results from RIA analysis or inhibition data (see Halfman, 1981¹⁰²). The binding or inhibition is then characterized in terms of an IC₅₀ value:

$$y = \frac{\text{Range}}{1 + \left(\frac{x}{IC_{50}}\right)^s} + \text{background}$$

Here *Range* is the maximum *y* range, and *s* is a slope factor. The *x* axis represents the concentration of analyte. Data fitted to this equation are usually displayed with a logarithmically scaled *x* axis. For convenience, the equation as defined attempts to provide initial estimates for the parameter values.

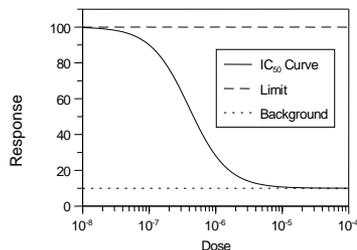

Note Estimating starting values is difficult when there are 4 unknown parameters. The method used within this definition works most of the time, but can give occasional poor fits. Under such circumstances, enter initial estimates explicitly in the non-linear fit window.

The curve displays inhibition with increasing dose. It is also possible to use the same equation for increasing response with dose. In this case, the *y* range is a negative value, and the “Background” is the maximum response.

T IC50 full 4 param - show background

This transformation displays the calculated background value.

T IC50 full 4 param - show IC50

This transformation displays the calculated background value.

T IC50 full 4 param - show limit

This transformation displays the limiting value of the data.

Ligand Binding

This page contains various equations describing the binding of ligands to surfaces, receptors, enzymes etc.


¹⁰² Halfman, C.J. (1981) *Methods in Enzymology* **74**, 481-508

Σ 1 Site

This is the simplest form of binding equation. The equation is


$$y = \frac{[L] \cdot Cap}{K_d + [L]}$$

Here, y is either the amount bound or is some factor proportional to it (e.g. radioactivity, absorbance etc.), or the number of moles bound per mole of receptor. The capacity for binding ligand, Cap , is either an amount bound in terms of moles of ligand, or a stoichiometric quantity in terms of moles of ligand per mole receptor depending upon the definition of y . The dissociation constant of ligand from the receptor is K_d . The concentration of free ligand, $[L]$, is the x data axis.

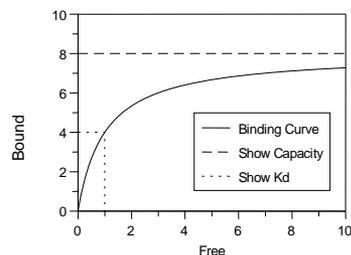

T 1 Site – Scatchard Plot

The Scatchard plot is the usual way to linearize binding data. This display equation takes the Cap and K_d values calculated by fitting to a single-site binding curve, and draws the line for the Scatchard plot. Use this display curve with the Scatchard data transformation.

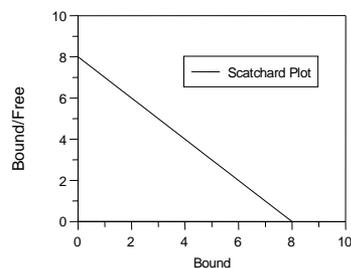

T 1 Site – Show Capacity


This transformation shows the capacity on the plot.

T 1 Site – Show Kd (1 site)

This transformation shows the dissociation constant, K_d . For best results set the curve resolution to **High** when using this transformation.

Σ 1 Site (with background)

Often binding experiments have an inherent background y value. This value can usually be determined separately by control experiments. However, such controls sometimes do not accurately reflect the real experiment, and it can be preferable to treat this background as a separate parameter to be obtained from the data fitting. The equation is:

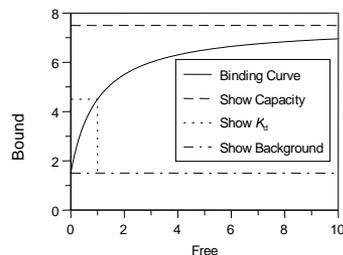

$$y = \frac{[L] \cdot Cap}{K_d + [L]} + \text{background}$$

Here, the parameters are the same as for the single site ligand binding except for the inclusion of a background or offset value.

T 1 Site (with background) - Scatchard

Use this transformation to display Scatchard plots for binding data containing a background value.

T 1 Site (with background) - Show background

This transformation shows the background value.

T 1 Site (with background) - Show Capacity

This transformation shows the capacity.


T 1 Site (with background) - Show Kd

This transformation shows the dissociation constant, K_d . For best results set the curve resolution to **High** when using this transformation.

Σ 2 Site

Often there is more than one type of site involved in binding the ligand. This equation describes the case of two sites characterized by individual K_d and Capacity values. The equation is:

$$y = \frac{[L] \cdot Cap_{(1)}}{K_{d(1)} + [L]} + \frac{[L] \cdot Cap_{(2)}}{K_{d(2)} + [L]}$$

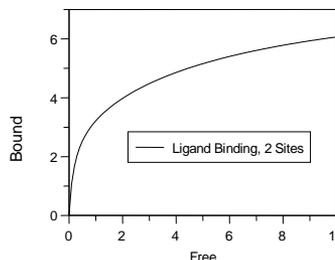

T 2 Site - Show Capacity (total)

This transformation shows the total binding capacity for ligand.

T 2 Site - Show Capacity 1

This transformation shows the capacity of the first binding site.

T 2 Site - Show Curve 1

This transformation shows the binding curve for the first binding site in isolation.

T 2 Site - Show Curve 2

This transformation shows the binding curve for the second binding site in isolation.

Σ 2 Site (with background)

This equation describes the situation where two independent binding sites are present, but the y value at zero ligand concentration is not zero. The equation for this is:

$$y = \frac{[L] \cdot Cap_{(1)}}{K_{d(1)} + [L]} + \frac{[L] \cdot Cap_{(2)}}{K_{d(2)} + [L]} + \text{background}$$

T 2 Site (with background) - Show background

This transformation shows the background binding.

T 2 Site (with background) - Show Capacity 1

This transformation shows the capacity of the first binding site. This value is on top of the background.

T 2 Site (with background) - Show Capacity (total)

This transformation shows the total binding capacity for ligand, on top of the background value.

T 2 Site (with background) - Show Curve 1

This transformation shows the binding curve for the first binding site in isolation.

T 2 Site (with background) - Show Curve 2

This transformation shows the binding curve for the second binding site in isolation.

Σ Adair (2 sites)

The Adair equation describes the binding of ligands to multiple sites on a receptor when the dissociation constants for binding of successive ligands are to be individually defined from the data. The equation for n sites is:

$$B = \frac{Cap \cdot \{L/K_{d1} + 2L^2/(K_{d1}K_{d2}) + \dots + nL^n/(K_{d1}K_{d2} \dots K_{dn})\}}{n \cdot \{1 + L/K_{d1} + L^2/(K_{d1}K_{d2}) + \dots + L^n/(K_{d1}K_{d2} \dots K_{dn})\}}$$

In this equation, B is the amount bound, L the concentration of free ligand, Cap the maximum capacity of the system to bind ligand, and K_{d1} , K_{d2} etc. the successive dissociation constants for ligand.

Σ Adair (3 sites)


Use this equation for 3 binding sites.

Σ Adair (4 sites)

Use this equation for 4 binding sites.

Σ Bound vs Total

The previous equations require the unbound ligand concentration to be known. However, although the total ligand concentration is always known (this is what you have added), the proportion bound and free is not always certain. This is particularly the case if the y axis is an observation that is not directly related to the bound concentration, for example fluorescence intensity or absorbance. Here the equation relates the y axis to the total ligand added, t .

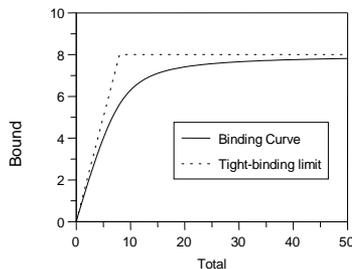

$$y = -\frac{-(K_d + t + Cap) + \sqrt{(K_d + t + Cap)^2 - 4 \cdot t \cdot Cap}}{2}$$

T Bound vs Total - tight binding limit

This transformation shows the maximum possible amount bound (or change in absorbance etc.) for a given total amount of ligand added if the binding were infinitely tight.

Σ Cooperative

A generalized equation describing cooperative binding (or multiple binding sites of differing affinities) is:

$$y = \frac{[L]^n \cdot Cap}{K^n + [L]^n}$$

Here, n is a measure of the co-operativity. For normal binding to a single site, n should not differ significantly from 1.0.

Σ Cooperative (with background)

This is the same as the equation above, but includes an unknown background level in the y axis observations, i.e.

$$y = \frac{[L]^n \cdot Cap}{K^n + [L]^n} + \text{background}$$


pH

These equations are used when the observed value varies with pH. They usually involve the determination of one or more pK_a values from the experimental data.

Σ pKa, (min = 0)

This equation is used when the observed parameter varies with pH from 0 at low pH to a limiting value at high pH. Examples include various enzymatic activities. The equation is:

$$y = \frac{\text{Limit} \cdot 10^{(\text{pH} - \text{pK}_a)}}{10^{(\text{pH} - \text{pK}_a)} + 1}$$

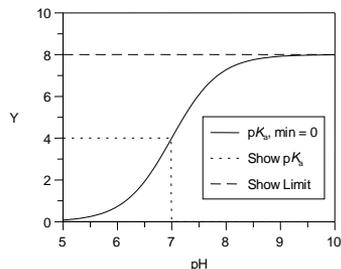

T pKa, (min = 0), Linearized Plot

This transformation displays the line for the linearized titration curve. The linearized equation is:

$$y = -\frac{K_a \cdot y}{[\text{H}^+]} + \text{Limit}$$

T pKa, (min = 0), Show Limit


This transformation displays the limit value of the y axis at high pH.

T pKa, (min = 0), Show pKa

This transformation is used to represent the position of the pK_a on the curve where there is a single ionization and the starting value is 0. For best results set the curve resolution to **High**.

Σ pKa, (min = 0), High to Low

This equation is used when the observed parameter varies with pH from a limiting value at low pH to 0 at high pH.

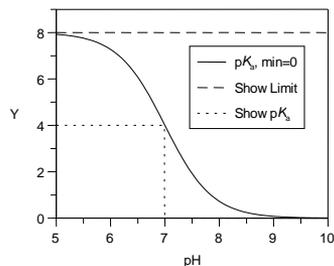

T pKa, (min = 0), High to Low, Show Limit


This transformation displays the limit value of the y axis at low pH.

T pKa, (min = 0), High to Low, Show pKa

This transformation is used to represent the position of the pK_a on the above curve. For best results set the curve resolution to **High**.

Σ pKa, Double

This equation is used when the observed parameter varies from one value at low pH to another at high pH, with the variation depending upon the ionization of two groups. The pK_a values of these groups determine the shape of the curve.

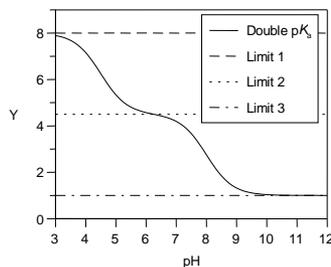

$$y = Limit_1 + \frac{(Limit_2 - Limit_1) \cdot 10^{(pH - pK_{a(1)})}}{10^{(pH - pK_{a(1)})} + 1} + \frac{(Limit_3 - Limit_2) \cdot 10^{(pH - pK_{a(2)})}}{10^{(pH - pK_{a(2)})} + 1}$$

T pKa, Double, Show Curve 1

This transformation shows the first of the two ionization curves in isolation.

T pKa, Double, Show Curve 2

This transformation shows the second of the two ionization curves in isolation.

T pKa, Double, Show Limit 1

This transformation displays the lower limit of the double curve.

T pKa, Double, Show Limit 2


This transformation displays the middle limit of the double curve.

T pKa, Double, Show Limit 3

This transformation displays the upper limit of the double curve.

Σ pKa, Double Bell

Sometimes two ionizations determine the value of the observed parameter. If the first ionization results in increasing activity and the second destroys activity this results in a bell-shaped curve. The equation assumes that activity is zero at low and high

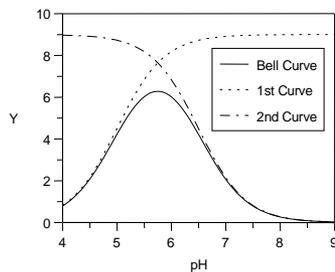

T pKa, Double Bell, Curve 1

This transformation is used to display the effect of the first ionization in the bell-shaped curve.

T pKa, Double Bell, Curve 2

This transformation is used to display the effect of the second ionization in the bell-shaped curve.

T pKa, Double Bell, Limit


This transformation shows the upper limit of the bell-shaped activity profile.

Σ pKa, Single

This equation is used when the observed parameter varies from one value at low pH to another at high pH, with the variation depending upon the ionization of a group. The pK_a value of this group determines the inflection point of the curve. The equation is:

$$y = \text{Limit}_1 + \frac{(\text{Limit}_2 - \text{Limit}_1) \cdot 10^{(\text{pH} - pK_{a(1)})}}{10^{(\text{pH} - pK_{a(1)})} + 1}$$

In this equation, Limit_1 and Limit_2 are the limits of the titration curve at low and high pH.

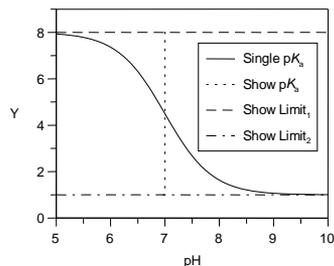

T pKa, Single, Linearized plot

This transformation displays the line for the linearized titration curve.

T pKa, Single, Show Limit 1

This transformation displays the limiting y value at one extreme of pH.

T pKa, Single, Show Limit 2

This transformation displays the limiting y value at the other extreme of pH.

T pKa, Single, Show pKa

This transformation is used to represent the position of the pK_a on the general ionization curve. For best results set the curve resolution to **High**.

Rate


This set of equations fits data describing a first order rate equation.

Σ 1st order

The equation for a process where the observed value increases with time is:


$$A_t = A_\infty(1 - e^{-kt})$$

Here, A_t is the amount of A at time t , A_∞ is the maximum amount of A formed, and k is the rate constant that describes the process.

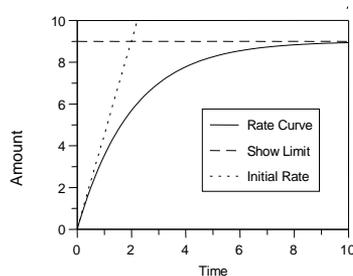

T 1st order, decay curve

This transformation converts the data into an exponential decay curve. To linearize the curve, display the y axis with a logarithmic scale.

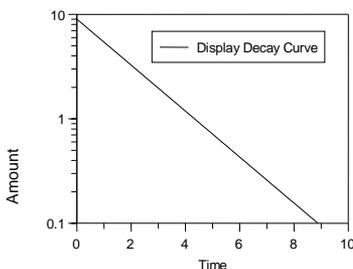

T 1st order, show limit

This transformation displays the limiting value, A_∞ .


T 1st order, show initial

The initial part of the full curve (approximately the first 10%) is linear, with the slope being $k.A_\infty$. This transformation shows the initial linear portion.

Σ 1st order (offset)

If the observed y value has a background or offset at time $t=0$ (due to experimental constraints, for example a starting absorbance may not be zero) the equation becomes:

$$A_t = A_\infty(1 - e^{-kt}) + \text{offset}$$

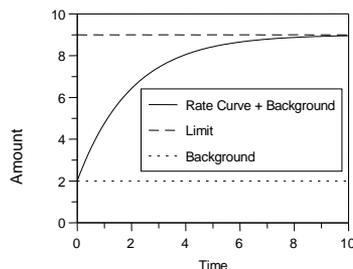

T 1st order (offset), decay curve

This transformation plots the rate equation as an exponential decay. Use it with 'Decay Curve (+Offset)' transformed data.

T 1st order (offset), show limit

This transformation displays the limiting value for the equation containing an offset.


T 1st order (offset), show offset

Use this transformation to show the offset value.

Σ 1st order rate versus [S]

To find enzyme kinetics parameters, it is often necessary to carry out a series of first order rate reactions at different substrate concentrations. The rate constant versus [S] is then analyzed using the Michaelis Menten equation. Using this 3D equation allows these procedures to be combined, giving:

$$A_t = A_\infty \left(1 - \exp \left(- \left(\frac{V_{max} \cdot [S]}{K_m + [S]} \right) \cdot t \right) \right)$$

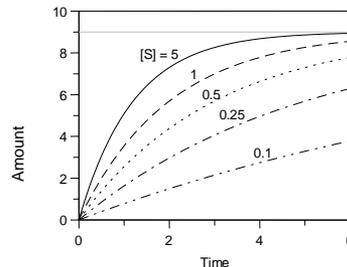

Σ 1st order rate versus [S] (offset)

This equation is the same as that above, but also includes an unknown background value:

$$A_t = A_\infty \left(1 - \exp \left(- \left(\frac{V_{max} \cdot [S]}{K_m + [S]} \right) \cdot t \right) \right) + \text{background}$$

Standard Templates


Data fitting templates provide a quick way of analyzing data and plotting graphs. In general, you are recommended to construct your own specific templates that are best suited to the data you wish to analyze. However, GraFit comes with a few standard templates that can be used for common data fitting situations.

About the Standard Templates

All the standard data fitting templates are defined so that as the data values are entered, the results are calculated and data plotted in the background.

Accessing Templates

To access one of the data fitting templates, select the **File ▶ New** command. Choose the template from the list shown.

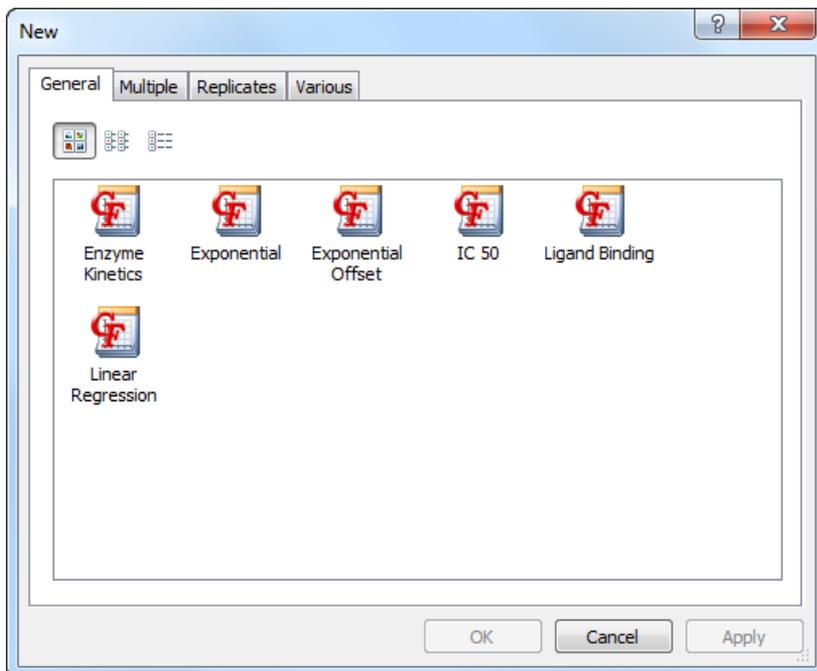

Template Layout

All the standard templates are formatted so that the data are entered into a data table that is on the left of the screen, and the fitted data and results are updated in the graph on the right. This graph is plotted as the data are being typed into the program.

The screenshot displays the GraFit software interface with the following components:

- Menu Bar:** Home, Text, View, Table, Options.
- Toolbar:** Import Data, Add Stock Item, Store Stock Item, Paste, Undo, Redo, Cut, Copy, Clear, Select All, Edit Text, Object.
- Items Panel:** Contains icons for Notes, Data, Graph, Non-linear Fit, Enzyme Kinetics, Lineweaver Burk, and Text Strings.
- Data Table:**

Type	1	2
Column Name	[Substrate]	Rate
1		
2		
3		
4		
5		
6		
7		
8		
9		
10		
11		
12		
13		
14		
15		
16		
17		
18		
19		
20		
- Graph:** Titled "Enzyme Kinetics Data". The main plot shows Rate vs. [Substrate] with axes ranging from 0 to 10. An inset plot shows 1/Rate vs. 1/[Substrate] with axes ranging from 0 to 10. The text "Results not fitted" is displayed at the bottom of the graph.
- Notes Panel:**

Instructions

 1. Enter substrate, rate data into the data table
 2. The data are fitted and the graph plotted automatically

The x axis name is linked to the first data column name; the y axis name is taken from the second data column name. The title of the graph is set from the "Graph Title" entry in the Text Strings data table.

Screen Size

To ensure compatibility with the majority of systems, the standard templates are formatted so that they display correctly when GraFit is run as a maximized window at a screen resolution of 1024 × 768.

Using Your Own Templates

In general, we recommend that you build your own templates for common data fitting situations. The advantages of building your own templates are summarized as follows:

- Specific data analyses can be created.
- Graphs can incorporate user-specific features such as corporate logos or the name of the user.
- Graphs can be customized to use specific fonts and/or styles. Users with color printers may wish to incorporate color into their graphs; those with monochrome printers will probably wish to create graphs in black and white.

- Users with high-resolution displays can make best use of the available screen resolution when laying out the various windows.

See Chapter 4 for more information about working with template files.

Description of the Standard Templates

The standard templates generally come in three related versions:

- **General Tab.** These analyses are appropriate for a single set of x / y data
- **Multiple Tab.** These analyses allow a number of data sets (typically up to 8) to be analyzed together and comparison plots made. The layout of the data in these templates is of the form $\{x, y_1, y_2, y_3, \dots\}$ and so should be used when the x data values are common to a set of different y values.
- **Multiple Separate X Tab.** These analyses allow a number of data sets (typically up to 8) to be analyzed together and comparison plots made. The layout of the data in these templates is of the form $\{x_1, y_1, x_2, y_2, x_3, y_3, \dots\}$ and so should be used where separate x and y values are present for each data set.
- **Replicates Tab.** These analyses are used where there are replicate data sets (up to 8 replicates for each data point). The data are fitted to the mean of the replicates, with error bars drawn to illustrate the spread of the data.

The descriptions that follow apply to the **General** case only, but can be extrapolated to the **Multiple** and **Replicates** cases. The files themselves contain instructions on how to use them and should be consulted for more specific information.


Enzyme Kinetics

This template analyzes enzyme kinetic data to the basic Michaelis-Menten enzyme kinetics equation:

$$v = \frac{V_{max}[S]}{K_m + [S]}$$

A graph is made of the fitted data together with an inset Lineweaver Burk plot, as shown below.

Enzyme Kinetics Data

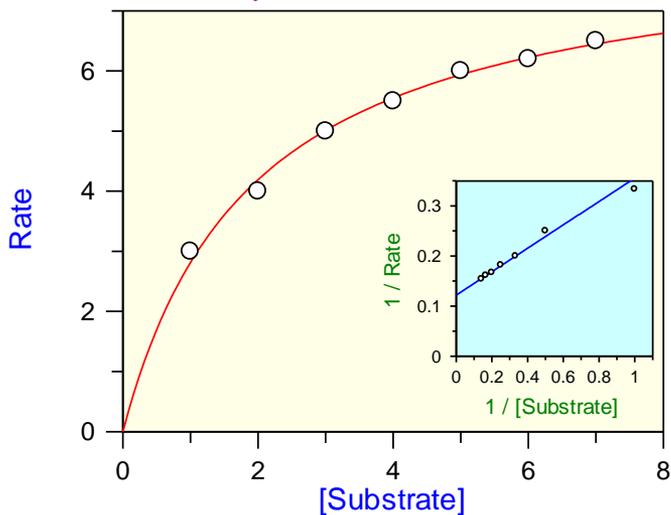

Parameter	Value	Std. Error
Vmax	8.2186	0.2298
Km	1.9235	0.1621


IC 50

This fits data to a full 4-parameter equation. The binding or inhibition is then characterized in terms of an IC_{50} value using the following equation:

$$y = \frac{\text{Range}}{1 + \left(\frac{x}{IC_{50}}\right)^s} + \text{background}$$

A graph is made of the fitted data, as shown below.

IC₅₀ Data

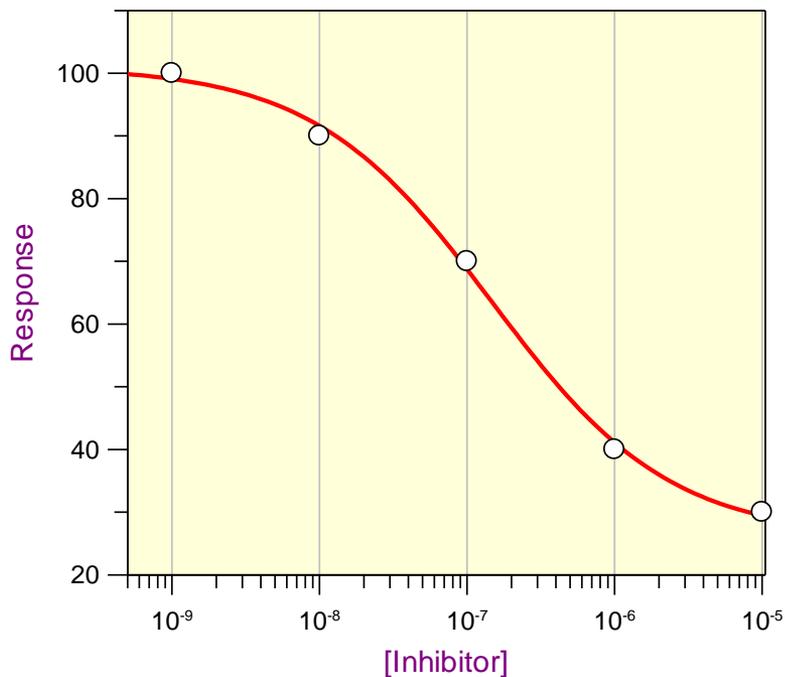

Parameter	Value	Std. Error
Y Range	75.0926	6.9937
IC 50	1.47142e-007	3.99433e-008
Slope factor	0.7194	0.1517
Background	26.0055	4.5106


IC 50 Back Corrected

This fits data to a version of the IC₅₀ equation that pre-supposes that the data have already been background corrected. The binding or inhibition is then characterized in terms of an IC₅₀ value using the following equation:

$$y = \frac{\text{Range}}{1 + \left(\frac{x}{IC_{50}}\right)^s}$$

A graph is made of the fitted data, as shown below.

IC50 Background Corrected Data

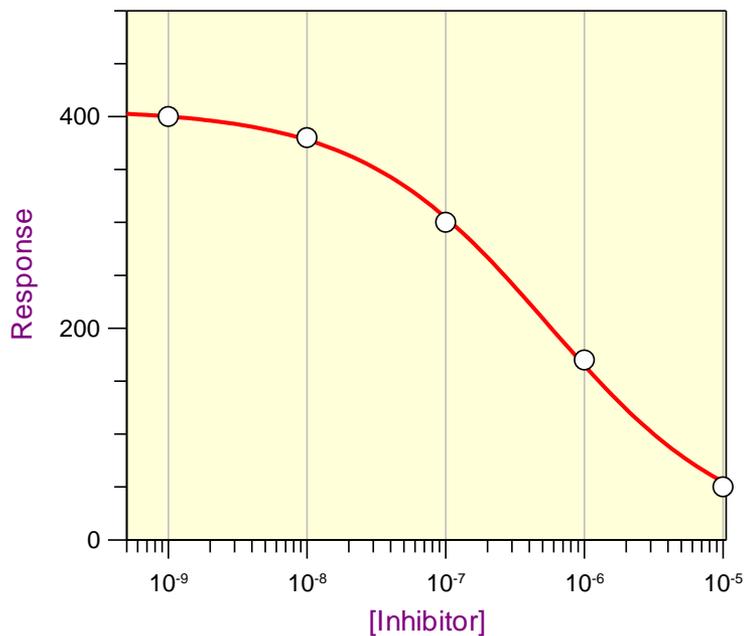

Parameter	Value	Std. Error
Y Range	407.1703	6.5733
IC 50	5.46427e-007	5.51388e-008
Slope factor	0.6397	0.0355


Ligand Binding

This template analyzes ligand binding data to single-site saturation equation:

$$y = \frac{[L] \cdot Cap}{K_d + [L]}$$

A graph is made of the fitted data together with an inset Scatchard plot, as shown below.


Single-site Ligand Binding

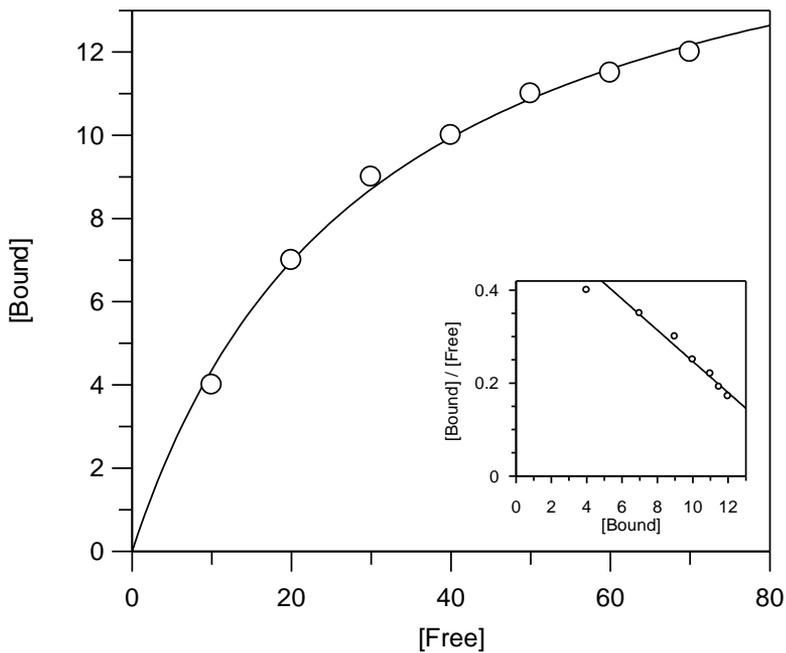

Parameter	Value	Std. Error
Capacity	17.3459	0.6126
Kd value	29.7734	2.5451

Linear Regression

This template allows data to be fitted using linear regression, and the results plotted.

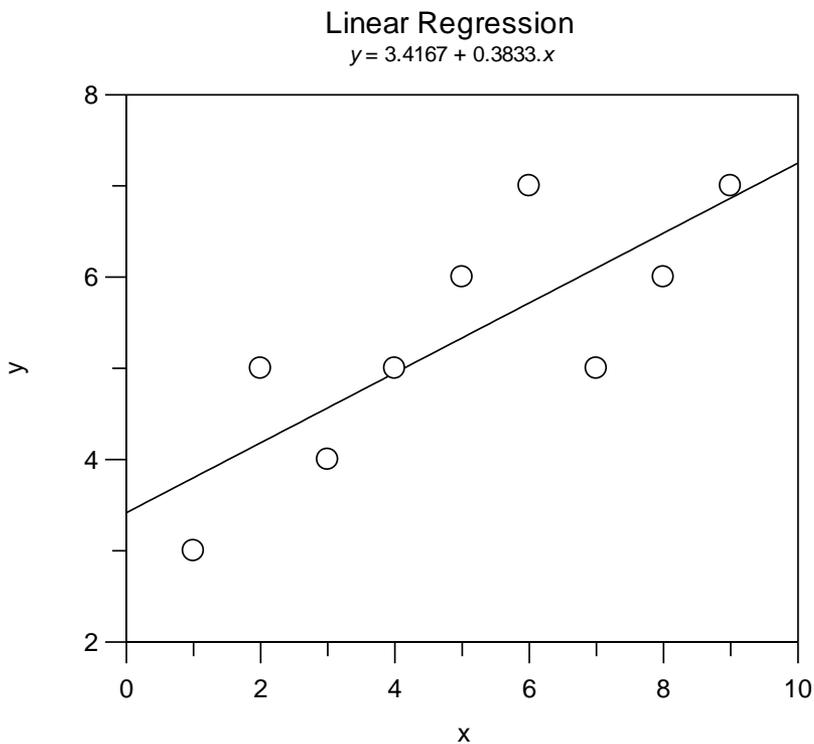

Correlation coefficient: 0.7936

Parameter	Value	Std. Error
a (intercept)	3.4167	0.6251
b (gradient)	0.3833	0.1111


Exponential

This template fits data to a single exponential decay function, using the following equation:

$$y = A_0 \cdot e^{-kt}$$

A graph is created showing the fit.

Exponential Decay

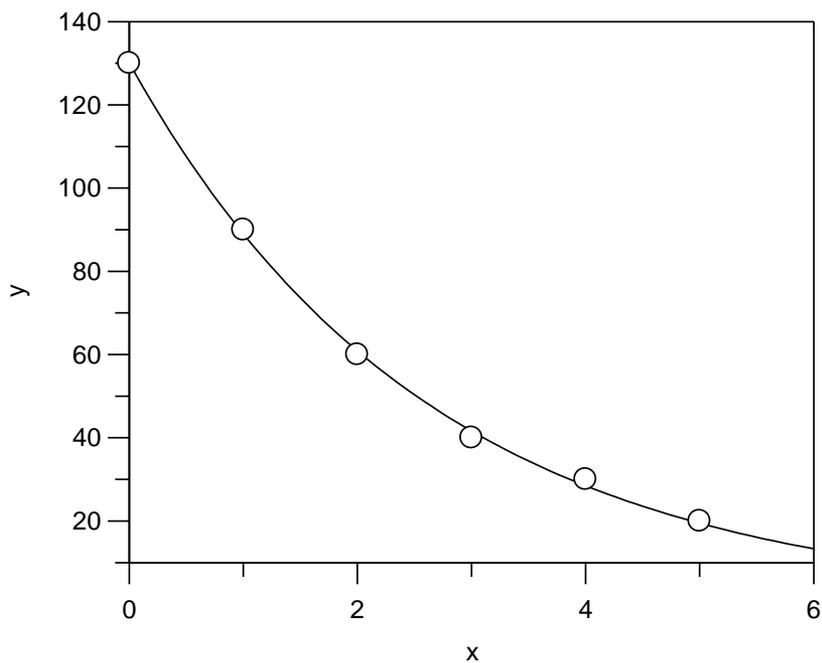


Parameter	Value	Std. Error
Initial value	130.1490	1.2017
Rate constant	0.3799	0.0066

Exponential Offset

This template fits data to a single exponential decay function that goes to a background, non-zero value, using the following equation:

$$y = A_0 \cdot e^{-kt} + \text{background}$$

A graph is created showing the fit.

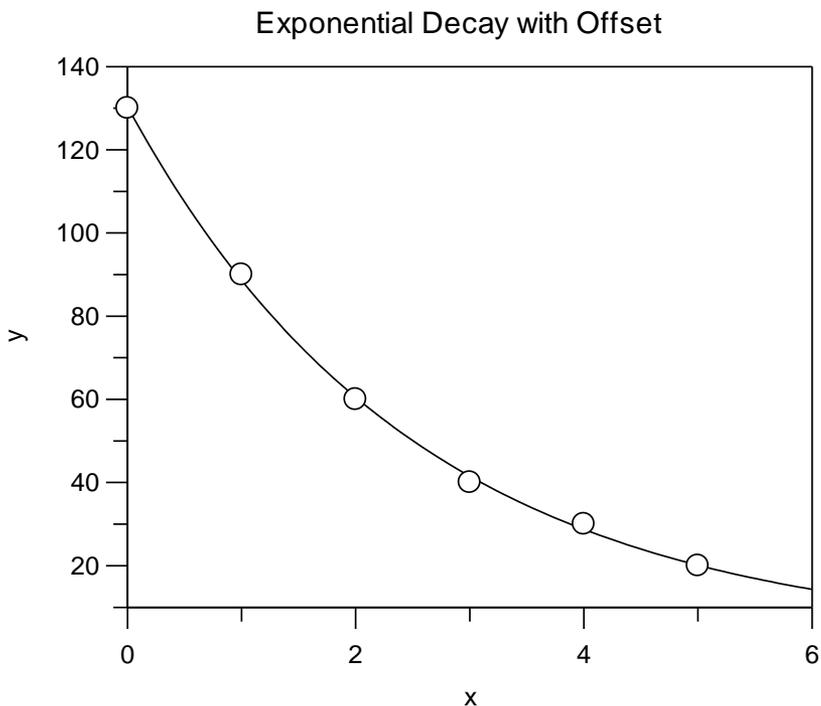

Parameter	Value	Std. Error
Initial value	128.1367	3.3007
Rate constant	0.3947	0.0243
Offset	2.3126	3.4885


Exponential Decay Various Models

This template fits a single data set to a variety of exponential decay models:

- Single Exponential
- Single Exponential with Offset
- Double Exponential
- Double Exponential with Offset

By comparing the fits obtained using the various models, it is possible to assess whether the data justify fitting to models that are more complex than a simple single exponential.

Exponential Decay, various models

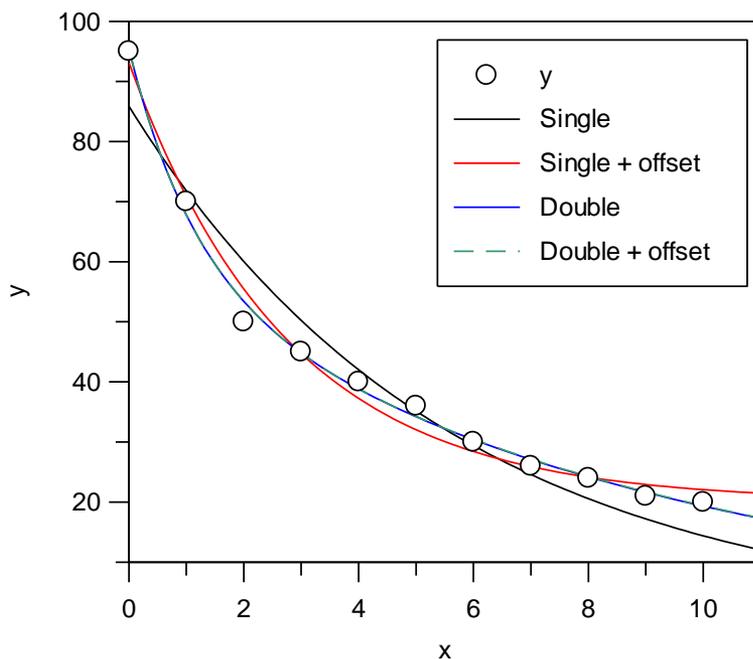
**Ligand Binding Various Models**

This template fits a single data set to a variety of Ligand binding models:

- Single Site
- Single Site with Background
- Two Sites

By comparing the fits obtained using the various models, it is possible to assess whether the data justify fitting to models that are more complex than a simple single site model.

Ligand Binding, Various Models

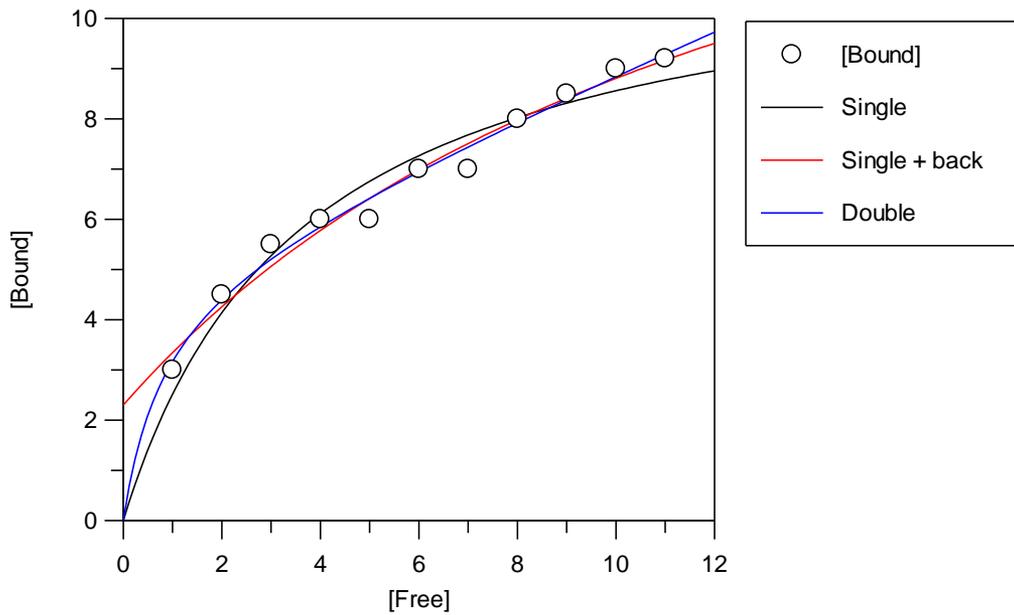

Passing Data between Applications


Beyond Embedding and Linking

In Chapter 18 we saw how to insert information from one application into another. This means that it is possible to incorporate a GraFit analysis into another document, for example into Microsoft Word or Microsoft Excel, in such way that all the information that is needed to create the displayed graph is still present. Double-clicking on the embedded graph image will launch GraFit, allowing the displayed graph to be customized or the data modified.

This form of data transfer works well if all the information used to create the graph is held exclusively within GraFit. However, it does not permit proper two-way transfer of information. Embedded GraFit objects allow more sophisticated data communication using a process known as OLE automation; this means that the GraFit object can be controlled programmatically from within an application that supports this.

Security Warnings When Using OLE Automation

Microsoft Excel 2007 and later is very strict about the use of “active content” in order to prevent files from unknown sources being used for malicious purposes. When you open a data file containing VBA program code, as used by the OLE automation routines, by default Excel will prevent this from running until you have confirmed that it is OK to do so.

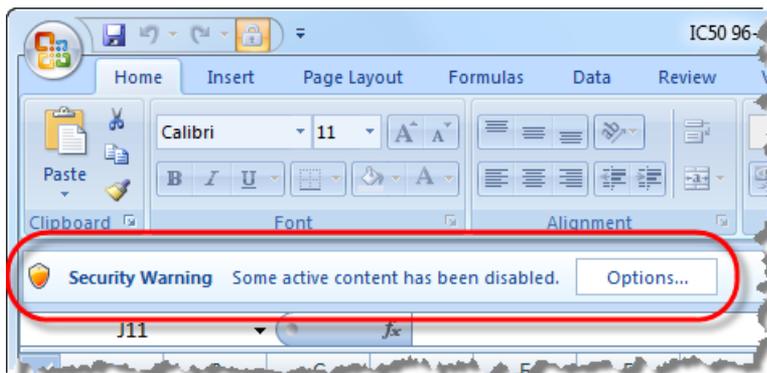

Security warning shown by Excel 2007 after opening a file containing macros

Clicking on the **Options** button results in the following dialog box.



Enabling macros in Excel 2007 and later. Unless you enable the macros it is not possible for Excel and GraFit to communicate

For the automation to function correctly this one time only click the **Enable this content** button, then click **OK**. If you do this then the next time you load the file you will again see the same security warning.

Using Content that has been Digitally Signed

The dialog box shown above indicates that the file in question has been digitally signed, in this instance by Erithacus Software Limited. If you wish to trust all macros signed by Erithacus Software click the **Trust all documents from this publisher** button then click **OK**. If you do this then you will no longer get warnings about files signed by Erithacus Software.


*Note We recommend that you select the **Trust all documents from this publisher** for files signed by Erithacus Software Limited. This will ensure that all our sample files will load without generating further warnings.*

When producing Excel files it is possible to add your own personal digital signature. See <http://office.microsoft.com/en-us/excel/HA012317811033.aspx> for more information.


Example: Passing Data to and from Microsoft Excel

One of the standard GraFit templates is the Enzyme Kinetics template (see Chapter 23). Let us assume that we want to take data from Microsoft Excel and fit these data using this GraFit template.


The easiest way of doing this is to carry out the following procedures, assuming that the data are present in an Excel spreadsheet:

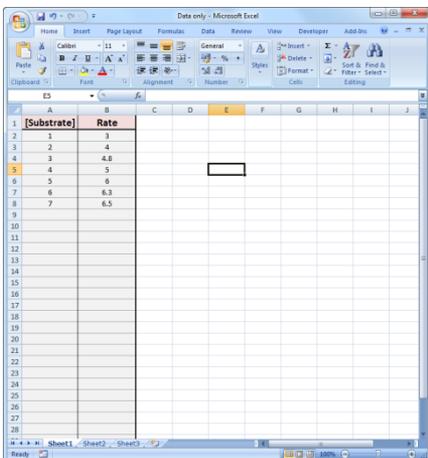

1. The data to be fitted are present in a Microsoft Excel file

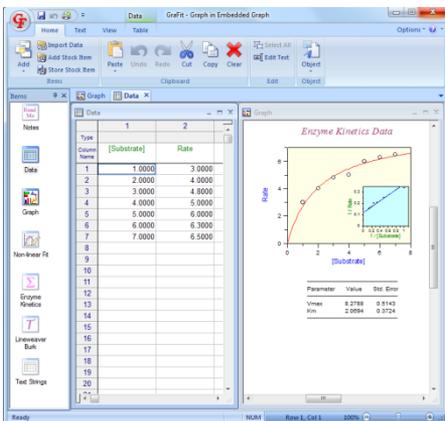

2. These data are copied and pasted into the standard GraFit enzyme kinetics template

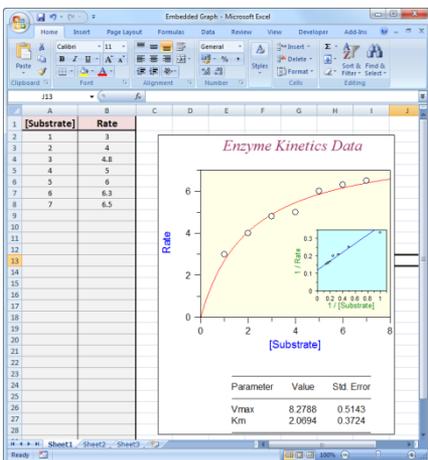

3. The GraFit graph is embedded back into the original Microsoft Excel file.

While this procedure produces the correct fit for the current set of data, it suffers from a number of limitations:


- If the data in the Excel spreadsheet change this is not reflected in the embedded graph. To ensure that the graph reflects the current set of data, these values would need to be copied and re-pasted into GraFit.
- Although the embedded graph displays the fitted parameters, the parameter values are not directly available within the Excel file. If you wished to use these fitted parameters, perhaps to perform further data manipulations within Excel, you would need to open GraFit, select the values and paste them into Excel. Once again, if the data or the fit changed then you would need to repeat this operation to get the updated values.

The following sections describe two approaches to achieving dynamic data transfer between Excel and GraFit.


Transferring Data using the “GraFit Fitter” Excel Template

The Excel file “GraFit Fitter” makes connecting Excel and GraFit a simple process that requires no programming skills. This file has all the “wiring” necessary to drive the data transfer, even allowing complex analysis scenarios with multiple data sets and GraFit objects; all you need to do is define where your data are located, where these data should be transferred to and where you wish the results to be returned.

Using “GraFit Fitter” to Fit the Example File

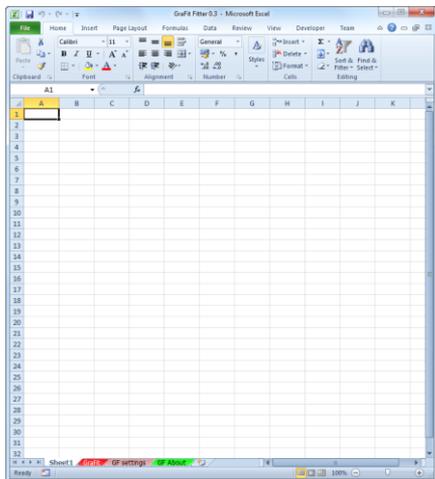


1. Open the “GraFit Fitter” file (this file is pre-installed into your “GraFit Data\Sample Data” directory).

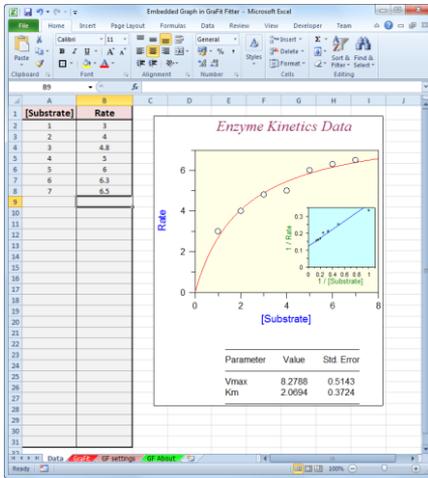


GraFit analyses		Data #1					
Page	Name	Page	Top Left	Bottom Right	Data Item	Start Column	Page
Data	GraFit 1	Data	A1	B31	Data		1

Results #1					
Fit Item	Num fits	Num Param	Page	Top Left	Fit Item
Non-linear	1	2	Data	D30	

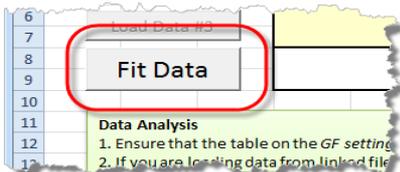


2. Incorporate the data to be fitted into this file, either on the blank initial worksheet page (as here) or by adding additional worksheets. Embed the GraFit graph into the file.
3. Setting up the analysis involves making entries on the "GF settings" page to define the way the data transfer is done
4. On the "GF settings" page, specify that the GraFit data fitting object is called (in this case) "GraFit 1" and is on the "Data" page; that the data to be fitted is also found on this page and occupies the cells A1 to B31. The data is inserted into the GraFit object's "Data" data table, starting at column number 1. Note that the top line of this data region is used to provide the column headings; the remainder provides the data to be fitted.
5. On the "GF settings" page, specify that the results are held in the GraFit item called "Non-linear Fit", that there is a single fit and that the analysis involves 2 parameters. These results will be transferred back to Excel and placed into worksheet page "Data", starting at cell position D30.
6. To fit the data and pass the results back to Excel, turn to the "GraFit" worksheet page and click the "Fit Data" button

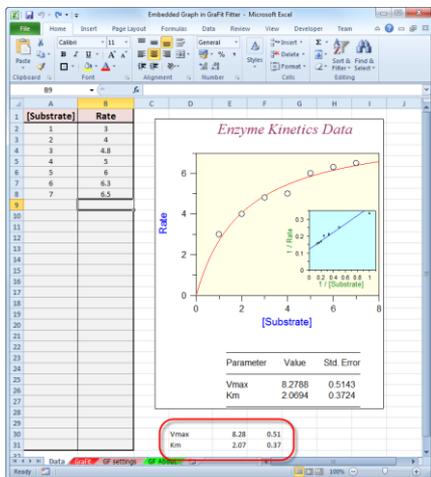

- The data are automatically transferred to the embedded GraFit analysis, fitted, then the fitted values passed back to Excel.

Sheet name	Region for Fit invalidate
Data	A1:B31

- As a final refinement, it is possible to specify that when the underlying data change in Excel, the embedded GraFit analysis should be marked as out of date. To do this, on the "GF settings" page specify one or more data regions where changes are monitored. In this case, on worksheet "Data" the region where the raw data are located is A1:B31.

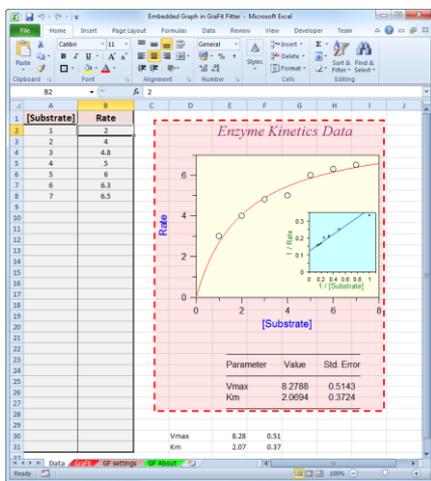

- Any changes to the raw data result in the embedded GraFit analysis being drawn with a dotted red outline and red background to show that it is out of date. To update the analysis turn to the "GraFit" worksheet page and click the "Fit Data" button again.

"GraFit Fitter" Reference

The "GraFit Fitter" file is installed into your "GraFit Data\Sample Data" directory. The file is saved as an Excel Template, which means that double-clicking on the file icon will open the information into a new unnamed file.

The file contains 4 worksheet pages. The first (Sheet 1) is blank and can be used to enter your data for analysis. You can also add any number of additional worksheets as required.



Worksheets in the "GraFit Fitter" file

Setting up Your Data Analysis

The first task is to set up this file so that it holds the data to be analyzed. These data can be in any required format, and may be included on the blank "Sheet 1" worksheet page or on any additional page. However, do not alter the three worksheet pages that are specific to the "GraFit Fitter" file—we will see how these are used shortly.

Ensure that Suitably Arranged Data are Available

You can use Excel to perform any desired rearrangement of data, but it is essential to ensure that in some defined part of the file there exists a region of cells that holds the data that will be passed to GraFit. This needs to be a contiguous region where the data is present in vertical columns, with the top row being used to generate the names of the columns when this is passed to GraFit. The region can be any width and any height.

	A	B
1	[Substrate]	Rate
2	1	3
3	2	4
4	3	4.8
5	4	5
6	5	6
7	6	6.3
8	7	6.5
9		
10		
11		
12		
13		
14		
15		
16		

This region spans cells A1 (top left) to B8 (bottom right)

Data values will be transferred between a designated region on one of the Excel worksheet pages...


	1	2
Type		
Column Name	[Substrate]	Rate
1	1.0000	3.0000
2	2.0000	4.0000
3	3.0000	4.8000
4	4.0000	5.0000
5	5.0000	6.0000
6	6.0000	6.3000
7	7.0000	6.5000
8		
9		
10		

... and a named Data Item within the embedded GraFit analysis file.

Hint If you have a pre-existing Excel file that you use to analyze data, you can copy and paste the relevant worksheet pages into the "GraFit Fitter" file. If your file does not already hold the data in an appropriate column format as shown above, use the ability of Excel to produce a tabulated set of information on a separate worksheet page that is just used to hold the intermediate values. See the example files for a case where this has been done.

Embed One or More GraFit Files into the Excel File

After creating a GraFit file that performs the data analysis you require, this file is embedded within the GraFit Fitter Excel file. To do this, set up your GraFit file to perform the data fitting, select all items in one of the graph window and use the **Copy All** command. Paste the graph into Excel, either on "Sheet 1" or any other worksheet page. Once pasted, select the graph and give it a name by which to can refer to it within the Excel file

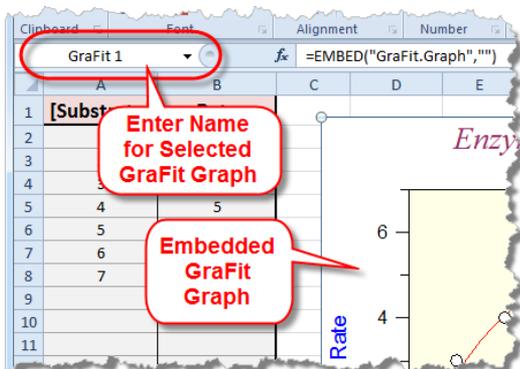

Naming an embedding GraFit item within the Excel file

It is possible to embed multiple GraFit analyses within a single Excel file if required.

Note During data fitting, when passing data to an embedded GraFit file there is a short delay as each file is opened. This means it can be inefficient to use very large numbers of embedded files if this is not necessary. For example, data from a 96-well plate assay might be analyzed by 96 different embedded files, or by 12 embedded files that each processes 8 sets of data. The latter would run far more quickly.

GF About Page

This page records the version number of the “GraFit Fitter” file. It is present for housekeeping purposes only.

GF Settings Page

The entries on this page are used to control the data transfer between GraFit and Excel.

GraFit Analyses Section

Fill in the GraFit analyses section so that there is one entry for each embedded GraFit file. These embedded files can be on any worksheet page, with the page name and the name of the GraFit graph being entered into the table

Page	Name	Page
Results	GraFit 1	Wc
Results	GraFit 2	Wc
Results	GraFit 3	Wc
Results	GraFit 4	Wc
Results	GraFit 5	Wc
Results	GraFit 6	Wc
Results	GraFit 7	Wc
Results	GraFit 8	Wc

In this example are 8 separate analyses

A blank entry marks the end of the analysis list

Setting up the GraFit analyses table on the GF Settings page. The number of analyses performed depends on how many entries are present in the analysis list

Data Section

For each of the embedded GraFit graphs it is possible to transfer one or two regions of data from Excel into a named Data Item within GraFit. The location of the source data within Excel is specified in terms of the worksheet page and region of cells that hold the information: the coordinate of the top left and bottom right cells needs to be specified. The data are transferred to GraFit and placed in the Data Item that is specified; the data are transferred starting with the designated Start Column. For example, if the Start Column is 1 then the data are copied into the first column of the Data Item. The top row of data within Excel has special meaning in that these data are used for the column names within GraFit (see above).

Data #1					
Source (Excel)			Destination (GraFit)		
Page	Top Left	Bottom Right	Data Item	Start Column	Page
Working	A2	I12	Data		1
Working	A15	I25	Data		1
Working	A28	I38	Data		1
Working	A41	I51	Data		1
Working	A54	I64	Data		1
Working	A67	I77	Data		1
Working	A80	I90	Data		1
Working	A93	I103	Data		1

Specifying the Data locations for the first data region on the GF settings page

In many cases it will only be necessary to transfer a single region of data, using the “Data #1” part of the table. If required, however, a second data region can be specified by completing the “Data #2” entries; should these not be required then leave the entries blank.

Results Section

When the data are copied into the GraFit file, GraFit will automatically carry out the data fitting procedures that are built into the file. Once these calculations are complete, the results can then be transferred back to Excel. The Results section allows one to specify the name of the data fitting item, the number of fitting results (1 or more) and the number of parameters that are involved. These are then copied back into Excel on the worksheet page that is entered, at a location that starts at the cell specified (“Top Left”).

Results #1					
Source (GraFit)			Destination (Excel)		
Fit Item	Num fits	Num Param	Page	Top Left	Page
Fit Full	8	4	Output	A2	Fit
Fit Full	8	4	Output	A11	Fit
Fit Full	8	4	Output	A20	Fit
Fit Full	8	4	Output	A29	Fit
Fit Full	8	4	Output	A38	Fit
Fit Full	8	4	Output	A47	Fit
Fit Full	8	4	Output	A56	Fit
Fit Full	8	4	Output	A65	Fit

Specifying the Results locations for the first results region on the GF settings page

The results are copied using a separate row for each parameter, with columns containing alternating fitted parameter values and fitted standard error values for each fitting result in term. So for example, if there are 3 fitting results and 2 different parameters the following values will be copied back to Excel:


Param name 1	Param 1, fit 1	Error 1, fit 1	Param 1, fit 2	Error 1, fit 2	Param 1, fit 3	Error 1, fit 3
Param name 2	Param 2, fit 1	Error 2, fit 1	Param 2, fit 2	Error 2, fit 2	Param 2, fit 3	Error 2, fit 3

The first column lists the names of the parameter values.

In many cases it will only be necessary to transfer a single set of results, using the “Results #1” part of the table. If required, however, up to three sets of results can be copied by completing the “Results #2” and “Results #3” entries; should these not be required then leave the entries blank.

Error Checking

The final column of the table provides some error checking to confirm that the entries are valid. Ensure that all entries read “OK”.

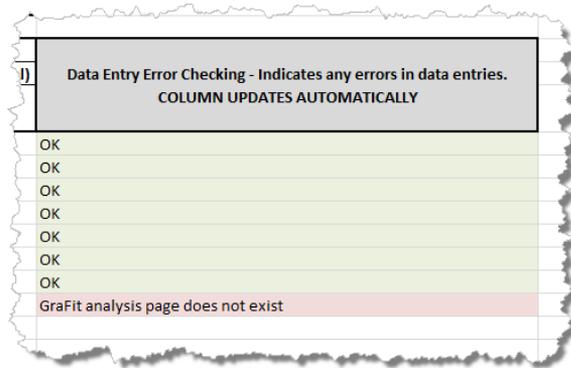


Specifying the Results locations for the first results region on the GF settings page

Ensuring that the Fit is marked as Invalid when Data are Altered

The final set of entries on the GF settings page is used to ensure that if the raw data within the Excel file changes, perhaps because new data values have been pasted in, the GraFit fits are flagged to show that they are out of date and so need refitting.

Up to 4 regions within Excel can be watched for changes; specify the name of the worksheet and the list of cells to be monitored. The region(s) specified can include regions (for example C4:F27) and combine regions and individual cells with commas (for example B3, D3:F12, H1:H24)

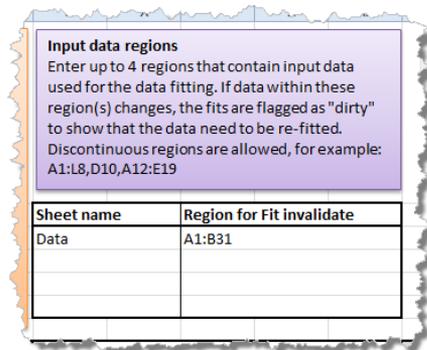

Specifying which regions of the Excel file should be monitored to see when the GraFit fit should be flagged as “dirty”

GraFit Page


The GraFit worksheet page is used to load any linked data fits and also to begin the data fitting.

	Data Page	Data file location
Load Data #1	Plate 1 Raw	C:\Users\Robin\Desktop\Test Data\Plate 1.csv
Load Data #2	Plate 2 Raw	C:\Users\Robin\Desktop\Test Data\Plate 2.csv
Load Data #3		(Not loaded)
Fit Data		(Not fitted)

The GraFit worksheet page contains the **Fit Data** button and also allows loading linked data

Linked Data Files

Data can be loaded into the Excel file in a number of ways. It can be typed in manually and it can be pasted from other sources. However, a very convenient way of including data is to use the **Data ► Get External Data From Text** command.

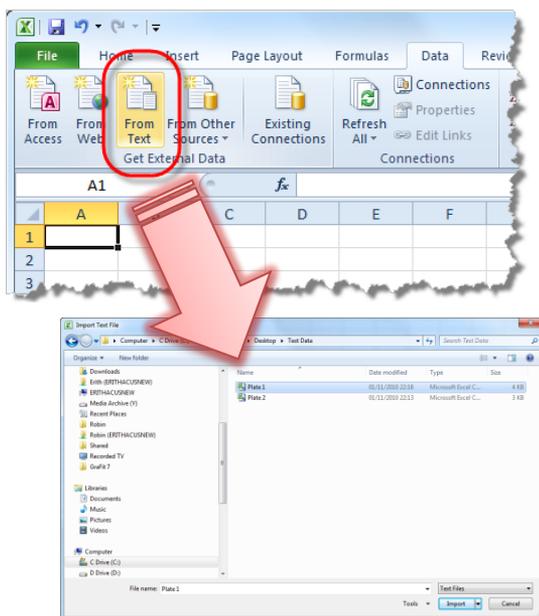

Using the **Data ► Get External Data From Text** command (Excel 2010 shown though equivalent commands are available in Excel 2007 and 2003).

The Excel text import Wizard allows the data to be inserted into a worksheet page in a very flexible manner (see the Excel documentation for more information).

GraFit can use this mechanism to load up to three user-specified files into the GraFit Fitter file. To do this:


1. Insert a new blank worksheet page for each of the datasets you wish to link to.
2. Create the file links using the **Data ► Get External Data From Text** command.
3. Format the worksheet in any way you wish.
4. For each linked data file, enter the name of the worksheet on which the data are loaded into the Data Page cell on the GraFit page. In the example above there are two linked data files, on pages called "Plate 1 Raw" and "Plate 2 Raw".

- Use the **Load Data #n** button for each of the linked data files in turn to load their contents.

Hint Data from laboratory instruments such as plate readers is routinely saved in a form that allows it to be read in this manner.

Fitting the Data

Once the GF Settings page has been fully set up and the data loaded, click the **Fit Data** button to start the data fitting.

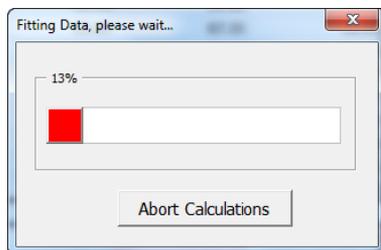


*During data fitting a progress dialog box is shown. Click **Abort Calculations** if you wish to cancel the data fitting.*

Transferring Data using Visual Basic (Advanced!)

It is possible to gain more complete control over the data transfer process by programmatically defining how data is passed between Excel and GraFit. This is done within Excel using Visual Basic for Applications (VBA). VBA is a programming language that allows complete control over the way that information is processed.

Note This is for advanced users only. A full description of how VBA works is beyond the scope of this manual and you should consult reference material on VBA if you are unfamiliar with how to use this application programming environment.

During installation of GraFit, a sample Excel file called *Inset data Kinetics* is placed in the GraFit *Sample Data* directory. This file contains the same data and embedded graph as above, but this file contains some VBA programming that links together the GraFit graph and the raw data. In addition, the fitted parameter values are brought back into Excel and presented as a table.

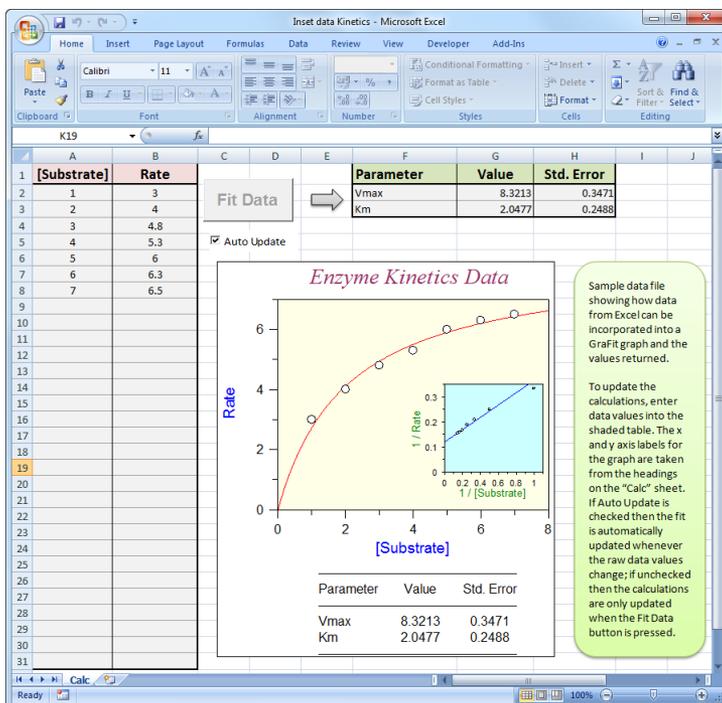


The sample Excel data file "Inset data Kinetics" shows how it is possible to link together data in Excel and an embedded GraFit graph


This file is set up so that whenever the data in the Excel sheet changes, the embedded GraFit graph also updates in real time. For large data sets, or when using a less powerful computer, live update might be a little slow. The example file therefore contains a check box, *Auto Update*, that controls whether updates are made continuously. If *Auto Update* is not checked, the embedded GraFit graph is only updated when the **Fit Data** button is clicked.

Getting Data Indirectly


In the example above, data values are typed directly into the data table and then act as the source of data used by GraFit. More generally, however, you will want to use data that has been manipulated within Excel. For example, the data might come from another program or datasheet entirely, or be processed within Excel to perform averaging or background subtraction etc. The sample file *Linked data IC50* shows how it is possible to take data from a separate Excel worksheet page and perform the equivalent analysis to that of the GraFit template file *IC 50*.

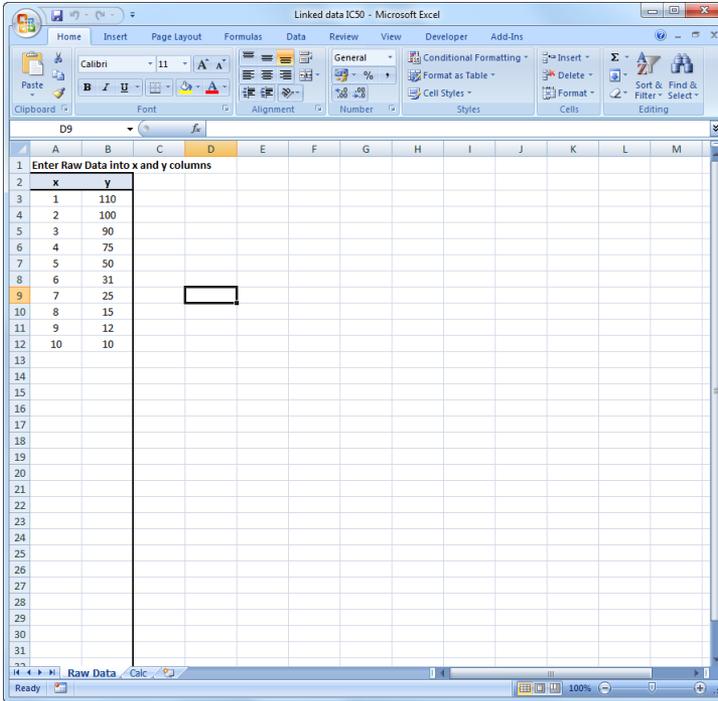

Data are present within the "Raw Data" sheet of the file

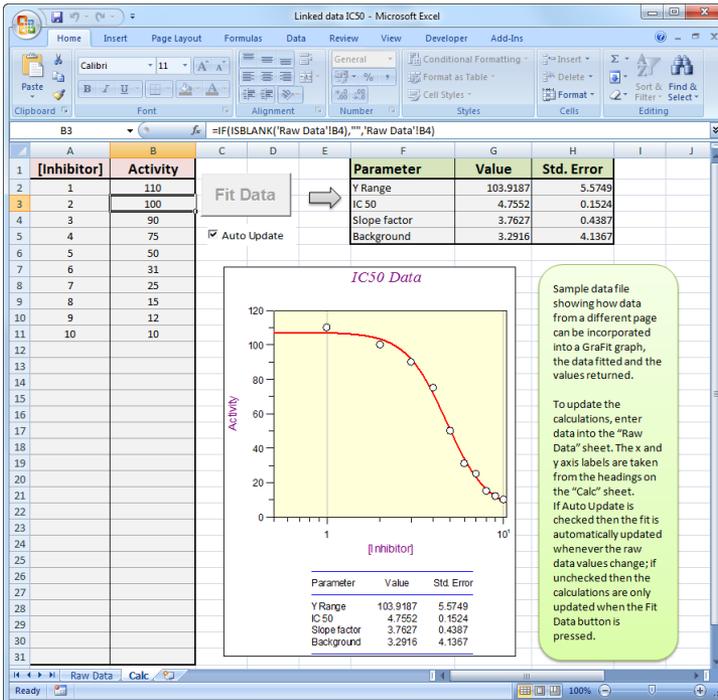

The values in the Raw Data worksheet are transferred to the Calc worksheet for analysis

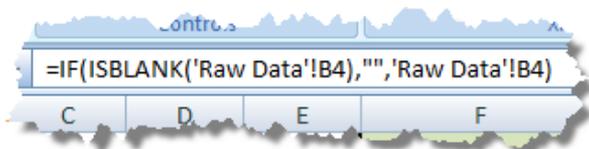

Detailed view of how the data are linked together between the worksheets

How is Data Transfer Done?

The data transfer is achieved using VBA. To illustrate, the main data fitting subroutine, *FitData()*, in the sample above is reproduced below:

```

Linked data IC50.xlsm - Module1 (Code)
(General) FitData
Public Const arr_size = 30
Public Const output_range = "F2:H5"
Public mydata() As Variant
Public Const CalcSheet$ = "Calc"
Public Const GFObjct$ = "GraFit 1"
Public Const GFData$ = "Data"
Public Const GFFit$ = "Non-linear Fit"

Public Sub FitData()

Rem Get the data from the worksheet
Dim a() As Variant
Dim b() As Variant

a = Worksheets(CalcSheet$).Range("A2:A" + Format(1 + arr_size)).Value
b = Worksheets(CalcSheet$).Range("B2:B" + Format(1 + arr_size)).Value

Rem We need to put the data into a 1-dimensional array
Dim c(1 To arr_size) As Variant
Dim d(1 To arr_size) As Variant

For i = 1 To arr_size
    Let c(i) = b(i, 1)
    Let d(i) = a(i, 1)
Next

Rem Pass the data to GraFit
Call GF_SetColumnData(CalcSheet$, GFObjct$, GFData$, 1, d)
Call GF_SetColumnData(CalcSheet$, GFObjct$, GFData$, 2, c)

Rem Update the names of the columns within GraFit
Call GF_SetColumnName(CalcSheet$, GFObjct$, GFData$, 1, _
Worksheets(CalcSheet$).Range("A1").Value)
Call GF_SetColumnName(CalcSheet$, GFObjct$, GFData$, 2, _
Worksheets(CalcSheet$).Range("B1").Value)

Rem Wait until calculation done
Call GF_WaitForCalculations(CalcSheet$, GFObjct$)

End Sub

```

The *FitData()* subroutine from the *Linked data IC50* sample file

This subroutine gathers together data from Excel and then passes the information to GraFit using a set of *GF_XXXXXX* subroutine calls. These subroutines are helper functions that make using GraFit automation easier and are collectively found in the sample data file in *GF_module*. The VBA program code for these is present in the *Sample Data* directory, in a file called *GF_module.bas*.

Note If you are planning use the automation features, we recommend that you edit the sample Excel files to see how this has been implemented.

Sample Data Files Installed


When GraFit is installed, some sample Excel data files are copied into your GraFit *Sample Data* directory. These files illustrate some of the ways in which you can transfer and analyze data.

Note The sample Excel files were created with Microsoft Excel 2010 but saved in Excel 97-2003 compatible format to be compatible with Excel 2003 or later.

Excel File	Description
Inset Data Kinetics	This file is described above
Kinetics (GraFit Fitter)	This file performs the equivalent fitting task to the “Inset Data Kinetics” file, but is built using the GraFit Fitter template
Linked Data IC50	This file is described above
IC50 (GraFit Fitter)	This file performs the equivalent fitting task to the “Linked Data IC50” file, but is built using the GraFit Fitter template

Excel File	Description
------------	-------------

IC50 96-well plate	This file shows how it is possible to carry out more extensive rearrangement data, including reorganizing data from a standard 96-well format.
--------------------	--

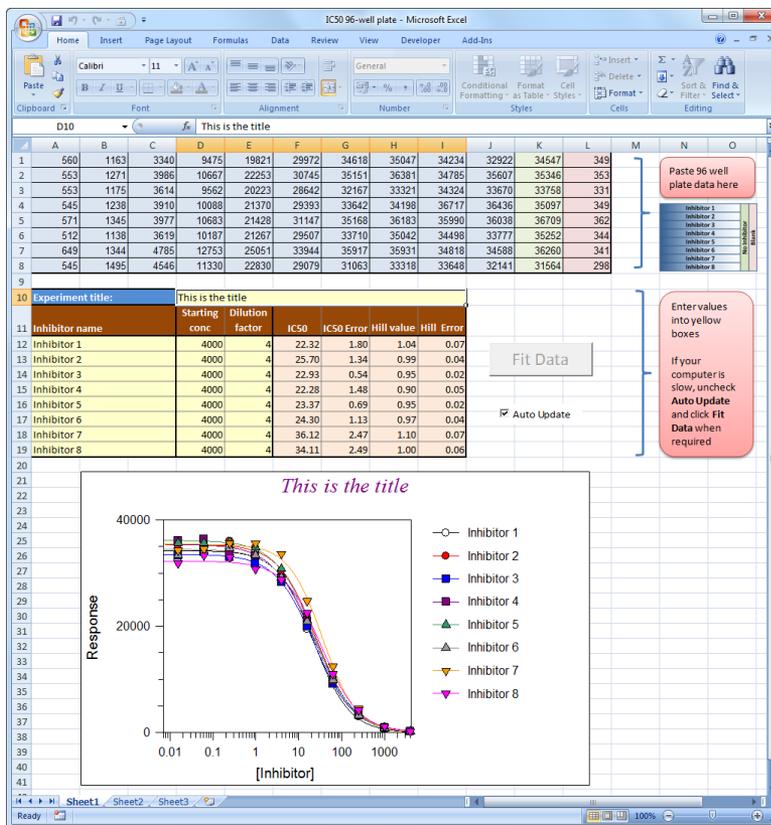


Note The online tutorials include a detailed description of how this file was generated.

IC50 96-well plate (GraFit Fitter)	This file performs the equivalent fitting task to the "IC50 96-well plate" file, but is built using the GraFit Fitter template
------------------------------------	--

IC50 96-well plate (GraFit Fitter) - Linked	This is a version of the "IC50 96-well plate (GraFit Fitter)" file that gets its data from a linked file.
---	---

Unprotecting the Files

The sample files (though not ones built from the GraFit Fitter template) are saved in a way that prevents entry of data in inappropriate parts of the spreadsheet but allows correct data entry—the relevant sheets are Protected. If you wish to edit these files it is first necessary to unlock them, which can be done using the **Format ► Unprotect Sheet** command.

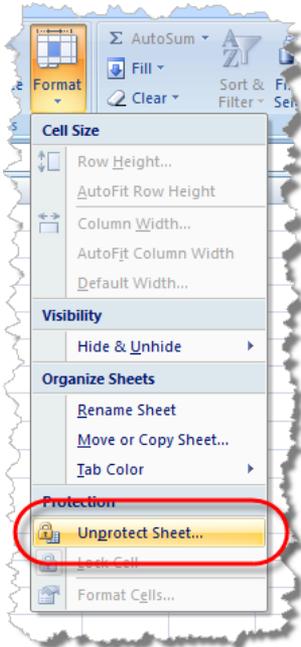

Use the **Format ► Unprotect Sheet** command to unlock the sample Excel files and allow them to be customized.

In Excel 2003 the equivalent command is **Tools ► Protection ► Unprotect Sheet**.

GraFit Automation Reference (Advanced!)



Using OLE Automation to control GraFit involves making calls to GraFit functions and subroutines that are exposed to programs that are automation-enabled. Not all programs are able to carry out automation and the examples shown in this manual use Microsoft Excel.

Using GraFit Automation with Microsoft Excel

To ensure that you are able to access the GraFit programming interface it is necessary to include the GraFit Type Library, *GraFit.tlb*, which is copied to the GraFit program directory during installation of the program. Use of this requires the following preliminary steps:

1. Ensure that the Developer Tab is Visible (Excel 2007)

Excel version 2007 does not by default show the programming tools ("Developer tab"). To reveal this, first select Excel Options from the Office button.

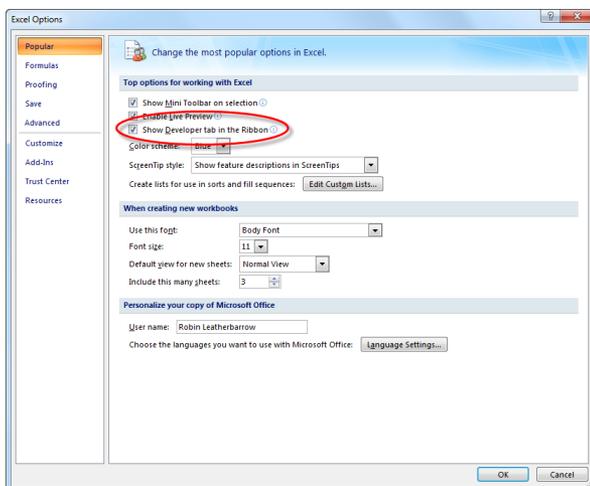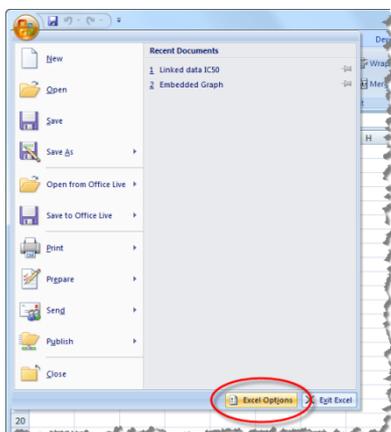


First, select Excel Options...

...then ensure Show Developer tab is checked

2. Launch Visual Basic


In Excel 2007, select the Developer tab and click **Launch Visual Basic**.

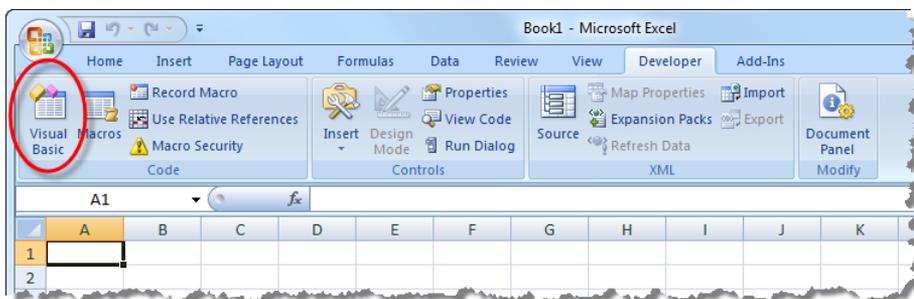

In Excel 2003 use the **Tools ► Macro ► Visual Basic Editor** command.

3. Ensure the GraFit Type Library is Selected

In the Visual Basic Editor, select the **Tools ► References** command to open the References dialog box.

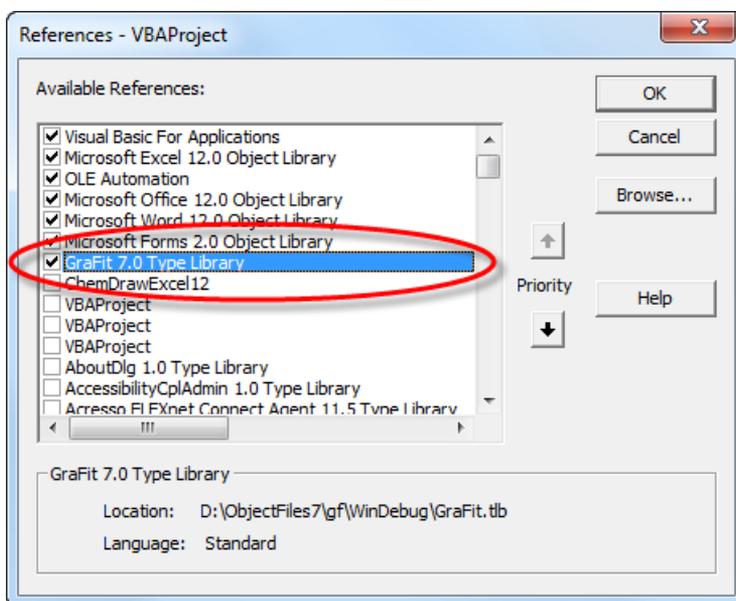

Ensure that the GraFit Type Library is loaded and selected if you wish to access the GraFit automation routines

The GraFit Type Library should have been loaded during installation. If it is not present, click the **Browse...** button to select the *GraFit.tlb* file that is present in the GraFit program files directory. To use the GraFit automation routines it is necessary to make sure that this type library is checked as shown above. It is necessary to select type libraries for each new VBA project that you create.

GraFit Automation Routines

The following classes are defined by GraFit. The main interface is the **Graph** interface, which represents the embedded graph object itself. To access the underlying fits and data, the **DataTable** interface holds members controlling access to a particular data table, the **DataTables** interface represents a collection of all the data tables in the GraFit object. The **FitResult** interface holds members controlling access to a particular data fitting object, the **FitResults** interface represents a collection of all the data fitting items within the GraFit object.

Class GraFit.Graph

GraFit.Graph Member	Description
Sub AbortAllCalculations ()	Abort all current fitting calculations
Function GetDataTables () As Object	Get the interface to the Data Tables
Function GetFitResults () As Object	Get the interface to the Fitting Results
Function GetVersionBuild () As Long	Get the Version Build number (o.o.o.#)
Function GetVersionMajor () As Long	Get the Major Version number (#.o.o.o)
Function GetVersionMinor () As Long	Get the Minor Version number (o.#.o.o)

GraFit.Graph Member	Description
Function GetVersionMod() As Long	Get the Version modification (o.o.#.o)
Function IsBusy() As Boolean	Test whether GraFit is currently calculating results

Class GraFit.DataTables

GraFit.DataTables Member	Description
Function GetCount() As Long	Get the number of data table items
Function GetTableIndex (index As Long) As Object	Get data table item by entry number
Function GetTableName (pStr As String) As Object	Get data table item by item window name

Class GraFit.DataTable

GraFit.DataTable Member	Description
Function SetColumn (col As Long, _array) As Boolean	Set the data values for a whole column from an array of data
Function SetColumnAssociatedData (column As Long, value1 As Double, value2 As Double, value3 As Double, value4 As Double) As Boolean	Set values for the associated data for a column (x[2], x[3], x[4], x[5])
Function SetColumnName (column As Long, name As String) As Boolean	Set the column name
Function SetData (col As Long, row As Long, data) As Boolean	Set any data value within the table
Function SetDataDouble (col As Long, row As Long, value As Double) As Boolean	Set a numeric data value within the table

Class GraFit.FitResults

GraFit.FitResults Member	Description
Function GetCount() As Long	Get the number of fitting item items (linear + nonlinear)
Function GetResultIndex (index As Long) As Object	Get fitting item by entry number
Function GetResultName (pStr As String) As Object	Get fitting item by item window name

Class GraFit.FitResult

GraFit.FitResult Member	Description
Function ChiSq (fitNo As Long) As Double	Get the chisq value for fit[fit number]
Function Constant (index As Long) As Variant	Get the value for constant[index]
Function ConstantName (index As Long) As String	Get the name for constant[index]
Function CorrelationCoeff (fitNo As Long) As Double	Get the correlation coefficient for fit[fit number]
Function DataSource () As String	Get the associated data table name
Function EquationName () As String	Get the associated equation name
Function GetDataFormat () As Long	Get the data format
Function GetErrColumnName () As String	Get the name of column holding the error values
Function GetFitStatus () As Long	Get the fitting status
Function GetWeighting () As String	Get the weighting type
Function GetXColumnName () As String	Get the name of column holding the x data values
Function GetYCalc (fitNo As Long, x As Double)	Get the calculated y value, for fit[fit number], [x value]
Function GetYColumnName () As String	Get the name of column holding the y data values
Function IsLinearFit () As Long	Is this a linear fit?
Function IsRobust () As Long	Is robust weighting used?
Function LimitXRange () As Long	Is there a limited x range?
Function LinearFitOrder () As Long	Get the order of the linear fit
Function NumConstants () As Long	Get the number of constants
Function NumData (fitNo As Long) As Long	Get the number of data points for fit[fit number]
Function NumFit () As Long	Get the number of fits
Function NumParameters () As Long	Get the number of parameters
Function Parameter (index As Long, fitNo As Long) As Variant	Get the value for parameter[index, fit Number]
Function ParamName (index As Long) As String	Get the name for parameter[index]

GraFit.FitResult Member	Description
Function RedChiSq (fitNo As Long) As Double	Get the reduced chisq value for fit[fit number]
Function StdError (index As Long, fitNo As Long) As Variant	Get the standard error value for parameter[index, fit Number]
Function XRangeLimitFrom () As Double	Get the lower limit of the x range
Function XRangeLimitTo () As Double	Get the upper limit of the x range

GraFit Helper Routines

The file *GF_module.bas*, found in the GraFit *Sample Data* directory, has a number of useful subroutines and functions that can be used when using automation. Add this file to your VBA project **File ► Import File...** command from the Visual Basic Editor. The following routines are present in the file:

```
Public Function GF_Parameter(ByRef strCalcSheet As String, _
    ByRef strGFOObject As String, _
    ByRef strGFFit As String, _
    ByVal i As Integer, _
    ByVal fitNo As Integer) As Variant
'-----
' Returns the value of a fitted parameter
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOObject      Name of the GraFit object
' strGFFit          Name of the fitting item within the GraFit object
' i                 The parameter number to return (start at 1)
' fitNo             The index (start at 1) of the fit in the object
'-----
Public Function GF_ParameterName(ByRef strCalcSheet As String, _
    ByRef strGFOObject As String, _
    ByRef strGFFit As String, _
    ByVal i As Integer) As String
'-----
' Returns the name of a parameter
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOObject      Name of the GraFit object
' strGFFit          Name of the fitting item within the GraFit object
' i                 The parameter number to return (start at 1)
'-----
Public Function GF_StandardError(ByRef strCalcSheet As String, _
    ByRef strGFOObject As String, _
    ByRef strGFFit As String, _
    ByVal i As Integer, ByVal fitNo As Integer) As Variant
'-----
' Returns the standard error value of a fitted parameter
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOObject      Name of the GraFit object
' strGFFit          Name of the fitting item within the GraFit object
' i                 The parameter number to return (start at 1)
' fitNo             The index (start at 1) of the fit in the object
'-----
```

```

Public Sub GF_SetColumnData(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String, _
    ByRef strGFData As String, _
    ByVal iColumn As Integer, _
    ByRef data() As Variant)
    '-----
    ' Fills a column of data in an embedded GraFit object
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    ' strGFData         Name of the data table item within the GraFit object
    ' iColumn           The number of the data column (start at 1)
    ' data()            The array of data (start at 1) to be inserted
    '-----

```

```

Public Sub GF_WaitForCalculations(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String)
    '-----
    ' Waits until embedded GraFit object has finished calculating
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    '-----

```

```

Public Sub GF_SetColumnName(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String, _
    ByRef strGFData As String, _
    ByVal iColumn As Integer, _
    ByRef strColumnName As String)
    '-----
    ' Sets the name of a data column in an embedded GraFit object
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    ' strGFData         Name of the data table item within the GraFit object
    ' iColumn           The number of the data column (start at 1)
    ' strColumnName     The new name of the column
    '-----

```

```

Sub GF_SetObjectDirty(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String)
    '-----
    ' Show an embedded GraFit object to indicate it is "dirty"
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    '-----

```

```

Sub GF_SetObjectClean(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String)
    '-----
    ' Show an embedded GraFit object to indicate it is "clean"
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    '-----

```

```

Public Sub GF_SetColumnAssociatedData(ByRef strCalcSheet As String, _
    ByRef strGFOobject As String, _
    ByRef strGFData As String, _
    ByVal iColumn As Integer, _
    ByVal data1 As Double, _
    ByVal data2 As Double, _
    ByVal data3 As Double, _
    ByVal data4 As Double)
-----
' Sets the data associated with a data column in an embedded GraFit object
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOobject      Name of the GraFit object
' strGFData         Name of the data table item within the GraFit object
' iColumn           The number of the data column (start at 1)
' data1             Data value # 1
' data2             Data value # 2
' data3             Data value # 3
' data4             Data value # 4
-----

Public Sub GF_SetColumnDataValue(ByRef strCalcSheet As String, _
    ByRef strGFOobject As String, _
    ByRef strGFData As String, _
    ByVal col As Long, _
    ByVal row As Long, _
    ByRef dataValue As Variant)
-----
' Sets a data value within a data table in an embedded GraFit object
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOobject      Name of the GraFit object
' strGFData         Name of the data table item within the GraFit object
' col               The number of the data column (start at 1)
' row               The number of the data row (start at 1)
' dataValue         The data value to be inserted
-----

Public Sub GF_EnsureInstanceClosed(ByRef strCalcSheet As String, _
    ByRef strGFOobject As String)
-----
' Ensures that this embedded GraFit object is fully closed. Activated
' objects run invisibly in the background. If you have more than one
' or two embedded objects then it is recommended you close them down
' when not needed in order to prevent memory becoming used up
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOobject      Name of the GraFit object
-----

Public Function GF_IsGraFitRunning(ByRef strCalcSheet As String, _
    ByRef strGFOobject As String) As Boolean
-----
' Returns True if this GraFit object is currently running
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOobject      Name of the GraFit object
-----

```

```

Public Sub GF_SetDataBlock(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String, _
    ByRef strGFData As String, _
    ByRef data() As Variant, _
    ByVal iColumnStart As Integer, _
    ByVal iColumnStep As Integer)
    '-----
    ' Sets data within a table in an embedded GraFit object
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    ' strGFData         Name of the data table item within the GraFit object
    ' data()            The 2D array of data to be inserted. Data will be provided
    '                  by the following route:
    '                  Dim d() As Variant
    '                  d = Worksheets(SHEETNAME).Range(RANGENAME).Value
    ' iColumnStart      The number of the data column to start setting data
    ' iColumnStep       How many GraFit columns to increment between steps
    '                  if iColumnStep=1 then the data block is contiguous
    '                  if iColumnStep=2 then every other column is set
    '                  (useful for setting xYxYxY.... data)
    '-----
Public Sub GF_SetColumnNameBlockVert(ByRef strCalcSheet As String, _
    ByRef strGFOBJECT As String, _
    ByRef strGFData As String, _
    ByRef data() As Variant, _
    ByVal iColumnStart As Integer, _
    ByVal iColumnStep As Integer)
    '-----
    ' Sets column headings within a table in an embedded GraFit object
    '
    ' strCalcSheet      Name of sheet holding the embedded GraFit object
    ' strGFOBJECT       Name of the GraFit object
    ' strGFData         Name of the data table item within the GraFit object
    ' data()            The 2D array of column names to be inserted. Data will be
    '                  provided by the following route:
    '                  Dim d() As Variant
    '                  d = Worksheets(SHEETNAME).Range(RANGENAME).Value
    '                  The data names are assumed to be in the leftmost column
    '                  of the specified data range
    ' iColumnStart      The number of the data column to start setting data
    ' iColumnStep       How many GraFit columns to increment between steps
    '                  if iColumnStep=1 then the data block is contiguous
    '                  if iColumnStep=2 then every other column is set
    '                  (useful for setting xYxYxY.... data)
    '-----

```

```

Public Sub GF_SetColumnNameBlockHoriz(ByRef strCalcSheet As String, _
    ByRef strGFOobject As String, _
    ByRef strGFData As String, _
    ByRef data() As Variant, _
    ByVal iColumnStart As Integer, _
    ByVal iColumnStep As Integer)
'-----
' Sets column headings within a table in an embedded GraFit object
'
' strCalcSheet      Name of sheet holding the embedded GraFit object
' strGFOobject      Name of the GraFit object
' strGFData         Name of the data table item within the GraFit object
' data()            The 2D array of column names to be inserted. Data will be
'                  provided by the following route:
'                  Dim d() As Variant
'                  d = Worksheets(SHEETNAME).Range(RANGENAME).Value
'                  The data names are assumed to be in the topmost row
'                  of the specified data range
' iColumnStart      The number of the data column to start setting data
' iColumnStep       How many GraFit columns to increment between steps
'                  if iColumnStep=1 then the data block is contiguous
'                  if iColumnStep=2 then every other column is set
'                  (useful for setting xYxYxY.... data)
'-----

```

Hard-coded Equations

Overview of Hard-coded Equations

Hard-coded equation definitions use specialized pre-compiled equation libraries to increase data fitting speed and/or provide greater flexibility in generating initial estimates. The GraFit program provides the capability to use hard-coded equations, though none are supplied as standard. At present some equipment manufacturers use these to allow GraFit to fit very large quantities of data more quickly than would be possible with the standard equation definition items. Hard-coded equation items simply substitute for the normal equation definitions when fitting data by non-linear regression (see Chapter 9). It should be stressed, however, that for the vast majority of applications the normal user-defined equation definition is quite sufficient.

What are Hard-coded Equations?

Hard-coded equations are pre-compiled code library files and have a .GFL file extension. They are normally supplied for specialist applications; it is not possible to generate them yourself. If you believe that you have an application that might benefit from use of a hard-coded equation library, please contact Erithacus Software technical support for further information.

Advantages

- Pre-compiled code results in faster analysis.
- The partial derivative calculations that are intrinsic to non-linear regression analysis can be performed explicitly rather than by numerical approximation, resulting in greater speed and better precision.

Disadvantages


- The definitions cannot be modified.
- It is not possible to write your own hard-coded equations.

Enabling Hard-coded Equations

The menu commands to insert hard-coded equations are not present by default, and must be explicitly enabled. The **Options** command¹⁰³ on the **View** tab, Options page controls this; to

¹⁰³ The **View ▶ Options** command in legacy mode

use hard-coded equations ensure that the **Enable hard-coded equation definitions** box is checked.

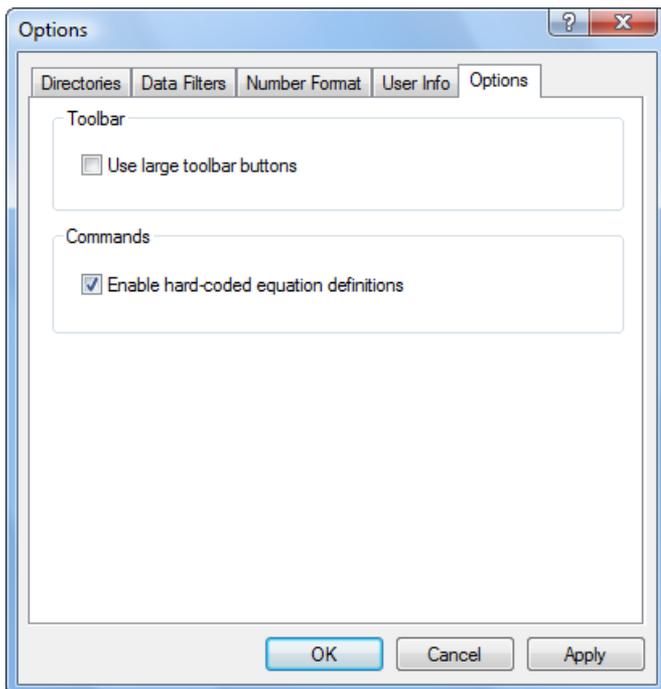


Inserting a Hard-coded Equation Item

Hard-coded equation items are added using the **Add ► Hard-coded Equation** command. The hard-coded equation window itself contains several panes that are used to define the equation. Each of the panes in this item window is discussed below.

General Pane

The General pane is used to define the location of the hard-coded equation library file, and shows a description of the equation as provided by the library supplier.

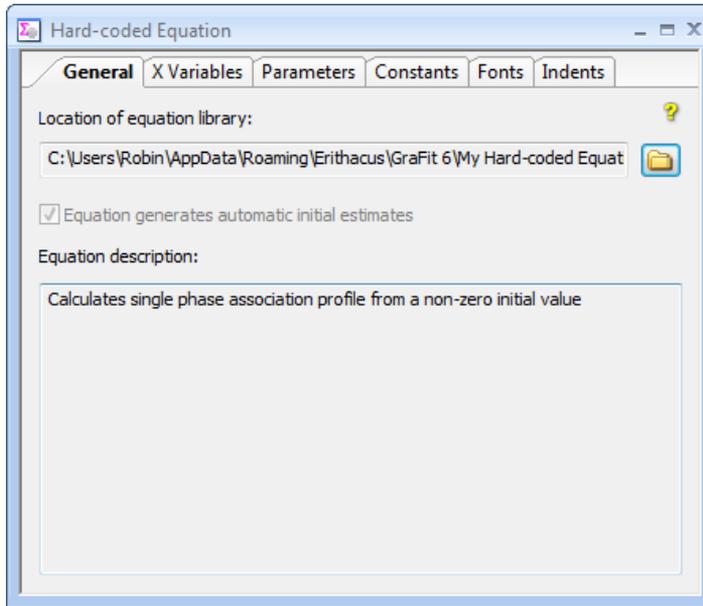


The General pane of the Hard-coded Equation item

Locating the Equation Library File

Click the browse button to locate the equation definition library file. Once loaded, the full path of this file is shown, and a description of the equation is displayed.

X Variables Pane


The X Variables pane is used to show the names and number of the independent (x) variables that are used in the equation. These names are pre-defined by the equation library vendor, and cannot be edited.

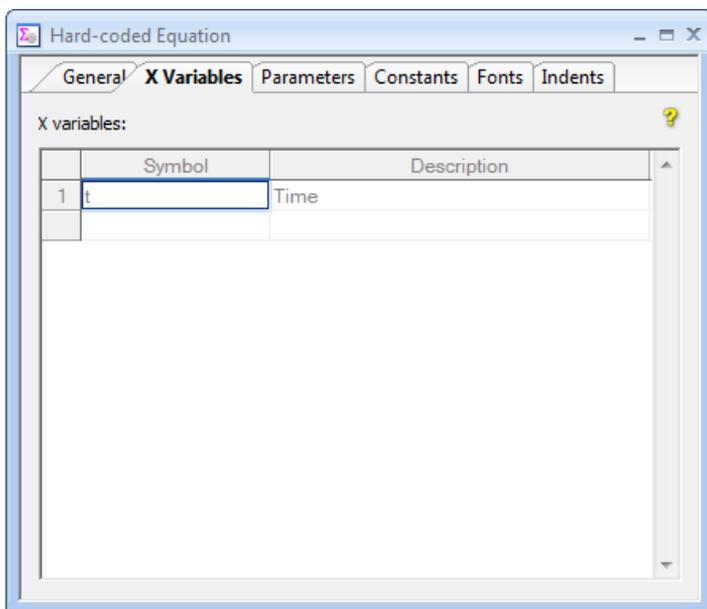

The X Variables pane of the Hard-coded Equation item

Parameters Pane


The Parameters pane is used to show the parameters in the equation definition. These names are pre-defined by the equation library vendor, and cannot be edited.

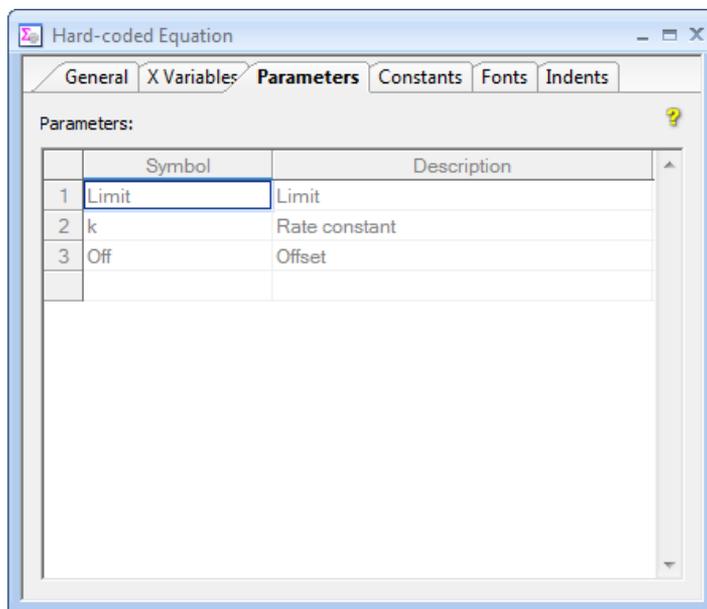

The Parameters pane of the Hard-coded Equation item

Constants Pane

The Constants pane is used to define constants in the equation definition. Constants are used to allow the entry of a value into an equation definition that is fixed in any single experiment, but varies from experiment to experiment. The presence of constants is optional, but when present an equation definition can include any number of constants. When constants are included in an equation used for non-linear regression, their values are entered into the non-linear fit item window as described in Chapter 6.

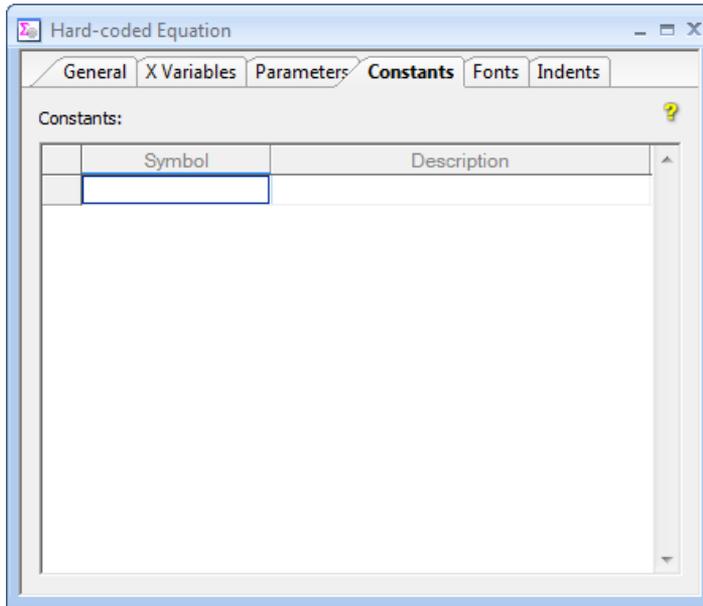

The Constants pane of the Hard-coded Equation item

GraFit Symbols Font

GraFit Symbols Font

GraFit is supplied with a font that can be used to represent the various graph drawing symbols used by the program. This can be useful when annotating graphs, for example, when labeling axes. The font is called **GraFit Symbols** and can be used just like any other Windows font.

The various symbols are associated with the characters as shown in the following table.

A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
○	□	△	▽	◇	◊	+	×	*	☆	÷	×	◆	⊕	☺	☹
a	b	c	d	e	f	g	h	i	j	k	l	m	n	o	p
●	■	▲	▼	◆	●	+	×	*	★	÷	×	◆	⊕	☺	☹
S	T	U	V	W	X	Y	Z	1	2	3	4	5			
▨	▩	▧	▦	▥	▤	▣	▢	□	■	▟	▞	▝	▜		

Note The GraFit symbols font does not contain any hinting information to assist reproduction on low-resolution displays. For this reason, it will not display well on screen, although it will print correctly.

Index

- - .BMP file, 154
 - .CSV file, 66, 67
 - .DIF file, 66, 67
 - .EMF file, 154
 - .GFC file, 49
 - .GFF file, 45
 - .GFL file, 325
 - .GFR file, 49
 - .TIF file, 154
 - .TXT file, 66, 67
 - .WMF file, 154
-
- _grnd_ function, 84, 119
 - _rnd_ function, 84, 119
- ## A
- Academic Licenses, 22
 - Activation
 - Software, 18
 - Activation FAQ, 31
 - Adair Equation, 277
 - Add command, 39
 - Add Data command, 186, 242
 - Add Data Table command, 53
 - Add menu, 39
 - Equation command, 109, 247
 - Graph command, 133, 139, 240, 243
 - Hard-coded Equation command, 326
 - Linear Fit command, 103
 - Non-linear Fit command, 89, 238, 248
 - Notepad command, 131
 - Transformation command, 125
 - Add Stock Item command, 42, 49, 238, 242
 - Adobe Acrobat, 18
 - Align command, 145
 - Allosteric Data, 246
 - Allosteric Kinetics, 270
 - Arithmetic operators, 117
 - Arrows, 152
 - ASCII format, 66, 67
 - Automation Reference, 315
 - Axis break symbol, 196
 - Axis manager, 190, 192, 193
 - Axis offset, 169
 - Axis scaling, 195
- ## B
- Bar charts, 136, 199
 - Bezier curves, 165
 - Bitmap file, 154
 - Bring to Front command, 146
- ## C
- Calculate menu

- Fit Data command, 103, 240
 - Calculated values, 78
 - Calibration Curve, 254
 - Changing Computers, 21
 - Chart Wizard, 199
 - Chi², 96
 - Chi-squared, 227, 230
 - Clear command, 39, 61, 62
 - Clipboard
 - Windows, 213, 214, 215
 - Clipboard, 211, 212
 - Column
 - Combinations, 71
 - Column charts, 133, 135, 199
 - Column Combinations, 71
 - Column constant, 83
 - Column Format command, 63
 - Column type, 67
 - A - B, 68, 71
 - A . B, 71
 - A / B, 71, 245
 - A + B, 71
 - Calculated Function, 82, 252
 - Column Summary, 72
 - Constants Entry, 78
 - Derivative, 85
 - Difference, 82, 243
 - F Test, 73
 - General, 69
 - Information, 77
 - Mean, 71
 - No. Data Points, 72
 - Parameter Entry, 77
 - Results, 76, 252
 - Series, 85
 - Smooth, 86
 - Standard Deviation, 72
 - Standard Error of the Mean, 72
 - Subset, 86
 - x Data, 69
 - x Values, 79, 254, 256
 - x(n) Data, 69
 - y Data, 70
 - y Data with Associated Value, 70, 258
 - y Values, 81, 254
 - Column Width command, 54
 - Columns
 - Copying and pasting, 59
 - Deleting, 62
 - Selecting, 59
 - Comment statements, 120
 - Competitive Inhibition, 260, 262
 - Compile command, 109, 125
 - Constants, 93, 114, 162, 185, 329
 - Linking to Data Table Entries, 93
 - Context Help, 43
 - Cooperative binding, 278
 - Copy command, 59, 60, 212, 215
 - Copy To command, 66
 - CSV format, 66, 67
 - Cursor shapes, 54, 144
 - Customize button, 36
 - Customizing Menus and Toolbars, 37
 - Cut and paste, 211
- ## D
- Data
 - Association, 194
 - Copying, 59
 - Deleting, 61
 - Export, 66
 - Fitting several data sets simultaneously, 256
 - Import, 66
 - Import / Export, 66
 - Inserting values, 62
 - Layout, 90, 104
 - Missing, 55
 - Numbers, 55
 - Overlaying, 189

- Pasting, 59
- Rescaling, 82
- Selecting, 58
- Smoothing, 86
- Sorting, 63
- Text, 55
- Types, 55

Data entry position, 53, 54, 57, 59, 66, 68, 95

Data Filters, 49

Data fitting

- Basics, 217
- Multiple Non-linear Regression, 224

Data Fitting, 75

Data function, 83

Data table

- Appearance, 63
- Locking columns, 87

Data table item, 53

Deactivation of GraFit, 21, 28

Definition, 115

Degrees of freedom, 227

Delete command, 62

Denaturation Curves, 259

Derivative, 85

DIF format, 66, 67

Difference calculation, 67

Digitally Signatures, 298

Directories, 47, 50

Dixon Plot, 261

Dose-response Curves, 273

Double exponential decay, 272

Drag and drop, 39

Draw menu

- Results Listing command, 98, 107, 240

Drawing Objects, 140

E

Edit Current Cell command, 54

Edit menu

- Clear command, 39, 61
- Copy command, 59, 212, 215
- Copy To command, 66
- Edit Current Cell command, 54
- Hide Object command, 167, 194
- Insert New Object command, 215
- Insert Object command, 214
- Paste command, 59, 212, 213, 214
- Paste From command, 49, 66
- Paste Link command, 212
- Paste Special command, 213, 214
- Properties command, 146, 149, 165, 194
- Rename command, 40
- Reset Original Size, 215
- Select All command, 215

Editing item names, 40

Email support, 15

Embedded object, 40

Embedded Objects, 211

Embedding and Linking, 297

Embedding objects, 211

Enhanced metafile, 154

Enzyme inhibition

- Competitive, 267
- Mixed, 267
- Non-competitive, 268, 269

Enzyme Inhibition, 260

Enzyme kinetics, 270, 287

- Ping Pong, 266
- Ternary Complex, 266

Equation definition, 109, 115, 116

Equation definition syntax, 112

Equation menu

- Compile command, 109, 124

Equations, 265

- Compiling, 124
- Defining, 246

- Error Messages, 124
- Functions, 118
- Multi-line Definitions, 119
- Multiple Regression, 224
- Temporary variables, 119
- Erithacus Software
 - Contacting, 15
- Error bars, 181
- Error Message, 129
- Error Messages, 124, 129
- Errors, 221, 227, 228, 229, 231, 233
- Estimates
 - Interactive, 95
- Estimatrix, 81, 94, 96, 238, 239
 - Linking to Data Table Entries, 95
- Evaluation Activation Code, 20
- Evaluation Mode, 19
- Excel
 - Enable Macros, 297
 - External Data, 308
 - Passing Data, 298
 - Security Warning, 297
- Exponential decay, 272, 292
- Exponential decay with Offset, 293
- Export Graph command, 154
- F**
- F test, 73, 235
 - Nested models, 236
- Field Code Command, 149
- Field Codes, 147
- File
 - .GFF, 45
- File menu, 35, 38
 - Add Stock Item command, 42, 49, 238, 265
 - Close command, 49
 - Export Graph command, 154
 - Import Data command, 49, 66
 - New command, 42, 50, 51, 285
 - Open command, 45, 238, 246
 - Page Setup command, 41, 139
 - Print command, 40, 97, 106
 - Print Preview command, 40, 97, 101, 102, 106, 240
 - Save As command, 48, 50
 - Save command, 48
 - Store Stock Item command, 42, 49
- Files
 - Recently-Used, 47
 - Creating new, 49
 - Exporting graphs, 154
 - Importing Data, 49
 - Opening, 45
 - Pinning to the File List, 48
 - Saving, 48
 - Saving Legacy Versions, 48
 - Template, 49
- Filters
 - Import / Export, 49
- Fit Data command, 40, 89, 240
- Fitting criteria, 96
- Font selection, 100, 101
- Font settings, 107, 131, 153, 154
- Footers, 42
- For statement, 121
- Formatting numbers, 63, 99, 172
- Formatting Options, 101
- Frequently asked Questions, 31
- Functions, 118
- G**
- GF_module.bas, 313, 320
- Global Analysis, 223
- Goto Statement, 121
- GraFit
 - Installing, 17
 - Running, 35
 - Technical Support, 15

- Templates, 42
 - Uninstalling, 21
 - GraFit (Academic), 22
 - GraFit (Student), 22
 - GraFit Editions, 22, 23
 - GraFit Fitter Excel Template, 300
 - GraFit Fitter Reference, 302
 - GraFit Items, 39
 - GraFit Pro, 22
 - GraFit Pro (Academic), 22
 - GraFit Workspace, 38, 45
 - Graph
 - Axis associations, 186
 - X ticks, 168
 - Graph Magnification, 142
 - Graph menu
 - Add Data command, 186, 242
 - Axis manager, 190
 - New Graph command, 133
 - Preferences command, 140, 147
 - Graph types, 133
 - Graphs, 139
 - Adding data, 186
 - Aligning objects, 145
 - Arrow properties, 152
 - Axis association, 183
 - Axis scale, 172
 - Axis text, 174
 - Chart bar style, 203
 - Colors, 146
 - Creating, 133
 - Curve resolution, 163
 - Data point style, 204
 - Data point styles, 161
 - Data/curve properties, 176
 - Default text properties, 147
 - Duplicating objects, 143
 - Error plots, 159, 181, 200, 208
 - Exporting, 154
 - Fill properties, 150
 - Font properties, 153
 - Grid settings, 169
 - Grouping objects, 145
 - Histogram bar style, 210
 - Initial scaling, 208
 - Initial scaling, 159
 - Join style, 205
 - Join styles, 164
 - Legend, 157, 200, 207
 - Legend box, 174, 186
 - Line properties, 151
 - Location, 172
 - Main graph region, 166
 - Multiple axes, 189
 - Object locations, 152
 - Plotting equations, 162
 - Rotation properties, 151
 - Scaling, 201
 - Selectable regions, 166
 - Selecting objects, 142
 - Sizing objects, 143
 - Split axes, 189
 - Text, 147
 - Text objects, 143
 - Tick styles, 160, 202, 209
 - Transformations, 159
 - Visibility, 167
 - X scale, 170
 - Y scale, 172
 - Y ticks, 171
 - Z ordering, 146
 - Grid Lines command, 140
 - Grid Magnification, 65
 - Group command, 145
- ## H
- Hard-coded Equation definition, 325
 - Hard-coded Equations
 - Enabling, 325
 - Inserting, 326
 - Library Files, 327
 - Hardware Requirements, 17

Headers, 42
Help commands, 42
Help menu
 Contents command, 43
Henderson-Hasselbalch equation, 237
Hide Object command, 167, 194
Hill Plot, 270
Histogram, 136, 207
Histogram Wizard, 207

I

IC50 equations, 273, 288, 289
Icon
 Changing, 40
Icon command, 40
Ignore Data command, 56
Import Data command, 49, 66
Import/Export Filters, 49
Indents, 100, 102, 107
Initial estimates, 94
Initial Estimates, 121, 123, 220, 238
Insert command, 62
Insert New Object command, 110, 215
Installation, 18
Interactive Estimator, 95, 238
Internet, 15
Interpolate between X limits, 80
Items
 Loading and Saving, 49

J

Junk email, 26

L

Legend Box, 174
 Margin settings, 175

Licence Control Panel, 25
License
 Renewal, 30
License Administrator, 25
Licensing FAQ, 31
Ligand binding, 290
Ligand binding equations, 274
Limiting x data range, 93, 105
Linear fitting, 89, 103
Linear regression, 230, 291
 Basics, 220
Lineweaver Burk plot, 222, 287
Linked object, 40
Linked Objects, 212
Linking, 214
Linking objects, 211
Lock Column command, 87
Logarithmic scaling, 170
Logical operators, 117

M

Magnification, 65, 141
 Graph, 142
 Grid, 65
Marquart algorithm, 232
Mask curve option, 161
Mean, 71, 236
Menu and Toolbars, 35
Metafile, 154
Michaelis-Menten enzyme kinetics, 287
Michaelis-Menten Equation, 270
Microsoft Excel, 298, 315
Mixed Inhibition, 262
Move Back command, 146
Move Forward command, 146
Multiple axes, 189

- Multiple Axes, 194
- Multiple Regression, 222, 251
 - Advantages, 224
 - Disadvantages, 224
 - Plotting results, 226
 - Uses, 223
- Multiple x variables, 91, 113
- Multi-user license, 21, 24, 26

- N**
- Nested models, 236
- New Graph command, 133
- NMR Titration Curve, 237
- Non-competitive Inhibition, 262
- Non-linear fitting, 89
- Non-linear regression, 232
- Non-linear regression Basics, 220
- Notepad item, 131
- Number formatting, 63, 98, 99, 106, 172, 173

- O**
- Object embedding and linking, 211
- Object menu, 146
 - Align command, 145
 - Group command, 145
 - Ungroup command, 145
- OLE automation, 297
- OLE objects, 111, 139, 140, 211, 215
- Options command, 46, 47, 49, 50, 64, 99, 173, 325
- Options menu
 - Style command, 43
- Outliers, 229, 233
- Outlines command, 140

- P**
- Page Settings, 41
- Page Setup command, 140
- Paired Data Sets, 256
- Parameters, 113, 328
- Passing Data, 297
 - Automation, 309
 - Sample Files, 313
- Paste command, 59, 60, 66, 212, 214, 215
- Paste From command, 49, 66
- Paste Special command, 214, 215
- PDF files, 18
- Ping Pong Kinetics, 266
- pKa determination, 279
- Polynomial regression
 - Basics, 220
- Preferences command, 68, 140, 147
- Print preview, 40
- Print Preview, 97, 98, 240
- Printer
 - Changing, 42
- Printing, 40
- Properties command, 68, 146, 147, 150, 165, 194, 195
- Proxy Settings, 19

- Q**
- Quick Access Toolbar, 35, 36

- R**
- Random numbers, 84
- Rate equations, 281
- Recently-Used Files, 47
- Reduced chi-squared, 227
- Regression analysis, 90, 220, 227
 - Assumptions, 228
 - Basics, 219
 - Linear, 230
 - Non-linear, 232

Reset Original Size command, 215
Residual plots, 228, 234
Residuals, 235, 242
Residuals plot, 128
Results, 97, 105
Results Listing command, 98, 107, 240
Ribbon Bar, 35
Robust weighting, 90, 229, 233, 245
Rotation of objects, 151
Row constant, 84
Row Height command, 54
Rows
 Selecting, 59

S

Savitzky-Golay smoothing, 85, 86
Scaling power, 174
Scatchard plot, 290
Screen resolution, 286
Security Warning
 Excel, 297
Select All command, 215
Selection Tool, 141
Send to Back command, 146
SetBlank function, 83, 84
Single-user license, 24
Size Column to Fit command, 54
Software Activation, 18
Software Updates, 15, 21
softWORKZ, 25
Sort command, 63
Spline curves, 165
Split axes, 189
Splitter box, 53
Standard Curve, 254

Standard deviation, 72, 236
Standard Error of the Mean, 72
Statistics, 71
Status Bar, 89, 103, 109, 126
Stock Items, 42
Store Stock Item command, 49
Student Licenses, 22
Symbol syntax, 112
Symbols font, 331
Syntax
 Symbol names, 112
System Requirements, 17

T

Table menu
 Column Format command, 63
 Column Properties command, 68
 Column Width command, 54
 Delete command, 62
 Ignore Data command, 56
 Insert command, 62
 Lock Column command, 87
 Preferences command, 63
 Row Height command, 54
 Size Column to Fit command, 54
 Sort Column(s) command, 63
Technical Support. *See* GraFit, Technical Support
Template files, 47, 49
Templates, 42, 51, 285
 Standard, 285
Temporary variables, 119
Text Alignment, 173
Text format, 66, 67
Text menu, 147
Text Orientation, 173
Text settings, 173
Theory, 227

TIFF file, 154

Transformation menu

 Compile command, 125, 129

Transformations, 125, 159, 179, 180

 Compiling, 129

 Curves, 128

 Data points, 126

 Defining, 129

 Error Messages, 129

Type Library, 315

U

Uncompetitive Inhibition, 262

Ungroup command, 145

Uninstalling GraFit, 21

Users

 Adding, 28

 deleting, 28

V

ValidData function, 83, 84

VBA, 309

View menu

 Grid Lines command, 140

 Options command, 46

 Show Outlines command, 140

 Zoom command, 65, 141

Visibility

 Graph objects, 167

Visual Basic, 316

Visual Basic for Applications, 309

W

Weighting, 90, 104, 221, 229, 231

 Checking, 242

 Determining weighting type, 234

 Explicit, 246

Window Display Styles, 43

Window menu

 Icon command, 40

Windows 2000, 46, 66, 155

Windows 7, 46, 66, 154

Windows metafile, 154

Windows Vista, 45, 66, 154

Windows XP, 46, 66, 155

X

X Variables, 111, 113, 327

X/Y Scatter Graph, 134, 157

XY graph wizard, 157, 176, 177, 178, 185, 186

Z

Zoom command, 65, 141

Zoom Slider, 65, 141

