


# iCom Bluetooth Pairing Guides Motorola

PHONAK

life is on

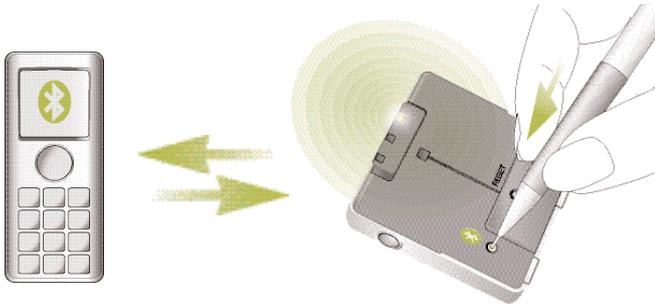

## Pairing and connecting iCom with your Motorola mobile phone

**1. Charge and turn on both mobile phone and iCom.  
Place them next to each other.**

### 2. Prepare your mobile phone for pairing

Find the connectivity settings of your phone. With most Motorola models one of the following navigation options applies.

Menu ▷ Connection ▷ Bluetooth link

or

Menu ▷ Bluetooth link ▷ Paired Devices ▷ Option

or

Settings ▷ Connection ▷ Bluetooth Link ▷ Hands-free

or

Control Panel ▷ Connections tab ▷ Bluetooth ▷ Devices tab

### 3. Prepare your iCom for pairing

Initiate the pairing process on your iCom by pressing the Bluetooth pairing button for about 2 seconds until the Bluetooth indicator flashes blue. Your iCom will remain in pairing mode for two minutes or until the pairing is completed.

### 4. Pair iCom and mobile phone

On your mobile phone, still in the Bluetooth menu, select one of the following (depending on the model):

[Add Audio Device]

or

[New Paired Device]

or

Look for Devices

or

Add

(If you have not activated the Bluetooth functionality, you may be asked if you wish to turn on Bluetooth. Doing so will start the search process for Bluetooth devices.)

Your mobile phone should present a list of Bluetooth devices discovered.

Select "iCom" from this list. Some phone models will ask **Bond with iCom?** If the phone asks to bond with iCom, confirm with "yes".

Enter passcode "0000" (four zeros).

iCom will be added to the Bluetooth devices list.

### 5. Connect iCom with mobile phone

Typically, iCom will be automatically connected after pairing.

On some phones you can configure automatic connection in the access settings of the device details. Choosing so allows iCom to connect automatically without confirmation. We recommend that you enable this feature if available.

If you have difficulties navigating to the Bluetooth settings, please refer to the mobile phone's user manual or ask the mobile phone dealer.