

Map of Medicine

EMIS installation guide

Contents


Purpose of the document	1
1. Integration steps	1
2. Activate Map of Medicine	2
3. Setup API Passwords	2
Add users to a product	3
Auto Login	4
After Installation - Using a Partner Product for the First Time	4

Purpose of the document

This document is intended for any EMIS Web Organisation attempting to connect to Map of Medicine. Please follow the steps below in order to activate Map of Medicine for use with your clinical system.


1. Integration steps

- Go to: System tools > EMAS Manager

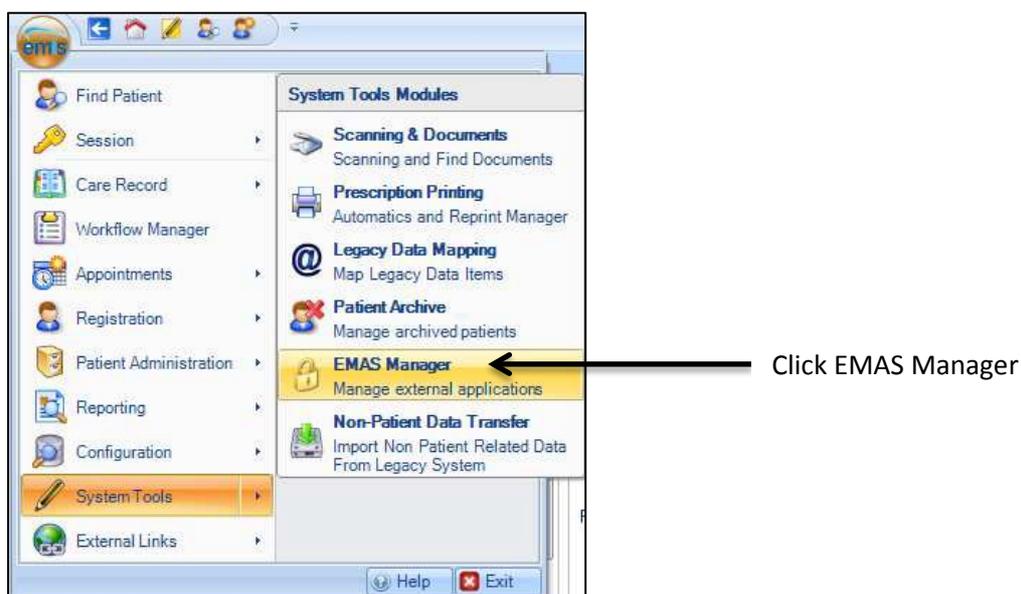

2. Activate Map of Medicine

- Select Map of Medicine from the Partner API list and then click Activate Application

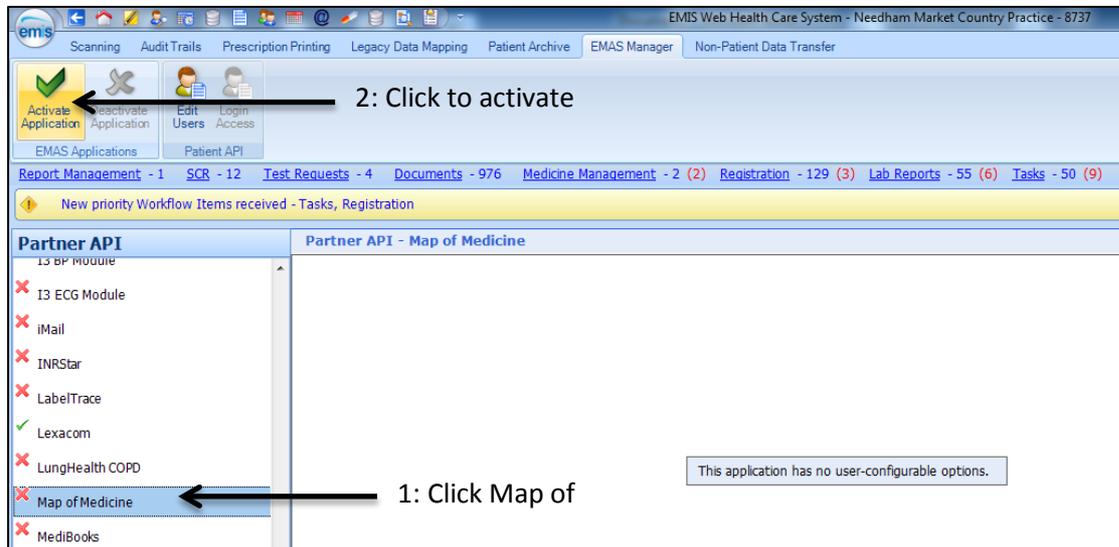


3. Setup API Passwords

You will need to specify which users in your practice you would like to access Map of Medicine through EMIS.

- Go to: System tools > EMAS Manager > Partner API > Edit Users (fig 1)
- Tick each user who requires access to (fig 2)




(Fig 1)

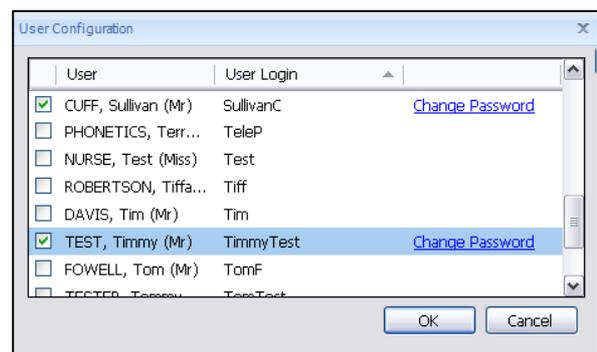

(Fig 2)


Note: If you are adding a new user that doesn't have a tick next to their name already, you will be prompted with the following screen

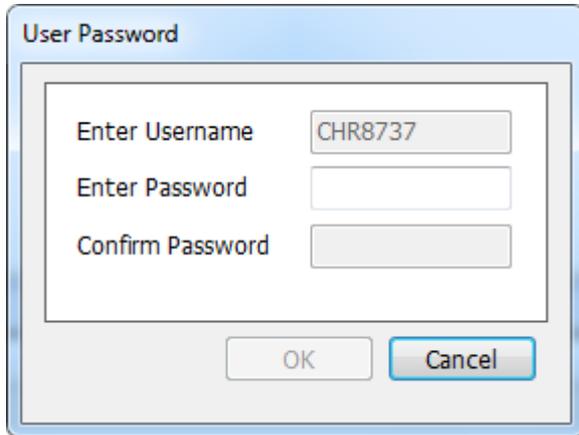


- You will need to add a password. In this case it is the user's manual EMIS login password (i.e not smartcard)

Add users to a product

- With the product still highlighted, click on *Login Access* (Fig 3)
- Tick *Auto Login* for each user that needs to use the Auto Login facility (see below for more details on the Auto Login facility).
- Tick *Allow Login* for each user that needs to have access to the product.
- Click on *OK* to save changes.

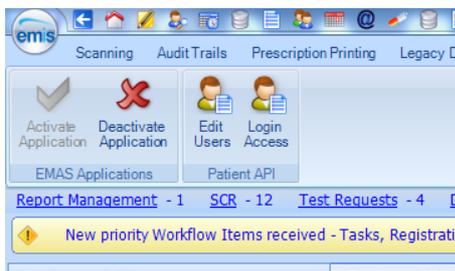

(Fig 3)

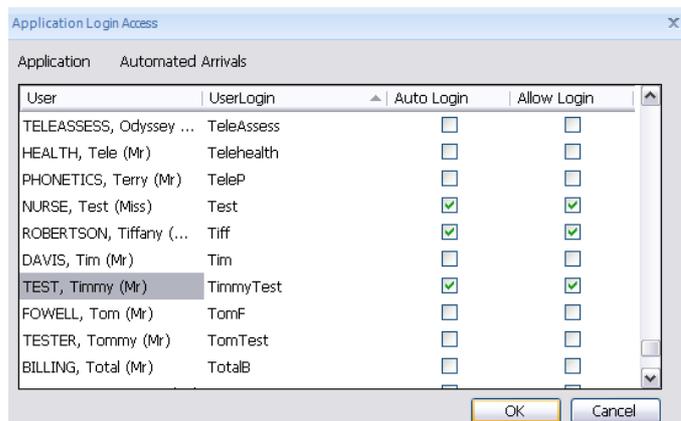

(Fig 4)

Auto Login

Auto Login allows authenticated EMIS Web users to connect to Partner products without having to re-enter their credentials. If Auto Login is enabled and the user is already logged into an EMIS Web client session, the API will automatically use their EMIS Web credentials to authenticate the login of a Partner product.

After Installation - Using Map of Medicine for the First Time

The user will be prompted to enter their EMIS Web username and password when a Partner product connects to their machine for the first time (*Fig5.*).



The screenshot shows a dialog box titled "emisWeb" with a sub-header "Activate New Partner Product". The text inside the dialog reads: "This is the first time this product has been used on this machine. Please enter your username and password to activate the product". There are three input fields: "Username:" with the value "Timmy Test", "Password:" with masked characters "xxxxxxxx", and "CDB:" with the value "50003". Below the fields are two buttons: "Activate Product" and "Cancel". To the right of the buttons are four colored squares (cyan, blue, cyan, blue). At the bottom of the dialog, there is a copyright notice: "Copyright © Egton Medical Information Systems 1988 - 2007. All rights reserved."

(Fig 5)