


Analytical Sample Preparation

2014 Catalog



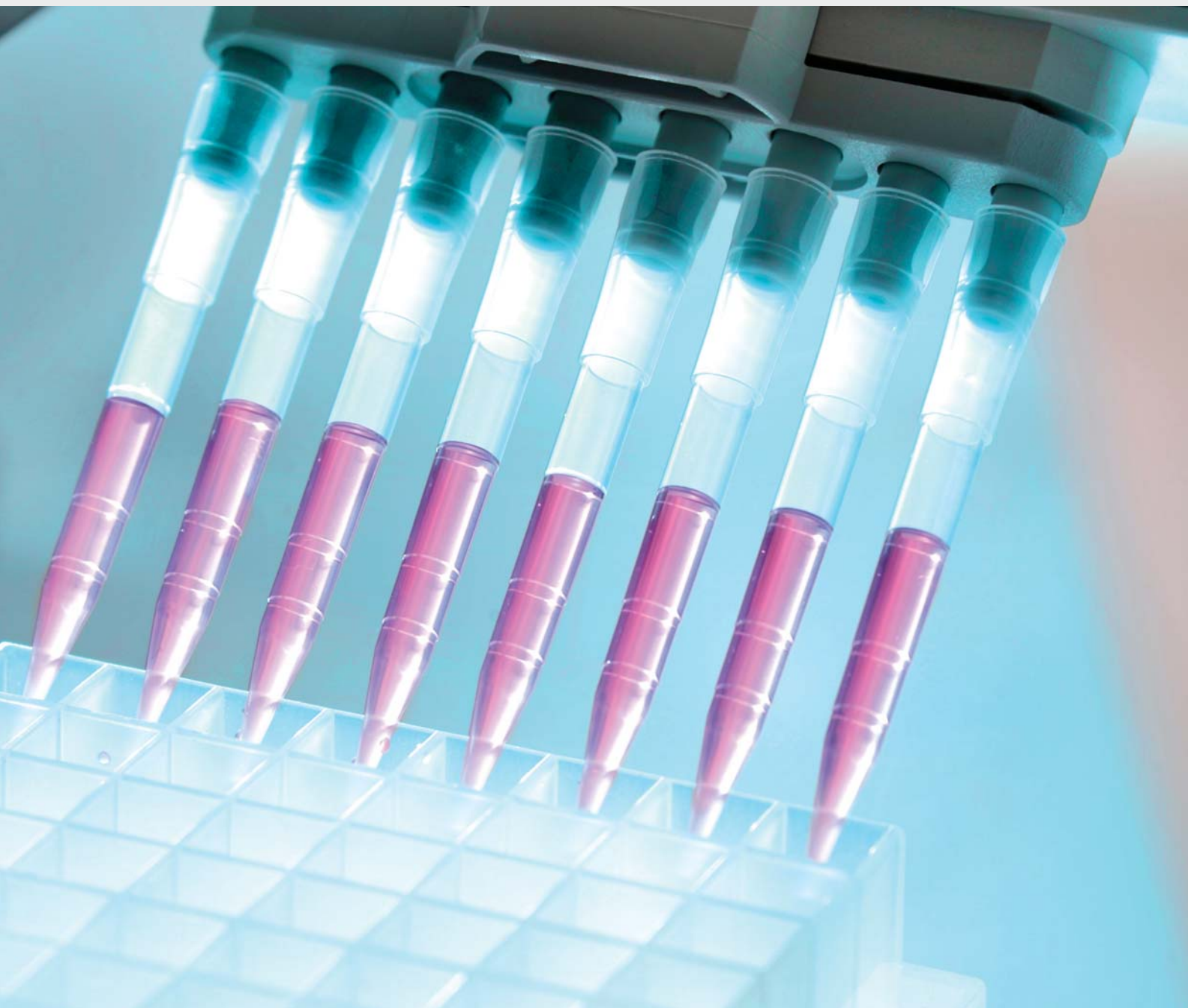

Analytical Sample Preparation
2014 Catalog

Contents

Biotage® 1-Point Support™	3
Sample Preparation Format Options	9
ISOLUTE® SLE+ Supported Liquid Extraction Columns and 96-Well Plates	13
EVOLUTE® Sample Preparation Products	15
ISOLUTE® SPE Columns and 96-Well Plates	19
ISOLUTE® Myco SPE Columns	37
ISOLUTE® Bulk SPE Sorbents	39
ISOLUTE® Protein Precipitation Products	41
ISOLUTE® QuEChERS	43
Sample Preparation Accessories	45
Manual Sample Processing Products and Accessories	49
Biotage® PRESSURE+ 48 and PRESSURE+ 96	
Biotage® VacMaster™-10 and -20	
Biotage® VacMaster™-96	
Sample Processing Accessories	
Gravity Rack	
Automated Sample Processing Products	55
RapidTrace®+	
Evaporation	57
Biotage® SPE Dry 96 and 96 Dual	
TurboVap®	
Index	60
Part Number Index	62



Biotage® 1-Point Support™

The Answer to All Your Questions

Method Development and Troubleshooting Advice

Biotage is ready to assist you with our team of Analytical Chemists who have many years of experience in providing practical, theoretical and technical knowledge relating to sample preparation techniques.

If you need help choosing the best sample preparation method for your application, or want to solve an existing sample preparation problem, please contact your local Biotage 1-Point Support™ team who will ensure you receive fast advice from our sample preparation experts.

Technical Information

Throughout this catalog, we highlight chemistry data sheets and application notes that provide detailed information on sample preparation techniques and SPE sorbents.

Sample Preparation Application Database

Biotage's Literature Library is a vital resource for analytical chemists looking to optimize sample preparation procedures. As sample preparation experts, Biotage are dedicated to finding the ideal solution to your sample preparation requirements with a regularly updated library containing application notes, technical notes and scientific presentations. The Biotage Literature Library has been optimized to provide the most appropriate, efficient and effective sample preparation solutions, taking into consideration matrix effects, analyte structure/functionality and relevant regulatory requirements.

The Biotage Literature Library is fully searchable by keyword, analyte, matrix, analytical technique, format, product type and industry type. This search feature allows you to obtain relevant application notes suitable for your needs in an efficient and dynamic manner. Visit the Biotage Literature Library www.biotage.com/applications.

ISOLUTE® SLE+ User Guide


Supported Liquid Extraction (SLE) is a quick, easy and cost effective sample preparation technique for analytical chemistry. This guide describes how to develop supported liquid extraction methods using ISOLUTE® SLE+ products, along with hints and tips for optimizing performance and extending the range of analytes that can be extracted.

EVOLUTE® User Guide

The EVOLUTE® User Guide contains a wealth of technical information and practical tips to make method development using EVOLUTE SPE products fast and simple. The guide describes generic methodologies and method optimization strategies for use with all the EVOLUTE sorbents, and introduces the streamlined Load-Wash-Elute procedure for EVOLUTE EXPRESS products.

SPE QuickStart Guide

Use our QuickStart Guide for an introduction to SPE.

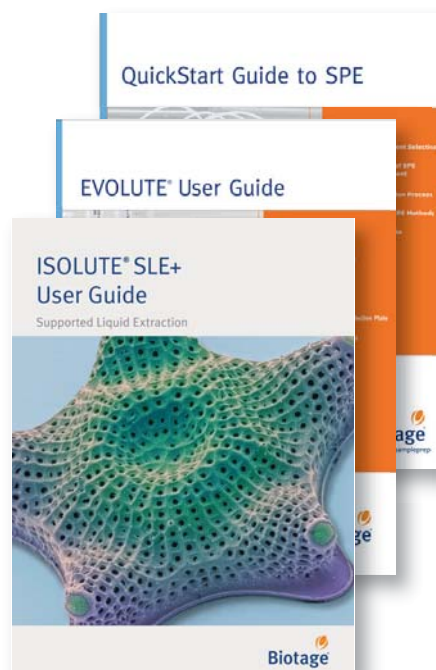

To request a copy of any of the above guides please visit www.biotage.com

Contact Biotage® 1-Point Support™

www.biotage.com

The Biotage web site offers our customers easy access to current information on new products, applications, and events.

Europe

Service and Support

Telephone: +46 18 56 59 11

E-mail: eu-1-pointsupport@biotage.com

North & Latin America

Service and Support

Telephone: +1 800 446 4752

Outside the USA: +1 704 654 4900

E-mail: us-1-pointsupport@biotage.com

Japan

Service and Support

Telephone: +81 3 5627 3123

E-mail: jp-1-pointsupport@biotage.com

China

Service and Support

Telephone: +86 21 2898 6655

E-mail: cn-1-pointsupport@biotage.com

To locate a distributor

please visit our web site at
www.biotage.com



How to Place Your Order

In this catalog product information is displayed with easy-to-follow application details and supporting documentation. Ordering information is found following each product listing.

How to Order

You can now order through our web site. Just log on to www.biotage.com and place your order. Alternatively you may place your order by telephone, email or fax.

Europe

Main Office: +46 18 56 59 00
Toll Free: +800 18 56 57 10
Fax: +46 18 59 19 22
Order Tel: +46 18 56 57 10
Order Fax: +46 18 56 57 05
order@biotage.com

North & Latin America

Main Office: + 1 704 654 4900
Toll Free: +1 800 446 4752
Fax: +1 704 654 4917
Order Tel: +1 704 654 4900
Order Fax: +1 434 296 8217
ordermailbox@biotage.com

When placing your order please have available:

- Your purchase order number
- Biotage part number(s)
- Product description(s)
- Shipping address
- Billing address
- Contact person, including telephone number
- Product user name and department

Japan

Tel: +81 3 5627 3123
Fax: +81 3 5627 3121
jp_order@biotage.com

China

Tel: +86 21 2898 6655
Fax: +86 21 2898 6153
cn_order@biotage.com
www.biotage.cn

Distributors

Please visit our Web site at
www.biotage.com for contact details.

Orders can also be placed using your VISA or MasterCard account (and American Express **IN THE US ONLY**).



Biotage is committed to developing and manufacturing sample preparation products of the highest quality. State-of-the-art manufacturing techniques are supported by a comprehensive Quality Control (QC) testing program documented under our ISO9001:2008 registered Quality Management System.

All the components used to make our sample preparation products (tubes, frits, 96-well plates and SPE sorbents) are rigorously cleaned and QC tested to ensure they meet our demanding purity specifications.

We use sophisticated instrumental techniques to confirm the physical and chemical nature of every batch of SPE sorbent – ensuring reproducible performance in your application.

Every ISOLUTE® or EVOLUTE® SPE product is accompanied by a detailed Quality Assurance (QA) report for your reference. This page explains the importance of the information it contains, and the impact this has on every sample preparation procedure you perform.

Reproducible sorbent mass packed into SPE columns:

- Column capacity and analyte elution volume requirements are constant
- Consistent column-to-column and batch-to-batch recoveries

Well controlled particle size distribution with minimal fines from column-to-column and batch-to-batch:


- Reliable sample processing using manifolds and automated devices
- Consistent flow through columns, from column-to-column and batch-to-batch
- Low back pressure – automation friendly
- Gravity loading of samples with some column configurations
- No fines in final extract to plug injectors or absorb analytes when sample is reconstituted in another solvent
- No channeling in sorbent beds minimizes sorbent mass requirements, reducing elution volumes and costs
- Minimal drying time, reproducible from batch-to-batch
- Large sorbent mass columns with good flow characteristics



Quality Assurance Report

Part No. 221-0010-A
Description ISOLUTE C18(EC) 100mg 1ml SPE Columns
Lot No. 10221002GP **Sorbent Batch No.** 221-0-02

This product has been manufactured, tested and assembled under the control of an ISO 9001:2008 registered quality system and has been subjected to the following QC tests:

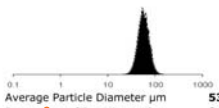Base Material Analysis Pore Diameter (Å) 50 Surface Area (m ² /g) 491 Surface pH N/A Elemental Analysis Carbon Loading (%C) 18.5 Bonded Phase Coverage (µmol/m ²) 1.6 UV Abs. of MeOH Extract @ 254 nm <0.1 Extraction Residue (%) <0.1 Particle Size Analysis C18(EC) BATCH No 221-0-02  Average Particle Diameter µm 53 Percentage of Particles <10µm 0.01 Percentage of Particles >90 2.73 Percentage of Particles 30-90 95.48 This ISOLUTE product has passed Biotage's rigorous QC tests.	Material HPLC Analysis Mobile Phase: 60%MeOH / 40%H2O Uracil Acetophenone Nitrobenzene Methyl Benzoate Toluene Capacity Factor K' Selectivity, Alpha 4,3 End Capping Test Exchange Capacity Special Recovery Test Packed Column and Component Test Column Flow Resistance Column Purity Test (GC) Frit Purity Test (GC) Material Weight Check N/A = Not Applicable Biotage welcomes feedback on all quality issues.
---	--

G. Phillips
Quality Control Manager

Europe Service and Support Telephone: +46 18 56 59 11 E-mail: EU-1-PointSupport@biotage.com	North America Service and Support Telephone: +1 800 446 4752 press (3) at the auto attendant e-mail: US-1-PointSupport@biotage.com	Japan Service and Support Telephone: +81 3 5627 3124 E-mail: JP-1-PointSupport@biotage.com
---	---	--

Europe and ROW
 If your country is not listed here, contact your local distributor. Please visit our website at www.biotage.com for your local distributor contact details.

© 2010 Biotage. All rights reserved. ISOLUTE is a registered trademark of Biotage AB.



Reproducible, optimized chemistry for ISOLUTE® and EVOLUTE® sorbents:

- Minimize method development time
- Eliminate need for time-consuming method changes when different batches of sorbent are used

High and reproducible capacities of ISOLUTE® ion exchange sorbents:

- Reproducible high recoveries in ion exchange SPE without using columns with excessively large sorbent mass
- Saves time, money and increases analyte concentration in final extract

High purity sorbents, frits, columns and 96-Well SPE plates:

- The purity of all the components used to manufacture ISOLUTE® and EVOLUTE® SPE products are monitored to ensure compatibility with the most demanding applications






Sample Preparation Format Options


Sample Preparation Format Options

True-to-Scale Columns

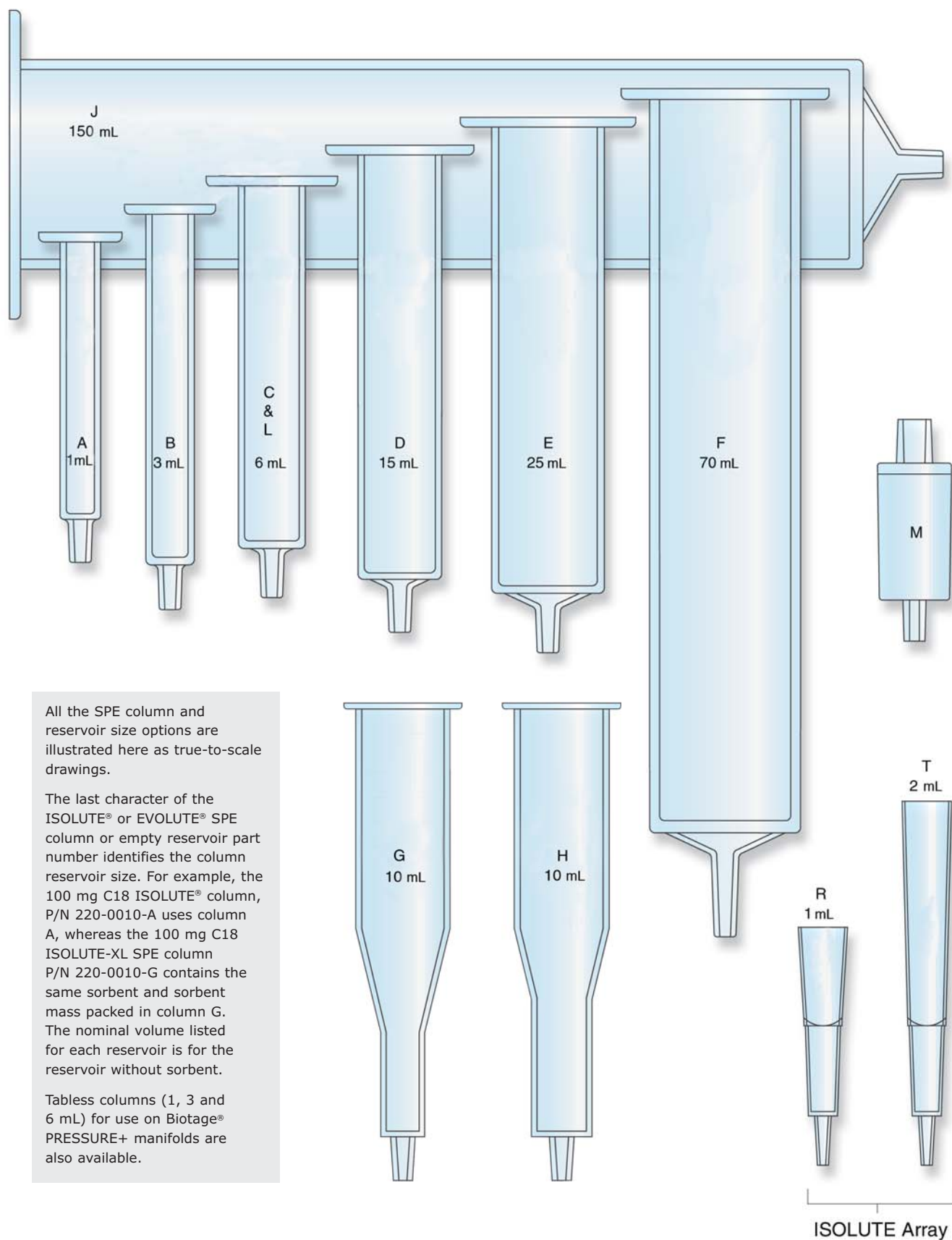


Sample Preparation Format Options

True-to-Scale 96-Well Plate (Fixed Well and Array Plate)

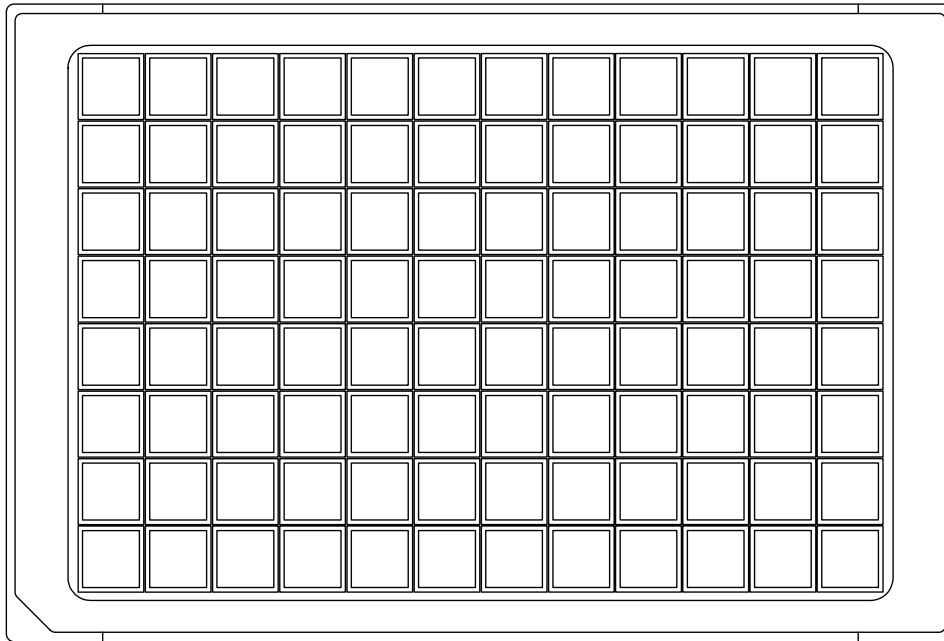

Biotage supplies two 96-Well plate formats, the fixed well plate and the modular Array. As with individual columns, the last character(s) of the part number identifies the format/column type. The suffix for fixed well plates is -P01, while the suffix for Array wells is -R (1 mL) and -T (2 mL). Pre-assembled Array plates have the suffixes -RP and -TP.

Scale diagrams (60% of actual size) of the fixed well and pre-assembled Array plate are shown below. Note the 'skirt' heights (distance from well outlet to sealing edge) of the two plates are different. To eliminate cross talk when processing plates, well outlets should penetrate the collection plate correctly. To accommodate this skirt height difference, the VacMaster-96 is supplied complete with a spacer for use with the Array format. This spacer is not required for processing the fixed well format.

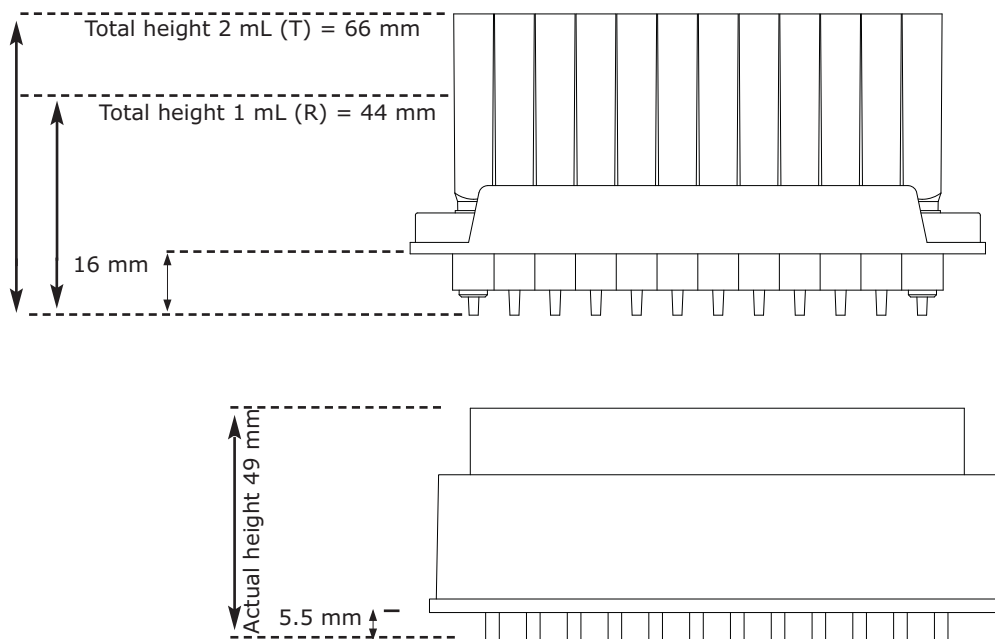

Side elevation (60% of actual size) of populated Array plate (top) and fixed well plate (bottom).





ISOLUTE® SLE+ Supported Liquid Extraction Columns and 96-Well Plates

Simple Load-Wait-Elute Methodology

Supported Liquid Extraction

Achieving Simplicity and Success in Sample Preparation

ISOLUTE® SLE+

Plates and Columns for Supported Liquid Extraction

Improve productivity and maximize analyte recovery using a Simple Load-Wait-Elute methodology.

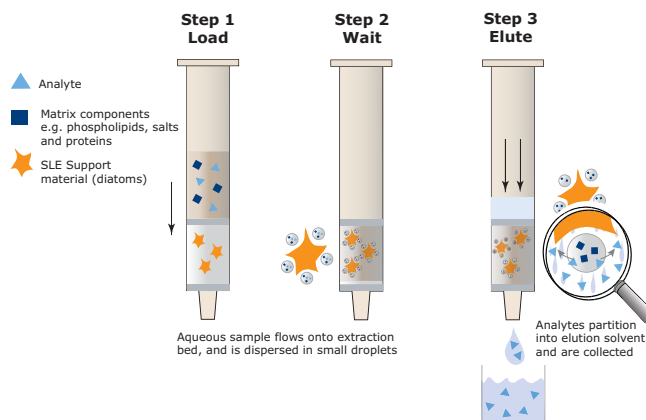


Figure 1. Supported Liquid Extraction Mechanism.

Chemical Description:

Modified form of diatomaceous earth.

Columns

Part Number	Description	Qty.
820-0055-B	ISOLUTE SLE+ 400 µL Sample Volume	50
820-0055-BG	ISOLUTE SLE+ 400 µL Sample Volume (Tabless)	50
820-0140-C	ISOLUTE SLE+ 1 mL Sample Volume	30
820-0140-CG	ISOLUTE SLE+ 1 mL Sample Volume (Tabless)	30
820-0290-D	ISOLUTE SLE+ 2 mL Sample Volume	20
820-0690-E	ISOLUTE SLE+ 5 mL Sample Volume	20
820-1420-F	ISOLUTE SLE+ 10 mL Sample Volume	16


Fixed Well Plates

Part Number	Description	Qty.
820-0200-P01	ISOLUTE SLE+ 200 µL Supported Liquid Extraction Plate	1
820-0400-P01	ISOLUTE SLE+ 400 µL Supported Liquid Extraction Plate	1

Array Wells

Part Number	Description	Qty.
820-0200-T	ISOLUTE SLE+ 200 µL Array Wells	100
820-0400-T	ISOLUTE SLE+ 400 µL Array Wells	100

ISOLUTE® SLE+ Supported Liquid Extraction Method Selection

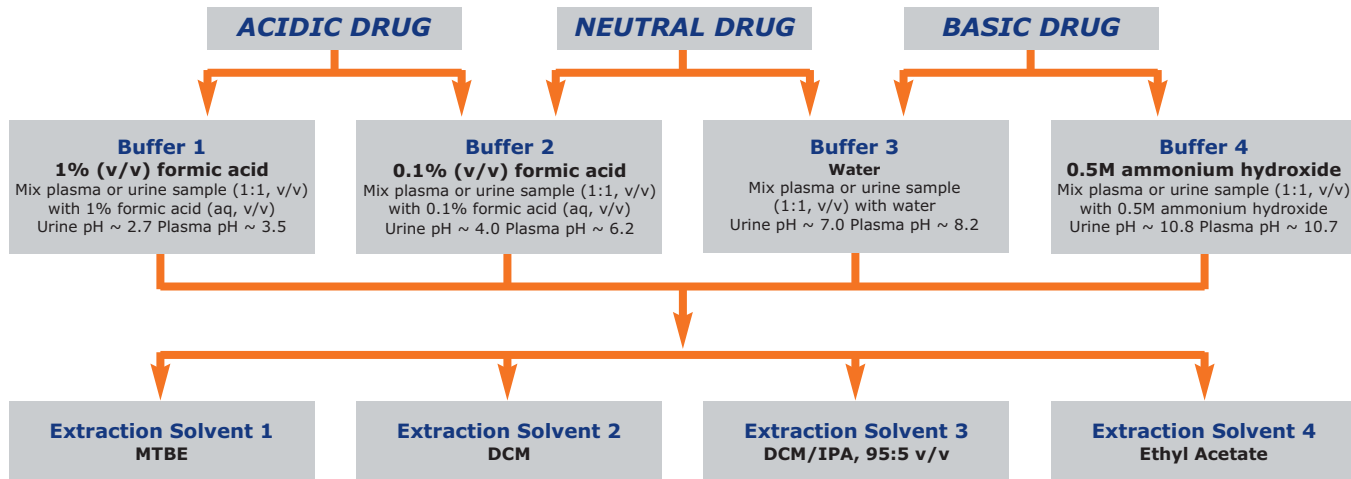

ISOLUTE® SLE+ User Guide

For further information please visit www.biotage.com to request a copy of our ISOLUTE® SLE+ User Guide






EVOLUTE® Sample Preparation Products

Advanced Polymeric Sorbents

EVOLUTE® ABN

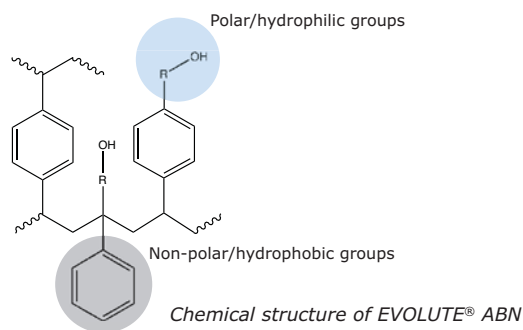


Chemical Description:

Water wettable polystyrene-divinylbenzene incorporating non-ionizable hydroxyl groups.

Average particle size	30 µm, 50 µm
Pore diameter	40 Å
Sorbent Type	Wettable non-polar sorbent with no secondary interactions

Application: EVOLUTE® ABN can be used to extract a diverse range of acidic, neutral and basic analytes from biological fluids and other aqueous matrices. Performance is not affected by drying the sorbent.

EVOLUTE® CX

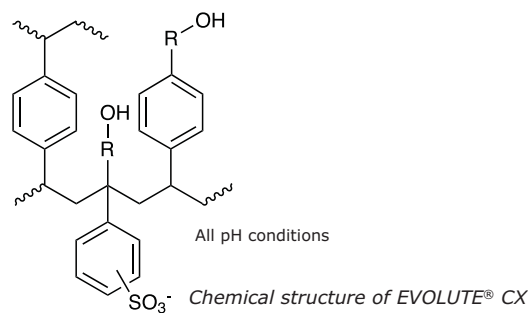

Chemical Description:

Sulfonic acid modified polystyrene-divinylbenzene incorporating non-ionizable hydroxyl groups.

Average particle size	30 µm, 50 µm
Pore diameter	40 Å
Sorbent Type	Mixed-mode non-polar/strong cation exchange
Exchange capacity	0.5 mmol/g

Application: EVOLUTE® CX can be used to extract basic analytes from biological fluids and other aqueous matrices. Performance is not affected by drying the sorbent.

EVOLUTE ABN 30 µm SPE Columns

Part Number	Description	Qty.
600-0001-A	EVOLUTE ABN 10 mg/1 mL	100
600-0002-A	EVOLUTE ABN 25 mg/1 mL	100
600-0002-AG	EVOLUTE ABN 25 mg/1 mL (Tabless)*	100
600-0002-B	EVOLUTE ABN 25 mg /3 mL	50
600-0002-H	EVOLUTE ABN 25 mg/10 mL XL	50

EVOLUTE ABN 50 µm SPE Columns

610-0005-B	EVOLUTE ABN 50 mg/3 mL	50
610-0005-BG	EVOLUTE ABN 50 mg/3 mL (Tabless)*	50
610-0010-B	EVOLUTE ABN 100 mg/3 mL	50
610-0010-H	EVOLUTE ABN 100 mg/10 mL XL	50
610-0020-B	EVOLUTE ABN 200 mg/3 mL	50
610-0020-C	EVOLUTE ABN 200 mg/6 mL	30
610-0050-C	EVOLUTE ABN 500 mg/6 mL	30

EVOLUTE EXPRESS ABN 30 µm Fixed Well Plates

600-0010-PX01	EVOLUTE EXPRESS ABN 10 mg Fixed Well Plate	1
600-0030-PX01	EVOLUTE EXPRESS ABN 30 mg Fixed Well Plate	1

EVOLUTE ABN 30 µm Pre-Assembled 96-Well Plates

600-0010-RP	EVOLUTE Array ABN 10 mg/1 mL Plate	1
600-0025-RP	EVOLUTE Array ABN 25 mg/1 mL Plate	1

EVOLUTE ABN 30 µm Array Loose Wells

600-0010-R	EVOLUTE Array ABN 10 mg Wells	100
600-0025-R	EVOLUTE Array ABN 25 mg Wells	100

EVOLUTE CX 30 µm SPE Columns

Part Number	Description	Qty.
601-0001-A	EVOLUTE CX 10 mg/1 mL	100
601-0001-AG	EVOLUTE CX 10 mg/1 mL (Tabless)*	100
601-0002-A	EVOLUTE CX 25 mg/1 mL	100
601-0002-AG	EVOLUTE CX 25 mg/1 mL (Tabless)*	100

EVOLUTE CX 50 µm SPE Columns

611-0005-B	EVOLUTE CX 50 mg/3 mL	50
611-0005-BG	EVOLUTE CX 50 mg/3 mL (Tabless)*	50
611-0010-B	EVOLUTE CX 100 mg/3 mL	50
611-0010-BG	EVOLUTE CX 100 mg/3 mL (Tabless)*	50
611-0010-H	EVOLUTE CX 100 mg/10 mL XL	50
611-0020-B	EVOLUTE CX 200 mg/3 mL	50
611-0020-C	EVOLUTE CX 200 mg/6 mL	30
611-0050-C	EVOLUTE CX 500 mg/6 mL	30
611-0050-CG	EVOLUTE CX 500 mg/6 mL (Tabless)*	30

EVOLUTE EXPRESS CX 30 µm Fixed Well Plates

601-0010-PX01	EVOLUTE EXPRESS CX 10 mg Fixed Well Plate	1
601-0030-PX01	EVOLUTE EXPRESS CX 30 mg Fixed Well Plate	1

EVOLUTE CX 30 µm Pre-Assembled 96-Well Plates

601-0010-RP	EVOLUTE Array CX 10 mg/1 mL Plate	1
601-0025-RP	EVOLUTE Array CX 25 mg/1 mL Plate	1


EVOLUTE CX 30 µm Array Loose Wells

601-0010-R	EVOLUTE Array CX 10 mg/1 mL Wells	100
601-0025-R	EVOLUTE Array CX 25 mg/1 mL Wells	100

EVOLUTE EXPRESS products can be used with the Load-Wash-Elute procedure. See www.biotage.com for more information.

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

EVOLUTE® AX

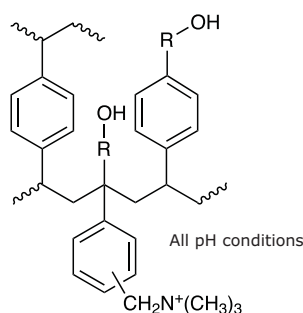

Chemical structure of EVOLUTE® AX

Chemical Description:

Quaternary amine modified polystyrene-divinylbenzene incorporating non-ionizable hydroxyl groups.

Average particle size	30 µm, 50 µm
Pore diameter	40 Å
Sorbent Type	Mixed-mode non-polar/strong anion exchange
Exchange capacity	0.7 mmol/g

Application: EVOLUTE® AX can be used to extract acidic analytes from biological fluids and other aqueous matrices. Performance is not affected by drying the sorbent.

EVOLUTE AX 30 µm SPE Columns

Part Number	Description	Qty.
603-0002-A	EVOLUTE AX 25 mg/1 mL	100
603-0002-AG	EVOLUTE AX 25 mg/1 mL (Tabless)*	100

EVOLUTE AX 50 µm SPE Columns

613-0005-B	EVOLUTE AX 50 mg/3 mL	50
613-0010-B	EVOLUTE AX 100 mg/3 mL	50
613-0010-BG	EVOLUTE AX 100 mg/3 mL (Tabless)*	50
613-0010-H	EVOLUTE AX 100 mg/10 mL XL	50
613-0020-C	EVOLUTE AX 200 mg/6 mL	30
613-0050-C	EVOLUTE AX 500 mg/6 mL	30
613-0050-CG	EVOLUTE AX 500 mg/6 mL (Tabless)*	30


EVOLUTE EXPRESS AX 30 µm Fixed Well Plates

603-0010-PX01	EVOLUTE EXPRESS AX 10 mg Fixed Well Plate	1
603-0030-PX01	EVOLUTE EXPRESS AX 30 mg Fixed Well Plate	1

EVOLUTE AX 30 µm Array Loose Wells

603-0025-R	EVOLUTE Array AX 25 mg/1 mL Wells	100
------------	-----------------------------------	-----

EVOLUTE® WCX

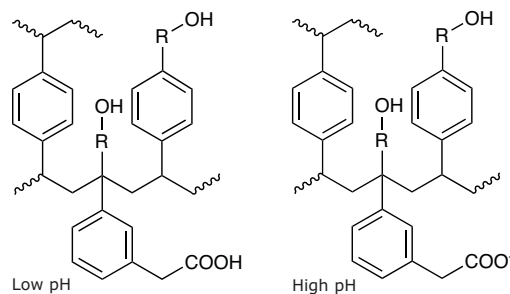

Chemical structure of EVOLUTE® WCX

Chemical Description:

Carboxylic acid modified polystyrene-divinylbenzene incorporating non-ionizable hydroxyl groups.

Average particle size	30 µm, 50 µm
Pore diameter	40 Å
Sorbent Type	Mixed-mode non-polar/weak cation exchange (pK _a ~5)
Exchange capacity	0.4 mmol/g

Application: EVOLUTE® WCX can be used to extract strongly basic analytes (e.g. quaternary amines) from biological fluids and other aqueous matrices. Performance is not affected by drying the sorbent.

EVOLUTE WCX 30 µm SPE Columns

Part Number	Description	Qty.
602-0002-A	EVOLUTE WCX 25 mg/1 mL	100

EVOLUTE WCX 50 µm SPE Columns

612-0005-B	EVOLUTE WCX 50 mg/3 mL	50
612-0010-B	EVOLUTE WCX 100 mg/3 mL	50
612-0010-H	EVOLUTE WCX 100 mg/3 mL XL	50
612-0020-C	EVOLUTE WCX 200 mg/6 mL	30
612-0050-C	EVOLUTE WCX 500 mg/6 mL	30
612-0050-CG	EVOLUTE WCX 500 mg/6 mL (Tabless)*	30

EVOLUTE EXPRESS WCX 30 µm Fixed Well Plates

602-0010-PX01	EVOLUTE EXPRESS WCX 10 mg Fixed Well Plate	1
602-0030-PX01	EVOLUTE EXPRESS WCX 30 mg Fixed Well Plate	1


EVOLUTE WCX 30 µm Array Loose Wells

602-0025-R	EVOLUTE Array WCX 25 mg/1 mL Wells	100
------------	------------------------------------	-----

EVOLUTE EXPRESS products can be used with the Load-Wash-Elute procedure. See www.biotage.com for more information.

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

EVOLUTE® WAX

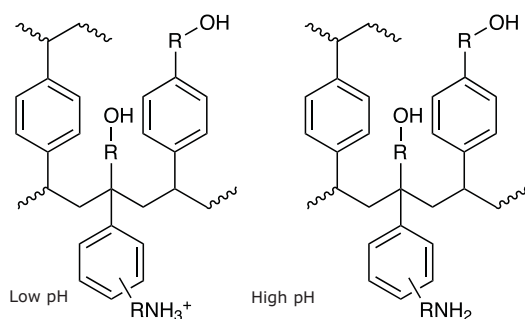

Chemical structure of EVOLUTE® WAX

Chemical Description:

Primary-secondary amine modified polystyrene-divinylbenzene incorporating non-ionizable hydroxyl groups.

Average particle size	30 μm , 50 μm
Pore diameter	40 Å
Sorbent Type	Mixed-mode non-polar/weak cation exchange ($\text{pK}_a \sim 5$)
Exchange capacity	0.3 mmol/g, 0.7 mmol/g

Application: EVOLUTE WAX can be used to extract strongly acidic analytes (e.g. sulfonic acids) from biological fluids and other aqueous matrices. Performance is not affected by drying the sorbent.

EVOLUTE WAX 30 μm SPE Columns

Part Number	Description	Qty.
604-0002-A	EVOLUTE WAX 25 mg/1 mL	100

EVOLUTE WAX 50 μm SPE Columns

614-0005-B	EVOLUTE WAX 50 mg/3 mL	50
614-0010-B	EVOLUTE WAX 100 mg/3 mL	50
614-0010-H	EVOLUTE WAX 100 mg/10 mL XL	50
614-0020-C	EVOLUTE WAX 200 mg/6 mL	30
614-0050-C	EVOLUTE WAX 500 mg/6 mL	30

EVOLUTE EXPRESS WAX 30 μm Fixed Well Plates

Part Number	Description	Qty.
604-0010-PX01	EVOLUTE EXPRESS WAX 10 mg Fixed Well Plate	1
604-0030-PX01	EVOLUTE EXPRESS WAX 30 mg Fixed Well Plate	1

EVOLUTE WAX 30 μm Array Loose Wells

Part Number	Description	Qty.
604-0025-R	EVOLUTE Array WAX 25 mg/1mL Wells	100

EVOLUTE® EXPRESS Sorbent Selection Plate


The EVOLUTE® EXPRESS Sorbent Selection plate contains the four mixed-mode EVOLUTE sorbents CX, WCX, AX and WAX on one plate.

	EVOLUTE CX			EVOLUTE WAX			EVOLUTE AX			EVOLUTE WCX		
	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												
	Method 1						Method 2					

Part Number	Description	Qty.
650-0010-PX01	EVOLUTE EXPRESS 10 mg Sorbent Selection Plate	1
650-0030-PX01	EVOLUTE EXPRESS 30 mg Sorbent Selection Plate	1

EVOLUTE® User Guide

For further information please visit www.biotage.com to request a copy of our EVOLUTE User Guide.

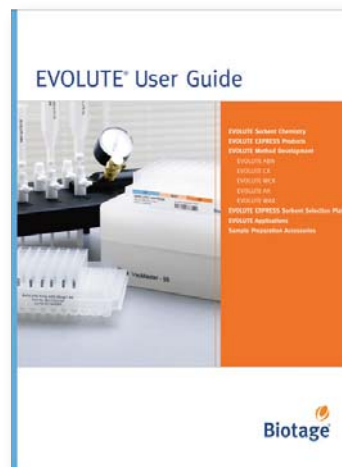


EVOLUTE EXPRESS products can be used with the Load-Wash-Elute procedure. See www.biotage.com for more information.



ISOLUTE[®] SPE Columns and 96-Well Plates

Sorbents and Formats for a Wide Range
of Sample Preparation Applications

ISOLUTE® ENV+

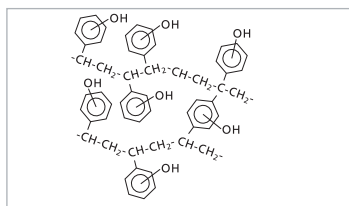


Chemical structure of ISOLUTE ENV+, a hydroxylated polystyrene divinylbenzene co-polymer

Average particle size	110 µm
Pore diameter	800 Å
Sorbent Type	Non-Polar

Application: Extraction of very polar compounds that are not retained by C8 and C18 non-polar silica based sorbents.

ISOLUTE® 101

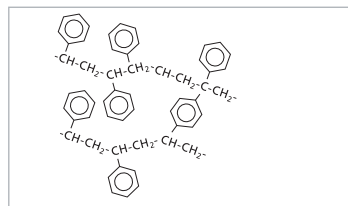

Chemical structure of ISOLUTE 101, an unmodified polystyrene divinylbenzene co-polymer

Average particle size	65 µm
Pore diameter	100 Å
Sorbent Type	Non-Polar

Application: Extraction of polar compounds from aqueous sample matrices.

ISOLUTE ENV+ SPE Columns

Part Number	Description	Qty.
915-0001-A	ISOLUTE ENV+ 10 mg/1 mL	100
915-0002-A	ISOLUTE ENV+ 25 mg/1 mL	100
915-0002-G	ISOLUTE ENV+ 25 mg/10 mL	50
915-0005-A	ISOLUTE ENV+ 50 mg/1 mL	100
915-0005-AG	ISOLUTE ENV+ 50 mg/1 mL (Tabless)*	100
915-0005-B	ISOLUTE ENV+ 50 mg/3 mL	50
915-0005-G	ISOLUTE ENV+ 50 mg/10 mL	50
915-0010-A	ISOLUTE ENV+ 100 mg/1 mL	100
915-0010-B	ISOLUTE ENV+ 100 mg/3 mL	50
915-0010-BG	ISOLUTE ENV+ 100 mg/3 mL (Tabless)*	50
915-0010-C	ISOLUTE ENV+ 100 mg/6 mL	30
915-0010-H	ISOLUTE ENV+ 100 mg/10 mL	50
915-0020-B	ISOLUTE ENV+ 200 mg/3 mL	50
915-0020-C	ISOLUTE ENV+ 200 mg/6 mL	30
915-0020-CD	ISOLUTE ENV+ 200 mg/6 mL (Depth filter)	30
915-0020-CG	ISOLUTE ENV+ 200 mg/6 mL (Tabless)*	30
915-0050-B	ISOLUTE ENV+ 500 mg/3 mL	50
915-0050-BG	ISOLUTE ENV+ 500 mg/3 mL (Tabless)*	50
915-0050-C	ISOLUTE ENV+ 500 mg/6 mL	30
915-0050-D	ISOLUTE ENV+ 500 mg/15 mL	20
915-0050-L	ISOLUTE ENV+ 500 mg/6 mL (Glass)	30
915-0100-C	ISOLUTE ENV+ 1 g/6 mL	30
915-0100-CG	ISOLUTE ENV+ 1 g/6 mL (Tabless)*	30
915-0100-E	ISOLUTE ENV+ 1 g/15 mL	20

ISOLUTE ENV+ Fixed Well Plates

Part Number	Description	Qty.
915-0010-P01	ISOLUTE-96 ENV+ 10 mg plate	1
915-0025-P01	ISOLUTE-96 ENV+ 25 mg plate	1
915-0050-P01	ISOLUTE-96 ENV+ 50 mg plate	1

ISOLUTE Array ENV+ Loose Wells

Part Number	Description	Qty.
915-0010-R	ISOLUTE Array ENV+ 10 mg/1 mL wells	100
915-0025-R	ISOLUTE Array ENV+ 25 mg/1 mL wells	100

Pre-assembled ISOLUTE Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® ENV+

TN109: Method Development in Solid Phase Extraction using ISOLUTE ENV+ SPE Columns for the Extraction of Aqueous Samples

ISOLUTE 101 SPE Columns


Part Number	Description	Qty.
101-0010-B	ISOLUTE 101 100 mg/3 mL	50
101-0020-B	ISOLUTE 101 200 mg/3 mL	50
101-0020-C	ISOLUTE 101 200 mg/6 mL	30
101-0050-B	ISOLUTE 101 500 mg/3 mL	50
101-0050-C	ISOLUTE 101 500 mg/6 mL	30

Support Documents for ISOLUTE® 101

TN119: Method Development in Solid Phase Extraction using ISOLUTE 101 SPE Columns for the Extraction of Aqueous Samples

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

ISOLUTE® C18

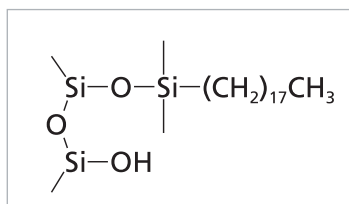

Chemical structure of C18 Octadecyl (non-endcapped) silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Most commonly used C18 sorbent for the extraction of acidic, neutral and basic compounds from aqueous matrices. Secondary silanol or ionic interactions can be used to enhance extract purity and method robustness for basic compounds.

ISOLUTE C18 SPE Columns

Part Number	Description	Qty.
220-0002-A	ISOLUTE C18 25 mg/1 mL	100
220-0002-AG	ISOLUTE C18 25 mg/1 mL (Tabless)*	100
220-0005-A	ISOLUTE C18 50 mg/1 mL	100
220-0005-G	ISOLUTE C18 50 mg/10 mL	50
220-0010-A	ISOLUTE C18 100 mg/1 mL	100
220-0010-B	ISOLUTE C18 100 mg/3 mL	50
220-0010-BG	ISOLUTE C18 100 mg/3 mL (Tabless)*	50
220-0010-G	ISOLUTE C18 100 mg/10 mL	50
220-0020-B	ISOLUTE C18 200 mg/3 mL	50
220-0020-C	ISOLUTE C18 200 mg/6 mL	30
220-0020-H	ISOLUTE C18 200 mg/10 mL	50
220-0050-B	ISOLUTE C18 500 mg/3 mL	50
220-0050-BG	ISOLUTE C18 500 mg/3 mL (Tabless)*	50
220-0050-C	ISOLUTE C18 500 mg/6 mL	30
220-0050-H	ISOLUTE C18 500 mg/10 mL	50
220-0050-L	ISOLUTE C18 500 mg/6 mL (Glass)	30
220-0100-B	ISOLUTE C18 1 g/3 mL	50
220-0100-C	ISOLUTE C18 1 g/6 mL	30
220-0100-CG	ISOLUTE C18 1 g/6 mL (Tabless)*	30
220-0200-C	ISOLUTE C18 2 g/6 mL	30
220-0200-D	ISOLUTE C18 2 g/15 mL	20
220-0500-E	ISOLUTE C18 5 g/25 mL	20
220-1000-F	ISOLUTE C18 10 g/70 mL	16

ISOLUTE-96 C18 Fixed Well Plates

Part Number	Description	Qty.
220-0025-P01	ISOLUTE-96 C18 25 mg plate	1
220-0050-P01	ISOLUTE-96 C18 50 mg plate	1
220-0100-P01	ISOLUTE-96 C18 100 mg plate	1

ISOLUTE Array C18 Loose Wells

Part number	Description	Qty.
220-0025-R	ISOLUTE Array C18 25 mg/1 mL wells	100
220-0050-R	ISOLUTE Array C18 50 mg/1 mL wells	100
220-0100-T	ISOLUTE Array C18 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.


Support Documents for ISOLUTE® C18

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-endcapped Sorbents

TN126: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Drugs from Biological Fluid Samples

ISOLUTE® C18(EC)

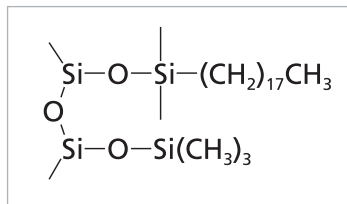

Chemical structure of C18 Octadecyl (endcapped) silane and trimethyl silyl group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction of polar compounds from aqueous sample matrices.

ISOLUTE C18(EC) SPE Columns

Part Number	Description	Qty.
221-0002-A	ISOLUTE C18(EC) 25 mg/1 mL	100
221-0005-A	ISOLUTE C18(EC) 50 mg/1 mL	100
221-0005-G	ISOLUTE C18(EC) 50 mg/10 mL	50
221-0010-A	ISOLUTE C18(EC) 100 mg/1 mL	100
221-0010-B	ISOLUTE C18(EC) 100 mg/3 mL	50
221-0010-C	ISOLUTE C18(EC) 100 mg/6 mL	30
221-0010-G	ISOLUTE C18(EC) 100 mg/10 mL	50
221-0020-B	ISOLUTE C18(EC) 200 mg/3 mL	50
221-0020-C	ISOLUTE C18(EC) 200 mg/6 mL	30
221-0020-H	ISOLUTE C18(EC) 200 mg/10 mL	50
221-0050-B	ISOLUTE C18(EC) 500 mg/3 mL	50
221-0050-BS	ISOLUTE C18(EC) 500 mg/3 mL**	50
221-0050-C	ISOLUTE C18(EC) 500 mg/6 mL	30
221-0050-CD	ISOLUTE C18(EC) 500 mg/6 mL (Depth filter)	30
221-0050-H	ISOLUTE C18(EC) 500 mg/10 mL	50
221-0100-B	ISOLUTE C18(EC) 1 g/3 mL	50
221-0100-C	ISOLUTE C18(EC) 1 g/6 mL	30
221-0200-C	ISOLUTE C18(EC) 2 g/6 mL	30
221-0200-D	ISOLUTE C18(EC) 2 g/15 mL	20
221-0500-E	ISOLUTE C18(EC) 5 g/25 mL	20
221-1000-F	ISOLUTE C18(EC) 10 g/70 mL	16

ISOLUTE-96 C18(EC) Fixed Well Plates

Part Number	Description	Qty.
221-0025-P01	ISOLUTE-96 C18(EC) 25 mg plate	1
221-0050-P01	ISOLUTE-96 C18(EC) 50 mg plate	1
221-0100-P01	ISOLUTE-96 C18(EC) 100 mg plate	1

ISOLUTE Array C18(EC) Loose Wells

Part Number	Description	Qty.
221-0025-R	ISOLUTE Array C18(EC) 25 mg/1 mL wells	100
221-0050-R	ISOLUTE Array C18(EC) 50 mg/1 mL wells	100
221-0100-T	ISOLUTE Array C18(EC) 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.


Support Documents ISOLUTE C18 (EC)

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

**Column contains stainless steel frits.

ISOLUTE® MFC18

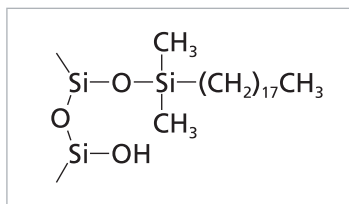

Chemical structure of monofunctional C18 silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	125 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using both non-polar interactions (for acidic, neutral and basic compounds) and readily accessible secondary silanol or ionic interactions (for basic compounds only). Not the ideal choice when working at extreme pH, C18 is the preferred option for these methods.

ISOLUTE MFC18 SPE Columns

Part Number	Description	Qty.
240-0005-A	ISOLUTE MFC18 50 mg/1 mL	100
240-0005-G	ISOLUTE MFC18 50 mg/10 mL	50
240-0010-A	ISOLUTE MFC18 100 mg/1 mL	100
240-0010-B	ISOLUTE MFC18 100 mg/3 mL	50
240-0010-G	ISOLUTE MFC18 100 mg/10 mL	50
240-0020-B	ISOLUTE MFC18 200 mg/3 mL	50
240-0020-H	ISOLUTE MFC18 200 mg/10 mL	50
240-0050-B	ISOLUTE MFC18 500 mg/3 mL	50
240-0050-C	ISOLUTE MFC18 500 mg/6 mL	30
240-0050-H	ISOLUTE MFC18 500 mg/10 mL	50
240-0100-C	ISOLUTE MFC18 1 g/6 mL	30
240-0200-D	ISOLUTE MFC18 2 g/15 mL	20
240-0500-E	ISOLUTE MFC18 5 g/25 mL	20

ISOLUTE-96 MFC18 Fixed Well Plates


Part Number	Description	Qty.
240-0025-P01	ISOLUTE-96 MFC18 25 mg plate	1

Support Documents for ISOLUTE® MFC18

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

ISOLUTE® C8

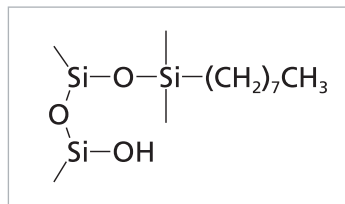

Chemical structure of C8 silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non Polar

Application: Most commonly used C8 sorbent for the extraction of acidic, neutral and basic compounds from aqueous matrices. Secondary silanol or ionic interactions can be used to enhance extract purity and method robustness for basic compounds.

ISOLUTE C8 SPE Columns

Part Number	Description	Qty.
290-0002-A	ISOLUTE C8 25 mg/1 mL	100
290-0005-A	ISOLUTE C8 50 mg/1 mL	100
290-0010-A	ISOLUTE C8 100 mg/1 mL	100
290-0010-B	ISOLUTE C8 100 mg/3 mL	50
290-0010-C	ISOLUTE C8 100 mg/6 mL	30
290-0020-B	ISOLUTE C8 200 mg/3 mL	50
290-0020-H	ISOLUTE C8 200 mg/10 mL	50
290-0050-B	ISOLUTE C8 500 mg/3 mL	50
290-0050-C	ISOLUTE C8 500 mg/6 mL	30
290-0100-C	ISOLUTE C8 1 g/6 mL	30
290-0200-D	ISOLUTE C8 2 g/15 mL	20

ISOLUTE-96 C8 Fixed Well Plates

Part Number	Description	Qty.
290-0025-P01	ISOLUTE-96 C8 25 mg plate	1
290-0050-P01	ISOLUTE-96 C8 50 mg plate	1
290-0100-P01	ISOLUTE-96 C8 100 mg plate	1

ISOLUTE Array C8 Loose Wells

Part Number	Description	Qty.
290-0025-R	ISOLUTE Array C8 25 mg/1 mL wells	100
290-0050-R	ISOLUTE Array C8 50 mg/1 mL wells	100
290-0100-T	ISOLUTE Array C8 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.


Support Documents for ISOLUTE® C8

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

TN126: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Biological Fluid Samples

ISOLUTE® C8(EC)

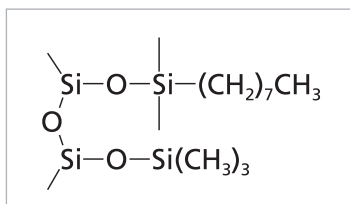

Chemical structure of C8 silane and trimethyl silyl group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using non-polar interactions (for acidic, neutral and basic compounds).

ISOLUTE C8(EC) SPE Columns

Part Number	Description	Qty.
291-0010-A	ISOLUTE C8(EC) 100 mg/1 mL	100
291-0010-B	ISOLUTE C8(EC) 100 mg/3 mL	50
291-0010-G	ISOLUTE C8(EC) 100 mg/10 mL	50
291-0020-B	ISOLUTE C8(EC) 200 mg/3 mL	50
291-0020-H	ISOLUTE C8(EC) 200 mg/10 mL	50
291-0050-B	ISOLUTE C8(EC) 500 mg/3 mL	50
291-0050-C	ISOLUTE C8(EC) 500 mg/6 mL	30
291-0100-C	ISOLUTE C8(EC) 1 g/6 mL	30
291-0200-D	ISOLUTE C8(EC) 2 g/15 mL	20
291-0500-E	ISOLUTE C8(EC) 5 g/25 mL	20

ISOLUTE-96 C8(EC) Fixed Well Plates

Part Number	Description	Qty.
291-0025-P01	ISOLUTE-96 C8(EC) 25 mg plate	1
291-0050-P01	ISOLUTE-96 C8(EC) 50 mg plate	1
291-0100-P01	ISOLUTE-96 C8(EC) 100 mg plate	1

ISOLUTE Array C8(EC) Loose Wells


Part Number	Description	Qty.
291-0025-R	ISOLUTE Array C8(EC) 25 mg/1 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® C8(EC)

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® C4

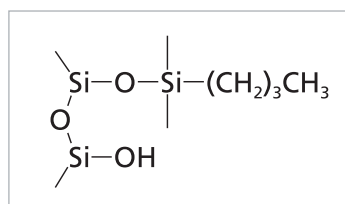

Chemical structure of C4 silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using non-polar interactions (for acidic, neutral and basic compounds). Secondary silanol or ionic interactions can be used to enhance extract purity and method robustness for basic compounds.

ISOLUTE C4 SPE Columns

Part Number	Description	Qty.
390-0010-A	ISOLUTE C4 100 mg/1 mL	100
390-0010-B	ISOLUTE C4 100 mg/3 mL	50
390-0020-B	ISOLUTE C4 200 mg/3 mL	50
390-0050-B	ISOLUTE C4 500 mg/3 mL	50
390-0050-C	ISOLUTE C4 500 mg/6 mL	30
390-0100-C	ISOLUTE C4 1 g/6 mL	30

ISOLUTE-96 C4 Fixed Well Plates

Part Number	Description	Qty.
390-0025-P01	ISOLUTE-96 C4 25 mg plate	1

ISOLUTE Array C4 Loose Wells

Part Number	Description	Qty.
390-0025-R	ISOLUTE Array C4 25 mg/1 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® C4

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

ISOLUTE® C2

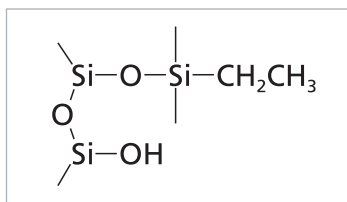

Chemical structure of C2 silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non Polar

Application: Most commonly used C2 sorbent for the extraction of acidic, neutral and basic compounds from aqueous matrices. Secondary silanol or ionic interactions can be used to enhance extract purity and method robustness for basic compounds.

ISOLUTE C2 SPE Columns

Part Number	Description	Qty.
320-0002-A	ISOLUTE C2 25 mg/1 mL	100
320-0005-A	ISOLUTE C2 50 mg/1 mL	100
320-0010-A	ISOLUTE C2 100 mg/1 mL	100
320-0010-B	ISOLUTE C2 100 mg/3 mL	50
320-0010-G	ISOLUTE C2 100 mg/10 mL	50
320-0020-B	ISOLUTE C2 200 mg/3 mL	50
320-0050-B	ISOLUTE C2 500 mg/3 mL	50
320-0050-C	ISOLUTE C2 500 mg/6 mL	30
320-0100-C	ISOLUTE C2 1 g/6 mL	30
320-0200-D	ISOLUTE C2 2 g/15 mL	20

ISOLUTE-96 C2 Fixed Well Plates

Part Number	Description	Qty.
320-0025-P01	ISOLUTE-96 C2 25 mg plate	1
320-0050-P01	ISOLUTE-96 C2 50 mg plate	1
320-0100-P01	ISOLUTE-96 C2 100 mg plate	1

ISOLUTE Array C2 Loose Wells

Part Number	Description	Qty.
320-0025-R	ISOLUTE Array C2 25 mg/1 mL wells	100
320-0050-R	ISOLUTE Array C2 50 mg/1 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.


Support Documents for ISOLUTE® C2

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

TN126: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Biological Fluid Samples

ISOLUTE® C2(EC)

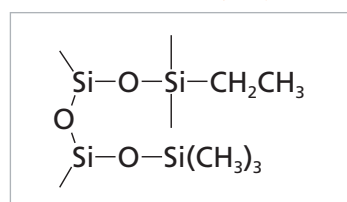

Chemical structure of C2 silane and trimethyl silyl group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using non-polar interactions (for acidic, neutral and basic compounds).

ISOLUTE C2(EC) SPE Columns

Part Number	Description	Qty.
321-0005-A	ISOLUTE C2(EC) 50 mg/1 mL	100
321-0010-A	ISOLUTE C2(EC) 100 mg/1 mL	100
321-0010-B	ISOLUTE C2(EC) 100 mg/3 mL	50
321-0010-G	ISOLUTE C2(EC) 100 mg/10 mL	50
321-0020-B	ISOLUTE C2(EC) 200 mg/3 mL	50
321-0050-B	ISOLUTE C2(EC) 500 mg/3 mL	50
321-0050-C	ISOLUTE C2(EC) 500 mg/6 mL	30
321-0050-H	ISOLUTE C2(EC) 500 mg/10 mL	50
321-0100-C	ISOLUTE C2(EC) 1 g/6 mL	30

ISOLUTE-96 C2(EC) Fixed Well Plates

Part Number	Description	Qty.
321-0025-P01	ISOLUTE-96 C2(EC) 25 mg plate	1

ISOLUTE Array C2(EC) Loose Wells


Part Number	Description	Qty.
321-0025-R	ISOLUTE Array C2(EC) 25 mg/1 mL wells	100
321-0050-R	ISOLUTE Array C2(EC) 50 mg/1 mL wells	100
321-0100-T	ISOLUTE Array C2(EC) 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® C2(EC)

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® PH

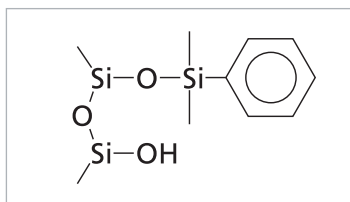

Chemical structure of phenyl silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using both non-polar interactions (for acidic, neutral and basic compounds) and secondary silanol or ionic interactions (for basic compounds only). This sorbent exhibits a different selectivity compared with C18 and C8 phases when both aromatic and non-aromatic compounds are being extracted.

ISOLUTE PH SPE Columns

Part Number	Description	Qty.
360-0002-A	ISOLUTE PH 25 mg/1 mL	100
360-0005-A	ISOLUTE PH 50 mg/1 mL	100
360-0010-A	ISOLUTE PH 100 mg/1 mL	100
360-0010-B	ISOLUTE PH 100 mg/3 mL	50
360-0050-B	ISOLUTE PH 500 mg/3 mL	50
360-0050-C	ISOLUTE PH 500 mg/6 mL	30
360-0100-C	ISOLUTE PH 1 g/6 mL	30

ISOLUTE-96 PH Fixed Well Plates

Part Number	Description	Qty.
360-0025-P01	ISOLUTE-96 PH 25 mg plate	1
360-0050-P01	ISOLUTE-96 PH 50 mg plate	1
360-0100-P01	ISOLUTE-96 PH 100 mg plate	1

ISOLUTE Array PH Loose Wells

Part Number	Description	Qty.
360-0025-R	ISOLUTE Array PH 25 mg/1 mL wells	100
360-0050-R	ISOLUTE Array PH 50 mg/1 mL wells	100
360-0100-T	ISOLUTE Array PH 100 mg/2 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® PH

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

ISOLUTE® CN

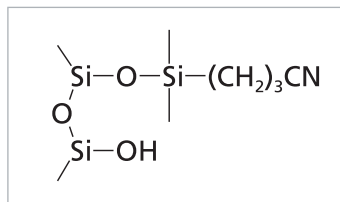

Chemical structure of cyanopropyl silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Non-Polar

Application: Extraction from aqueous matrix using both non-polar interactions (for acidic, neutral and basic compounds) and secondary silanol or ionic interactions (for basic compounds only). Can also be used in polar SPE mode.

ISOLUTE CN SPE Columns

Part Number	Description	Qty.
420-0010-A	ISOLUTE CN 100 mg/1 mL	100
420-0010-B	ISOLUTE CN 100 mg/3 mL	50
420-0020-B	ISOLUTE CN 200 mg/3 mL	50
420-0050-B	ISOLUTE CN 500 mg/3 mL	50
420-0050-C	ISOLUTE CN 500 mg/6 mL	30
420-0100-C	ISOLUTE CN 1 g/6 mL	30

ISOLUTE-96 CN Fixed Well Plates

Part Number	Description	Qty.
420-0025-P01	ISOLUTE-96 CN 25 mg plate	1
420-0100-P01	ISOLUTE-96 CN 100 mg plate	1

ISOLUTE Array CN Loose Wells

Part Number	Description	Qty.
420-0025-R	ISOLUTE Array CN 25 mg/1 mL wells	100
420-0100-T	ISOLUTE Array CN 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.


Support Documents for ISOLUTE® CN

TN101: Method Development in Solid Phase Extraction using Non-polar ISOLUTE® SPE Columns for the Extraction of Aqueous Samples

TN102: Method Development in Solid Phase Extraction using Polar ISOLUTE® SPE Columns for the Extraction of Non-aqueous Samples

TN112: General Approach to the Extraction of Basic Drugs from Biological Fluids using Non-polar Non-encapped Sorbents

ISOLUTE® HCX

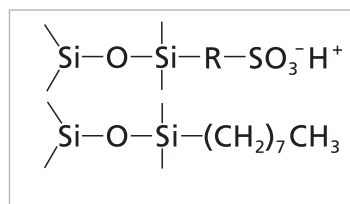

ISOLUTE® HCX combines C8 (Octyl) and sulfonic acid functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: The first choice sorbent for extracting drugs of abuse from biological fluid samples.

ISOLUTE HCX SPE Columns

Part Number	Description	Qty.
902-0002-A	ISOLUTE HCX 25 mg/1 mL	100
902-0002-AG	ISOLUTE HCX 25 mg/1 mL (Tabless)*	100
902-0005-A	ISOLUTE HCX 50 mg/1 mL	100
902-0010-A	ISOLUTE HCX 100 mg/1 mL	100
902-0010-B	ISOLUTE HCX 100 mg/3 mL	50
902-0013-A	ISOLUTE HCX 130 mg/1 mL	100
902-0013-B	ISOLUTE HCX 130 mg/3 mL	50
902-0013-BG	ISOLUTE HCX 130 mg/3 mL (Tabless)*	50
902-0013-C	ISOLUTE HCX 130 mg/6 mL	30
902-0013-CG	ISOLUTE HCX 130 mg/6 mL (Tabless)*	30
902-0013-H	ISOLUTE HCX 130 mg/10 mL	50
902-0020-B	ISOLUTE HCX 200 mg/3 mL	50
902-0020-H	ISOLUTE HCX 200 mg/10 mL	50
902-0030-B	ISOLUTE HCX 300 mg/3 mL	50
902-0030-C	ISOLUTE HCX 300 mg/6 mL	30
902-0030-CG	ISOLUTE HCX 300 mg/6 mL (Tabless)*	30
902-0030-H	ISOLUTE HCX 300 mg/10 mL	30
902-0050-C	ISOLUTE HCX 500 mg/6 mL	30

ISOLUTE-96 HCX Fixed Well Plates

Part Number	Description	Qty.
902-0025-P01	ISOLUTE-96 HCX 25 mg plate	1
902-0050-P01	ISOLUTE-96 HCX 50 mg plate	1
902-0100-P01	ISOLUTE-96 HCX 100 mg plate	1

ISOLUTE Array HCX Loose Wells

Part Number	Description	Qty.
902-0025-R	ISOLUTE Array HCX 25 mg/1 mL wells	100
902-0050-R	ISOLUTE Array HCX 50 mg/1 mL wells	100
902-0050-T	ISOLUTE Array HCX 50 mg/2 mL wells	100
902-0100-T	ISOLUTE Array HCX 100 mg/2 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® HCX

TN116: Generic Method for the Extraction of Basic Drugs from Biological Fluids using ISOLUTE® Mixed-mode SPE Columns and 96-well Plates

TN125: Method Development in Solid Phase Extraction using ISOLUTE® HCX for the Extraction of Drugs from Biological Fluid Samples

ISOLUTE® HCX-3

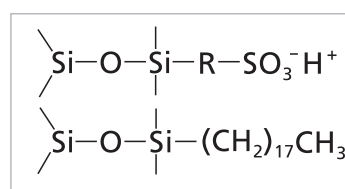

ISOLUTE® HCX-3 combines C18 (Octadecyl) and sulfonic acid functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: Extraction of basic analytes from aqueous matrix using dual non-polar and strong cation exchange interactions. A good alternative to HCX for basic compounds that require more retentive non-polar character from the mixed-mode sorbent.

ISOLUTE HCX-3 SPE Columns

Part Number	Description	Qty.
905-0002-A	ISOLUTE HCX-3 25 mg/1 mL	100
905-0005-A	ISOLUTE HCX-3 50 mg/1 mL	100
905-0010-A	ISOLUTE HCX-3 100 mg/1 mL	100

ISOLUTE-96 HCX-3 Fixed Well Plates

Part Number	Description	Qty.
905-0025-P01	ISOLUTE-96 HCX-3 25 mg plate	1
905-0100-P01	ISOLUTE-96 HCX-3 100 mg plate	1

ISOLUTE Array HCX-3 Loose Wells


Part Number	Description	Qty.
905-0025-R	ISOLUTE Array HCX-3 25 mg/1 mL wells	100
905-0100-T	ISOLUTE Array HCX-3 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® HCX-3

TN116: Generic Method for the Extraction of Basic Drugs from Biological Fluids using ISOLUTE® Mixed-mode SPE Columns and 96-Well Plates

ISOLUTE® HCX-5

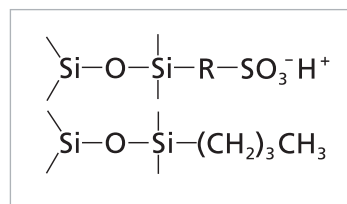

ISOLUTE® HCX-5 combines C4 (Butyl) and sulfonic acid functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: HCX-5 provides the cleanest extract of all the mixed-mode sorbents. Ideal choice where the basic analyte to be extracted has sufficient non-polar character to be well retained by the C4 non-polar component of the mixed-mode sorbent.

ISOLUTE HCX-5 SPE Columns

Part Number	Description	Qty.
906-0002-A	ISOLUTE HCX-5 25 mg/1 mL	100
906-0005-A	ISOLUTE HCX-5 50 mg/1 mL	100
906-0010-A	ISOLUTE HCX-5 100 mg/1 mL	100
906-0010-G	ISOLUTE HCX-5 100 mg/10 mL	50
906-0013-H	ISOLUTE HCX-5 130 mg/10 mL	50

ISOLUTE-96 HCX-5 Fixed Well Plates

Part Number	Description	Qty.
906-0025-P01	ISOLUTE-96 HCX-5 25 mg plate	1
906-0100-P01	ISOLUTE-96 HCX-5 100 mg plate	1

ISOLUTE Array HCX-5 Loose Wells


Part Number	Description	Qty.
906-0025-R	ISOLUTE Array HCX-5 25 mg/1 mL wells	100
906-0100-T	ISOLUTE Array HCX-5 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® HCX-5

TN116: Generic Method for the Extraction of Basic Drugs from Biological Fluids using ISOLUTE® Mixed-mode SPE Columns and 96-well Plates

ISOLUTE® HCX-Q

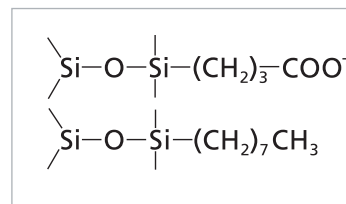

ISOLUTE® HCX-Q combines C8 (Octyl) and carboxylic acid functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: Extraction of quaternary amine and polybasic analytes from aqueous matrix using dual non-polar and weak cation exchange interactions.

ISOLUTE HCX-Q SPE Columns

Part Number	Description	Qty.
986-0002-A	ISOLUTE HCX-Q 25 mg/1 mL	100
986-0005-A	ISOLUTE HCX-Q 50 mg/1 mL	100
986-0010-A	ISOLUTE HCX-Q 100 mg/1 mL	100

ISOLUTE-96 HCX-Q Fixed Well Plates

Part Number	Description	Qty.
986-0025-P01	ISOLUTE-96 HCX-Q 25 mg plate	1
986-0100-P01	ISOLUTE-96 HCX-Q 100 mg plate	1

ISOLUTE Array HCX-Q Loose Wells


Part Number	Description	Qty.
986-0025-R	ISOLUTE Array HCX-Q 25 mg/1 mL wells	100
986-0100-T	ISOLUTE Array HCX-Q 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents ISOLUTE® HCX-Q

TN129: Generic Method for the Extraction of Quaternary Amine and Polybasic Drugs from Biological Fluids using ISOLUTE® HCX-Q SPE Columns and 96-well Plates

ISOLUTE® HAX

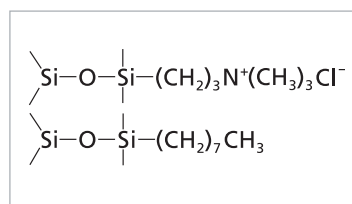

ISOLUTE® HAX combines C8 (Octyl) and quaternary amine functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: Extraction of acidic analytes from aqueous matrix using dual non-polar and strong anion exchange interactions. Suitable for a broad range of acidic compounds.

ISOLUTE HAX SPE Columns

Part Number	Description	Qty.
903-0005-A	ISOLUTE HAX 50 mg/1 mL	100
903-0010-A	ISOLUTE HAX 100 mg/1 mL	100
903-0020-B	ISOLUTE HAX 200 mg/3 mL	50
903-0020-C	ISOLUTE HAX 200 mg/6 mL	30
903-0020-H	ISOLUTE HAX 200 mg/10 mL	50

ISOLUTE-96 HAX Fixed Well Plates

Part Number	Description	Qty.
903-0025-P01	ISOLUTE-96 HAX 25 mg plate	1
903-0100-P01	ISOLUTE-96 HAX 100 mg plate	1

ISOLUTE Array HAX Loose Wells


Part Number	Description	Qty.
903-0025-R	ISOLUTE Array HAX 25 mg/1 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® HAX

TN127: Method Development in Solid Phase Extraction using ISOLUTE® HAX for the Extraction of Drugs from Biological Fluid Samples

ISOLUTE® Multimode

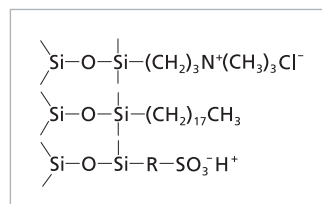

ISOLUTE® Multimode combines C18 (Octadecyl), quaternary amine and sulfonic acid functionalities

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Mixed-Mode

Application: Isolation of small neutral highly water soluble species from complex mixtures.

ISOLUTE Multimode SPE Columns


Part Number	Description	Qty.
904-0010-A	ISOLUTE Multimode 100 mg/1 mL	100
904-0030-B	ISOLUTE Multimode 300 mg/3 mL	50
904-0030-C	ISOLUTE Multimode 300 mg/6 mL	30
904-0050-B	ISOLUTE Multimode 500 mg/3 mL	50
904-0100-C	ISOLUTE Multimode 1 g/6 mL	30

Support Documents for ISOLUTE® Multimode

IST1022: Extraction of Aflatoxins from Cereals

IST1076: Extraction of Acrylamide from Cooked Foodstuffs

ISOLUTE® NH2

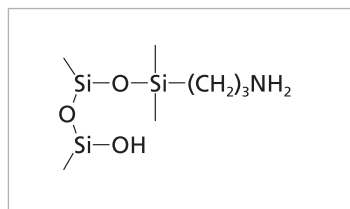

Chemical structure of NH2
aminopropyl silane covalently
bonded to the surface of a silica
particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Weak anion exchange (pKa 9.8) or polar
Exchange capacity	0.6 meq/g

Application: Extraction of strong acids and polyacidic compounds from aqueous sample matrix. Analyte elution can be performed by neutralizing charge on the sorbent. Sorbent supplied as the free base.

Alternatively, ISOLUTE NH2 can be used for extraction of polar compounds from a non-polar matrix using hydrogen bonding retention mechanism. Less retentive than SI.

ISOLUTE NH2 SPE Columns

Part Number	Description	Qty.
470-0002-A	ISOLUTE NH2 25 mg/1 mL	100
470-0005-A	ISOLUTE NH2 50 mg/1 mL	100
470-0005-G	ISOLUTE NH2 50 mg/10 mL	50
470-0010-A	ISOLUTE NH2 100 mg/1 mL	100
470-0010-B	ISOLUTE NH2 100 mg/3 mL	50
470-0010-G	ISOLUTE NH2 100 mg/10 mL	50
470-0020-B	ISOLUTE NH2 200 mg/3 mL	50
470-0050-B	ISOLUTE NH2 500 mg/3 mL	50
470-0050-C	ISOLUTE NH2 500 mg/6 mL	30
470-0050-H	ISOLUTE NH2 500 mg/10 mL	50
470-0100-C	ISOLUTE NH2 1 g/6 mL	30
470-0200-C	ISOLUTE NH2 2 g/6 mL	30
470-0200-D	ISOLUTE NH2 2 g/15 mL	20

ISOLUTE-96 NH2 Fixed Well Plates

Part Number	Description	Qty.
470-0050-P01	ISOLUTE-96 NH2 50 mg plate	1
470-0100-P01	ISOLUTE-96 NH2 100 mg plate	1

ISOLUTE Array NH2 Loose Wells

Part Number	Description	Qty.
470-0025-R	ISOLUTE Array NH2 25 mg/1 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® NH2

TN102: Method Development in Solid Phase Extraction using Polar ISOLUTE® SPE Columns for the Extraction of Non-aqueous Samples

TN104: Method Development in Solid Phase Extraction using ISOLUTE® NH2 SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® PSA

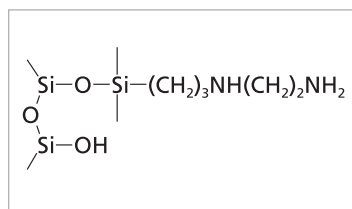

Chemical structure of PSA
ethylenediamine-n-propyl
silane covalently bonded to
the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Weak anion exchange (pKa 10.1 and 10.9) or polar
Exchange capacity	0.4 meq/g

Application: ISOLUTE PSA can be used to extract strong acids and polyacidic compound from aqueous sample matrix. Analyte elution can be achieved by neutralizing the charge on the sorbent. Sorbent will complex with certain metal ions. Supplied as free base.

Alternatively, ISOLUTE PSA can be used for extraction of polar compounds from a non-polar matrix using a hydrogen bonding retention mechanism. Less retentive than silica when used in this mode.

ISOLUTE PSA SPE Columns

Part Number	Description	Qty.
480-0010-A	ISOLUTE PSA 100 mg/1 mL	100
480-0010-B	ISOLUTE PSA 100 mg/3 mL	50
480-0020-B	ISOLUTE PSA 200 mg/3 mL	50
480-0050-B	ISOLUTE PSA 500 mg/3 mL	50
480-0050-C	ISOLUTE PSA 500 mg/6 mL	30
480-0100-C	ISOLUTE PSA 1 g/6 mL	30

ISOLUTE-96 PSA Fixed Well Plates

Part Number	Description	Qty.
480-0025-P01	ISOLUTE-96 PSA 25 mg plate	1
480-0100-P01	ISOLUTE-96 PSA 100 mg plate	1

ISOLUTE Array PSA Loose Wells

Part Number	Description	Qty.
480-0025-R	ISOLUTE Array PSA 25 mg/1 mL wells	100
480-0100-T	ISOLUTE Array PSA 100 mg/2 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® PSA

TN102: Method Development in Solid Phase Extraction using Polar ISOLUTE® SPE Columns for the Extraction of Non-aqueous Samples

TN105: Method Development in Solid Phase Extraction using ISOLUTE® PSA SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® SAX

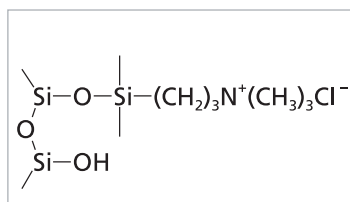

Chemical structure of SAX quaternary amine silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Strong anion exchange
Exchange capacity	0.6 meq/g

Application: Extraction of acidic analytes from aqueous sample matrix. Supplied with chloride counter ion.

ISOLUTE SAX SPE Columns

Part Number	Description	Qty.
500-0002-A	ISOLUTE SAX 25 mg/1 mL	100
500-0005-A	ISOLUTE SAX 50 mg/1 mL	100
500-0005-AG	ISOLUTE SAX 50 mg/1 mL (Tabless)*	100
500-0005-G	ISOLUTE SAX 50 mg/10 mL	50
500-0010-A	ISOLUTE SAX 100 mg/1 mL	100
500-0010-B	ISOLUTE SAX 100 mg/3 mL	50
500-0010-C	ISOLUTE SAX 100 mg/6 mL	30
500-0010-G	ISOLUTE SAX 100 mg/10 mL	50
500-0020-B	ISOLUTE SAX 200 mg/3 mL	50
500-0020-H	ISOLUTE SAX 200 mg/10 mL	50
500-0050-B	ISOLUTE SAX 500 mg/3 mL	50
500-0050-BG	ISOLUTE SAX 500 mg/3 mL (Tabless)*	50
500-0050-C	ISOLUTE SAX 500 mg/6 mL	30
500-0050-H	ISOLUTE SAX 500 mg/10 mL	50
500-0100-B	ISOLUTE SAX 1 g/3 mL	50
500-0100-C	ISOLUTE SAX 1 g/6 mL	30
500-0100-CG	ISOLUTE SAX 1 g/6 mL (Tabless)*	30
500-0200-D	ISOLUTE SAX 2 g/15 mL	20

ISOLUTE-96 SAX Fixed Well Plates

Part Number	Description	Qty.
500-0025-P01	ISOLUTE-96 SAX 25 mg plate	1
500-0050-P01	ISOLUTE-96 SAX 50 mg plate	1
500-0100-P01	ISOLUTE-96 SAX 100 mg plate	1

ISOLUTE Array SAX Loose Wells

Part Number	Description	Qty.
500-0050-R	ISOLUTE Array SAX 50 mg/1 mL wells	100
500-0100-T	ISOLUTE Array SAX 100 mg/2 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® SAX

TN103: Method Development in Solid Phase Extraction using ISOLUTE® PE-AX and SAX SPE Columns for the Extraction of Aqueous Samples

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

ISOLUTE® PE-AX

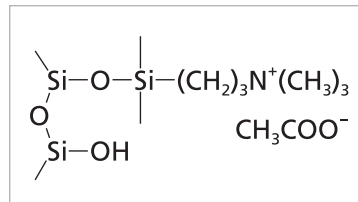

Chemical structure of PE-AX quaternary amine silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Strong anion exchange
Exchange capacity	0.6 meq/g

Application: Extraction of acidic analytes from aqueous sample matrix. Supplied with acetate counter ion for more efficient extraction of acidic analytes including those with polar/water soluble characteristics.

ISOLUTE PE-AX SPE Columns

Part Number	Description	Qty.
503-0002-A	ISOLUTE PE-AX 25 mg/1 mL	100
503-0005-A	ISOLUTE PE-AX 50 mg/1 mL	100
503-0010-A	ISOLUTE PE-AX 100 mg/1 mL	100
503-0010-B	ISOLUTE PE-AX 100 mg/3 mL	50
503-0010-C	ISOLUTE PE-AX 100 mg/6 mL	30
503-0020-B	ISOLUTE PE-AX 200 mg/3 mL	50
503-0050-B	ISOLUTE PE-AX 500 mg/3 mL	50
503-0050-C	ISOLUTE PE-AX 500 mg/6 mL	30
503-0100-C	ISOLUTE PE-AX 1 g/6 mL	30
503-0200-D	ISOLUTE PE-AX 2 g/15 mL	20

ISOLUTE-96 PE-AX Fixed Well Plates

Part Number	Description	Qty.
503-0025-P01	ISOLUTE-96 PE-AX 25 mg plate	1
503-0100-P01	ISOLUTE-96 PE-AX 100 mg plate	1

ISOLUTE Array PE-AX Loose Wells


Part Number	Description	Qty.
503-0050-R	ISOLUTE Array PE-AX 50 mg/1 mL wells	100
503-0100-T	ISOLUTE Array PE-AX 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® PE-AX

TN103: Method Development in Solid Phase Extraction using ISOLUTE® PE-AX and SAX SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® CBA

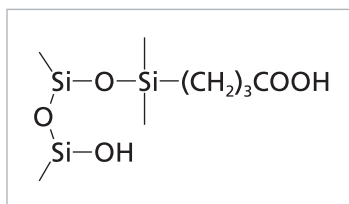

Chemical structure of CBA silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Weak anion exchange (pKa 4.8)
Exchange capacity	0.6 meq/g

Application: Extraction of strong bases and polybasic compounds from aqueous sample matrix. Analyte elution can be performed by neutralizing charge on the sorbent.

ISOLUTE CBA SPE Columns

Part Number	Description	Qty.
520-0002-A	ISOLUTE CBA 25 mg/1 mL	100
520-0005-A	ISOLUTE CBA 50 mg/1 mL	100
520-0010-A	ISOLUTE CBA 100 mg/1 mL	100
520-0010-B	ISOLUTE CBA 100 mg/3 mL	50
520-0020-B	ISOLUTE CBA 200 mg/3 mL	50
520-0050-B	ISOLUTE CBA 500 mg/3 mL	50
520-0050-C	ISOLUTE CBA 500 mg/6 mL	30
520-0100-C	ISOLUTE CBA 1 g/6 mL	30
520-0200-D	ISOLUTE CBA 2 g/15 mL	20

ISOLUTE-96 CBA Fixed Well Plates

Part Number	Description	Qty.
520-0025-P01	ISOLUTE-96 CBA 25 mg plate	1
520-0050-P01	ISOLUTE-96 CBA 50 mg plate	1
520-0100-P01	ISOLUTE-96 CBA 100 mg plate	1

ISOLUTE Array CBA Loose Wells


Part Number	Description	Qty.
520-0025-R	ISOLUTE Array CBA 25 mg/1 mL wells	100
520-0050-R	ISOLUTE Array CBA 50 mg/1 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® CBA

TN108: Method Development in Solid Phase Extraction using ISOLUTE® CBA SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® SCX

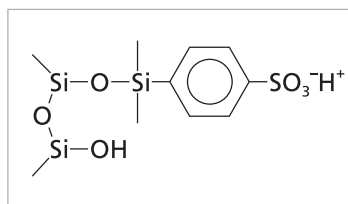

Chemical structure of benzenesulfonic acid functional group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Strong cation exchange
Exchange capacity	0.4 meq/g

Application: Extraction of basic analytes from aqueous or partially aqueous sample matrix. Supplied in the protonated form.

ISOLUTE SCX SPE Columns

Part Number	Description	Qty.
530-0005-A	ISOLUTE SCX 50 mg/1 mL	100
530-0010-A	ISOLUTE SCX 100 mg/1 mL	100
530-0010-B	ISOLUTE SCX 100 mg/3 mL	50
530-0020-B	ISOLUTE SCX 200 mg/3 mL	50
530-0050-B	ISOLUTE SCX 500 mg/3 mL	50
530-0050-C	ISOLUTE SCX 500 mg/6 mL	30
530-0050-H	ISOLUTE SCX 500 mg/10 mL	50
530-0100-B	ISOLUTE SCX 1 g/3 mL	50
530-0100-C	ISOLUTE SCX 1 g/6 mL	30
530-0200-D	ISOLUTE SCX 2 g/15 mL	20

ISOLUTE-96 SCX Fixed Well Plates

Part Number	Description	Qty.
530-0025-P01	ISOLUTE-96 SCX 25 mg plate	1
530-0050-P01	ISOLUTE-96 SCX 50 mg plate	1
530-0100-P01	ISOLUTE-96 SCX 100 mg plate	1

ISOLUTE Array SCX Loose Wells


Part Number	Description	Qty.
530-0050-R	ISOLUTE Array SCX 50 mg/1 mL wells	100
530-0100-T	ISOLUTE Array SCX 100 mg/2 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® SCX

TN106: Method Development in Solid Phase Extraction using ISOLUTE® SCX and SCX-3 SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® SCX-2

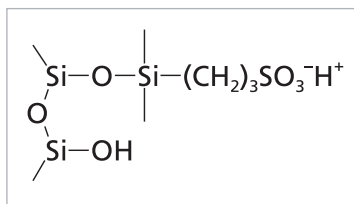

Chemical structure of propylsulfonic acid functional group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Strong cation exchanger
Exchange capacity	0.6 meq/g

Application: Extraction of basic analytes from aqueous or partially aqueous sample matrix. Supplied in the protonated form. Sorbent of choice if strong cation exchange is primary retention mechanism, or, if aqueous conditions required for analyte elution. For methods listing PRS as the SPE sorbent, SCX-2 is a direct replacement.

ISOLUTE SCX-2 SPE Columns

Part Number	Description	Qty.
532-0002-A	ISOLUTE SCX-2 25 mg/1 mL	100
532-0005-A	ISOLUTE SCX-2 50 mg/1 mL	100
532-0010-A	ISOLUTE SCX-2 100 mg/1 mL	100
532-0010-B	ISOLUTE SCX-2 100 mg/3 mL	50
532-0020-B	ISOLUTE SCX-2 200 mg/3 mL	50
532-0050-B	ISOLUTE SCX-2 500 mg/3 mL	50
532-0050-C	ISOLUTE SCX-2 500 mg/6 mL	30
532-0050-H	ISOLUTE SCX-2 500 mg/10 mL	50
532-0100-B	ISOLUTE SCX-2 1 g/3 mL	50
532-0100-C	ISOLUTE SCX-2 1 g/6 mL	30
532-0100-CG	ISOLUTE SCX-2 1 g/6 mL (Tabless)*	30
532-0200-D	ISOLUTE SCX-2 2 g/15 mL	20

ISOLUTE-96 SCX-2 Fixed Well Plates

Part Number	Description	Qty.
532-0025-P01	ISOLUTE-96 SCX-2 25 mg plate	1
532-0050-P01	ISOLUTE-96 SCX-2 50 mg plate	1
532-0100-P01	ISOLUTE-96 SCX-2 100 mg plate	1

ISOLUTE Array SCX-2 Loose Wells


Part Number	Description	Qty.
532-0025-R	ISOLUTE Array SCX-2 25 mg/1 mL wells	100

Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® SCX-2

TN107: Method Development in Solid Phase Extraction using ISOLUTE® SCX-2 SPE Columns for the Extraction of Aqueous Samples

ISOLUTE® SCX-3

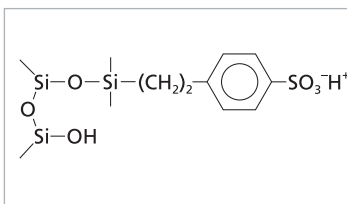

Chemical structure of ethylbenzene sulfonic acid functional group covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Strong cation exchanger
Exchange capacity	0.6 meq/g

Application: Extraction of basic analytes from aqueous or partially aqueous sample matrix. Supplied in the protonated form. Sorbent of choice where significant non-polar secondary interactions are required.

ISOLUTE SCX-3 Columns

Part Number	Description	Qty.
533-0002-A	ISOLUTE SCX-3 25 mg/1 mL	100
533-0005-A	ISOLUTE SCX-3 50 mg/1 mL	100
533-0010-A	ISOLUTE SCX-3 100 mg/1 mL	100
533-0010-B	ISOLUTE SCX-3 100 mg/3 mL	50
533-0020-B	ISOLUTE SCX-3 200 mg/3 mL	50
533-0050-B	ISOLUTE SCX-3 500 mg/3 mL	50
533-0050-C	ISOLUTE SCX-3 500 mg/6 mL	30
533-0050-H	ISOLUTE SCX-3 500 mg/10 mL	50
533-0100-B	ISOLUTE SCX-3 1 g/3 mL	50
533-0100-C	ISOLUTE SCX-3 1 g/6 mL	30
533-0200-D	ISOLUTE SCX-3 2 g/15 mL	20

ISOLUTE-96 SCX-3 Fixed Well Plates

Part Number	Description	Qty.
533-0025-P01	ISOLUTE-96 SCX-3 25 mg plate	1
533-0050-P01	ISOLUTE-96 SCX-3 50 mg plate	1
533-0100-P01	ISOLUTE-96 SCX-3 100 mg plate	1

ISOLUTE Array SCX-3 Loose Wells

Part Number	Description	Qty.
533-0025-R	ISOLUTE Array SCX-3 25 mg/1 mL wells	100
533-0050-R	ISOLUTE Array SCX-3 50 mg/1 mL wells	100
533-0100-T	ISOLUTE Array SCX-3 100 mg/2 mL wells	100


Pre-assembled ISOLUTE® Array plates are available. To order, add the suffix P to the equivalent loose well part number.

Support Documents for ISOLUTE® SCX-3

TN106: Method Development in Solid Phase Extraction using ISOLUTE® SCX and SCX-3 SPE Columns for the Extraction of Aqueous Samples

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

ISOLUTE® SI

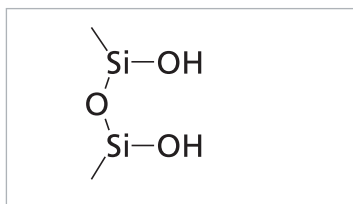

Chemical structure of silanol groups on the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Nominal moisture content	7%
Sorbent Type	Polar

Application: Extraction of polar compounds from a non-polar sample matrix using hydrogen bonding retention mechanism.


ISOLUTE SI SPE Columns

Part Number	Description	Qty.
460-0002-A	ISOLUTE SI 25 mg/1 mL	100
460-0010-A	ISOLUTE SI 100 mg/1 mL	100
460-0010-B	ISOLUTE SI 100 mg/3 mL	50
460-0010-G	ISOLUTE SI 100 mg/10 mL	50
460-0020-B	ISOLUTE SI 200 mg/3 mL	50
460-0020-H	ISOLUTE SI 200 mg/10 mL	50
460-0050-B	ISOLUTE SI 500 mg/3 mL	50
460-0050-C	ISOLUTE SI 500 mg/6 mL	30
460-0050-H	ISOLUTE SI 500 mg/10 mL	50
460-0100-B	ISOLUTE SI 1 g/3 mL	50
460-0100-C	ISOLUTE SI 1 g/6 mL	30
460-0200-C	ISOLUTE SI 2 g/6 mL	30
460-0200-L	ISOLUTE SI 2 g/6 mL (Glass)	30
460-0200-D	ISOLUTE SI 2 g/15 mL	20
460-0500-E	ISOLUTE SI 5 g/25 mL	20
460-1000-F	ISOLUTE SI 10 g/70 mL	16

Support Documents for ISOLUTE® SI

TN102: Method Development in Solid Phase Extraction using Polar ISOLUTE® SPE Columns for the Extraction of Non-aqueous Samples

ISOLUTE® DIOL

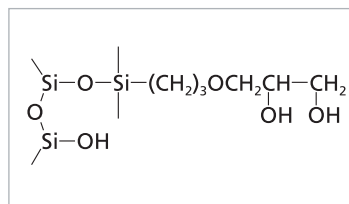

Chemical structure of DIOL silane covalently bonded to the surface of a silica particle

Average particle size	50 µm
Pore diameter	60 Å
Sorbent Type	Polar

Application: Extraction of polar compounds from a non-polar sample matrix using hydrogen bonding retention mechanism. Less retentive sorbent than SI. Can also be used as a non-polar sorbent with aqueous sample matrices.

ISOLUTE DIOL SPE Columns

Part Number	Description	Qty.
430-0010-A	ISOLUTE DIOL 100 mg/1 mL	100
430-0010-B	ISOLUTE DIOL 100 mg/3 mL	50
430-0010-G	ISOLUTE DIOL 100 mg/10 mL	50
430-0020-B	ISOLUTE DIOL 200 mg/3 mL	50
430-0050-B	ISOLUTE DIOL 500 mg/3 mL	50
430-0050-C	ISOLUTE DIOL 500 mg/6 mL	30
430-0050-H	ISOLUTE DIOL 500 mg/10 mL	50
430-0100-C	ISOLUTE DIOL 1 g/6 mL	30

Support Documents for ISOLUTE® DIOL

TN102: Method Development in Solid Phase Extraction using Polar ISOLUTE® SPE Columns for the Extraction of Non-aqueous Samples

ISOLUTE® FL

Average particle size	150–250 µm
Pore diameter	60 Å
Nominal moisture content	<2%
Surface pH	8
Sorbent Type	Polar

Application: Extraction of polar compounds from a non-polar sample matrix. Alternative to silica based polar sorbents. Minimal retention of basic compounds. Activated for separation of chlorinated pesticides.

ISOLUTE FL SPE Columns

Part Number	Description	Qty.
712-0010-A	ISOLUTE FL 100 mg/1 mL	100
712-0020-A	ISOLUTE FL 200 mg/1 mL	100
712-0020-B	ISOLUTE FL 200 mg/3 mL	50
712-0020-H	ISOLUTE FL 200 mg/10 mL	50
712-0050-B	ISOLUTE FL 500 mg/3 mL	50
712-0050-BG	ISOLUTE FL 500 mg/3 mL (Tabless)*	50
712-0050-C	ISOLUTE FL 500 mg/6 mL	30
712-0050-H	ISOLUTE FL 500 mg/10 mL	50
712-0050-L	ISOLUTE FL 500 mg/6 mL (Glass)	30
712-0100-B	ISOLUTE FL 1 g/3 mL	50
712-0100-C	ISOLUTE FL 1 g/6 mL	30
712-0200-D	ISOLUTE FL 2 g/15 mL	20
712-0500-E	ISOLUTE FL 5 g/25 mL	20
712-2000-F	ISOLUTE FL 20 g/70 mL	16

Support Documents for ISOLUTE FL

IST1080 Multi-residue Extraction and Clean-up of Pesticides from Fruits and Vegetables

ISOLUTE® HM-N

Columns for Supported Liquid Extraction

Part Number	Description	Qty.
800-0040-BM	ISOLUTE HM-N (300 µL sample)	100
800-0100-CM	ISOLUTE HM-N (1 mL sample)	100
800-0220-DM	ISOLUTE HM-N (3 mL sample)	100
800-0350-EM	ISOLUTE HM-N (5 mL sample)	100
800-0700-FM	ISOLUTE HM-N (10 mL sample)	50
800-1300-FM	ISOLUTE HM-N (20 mL sample)	50

*Tabless columns for use with Pressure+ Positive Pressure Manifolds and other automated SPE systems. Other tabless columns are available, contact Biotage for details.

ISOLUTE® Alumina

High activity, 50–200 µm particle size range alumina, available in acidic, neutral and basic surface pH options. The surface of the alumina can absorb molecules by interaction with the aluminum metal center, hydrogen bonding with surface hydroxyl groups, or by ion exchange if the surface carries a charge. The extent of these different interactions can be enhanced by control of the surface pH by treatment with acidic, basic or neutral solutions.

ISOLUTE AL-A

Acid washing results in a surface with decreased capacity for basic compounds. Compounds are retained by ion exchange with the positively charged surface or by specific interactions with the metal center. ISOLUTE® AL-A has a nominal moisture content of <0.1% (Brockman Activity I/Super 1 grade), and surface pH of 4.5.

ISOLUTE AL-N

Neutral surface allows interactions of the aluminum metal center with compounds whose heteroatoms are electronegative (e.g. N, O, P & S) or whose highly aromatic structure makes them “electron rich”. The adsorbent can be useful for retaining amines and aromatic compounds from either aqueous or non-aqueous solvents. ISOLUTE AL-N has a nominal moisture content of <0.1% (Brockman Activity I/Super 1 grade), and surface pH of 7.5.

ISOLUTE AL-B

Washing this material with a basic solution results in a net negative charge. Cationic compounds are retained on the negatively charged surface or by specific interactions with the metal center. ISOLUTE AL-B has a nominal moisture content of <0.1% (Brockman Activity I/Super 1 grade), and surface pH of 10.

Column configuration	Qty.	Part Numbers		
		AL-A	AL-N	AL-B
100 mg/1 mL	100	713-0010-A	714-0010-A	715-0010-A
200 mg/1 mL	100	713-0020-A	714-0020-A	715-0020-A
200 mg/3 mL	50	713-0020-B	714-0020-B	715-0020-B
500 mg/3 mL	50	713-0050-B	714-0050-B	715-0050-B
500 mg/6 mL	30	713-0050-C	714-0050-C	715-0050-C
1 g/3 mL	50	713-0100-B	714-0100-B	715-0100-B
1 g/6 mL	30	713-0100-C	714-0100-C	715-0100-C
2 g/6 mL	30	—	714-0200-C	—
2 g/15 mL	20	713-0200-D	714-0200-D	715-0200-D
5 g/25 mL	20	713-0500-E	714-0500-E	715-0500-E
10 g/70 mL	16	—	714-1000-F	715-1000-F

ISOLUTE® EPH SPE Columns

Fractionate extractable petroleum hydrocarbons in soil extracts prior to GC analysis. ISOLUTE EPH columns are optimized to ensure no breakthrough of lower MW aromatic hydrocarbons (PAHs) in the aliphatic fraction. Application notes with methodology optimized for automated SPE processing are available. Request application note AN703 (ASPEC) or AN704 (RapidTrace) for automated methods.



Part Number	Description	Qty.
928-0145-B	ISOLUTE EPH 1.45 g/3 mL Columns for automated method	50
928-0500-E	ISOLUTE EPH 5 g/25 mL Columns for manual method	20

ISOLUTE® Layered SPE Columns

Simultaneously extract analytes with a broad range of polarity characteristics from aqueous samples.

Part Number	Description	Qty.
933-0050-B	ISOLUTE C2/C18(EC) 500 mg/3 mL	50
933-0100-C	ISOLUTE C2/C18(EC) 1 g/6 mL	30
934-0040-C	ISOLUTE C8/ENV+ 400 mg/ 6 mL	30
935-0040-C	ISOLUTE C18/ENV+ 400 mg/ 6 mL	30
941-0020-B	ISOLUTE C2/ENV+ 200 mg/3 mL	50

ISOLUTE® PAH SPE Columns

Extract PAHs from water samples containing polar interferences such as humic acids. The layered column removes humic acids from the final extract, so that they cannot interfere in the final analysis. Download application note IST1025A for more information.

Part Number	Description	Qty.
927-0075-B	ISOLUTE PAH 750 mg/3 mL	50
927-0150-C	ISOLUTE PAH 1.5 g/6 mL	30
927-0150-CD	ISOLUTE PAH 1.5 g/6 mL (Depth filter)	30

ISOLUTE® SAX/PSA SPE Columns

Remove polar interferences and pigments from plant extracts prior to pesticide residue analysis. Download application note IST1027A for more information.

Part Number	Description	Qty.
924-0050-C	ISOLUTE SAX/PSA 500 mg/6 mL	30
924-0100-C	ISOLUTE SAX/PSA 1 g/6 mL	30

ISOLUTE® TPH SPE Columns

Extract total petroleum hydrocarbons from water samples and eliminate tedious shaking and emulsion problems common in liquid-liquid extraction methods. Manual and automated methodologies are available, download application notes IST1042 and IST1018 respectively for more information.

Part Number	Description	Qty.
752-0100-C	ISOLUTE TPH 1 g/6 mL	30

ISOLUTE® O&G SPE columns

Extract total Oil and Grease from water samples and fractionate into HEM and SGT-HEM fractions in one step. Application note IST1005 meets US EPA 1664 guidelines.



Part Number	Description	Qty.
753-0100-CD	ISOLUTE O&G 1 g/6 mL	30
753-0300-FD	ISOLUTE O&G 3 g/70 mL	16

ISOLUTE® Na₂SO₄/FL SPE columns

Remove halogenated interferences from mineral oil extracts. Application note IST1077 complies with the ISO9377-2:200 method.

Part Number	Description	Qty.
976-0400-C	ISOLUTE Na ₂ SO ₄ /FL 4 g/6 mL	30



ISOLUTE[®] MYCO SPE Columns

Solid Phase Extraction Columns for Mycotoxin Analysis

ISOLUTE® Myco SPE columns offer simple and efficient multiple mycotoxin sample preparation from a wide range of matrices, ideally suited for selective and fast LC-MS/MS analysis. ISOLUTE Myco SPE columns contain a novel polymer-based sorbent designed specifically to be selective enough to isolate a wide variety of different mycotoxins.

ISOLUTE Myco SPE columns are available in tabless (flangeless) format for use on both manual and automated sample processing manifolds and systems.

- » Matrix specific application notes: Each method confidently meets the maximum residue limits (MRL) set down by EU and US regulations in terms of detection limits, %RSD and recoveries, in most cases far exceeding these requirements.
- » ISOLUTE Myco clean up methods are quick, simple to use and require no offline steps; making the technique suitable for automation, especially useful for high throughput testing laboratories.
- » Unlike immunoaffinity columns, ISOLUTE Myco columns are cost effective, easy to use and have no special storage requirements.




ISOLUTE Myco SPE Columns

Part Number	Description	Qty.
150-0006-BG	ISOLUTE Myco 60 mg/3 mL (Tabless)	50

Matrix Specific Application Notes

Matrix specific application notes now available

- » AN781 Extraction of Patulin from Apple Juice Using ISOLUTE® Myco prior to LC-MS/MS Analysis
- » AN782 Extraction of Multiple Mycotoxins From Grain Using ISOLUTE® Myco prior to LC-MS/MS Analysis
- » AN783 Extraction of Deoxynivalenol From Grain Using ISOLUTE® Myco prior to LC-MS/MS Analysis
- » AN784 Extraction of Multiple Mycotoxins From Nuts Using ISOLUTE® Myco prior to LC-MS/MS Analysis

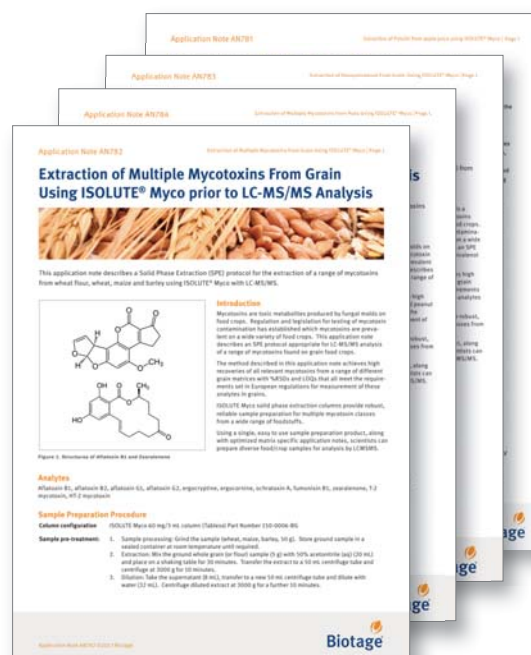



ISOLUTE[®] Bulk Sorbents

The same high quality sorbents used to manufacture ISOLUTE® SPE columns and 96-well plates are also available in bulk. The tightly controlled physical and chemical properties of these sorbents which provide reproducible performance in SPE procedures also ensure their suitability when used in applications that require loose material.

Non-Polar Sorbents

Part Number	Description	Qty.
9220-0100	ISOLUTE C18, bulk	100 g
9221-0100	ISOLUTE C18(EC)*, bulk	100 g
9290-0100	ISOLUTE C8, bulk	100 g
9291-0100	ISOLUTE C8(EC)*, bulk	100 g
9320-0100	ISOLUTE C2, bulk	100 g
9321-0100	ISOLUTE C2(EC)*, bulk	100 g

*(EC) - endcapped, a chemical process to reduce the concentration of silica surface silanol groups that provide polar and weak cation exchange secondary interactions. For more details see the QuickStart Guide to SPE.

Polar Sorbents

Part Number	Description	Qty.
9460-0100	ISOLUTE SI, bulk	100 g
9470-0100	ISOLUTE NH2, bulk	100 g
9480-0100	ISOLUTE PSA, bulk	100 g
9712-0100	ISOLUTE FL, bulk	100 g

Ion Exchange Sorbents

Part Number	Description	Qty.
9470-0100	ISOLUTE NH2, bulk	100 g
9480-0100	ISOLUTE PSA, bulk	100 g
9500-0100	ISOLUTE SAX, bulk	100 g
9503-0100	ISOLUTE PE-AX, bulk	100 g
9520-0100	ISOLUTE CBA, bulk	100 g
9532-0100	ISOLUTE SCX-2, bulk	100 g
9533-0100	ISOLUTE SCX-3, bulk	100 g

ISOLUTE® sorbents are available in additional quantities to those detailed in the above table.

ISOLUTE HM-N, a modified form of diatomaceous earth, is also available in bulk. The material is used for applications requiring an inert support as part of the sample preparation process (e.g. SFE). For details of pre-packed columns containing ISOLUTE HM-N, see [page 34](#).

Inert Support (modified diatomaceous earth)

Part Number	Description	Qty.
9800-1000	ISOLUTE HM-N, bulk	1 kg
9800-5000	ISOLUTE HM-N, bulk	5 kg

ISOLUTE® HM-N is available in additional quantities to those detailed in the above table. Please contact your local Biotage representative.





Protein Precipitation Products

Fast, Simple Protein Removal

Protein Precipitation Plates

The optimized filtration system in ISOLUTE® PPT+ plates provides an easy to automate solution for efficient protein removal from biological fluid samples.



96-Well Plates

Part Number	Description	Qty.
120-2040-P01	ISOLUTE PPT+ fixed well plate, 2 mL	1

ISOLUTE Array Wells

Part Number	Description	Qty.
120-2040-R	ISOLUTE Array PPT+ wells, 1 mL	100
120-2040-T	ISOLUTE Array PPT+ wells, 2 mL	100

ISOLUTE Array Accessories

Part Number	Description	Qty.
120-1000-P01	ISOLUTE Array base plate	1
120-1200	Strip of 8 base plate sealing plugs*	50
120-1201	Luer adaptors (to fit any standard sample processing manifold)	25
120-1202	Well removing tool	1

*Required when processing a partially populated ISOLUTE Array PPT+ plate.

Support Documents for ISOLUTE® PPT+

TN130: Sample Preparation using ISOLUTE PPT+ Protein Precipitation Plates.

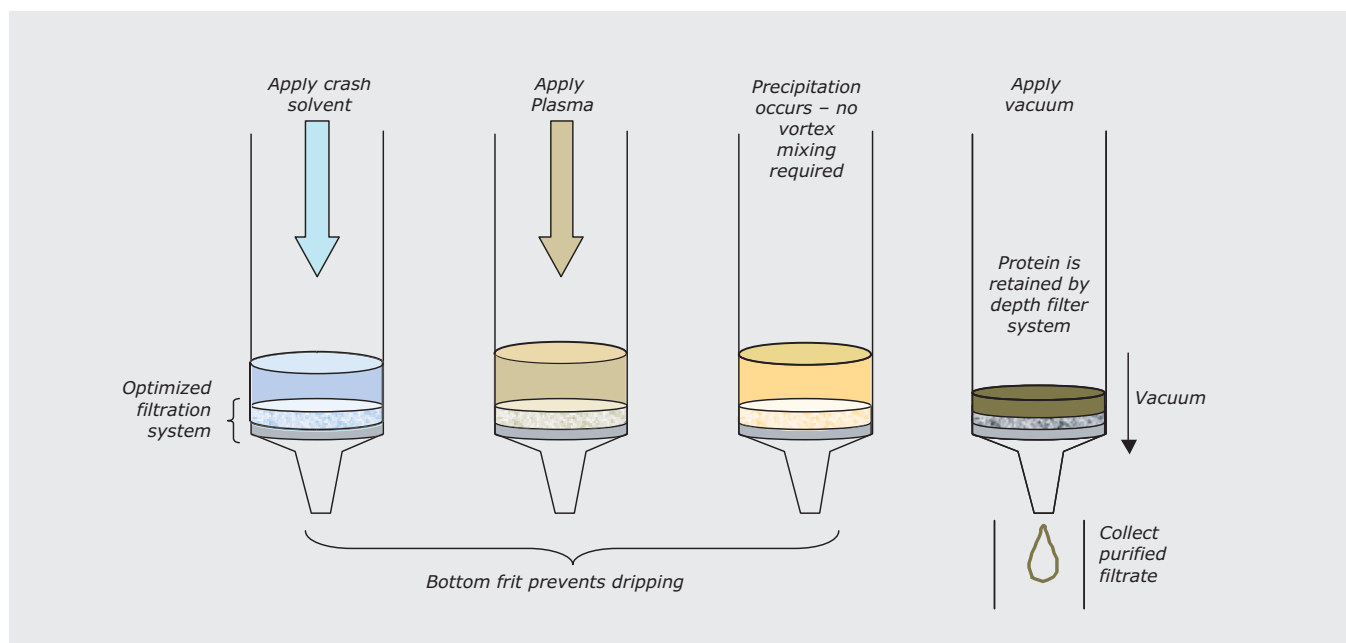


Figure 1. Protein precipitation procedure using ISOLUTE® PPT+ plates.



ISOLUTE® QuEChERS

Fast and Efficient Clean up of
Complex Food Samples

ISOLUTE® QuEChERS

ISOLUTE® QuEChERS

ISOLUTE® QuEChERS (**Q**uick, **E**asy, **C**heap, **E**ffective, **R**ugged and **S**afe) products provide simple clean up of complex samples using salt assisted extraction and partitioning, followed by dispersive SPE. The pre-weighed extraction and clean up tubes conform to AOAC and EN methodologies, and include options for waxed and highly pigmented fruit and vegetable samples.



Extraction Tubes

Part Number	Description	Tube	Pack Size	MgSO ₄	Na Acetate	Na Citrate	Na Citrate sesqui-hydrate	NaCl
Q0010-15V	15 g QuEChERS AOAC 15 mL Extraction Tube	15 mL	25	6 g	1.5 g			
Q0020-15V	10 g QuEChERS EN 15 mL Extraction Tube	15 mL	25	4 g		1 g	0.5 g	1 g

Clean up Tubes

Part Number	Description	Tube	Pack Size	PSA	MgSO ₄ (purest)	C18(EC)	GCB
Q0030-15V	AOAC Fruit and Vegetables Clean up Tube	15 mL	25	400 mg	1200 mg		
Q0035-15V	EN Fruit and Vegetable Clean up Tube	15 mL	25	150 mg	900 mg		
Q0050-15V	AOAC Waxed Fruit and Vegetables Clean up Tube	15 mL	25	400 mg	1200 mg	400 mg	
Q0060-15V	EN Waxed Fruit and Vegetables Clean up Tube	15 mL	25	150 mg	900 mg	150 mg	
Q0070-15V	AOAC Pigmented Fruit and Vegetables Clean up Tube	15 mL	25	400 mg	1200 mg		400 mg
Q0080-15V	EN Pigmented Fruit and Vegetables Clean up Tube	15 mL	25	150 mg	900 mg		15 mg
Q0090-15V	EN Highly Pigmented Fruit and Vegetables Clean up Tube	15 mL	25	150 mg	900 mg		45 mg

Centrifuge Tubes

Part Number	Description	Tube	Pack Size
Q0000-50V	50 mL Centrifuge Tube with Rack	50 mL	25



Sample Preparation Accessories


Accessories

For SPE and Filtration

ISOLUTE® Column Adaptors

Column adaptors attach empty reservoirs, filtration columns and depth filters to SPE columns.

Part Number	Description		Qty.
120-1100	PTFE Adaptor - 1, 3 & 6 mL and XL (columns A, B, C, G & H)	1	10
120-1101	PE Adaptor - 1, 3 & 6 mL (columns A, B & C)	2	10
120-1102	PE Adaptor - 15 & 25 mL (columns D & E)	3	10
120-1103	PE Adaptor - 70 mL (column F)	4	10
120-1104	PTFE Adaptor - 6 mL Glass columns	5	10

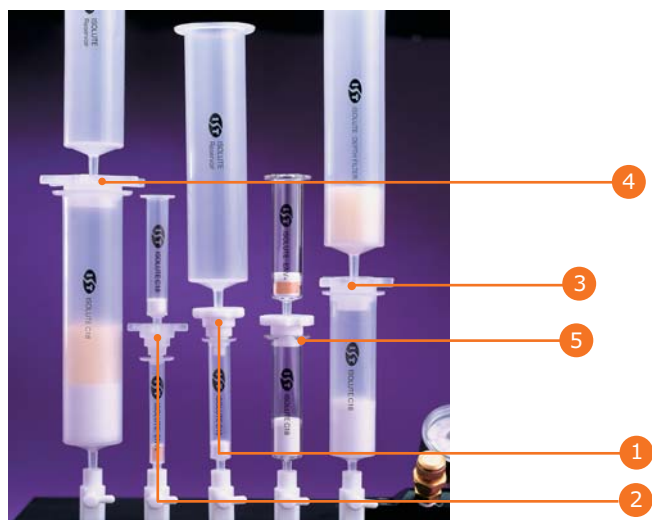

ISOLUTE® Depth Filter Reservoirs

Reservoirs with pre-fitted depth filters remove particulate matter from aqueous samples, and prevent column blockage. Also allow for separate analysis of particulate compounds.



Part Number	Description		Qty.
120-1003-CD	ISOLUTE depth filter, 6 mL	1	100
120-1009-FD	ISOLUTE depth filter, 70 mL	2	50

ISOLUTE® Column Caps

Seal prepared columns, samples loaded columns or immunoaffinity columns for transportation or storage.

Part Number	Description	Column Type	Qty.
1201-0120	Bottom Luer cap, fits all columns	All	100
1201-0121-A	Top cap, 1 mL column	A	100
1201-0122-B	Top cap, 3 mL column	B	100
1201-0123-C	Top cap, 6 mL column	C	100
1201-0125-H	Top cap, 10 mL - XL column	G & H	100
1201-0126-D	Top cap, 15 mL column	D	100
1201-0127-E	Top cap, 25 mL column	E	100
1201-0128-F	Top cap, 70 mL column	F	100



ISOLUTE® Empty Reservoirs

Stack above SPE columns to increase the reservoir volume.

Part Number	Description	Qty.
120-1001-A	ISOLUTE Reservoir 1 mL	100
120-1002-B	ISOLUTE Reservoir 3 mL	100
120-1003-C	ISOLUTE Reservoir 6 mL	100
120-1006-D	ISOLUTE Reservoir 15 mL	100
120-1007-E	ISOLUTE Reservoir 25 mL	100
120-1009-F	ISOLUTE Reservoir 70 mL	50



ISOLUTE® Frits

20 µm porosity sintered polyethylene frits to fit all SPE column sizes.

Part Number	Description	Column Type	Qty.
120-1031-A	ISOLUTE Frit, 6 mm (1/4") diameter – 20 µm PE	A (1 mL)	100
120-1032-G	ISOLUTE Frit, 6 mm (1/4") diameter – 20 µm PE	G (XL)	100
120-1033-B	ISOLUTE Frit, 9 mm (3/8") diameter – 20 µm PE	B (3 mL)	100
120-1034-H	ISOLUTE Frit, 9 mm (3/8") diameter – 20 µm PE	H (XL)	100
120-1035-C	ISOLUTE Frit, 13 mm (1/2") diameter – 20 µm PE	C (6 mL)	100
120-1036-D	ISOLUTE Frit, 16 mm (5/8") diameter – 20 µm PE	D (15 mL)	100
120-1037-E	ISOLUTE Frit, 20 mm (13/16") diameter – 20 µm PE	E (25 mL)	100
120-1038-F	ISOLUTE Frit, 27 mm (11/16") diameter – 20 µm PE	F (70 mL)	100



ISOLUTE® Filtration Columns

Pre-fitted with two 20 µm polyethylene frits. Use as standalone filters or stacked above an SPE column.

Part Number	Description	Qty.
120-1021-A	ISOLUTE double fritted filtration column, 1 mL–20 µm PE	100
120-1022-B	ISOLUTE double fritted filtration column, 3 mL–20 µm PE	100
120-1023-C	ISOLUTE double fritted filtration column, 6 mL–20 µm PE	100
120-1024-H	ISOLUTE double fritted filtration column, 10 mL–20 µm PE	50
120-1025-D	ISOLUTE double fritted filtration column, 15 mL–20 µm PE	100
120-1026-E	ISOLUTE double fritted filtration column, 25 mL–20 µm PE	100
120-1028-F	ISOLUTE double fritted filtration column, 70 mL–20 µm PE	50

ISOLUTE® Sodium Sulfate Drying Cartridges

Dry organic solvents with this easy to use format (contains 2.5 g of high purity sodium sulfate per cartridge). Stack beneath SPE columns during elution for efficient in-line solvent drying.

Part Number	Description	Qty.
802-0250-M	ISOLUTE Sodium Sulfate Drying Cartridge	50



Accessories

For SPE and Filtration

ISOLUTE® Filtration Plates

Each well contains a 20 µm polyethylene frit for high throughput filtration applications.

Part Number	Description	Qty.
120-1022-P01	Filtration plate, fixed well, 2 mL–20 µm	1



Plate Sealing Accessories

Seal wells and Luer outlets of ISOLUTE-96 fixed well plates. Piercable sealing cap can also be used to seal collection plates.

Part Number	Description	Qty.
121-5204	Piercable Sealing Cap	50
121-5205	Luer Cap Mat	25



Deep Well Collection Plates

1 or 2 mL volume per well, constructed from high purity, solvent resistant polypropylene. Seal with piercable sealing cap.

Part Number	Description	Qty.
121-5202	Collection plate, 1 mL	50
121-5203	Collection plate, 2 mL	50



Array Accessories

Populate the reusable base plate to transfer individual Array wells into microplate format. Seal any unused positions with base plate plugs. Align wells and depopulate base plate after use with the well removing tool. Use Luer adaptors to process individual Array wells using a VacMaster-10 or -20 vacuum manifold.

Part Number	Description	Qty.
120-1000-P01	ISOLUTE Array base plate	1
120-1200	ISOLUTE base plate sealing strips (strips of 8)	50
120-1201	Luer Adaptors (to fit any vacuum manifold)	25
120-1202	Well removing tool	1

Manual Sample Processing Products and Accessories

Biotage® PRESSURE+ 48 and PRESSURE+ 96

Positive Pressure Manifolds

Biotage® PRESSURE+ 96

Part Number	Description
PPM-96	PRESSURE+ 96 Positive Pressure Manifold 96 position. Includes 1 x 1 mL 96 well collection plate, 1 x 2 mL 96 well collection plate, 1 x 10 mL 24 well collection plate (waste plate), 96-column sealing gasket, gas supply adaptor kit (6' of 1/8" i.d. polyethylene tubing and 1/8" and 1/4" NPT connectors), User Manual CD-ROM.

Accessories

PPM-A96-GSKT	PRESSURE+ 96 Sealing Gasket 96 position
PPM-A96-1024	PRESSURE+ 96 Collection Tray 10 mL 24 well
PPM-GA	PRESSURE+ Gas Supply Adaptor for all models
121-5202	96-well collection plate, 1 mL, pack of 50
121-5203	96-well collection plate, 2 mL, pack of 50



Biotage® PRESSURE+ 48

Part Number	Description
PPM-48	PRESSURE+ 48 Positive Pressure Manifold 48 position. Includes Sealing gasket 48 position, gas supply adaptor kit (6' of 1/8" i.d. polyethylene tubing and 1/8" and 1/4" NPT connectors), waste bin rack, waste bin inserts, User Manual CD-ROM. Racks must be ordered separately.

Racks

PPM-A48-1RCK	PRESSURE+ 48 SPE Column Rack 1 mL
PPM-A48-3RCK	PRESSURE+ 48 SPE Column Rack 3 mL
PPM-A48-6RCK	PRESSURE+ 48 SPE Column Rack 6 mL
PPM-A48-1275	PRESSURE+ 48 Collection Rack 12 X 75 mm
PPM-A48-13100	PRESSURE+ 48 Collection Rack 13 x 100 mm
PPM-A48-1232	PRESSURE+ 48 Sample Vial Rack 12 x 32 mm
PPM-A48-16100	PRESSURE+ 48 Collection Rack 16 x 100 mm

Accessories

PPM-A48-GSKT	PRESSURE+ 48 Sealing Gasket 48 position
PPM-A48-WST	PRESSURE+ 48 Waste Bin Inserts
PPM-GA	PRESSURE+ Gas Supply Adaptor for all models

Service Options

Part Number	Description
SER-PPM-IN	PRESSURE+ Manifold Service Installation
SER-PPM-PL	PRESSURE+ Manifold Platinum Service Contract - 1 year
SER-PPM-SLVR	PRESSURE+ Manifold Silver Repair Depot Service Contract - 1 year
SER-PPM-FYMP	PRESSURE+ Manifold 9 month extended warranty (1 year total)

Biotage® VacMaster™ -10 and -20

10 or 20 position vacuum manifolds for SPE, SLE and filtration applications.



VacMaster-10 Sample Processing Manifold

Part Number	Description
121-1010	VacMaster-10 Sample Processing Manifold (with 10 mm rack)
121-1012	VacMaster-10 Sample Processing Manifold (with 12 mm rack)
121-1016	VacMaster-10 Sample Processing Manifold (with 16 mm rack)

VacMaster-10 Replacement Parts

Part Number	Description	Qty.
121-1039	Replacement tank with fittings	1
121-1045	Replacement lid with fittings	1
121-1030	Silicone Lid Gasket	1
121-1132	Replacement rack, 10 mm	1
121-1133	Replacement rack, 12 mm	1
121-1134	Replacement rack, 16 mm	1

VacMaster-20 Sample Processing Manifold

Part Number	Description
121-2010	Sample Processing Manifold (with 10 mm rack)
121-2012	Sample Processing Manifold (with 12 mm rack)
121-2016	Sample Processing Manifold (with 16 mm rack)

VacMaster-20 Replacement Parts

Part Number	Description	Qty.
121-2068	Replacement tank with fittings	1
121-2075	Replacement lid with fittings	1
121-2059	Silicone Lid Gasket	1
121-2161	Replacement rack, 10 mm	1
121-2162	Replacement rack, 12 mm	1
121-2163	Replacement rack, 16 mm	1

VacMaster Accessories

Part Number	Description	Qty.
121-0010	PTFE 'T' valve (for connecting 2 tanks)	1

PTFE Stopcock Option and Spare Parts

Part Number	Description	Qty.
121-0009	Universal PTFE stopcock 1	10
121-0001	PTFE stopcock/needle unit 2	10
121-0002	PTFE needle 3	10
121-0003	Stainless steel needle 4	20
121-0004	Stainless steel needle retainer	10
121-0005	Port sealing plugs	30
121-0009-S	PTFE stopcock - pressure positive	10



Biotage® VacMaster™-96

Sample Processing Manifold

Biotage® VacMaster™-96

Vacuum manifold for processing 96-well SPE, SLE or filtration plates.

Part Number	Description	Qty.
121-9600	VacMaster-96 Sample Processing Manifold (without vacuum control)	1

Vacuum Control Options

121-9601	VacMaster VCU-1 Vacuum Control Unit	1
121-9602	VacMaster VCU-2 Vacuum Control and Generation Unit	1


96 Well Collection Plates - Deep Well

121-5202	Collection plate, 1 mL	50
121-5203	Polypropylene, Square well - 2 mL	50



Replacement Parts

Part Number	Description		Qty.
121-9612	VacMaster-96 replacement gasket	1	1
121-9613	VacMaster-96 replacement o-ring	2	1
121-9614	Collection plate spacer (2 mm) for deep well collection plates	3	1
121-9610	Array insert (6 mm), acetal	4	1
121-9615	Collection plate spacer (29 mm) for shallow well collection plates	5	1
121-9611	VacMaster-96 insert (12 mm) for "shallow skirt" plates	6	1

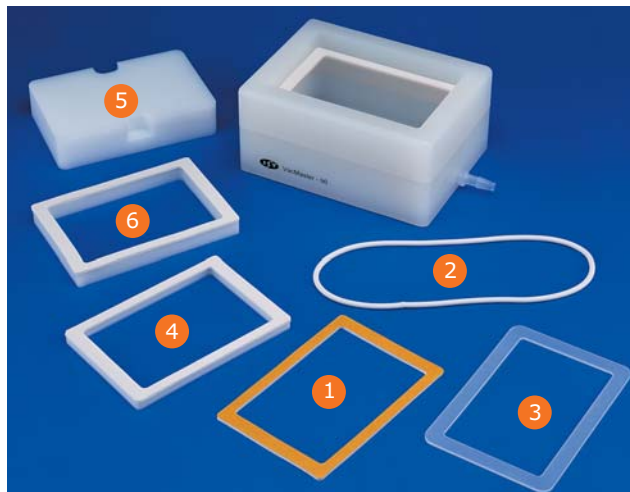

Sample Processing Accessories

Biotage® VacMaster™ Drying Adaptor

Connect to laboratory air or nitrogen supply to dry 10 or 20 SPE columns simultaneously.

Part Number	Description	Qty.
124-1001	VacMaster-10 Drying Adaptor for 1,3 & 6 mL columns	1
124-2001	VacMaster-20 Drying Adaptor for 1,3 & 6 mL columns	1



Biotage® VacMaster™ Trap Kit

Connect between VacMaster (-10, -20, or -96) and vacuum source to collect waste. 1 L volume.

Part Number	Description	Qty.
121-2095	VacMaster trap kit (complete)	1



Biotage® VacMaster™ Large Volume Extraction (LVE) Kit

For unattended loading of large volume samples. Inert PTFE tubing prevents sample contamination.

Part Number	Description	Qty.
121-2090	VacMaster LVE kit (complete) for 1, 3 and 6 mL columns	1



Gravity Rack

Gravity Rack

Process up to 20 samples simultaneously under gravity.


Part Number	Description	Qty.
123-2016	Gravity Rack with 16 mm collection tube rack	1
123-2019	Gravity Rack with 19 mm collection tube rack	1
Accessories		
121-0009	Universal PTFE stopcock	10
121-0001	PTFE stopcock/needle unit	10
121-0002	PTFE needle unit	10
121-0003	Stainless steel needle	20
121-0004	Stainless steel needle retainer	10



Automated Sample Processing Products

RapidTrace®+

Automated High Throughput Solid Phase Extraction

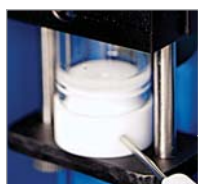

Closed fluid path system with segregated waste lines maintain integrity of sample.



Eight solvent/reagent lines for maximum methods flexibility.



Magnetically encoded racks of routine samples.



Reagent mixing chamber allows automated online reagent blending in the method.



Positive pressure, stepper motor syringe pump delivers precise flow rates.




Compatible 1, 3 and 6 mL version SPE columns.

Part Number	Description
C50000	RapidTrace+ Workstation, 1 mL and 3 mL (10 columns)
C125713	RapidTrace+ Workstation, 6 mL (5 columns). 3 mL and 1 mL operation possible by purchasing appropriate plunger (see accessories brochure)
C52006	RapidTrace Start-Up Kit with Software
C52689	RapidTrace Notebook Controller

RapidTrace®+ Brochure

For further information please visit www.biotage.com to request a copy of our RapidTrace brochure.

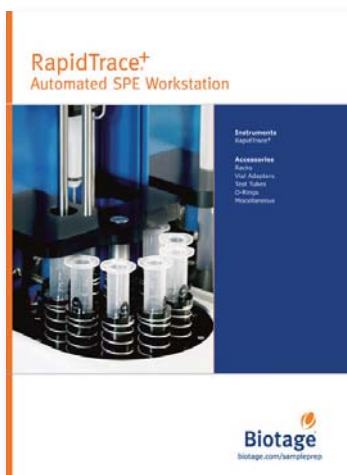

Part Number	Description	Qty.
Racks (Each RapidTrace+ requires one rack)		
C50974	13 x 100 mm Sample Tube 12 x 75 Fraction Tube	1
C50976	13 x 100 mm Sample Tube 12 x 75 Fraction Tube	Pk/5
C58309	13 x 100 mm Sample Tube 13 x 100 mm Fraction Tube	1
C56786	16 x 100 mm Sample Tube 16 x 100 mm Fraction Tube	1
C56536	13 x 100 mm Sample Tube 12 x 75 mm Fraction Tube (9 Position Chilled Rack)	1
C133968	40 mL Sample Flask (5 position) 16 x 100 Fraction Tube (10 position)	1

Evaporation

Biotage® SPE Dry 96 & 96 Dual

Sample Concentrator System

Biotage® SPE Dry 96 Sample Concentrator System

Part Number	Description
SD-9600-DHS-NA	SPE Dry 96 Sample Concentrator System, 100/120V
SD-9600-DHS-EU	SPE Dry 96 Sample Concentrator System, 220/240V
SD-9600-DHS-T-NA	SPE Dry 96 Sample Concentrator System, with PTFE coated needles (Top Head Only), 100/120V
SD-9600-DHS-T-EU	SPE Dry 96 Sample Concentrator System, with PTFE coated needles (Top Head Only), 220/240V


Biotage® SPE Dry 96 Dry Dual Sample Concentrator System

Part Number	Description
SD2-9600-DHS-NA	SPE Dry 96 Dual Sample Concentrator System, 100/120V
SD2-9600-DHS-EU	SPE Dry 96 Dual Sample Concentrator System, 220/240V
SD2-9600-DHS-T-NA	SPE Dry 96 Dual Sample Concentrator System, with PTFE coated needles (Top Head Only), 100/120V
SD2-9600-DHS-T-EU	SPE Dry 96 Dual Sample Concentrator System, with PTFE coated needles (Top Head Only), 220/240V



Evaporation Brochure

For further information please visit www.biotage.com to request a copy of our Evaporation brochure.

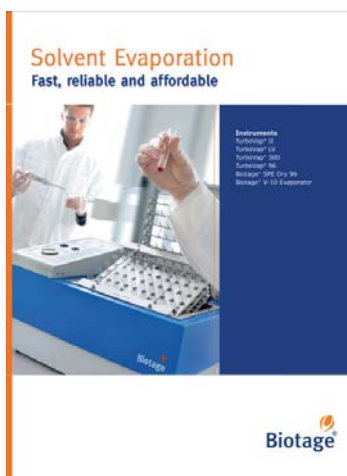

TurboVap® Product Family



TurboVap 500



TurboVap 96



TurboVap II



TurboVap LV

Part Number	Description
-------------	-------------

TurboVap 500

C103202	TurboVap 500, 100/120V, 50/60Hz
C103203	TurboVap 500, 220/240V, 50/60Hz

TurboVap 96

C103263	TurboVap 96, 100/120V
C103264	TurboVap 96, 220/240V

TurboVap II

C103186	TurboVap II, 200 mL with 0.5 mL Endpoint Stem, 100/120V
C103187	TurboVap II, 200 mL with 1.0 mL Endpoint Stem, 100/120V
C103188	TurboVap II, 50 mL with 0.5 mL Endpoint Stem, 100/120V
C103189	TurboVap II, 50 mL with 1.0 mL Endpoint Stem, 100/120V
C103190	TurboVap II, 200 mL with 0.5 mL Endpoint Stem, 220/240V
C103192	TurboVap II, 200 mL with 1.0 mL Endpoint Stem, 220/240V
C103193	TurboVap II, 50 mL with 0.5 mL Endpoint Stem, 220/240V
C103194	TurboVap II, 50 mL with 1.0 mL Endpoint Stem, 220/240V
C103196	TurboVap II, 40 mL Dionex ASE Vial Compatible, 100/120V
C103197	TurboVap II, 40 mL Dionex ASE Vial Compatible, 220/240V

TurboVap LV

C103198	TurboVap LV without Racks, 100/120V
C103199	TurboVap LV without Racks, 220/240V
C103200	TurboVap LV, ASE, 40/60 mL vial, 100/120V
C103201	TurboVap LV, ASE, 40/60 mL vial, 220/240V
C112175	TurboVap LV, PTFE, 50 position, 100/120V
C133718	TurboVap LV, PTFE, 50 position, 220/240V
C55941	TurboVap LV, PTFE-Coated Upper Manifold Replacement Kit
C109301	TurboVap LV, Assy-Water Bath, PTFE Coated

Tube Racks

C48950	Tube Rack (for 10 x 75 mm tubes P/N C48985) 250 tubes included
C44577	Tube Rack (for 12 x 75 mm tubes P/N C44651) 250 tubes included
C45286	Tube Rack (for 13 x 100 mm tubes P/N C40707) 250 tubes included
C44139	Tube Rack (for 16 x 100 mm tubes P/N C40708) 250 tubes included
C45269	Tube Rack (for 16 x 125 mm tubes P/N C45273) 250 tubes included
C44282	Tube Rack (for 20 x 150 mm tubes P/N C40709) 250 tubes included
C44283	Tube Rack (16 x 100 mm-Prelude compatible)
C44880	Tube Rack (for 15 mL Conical Centrifuge Tubes P/N C44941) 125 tubes included
C47818	Tube Rack (15 mL Autotrace tubes with conical bottom, suitable for tubes C47816)
C47820	Tube Rack (10 mL conical tubes, suitable for tubes C47811)
C48928	Tube Rack (1.5-2.0 mL flip cap microcentrifuge tubes)
C48929	Tube Rack (1.5-2.0 mL screw cap microcentrifuge tubes)
C59546	Tube Rack (GC vials)
C61345	Tube Rack (ASE 40/60 mL IChem Vials)

For further information please visit www.biotage.com to request a copy of our Evaporation brochure.

Index

101 (Polystyrene-divinylbenzene)	
SPE columns	20
1-Point-Support™	5
96-well accessories	48
96-well filtration plates	48

A

Accessories, SPE and filtration	45–48
Adaptors, polyethylene, for SPE columns	46
Adaptors, PTFE, for SPE columns	46
AL (Alumina) SPE columns	34
Aminopropyl, see NH ₂	
Applications	4
Aromatic sulfonic acid, see SCX, SCX-3	
Array – SPE wells and 96-well plates	11, 48
Automated Sample Processing Products	55

B

Benzenesulfonic acid, see SCX, SCX-3	
Bulk sorbents, ISOLUTE®	39
Butyl, see C4	

C

C2, C2(EC) (Ethyl)	
96-well SPE plates	24
C4 (Butyl)	
96-well SPE plates	23
SPE columns	23
C8, C8(EC) (Octyl)	
96-well SPE plates	22–23
Bulk sorbent	40
SPE columns	22–23
C18, C18(EC) (n-Octadecyl)	
96-well SPE plates	21–22
Bulk sorbent	40
SPE columns	21–22
Caps, SPE column	46
Carboxypropyl acid, see CBA	
CBA (Carboxypropyl acid)	
96-well SPE plates	31
Bulk sorbent	40
SPE columns	31

CN, CN(EC) (Cyanopropyl)	
96-well SPE plates	25
SPE columns	25
Collection plates, 96-well	48
Column adaptors	46
Cyanopropyl, see CN, CN(EC)	
Cyclohexyl, see CH(EC)	

D

Deep well collection plates	48
Depth filter reservoirs	46
DIOL (2,3-Dihydroxypropoxypropyl)	
SPE columns	33
Drying adaptor, VacMaster	53

E

Empty reservoirs	47
ENV+ (Hydroxylated polystyrene- divinylbenzene)	
96-well SPE plates	20
SPE columns	20
EPH SPE columns	35
Ethyl, see C2, C2(EC)	
Ethylbenzene sulfonic acid, see SCX-3	
Evaporation	57
EVOLUTE®	4, 15–18

F

Filtration columns	47
Filtration plates	48
Florisil SPE columns	34
Frits, polyethylene	47

G

Gravity Rack for sample processing under gravity	52
---	----

H

HAX	
96-well SPE plates	28
SPE columns	28
HGX	
96-well SPE plates	26
SPE columns	26
HGX-3	
96-well SPE plates	26
SPE columns	26

HGX-5	
96-well SPE plates	27
SPE columns	27
HGX-Q	
96-well SPE plates	27
SPE columns	27
Hexyl, see C6	
HM-N (diatomaceous earth)	
Bulk material	40
Packed columns	34

I

Ion-exchange sorbents	40
Ion-exchange SPE columns and 96-well plates	29–32
ISOLUTE®-96 SPE plates 11–14, 19–32	
ISOLUTE Array wells and SPE plates	11, 19–32
ISOLUTE Bulk sorbents	40
ISOLUTE EPH	35
ISOLUTE HM-N (diatomaceous earth)	
Bulk sorbent	40
Packed columns	34
ISOLUTE Myco SPE columns	37
ISOLUTE SLE+	4, 13
ISOLUTE SPE columns see individual sorbent chemistries. <i>For example,</i> <i>for ISOLUTE C18, see C18</i>	

L

Large volume extraction (LVE) kits	53
Layered columns	35

M

Manual sample processing products and accessories	49–54
Method development assistance	4
MFC18 (n-Octadecyl- monofunctional)	
96-well SPE plates	22
SPE columns	22
Mixed-mode SPE columns	26–28
Multimode	28

N

Needles – PTFE and stainless steel	51
NH2 (Aminopropyl)	
<i>96-well SPE plates</i>	29
<i>Bulk sorbent</i>	40
<i>SPE columns</i>	29
Non-polar sorbents	40
Non-polar ISOLUTE® SPE columns	
and 96-well SPE plates	20–25

O

Octadecyl,	
see C18, C18(EC), MFC18	
Octyl, see C8, C8(EC)	
Oil & Grease column – (O&G)	36

P

Part numbering system,	
SPE column	10
Part Number Index	62
PAH SPE columns	35
PE-AX (Pre-equilibrated anion exchanger)	
<i>96-well SPE plates</i>	30
<i>SPE columns</i>	30
PH, PH(EC) (Phenyl)	
<i>96-well SPE plates</i>	25
<i>SPE columns</i>	25
Phenyl, see PH, PH(EC)	
Plate sealing accessories	48
Polar sorbents	40
Polar SPE columns and	
96-well SPE plates	33, 34
Polyethylene frits	47
Polypropylene reservoirs	47
Polystyrene-divinylbenzene	
see EVOLUTE® ABN, 101 and ENV+	
Primary secondary amine, see PSA	
PRESSURE+ 48 and PRESSURE+ 96	
Positive pressure manifolds	50
Protein precipitation plates	41
Propylsulfonic acid, see SCX-2 (PRS)	
PSA (Primary secondary amine)	
<i>96-well SPE plates</i>	29
<i>Bulk sorbent</i>	40
<i>SPE columns</i>	29
PTFE adaptors	46
PTFE needles: VacMaster	51

PTFE stopcocks: VacMaster	51
PTFE stopcock / needles:	
VacMaster	51
PTFE stopcocks	
(for universal use)	51

Q

Quality issues	7
Quaternary amine,	
see PE-AX and SAX	
QuEChERS	43
QuickStart Guide to SPE	4

R

Rack, Gravity	54
Rack TurboVap®	59
RapidTrace®+	56
Rack for VacMaster	51
Reservoirs	47

S

Sample Preparation	
Format Options	10–11
SAX (Quaternary amine)	
<i>96-well SPE plates</i>	30
<i>SPE columns</i>	30
SAX / PSA SPE columns	35
SCX (Benzenesulfonic acid)	
<i>96-well SPE plates</i>	31
<i>SPE columns</i>	31
SCX-2 (Propylsulfonic acid)	
<i>96-well SPE plates</i>	32
<i>Bulk sorbent</i>	40
<i>SPE columns</i>	32
SCX-3 (Ethylbenzene sulfonic acid)	
<i>96-well SPE plates</i>	32
<i>SPE columns</i>	32
SI (Silica)	
<i>Bulk sorbent</i>	40
<i>SPE columns</i>	33
SLE+ plates and columns	14
Sodium sulfate	
drying cartridge	47
Solvent evaporation	57–59
Solvent trap	53
Sorbents, Bulk ISOLUTE	40
SPE Dry 96 sample	
concentrator systems	58

Stainless steel needles:	
VacMaster	51
Stopcock needles, PTFE	51
Stopcocks, PTFE for universal use	51
Strong anion exchange,	
see SAX, PE-AX	
Strong cation exchange,	
see SCX, -2, -3	
Supported Liquid Extraction (SLE),	
see HM-N columns and SLE+ plates	
Styrene divinyl benzene, SDVB	
see EVOLUTE ABN, ENV+ and 101	

T

Tabless columns	16, 17, 20, 21, 26, 30, 32, 38
Technical support	4
Trap Kit	
VacMaster	53
TPH SPE columns	35
TurboVap	59
Tube Rack TurboVap	59

U

Universal PTFE stopcock	51
-------------------------	----

V

VacMaster -10,-20	51
VacMaster-96	52
VacMaster accessories	51–53
VacMaster drying adaptor	53
VacMaster LVE kit	53
VacMaster sample	
processing manifolds	51–53
Vacuum control units	
VCU-1	52
VCU-2 (with integral vacuum generator)	52
Vacuum manifold,	
see VacMaster and VacMaster-96	

W

Weak anion exchange,	
see NH2, PSA and EVOLUTE	
Weak cation exchange,	
see CBA XL SPE columns and EVOLUTE	

Part Number Index

Part Number	Page	120-2040-T	42	220-0005-A	21	221-0100-T	21
101-0010-B	20	121-0001	51, 54	220-0005-G	21	221-0200-C	21
101-0020-B	20	121-0002	51, 54	220-0010-A	21	221-0200-D	21
101-0020-C	20	121-0003	51, 54	220-0010-B	21	221-0500-E	21
101-0050-B	20	121-0004	51, 54	220-0010-BG	21	221-1000-F	21
101-0050-C	20	121-0005	51	220-0010-G	21	240-0005-A	22
120-1000-P01	42, 48	121-0009	51, 54	220-0020-B	21	240-0005-G	22
120-1001-A	47	121-0009-S	51	220-0020-C	21	240-0010-A	22
120-1002-B	47	121-0010	51	220-0020-H	21	240-0010-B	22
120-1003-C	47	121-1010	51	220-0025-P01	21	240-0010-G	22
120-1003-CD	46	121-1012	51	220-0025-R	21	240-0020-B	22
120-1006-D	47	121-1016	51	220-0050-B	21	240-0020-H	22
120-1007-E	47	121-1030	51	220-0050-BG	21	240-0025-P01	22
120-1009-F	47	121-1039	51	220-0050-C	21	240-0050-B	22
120-1009-FD	46	121-1045	51	220-0050-H	21	240-0050-C	22
1201-0120	46	121-1132	51	220-0050-L	21	240-0050-H	22
1201-0121-A	46	121-1133	51	220-0050-P01	21	240-0100-C	22
1201-0122-B	46	121-1134	51	220-0050-R	21	240-0200-D	22
1201-0123-C	46	121-2010	51	220-0100-B	21	240-0500-E	22
1201-0125-H	46	121-2012	51	220-0100-C	21	290-0002-A	22
1201-0126-D	46	121-2016	51	220-0100-CG	21	290-0005-A	22
1201-0127-E	46	121-2059	51	220-0100-P01	21	290-0010-A	22
1201-0128-F	46	121-2068	51	220-0100-T	21	290-0010-B	22
120-1021-A	47	121-2075	51	220-0200-C	21	290-0010-C	22
120-1022-B	47	121-2090	53	220-0200-D	21	290-0020-B	22
120-1022-P01	48	121-2095	53	220-0500-E	21	290-0020-H	22
120-1023-C	47	121-2161	51	220-1000-F	21	290-0025-P01	22
120-1024-H	47	121-2162	51	221-0002-A	21	290-0025-R	22
120-1025-D	47	121-2163	51	221-0005-A	21	290-0050-B	22
120-1026-E	47	121-5202	48, 50, 52	221-0005-G	21	290-0050-C	22
120-1028-F	47	121-5203	48, 50, 52	221-0010-A	21	290-0050-P01	22
120-1031-A	47	121-5204	48	221-0010-B	21	290-0050-R	22
120-1032-G	47	121-5205	48	221-0010-C	21	290-0100-C	22
120-1033-B	47	121-9600	52	221-0010-G	21	290-0100-P01	22
120-1034-H	47	121-9601	52	221-0020-B	21	290-0100-T	22
120-1035-C	47	121-9602	52	221-0020-C	21	290-0200-D	22
120-1036-D	47	121-9610	52	221-0020-H	21	291-0010-A	23
120-1037-E	47	121-9611	52	221-0025-P01	21	291-0010-B	23
120-1038-F	47	121-9612	52	221-0025-R	21	291-0010-G	23
120-1100	46	121-9613	52	221-0050-B	21	291-0020-B	23
120-1101	46	121-9614	52	221-0050-BS	21	291-0020-H	23
120-1102	46	121-9615	52	221-0050-C	21	291-0025-P01	23
120-1103	46	123-2016	54	221-0050-CD	21	291-0025-R	23
120-1104	46	123-2019	54	221-0050-H	21	291-0050-B	23
120-1200	42, 48	124-1001	53	221-0050-P01	21	291-0050-C	23
120-1201	42, 48	124-2001	53	221-0050-R	21	291-0050-P01	23
120-1202	42, 48	150-0006-BG	38	221-0100-B	21	291-0100-C	23
120-2040-P01	42	220-0002-A	21	221-0100-C	21	291-0100-P01	23
120-2040-R	42	220-0002-AG	21	221-0100-P01	21	291-0200-D	23

291-0500-E	23	390-0100-C	23	470-0200-C	29	520-0002-A	31
320-0002-A	24	420-0010-A	25	470-0200-D	29	520-0005-A	31
320-0005-A	24	420-0010-B	25	480-0010-A	29	520-0010-A	31
320-0010-A	24	420-0020-B	25	480-0010-B	29	520-0010-B	31
320-0010-B	24	420-0025-P01	25	480-0020-B	29	520-0020-B	31
320-0010-G	24	420-0025-R	25	480-0025-P01	29	520-0025-P01	31
320-0020-B	24	420-0050-B	25	480-0025-R	29	520-0025-R	31
320-0025-P01	24	420-0050-C	25	480-0050-B	29	520-0050-B	31
320-0025-R	24	420-0100-C	25	480-0050-C	29	520-0050-C	31
320-0050-B	24	420-0100-P01	25	480-0100-C	29	520-0050-P01	31
320-0050-C	24	420-0100-T	25	480-0100-P01	29	520-0050-R	31
320-0050-P01	24	430-0010-A	33	480-0100-T	29	520-0100-C	31
320-0050-R	24	430-0010-B	33	500-0002-A	30	520-0100-P01	31
320-0100-C	24	430-0010-G	33	500-0005-A	30	520-0200-D	31
320-0100-P01	24	430-0020-B	33	500-0005-AG	30	530-0005-A	31
320-0200-D	24	430-0050-B	33	500-0005-G	30	530-0010-A	31
321-0005-A	24	430-0050-C	33	500-0010-A	30	530-0010-B	31
321-0010-A	24	430-0050-H	33	500-0010-B	30	530-0020-B	31
321-0010-B	24	430-0100-C	33	500-0010-C	30	530-0025-P01	31
321-0010-G	24	460-0002-A	33	500-0010-G	30	530-0050-B	31
321-0020-B	24	460-0010-A	33	500-0020-B	30	530-0050-C	31
321-0025-P01	24	460-0010-B	33	500-0020-H	30	530-0050-H	31
321-0025-R	24	460-0010-G	33	500-0025-P01	30	530-0050-P01	31
321-0050-B	24	460-0020-B	33	500-0050-B	30	530-0050-R	31
321-0050-C	24	460-0020-H	33	500-0050-BG	30	530-0100-B	31
321-0050-H	24	460-0050-B	33	500-0050-C	30	530-0100-C	31
321-0050-R	24	460-0050-C	33	500-0050-H	30	530-0100-P01	31
321-0100-C	24	460-0050-H	33	500-0050-P01	30	530-0100-T	31
321-0100-T	24	460-0100-B	33	500-0050-R	30	530-0200-D	31
360-0002-A	25	460-0100-C	33	500-0100-B	30	532-0002-A	32
360-0005-A	25	460-0200-C	33	500-0100-C	30	532-0005-A	32
360-0010-A	25	460-0200-D	33	500-0100-CG	30	532-0010-A	32
360-0010-B	25	460-0200-L	33	500-0100-P01	30	532-0010-B	32
360-0025-P01	25	460-0500-E	33	500-0100-T	30	532-0020-B	32
360-0025-R	25	460-1000-F	33	500-0200-D	30	532-0025-P01	32
360-0050-B	25	470-0002-A	29	503-0002-A	30	532-0025-R	32
360-0050-C	25	470-0005-A	29	503-0005-A	30	532-0050-B	32
360-0050-P01	25	470-0005-G	29	503-0010-A	30	532-0050-C	32
360-0050-R	25	470-0010-A	29	503-0010-B	30	532-0050-H	32
360-0100-C	25	470-0010-B	29	503-0010-C	30	532-0050-P01	32
360-0100-P01	25	470-0010-G	29	503-0020-B	30	532-0100-B	32
360-0100-T	25	470-0020-B	29	503-0025-P01	30	532-0100-C	32
390-0010-A	23	470-0025-R	29	503-0050-B	30	532-0100-CG	32
390-0010-B	23	470-0050-B	29	503-0050-C	30	532-0100-P01	32
390-0020-B	23	470-0050-C	29	503-0050-R	30	532-0200-D	32
390-0025-P01	23	470-0050-H	29	503-0100-C	30	533-0002-A	32
390-0025-R	23	470-0050-P01	29	503-0100-P01	30	533-0005-A	32
390-0050-B	23	470-0100-C	29	503-0100-T	30	533-0010-A	32
390-0050-C	23	470-0100-P01	29	503-0200-D	30	533-0010-B	32

533-0020-B	32	610-0010-B	16	713-0020-A	34	820-1420-F	14
533-0025-P01	32	610-0010-H	16	713-0020-B	34	902-0002-A	26
533-0025-R	32	610-0020-B	16	713-0050-B	34	902-0002-AG	26
533-0050-B	32	610-0020-C	16	713-0050-C	34	902-0005-A	26
533-0050-C	32	610-0050-C	16	713-0100-B	34	902-0010-A	26
533-0050-H	32	611-0005-B	16	713-0100-C	34	902-0010-B	26
533-0050-P01	32	611-0005-BG	16	713-0200-D	34	902-0013-A	26
533-0050-R	32	611-0010-B	16	713-0500-E	34	902-0013-B	26
533-0100-B	32	611-0010-BG	16	714-0010-A	34	902-0013-BG	26
533-0100-C	32	611-0010-H	16	714-0020-A	34	902-0013-C	26
533-0100-P01	32	611-0020-B	16	714-0020-B	34	902-0013-CG	26
533-0100-T	32	611-0020-C	16	714-0050-B	34	902-0013-H	26
533-0200-D	32	611-0050-C	16	714-0050-C	34	902-0020-B	26
600-0001-A	16	611-0050-CG	16	714-0100-B	34	902-0020-H	26
600-0002-A	16	612-0005-B	17	714-0100-C	34	902-0025-P01	26
600-0002-AG	16	612-0010-B	17	714-0200-C	34	902-0025-R	26
600-0002-B	16	612-0010-H	17	714-0200-D	34	902-0030-B	26
600-0002-H	16	612-0020-C	17	714-0500-E	34	902-0030-C	26
600-0010-PX01	16	612-0050-C	17	714-1000-F	34	902-0030-CG	26
600-0010-R	16	612-0050-CG	17	715-0010-A	34	902-0030-H	26
600-0010-RP	16	613-0005-B	17	715-0020-A	34	902-0050-C	26
600-0025-R	16	613-0010-B	17	715-0020-B	34	902-0050-P01	26
600-0025-RP	16	613-0010-BG	17	715-0050-B	34	902-0050-R	26
600-0030-PX01	16	613-0010-H	17	715-0050-C	34	902-0050-T	26
601-0001-A	16	613-0020-C	17	715-0100-B	34	902-0100-P01	26
601-0001-AG	16	613-0050-C	17	715-0100-C	34	902-0100-T	26
601-0002-A	16	613-0050-CG	17	715-0200-D	34	903-0005-A	28
601-0002-AG	16	614-0005-B	18	715-0500-E	34	903-0010-A	28
601-0010-PX01	16	614-0010-B	18	715-1000-F	34	903-0020-B	28
601-0010-R	16	614-0010-H	18	752-0100-C	35	903-0020-C	28
601-0010-RP	16	614-0020-C	18	753-0100-CD	36	903-0020-H	28
601-0025-R	16	614-0050-C	18	753-0300-FD	36	903-0025-P01	28
601-0025-RP	16	650-0010-PX01	18	800-0040-BM	34	903-0025-R	28
601-0030-PX01	16	650-0030-PX01	18	800-0100-CM	34	903-0100-P01	28
602-0002-A	17	712-0010-A	34	800-0220-DM	34	904-0010-A	28
602-0010-PX01	17	712-0020-A	34	800-0350-EM	34	904-0030-B	28
602-0025-R	17	712-0020-B	34	800-0700-FM	34	904-0030-C	28
602-0030-PX01	17	712-0020-H	34	800-1300-FM	34	904-0050-B	28
603-0002-A	17	712-0050-B	34	802-0250-M	47	904-0100-C	28
603-0002-AG	17	712-0050-BG	34	820-0055-B	14	905-0002-A	26
603-0010-PX01	17	712-0050-C	34	820-0055-BG	14	905-0005-A	26
603-0025-R	17	712-0050-H	34	820-0140-C	14	905-0010-A	26
603-0030-PX01	17	712-0050-L	34	820-0140-CG	14	905-0025-P01	26
604-0002-A	18	712-0100-B	34	820-0200-P01	14	905-0025-R	26
604-0010-PX01	18	712-0100-C	34	820-0200-T	14	905-0100-P01	26
604-0025-R	18	712-0200-D	34	820-0290-D	14	905-0100-T	26
604-0030-PX01	18	712-0500-E	34	820-0400-P01	14	906-0002-A	27
610-0005-B	16	712-2000-F	34	820-0400-T	14	906-0005-A	27
610-0005-BG	16	713-0010-A	34	820-0690-E	14	906-0010-A	27

906-0010-G	27	933-0100-C	35	C44577	59	SD-9600-DHS-NA	58
906-0013-H	27	934-0040-C	35	C44880	59	SD-9600-DHS-T-EU	58
906-0025-P01	27	935-0040-C	35	C45269	59	SD-9600-DHS-T-NA	58
906-0025-R	27	941-0020-B	35	C45286	59	SER-PPM-FYMP	50
906-0100-P01	27	9460-0100	40	C47818	59	SER-PPM-IN	50
906-0100-T	27	9470-0100	40	C47820	59	SER-PPM-PL	50
915-0001-A	20	9480-0100	40	C48928	59	SER-PPM-SLVR	50
915-0002-A	20	9500-0100	40	C48929	59		
915-0002-G	20	9503-0100	40	C48950	59		
915-0005-A	20	9520-0100	40	C50000	56		
915-0005-AG	20	9532-0100	40	C50974	56		
915-0005-B	20	9533-0100	40	C50976	56		
915-0005-G	20	9712-0100	40	C52006	56		
915-0010-A	20	976-0400-C	36	C52689	56		
915-0010-B	20	9800-1000	40	C55941	59		
915-0010-BG	20	9800-5000	40	C56536	56		
915-0010-C	20	986-0002-A	27	C56786	56		
915-0010-H	20	986-0005-A	27	C58309	56		
915-0010-P01	20	986-0010-A	27	C59546	59		
915-0010-R	20	986-0025-P01	27	C61345	59		
915-0020-B	20	986-0025-R	27	PPM-48	50		
915-0020-C	20	986-0100-P01	27	PPM-96	50		
915-0020-CD	20	986-0100-T	27	PPM-A48-1232	50		
915-0020-CG	20	C103186	59	PPM-A48-1275	50		
915-0025-P01	20	C103187	59	PPM-A48-13100	50		
915-0025-R	20	C103188	59	PPM-A48-16100	50		
915-0050-B	20	C103189	59	PPM-A48-1RCK	50		
915-0050-BG	20	C103190	59	PPM-A48-3RCK	50		
915-0050-C	20	C103192	59	PPM-A48-6RCK	50		
915-0050-D	20	C103193	59	PPM-A48-GSKT	50		
915-0050-L	20	C103194	59	PPM-A48-WST	50		
915-0050-P01	20	C103196	59	PPM-A96-1024	50		
915-0100-C	20	C103197	59	PPM-A96-GSKT	50		
915-0100-CG	20	C103198	59	PPM-GA	50		
915-0100-E	20	C103199	59	Q0000-50V	44		
9220-0100	40	C103200	59	Q0010-15V	44		
9221-0100	40	C103201	59	Q0020-15V	44		
924-0050-C	35	C103202	59	Q0030-15V	44		
924-0100-C	35	C103203	59	Q0035-15V	44		
927-0075-B	35	C103263	59	Q0050-15V	44		
927-0150-C	35	C103264	59	Q0060-15V	44		
927-0150-CD	35	C109301	59	Q0070-15V	44		
928-0145-B	35	C112175	59	Q0080-15V	44		
928-0500-E	35	C125713	56	Q0090-15V	44		
9290-0100	40	C133718	59	SD2-9600-DHS-EU	58		
9291-0100	40	C133968	56	SD2-9600-DHS-NA	58		
9320-0100	40	C44139	59	SD2-9600-DHS-T-EU	58		
9321-0100	40	C44282	59	SD2-9600-DHS-T-NA	58		
933-0050-B	35	C44283	59	SD-9600-DHS-EU	58		

Limited Warranty

Biotage warrants that Biotage Consumables, including all FLASH®, FLASH+®, ISOLUTE®, EVOLUTE® and Samplet® purification cartridges, scavenger resins, solid-bound reagents and microwave vials, will be of good quality and workmanship, and will meet the applicable product specification. This warranty applies only to the initial test performed at the Customer's facility upon the initial start-up of the consumables and expires at the time the user applies an actual sample to the consumable. If the cartridge is packed with media provided by the Customer, the Biotage warranty applies only to the plastic tube, frits, and labor required for packing and testing the cartridge. Biotage will not be liable for any damage to media provided by the Customer that may be caused when Biotage packs such media in accordance with Biotage's standard operating procedures.

Should any Biotage consumable fail to meet the limited warranty above after being tested in accordance with the applicable Biotage standard operating procedures, Biotage will provide, at its sole option, either a replacement cartridge or reaction vial. If a cartridge, it will be packed with the original media, or packed with new media, at no cost to the Purchaser. If such failed cartridge is packed with media provided by the Purchaser, Biotage will make a reasonable effort to re-pack the original media, or pack a replacement cartridge with new media provided by the Purchaser.

Returns of any Consumable must be authorized in advance. Please contact Biotage for a Return Authorization (RA) number and shipping instructions. All claims must be made within thirty (30) days of shipment from Biotage, or after the initial test of the cartridge at the Customer's facility, whichever is first.

BIOTAGE MAKES NO OTHER WARRANTY, EXPRESSED OR IMPLIED, INCLUDING MERCHANTABILITY OR SUITABILITY FOR ANY PURPOSE.

BIOTAGE WILL NOT BE LIABLE IN ANY EVENT FOR SPECIAL, INDIRECT OR CONSEQUENTIAL DAMAGES WHETHER ARISING IN TORT, UNDER ANY WARRANTY OR OTHERWISE.

Trademark Acknowledgements

©2014 Biotage

The following trademarks are owned by Biotage AB: Advancer, Advancer 350, Advancer Kilobatch, AFFINILUTE, Biotage, Biotage ZIP, Endeavor, EVOLUTE, EVOLUTE EXPRESS, ExploraSep, Extrahera,  Firefly design, FLASH+, FlashMaster, FlashVac, Flash 75, Flash 150, Flash 400, Horizon, HPFC, HP-SIL, HP-Sphere, Initiator, Initiator Peptide Workstation, Initiator+, Initiator+ Alstra, Initiator+ Robot 60, Initiator+ Robot 8, Initiator+ SP Wave, Isolera, Isolera Dalton, Isolera Dalton Mass Detector, Isolera Dalton Nanolink, Isolera Dalton System, Isolera Four, Isolera LS, Isolera One, Isolera Prime, Isolera Spektra, Isolera Spektra Four, Isolera Spektra LS, Isolera Spektra One, ISOLUTE, ISOLUTE Myco, ISOLUTE QuEChERS, IST,  IST design, KILOPREP, KP-C18-HS, KP-C18-WP, KP-C4-WP, KP-NH, KP-Sil, KP-Sphere, MIP Rule of 6, MIP[4]Process, MIP[4]Proteins, MIP4SPE, PathFinder, PRESSURE+, PRESSURE+ 48, PRESSURE+ 96, RapidTrace, RapidTrace+, RENSA, Resolux, Robot 60, Robot 8, Samplet, SIM, SNAP, SNAP Ultra, SNAP XL, SP Wave, SP1, SP4, SPE Dry, SPE Dry 96, SPE Dry 96 Dual, SPx, Syro Wave, TurboVap, Universal Phase Separator, V-10, VacMaster, ZIF, ZIF-SIM, ZIP-Sphere, 1-Point Support.

Other product and company names mentioned herein may be trademarks or registered trademarks and/or service marks of their respective owners, and are used only for explanation and to the owners' benefit, without intent to infringe.

Notes



EVOLUTE® EXPRESS



ISOLUTE® SLE+



Biotage® PRESSURE+ 96 & 48

Your Complete Partner for Effective Chemistry

Biotage is a worldwide supplier of instruments and accessories designed to facilitate the work of laboratory chemists. With our deep knowledge of the industry, academic contacts and in-house R&D teams, we can deliver the best solutions to your challenges. We take great pride in our flexibility and ability to meet our customer's individual needs. With strong foundations in both analytical and organic chemistry, we can offer the widest range of solutions available on the market.

EUROPE

Main Office: +46 18 565900
Toll Free: +800 18 565710
Fax: +46 18 591922
Order Tel: +46 18 565710
Order Fax: +46 18 565705
order@biotage.com
Support Tel: +46 18 56 59 11
Support Fax: + 46 18 56 57 11
eu-1-pointsupport@biotage.com

NORTH & LATIN AMERICA

Main Office: +1 704 654 4900
Toll Free: +1 800 446 4752
Fax: +1 704 654 4917
Order Tel: +1 704 654 4900
Order Fax: +1 434 296 8217
ordermailbox@biotage.com
Support Tel: +1 800 446 4752
Outside US: +1 704 654 4900
us-1-pointsupport@biotage.com

JAPAN

Tel: +81 3 5627 3123
Fax: +81 3 5627 3121
jp_order@biotage.com
jp-1-pointsupport@biotage.com

CHINA

Tel: +86 21 2898 6655
Fax: +86 21 2898 6153
cn_order@biotage.com
cn-1-pointsupport@biotage.com

To locate a distributor,
please visit our website at
www.biotage.com