



Revision No.: ZJ0004
Issue Date: Jul 1st, 2012

Enterovirus Real Time RT-PCR Kit

User Manual

For In Vitro Diagnostic Use Only

REF QR-0206-01

For use with LightCycler1.0/2.0 Instrument



Obelis S.A.
Boulevard Général Wahis 53
1030 Brussels, BELGIUM
Tel: +(32) 2.732.59.54
Fax: +(32) 2.732.60.03
E-Mail : mail@obelis.net



Shanghai ZJ Bio-Tech Co., Ltd.
www.liferiver.com.cn Tel: +86-21-34680596
trade@liferiver.com.cn Fax: +86-21-34680595
2nd floor, No.15 Building, No.188 Xijunhuan Road,
Pujiang Hi-tech Park, Shanghai, China



1. Intended Use

By using real time PCR systems, Enterovirus real time PCR kit is used for the detection of Coxsackie Virus in samples like nasal and pharyngeal secretions, sputum, provoked sputum, stool, C.S.F and etc.

2. Principle of Real-Time PCR

The principle of the real-time detection is based on the fluorogenic 5' nuclease assay. During the PCR reaction, the DNA polymerase cleaves the probe at the 5' end and separates the reporter dye from the quencher dye only when the probe hybridizes to the target DNA. This cleavage results in the fluorescent signal generated by the cleaved reporter dye, which is monitored real-time by the PCR detection system. The PCR cycle at which an increase in the fluorescence signal is detected initially is proportional to the amount of the specific PCR product. Monitoring the fluorescence intensities in real-time allows the detection of the accumulating product without having to re-open the reaction tube after the amplification.

3. Product Description

Enteroviruses are small viruses that are made of ribonucleic acid (RNA) and protein. In addition to the three different polioviruses, there are 61 non-polio enteroviruses that can cause disease in humans: 29 Coxsackieviruses (23 Coxsackie A viruses and 6 Coxsackie B viruses), 28 echoviruses, and 4 other enteroviruses.

Enteroviruses can be found in the respiratory secretions (e.g., saliva, sputum, or nasal mucus) and stool of an infected person. Other persons may become infected by direct contact with secretions from an infected person or by contact with contaminated surfaces or objects, such as a drinking glass or telephone. Parents, teachers, and child care center workers may also become infected by contamination of the hands with stool from an infected infant or toddler during diaper changes.

The Enterovirus real time RT-PCR kit contains a specific ready-to-use system for the detection of the Enterovirus using RT-PCR (Reverse Transcription Polymerase Chain Reaction) in the real-time PCR system. The master contains a Super Mix for the specific amplification of the Enterovirus RNA including CA16, CA4, CA5, CA7, CA9, CA10, CB2, CB5, CB13, echoviruses and EV71. The reaction is done in one step real time RT-PCR. The first step is a reverse transcription (RT), during which the Enterovirus RNA is transcribed into cDNA. Afterwards, a thermostable DNA polymerase is used to amplify the specific gene fragments by means of PCR (polymerase chain reaction). Fluorescence is emitted and measured by the real time systems' optical unit during the PCR. The detection of amplified Enterovirus DNA fragment is performed in fluorimeter **channel 530nm** with the fluorescent quencher BHQ1. DNA extraction buffer is available in the kit and bronchial lavage sample or lung section sample is used for the extraction of the DNA. In addition, the kit contains a system to identify possible PCR inhibition by measuring the **560nm** fluorescence of the internal control (IC). An external positive control (1×10⁷ copies/ml) contained, allows the determination of the gene load. For further information, please refer to section 9.2 Quantitation.

4. Kit Contents

Ref.	Type of reagent	Presentation	25rxns
1	EV Super Mix	1 vial, 350µl	
2	RT-PCR Enzyme Mix	1 vial, 28µl	
3	Molecular Grade Water	1 vial, 400µl	
4	Internal Control	1 vial, 30µl	
5	EV Positive Control (1×10 ⁷ copies/ml)	1 vial, 30µl	

Analysis sensitivity: 5×10³ copies/ml

LOQ: 1×10⁴~1×10⁶ copies/ml

Note: Analysis sensitivity depends on the sample volume, elution volume, nucleic acid extraction methods and other factors. If you use the RNA extraction kits recommended, the analysis sensitivity is the same as it declares. However, when the sample volume is dozens or even hundreds of times greater than elution volume by some concentrating method, it can be much higher.

5. Storage

- All reagents should be stored at -20°C. Storage at +4°C is not recommended.
- All reagents can be used until the expiration date indicated on the kit label.
- Repeated thawing and freezing (> 3x) should be avoided.
- Cool all reagents during the working steps.
- Super Mix and Reaction Mix should be stored in the dark.

6. Additionally Required Materials and Devices

- Biological cabinet
- Vortex mixer
- Cryo-container
- Sterile filter tips for micro pipets
- Disposable gloves, powderless
- Refrigerator and Freezer
- Desktop microcentrifuge for "ependorf" type tubes (RCF max. 16,000 x g)
- Real time PCR system
- Real time PCR reaction tubes/plates
- Pipets (0.5µl – 1000µl)
- Sterile microtubes
- Biohazard waste container
- Tube racks

7. ⚠ Warnings and Precaution

- Carefully read this instruction before starting the procedure.
- For in vitro diagnostic use only.
- This assay needs to be carried out by skilled personnel.
- Clinical samples should be regarded as potentially infectious materials and should be prepared in a laminar flow hood.
- This assay needs to be run according to Good Laboratory Practice.
- Do not use the kit after its expiration date.
- Avoid repeated thawing and freezing of the reagents, this may reduce the sensitivity of the test.
- Once the reagents have been thawed, vortex and centrifuge briefly the tubes before use.
- Prepare quickly the Reaction mix on ice or in the cooling block.
- Set up two separate working areas: 1) Isolation of the RNA/ DNA and 2) Amplification/ detection of amplification products.
- Pipets, vials and other working materials should not circulate among working units.
- Use always sterile pipette tips with filters.
- Wear separate coats and gloves in each area.
- Do not pipette by mouth. Do not eat, drink, smoke in laboratory.

- Avoid aerosols.

8. Sample Collection, Storage and transport

- Collected samples in sterile tubes.
- Specimens can be extracted immediately or frozen at -20°C to -80°C.

9. Procedure

9.1 RNA-Extraction

Different brands of RNA extraction kits are available. You may use your own extraction systems or the commercial kit based on the yield. For RNA extraction kit, please comply with manufacturer's instructions. The recommended extraction kit is as follows:

Nucleic Acid Isolation Kit	Cat. Number	Manufacturer
RNA Isolation Kit	ME-0010/ME-0012	ZJ Biotech
QIAamp Viral RNA Mini Extraction Kit (50)	52904	QIAGEN


9.2 Internal Control

It is necessary to add internal control (IC) in the reaction mix. Internal control (IC) allows the user to determine and control the possibility of RT-PCR inhibition.

Add the internal control (IC) 1µl/rxn and the result will be shown in the **channel 560nm**.

9.3 RT-PCR Protocol

The Master Mix volume for each reaction should be pipetted as follows:

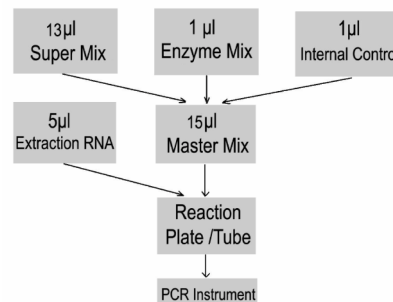

※PCR system without channel 560nm may be treated with 1µl Molecular Grade Water instead of 1µl IC.

- 1) The volumes of Super Mix and Enzyme Mix per reaction multiply with the number of samples, which includes the number of controls, standards, and sample prepared. Molecular Grade Water is used as the negative control. For reasons of unprecise pipetting, always add an extra virtual sample. Mix completely then spin down briefly in a centrifuge.
- 2) Pipet **15µl Master Mix** with micropipets of sterile filter tips to each real time PCR reaction plate/tubes. Separately add **5µl RNA** sample, positive and negative controls to different reaction plate/tubes. Immediately close the plate/tubes to avoid contamination.
- 3) Spin down briefly in order to collect the Master Mix in the bottom of the reaction tubes.
- 4) Perform the following protocol in the instrument:

45°C for 10min	1cycle	Selection of fluorescence channels	
95°C for 15min	1cycle	530nm	Target Nucleic Acid
95°C for 5sec, 60°C for 30sec	40cycles	560nm	IC
(Fluorescence measured at 60°C)			

10. Threshold setting: Choose **Arithmetic** as back ground and **none** as Noise Band method, then adjust the Noise band just above the maximum level of molecular grade water, and adjust the threshold just under the minimum of the positive control.

11. Calibration for quantitative detection: Input each concentration of standard controls at the end of run, and a standard curve will be automatically formed.

12. Quality control:

Negative control, positive control and internal control must be performed correctly, otherwise the sample results is invalid.

Control	Channel	Crossing point value	
		530nm	560nm
Molecular Grade Water		Blank	25~35
Positive Control(qualitative assay)		≤35	—
QS (quantitative detection)		Correlation coefficient of QS curve ≤ -0.98	

13. Data Analysis and Interpretation

The following sample results are possible:

	Crossing point value		Result Analysis
	530nm	560nm	
1#	Blank	25~35	Below the detection limit or negative
2#	≤38	—	Positive; and the software displays the quantitative value
3#	38~40	25~35	Re-test; if it is still 38~40, report as 1#
4#	Blank	Blank	RT-PCR Inhibition; no diagnosis can be concluded.

For further questions or problem, please contact our technical support at trade@liferiver.com.cn