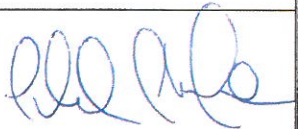
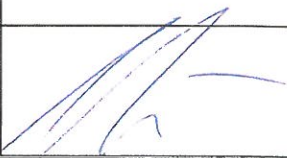


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 1 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

User Manual Cato[®] HL7 2.32 Interface

Creation / Last Modification <i>Erstellung / letzte Bearbeitung</i>		Review <i>Review</i>	
Created / Modified by: <i>Erstellt / bearbeitet von:</i>	Stefan Schuh	Reviewed by: <i>Review durch:</i>	Samuel Gonzalo
Date: <i>Datum:</i>	11 Mar 2014	Date: <i>Datum:</i>	11 Mar 2014
Signature: <i>Unterschrift:</i>		Signature: <i>Unterschrift:</i>	

Release <i>Freigabe</i>	
Released by: <i>Freigegeben von:</i>	Helmut Hörner
Date: <i>Datum:</i>	14 Mar 2014
Signature: <i>Unterschrift:</i>	

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 2 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Change History

Änderungshistorie

Rev. Rev.	ECO ECO	Date Datum	Auth. Autor	Short Description Kurzbeschreibung	Release date Freigabe- datum	Released by Freigeg- eben von
01	-	11 Mar 2014	SS	Created based MWS-SW-USRMAN-CATO-HL7-2.30-EN	14 Mar 2014	HH

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 3 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Table of Contents

1	Preface	5
2	Supported Message Types	6
2.1	ADT Messages	6
2.1.1	Description	6
2.1.2	Simplified Workflow	10
2.2	ORU messages	30
2.2.1	Short description.....	30
2.2.2	Workflow.....	30
2.3	BAR Messages	31
2.3.1	Description	31
2.3.2	Simplified Workflow	32
2.3.3	Supported Segments	33
2.4	PID – Patient Master Data	33
2.5	PV1 – Patient Visits (Ward Assignments)	35
2.6	ZBE – Patient Moving Information.....	37
2.7	MRG – Patient Merge Information	37
2.8	DG1 – Diagnose Information	38
2.9	OBX – Observation / Medical results	39
2.10	OBR – Observation request segment	40
2.11	IN1 – Insurance Segment	40
3	Preparing the cato® database for HL7	41
3.1	Synchronization.....	41
3.1.1	Description	41
3.1.2	Patient synchronization	42
3.1.3	Ward assignment synchronization.....	43
3.1.4	Case synchronization.....	43
3.1.5	Ward synchronization.....	44
3.1.6	Medical results	44
4	Configuration.....	45
4.1	Licenses	45
4.2	The Configuration and Log-Viewer Window	45
4.3	Changing the Configuration	46
4.4	HL7 Settings.....	48
4.4.1	Fundamentals.....	48

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 4 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2	System Settings	52
4.4.3	Database Settings.....	99
4.4.4	Client Settings	110
4.4.5	Connection Settings.....	113
5	Logging Output	125
5.1	Rejected HL7 Messages	125
5.2	Accepted HL7 Messages:	127

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 5 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

1 Preface

This document describes the configuration and operation of the cato[®] HL7 Interface. The interface imports patient relevant data via the open standard HL7 in the version 2.5. For the use with other versions please make sure that the compatibility is proven for all message segments and fields used by the cato[®] HL7 Interface as described in this document.

Further information regarding using the interface for processing RDE Order Messages can be found in Manual **User Manual Cato HL7 2.32 RDE Order Interface**

This document is written for IT-administrators or advanced users.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 6 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2 Supported Message Types

The cato[®] HL7 Interface supports the following message types.

2.1 ADT Messages

2.1.1 Description

Source: *HL7 manual version 2.5, chapter 3 – Patient administration*

2.1.1.1 ADT/A01 – Admit/visit notification

An A01 event is intended to be used for ‘Admitted’ patients only. An A01 event is sent as a result of a patient undergoing the admission process which assigns the patient to a bed. It signals the beginning of a patient’s stay in a healthcare facility. Normally, this information is entered in the primary Patient Administration system and broadcast to the nursing units and ancillary systems.

2.1.1.2 ADT/A02 – Transfer a patient

An A02 event is issued as a result of the patient changing his or her assigned physical location.

The fields included when this message is sent should be the fields pertinent to communicate this event. When other important fields change, it is recommended that the A08 (update patient information) event be used in addition.

2.1.1.3 ADT/A03 – Discharge/end visit

An A03 event signals the end of a patient’s stay in a healthcare facility. It signals that the patient’s status has changed to ‘discharged’ and that a discharge date has been recorded. The patient is no longer in the facility.

The patient’s location prior to discharge should be entered in *PV1-3 - Assigned Patient Location*.

2.1.1.4 ADT/A04 – Register a patient

An A04 event signals that the patient has arrived or checked in as a one-time, or recurring outpatient, and is not assigned to a bed. One example might be its use to signal the beginning of a visit to the Emergency Room (= Casualty, etc.) admissions. *PV1-44 - Admit Date/Time* is used for the visit start date/time.

2.1.1.5 ADT/A05 – Pre-admit a patient

An A05 event is sent when a patient undergoes the pre-admission process.

For example, a pre-admit may be performed prior to inpatient or outpatient surgery so that lab tests can be performed prior to the surgery. This event can also be used to pre-register a non-admitted patient.

2.1.1.6 ADT/A06 – Change an outpatient to an inpatient

An A06 event is sent when a patient who was present for a non-admitted visit is being admitted after an evaluation of the seriousness of the patient’s condition. This event changes a patient’s status from nonadmitted to admitted. The new patient location should appear in *PV1-3 - Assigned Patient Location*, while the old patient location (if different) should appear in *PV1-6 - Prior Patient Location*. The new patient class should appear in *PV1-2 - Patient Class*.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 7 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.1.7 ADT/A07 – Change an inpatient to an outpatient

An A07 event is sent when a patient who was admitted changes his/her status to ‘no longer admitted’ but is still being seen for this episode of care. This event changes a patient from an ‘admitted’ to a ‘nonadmitted’ status. The new patient location should appear in *PV1-3 - Assigned Patient Location*, while the old patient location (if different) should appear in *PV1-6 - Prior Patient Location*.

2.1.1.8 ADT/A08 – Update patient information

This trigger event is used when any patient information has changed but when no other trigger event has occurred. For example, an A08 event can be used to notify the receiving systems of a change of address or a name change.

2.1.1.9 ADT/A11 – Cancel admit / visit notification

For ‘admitted’ patients, the A11 event is sent when an A01 (admit/visit notification) event is cancelled, either because of an erroneous entry of the A01 event, or because of a decision not to admit the patient after all.

For ‘non-admitted’ patients, the A11 event is sent when an A04 (register a patient) event is cancelled, either because of an erroneous entry of the A04 event, or because of a decision not to check the patient in for the visit after all. To cancel an A05 (pre-admit a patient) event, use the A38 (cancel pre-admit), which was new for Version 2.3 of this Standard.

2.1.1.10 ADT/A12 – Cancel transfer

The A12 event is sent when an A02 (transfer a patient) event is cancelled, either because of erroneous entry of the A02 event or because of a decision not to transfer the patient after all. *PV1-3 - Assigned Patient Location* must show the location of the patient prior to the original transfer.

2.1.1.11 ADT/A13 – Cancel discharge / end visit

The A13 event is sent when an A03 (discharge/end visit) event is cancelled, either because of erroneous entry of the A03 event or because of a decision not to discharge or end the visit of the patient after all. *PV1-3-Assigned Patient Location* should reflect the location of the patient after the cancellation has been processed.

2.1.1.12 ADT/A21 – Patient goes on a leave of absence

An A21 event is sent to notify systems that an admitted patient has left the healthcare institution temporarily. It is used for systems in which a bed is still assigned to the patient, and it puts the current admitted patient activities on hold. For example, it is used to notify dietary services and laboratory systems when the patient goes home for the weekend.

2.1.1.13 ADT/A22 – Patient returns from a leave of absence

An A22 event is sent to notify systems that an admitted patient has returned to the healthcare institution after a temporary ‘leave of absence’. It is used for systems in which a bed is still assigned to the patient, and it takes their current admitted patient activities off of ‘hold’ status. For example, it is used to notify dietary services and laboratory systems when the patient returns from a weekend trip to his/her home.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 8 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.1.14 ADT/A28 – Add person or patient information

The purpose of this and the three following messages was to allow sites with multiple systems and respective master patient databases to communicate activity related to a person regardless of whether that person is currently a patient on each system. Each system has an interest in the database activity of the others in order to maintain data integrity across an institution. Though they are defined within the ADT message set, these messages differ in that they are not patient-specific.

2.1.1.15 ADT/A29 – Delete person information

An A29 event can be used to delete all demographic information related to a given person. This event ,undoes‘ an A28 (add person information) event. The information from the A28 event is deleted. This event is used, for example, when adding the information was performed in error, or when another record already exists for the person, or when one wants to purge the person from the database. When this event occurs, all visit and account level data for this person is also purged.

2.1.1.16 ADT/A31 – Update person information

An A31 event can be used to update person information on an MPI. It is similar to an A08 (update patient information) event, but an A08 (update patient information) event should be used to update patient information for a current episode. An A28 (add person information) or A31 can also be used for backloading MPI information for the person, or for backloading person and historical information.

2.1.1.17 ADT/A34 – Merge patient information

Event A34 has been retained for backward compatibility only. From V2.3.1 onwards, event A40 (Merge patient - patient identifier list) should be used instead.

2.1.1.18 ADT/A38 – Cancel Pre-Admit

The A38 event is sent when an A05 (pre-admit a patient) event is cancelled, either because of erroneous entry of the A05 event or because of a decision not to pre-admit the patient after all.

2.1.1.19 ADT/A40 – Merge patient information – PID

A merge has been done at the patient identifier list level. That is, two *PID-3 - Patient Identifier List* identifiers have been merged into one.

An A40 event is used to signal a merge of records for a patient that was incorrectly filed under two different identifiers. The ‘incorrect source identifier’ identified in the MRG segment (*MRG-1 - Prior Patient Identifier List*) is to be merged with the required ‘correct target identifier’ of the same ‘identifier type code’ component identified in the PID segment (*PID-3 - Patient Identifier List*).

2.1.1.20 ADT/A45 – Move visit information

A move has been done at the visit identifier level. That is, a *PV1-19 - Visit Number* or *PV1-50 - Alternate Visit ID* associated with one account identifier (*PID-18 - Patient Account Number*) has been moved to another account identifier.

An A45 event is used to signal a move of records identified by the *MRG-5 - Prior Visit Number* or the *MRG-6 - Prior Alternate Visit ID* from the ‘incorrect source account identifier’ identified in the MRG segment (*MRG-3 - Prior Patient Account Number*) to the ‘correct target account identifier’ identified in the PID segment (*PID-18 - Patient Account Number*).

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 9 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.1.1.21 ADT-A52 – Cancel leave of absence

The A52 event is sent when an A21 (patient goes on ‘leave of absence’) event is cancelled, either because of erroneous entry of the A21 event or because of a decision not to put the patient on ‘leave of absence’ after all.


2.1.1.22 ADT-A53 – Cancel return from a leave of absence

The A53 event is sent when an A22 (patient returns from ‘leave of absence’) event is cancelled, either because of erroneous entry of the A22 event or because of a decision not to return the patient from ‘leave of absence’ after all.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 10 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2 Simplified Workflow

2.1.2.1 ADT/A01 – Admit/visit notification

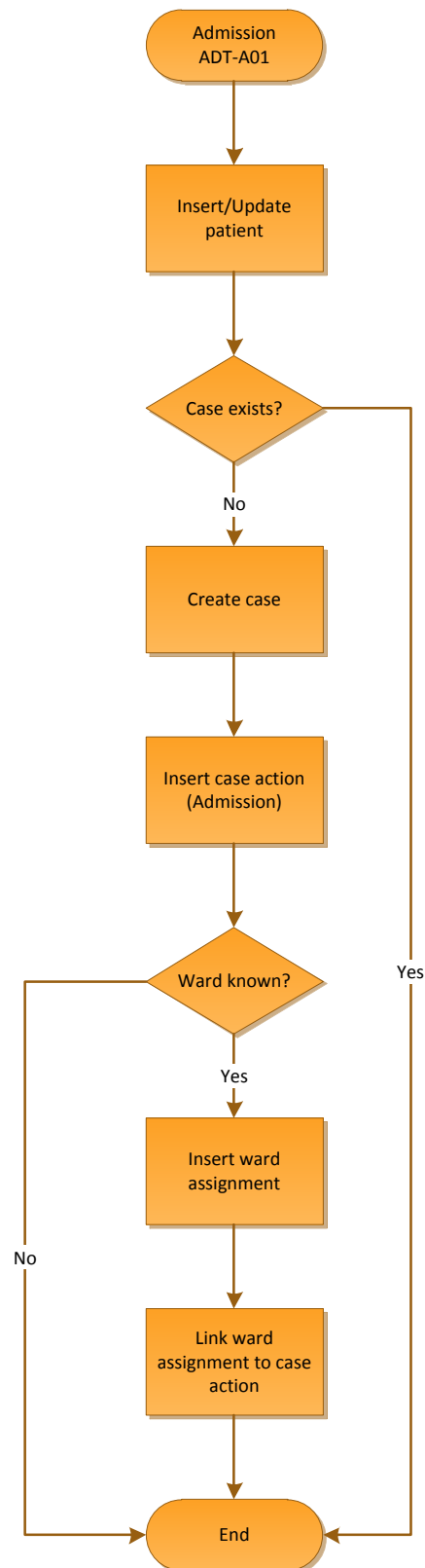

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 11 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.2 ADT/A02 – Transfer a patient

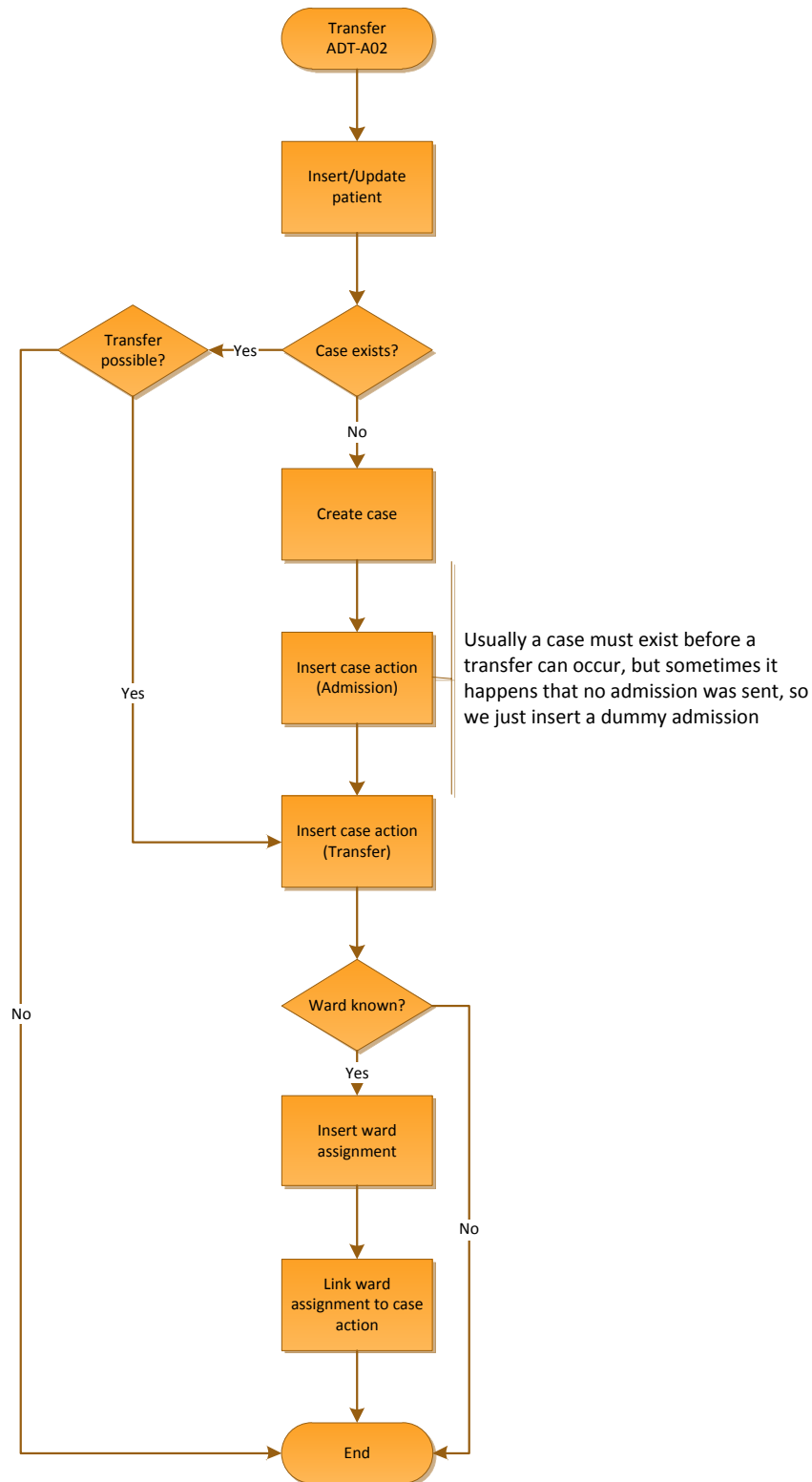

A transfer is possible if the last case action was an admission, a transfer or a return from a leave of absence

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 12 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.3 ADT/A03 – Discharge/end visit

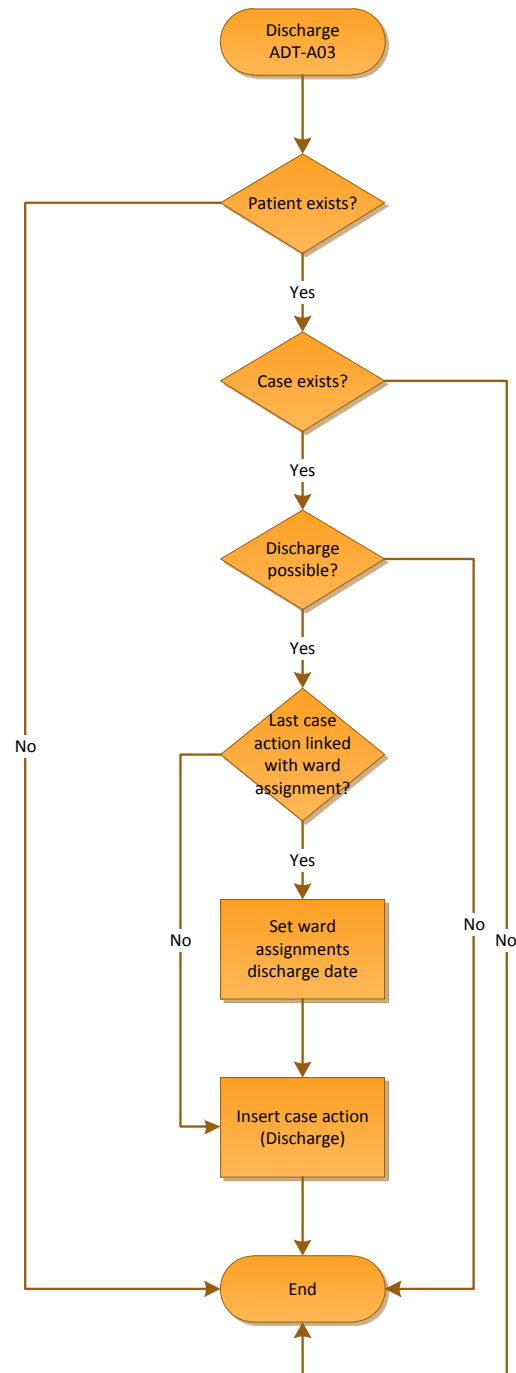

A discharge is possible if the last case action was an admission, a transfer or a return from a leave of absence

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 13 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.4 ADT/A04 – Register a patient

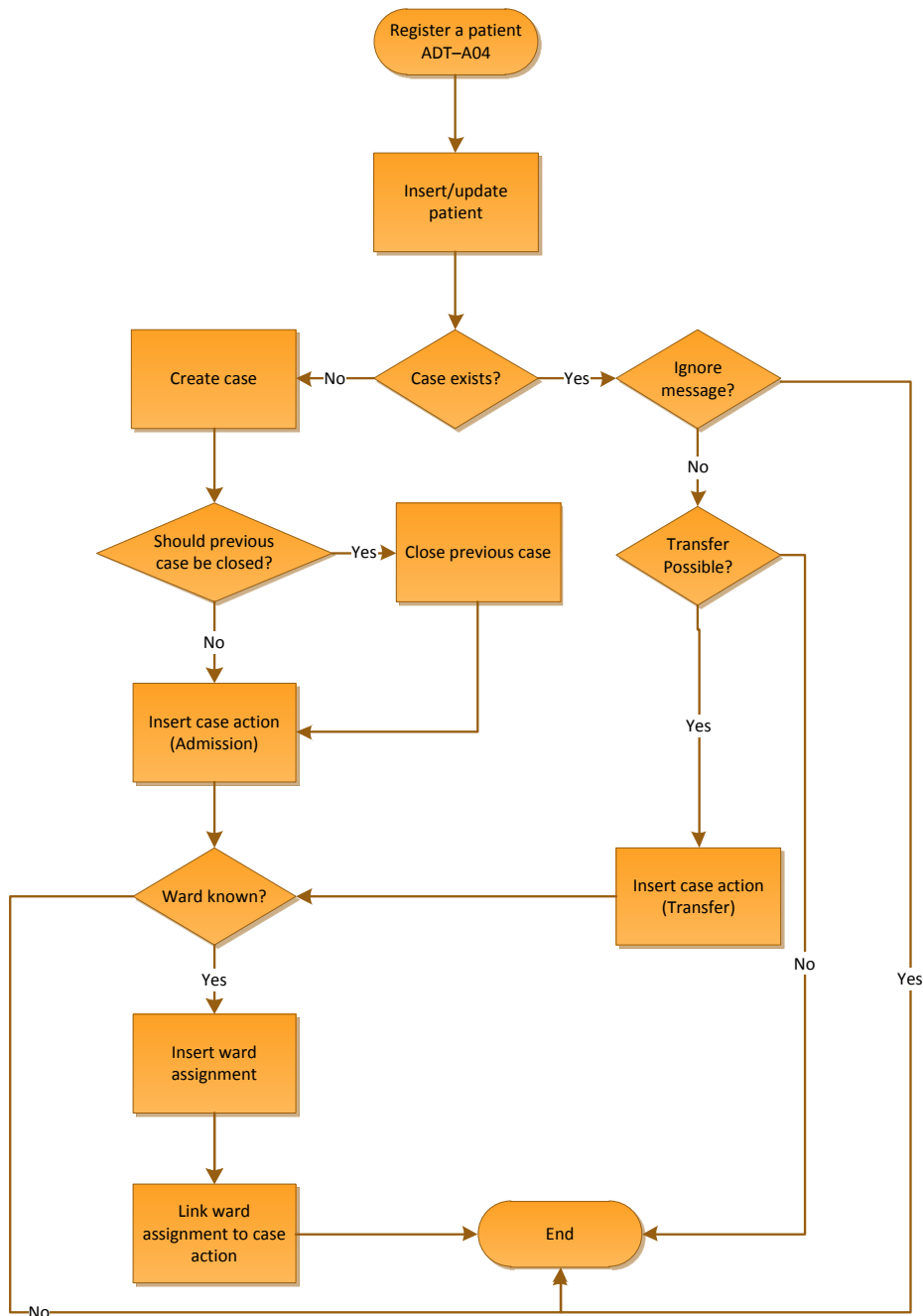


A transfer is possible if the last case action was an admission, a transfer or a return from a leave of absence.
The message can be ignored by setting the option ,A04-insert settings' to ,Ignore if a case exists'

A previous case will only be closed if both cases (the new and the previous) are outpatient cases, the previous case isn't already closed and the ward Ids are equal. The message can be activated by changing the option ,Close outpatient case on new admit' to true.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 14 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.5 ADT/A05 – Pre-admit a patient

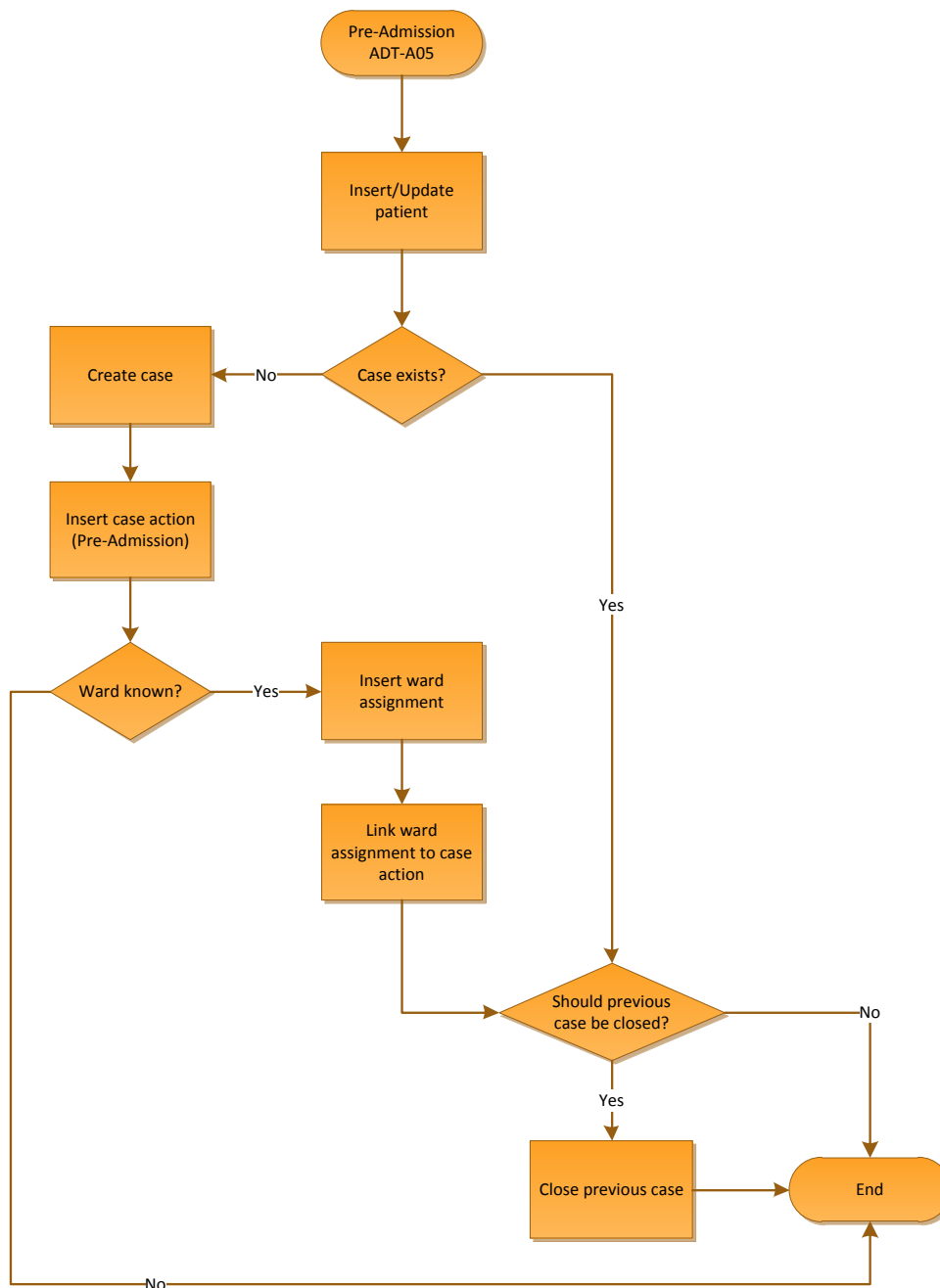

A previous case will only be closed if both cases (the new and the previous) are pre-admitted cases, the previous case isn't already closed and the ward ids are equal. The message can be activated by changing the option ,Close outpatient case during pre-admit' to true.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 15 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.1.2.6 ADT/A06 – Change an outpatient to an inpatient

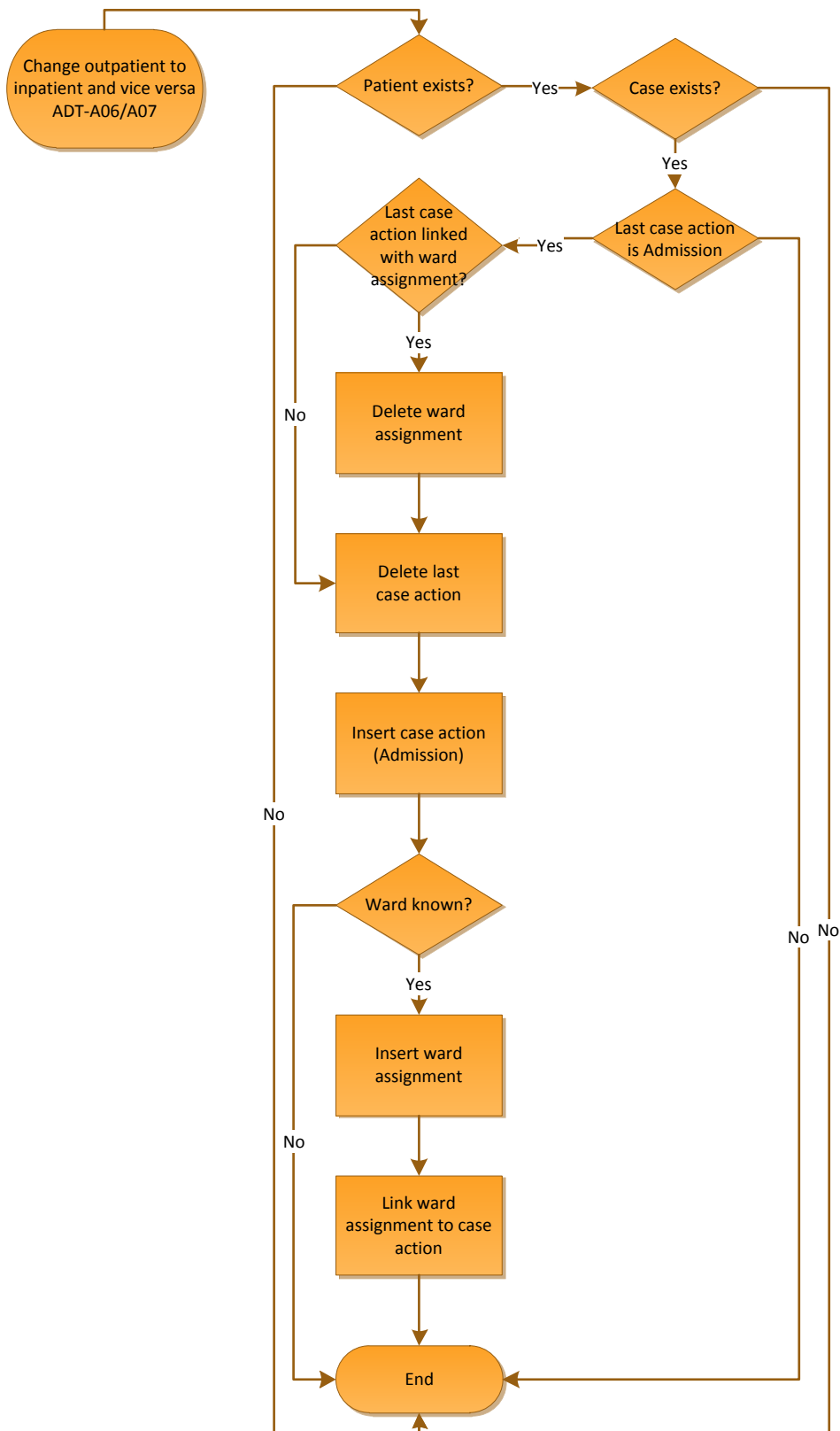


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 16 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.7 ADT/A07 – Change an inpatient to an outpatient


Workflow: See ADT-A06

2.1.2.8 ADT/A08 – Update patient information

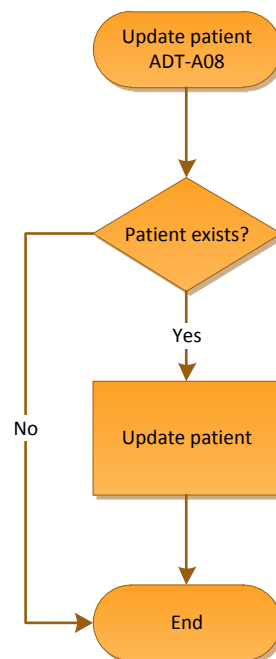


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 17 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.9 ADT/A11 – Cancel admit / visit notification

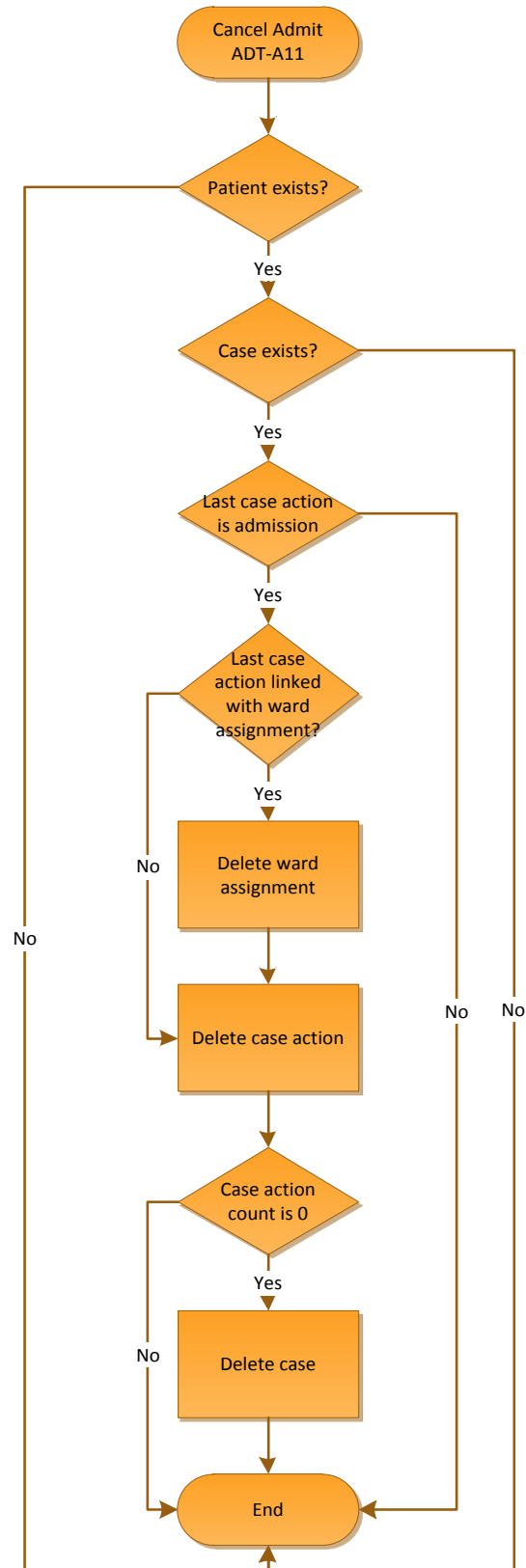


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 18 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.10 ADT/A12 – Cancel transfer

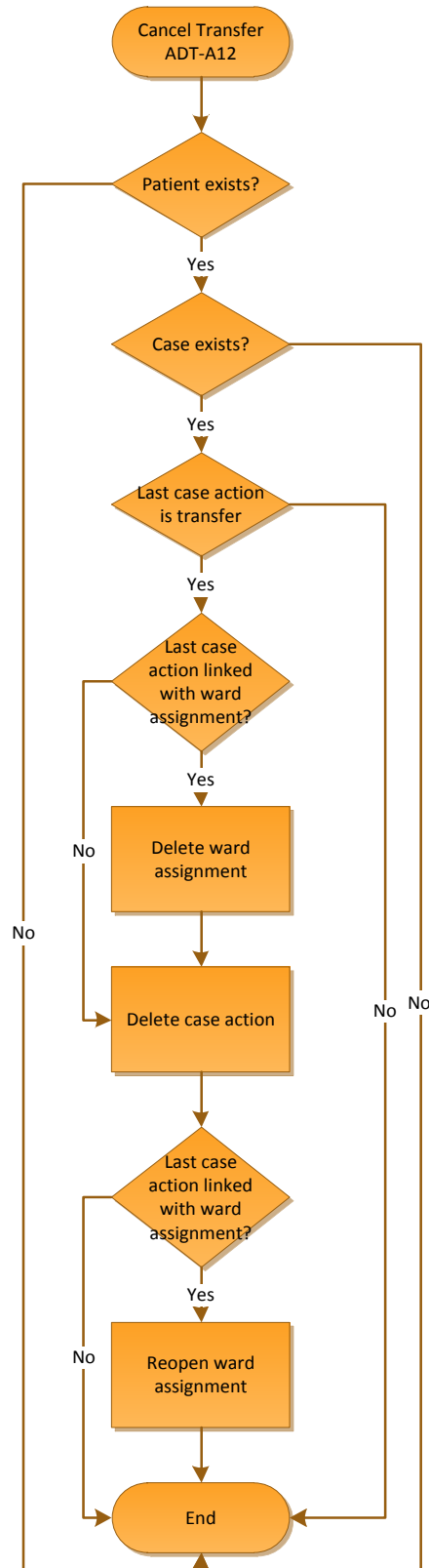


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 19 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.11 ADT/A13 – Cancel discharge / end visit

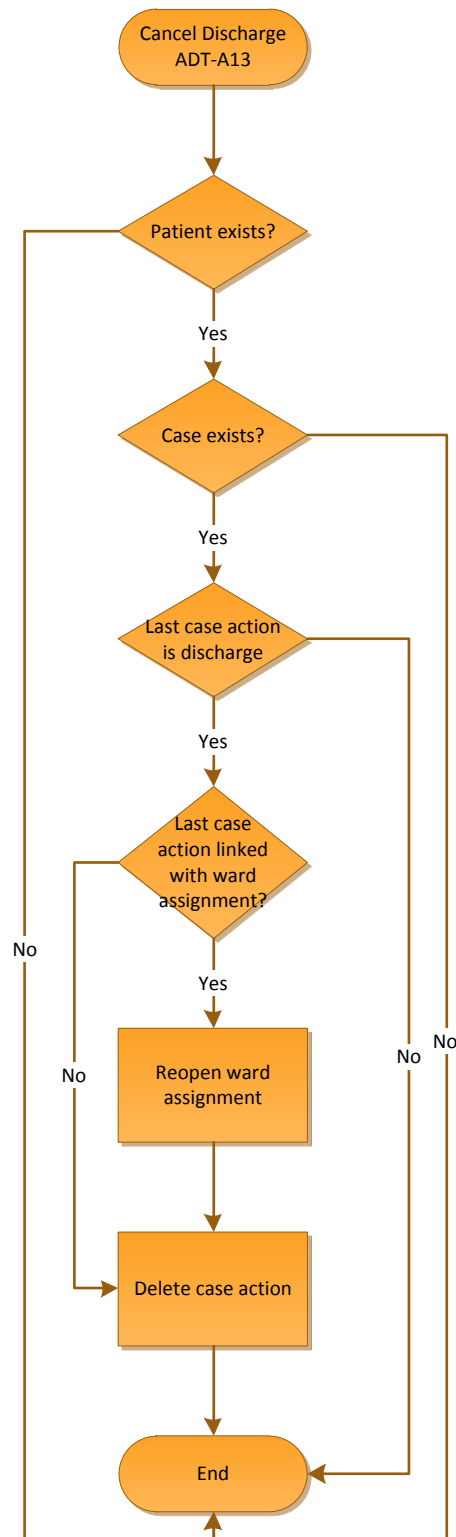

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 20 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.12 ADT/A21 – Patient goes on a leave of absence

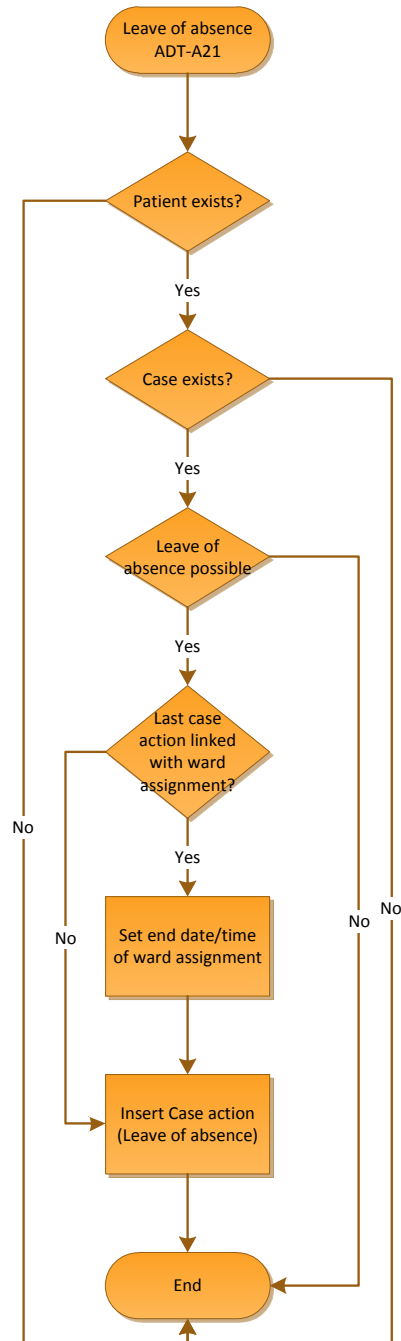

A leave of absence is possible if the last case action was an admission, a transfer or a return from a leave of absence.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 21 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\	User Manual Cato HL7 2.32
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.13 ADT/A22 – Patient returns from a leave of absence

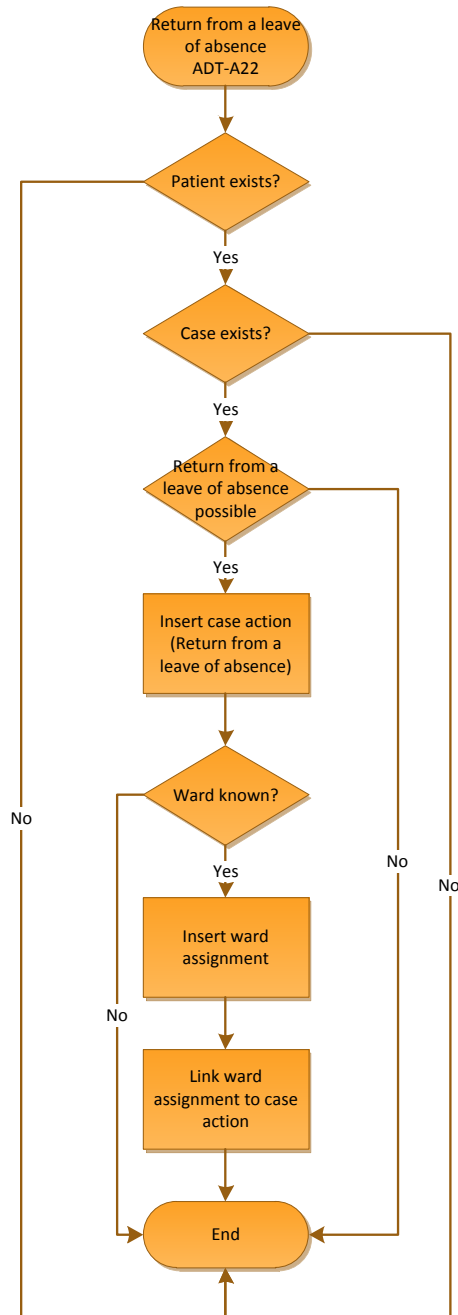


A return from a leave of absence is possible if the last case action was a leave of absence. The new case action and the ward assignment will be inserted with the ward of last presence (admission, transfer, return from a leave of absence) before the leave of absence. The ward can be unknown, too!

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 22 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.14 ADT/A28 – Add person or patient information

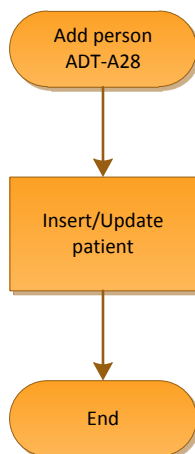


2.1.2.15 ADT/A29 – Delete person information

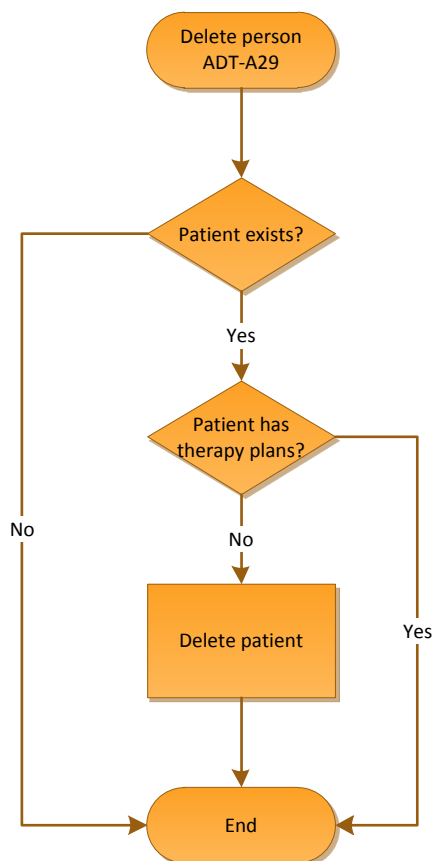


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 23 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.16 ADT/A31 – Update person information

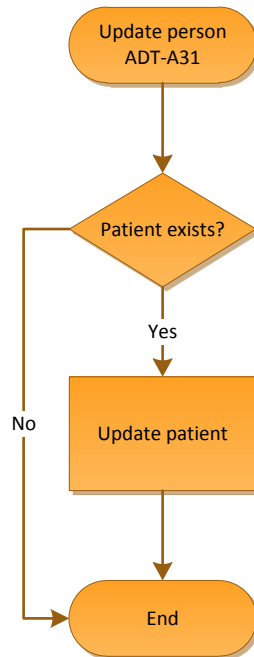

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 24 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.17 ADT/A34 – Merge patient information

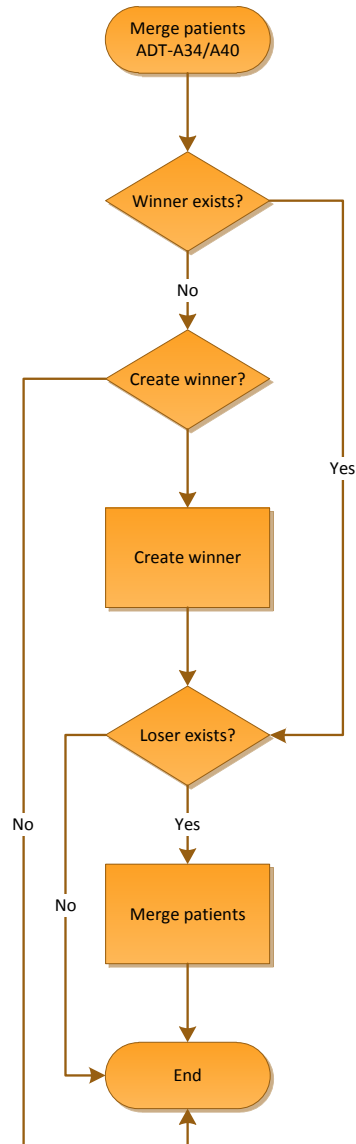

Creating the winner depends on the setting ,Create patients' which can be set to ,Always' , ,If the merge-slave exists' , ,Never'

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 25 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.18 ADT/A38 – Cancel Pre-Admit

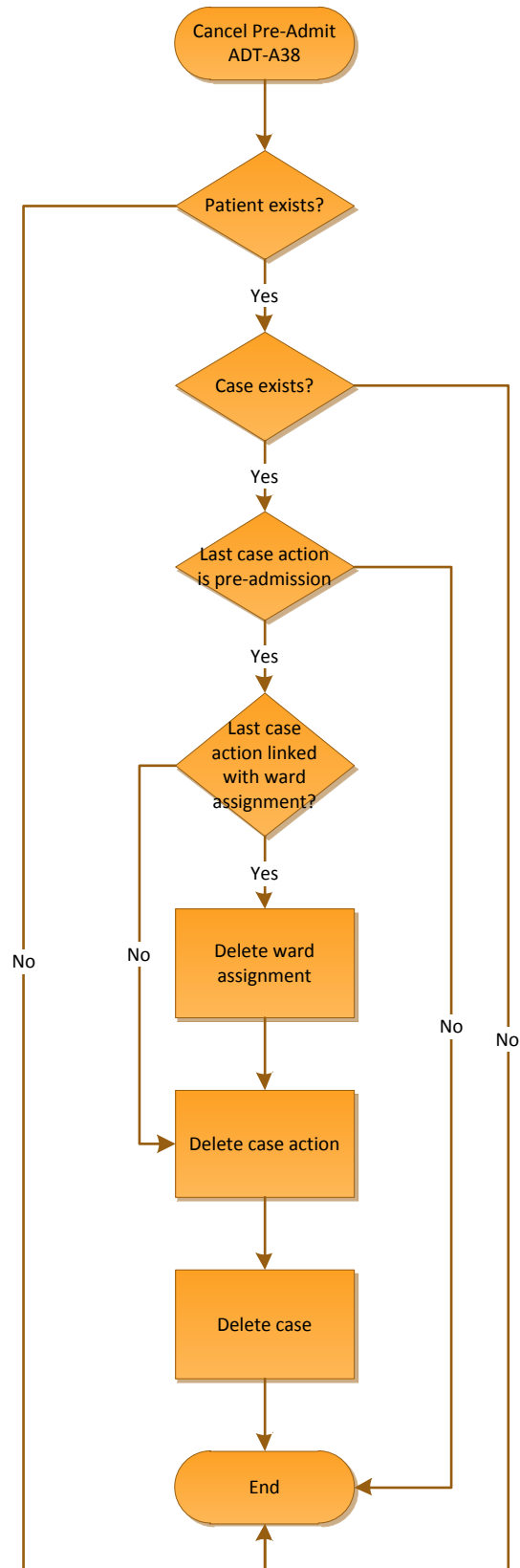


	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 26 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.1.2.19 ADT/A40 – Merge patient information - PID

Workflow: See ADT-A40

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 27 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.1.2.20 ADT/A45 – Move visit information

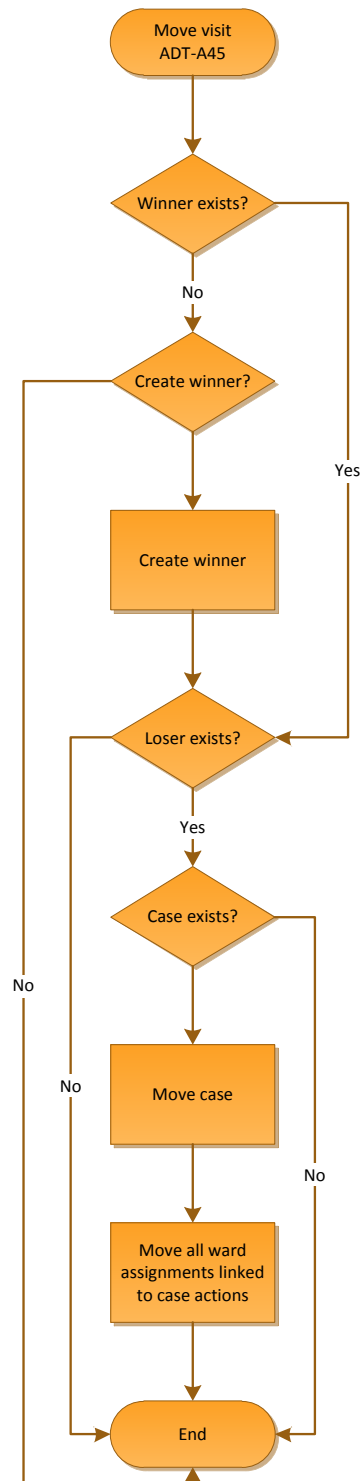

Creating the winner depends on the setting ,Create patients‘ which can be set to ,Always‘ , ,If the merge-slave exists‘ , ,Never‘


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 28 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.21 ADT-A52 – Cancel leave of absence

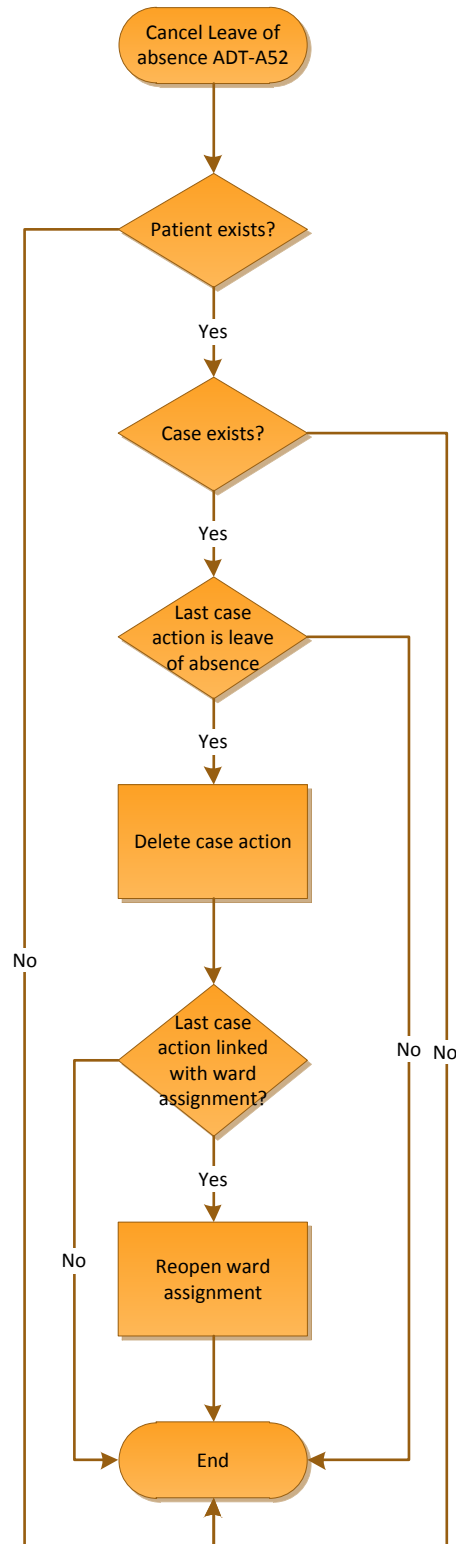

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 29 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.1.2.22 ADT-A53 – Cancel return from a leave of absence

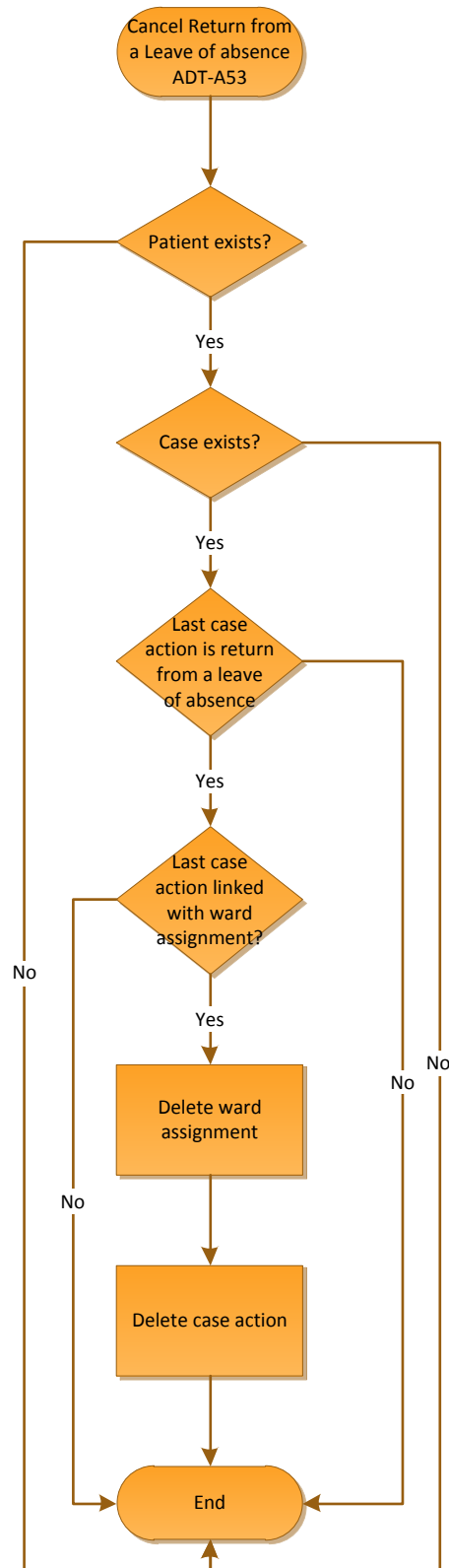

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 30 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.2 ORU messages


2.2.1 Short description

2.2.1.1 ORU/R01 – Unsolicited transmission of an observation message

The ORU message is for transmitting laboratory results to other systems

2.2.2 Workflow

2.2.2.1 ORU/R01 – Unsolicited transmission of an observation message

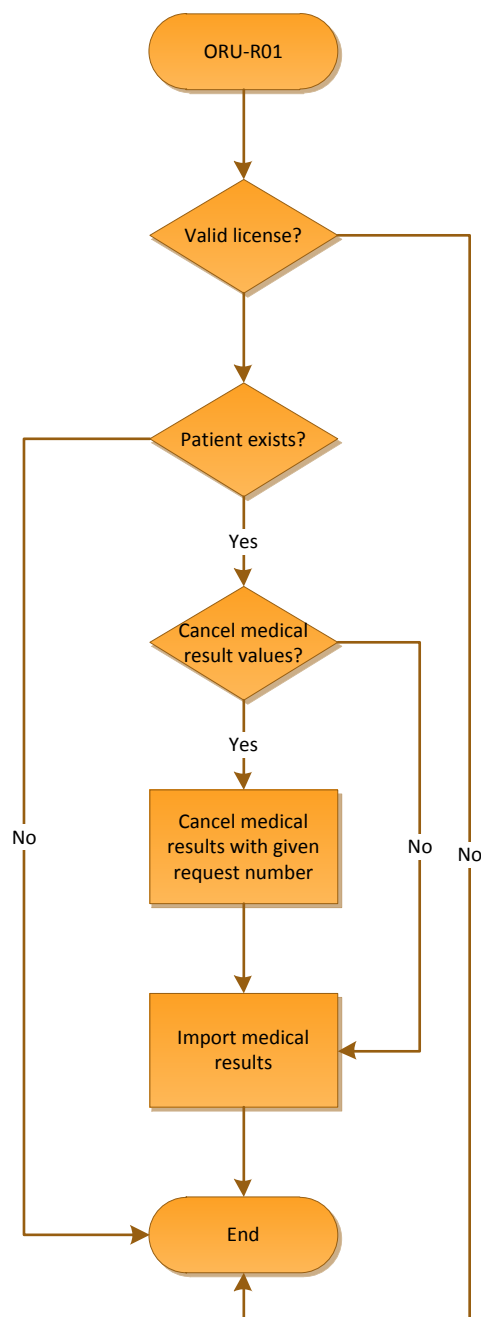

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 31 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.3 BAR Messages

2.3.1 Description

2.3.1.1 BAR^P01

This message is used in caton[®] for inserting or updating of diagnoses only.

2.3.1.2 BAR^P02

This message is used in caton[®] for cancelling of diagnoses only.

2.3.1.3 BAR^P05

This message is used in caton[®] for inserting or updating of diagnoses only.


2.3.1.4 BAR^P12

This message is used in caton[®] for inserting or updating of diagnoses only.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 32 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


2.3.2 Simplified Workflow

2.3.2.1 BAR^P01, BAR^P05, BAR^P12

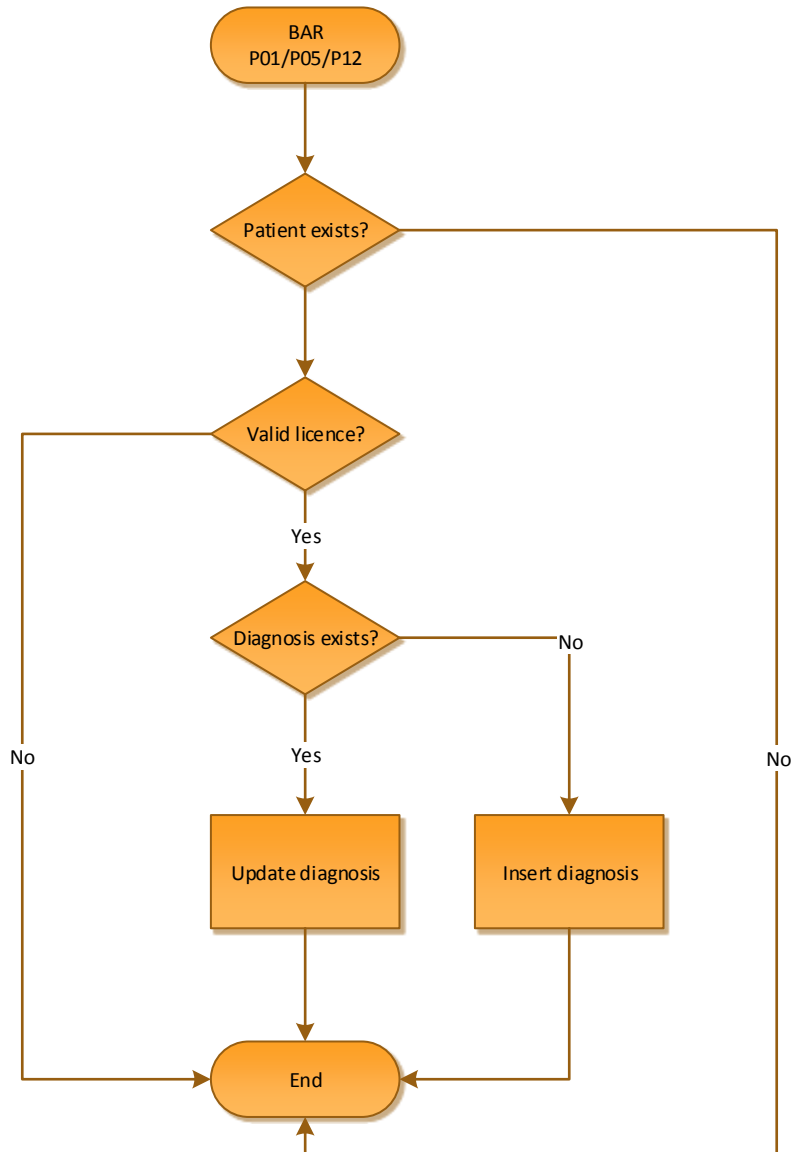

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 33 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.3.2.2 BAR^P02

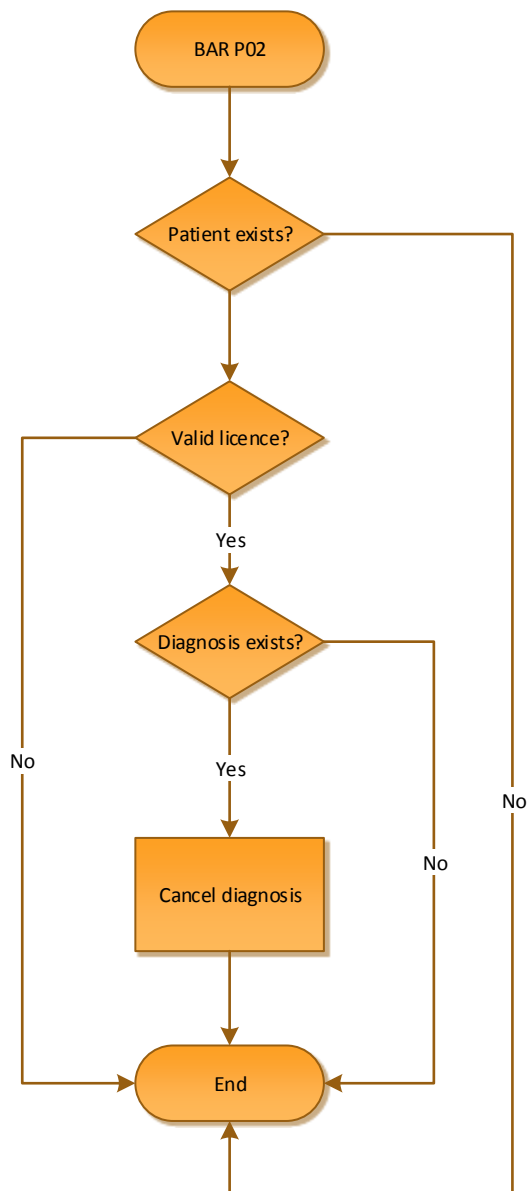

2.3.3 Supported Segments

The segments described in this chapter are interpreted by the cato® HL7 Interface. For the standard segments (MSH, EVN, ERR, MSA, and so on) not mentioned here please refer to a HL7 documentation.

Fields not listed here are not used by cato® and therefore not read.

2.4 PID – Patient Master Data

The PID segment is used as the primary means of communicating patient identification information. This segment contains permanent patient identifying and demographic information that, for the most part, is not likely to change frequently.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 34 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

- **Seq 2 – Patient ID**
Depending on [4.4.2.1.1.4 Import Configuration] this field is synchronized with the patient number or the patient UID registered in cato®.
NOTE: Used for backward compatibility only!
- **Seq 3– Patient Identifier List**
One single patient identifier or a list of unique patient identifiers. Depending on [4.4.2.1.1.4 Import Configuration] this field is synchronized with the patient number or the patient UID registered in cato®.
- **Seq 4 – Alternate Patient ID**
- **Seq 5 – Patient Name**
- **Seq 6 – Mother's Maiden Name**
- **Seq 7 – Date of Birth**
- **Seq 8 – Sex**
- **Seq 9 – Patient Alias**
- **Seq 11 – Patient Address**
- **Seq 13 – Phone Nr. 1**
- **Seq 14 – Phone Nr. 2**
Used only if Phone Nr. 1 is left blank
- **Seq 19 – SSNR (Social Security Number)**
- **Seq 21 – Mother's Identifier**
- **Seq 29 – Death Date**
- **Seq 30 – Death Indicator**

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 35 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.5 PV1 – Patient Visits (Ward Assignments)

The PV1 segment is used by Registration/Patient Administration applications to communicate information on a visit-specific basis. This segment can be used to send multiple-visit statistic records to the same patient account or single-visit records to more than one account.

- **Seq 2 – Patient Class**
Defines a patient as in- or outpatient. In cato® the patient class is defined by the ward which can be marked as inpatient or outpatient.
- **Seq 3 – Assigned Patient Location**
The current ward of the patient. Mainly used for synchronization of wards (See 4.4.2.1.3.1 Ward's ID). Please notice that not only the UID has to match but also the type of the ward (in-patient ward or outpatient ward) depending on Seq 2 – Patient Class.
- **Seq 6 – Prior Patient Location**
The ward the patient is moved from. Optional
- **Seq 7 – Attending doctor**
This field contains the attending physician information. Optional
- **Seq 8 – Referring doctor**
This field contains the referring physician information. Optional
- **Seq 11 – Temporary Location**

Seq 19 – Visit Number

This field can be mapped to either the **admission number** of a cato® ward assignment or to the **case number**. (See 0

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 36 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

- Database field for saving the stay number)
- **Seq 20 – Financial Class**
Insurance status, optional
- **Seq 24 – Contract Code**
Type of case, optional
- **Seq 39 – Servicing Facility**
Defines the healthcare facility with which the visit is associated.
- **Seq 44 – Admit Date Time**
Starting date time of the visit.
- **Seq 45 – Discharge Date Time**
Ending date and time of the visit. Optional.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 37 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

2.6 ZBE – Patient Moving Information

For a correct synchronization of the full history of the patient moves within the hospital some systems use the ZBE segment. With this segment each move of a patient has its own entity with a unique ID which allows a full synchronization with the ward assignments in caton®.

- **Seq 1 – Unique Move ID**
The ID representing the specific transfer of the patient
- **Seq 2 – Start Date Time**
Starting date time of the transfer / visit.
- **Seq 3 – End Date Time**
Ending date time of the transfer / visit.
- **Seq 4 – Action**
Possible values are:
 - INSERT
 - UPDATE
 - DELETE

2.7 MRG – Patient Merge Information

The MRG segment provides caton® with information necessary to initiate the merging of patient data. In caton® it is possible to merge two patients (ADT-A34, A40) or to shift cases/visits from one patient to another (ADT-A45)

- **Seq 1 – Prior Patient Identifier List**
The UID of the patient losing the merge process. The patient specified here will no longer exist in the caton® database, all its visits, medical results, diagnoses, therapies and so on are transferred (merged) with the winning patient (as specified in the PID segment)
- **Seq 5 – Prior Visit ID**
Used for the A45 Message to indicate the old (wrong) visit ID.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 38 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.8 DG1 – Diagnose Information

The DG1 Segment contains information about the patient's diagnoses (e.g. coded in ICD10).

- **Seq 2 – Coding Method**

This field is read for backward compatibility only! Seq 3 (Diagnose Code) should be used instead.

Seq 3 – Diagnose Code

Contains information about the diagnose code, and (for newer versions) the coding system and a description. Valid coding systems for cato® are ICD9 and ICD10, but it is also possible to specify a diagnose code without a coding system. (See 0

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 39 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

- Diagnosis type)
- **Seq 4 – Diagnose Description**
This field is read for backward compatibility only! Seq 3 (Diagnose Code) should be used instead.
- **Seq 5 – Diagnose Date Time**
- **Seq 6 – Diagnose Type**
For future use.
- **Seq 7 – Diagnose Priority**
For future use.

2.9 OBX – Observation / Medical results

The OBX segment is used to transfer a single patient's medical result.

- **Seq 2 – Value Type**
Type of the result value. (See 4.4.2.1.5.4 Valid observation value types (OBX-2))
- **Seq 3 – Observation Identifier**
The medical result UID. This field is used to synchronize with the possible medical results of the cato® database.
- **Seq 5 – Observation value**
Value of the medical result.
- **Seq 6 – Observation unit**
Used for comparing the configured sending system's unit with the actually sent unit. If these units differ the medical result value will NOT be imported!
- **Seq 7 – References range**
Defines the standard range for a particular medical result. (See 4.4.2.1.5.8 Import standard range)
- **Seq 11 – Result Status**
This field defines the status of the result. (See 4.4.2.1.5.5 Medical result status (preliminary) (OBX-11) and 4.4.2.1.5.6 Medical result status (invalid) (OBX-11))
- **Seq 14 – Date Time of the Observation**

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 40 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

2.10 OBR – Observation request segment

The OBR segment acts (in some message types) as header for a list of OBX segments.

- **Seq 2 – Placer Order Number**
Can be used as request number for medical results
- **Seq 3 – Filler Order Number**
Can be used as request number for medical results
- **Seq 7 – Observation Start date/time**
First fallback if no date/time field stated in OBX-segment
- **Seq 8 – Observation End date/time**
Second fallback if no date/time field stated in OBX-segment
- **Seq 14 – Receiving date/time**
Third fallback if no date/time field stated in OBX-segment

2.11 IN1 – Insurance Segment


- **Seq 2 – Insurance Plan ID**
Unique identifier for the insurance plan.
- **Seq 3 – Insurance Company ID**
Unique identifier for the insurance company.
- **Seq 4 – Insurance Company Name**
Insurance company name
- **Seq 5 – Insurance Company Address**
Address of the insurance company
- **Seq 12 – Plan Effective Date**
Can be used for cost unit validity
- **Seq 13 – Plan Expiration Date**
Insurance Card Validity; can also be used for cost unit validity
- **Seq 29 – Verification Date/Time**
Can be used for cost unit validity
- **Seq 36 – Policy Number**
Can be used as patient's social security number

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 41 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

3 Preparing the cato® database for HL7

3.1 Synchronization

The synchronization can be changed as seen below: (See 4.4.2.1.1.4 Import Configuration)

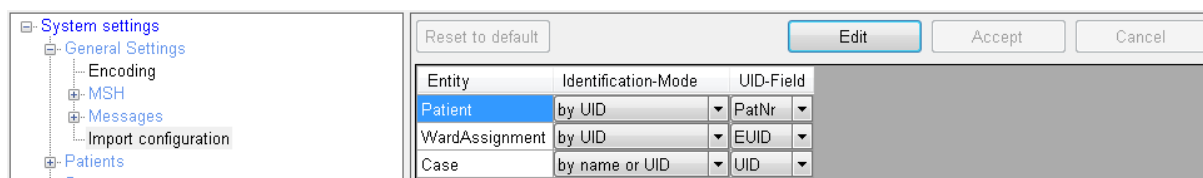

Entity	Identification-Mode	UID-Field
Patient	by UID	PatNr
WardAssignment	by UID	EUID
Case	by name or UID	UID

3.1.1 Description

3.1.1.1 Identification-Mode

Defines how an entity should be identified.

Possible modes:

- **By name**
A patient will be identified by name and birth date.
A ward assignment will be identified by ward and visit start date/time.
- **By UID**
The entity will be identified by the field specified by 'UID-Field'
- **By name or UID**
A combination of UID and name; at first the entity will be identified by UID, if the entity was not found then it will be identified by name
- **No selection**
The entity should not be identified.

NOTE: This mode is not recommended for most use cases.

3.1.1.2 Synchronization-Mode

Not used

3.1.1.3 UID-Field

The entity will be identified by UID-Field if the 'Identification-Mode' is set to 'By UID' or 'By name or UID'. Possible fields depend on the entity, so they will be described below.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 42 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

3.1.2 Patient synchronization

3.1.2.1 Identification-Mode

- **By name**
Patient is identified by first name, last name and birth date
- **By UID (Recommended)**
Patient is identified by a unique identifier defined in [3.1.2.2 UID-Fields]
- **By name or UID**
As defined in [3.1.1.1 Identification-Mode]
- **No selection**
As defined in [3.1.1.1 Identification-Mode]

3.1.2.2 UID-Fields

- **PatNr (Recommended)**
The patient's number.

3.1.2.3 Starting the interface connecting to an existing database

When starting the HL7 Interface in an environment where cato® was already used without the HL7 interface, the attention should be drawn on the patient synchronization, **to avoid the doubling of patients!**

There are two possibilities:

1. If patient numbers have been correctly **manually entered** by the users, there should be no problem starting the interface.
2. If no patient numbers exist, the following configuration modes are possible:
 - a. Update all patient numbers manually (the secure method) or
 - b. **Allow the synchronization by Name and Birth Date:**

This means: If the patient to be inserted or updated, could not be found by patient number (because it was not entered in cato®), the HL7 Interface tries to find the patient by Name (**first name and last name**) and the **birth date**.

Problems with this method of synchronization can occur, if multiple patients have the same firstname, lastname and birth date. Different spelling of names is also a problem. E.g. The patient last name is stored in the hospital information system as 'MacAllister', but in cato® as 'McAllister'.

If the decision is to use this faulty method of synchronization, please choose 'By Name or UID' as synchronization mode for the patient.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 43 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

3.1.3 Ward assignment synchronization

3.1.3.1 Identification-Mode

- **By name (Recommended)**
Ward assignment is identified by patient and visit start date.
- **By UID**
Ward assignment is identified by a unique identifier defined in [3.1.3.2 UID-Fields]
By name or UID
As defined in [3.1.1.1 Identification-Mode]
- **No selection**
As defined in [3.1.1.1 Identification-Mode]

3.1.3.2 UID-Fields

- **EUID**
Identifier as defined in ZBE segment
- **CaseNumber**
Value stored in case number
- **AdmissionNumber**
Value stored in admission number

3.1.4 Case synchronization

3.1.4.1 Identification-Mode

- **By name (Recommended)**
Case is identified by patient and case number.
To avoid problems with multiple case numbers, this mode is the preferred option for case identification
- **By UID**
Case is identified by a unique identifier defined in [3.1.4.2 UID-Fields]
Should only be used if case numbers are unique for all patients.
- **By name or UID**
As defined in [3.1.1.1 Identification-Mode]
- **No selection**
As defined in [3.1.1.1 Identification-Mode]

3.1.4.2 UID-Fields

- **UID**
Same as case number


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 44 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

3.1.5 Ward synchronization

Before a patient's visit can be successfully transferred, the unique Ids of the wards need to be stored in cato®. This can be done with the cato® client:

Ward name	Chief physician	Cost unit	Ward no.	Preparation site no.	Type	LTPN	Phone
1A - Testward	Dr. Queen		1A@13C		Ward		
2A - Testward	Dr. King		2A@34X		Outpatient ward		

The value of the ward no. depends on the configuration specified in the HL7 configuration:

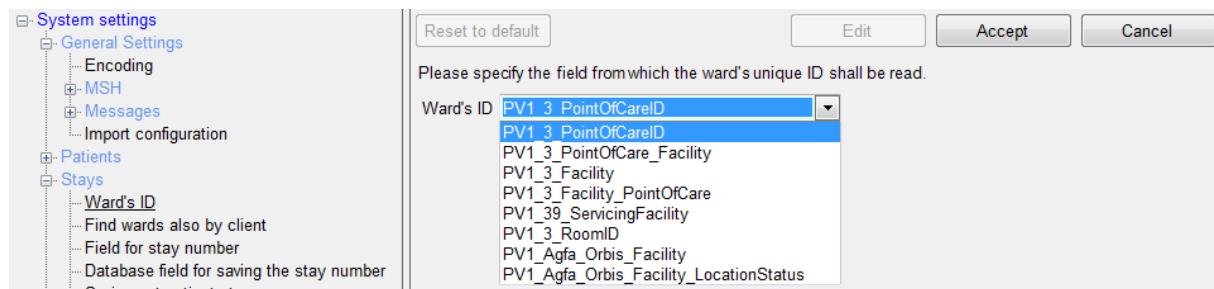

For a description please refer to 4.4.2.1.3.1 Ward's ID

3.1.6 Medical results

See 4.4.3.2.3.2 Medical result values

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 45 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4 Configuration

4.1 Licenses

NOTE: Before starting the service please make sure to have the appropriate licenses for HL7 added to the cato® database! Inserting a license can be done by using the license management module of the cato® client.

Three different types of licenses exist:

- The HL7 basic license for receiving patient data messages
- The HL7 medical result license for receiving patient's medical results (additional)
- The HL7 diagnoses license for receiving patient's diagnoses (e.g. ICD10) (additional)


4.2 The Configuration and Log-Viewer Window

For all settings of the interface, a front end is installed together with the service. After starting the cato® HL7 Configuration tool a new icon is shown in the taskbar:




A red icon indicates that the service is stopped; a green icon stands for a running service.

By right-clicking on this icon, the service can be started and stopped and the configuration tool can be accessed:

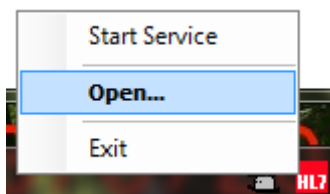


Click on “Open...” (or double-click the icon) to open the configuration dialog:

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 46 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

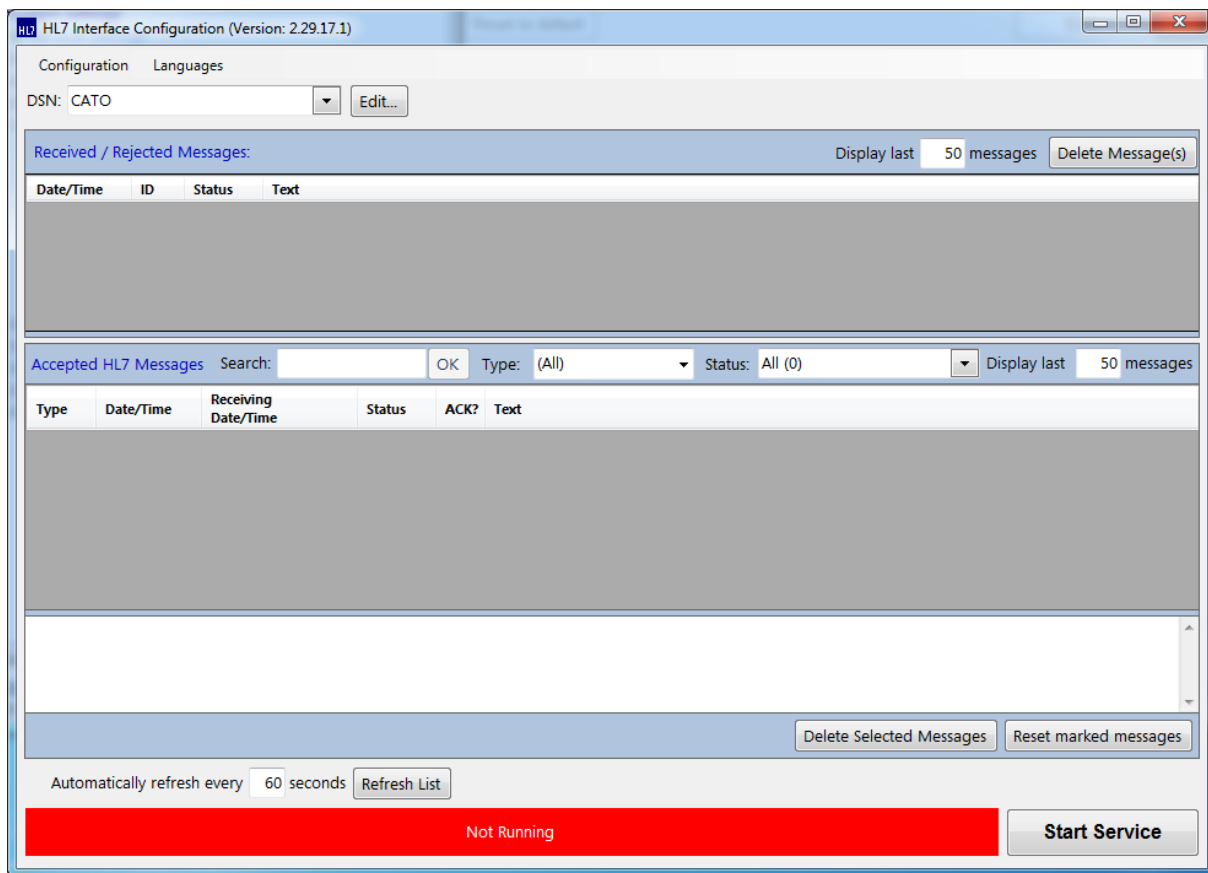


To select the preferred language click on menu item 'Languages':

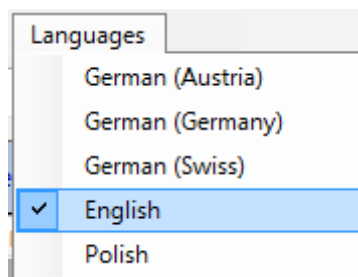


4.3 Changing the Configuration

To set up the configuration select 'Edit Configuration' from menu 'Configuration':

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 47 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

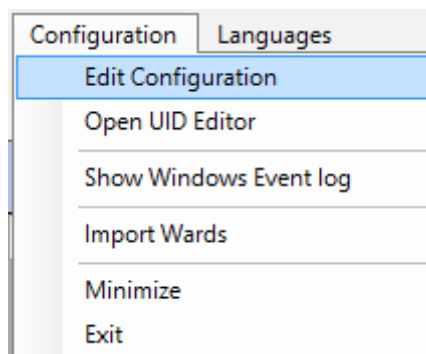

The HL7 Settings Window will be shown:

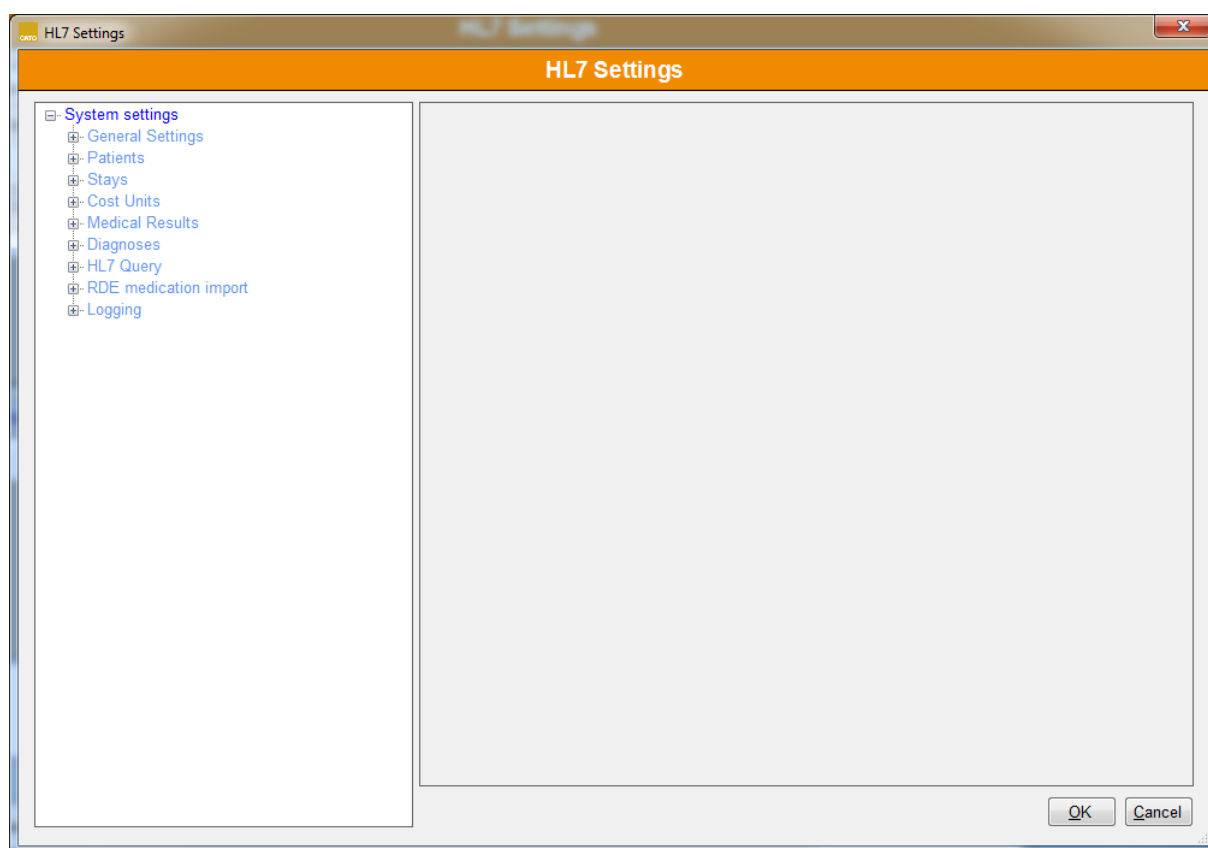

The following chapters describe the configuration options of the cato® HL7 Interface.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 48 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4 HL7 Settings

4.4.1 Fundamentals

HL7 settings are hierarchically structured and they consist of levels, sections and settings.

Example:



	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 49 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

NOTE: Levels are displayed in dark blue, sections in light blue. Settings are displayed in black (or orange, if they are inherited).

4.4.1.1 Levels

There are four different levels:

- **System settings** (Occurrence: Once)
 - **Database settings** (Occurrence: Multiple)
 - **Client settings** (Occurrence: Multiple)
 - **Connection settings** (Occurrence: Multiple)

Each level contains several sections.

NOTE: Levels appear in dark blue

4.4.1.2 Sections

A section contains logically related settings.

NOTE: Sections appear in light blue

4.4.1.3 Settings


A setting will be inherited from the parent-level unless the setting is overwritten.

NOTE: If you have multiple databases with multiple clients and connections the easiest way is to configure settings in the ‘System settings’ level and to overwrite specific options in underlying levels!


E.g.: The encoding for each connection in you configuration is ‘Western European’ and only one connection requires ‘Unicode (UTF-8)’ then set the encoding in the ‘System settings’ to ‘Western European’ and only in the specific connection you have to set it to ‘Unicode (UTF-8)’.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 50 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

When a setting is selected the value of the setting can be seen on the right side.

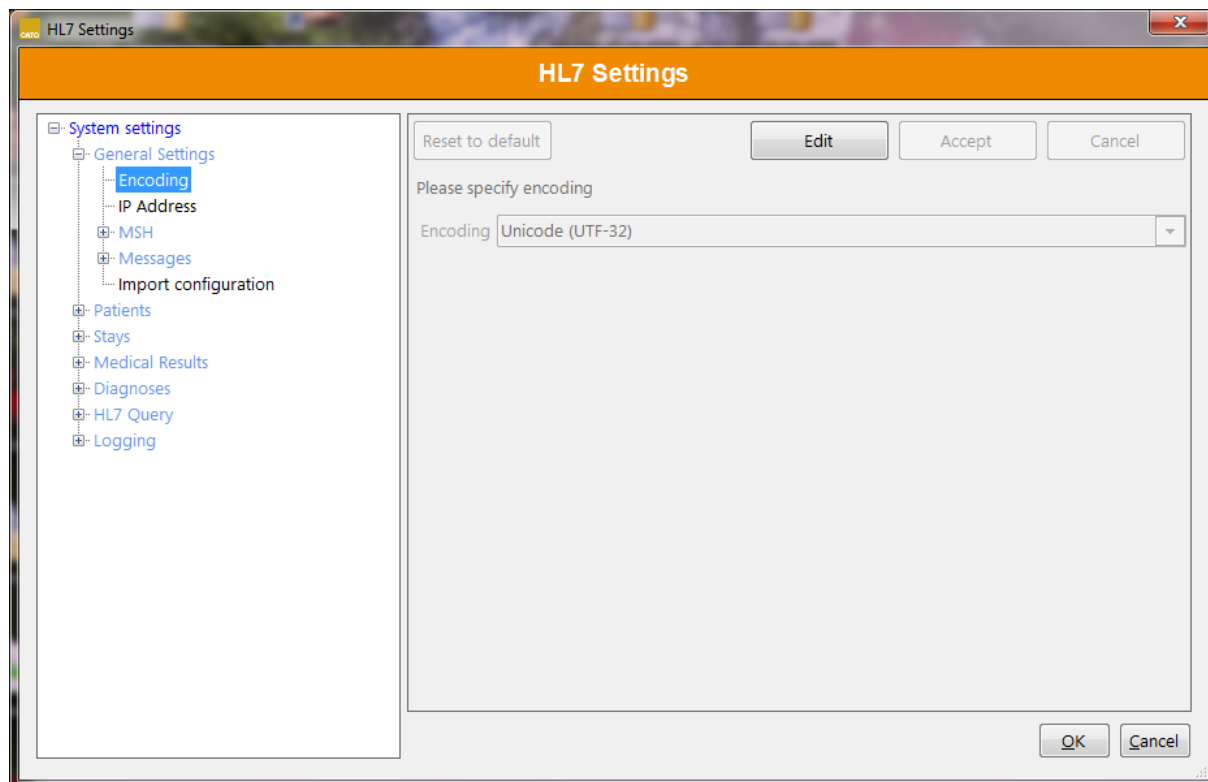

By clicking the button “Edit” the setting can be modified:

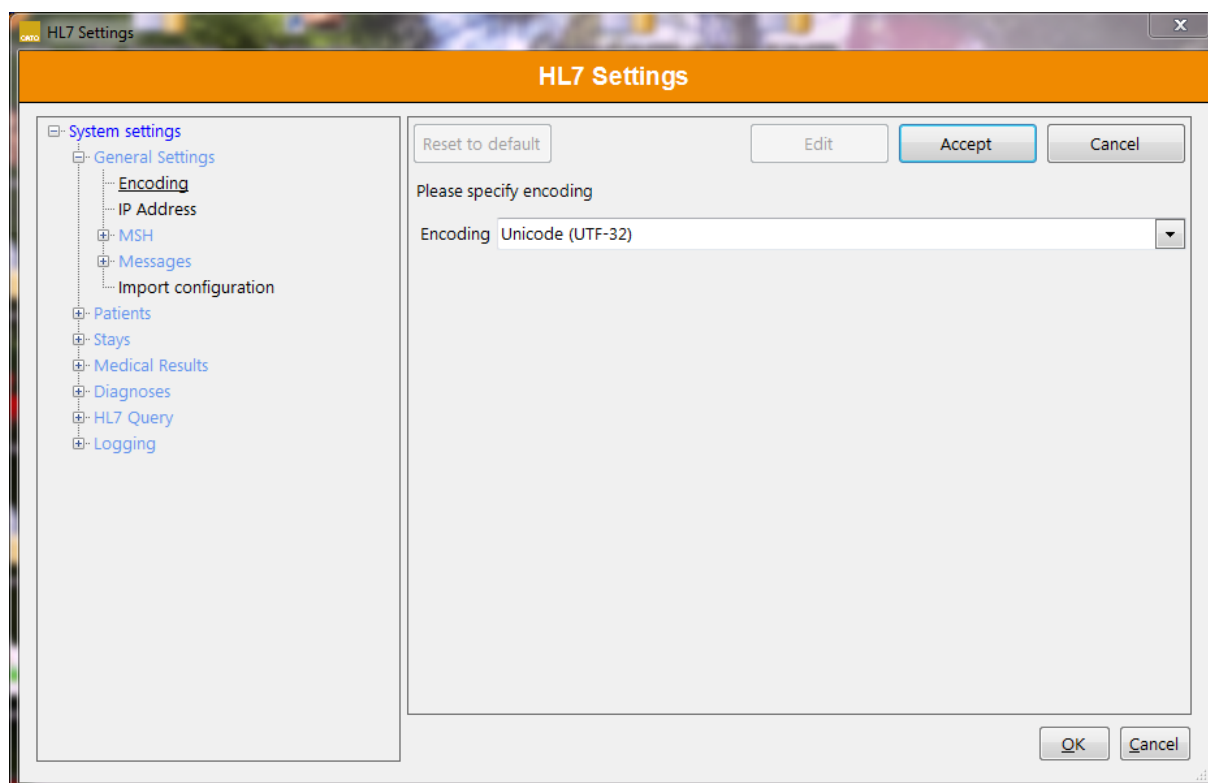


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 51 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

By clicking the button ‘Accept’ the setting will be changed.

By clicking the button ‘Cancel’ no changes will be applied.

NOTE: The button ‘Reset to default’ is only available when a setting is overwritten in a sub-level. By clicking this button the setting will be cleared and so inheritance will take effect again.

NOTE: Settings appear in black unless they are inherited; inherited settings appear in orange.

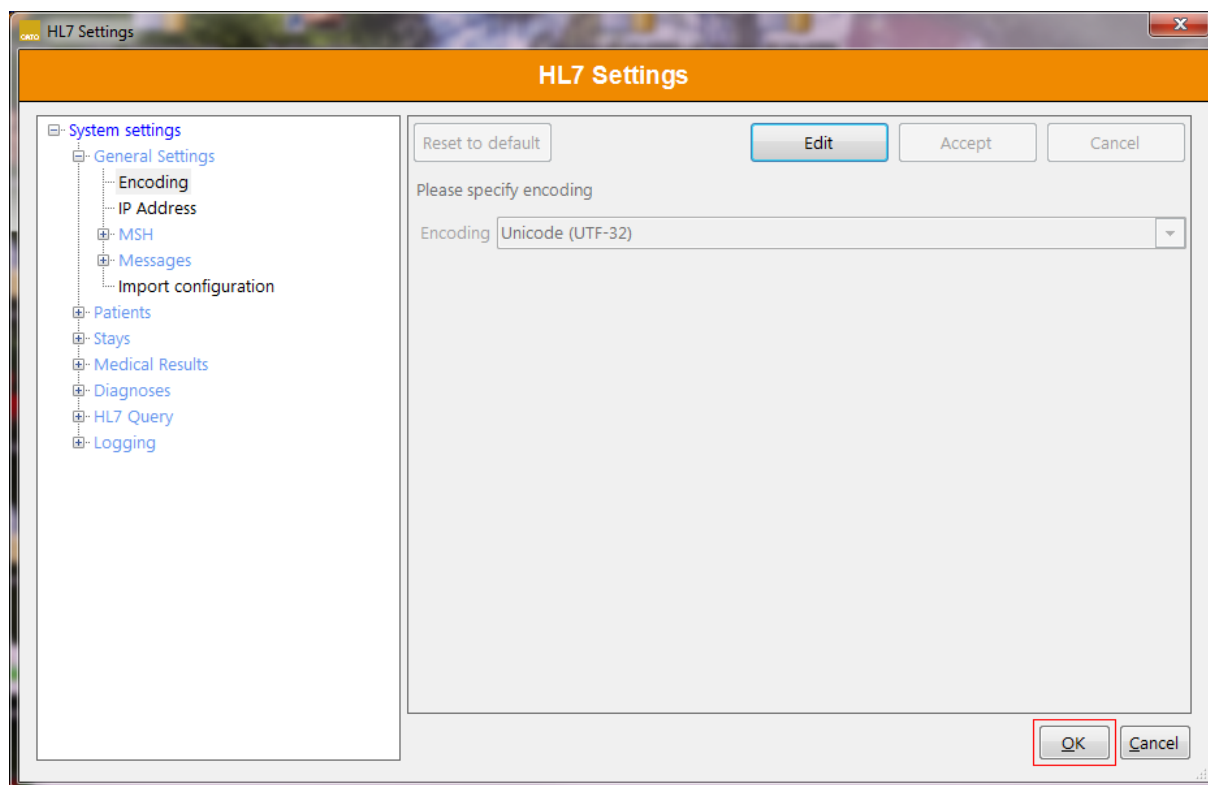

NOTE: All changes are not saved to disk unless the button OK in the bottom right part of the window is clicked!

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 52 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	


4.4.2 System Settings

4.4.2.1 Settings

4.4.2.1.1 General Settings (Section)

4.4.2.1.1.1 Encoding

Enter the encoding of the received messages

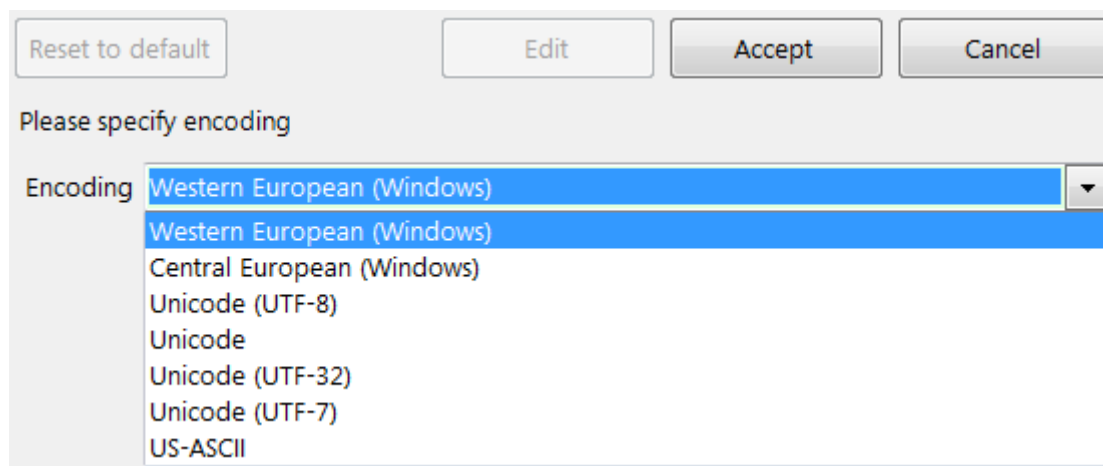

Reset to default Edit Accept Cancel

Please specify encoding

Encoding: Western European (Windows)

- Western European (Windows)
- Central European (Windows)
- Unicode (UTF-8)
- Unicode
- Unicode (UTF-32)
- Unicode (UTF-7)
- US-ASCII


Default: Western European (Windows)

NOTE: If problems with the received messages (such as an unreadable text or incorrect signs) occur this can signal a wrong encoding

4.4.2.1.1.2 MSH (Section)

4.4.2.1.1.2.1 HL7 Sending Application

Enter the Sending Application for HL7 Acknowledgements.

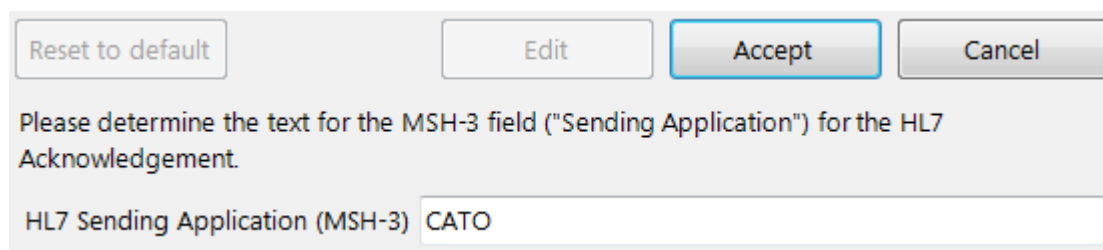

Reset to default Edit Accept Cancel

Please determine the text for the MSH-3 field ("Sending Application") for the HL7 Acknowledgement.

HL7 Sending Application (MSH-3) CATO

Default: CATO

4.4.2.1.1.2.2 HL7 Sending Facility

Enter the Sending Facility for HL7 Acknowledgements.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 53 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please determine the text for the MSH-4 field ("Sending Facility") for the HL7 Acknowledgement.

HL7 Sending Facility (MSH-4) CATO

Default: CATO

4.4.2.1.1.3 Messages (Section)

4.4.2.1.1.3.1 Message mapping

Specify the mapping of incoming messages:

- **No mapping**
Process message as it is
- **Neuchatel**
Special mapping for Neuchatel
- **ORU**
Maps ORU (without R01) to ORU^R01

Reset to default
Edit
Accept
Cancel

Please specify the message mapping

Message mapping No mapping ▼

- No mapping
- Neuchatel
- ORU

Default: No mapping

NOTE: Usually there is no need of a message mapping. If you need a special mapping please contact Cato Software Solutions.

4.4.2.1.1.3.2 General Settings

- **Create patients** (Default: yes/true)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Create as temporary messages** (Default: no/false)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
Only changeable if save cases is selected
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 54 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Please specify the general message settings.

☒ Create patients

☒ Update patients

☒ Import diagnoses

☒ Process ZBE updates

☐ Create as temporary messages

☐ Ignore consultation

☒ Save cases

☒ Save ward assignments

☐ Process message anew by default, when ZBE-UID non-existent

NOTE: Please be careful when deactivating one or more of the options! You avoid the stated action for ALL message types when deactivating a setting.

This means that NO patient will be created on any message received when removing the hook from ‘Create patients’ (and similar for the other settings) independent from specific message settings!

4.4.2.1.1.3.3 ADT-A01 (Section)

4.4.2.1.1.3.3.1 ADT-01

- **Create patients** (Default: If cases are inserted)
 - Always
 - **If cases are inserted**
Then a patient will only be created when the setting “Save cases” is set to yes/true
 - Never
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
Only selectable if save cases is chosen
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 55 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

- **Save visit start** (Default: Unchanged)
The begin of a visit can be changed if needed.
 - **Unchanged**
Process visit start as received
 - **Current timestamp**
Only used for some HL7 Query customers, sets the current date/time
 - **Begin of day**
Visit start time is set to 00:00 or the earliest possible time (depending of previous visits)
- **Save visit end** (Default: Unchanged)
The end of a visit can be changed if needed.
 - **Unchanged**
Process visit start as received
 - **Keep open**
End of visit will not be saved
 - **End of day**
Visit end is set to 00:00 of the next day (in relation to the received visit end)
 - **End of day, if time is not available**
Visit end is set to 00:00 of the next day (in relation to the received visit end), if no time is stated
- **Visit start field** (Default: ZBE_2, PV1_44)
Several fields with a specific order can be specified for the visit start. So the HL7 interface will step through each field until it finds a valid entry (valid -> when a date is stated).
Fallbacks can be specified - if for some reasons - fields are not filled properly.

Possible values are:

- **EVN_2**
Segment EVN, field 2 (Recorded date/time)
- **EVN_3**
Segment EVN, field 3 (Date/time planned)
- **EVN_6**
Segment EVN, field 3 (Date/time event occurred)
- **ZBE_2**
Segment ZBE, field 2 (Start date/time)
- **PV1_44**
Segment PV1, field 44 (Admission date/time)
- **MSH_7**
Segment MSH, field 7 (Message date/time)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 56 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

- **Visit end field** (Default: <Empty List>)
Several fields with a specific order can be specified for the visit end. Same behavior as visit start field

Possible values are

- **EVN_2**
Described above
 - **EVN_3**
Described above
 - **EVN_6**
Described above
 - **ZBE_2**
Described above
 - **ZBE_3**
Segment ZBE, field 3 (End date/time)
 - **PV1_45**
Segment PV1, field 45 (Discharge date/time)
 - **MSH_7**
Described above
- **ZBE-Update: Date field** (Default: ZBE_2)
Several fields for visit updates can be selected with a specific order, when a ZBE segment is stated.

Possible values are:

- **ZBE_2**
Described above
- **ZBE_3**
Described above

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 57 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Please specify the settings for message type ADT-A01

Create patients: If cases are inserted ▼

☒ Update patients

☒ Import diagnoses

☒ Process ZBE updates

☐ Ignore consultation

☒ Save cases

☒ Save ward assignments

☐ Process message anew by default, when ZBE-UID non-existent

Ward assignments

Save visit start: Unchanged ▼

Save visit end: Unchanged ▼

Visit start field

EVN_2	→	ZBE_2	↑
EVN_3		PV1_44	
EVN_6	←		↓
MSH_7			

Visit end field

PV1_45	→		↑
ZBE_3			
EVN_2	←		↓
EVN_3			
EVN_6			

ZBE-Update: Date field

ZBE_3	→	ZBE_2	↑
	←		↓

4.4.2.1.1.3.3.2 Admission following preadmission

Behavior in case an admission message (A01) is received after a preadmission one (A05):

- Set up admission
- Delete preadmission and set up admission
- Ignore admission

Please determine how admissions (A01) are to be processed after preadmissions (A05)

Admission following preadmission Set up admission ▼

- Set up admission
- Delete preadmission and set up admission
- Ignore admission

Default: Set up admission

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 58 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.4 ADT-A02

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: If cases are inserted)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, EVN_6, EVN_2)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.5 ADT-A03 (Section)

4.4.2.1.1.3.5.1 ADT-03

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: <Empty List>)
- **Visit end field** (Default: ZBE_3, PV1_45, EVN_6)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 59 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.5.2 Conclude last valid unit assignment

Allows configuring whether unit assignments should be concluded even if the unit assignment is not the last in the case when receiving an A03 message.

Conclude unit assignments, even if the unit assignment is not the last in the case.

☐ Conclude last valid unit assignment

Default: false

4.4.2.1.1.3.6 ADT-A04 (Section)

4.4.2.1.1.3.6.1 ADT-A04

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: If cases are inserted)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44, EVN_6, EVN_2)
- **Visit end field** (Default: < Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 60 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.3.6.2 A04-insert settings

Specify when an A04 message should be processed:

- **Always**
- **Ignore if a case exists**
Will be ignored when a case exists with given case number

Reset to default
Edit
Accept
Cancel

Please determine in which cases the message A04 should be imported.

A04-insert settings always

always
Ignore if a case exists

Default: always

NOTE: To avoid transfers on outpatient cases you have to change this setting to 'Ignore if a case exists'

4.4.2.1.1.3.6.3 Ambulatory daytime admissions

Defines whether an A04-message should be treated as daytime admission.

- **Yes** (ward assignments will be inserted as a day entry)
- **No** (ward assignments will be inserted with the given date/time)

Reset to default
Edit
Accept
Cancel

Please determine if A04 news should be treated as a daytime admission in ambulatory cases.

☐ Yes ☒ No

Default: No

4.4.2.1.1.3.6.4 Close outpatient case on new admit

Defines whether an existing open outpatient case will be closed if a new equivalent one is inserted.

- **Yes**
- **No**

Reset to default
Edit
Accept
Cancel

Please define whether an open case shall be closed when receiving an admit message (A04) for the same ward.

☐ Yes ☒ No

Default: No


NOTE: Equivalent means that both cases are linked to the same ward ID and both cases are outpatient cases.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 61 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.6.5 Processing after pre-admission

Determines how processing is to be conducted after preadmissions when an A04 message is received:

- Process notification
- Delete preadmission and process notification
- Ignore notification

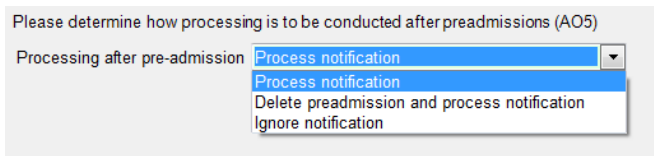

Default: Process notification

4.4.2.1.1.3.7 ADT-A05 (Section)

4.4.2.1.1.3.7.1 ADT-A05

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: If cases are inserted)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, EVN_3)
- **Visit end field** (Default: < Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 62 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.3.7.2 Close outpatient case during pre-admit

Defines whether an existing open outpatient pre-admitted case will be closed if a new equivalent one is inserted.

- **Yes**
- **No**

Reset to default
Edit
Accept
Cancel

Please determine if an open outpatient case may be closed through a pre-admit message with a new outpatient case for the same organizational unit.

☐ Yes
☒ No

Default: No

NOTE: Equivalent means that both cases are linked to the same ward ID and both cases are outpatient pre-admitted cases.

4.4.2.1.1.3.8 ADT-A06

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 63 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.3.9 ADT-A07

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.10 ADT-A08 (Section)

4.4.2.1.1.3.10.1 ADT-A08

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Save cases** (Default: no/false)
 - **Save ward assignments** (Default: no/false)
 - **Process A08 without ZBE segment**
Some systems also allow inserting of cases/ward assignment with A08. So this setting must be activated.
NOTE: This setting should not be activated without consultation of Cato Software Solutions!
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 64 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.10.2 *Processing after pre-admission*

See A04 (0)

4.4.2.1.1.3.11 ADT-A11

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.12 ADT-A12

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 65 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.13 ADT-A13

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.14 ADT-A21

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 66 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.15 ADT-A22

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.16 ADT-A28

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: no/false)
 - **Save ward assignments** (Default: no/false)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 67 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.17 ADT-A29

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: no/false)
 - **Save ward assignments** (Default: no/false)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.18 ADT-A31

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: no/false)
 - **Save ward assignments** (Default: no/false)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 68 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.19 ADT-A34

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
 - **Always**
 - **If the merge-slave exists**
Then the merge-winner will only created when the merge-loser (slave) exists
 - **Never**
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.20 ADT-A38

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 69 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.21 ADT-A40

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
 - **Always**
 - **If the merge-slave exists**
Then the merge-winner will only created when the merge-loser (slave) exists
 - **Never**
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 70 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.22 ADT-A45

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Always)
 - **Always**
 - **If the merge-slave exists**
Then the merge-winner will only created when the merge-loser (slave) exists
 - **Never**
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.23 ADT-A52

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 71 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.1.3.24 ADT-A53

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: Never)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: ZBE_2, PV1_44)
- **Visit end field** (Default: <Empty List>)
- **ZBE-Update: Date field** (Default: ZBE_2)

4.4.2.1.1.3.25 Skyline

Description: See ADT-A01, Defaults as described below

- **Create patients** (Default: always)
- **Update patients** (Default: yes/true)
- **Import diagnoses** (Default: yes/true)
- **Process ZBE updates** (Default: yes/true)
- **Ignore consultation** (Default: no/false)
- **Save cases** (Default: yes/true)
 - **Save ward assignments** (Default: yes/true)
- **Process message anew by default, when ZBE-UID non-existent** (Default: no/false)
- **Save visit start** (Default: Unchanged)
- **Save visit end** (Default: Unchanged)
- **Visit start field** (Default: PV1_44, ZBE_2, EVN_6, EVN_2)
- **Visit end field** (Default: PV1_45)
- **ZBE-Update: Date field** (Default: ZBE_2)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 72 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.3.26 Reply to already received messages with Application Accept (AA)

Specify if already received messages (Message Control ID is already stored in the database) should be accepted or rejected.

Please specify if you want to reply to already received messages with Application Accept (AA) instead of Application Reject (AR)

☐ Reply to already received messages with Application Accept (AA)

Default: no/false

NOTE: In some cases when a sending system receives an 'Application Reject' it resends the same message again. This would lead to a never-ending process of sending this message because the cato® HL7 Interface will always reply with an 'Application Reject' if the message was already received.

So changing this setting to yes/no would be reasonable.

4.4.2.1.1.3.27 Special option: ignore day stays

This option should only be activated if you wish to ignore stays which are communicated as day cases.

Please only activate this option if you wish to ignore stays which are communicated as day cases.

☐ Special option: ignore day stays

Default: false

4.4.2.1.1.4 Import Configuration

Defines how specific entities should be imported. (Further described in 3.1 Synchronization)

Entity	Identification-Mode	UID-Field
Patient	by UID	PatNr
WardAssignment	by name or UID	EUID
Case	by name	UID
Diagnosis	by name	UID

Default:


- **Patient**
Identification Mode: by UID, UID-Field: PatNr
- **WardAssignment**
Identification Mode: by name or UID, UID-Field: EUID
- **Case**
Identification Mode by name
- **Diagnosis**
Identification Mode by name

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 73 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.5 Shutdown HL7 service on critical errors

Defines whether the service should be shut down when a critical error occurs.

- Yes/true
- No/false

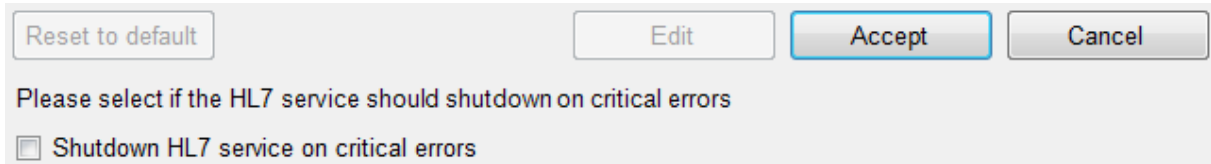

Default: No/false

NOTE: Please do not change this setting unless there is no other way to handle critical errors.

4.4.2.1.1.6 File forwarding

If you have a file-based installation and there is a need of duplicate (or multiply) messages you can select a source folder and file type and one or more destination folders.

Use authentication if needed.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 74 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please configure file forwarding

☒ Active

Source

Network path: a
Search

☐ Use authentication

User:

Password:

Test read access
Test write access

File type: *.hl7

Targets

New
Maintain
Delete

Path
b
c

NOTE: In this sample configuration all incoming files of type 'hl7' in the path a will be copied to path b and c.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 75 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.1.7 Settings for files with several messages

It is possible to configure HL7 service to process files containing multiple messages.

Default: not active

Please specify settings for files with several messages

☒ Active

Encoding: Western European (Windows) ▼

Incoming data

Network path: Search

☐ Use authentication

User:

Password:

Test read access Test write access

File type:

Outgoing data

Network path: Search

☐ Use authentication

User:

Password:

Test read access Test write access

File type:

In the incoming data section, a network path expecting the files with multiple messages can be configured.

In the outgoing data section, a network path where the split files will be saved can be configured.

The file type field expects a file mask in the form “*.hl7” for instance.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 76 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.2 Patients (Section)

4.4.2.1.2.1 Sex IDs (Section)

4.4.2.1.2.1.1 Sex identifier male

Specify comma separated sex identifiers for male patients.

Reset to default
Edit
Accept
Cancel

Please enter sex identifiers for male patients. Several comma-separated values can be entered.

Sex identifier male

Default: M

4.4.2.1.2.1.2 Sex identifier female

Specify comma separated sex identifiers for female patients.

Reset to default
Edit
Accept
Cancel

Please enter sex identifiers for female patients. Several comma-separated values can be entered.

Sex identifier female

Default: F;W

4.4.2.1.2.2 Select preferred identification field for patients

Specify which field contains the identifier for synchronizing patients:

- **PID<3>**
Segment PID, field 3 (Patient identifier list)
- **PID<2>**
Segment PID, field 2 (Patient ID)

Reset to default
Edit
Accept
Cancel

Please set the preferred identification field for patients

Select preferred identification field for patients

PID<3>
PID<2>

Default PID<3>

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 77 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.2.3 Integrate client name in patient number

Specify if you want to integrate the client into the patient number (separated by an underscore)

Reset to default
Edit
Accept
Cancel

Please set this option if you wish to integrate the client name in the patient number (Client_Patient number)

☐ Integrate client name in patient number

Default: no/false

NOTE: You have multiple clients in one database with a separate set of Patient IDs (so a Patient ID can occur multiple) you have to distinguish them by a synthetic Patient ID.

4.4.2.1.2.4 Field for social security number

Specifies the field containing the patient's social security number.

Reset to default
Edit
Accept
Cancel

Please define the field for importing the social security number of the patient

Field for social security number PID<19>

PID<19>
IN1<36>

Default: PID<19>

4.4.2.1.2.5 Accept 4-digit Social Security number

Specify if you want to accept 4-digit social security numbers.

Reset to default
Edit
Accept
Cancel

Please determine, whether the sending system transmits 4-digit Social Security Numbers.

☐ Accept 4-digit Social Security Number

Default: no/false

NOTE: In Austria the SSNR consists of a 4-digit number and then 6-digit coded birth date (ddmmyy) and some systems only send the first 4 digits of this SSNR, so the SSNR can be automatically put together.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 78 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.2.6 Disable change check

Specify if you want to disable the change check.

Reset to default
Edit
Accept
Cancel

Attention: When setting this option, it is possible that changes on more than one attribute (Name, Sex, Birthdate) are applied.

☐ Disable Change-Check

Default: no/false

NOTE: The cato® HL7 Interface checks if there are relevant changes on existing patients (changing name, sex, birth date; at least 2 changes must apply) so errors by the sending system can be avoided.

4.4.2.1.2.7 Ignore patient prefix

Set this option if the patient prefix (Segment PID, field 5, component 5) should be ignored.

Reset to default
Edit
Accept
Cancel

Select this option if you do not wish to adopt the patient prefix.

☐ Ignore patient prefix.

Default: no/false

4.4.2.1.2.8 Ignore patient suffix

Set this option if the patient prefix (Segment PID, field 5, component 4) should be ignored.

Reset to default
Edit
Accept
Cancel

Select this option if you do not wish to adopt the patient suffix.

☐ Ignore patient suffix

Default: no/false

4.4.2.1.2.9 Delete patients without therapy plan

Set this option if you want to delete a patient without a therapy plan when receiving an A03 message (Discharge a patient)

Reset to default
Edit
Accept
Cancel

Activate this option if you wish to delete patients without therapy plans upon reception of an A03 message.

☐ Delete patients without therapy plan

Default: no/false

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 79 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.2.10 Save patient number also in the barcode field

Set this option if the patient number shall also be saved in the barcode field.

Reset to default
Edit
Accept
Cancel

Please determine if the patient number shall also be saved in the barcode field.

☐ Save patient number also in the barcode field.

Default: no/false

4.4.2.1.2.11 Attach client to first name

Set this option if you want to attach the client to the patient's first name.

Reset to default
Edit
Accept
Cancel

Please determine if the client should be attached to the first name

☐ Attach client to first name

Default: no/false

NOTE: When you have multiple clients in one database with a separate set of Patient IDs (so Patient IDs can occur multiple) users can easier find and identify a patient when the client is attached to the first name.

4.4.2.1.2.12 Standard patient type for import

Specifies the standard patient type for import.

Reset to default
Edit
Accept
Cancel

Please determine the standard patient type for importing

Standard patient type for import Both

- Both
- Cyto
- TPN

Default: Both (Cyto and TPN)

4.4.2.1.2.13 Extract ZIP and city from County/Parish code

Set this option to true if you want to import the ZIP and City from PID-County/Parish code.

Reset to default
Edit
Accept
Cancel

Please specify if the ZIP and city should be extracted from County/Parish code.

☐ Extract ZIP and city from County/Parish code

Default: no/false

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 80 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.2.14 Patient's 2nd phone number from PID<14>

Set this option to true if you want to import the second phone number of the patient from PID<14>.

Reset to default	Edit	Accept	Cancel
Import patient's second phone number from PID<14>			
☐ Patient's 2. phone number from PID<14>			

4.4.2.1.2.15 Patients without a valid unit assignment to be entered as blocked

Specify whether a new patient without (valid) ward assignments will be entered as blocked.

- **Yes/true**
- **No/false**

Reset to default	Edit	Accept	Cancel
Should patients without a valid unit assignment be entered as blocked?			
☐ Patients without a valid unit assignment to be entered as blocked			

Default: yes/true

4.4.2.1.2.16 Save case# as patient barcode

Specify for which incoming message types the case# should also be stored in the patient barcode field.

- **A01, A02, A03, A04, A05, A06, A07, A08, A11, A12, A13, A19, A21, A22, A23, A28, A29, A31, A34, A38, A39, A40, A45, A52, A53, R01, Skyline**

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 81 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please choose the types of messages for which the case# will be saved as patient barcode.

A01 A02 A03 A04 A05 A06 A07 A08 A11 A12 A13 A19 A21 A22 A23 A28 A29 A31 A34 A38 A39 A40 A45 A52 A53 R01 Skyline	>> >>> << <<<	
---	------------------------	--

Default: <Empty list>

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 82 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.3 Stays (Section)

4.4.2.1.3.1 Ward's ID

Specify which field(s) contains the ward ID (See 0

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 83 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Ward synchronization):

- **PV1_3_PointOfCareID**
Segment PV1, field 3, component 1
- **PV1_3_PointOfCare_Facility**
Segment PV1, field 3, components 1 & 4 separated by @ e.g.: INT@MED1
- **PV1_3_PointOfCareID_PV1_39_ServicingFacility**
Segment PV1, field 3, component 1 & segment PV1, field 39 separated by @ e.g.: CYTO@S
- **PV1_3_PointOfCareID_AttendingDoctor**
Segment PV1, field 3, component 1 & segment PV1, field 7 separated by @ e.g.: A@B
- **PV1_7_AttendingDoctor**
Segment PV1, field 7
- **PV1_3_Facility**
Segment PV1, field 3, component 4
- **PV1_3_Facility_PointOfCare**
Segment PV1, field 3, components 4 & 1 separated by @ e.g.: MED1@INT
- **PV1_39_ServicingFacility**
Segment PV1, field 39
- **PV1_3_RoomID**
Segment PV1, field 3, component 2
- **PV1_Agfa_Orbis_Facility**
Special Setting for Agfa Orbis, please contact Cato Software Solutions if none of the above settings is appropriate
- **PV1_Agfa_Orbis_Facility_LocationStatus**
Special Setting for Agfa Orbis, please contact Cato Software Solutions if none of the above settings is appropriate

Reset to default
Edit
Accept
Cancel

Please specify the field from which the ward's unique ID shall be read.

Ward's ID PV1_3_Facility_PointOfCare

- PV1_3_Facility_PointOfCare
- PV1_3_PointOfCareID
- PV1_3_PointOfCare_Facility
- PV1_3_PointOfCareID_PV1_39_ServicingFacility
- PV1_3_PointOfCareID_AttendingDoctor
- PV1_7_AttendingDoctor
- PV1_3_Facility
- PV1_3_Facility_PointOfCare
- PV1_39_ServicingFacility
- PV1_3_RoomID
- PV1_Agfa_Orbis_Facility
- PV1_Agfa_Orbis_Facility_LocationStatus

Default: PV1_3_PointOfCareID

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 84 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.3.2 Find wards also by client

Specify if a ward should be also identified by the client.

Reset to default
Edit
Accept
Cancel

Please specify if a ward should be also identified by the client (coded in the hospital)

☐ Find wards also by client

Default: no/false

NOTE: This is useful if you have multiple wards with the same identifier but in different hospitals/facilities.

4.4.2.1.3.3 Field for stay number

Specify which field contains the stay number (case number, admission number)

- **Ignore**
- **PV1_19**
Segment PV1, Field 19 (Visit Number)
- **ZBE_UID**
Segment ZBE, Field 1
- **PID_4**
Segment PID, Field 4, Alternate Patient ID

Reset to default
Edit
Accept
Cancel

Please specify the field (of the message) from which the stay number shall be determined

Field for stay number PV1_19 ▼

Ignore
PV1_19
ZBE_UID
PID_4

Default: PV1_19

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 85 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.3.4 Database field for saving the stay number

Specify the database field where the stay number shall be saved:

- **Do not apply**
Ignore, no stay number will be saved
- **Admission no.**
Admission Number
- **Case number**

Reset to default
Edit
Accept
Cancel

Please specify the database field where the stay number shall be saved

Database field for saving the stay number Case number ▼

Do not apply
Admission no.
Case number

Default: Case number

4.4.2.1.3.5 Saving out-patient stays

Define which type of outpatient stays should be saved.

- **All out-patient stays**
- **No out-patient stays**
- **All except those irrelevant for billing**
Account types can be entered
- **None except those relevant for billing**
Account types can be entered

Account types can be entered separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Please specify whether and which out-patient stays of the patient shall be saved in the transfer table (patient's assigned wards).

Save type: All out-patient stays ▼

Account type:

Default: All out-patient stays

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 86 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.2.1.3.6 Parallel outpatient cases

Set this option if the sending system provides parallel outpatient cases.

Reset to default
Edit
Accept
Cancel

Activate this option if outpatient cases are handled parallel to inpatient cases in your hospital.

☐ Parallel outpatient cases

Default: no/false

NOTE: When setting this option, after a discharge (of an inpatient case) the last outpatient case will be reopened in the ward assignment table.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 87 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.3.7 Stay types (Section)

4.4.2.1.3.7.1 Optional tag for inpatient stays

You can specify optional tags for inpatient stays if the sending system provides other tags than the default values. They are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Set this field with the tag for inpatient stays, which will be transmitted by the external system in the PV1, Field - Patient-Class segment. (separated by ;)

Optional tag for inpatient stays I;S

Default: I;S

4.4.2.1.3.7.2 Optional tag for outpatient stays

You can specify optional tags for outpatient stays if the sending system provides other tags than the default values. They are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Set this field with the tag for outpatient stays, which will be transmitted by the external system in the PV1, Field - Patient- Class segment. (separated by ;)

Optional tag for outpatient stays O;A

Default: O;A

4.4.2.1.3.7.3 Optional tag for partly in-patient stays

You can specify optional tags for partly inpatient stays if the sending system provides other tags than the default values. They are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Set this field with the flag for partly in-patient stays, which will be transmitted in the PV1 segment Patient-Class by the third party system.

Optional flag for partly in-patient stays T;TS

Default: T;TS

4.4.2.1.3.7.4 Interpret missing stay type as

Some systems do not send a stay type but in fact it is e.g. an outpatient stay.

- **inpatient stay**
- **outpatient stay**
- **Partly in-patient**
- **invalid PV1 segment**

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 88 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please specify how a missing stay type shall be interpreted (PV1-2)

Interpret missing stay type as invalid PV1 segment ▼

- invalid PV1 segment
- inpatient stay
- outpatient stay
- Partly in-patient
- invalid PV1 segment

Default: invalid PV1 segment

4.4.2.1.3.8 Include private ID (PID-21) in ward's ID

Set this option if you want to include segment PID, field 21 in the ward id (separated by underscore)

Reset to default
Edit
Accept
Cancel

Please specify whether the private ID (PID-21) shall be included in the ward's UID.

☐ Include private ID (PID-21) in ward's ID

Default: no/false

4.4.2.1.3.9 Case-related editing

Set this option if you want to include the visit number to visit and case UID, separated by a dot.

Reset to default
Edit
Accept
Cancel

Please specify whether messages shall be processed case-related.

☒ Case-related editing

Default: yes/true

NOTE: This option must be set when the identifier stated in the ZBE-segment is not unique over all visits but only in the specific case.

4.4.2.1.3.10 Don't change existing ward assignments in case of unknown wards

Please specify whether an existing ward assignment should not be truncated when receiving an action with unknown ward.

- A01, A02, A03, A04, A05, A06, A07, A08, A11, A12, A13, A19, A21, A22, A23, A28, A29, A31, A34, A38, A39, A40, A45, A52, A53, R01, Skyline

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 89 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please choose the types of messages for which ward assignments will not be changed in case of unknown wards.

A01 A02 A03 A04 A05 A06 A07 A08 A11 A12 A13 A19 A21 A22 A23 A28 A29 A31 A34 A38 A39 A40 A45 A52 A53 R01 Skyline	>> >>> << <<<	
---	------------------------	--

Default: <Empty list>

Example: If there is an admission for case 12345 on ward ABC (May 4th, 2012, 10am) and then you receive a transfer to ward XYZ (which is not known) by May 5th, 11am usually the first ward assignment (ABC) will be cut and it ends by May 5th, 10:59am. With this setting activated the prior ward assignment will NOT be cut and remains open.

4.4.2.1.3.11 Evaluate PV1-10 for ward-Ids

Determines whether PV1-10 shall be evaluated for wards-Ids

Please determine whether PV1-10 shall be evaluated for ward-Ids

☐ Evaluate PV1-10 for ward-Ids

Default: false

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 90 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.4 Cost units

4.4.2.1.4.1 Settings for insured status import

Specifies the settings for insurance status import.

Please determine the settings for insured status import

☐ Active

Field:

Country:

☐ Create automatically, if non-existent

Default: Not active / Field = PV1<20>

NOTE: Only relevant for Austria, Germany, Switzerland

4.4.2.1.4.2 Setting for co-payment waived import

Allows configuring the settings for co-payment waived import

Please specify the settings for co-payment waived import

☐ Active

Field:

Default: Inactive / Field = IN1.15.3

4.4.2.1.4.3 Import validity of the insurance card

Set this option if you want to save the validity of the insurance card (Segment IN1, field 13). You can also specify a format for the date.

Reset to default Edit Accept Cancel

Please determine if the validity of the insurance card should be imported

☐ Activate

Format:

Default: not active

Only relevant for Austria, Germany, Switzerland

4.4.2.1.4.4 Import insurance card number from IN1<36>

Specify whether you want to import the insurance card number.

Reset to default Edit Accept Cancel

Import patient's insurance card number from IN1<36>

☐ Import insurance card number from IN1<36>

Default: no/false

Only relevant for Austria, Germany, Switzerland

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 91 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.5 Medical Results (Section)

NOTE: Medical results are only imported if a valid license is available!

4.4.2.1.5.1 Medical result ID

Define which field contains the medical result's unique ID

- **OBX_3_Identifier**
Segment OBX, field 3, component 1
- **OBX_3_IdentifierAndText**
Combination of segment OBX, field 3, component 1 & 2
- **OBR_4**
Segment OBR, field 4

Reset to default
Edit
Accept
Cancel

Please specify the field the medical result's unique ID shall be read from

Medical result ID OBX_3_Identifier

OBX_3_Identifier
OBX_3_IdentifierAndText
OBR_4

Default: OBX_3_Identifier

4.4.2.1.5.2 Medical result request number

Define which field contains the medical result request number.

- **OBR2**
Segment OBR, field 2
- **ORC3**
Segment ORC, field 3
- **OBR3**
Segment OBR, field 3

Reset to default
Edit
Accept
Cancel

Please define which field the medical result request number shall be retrieved from.

Medical result request number OBR2

OBR2
ORC3
OBR3

Default: OBR2

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 92 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.5.3 Deleting all medical results with a certain request number before updating

Set this option if you want to delete (in fact: cancel) all medical results with the given request number.

Reset to default
Edit
Accept
Cancel

Activate this option if you wish to delete all medical results with a received RequestNumber prior to the updating of medical results.

☐ Deleting all medical results with a certain request number before updating.

Default: no/false

NOTE: This option can be useful if e.g. 4 valid medical results were sent in a previous message and now 1 of these is invalid because of an error. Then all 4 will be cancelled and only the 3 valid results will be inserted again.

4.4.2.1.5.4 Valid observation value types (OBX-2)

You can specify valid observation values types, if the sending system provides other than the default entries. Entries are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Please specify valid observation value types (OBX-2)

Valid observation value types (OBX-2)

Default: ST;NM;

4.4.2.1.5.5 Medical result status (preliminary) (OBX-11)

You can specify an observation result status for temporary medical results, if the sending system provides this. Entries are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Please specify, which observation result status (OBX-11) shall be defined as invalid (separated by commas).

Observation result status (invalid) (OBX-11)

Default: <Empty>

4.4.2.1.5.6 Medical result status (invalid) (OBX-11)

You can specify an observation result status for invalid medical results, if the sending system provides this. Entries are separated by comma or semicolon.

Reset to default
Edit
Accept
Cancel

Please specify, which observation result status (OBX-11) shall be defined as temporary (separated by commas).

Observation result status (temporary) (OBX-11)

Default: <Empty>

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 93 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.5.7 Save request number as remark

You can define that the request number is saved as additional remark to the medical result. You can also specify a text before and a text after the request number.

Currently there is one variable which can be used in both text fields:


- **[DateTime]**

Simply insert [DateTime] into the text fields where you want to put it:

E.g:

Text before: Request number:


Text after: at [DateTime]

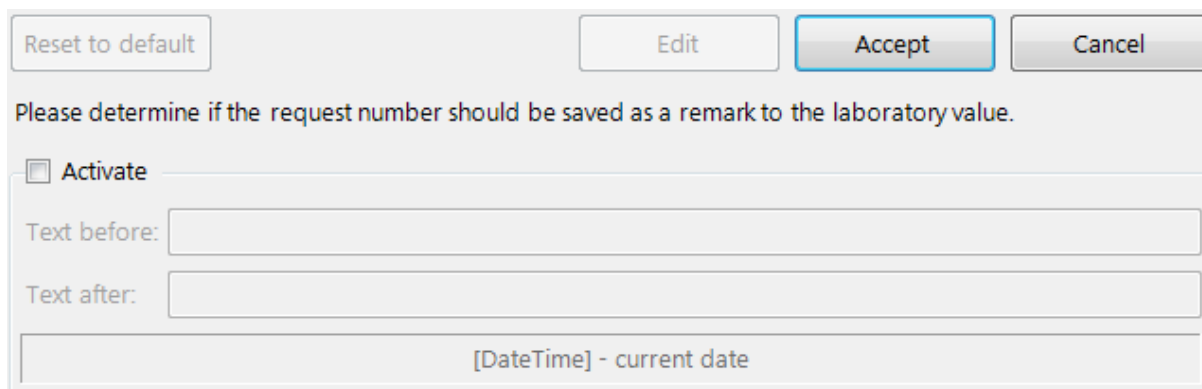

Default: Not active

4.4.2.1.5.8 Import standard range

By activating this setting standard ranges for each medical result value are imported.

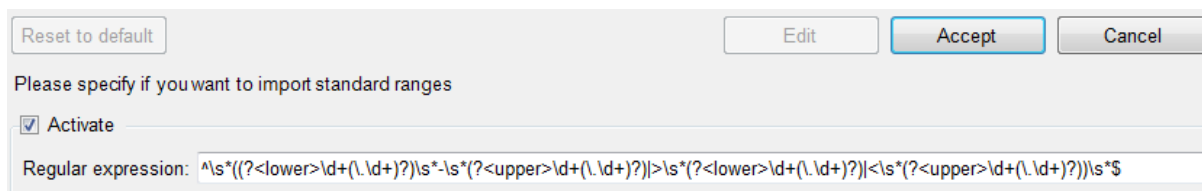

NOTE: The regular expression should not be changed without consultation of Cato Software Solutions.

Default: Active

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 94 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.5.9 Apply current time as medical result date

Allows configuring whether the current time should be used as a medical result date if no other time has been specified.

ATTENTION: This setting can lead to unintended problems and should only be used in case of necessity.

Default: False

4.4.2.1.6 Diagnoses (Section)

Note: Diagnoses are only imported if a valid license is available!

4.4.2.1.6.1 Import diagnoses

Set this option if you want to import diagnoses.

Reset to default
Edit
Accept
Cancel

Please specify whether diagnoses shall be imported or not

☒ Import diagnoses


Default: yes/true

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 95 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.6.2 Diagnosis type

You can specify which type of diagnoses is sent.

- ICD10
- ICD9

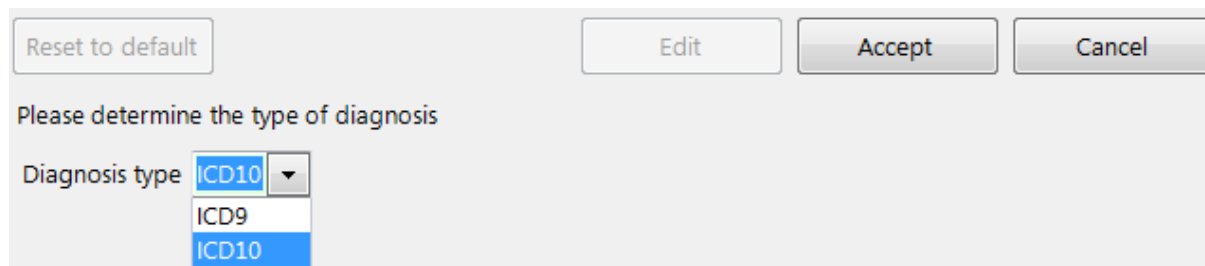

Reset to default Edit Accept Cancel

Please determine the type of diagnosis

Diagnosis type ICD10 ▼

ICD9

ICD10

Default: ICD10

4.4.2.1.7 HL7 Query (Section)

NOTE: For internal use

4.4.2.1.7.1 HL7 query – activation

NOTE: For internal use

4.4.2.1.8 RDE medication import

Medication import is a part of the HL7-interface but is described in the RDE-manual.

4.4.2.1.9 Logging (Section)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 96 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.9.1 General logging settings

You can specify general logging settings:

- **Write to Windows-Event log** (Default: yes/true)
- **Time before deleting valid messages in hours** (Default: 120)
When setting to 0, valid message will never be deleted (useful for a testing environment)
- **Time before deleting invalid message in hours** (Default: 240)
When setting to 0, invalid message will never be deleted (useful for a testing environment)
- **Send mails** (Default: no/false)
 - **SMTP-Server**
 - **SMTP-Authentication**
 - **Login**
 - **Password**
 - **Sender**
 - **Receiver**
 - **Log Level:**
 - **Errors only** (default)
 - **Errors and warnings**
 - **Any information**
 - **Shortest time between two mails (in minutes)**
 - **Allow HTML**



Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 97 of 128
Revision:	01	
Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
Template:	FO Other Documents (MWS-FO-OD)	
Template-Rev.:	02	

[Reset to default](#)[Edit](#)[Accept](#)[Cancel](#)

Please specify the general logging settings

☒ Write to Windows-Eventlog

Time before deleting valid messages: hours

Time before deleting invalid messages: hours

☒ Send mails

SMTP-Server [Test](#)

☐ SMTP Authentication

Login:

Password:

Sender:

Receiver:

Log Level: **Errors Only** ▼

Shortest time between two mails: minutes

☐ Allow HTML

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 98 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.2.1.9.2 Log incorrect wards

Set this option if you want to save a copy of a message with unknown ward to a network path.

You can use authentication if needed.

Reset to default
Edit
Accept
Cancel

Activate this option if messages with nonexistent wards shall be saved in a certain folder.

☐ Active

Network path:
Search

☐ Use authentication

User:
Password:

Test read access
Test write access

Default: Inactive

4.4.2.1.9.3 Log accepted messages

Set this option if you want to save a copy of an accepted message to a network path.

You can use authentication if needed.

Reset to default
Edit
Accept
Cancel

Activate this option if accepted messages shall be saved in the file system.

☐ Active

Network path:
Search

☐ Use authentication

User:
Password:

Test read access
Test write access

Default: Inactive

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 99 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	



4.4.3 Database Settings

For each connected database settings can be added.

4.4.3.1 Manage database settings

4.4.3.1.1 Add database settings

By right-clicking on 'System settings' new database settings can be added:

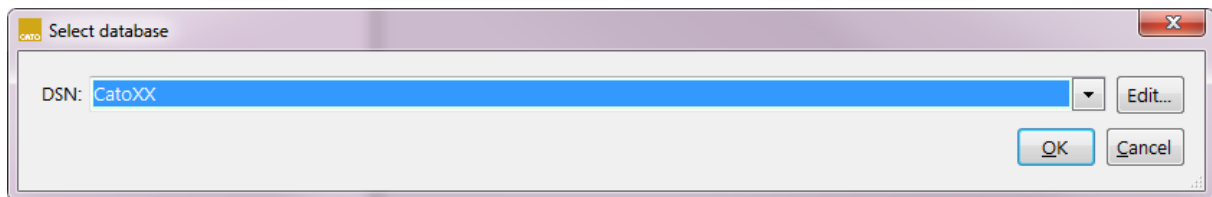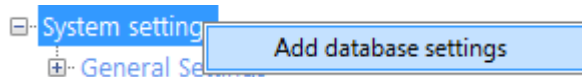



By clicking OK new database settings will be added:



 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 100 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.3.1.2 Edit database settings


By right-clicking on specific database settings, they can be modified:

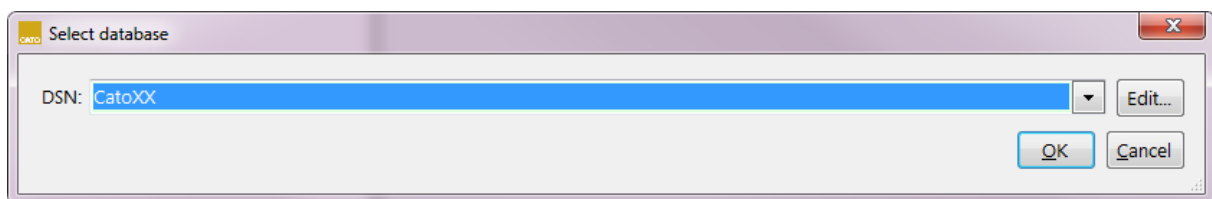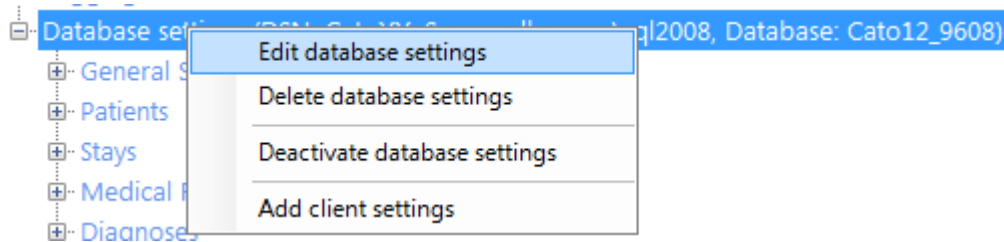

By clicking OK database settings are changed.

4.4.3.1.3 Delete database settings

By right-clicking on specific database settings they can be deleted:

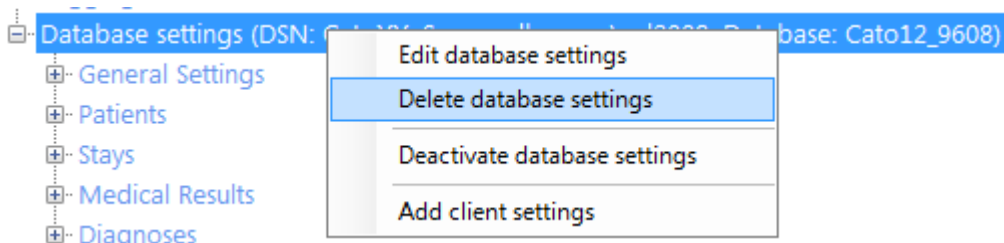


NOTE: By deleting database settings all existing sub-levels (client settings, connection settings) will also be deleted!


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 101 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.3.1.4 Deactivate database settings

By right-clicking on specific database settings they can be set inactive (when active):

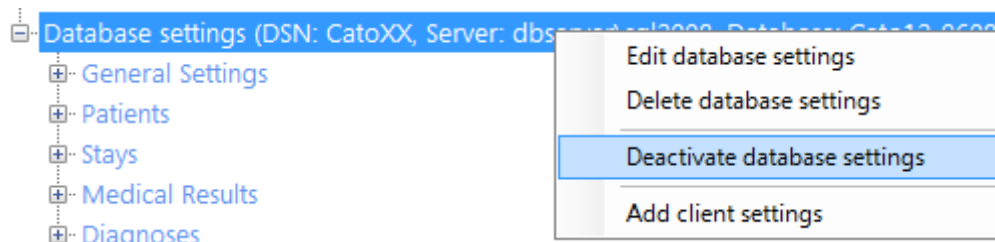


Now database settings (and all sub-levels) are inactive:

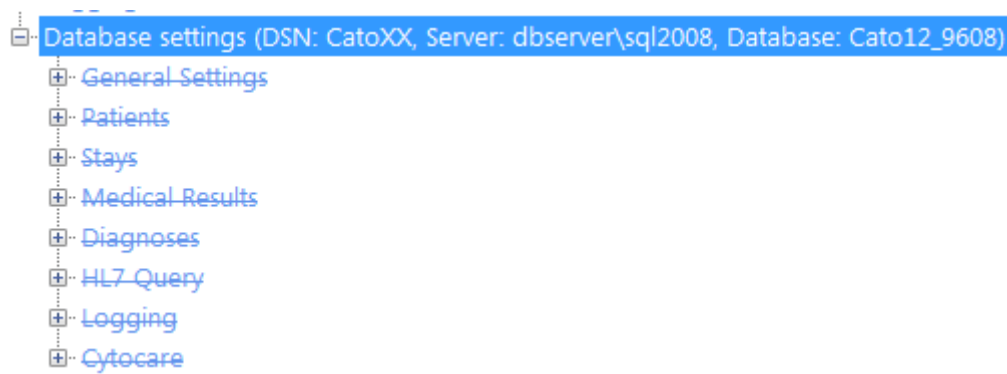

NOTE: This is useful for testing environments to easily deactivate database settings. After testing they can be reactivated.

4.4.3.1.5 Activate database settings

By right-clicking on specific database settings database settings can be set activate (when inactive):

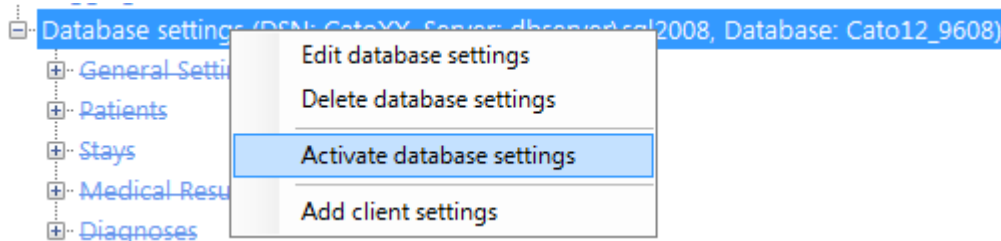

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 102 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.3.2 Settings


NOTE: Only new (added at this level) settings are listed.

4.4.3.2.1 General Settings (Section)

4.4.3.2.1.1 Database user

An import user must be specified. This user is listed in the cato® event log for all actions that happen by processing a message.

E.g.: Import patients, import medical results, and so on

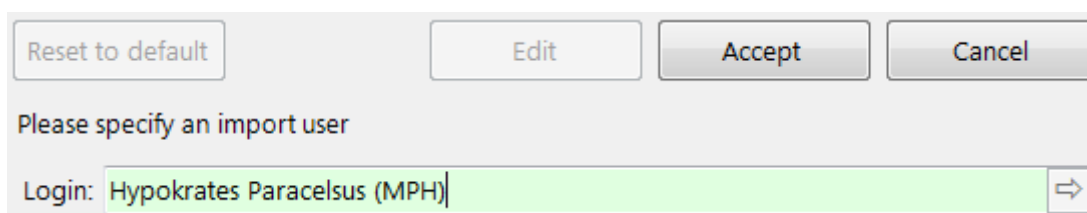


Default: <No selection>

Note: When the import user isn't set the HL7Service will not start! An error message will be written to the Windows event log.

4.4.3.2.1.2 Messages (Section)

4.4.3.2.1.2.1 Maximum number of messages to be processed

You can set how many messages will be processed at once.

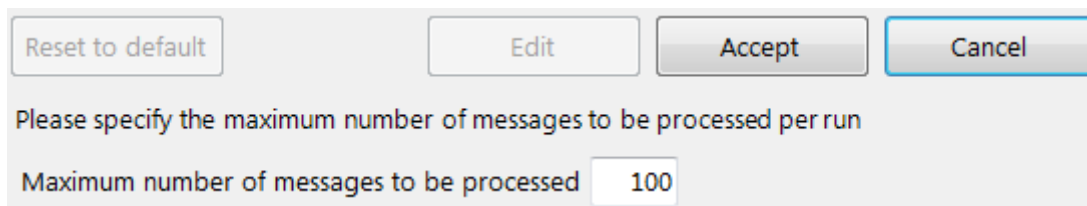

Default: 100

NOTE: By settings this option to 0 all available messages will be processed.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 103 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.3.2.2 Cost Units (Section)

4.4.3.2.2.1 Create cost units

If you need to automatically create cost units then you have to active this setting.

You can specify following settings:

- UID-Field
 - IN1_3 segment IN1, field 3 (Insurance Company ID) (default)
- Billing models
 - Outpatient
 - n.a. – none available (default)
 - all existing billing models
 - Alternative outpatient
 - n.a. – none available (default)
 - all existing billing models
 - Oncol. Office
 - n.a. – none available (default)
 - all existing billing models
 - Alternative oncol. Office
 - n.a. – none available (default)
 - all existing billing models
- Valid message types (Defines for which message type cost unit should be created)
 - A01, A02, A03, A04, A05, A06, A07, A08, A11, A12, A13, A19, A21, A22, A23, A28, A29, A31, A34, A38, A39, A40, A45, A52, A53, R01, Skyline

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 104 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please define settings for cost unit creation

☒ Active

UID-Field: IN1 3

Billing models

Outpatient: n. a. *

Alternative outpatient: n. a.

Oncol. office: n. a.

Alternative oncol. office: n. a.

Valid message types

A01	>>	
A02	>>>	
A03	<<	
A04	<<<	
A05		
A06		
A07		
A08		
A11		
A12		
A13		
A19		
A21		
A22		
A23		
A28		
A29		
A31		
A34		
A38		
A39		
A40		
A45		
A52		
A53		
R01		
Skyline		

NOTE: Only relevant for Austria, Germany, Switzerland

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 105 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.3.2.2 Import cost units

If you need to import cost units then you have to active this setting.

You can specify following settings:

- UID-Field
 - **IN1_3 – Segment IN1, field 3 (Insurance Company ID) (default)**
- Valid from-field
 - IN1_12 – Segment IN1, field 12 (Plan Effective Date) (default)
 - IN1_29 – Segment IN1, field 29 (Verification Date/Time)
 - EVN_2 – Segment EVN, field 2 (Recorded Date/Time)
 - PV1_44 – Segment PV1, field 44 (Admission Date/Time)
- Valid from-regex
 - A regular expression for parsing the provided date/time
- Valid to-field
 - None (default) – The cost unit will be stored with no end-date/time
 - IN1_13 – Segment IN1, field 13 (Plan Expiration Date)
 - PV1_45 – Segment PV1, field 45 (Discharge Date/Time)

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 106 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

Reset to default
Edit
Accept
Cancel

Please define settings for cost unit import

☒ Active

UID-Field: IN1 3

Valid from-field: IN1 12

Valid from-regex: ^(?<year>\d{4})(?<month>\d{2})(?<day>\d{2})

Valid to-field: None

Valid to-regex: ^(?<year>\d{4})(?<month>\d{2})(?<day>\d{2})

Valid message types

A01
A02
A03
A04
A05
A06
A07
A08
A11
A12
A13
A19
A21
A22
A23
A28
A29
A31
A34
A38
A39
A40
A45
A52
A53
R01
Skyline

>>
>>>
<<
<<<

NOTE: Only relevant for Austria, Germany, Switzerland

4.4.3.2.3 Medical Results (Section)

Note: Medical results are only imported if a valid license is available!

4.4.3.2.3.1 Database user

You can specify a user for importing medical results. If no user is specified the default database user specified in General Settings is used.

Reset to default
Edit
Accept
Cancel

Please specify an import user

Login:

Default: <No selection>

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 107 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.3.2.3.2 Medical result values

All medical results which shall be imported must be configured here:

Reset to default
Edit
Accept
Cancel

Please specify the medical result values to be imported

New
Edit
Delete

Active	UID send. system	Unit send. system	Med. result	Unit in Cato	Factor

By clicking the button 'New' you can enter a new medical result:

Reset to default
Edit
Accept
Cancel

Please specify the medical result values to be imported

☒ Active

UID in sending system:

Unit in sending system:

Medical result: 

Unit in Cato:

Conversion factor:

Add
Cancel

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 108 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

Several fields must be entered:

- **Active**
A medical result can entered but set to inactive, when unused temporarily
- **UID in sending system**
The unique identifier of the medical result how it is provided by the sending system
- **Unit in sending system**
The unit of the medical result how it is provided by the sending system
- **Medical result**
The medical result to be imported
- **Unit in Cato**
Can be entered manually but also by clicking the button 'Select' (Enabled when a medical result is selected)
- **Conversion factor**
Sometimes a unit provided by the sending system is different to the unit stored in the cato® database. You can specify a conversion factor.

NOTE: The incoming value from the sending system is multiplied with the conversion factor and this value will be stored as medical result value.

Conversion factor example 1:

Unit in sending system: milligram per liter

Unit in cato: milligram per milliliter

Conversion factor must be set to 0.001

Conversion factor example 2:

Unit in sending system: gram per deciliter

Unit in cato: milligram per deciliter

Conversion factor must be set to 1000

By pressing the button 'Add' the medical result will be added.

By pressing the button 'Cancel' no entry will be added.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 109 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

An example for leukocytes:

Reset to default
Edit
Accept
Cancel

Please specify the medical result values to be imported

New
Edit
Delete

Active	UID send. system	Unit send. system	Med. result	Unit in Cato	Factor
☒	LEU	Billion/L	Leukocytes	10^9/liter	1

The sending system provides leukocytes with identifier “LEU” and unit “Billion/L” (Billion per liter). In cat^o® the unit for leukocytes is stored as 10^9/liter, which means the same as Billion/L, so the factor is 1.

By pressing the button ‘Edit’ the settings for the selected medical result can be changed.

By pressing the button ‘Delete’ the selected medical result will be deleted.

NOTE: Entering units is an additional security feature! A medical result will not be saved when the units do not match:

- When the sending system provides another unit than configured
- When the unit in cat^o® has been modified

4.4.3.2.4 Diagnoses (Section)

Note: Diagnoses are only imported if a valid license is available!

4.4.3.2.4.1 Database user

You can specify a user for importing diagnoses. If no user is specified the default database user specified in General Settings is used.

Reset to default
Edit
Accept
Cancel

Please specify an import user

Login:

Default: <No selection>



 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 110 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.4 Client Settings

4.4.4.1 Manage client settings

4.4.4.1.1 Add client settings


By right-clicking on specific database settings new client settings can be added:

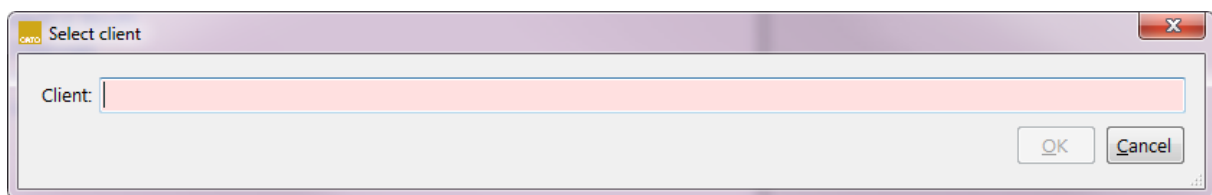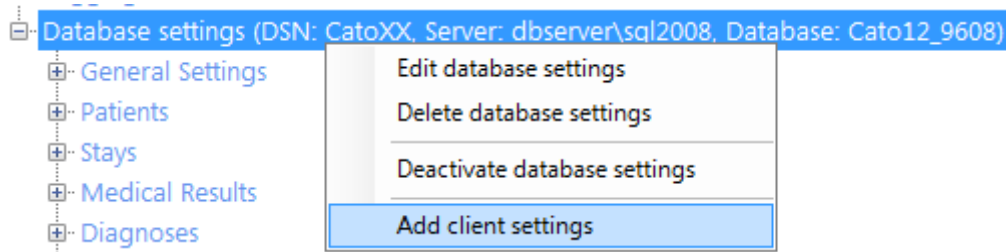

NOTE: If you don't work with clients, just enter "*", "default" or similar.



By clicking 'OK' new client settings will be added:

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 111 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

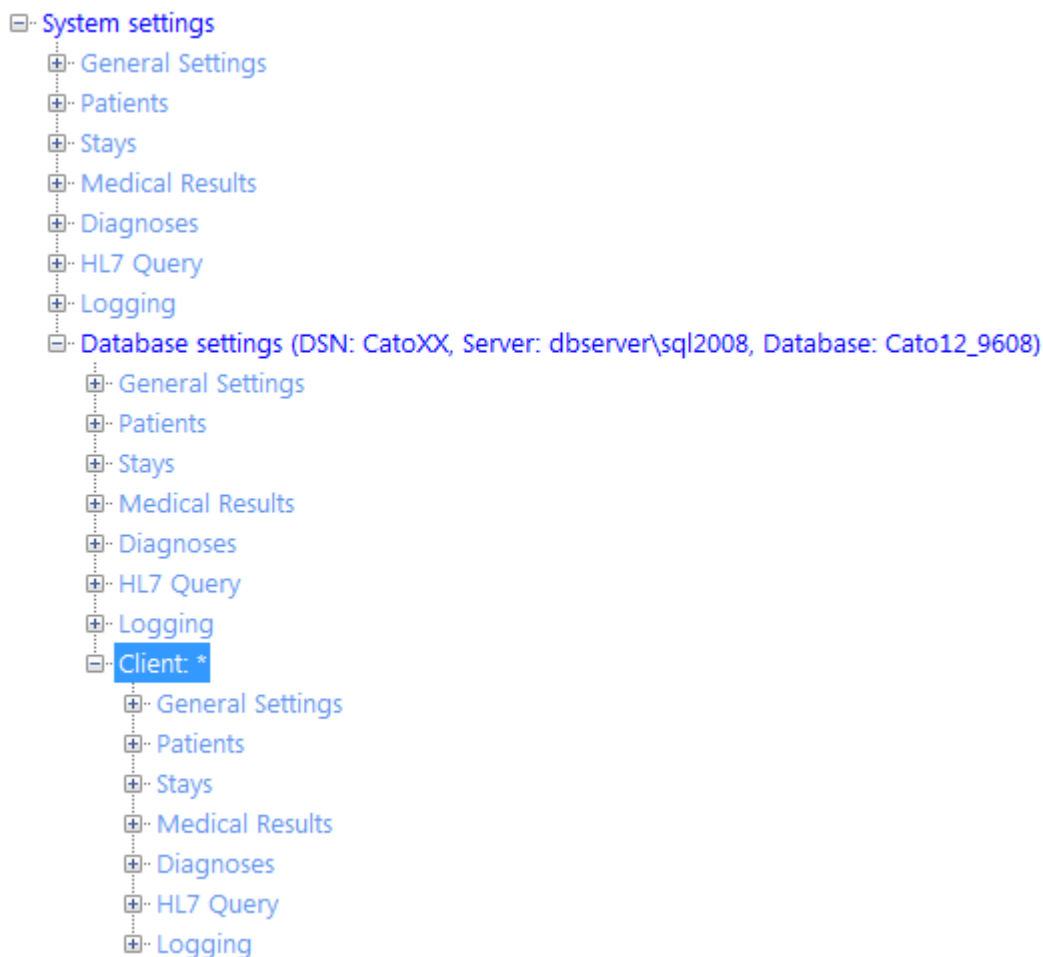

4.4.4.1.2 Edit client settings

By right-clicking on specific client settings they can be modified:

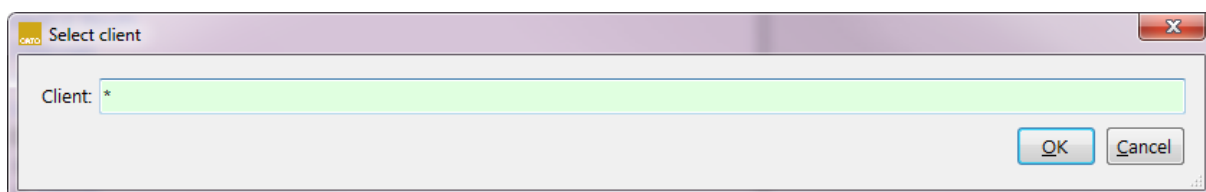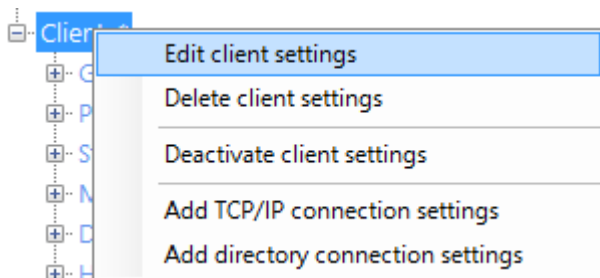


By clicking 'OK' client settings are changed.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 112 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.4.1.3 Delete client settings


By right-clicking on specific client settings they can be deleted:

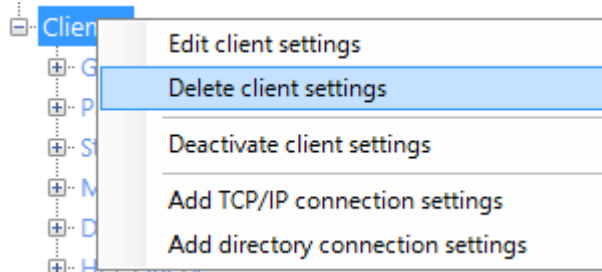


NOTE: By deleting client settings all existing sub-levels (connection settings) will also be deleted!

4.4.4.1.4 Deactivate client settings

By right-clicking on specific client settings they can set inactive (when active):

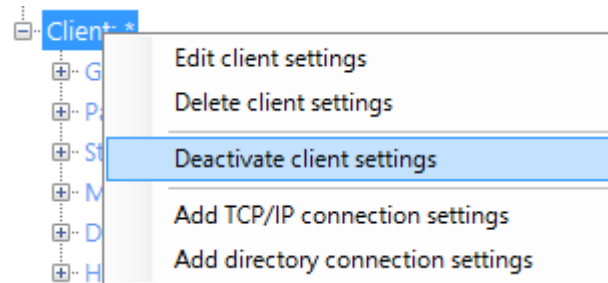


Now the client settings (and all sub-levels) are inactive:

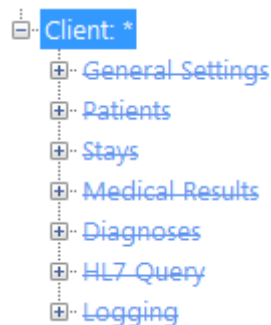

Note: This is useful for testing environments to easily deactivate client settings. After testing they can be reactivated.

4.4.4.1.5 Activate client settings

By right-clicking on specific client settings they can be set active (when inactive):

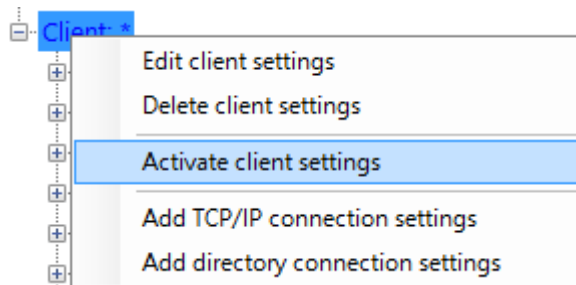

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 113 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.4.2 Settings

No additional settings!

4.4.5 Connection Settings

Two different connection settings are available:



- **TCP/IP connection settings**
- **Directory connection settings**

4.4.5.1 TCP/IP connection settings

4.4.5.1.1 Manage TCP/IP connection settings

4.4.5.1.1.1 Add TCP/IP connection settings

By right-clicking on specific client settings new TCP/IP connection settings can be added:

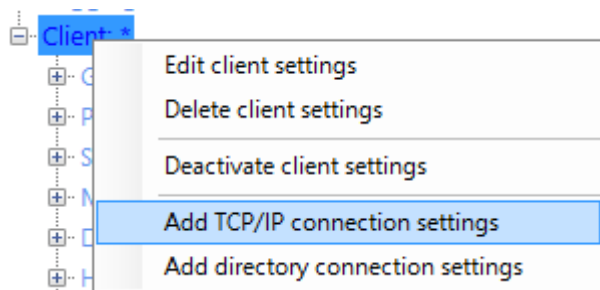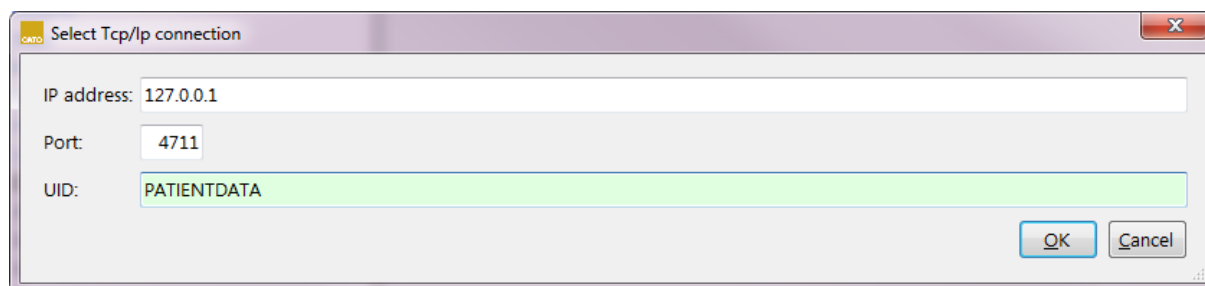Please specify a unique identifier (Unique for all connection settings)

By clicking 'OK' new TCP/IP connection settings will be added:

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 114 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

[-] System settings

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

[-] Database settings (DSN: CatoXX, Server: dbserver\sql2008, Database: Cato12_9608)

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

[-] Client: *

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

[-] IP: 127.0.0.1:4711

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses



+ HL7 Query

+ Logging

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 115 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.1.1.2 Edit TCP/IP connection settings


By right-clicking on specific TCP/IP connection settings they can be modified:

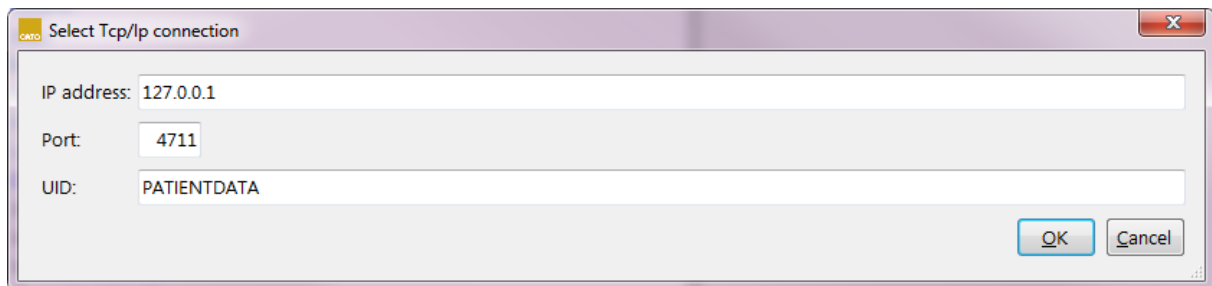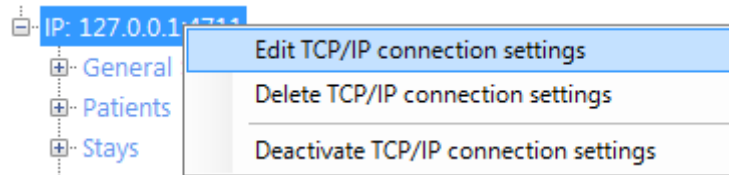

By clicking 'OK' TCP/IP connection settings are changed.

4.4.5.1.1.3 Delete TCP/IP connection settings


By right-clicking on specific TCP/IP connection settings they can be deleted:

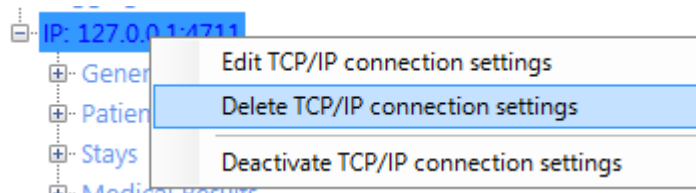


 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 116 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.1.1.4 Deactivate TCP/IP connection settings

By right-clicking on specific TCP/IP connection settings they can set inactive (when active):

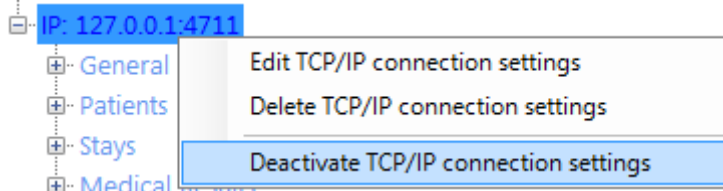


Now the TCP/IP connection settings are inactive:

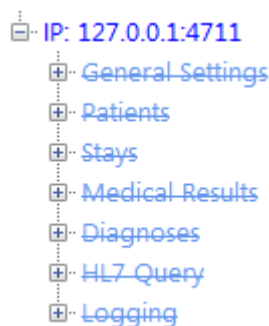

Note: This is useful for testing environments to easily deactivate TCP/IP connection settings. After testing they can be reactivated.

4.4.5.1.1.5 Activate TCP/IP connection settings

By right-clicking on specific TCP/IP connection settings they can be set active (when inactive):

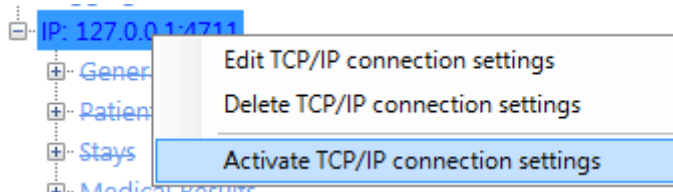

4.4.5.1.2 Settings

NOTE: Only new (added at this level) settings are listed.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 117 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.1.2.1 General Settings (Section)

4.4.5.1.2.1.1 KIR2008 – Response path

This is a special setting for a specific system where the response path can be specified.



4.4.5.1.2.1.2 LLTP Framing (Section)


LLTP framing specifies the start and end of a message in a stream.

NOTE: This setting shouldn't be changed unless the sending system provides other parameters.

4.4.5.1.2.1.2.1 Start String

This option specifies the start of a message in a stream.

Input: Hexadecimal, separated by comma or semicolon.

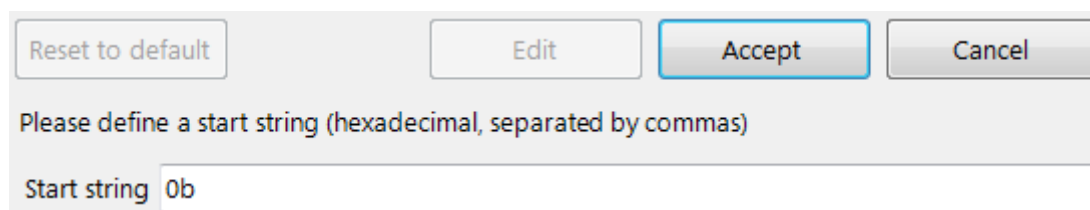


Default: 0b (Hexadecimal)

4.4.5.1.2.1.2.2 End String

This option specifies the end of a message in a stream.

Input: Hexadecimal, separated by comma or semicolon.

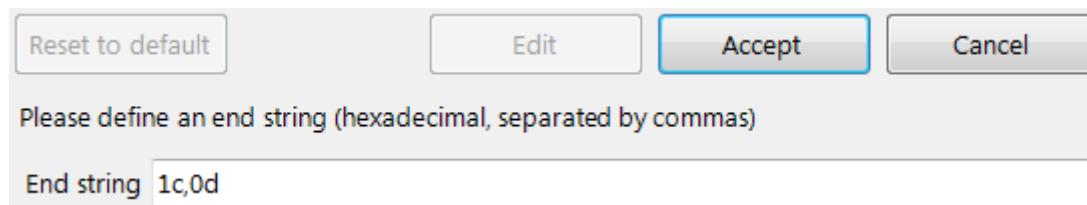

Default: 1c,0d (Hexadecimal)

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 118 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.1.2.1.3 Closing of socket after receiving a message

Specify if you want to close the socket after receiving a message.

Reset to default
Edit
Accept
Cancel

Please set this option if the link should be closed automatically after receiving a message.

☒ Closing of socket after receiving a message

Default: yes/true

Note: Some systems need the socket to be closed after receiving a message to shutdown communication. For a faster processing it is useful to test communication without closing the socket.



 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 119 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.2 Directory connection settings

4.4.5.2.1 Manage directory connection settings

4.4.5.2.1.1 Add directory connection settings

By right-clicking on specific client settings new directory connection settings can be added:

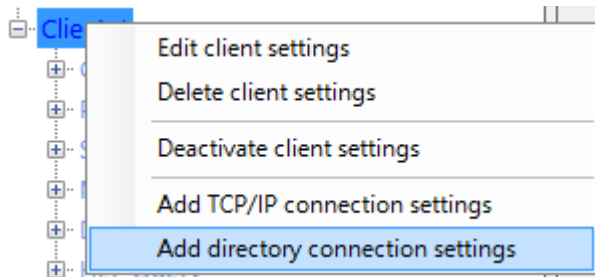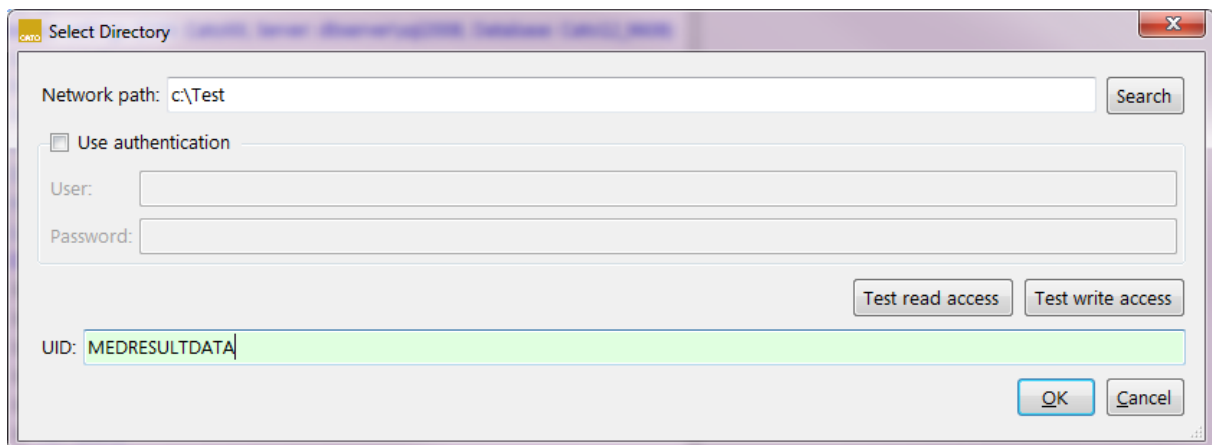Please specify a unique identifier (Unique for all connection settings)

By clicking 'OK' new directory connection settings will be added:

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 120 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

[-] System settings

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

[-] Database settings (DSN: CatoXX, Server: dbserver\sql2008, Database: Cato12_9608)

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

[-] Client: *

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses

+ HL7 Query

+ Logging

+ IP: 127.0.0.1:4711

[-] Path: c:\Test

+ General Settings

+ Patients

+ Stays

+ Medical Results

+ Diagnoses



+ HL7 Query

+ Logging

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 121 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.2.1.2 Edit directory connection settings


By right-clicking on specific directory connection settings they can be modified:

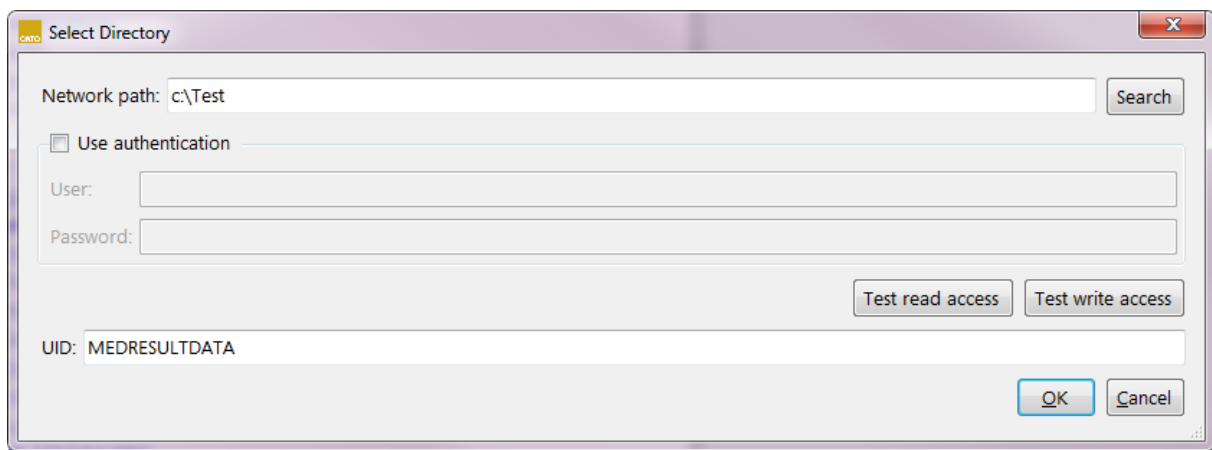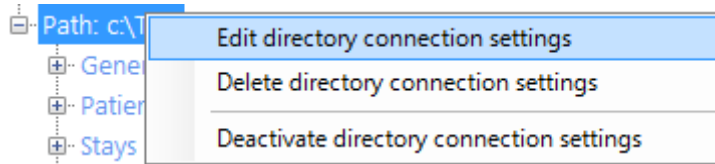

By clicking 'OK' directory connection settings are changed.

4.4.5.2.1.3 Delete directory connection settings


By right-clicking on specific directory connection settings they can be deleted:

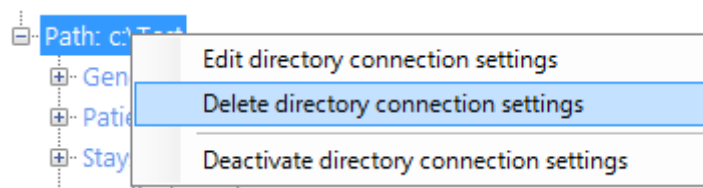


	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 122 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.5.2.1.4 Deactivate directory connection settings

By right-clicking on specific directory connection settings they can set inactive (when active):

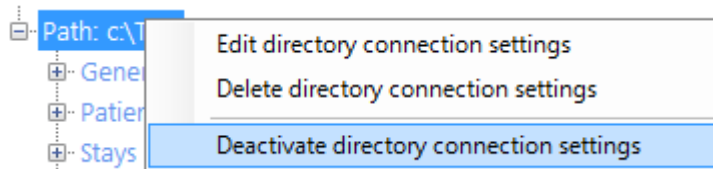


Now the directory connection settings are inactive:

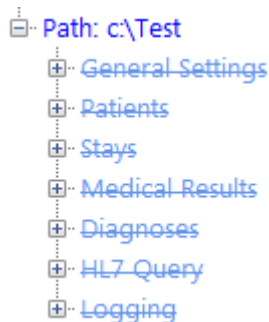

Note: This is useful for testing environments to easily deactivate directory connection settings. After testing they can be reactivated.

4.4.5.2.1.5 Activate directory connection settings

By right-clicking on specific directory connection settings they can be set active (when inactive):

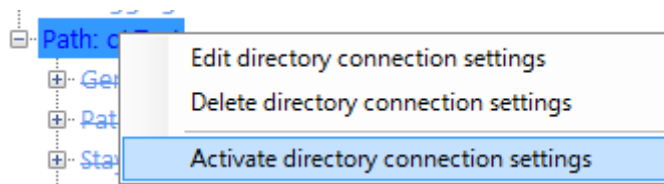

4.4.5.2.2 Settings


NOTE: Only new (added at this level) settings are listed.

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 123 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

4.4.5.2.2.1 General Settings (Section)


4.4.5.2.2.1.1 KIR2008 – Response path

This is a special setting for a specific system where the response path can be specified.

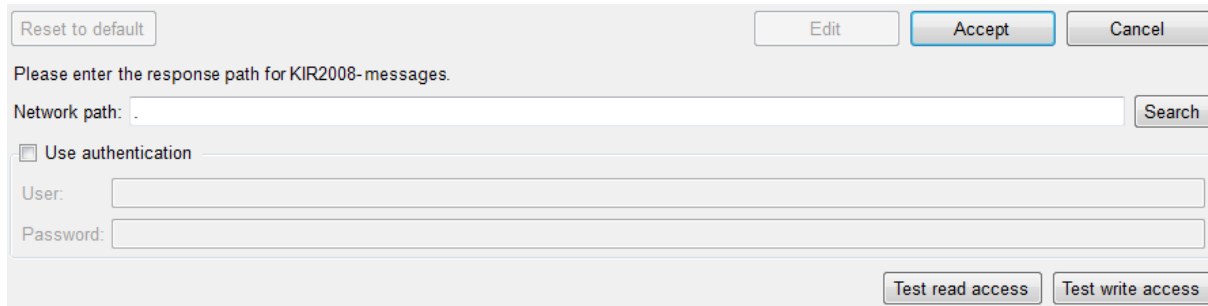

4.4.5.2.2.1.2 File type


A file type can be specified. E.g.: *.hl7, *.msg

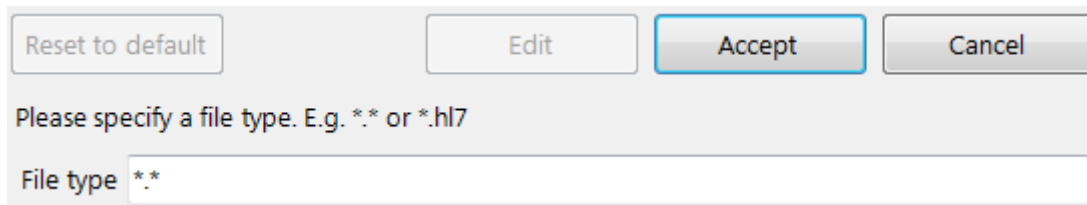

Default: *.* (All files)

4.4.5.2.2.1.3 Path for reply messages

A path for reply messages can be specified.

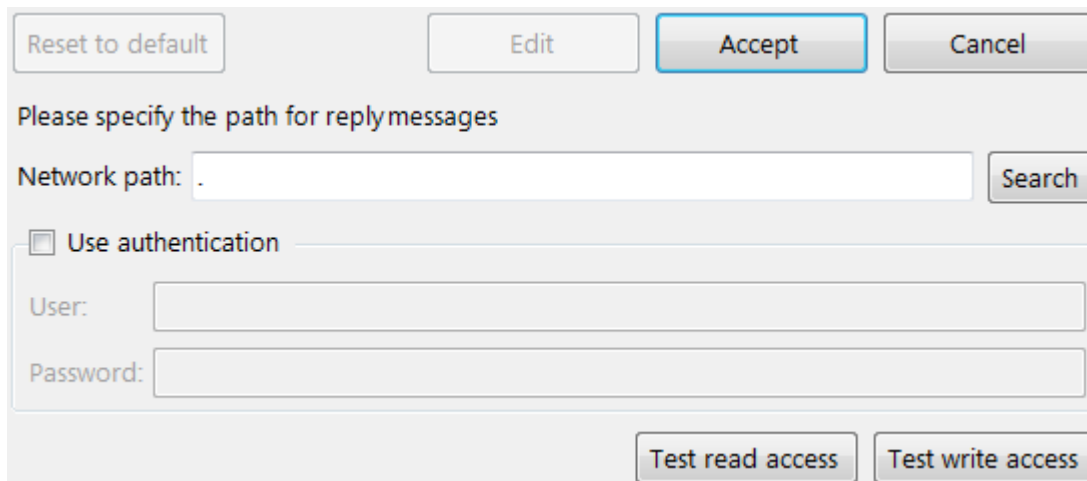

Default: .

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 124 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

4.4.5.2.2.1.4 Path for invalid messages

A path for invalid messages can be specified.

Reset to default
Edit
Accept
Cancel

Please specify the path for invalid messages

Network path: Search

☐ Use authentication

User:

Password:


Test read access
Test write access

Default: .

 Medication Workflow Solutions	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 125 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
	Template-Rev.:	02	

5 Logging Output

The output can be seen in the following tab after selected the appropriate DSN (Data source name):

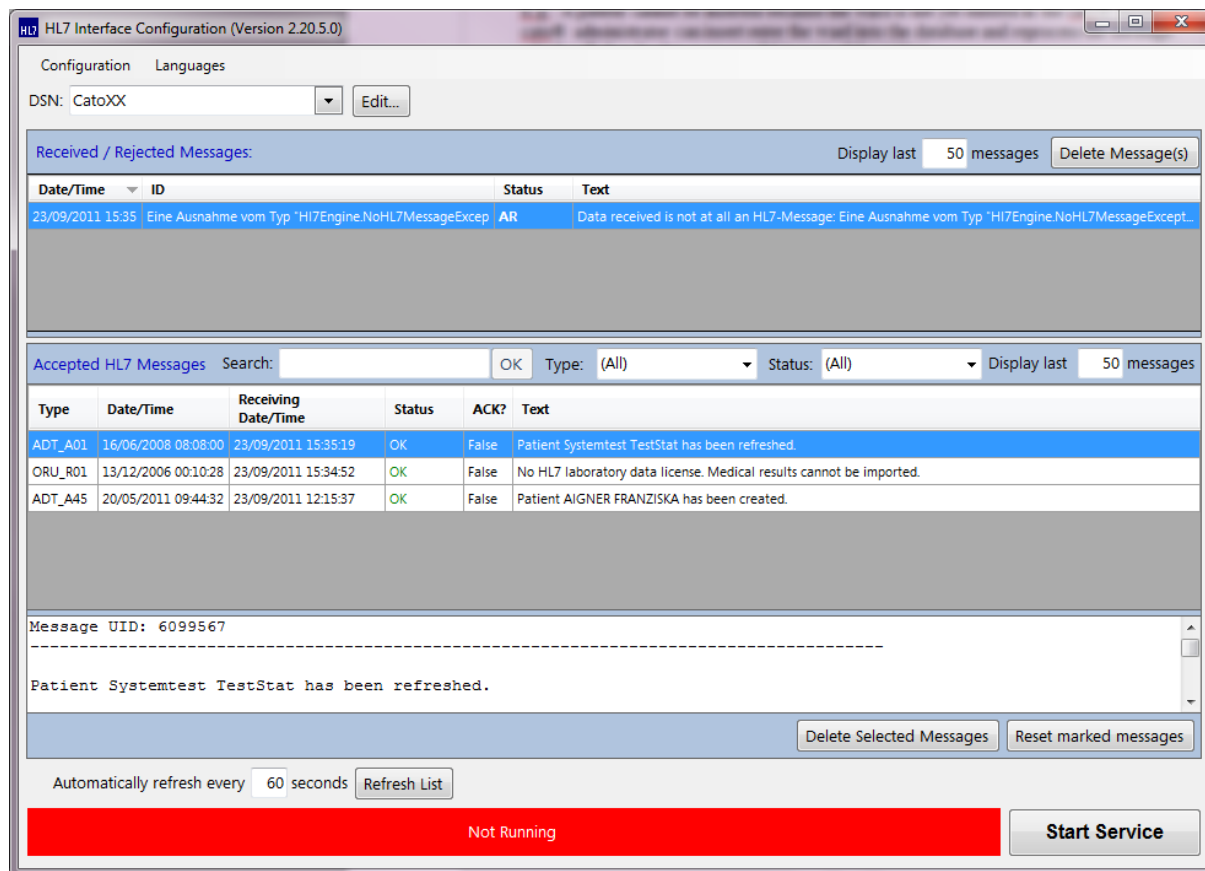

HL7 Interface Configuration (Version 2.20.5.0)

Configuration Languages

DSN: CatoXX Edit...

Received / Rejected Messages: Display last 50 messages Delete Message(s)

Date/Time	ID	Status	Text
23/09/2011 15:35	Eine Ausnahme vom Typ 'Hi7Engine.NoHL7MessageExcept	AR	Data received is not at all an HL7-Message: Eine Ausnahme vom Typ 'Hi7Engine.NoHL7MessageExcept...

Accepted HL7 Messages Search: OK Type: (All) Status: (All) Display last 50 messages

Type	Date/Time	Receiving Date/Time	Status	ACK?	Text
ADT_A01	16/06/2008 08:08:00	23/09/2011 15:35:19	OK	False	Patient Systemtest TestStat has been refreshed.
ORU_R01	13/12/2006 00:10:28	23/09/2011 15:34:52	OK	False	No HL7 laboratory data license. Medical results cannot be imported.
ADT_A45	20/05/2011 09:44:32	23/09/2011 12:15:37	OK	False	Patient AIGNER FRANZISKA has been created.

Message UID: 6099567

Patient Systemtest TestStat has been refreshed.

Delete Selected Messages Reset marked messages

Automatically refresh every 60 seconds Refresh List

Not Running Start Service

5.1 Rejected HL7 Messages

The upper list shows rejected HL7 messages which either contain not supported message types / events or which are syntactically incorrect. A Message is also rejected if it was already received, i.e. a message with the same message control ID was sent earlier.

A rejected message will be answered by an 'Application Error (AE)' or an 'Application Rejected (AR)' Acknowledgement. Please refer to [0

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 126 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

General logging settings] to set the interval after which the log entry is deleted.

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 127 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

5.2 Accepted HL7 Messages:

The lower list shows information about the processing of accepted HL7 messages.

A message can be in one of the following states:

- **OK**
Everything OK, information about the data imported is displayed on click
- **Error**
An error occurred. More information about the error is displayed on click of the erroneous message.
- **Ignored**
Everything OK, but no data has been imported due to the business logic.
E.g.: The patient was not found.
- **Case not found**
E.g.: When a discharge message (ADT-A03) is received and the case is not existent in caton® this special ignore status will be set.
- **Case event not possible**
E.g.: When a case is already closed after the discharge event (ADT-A03) an other transfer is sent, then this special ignore status will be set.

The list can be filtered by:

- **Message type**
- **Message status**
- **Text**

All error information are also logged to the windows event log or sent by mail if the appropriate options are selected (see 0

	Doc.-ID:	MWS-SW-USRMAN-CATO-HL7-2.32-EN	Page 128 of 128
	Revision:	01	
	Folder and File:	001 Product Documents\Cato HL7\2.32 EN\ User Manual Cato HL7 2.32	
	Template:	FO Other Documents (MWS-FO-OD)	
Medication Workflow Solutions	Template-Rev.:	02	

General logging settings)