Caeppas TEA





Life's obstacles surround us each and every day and we have the choice to neglect them or accept them as

opportunities. The Crossfire™ T6A State-of-the-Art frame technology and unsurpassed fabrication make all obstacles an opportunity!

The Crossfire T6A features an adjustable, open-frame design unlike any other performance rigid chairs available today. Where most adjustable rigid chairs sacrifice performance by changing the geometry of the frame relative to the ground, Top End's dual-rigid frame is designed for performance on our lower frame and positioning and comfort on the upper seat frame.

Together these two frames deliver peak performance in push efficiency and the ability to make necessary seating adjustments without

sacrificing performance over time.

The months spent engineering this design resulted in a super lightweight chair that can be adjusted in minutes precisely fit for life's opportunities.





Crossfire



The seat frame allows for ample rear seat-to-floor adjustment with a durable and eye-catching aluminum crescent link. The side frame offers 3.25" of rear wheel center-of-gravity adjustment.



Each custom frame is built to accommodate a short front fork to reduce overall caster trail and keep the user's center-of-gravity low to the ground for improved balance.






The back angle is adjusted simply by loosening a button-head screw and rotating an adjustment square located on the inside of this back bracket.



Two pressed bearings located in the fork allow for a more controlled rotation and help to eliminate flutter. The T6A also eliminates problematic caster angle adjustment with our seat frame above the side frame design.






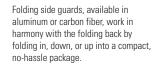




The durable folding back locks into place creating a compact package ideal for transfers and hasty airline cargo handlers!



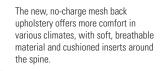



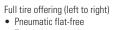

- Treaded pneumatic
- Primo
 Schwalbe®
- Knobby tireFlat-free insert





Top End® Crossfire™ T6A Active Performance Wheelchair



Frame	
Width	6.75" plus seat width w/0° camber
Seat widths	12" - 20"
Seat depths	14" - 20"
Front seat-to-floor*	16" - 21" (1/2" increments)
Rear seat-to-floor	14" - 20"
Axle positions	Active, Extended
Center-of-gravity range	3.25"
Front frame modifications	Tapered, V-front
Front frame angles	70°, 75°, 80°, 85°
Side & seat frame color choices	16 colors
Anodizing color options	Victory Red, Shadow Blue, Bold Black, Classic Silver
*Measurement taken from front to Wheels	p of seat rail to bottom of footrest
Wheel options	22", 24", 25", 26"
Tire options	Schwalbe®, Primo, Knobby, Kenda, Treaded pneumatic, Pneumatic flat free
Camber choices	0°, 3°, 6°, 9°
Spinergy wheel spoke color choice	s Red, Blue, White, Yellow, Black
Caster options	3", 4", 5"
Handrim options	Black or Silver Anodized, Natural F Titanium, Nylon Coated

Knee-to-heel* (range adjusts	in 1" increments)
16" 11.25" - 15.25"	41 + 31 (2)
17" 12.25" - 16.25"	
18" 13.25" - 17.25"	
19" 14.25" - 18.25"	
20" 15.25" - 19.25"	
21" 16.25" - 20.25"	
Positioning	
Arm type	T-Arm, Swingaway
Back type	Folding, lock-down, angle-adj. back
Back height	8"-11", 10"-14", 12"-16", 14"-18", 16"-20"
Adjustable tension back upholstery options	Breathable mesh, Ballistic nylon
Side guards	Fold-in, Folding, Removable
Weight capacity	250 lb.
Chair weight*	15.6 lb.
*16" x 16" with no rear wheels a	nd standard options

Anodizing Colors



Frame Colors



Invacare Corporation

www.invacare.com

Top End

www.topendwheelchair.com Sports and Recreation Products 4501 63rd Circle North Pinellas Park, Florida 33781-5914 (727) 522-8677 (800) 532-8677 Fax: (727) 522-1007

USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

Canada

570 Matheson Blvd. E. Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

©2008 Invacare Corporation.
All rights reserved. Trademarks are identified by the symbols TM and [®]. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.
Schwalbe is a registered trademark of RALF BOHLE GmbH Reichshof Germany.

Form No. 08-018

