

Invacare® Pronto® M41™ Power Wheelchair

POWER Custom  
REHAB  
Sports & Rec  
Seating

# Pronto M41™



Yes, you can.®





Yes, you can.®

The Invacare® Pronto® M41™ is the latest offering in the sleek and versatile Pronto® series of power wheelchairs. The M41™ features the latest technology in electronics and in-line motors. This combination provides increased efficiency and battery range, as well as an increased top speed. The M41 also features a center-wheel drive platform, simplified SureStep®, and six wheels on the ground providing stability and maneuverability similar to other Pronto® models. All of these features are packed into a simple, clean-looking, and easy-to-service compact frame that provides mobility with style for people on the go!



Pronto® M41™ with Fold-Down Seat



Pronto® M41™ with Semi-Recline Van Seat

## Features

- Compact size for a tight turning radius
- High efficiency in-line motors for greater speed and longer battery range
- SureStep® suspension with six wheels on the ground for smooth driving and optimum stability
- Charger connection conveniently located on joystick
- Stylish design
- 3-amp off-board charger
- Three seat options
  - Medium-back fold-down seat
  - High-back semi-recline van seat
  - High-back semi-recline van seat with solid seat pan
- Auto style positioning strap
- U1 batteries
- 300 lb. weight capacity


## Pronto® M41™ Specifications

Overall dimensions (measured with fold-down seat)			
Length	34"		
Width (including joystick)	24"		
Height	33"		
Seat-to-floor	19"		
Seat dimensions	(fold-down seat)	(semi-recline seat)	(semi-recline seat w/solid seat pan)
Width	18"	16" - 20"	16" - 20"
Depth	17"	16" - 18"	16" - 20"
Seat angle	92°	80° - 115°	90° - 115°
Weight capacity	300 lb.		
Wheel sizes			
Front/rear casters	6" x 2"		
Drive wheels (Flat-free)	10.5" x 3.5"		
Maximum speed	5 mph		
Motor	250 watts		
Incline capability	7° grade		
Turning radius			
Front	17.75"		
Rear	19.5"		
Weight	(fold-down seat)	(semi-recline seat)	(semi-recline seat w/solid seat pan)
with U1 batteries	161 lb.	175 lb.	180 lb.
w/out batteries and foot plate	110 lb.	124 lb.	129 lb.
Ground clearance	2"		
Battery	U1		
Battery charger	3 amp off-board		
Per-charge battery range	Up to 17 miles*		
Limited warranties			
Base frame	60 months		
Seat frame	36 months		
Electronics	12 months		
Motor/gearbox	12 months		


\* Travel range will vary depending on user weight, speed, battery charge and terrain.

## Colors Available

Candy Red

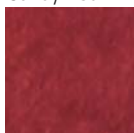

Deep Blue

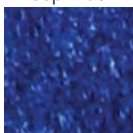

**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2008

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.

Specifications are subject to change without notification.

Form No. 06-033  
Rev. 04/08