

 \mathbf{T} he Invacare® Turn-Q™ LTM is a powered lateral turning mattress which automatically turns the user up to 30° bilaterally at preset intervals of 30, 60, or 90 minutes. The LTM reduces interface pressure below capillary closure for the treatment of stage I through stage IV pressure ulcers.

STANDARD FEATURES

- · Automatic turning and static flotation modes
- Bilaterally turns patient up to 30° (depends on patient weight and body mass)
- Turn intervals 30, 60, and 90 minutes
- Power-failure mode keeps the lower mattress inflated for up to 12 hours in the event of a power failure
- · Easy-to-understand, caregiver-friendly membrane control panel



INVACARE® TURN-Q™ LTM LATERAL TURNING MATTRESS



USER-FRIENDLY CONTROLS

- Power button with LED indicator displays whether the unit is in standby or operating mode
- · Auto-firm with LED indicator inflates the mattress to maximum firmness to facilitate patient handling
- · Patient weight adjustable from 70 350 lbs.
- Turning intervals of 30, 60, and 90 minutes, and a non-turning mode
- Auto-firm off will manually cancel the auto-firm mode; auto-firm mode will automatically return to normal operation if left unattended for 30 minutes or more




PORTABILITY

- Portable, lightweight mattress-replacement system fits most standard hospital beds
- Control unit easily attaches to top of footboard
- CPR quick-deflate handle activates rapid deflation of the mattress in approximately 30 seconds

PRESSURE RELIEF

Average pressure per pressure point for 10 subjects

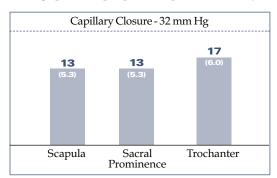

() Standard deviation

SPECIFICATIONS

POWER	115 VAC, 60 Hz, 15 watts
DIMENSIONS Standard mattress Power unit	35"W x 80"L x 11"H 20.5"W x 10.8"L x 5.6"H
WEIGHT Mattress Power unit Maximum patient weight	29 lbs. 23 lbs. 350 lbs.
SAFETY Code approval Fire	UL544 CAL 106 & 117
WARRANTY Power unit Soft goods	1-year limited 6-month limited
MATERIAL Power unit Therapy cover	Cast, machined aluminum 200-denier ripstop nylon with waterproof coating



Invacare Corporation www.invacare.com

USA One Invacare Way Elyria, Ohio 44036-4028 (800) 333-6900 Canada 5970 Chedworth Way Mississauga, Ontario L5R 3T9 Canada (800) 668-5324 Invacare, Yes, you can. and Turn-Q are trademarks of Invacare Corporation.

Specifications are subject to change without notification.

© 2000 Invacare Corporation Form No. 97-35 Rev. 05/00

^{*}Testing completed in an independent test lab